




Philips Avent
Niplette™

A non-surgical approach

For flat or inverted nipples
2 Niplettes and 2 Breast pads

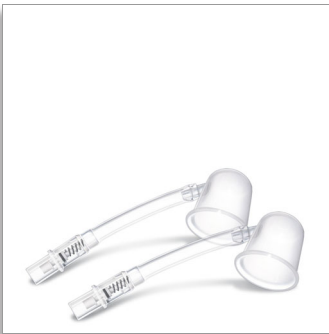

SCF152/02

Simple, effective solution for inverted nipples

Clinically proven results*

The Niplette™ helps moms with flat or inverted nipples to breastfeed. This unique and revolutionary instrument provides an easy, non surgical and long lasting solution. In a matter of weeks of continued wear, the nipple will stay erect.

Easy and comfortable

- Additional information
- A device for flat or inverted nipples
- A simple procedure

Non-surgical – gentle sucking action

- Clinically proven results

Simple and discreet - suitable for any nipple size

- Ideal use before or in the first 6 months of pregnancy

PHILIPS
AVENT

Highlights

For flat or inverted nipples




Inverted or non-protractile nipples affect up to 10% of women causing psychological distress and making breast feeding difficult for mother and baby. The suckling action of the baby should draw out the nipple. If not, the Niplette™ is a simple, comfortable solution that can help. The device makes it possible for women with flat or inverted nipples to comfortably breastfeed without the need for invasive surgery*. It consists of a transparent nipple mould with a sealing flange, attached to a valve and a syringe port.

When the use is recommended

Ideally, the Niplette should be used before pregnancy and be worn in 8-hour periods per day or night*. If breasts are not too sensitive it can also be used in the first six months of pregnancy to achieve a permanent correction or after the birth of the baby, for a few minutes before each feed. The Niplette will suck the nipple out enabling the baby to latch on easily and help establish breastfeeding during the first few days. The permanent cosmetic correction can then be effected once breastfeeding is over, if this is the case the Niplette can be applied again from time to time.

Gentle sucking action

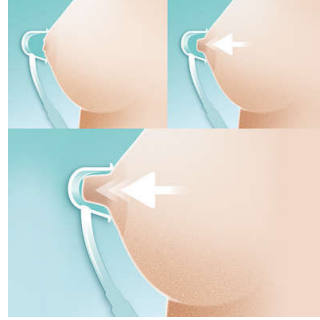

The cup is held over the nipple areola with one hand and air is withdrawn using a 5 ml syringe so that the nipple can be sucked into it. The user is in control of the suction and can pull on the nipple as firmly as comfortable. When the nipple has been pulled out the user after carefully separating the syringe from the valve can continue with their normal activities and wear the Niplette discreetly inside the bra. The Initial usage is encouraged as much as possible*.

Proven clinical success

Clinical trials* have been carried out on the Niplette proved success in achieving a permanent correction to flat or inverted nipples, whether the condition existed since puberty, or was as a result of breast reduction surgery. Provided they had not undergone breast surgery, pregnant women who had been concerned about breast feeding all went on to breast feed successfully. A permanent correction is normally achieved in between one and three months of continuous use.

Additional information

Niplette will not function once milk flow comes in fully. However by then baby will have become used to successfully latching on to the nipple.

Specifications

Development stages

- Stage: Before Pregnancy, 0 - 6 months, Breastfeeding initiation phase

Material

- Breast pads: Dermatologically tested
- Niplette: Polypropylene

What is included

- Disposable breast pad: 2 pcs

What is included

- Niplette: 2 pcs

Design

- Discreet design

Ease of use

- Easy concealing under clothing: Push syringe to release vacuum

Functions

- Solution for flat nipples: Use gentle suction

Easy to clean

- Clean in hot soapy water



Issue date 2021-03-09

Version: 18.3.1

EAN: 05 01290 90052 72

© 2021 Koninklijke Philips N.V.
All Rights reserved.

Specifications are subject to change without notice.
Trademarks are the property of Koninklijke Philips N.V.
or their respective owners.

www.philips.com

* McGeorge, Mr. D, FRCS (Plast), The Niplette: an instrument for the non-surgical correction of inverted nipples, British Journal of Plastic Surgery (1994) Vol 47, Pages 46-49

* This section contains consumer opinions on the product. Philips dissociates itself from the contents entered by consumers in this section and consequently any technical information and / or advice on the use of the product included therein are not intended as official Philips information.