




Philips Avent  
Nipplette™

**A non-surgical approach**

For flat or inverted nipples  
1 Nipplette and 2 Breast pads

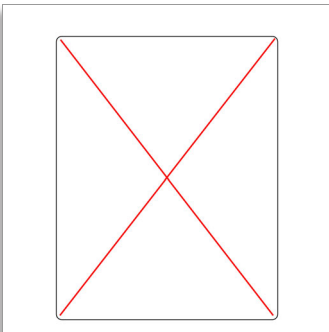

SCF152/01

## Simple, effective solution for inverted nipples

### Clinically proven results\*

The Nipplette™ helps moms with flat or inverted nipples to breastfeed. This unique and revolutionary instrument provides an easy, non surgical and long lasting solution. In a matter of weeks of continued wear, the nipple will stay erect.

#### **Easy and comfortable**

- Additional information
- A device for flat or inverted nipples
- A simple procedure

#### **Non-surgical – gentle sucking action**

- Clinically proven results

#### **Simple and discreet - suitable for any nipple size**

- Ideal use before or in the first 6 months of pregnancy

**PHILIPS**  
**AVENT**

# Specifications

## Easy to clean

- Clean in hot soapy water

## Design

- Discreet design

## Material

- Breast pads: Dermatologically tested
- Niplette: Polypropylene

## What is included

- Disposable breast pad: 2 pcs
- Niplette: 1 pcs

## Ease of use

- Easy concealing under clothing: Push syringe to release vacuum

## Functions

- Solution for flat nipples: Use gentle suction

## Development stages

- Stage: Before Pregnancy, 0 - 6 months, Breastfeeding initiation phase

# Highlights

## For flat or inverted nipples

Inverted or non-protractile nipples affect up to 10% of women causing psychological distress and making breast feeding difficult for mother and baby. The suckling action of the baby should draw out the nipple. If not, the Niplette™ is a simple, comfortable solution that can help. The device makes it possible for women with flat or inverted nipples to comfortably breastfeed without the need for invasive surgery\*. It consists of a transparent nipple mould with a sealing flange, attached to a valve and a syringe port.

## When the use is recommended

Ideally, the Niplette should be used before pregnancy and be worn in 8-hour periods per day or night\*. If breasts are not too sensitive it can also be used in the first six months of pregnancy to achieve a permanent correction or after the birth of the baby, for a few minutes before each feed. The Niplette will suck the nipple out enabling the baby to latch on easily and help establish breastfeeding during the first few days. The permanent cosmetic correction can then be effected once breastfeeding is over, if this is the case the Niplette can be applied again from time to time.

## Gentle sucking action

The cup is held over the nipple areola with one hand and air is withdrawn using a 5 ml syringe so that the nipple can be sucked into it. The user is in control of the suction and can pull on the nipple as firmly as comfortable. When the nipple has been pulled out the user after carefully separating the syringe from the valve can continue with their normal activities and wear the Niplette discreetly inside the bra. The Initial usage is encouraged as much as possible\*.

## Proven clinical success

Clinical trials\* have been carried out on the Niplette proved success in achieving a permanent correction to flat or inverted nipples, whether the condition existed since puberty, or was as a result of breast reduction surgery. Provided they had not undergone breast surgery, pregnant women who had been concerned about breast feeding all went on to breast feed successfully. A permanent correction is normally achieved in between one and three months of continuous use.

## Additional information

Niplette will not function once milk flow comes in fully. However by then baby will have become used to successfully latching on to the nipple.



Issue date 2023-10-05

Version: 18.3.1

© 2023 Koninklijke Philips N.V.  
All Rights reserved.

Specifications are subject to change without notice.  
Trademarks are the property of Koninklijke Philips N.V.  
or their respective owners.

[www.philips.com](http://www.philips.com)

\* McGeorge, Mr. D. FRCS (Plast), The Niplette: an instrument for the non-surgical correction of inverted nipples, British Journal of Plastic Surgery (1994) Vol 47, Pages 46-49