

HeartStart Event Review Pro User Guide

PHILIPS

About this edition

HeartStart Event Review Pro version 4.2

Publication date: November 2011

Document part number: 453564274562

The information in this document applies to the product version indicated above. This information is subject to change without notice.

Philips Healthcare shall not be liable for errors contained herein or for incidental or consequential damages in connection with the furnishing, performance, or use of the material.


Copyright © 2011

Philips Electronics North America Corp.

No part of this publication may be reproduced, transmitted, transcribed, stored in a retrieval system, or translated into any human or computer language in any form by any means without the consent of the copyright holder.

Unauthorized copying of this publication may not only infringe copyright but also reduce the ability of Philips Healthcare to provide accurate and up-to-date information to users and operators alike.

Medical device directive

Event Review Pro complies with the requirements of the Medical Device Directive 93/42/EEC and carries the  mark accordingly.

Trademarks

Philips HeartStart Configure, HeartStart FR3, HeartStart FRx, HeartStart HS1, HeartStart Defibrillator, and the HeartStart logo are either trademarks or registered trademarks of Koninklijke Philips Electronics N.V.

The Event Review Pro application uses *Bluetooth* wireless technology. The *Bluetooth* wordmark and logos are either trademarks or registered trademarks of *Bluetooth* SIG, Inc. The *Bluetooth* wordmark and logos are owned by *Bluetooth* SIG, Inc.

Acrobat, Adobe Reader, and PDF are either trademarks or registered trademarks of Adobe Systems Incorporated in the United States and/or other countries. Microsoft®Windows, Windows XP, Windows Vista, Windows 7, Windows Explorer, Microsoft .NET Framework, and Internet Explorer are either trademarks or registered trademarks of Microsoft® Corporation in the United States and/or other countries.

Actiontec is either a trademark or a registered trademark of Actiontec Electronics, Inc. Adobe, Acrobat, Acrobat Reader, and PDF are either trademarks or registered trademarks of Adobe Systems Incorporated in the United States and/or other countries. CompactFlash is either a trademark or a registered trademark of ScanDisk Corporation. IrDA and IR logo are trademarks or registered trademarks of Infrared Data Association. IVT is either a trademark or a registered trademark of Interactive Video Technologies, Inc. Kingston is either a trademark or a registered trademark of Kingston Technologies Company. Microsoft, Windows, Internet Explorer, Excel, Word, and Outlook are either trademarks or registered trademarks of Microsoft Corporation in the United States and/or other countries. Palm and HotSync are among the trademarks or registered trademarks owned by or licensed to Palm, Inc. Pentium is either a trademark or a registered trademark of Intel Corporation or its subsidiaries in the United States and other countries. Q-CPR is either a trademark or a registered trademark of Laerdal Medical. Toshiba is either a trademark or a registered trademark of TOSHIBA Corporation. Widcomm is either a trademark or a registered trademark of Broadcom Corporation.

Authorized EU representative



Philips Medizin Systeme
Boeblingen GmbH
Hewlett-Packard Strasse 2
71034 Boeblingen, Germany
(+49) 7031 463-2254

Device manufacturer



Philips Medical Systems
22100 Bothell Everett Highway
Bothell, WA 98021-8431, USA



Table of contents

Installing Event Review Pro	8
System requirements.....	9
Downloading the application.....	12
Checking for updates.....	12
The installation procedure.....	13
Activating the software.....	14
Activate over the Internet.....	15
Activate by email.....	15
Activate later.....	17
Registering the application.....	17
Uninstalling Event Review Pro.....	17
Product compatibility.....	18
Getting started	19
Intended use.....	19
What's new in this release.....	19
Getting assistance.....	20
Conventions used in this guide.....	21
Starting the application	23
Understanding the application window.....	23
Completing fields.....	25
Using the application tables and logs.....	25
Resizing panes and workspaces.....	27
Managing mouse behavior.....	27
Using a Tablet PC.....	27
Saving your work.....	28
Restoring confirmation messages.....	28
Using the Navigation pane.....	28
Working with cases	30
Understanding the Cases navigation pane.....	30
Creating a case (including ECG) with the Case Wizard.....	31
Creating a case.....	33
Adding ECGs.....	34

Adding ECGs table.....	35
Adding an ECG from a data card.....	36
Adding an ECG from a USB port.....	37
Downloading an ECG using an FR3 Bluetooth transmission.....	37
Downloading an ECG using an MRx Bluetooth transmission.....	38
Downloading an ECG from an infrared connection.....	39
Adding an ECG from a file.....	40
After you add an ECG.....	41
Adding multiple ECGs.....	41
Removing ECGs from a case.....	42
Displaying case details.....	43
Sorting and grouping cases.....	43
Hiding and displaying cases.....	44
Importing case files manually.....	44
Importing case files automatically.....	46
Saving cases.....	46
Adjusting the date and time of the defibrillator data.....	46
Erasing the data source.....	47
Printing case information.....	49
Exporting cases.....	49
Emailing a case.....	51
Deleting cases.....	52
Reviewing cases with duplicate ECGs.....	52
Adding case details.....	54
Identifying the case.....	55
Describing the conditions at the scene.....	60
Documenting events.....	60
Adding, describing, and removing events.....	61
Hiding and displaying event details.....	62
Sorting events.....	62
Changing the date and time of the event.....	62
Documenting the patient's outcome.....	63
Working with ECGs.....	65
Working in the ECG tab.....	65
Using the Events pane.....	68
Viewing waveforms.....	71
Working with waveforms.....	72
Magnifying waveforms.....	73
Displaying and hiding channels.....	73

Arranging sections of the ECG display.....	73
Changing the display time.....	74
Changing the time scale.....	75
Scrolling through the channel.....	75
Customizing channels.....	76
Changing the waveform offset.....	77
Changing the waveform scale.....	78
Reviewing channel information.....	79
Defining and viewing ECG selections.....	81
Creating ECG selections.....	81
Managing ECG selections.....	82
Naming, finding, and sorting ECG selections.....	83
Using selection templates.....	84
Copying ECGs to the Clipboard.....	85
Exporting waveform data.....	86
Managing notes on the waveform.....	86
Zooming in and out of the waveform.....	87
Customizing the ECG display.....	88
Reviewing vital trends data.....	89
Working with 12-lead ECGs.....	90
Configuring the 12-Lead view.....	91
Reviewing the 12-lead ECG.....	92
Magnifying 12-lead ECGs.....	93
Customizing Q-CPR episodes.....	93
Creating Q-CPR exclusions.....	94
Removing Q-CPR exclusions.....	95
Adding notes to the Ventilation channel.....	95
Using key command and mouse shortcuts.....	96
Key command shortcuts.....	96
Mouse shortcuts.....	98
Reviewing case details.....	100
Reviewing case events.....	100
Attaching (adding) files to a case.....	102
Working with reports.....	106
Generating reports.....	107
Generating a report for an open case.....	107
Generating a report from the database.....	107
Using the Reports toolbar.....	108
Working with case reports.....	109

Working with ECG reports.....	110
Working with Q-CPR reports.....	111
Working with Q-CPR Details reports.....	112
Q-CPR report data.....	112
Working with Vital Trends reports.....	114
Working with Response Times reports.....	114
Working with Utstein reports.....	116
Printing reports.....	117
Exporting reports.....	118
Emailing reports.....	119
Working with the System Log.....	120
Working with columns.....	121
Grouping and sorting entries.....	121
Filtering entries.....	122
Removing filters.....	125
Printing System Log entries.....	125
Using the Import Service.....	127
Managing Import Service inboxes.....	127
Configuring the Import Service.....	128
Managing Import Service archives.....	129
Using the HTTP Import Service.....	130
Customer support.....	132
Philips Response Center telephone numbers.....	133
Comments or suggestions?.....	135
Supported help.....	135
Unsupported help.....	136
Helping us help you.....	136
Working with defibrillators.....	137
Supported defibrillators.....	137
Selecting accessories for data transfer.....	138
Choosing data transfer adapters and readers.....	138
Quick reference to data cards.....	139
Choosing Bluetooth accessories.....	140
Reading FR3 series cards.....	140
Retrieving FR3 data using Bluetooth transmission.....	140
Setting up Bluetooth transmissions for the MRx and FR3.....	141
Bluetooth option prerequisites for the MRx and FR3.....	141
Pairing and testing the Bluetooth option with the computer: MRx and FR3.....	142
MRx.....	142

FR3.....	143
Reading HeartStart MRx cards.....	144
Sending the MRx Bluetooth transmission.....	144
Reading FR2 series cards.....	145
Reading XL cards.....	146
Transferring XL+ data.....	146
Determining the HS1 and FRx case date and time.....	147
Using infrared connections for the HS1 and FRx.....	148
Setting up an infrared adapter.....	148
Setting up the infrared connection.....	148
Understanding voice and system messages.....	151
Emailing device history data.....	151
Managing the database.....	153
Using the Event Review Pro database on a remote database server.....	153
Changing the database server.....	154
Using Microsoft SQL.....	155
Migrating cases from previous versions.....	157
CPC and OPC.....	160
Glasgow Coma score.....	162
Glossary.....	164
Index.....	171

Installing Event Review Pro

You can install HeartStart Event Review Pro 4.2 as a new installation or as an upgrade from Event Review Pro 4.0/4.1. If you install Event Review Pro 4.2 as a new installation, you can install and run the new version for 60 days while you purchase a license. Event Review Pro stops working if you do not activate the software within 60 days of installation.

You can install Event Review Pro on the same computer as Event Review Pro 3.5 or Event Review Pro 4.0/1. Event Review Pro does not remove Event Review Pro 3.5, but it does remove other, more recent versions. You can migrate case data from previous versions after installation.

The following table lists the upgrade path from earlier versions of Event Review Pro to Event Review Pro 4.2. Note that, if you are using a version of Event Review Pro before 3.5, you need to upgrade to 3.5 *before* upgrading to 4.2

From	Install and upgrade to	Download and install
Event Review Pro 3.0/3.2	Event Review Pro 3.5	Event Review Pro 4.2
Event Review Pro 3.5		Event Review Pro 4.2
Event Review Pro 4.0/4.1		Event Review Pro 4.2

- For full details about hardware and software requirements, see [System requirements on page 9](#).
- For help with downloading, see [Downloading the application on page 12](#).
- For help with installing, see [The installation procedure on page 13](#).
- For help setting up a remote database for use with Event Review Pro, see [Using the Event Review Pro database on a remote database server on page 153](#).
- For help with migrating data from an earlier version of the software to the current version, see [Migrating cases from previous versions on page 157](#).

System requirements

Event Review Pro operates on a computer running either the Microsoft® Windows XP™ or Microsoft® Windows 7™ operating system. Make sure that any hardware you choose is certified as Microsoft® - compatible.

Philips provides only the Event Review Pro software. The software and hardware listed in the following table are provided by the customer, unless noted.

This table has the following sections:

- Software requirements
- Hardware requirements
- Accessories

Software requirements table

Component	Requirement
Operating system	Microsoft Windows XP Professional SP 3, or Microsoft Windows XP tablet Edition SP3 Microsoft Windows 7 Microsoft Windows® Server® 2008 R2, 2008, 2003 R2, or 2003 for a remote server running the shared database
Database	Microsoft SQL Server 2008 R2, 2008, or 2005


Hardware requirements table

Component	Requirements
Processor speed	Minimum: 1 GHZ Recommended: 2 GHZ core duo or higher
Display	Minimum: 1024 x 768 Recommended: 1400 x 1050 or higher
Video memory	Minimum: 64 MB video memory Recommended: 256 MB video memory

Component	Requirements
Memory	Minimum: 1 GB Recommended: 2 GB or higher
Disk space	Minimum: 2 GB of disk space Recommended: 5 GB of disk space Hard-disk space requirements vary depending on usage and defibrillator type. Variables affecting disk-space requirements include the number of cases archived and the amount of audio information archived. For example, a 15-minute FR2 ECG with no audio is approximately 100 KB. The same ECG with audio can exceed 5 MB. An ECG and audio from a HeartStart MRx defibrillator can be as large as 256 MB.
Internet connection	Required to activate the application software, to use the Email feature, and to receive software updates.

Accessories table

Component	Requirement
Printer	To print reports Preferred: Color printer
PDF Reader	To view <i>HeartStart Event Review Pro user guide</i> Recommended: Adobe Reader, latest version. For more information, see the following Web site: http://get.adobe.com/reader
Backup and restore tool	To help prevent data loss and corruption
Connection to a time server	To synchronize your computer date and time Note: If you use a time server, do not alter your system clock manually.
Email application	To activate the software by email or to send cases using email, and to configure a MAPI-compliant email client
Spreadsheet application	To display Q-CPR data in charts; data is exported in .xls format only

Component	Requirement
Sound card and speaker	To play audio from the defibrillator
Card reader	<p>To read information from a defibrillator, see your Philips customer representative to purchase the appropriate card reader.</p> <p>See Selecting accessories for data transfer on page 138 for a summary of the major card reader types.</p>
IrDA support	<p>Required to read information from an HS1 or FRx defibrillator</p> <p>Your computer must support the following technology:</p> <ul style="list-style-type: none"> • IrDA functionality • An infrared transceiver or an infrared adapter <p>The transceiver will appear as a small, dark red or black window on the computer and device.</p>
<i>Bluetooth</i> adapter	To transfer HeartStart MRx and HeartStart FR3 data to your computer
<i>Bluetooth</i> stack	<p>You can use <i>Bluetooth</i> to transfer MRx or FR3 data from the defibrillator to your computer.</p> <p><i>Bluetooth</i> version: 2.0 or higher</p> <hr/> <p> You must have a supported adapter for Bluetooth transmission. Most Bluetooth adapters, such as IO Gear GBU421) are suitable. You can order an adapter from your Philips representative.</p> <hr/> <p>The <i>Bluetooth</i> protocol varies depending on your Philips defibrillator. It becomes more complicated as the <i>Bluetooth</i> stacks from different vendors implement different subsets of the protocols.</p> <ul style="list-style-type: none"> • If you are working with both the MRx and the FR3, the Widcomm stack currently can support both at the same time. • If you are working only with the FR3, the Windows built-in Microsoft stack is the preferred stack and provides the best performance. • If you are working only with the MRx, many other stacks(such as the Toshiba stack that comes with the Panasonic Toughbook) work well, in addition to the Widcomm stack. Consult the MRx documentation for information on other supported stacks. One notable exception is the Microsoft stack, which does not have the protocol implementation that can support the MRx.

Component	Requirement
	The <i>Bluetooth</i> stack requirement dictates the <i>Bluetooth</i> transceiver/dongle that you should use.

Downloading the application

You can download and use Event Review Pro on a trial basis for 60 days before purchasing and activating the software.

Use the authorization that you receive from your sales representative to download the application from the Internet. Save the application file to your computer. You can use the application on a trial basis and activate the software later.

To download the application

1. Start **Windows**.
2. Start **Internet Explorer**.
3. Navigate to the URL provided on your Proof of Purchase Certificate or by your sales representative.
4. On the left pane, select the product.
5. On the right pane, click the link in Software Download. Navigate to the download link.
6. Save the installation file to your computer.

For information on possible updates to the application, see [Checking for updates on page 12](#).

For instructions on installation, see [The installation procedure on page 13](#).

Checking for updates

You can check at any time to see if new updates have been released for Event Review Pro.



Before you install any software updates on the computer, make sure that you have Windows Administrator privileges for the computer.

To check for updates

1. On the Help menu, click **Check for Updates**.
The application connects to the Web site.
2. Click **Show Updates**.
The Web site displays a list of available updates. If a new update is available, follow the rest of this procedure to download and install it.
3. Click **Get Update for the Application**.
4. Click **Download** and follow the instructions to save the update.
5. Note the name and location of the saved download file and follow the instructions to install it.
6. When the file finishes installing, close the window.
7. Close the older version of the application.

The installation procedure

Once you have saved the installation file to your computer (see [Downloading the application on page 12](#)), follow these instructions.

To install Event Review Pro

1. In the folder where the installation file has been saved, double-click the file.
You may see a User Account Control box that asks you to authorize this program to make changes. Click **Yes**.
2. The Install Shield wizard checks for your database and for versions of the Microsoft .NET Framework. If it needs to install or upgrade either of these, it does so. Follow any instructions during this process (which can take some time); if the User Account Control box asks you to authorize this program to make changes, click **Yes**.
3. The Install Shield wizard dialog box is now ready to begin installation of Event Review Pro. Click **Install**.
4. On the **Welcome** screen, click **Next**.
5. On the next page, click **I Accept the Terms....** and click **Next**.
6. The next page displays the current readme file, containing any late-breaking information. When you have read this, click **Next**.

7. On the next page, click the edition of Event Review Pro that you want to install. The EMS edition includes case information specific to nonhospital settings; the Hospital edition includes case information, plus extra reports, that are specific to hospital settings. Click **Next**.
8. In the Custom Setup dialog box, you have the opportunity to save Event Review Pro in a different location from the default. If you need to check this new location for available space, click the **Space** button. If you want to choose a different location, click **Change** and navigate to the new location. When you are finished, or if you want to accept the default location, click **Next**.
9. Confirm that you want to continue install; click **Install**.
10. When the installation is finished, click **Finish**.
11. An entry for Event Review Pro now appears on the Start menu. See [Starting the application on page 23](#).

Activating the software

After you install the software, activate the product.

You can use the Internet or email to activate the software at any time after you install the application. Until you complete the activation, you will see the Philips HeartStart Activation Wizard window when the application starts. At the end of 60 days, however, if you have not yet activated, you cannot run the application on this computer.

You can view the serial number by clicking the Help menu and clicking **About** Event Review Pro. In addition, save the serial number in a safe place that you can access in case you need to re-install the software.

Email activation support is available in English only at:
activation.support@philips.com

Product support via email is available in English only at:
eventreview.support@philips.com

For complete information on customer support, see [Customer support on page 132](#).

For telephone support in English only, call the following numbers between 9:00 AM and 5:00 PM, Pacific Time:

- (800) 263-3342, inside the United States
- +1.978.659-3332, outside the United States

Activate over the Internet

You can activate the software by the Internet at any point when you first install the application, preferably within 60 days of installation.

At the end of the 60 days, however, if you have not yet activated, you cannot run the application on this computer

The activation wizard sends encoded information to the Philips server to validate that the serial number you enter has not already been used for the allowed number of installations and to activate your software. If you are not connected to the Internet, the wizard alerts you that there is no connection.

To activate the software by the Internet

1. Start the application.
The application opens HeartStart Activation Wizard.
2. Type the serial number that you received.
3. Click **I want to activate the software over the Internet**.
4. Click **Activate**.

The wizard starts the application.

Activate by email

You can request to activate the software by email at any time after you first install the application—preferably within 60 days of installation. When you use the activation wizard, it attaches an encrypted file when it sends the activation request. Customer Support sends you the activation code.

An email application must be installed on this computer, or you must have access to an installed email application. You can continue to use the software for up to 60 days before you complete the activation. At the end of the 60 days, however, if you have not yet activated, you cannot run the application on this computer.

If an email application is not installed on this computer, contact Customer Support.

Email activation support is available in English only at:
activation.support@philips.com

Product support via email is available in English only at:
eventreview.support@philips.com

For complete information on customer support, see [Customer support on page 132](#).

For telephone support in English only, call the following numbers between 9:00 AM and 5:00 PM, Pacific Time:

- (800) 263-3342, inside the United States
- +1.978.659-3332, outside the United States

The Activation Wizard online Help includes the steps on how to activate the software by email. Print the Help for future use.



The Activation Wizard online Help includes the steps on how to activate the software by email. Print the Help for future use.

To activate the software by email

1. Start the application.
The application opens HeartStart Activation Wizard.
2. Type the serial number that you received.
3. Click **I want to activate the software by email**.
4. Click **Next**.
5. Click **Request Activation Code** to open your email application with a pre-addressed email activation request.
6. Click **Send**.
Customer Support will send you an email with the activation code.

When you receive an email from Customer Support, complete the following:

1. Start the application.
The application opens the HeartStart Activation Wizard.
2. Verify or type the serial number that you received.
3. Click **I want to activate the software by email**.
4. Click **Next** to access the Activate by Email page.
5. Type the activation code in the **Enter your activation code here** field.
6. Click **Activate**.
The wizard validates your entry and displays a congratulations message.
7. Click **OK**.
The wizard starts the application.

Activate later

If you do not want to activate your copy of the application when you start it, you can activate the software at any time when you start the application. At the end of the 60 days, however, if you have not yet activated, you cannot run the application on this computer.

During the 60 days, you can use the application and save your work.

Until you activate the software, you will see the HeartStart Activation Wizard window each time that you start the application.

To use the application without activating the software

1. Start the application.
2. When HeartStart Activation Wizard opens, click **Skip**.

The wizard starts the application.

Registering the application

You can register your copy of the application using the Internet. Registration is optional and is not the same as the activation, which is required.

To register the software

1. Start Windows.
2. Start Internet Explorer.
3. Navigate to the following URL
http://www.philips.com/software_registration
4. Click a language.
5. On the language page, click the product.
6. Complete the registration form.

Uninstalling Event Review Pro

If you need to remove the application, in the Windows Control Panel, click **Programs** and then **Uninstall a Program** (Windows 7) or click **Add/Remove Programs** (Windows XP). Then click the name of the application.

If you have an active Internet connection at the time you remove the application, the removal process contacts the Philips activation server and releases your serial number for installation on another computer. The removal process does not remove the application database.



If you had Event Review Pro 4.0 or 4.1 installed, when you install Event Review Pro 4.2, the earlier version is automatically uninstalled. However, if you had Event Review Pro 3.5 installed and do not want to run both versions on the same machine, you must manually uninstall it.

Product compatibility

You can install HeartStart Event Review Pro 4.2 on the same computer as Event Review Pro 3.5 or Event Review Pro 4.0/4.1. Event Review Pro does not remove Event Review Pro 3.5, but it does remove Event Review Pro 4.0/4.1.

You can migrate the Event Review Pro 3.5 and Event Review Pro 4.0/4.1 data to Event Review Pro 4.2 after installing HeartStart Event Review Pro 4.2. For more information, see [Managing the database on page 153](#).

Getting started

In Event Review Pro, a "case" consists of all the information relating to a specific cardiac emergency. This includes information recorded by supported defibrillators as well as information entered by the device user. Event Review Pro creates a case by downloading ECG information, shock decisions, and any recorded audio from Philips defibrillators. You can add details to the case as well.

Intended use

HeartStart Event Review Pro captures, reports, and manages cardiac arrest information. Event Review Pro has tools for collecting and analyzing information, reviewing cases, identifying trends, and evaluating the emergency response. All information relating to a specific cardiac emergency is created in a case by downloading ECG information, shock decisions, and recorded audio from Philips and Laerdal automated external defibrillators (AED). Event Review Pro also allows users to annotate case details.

In addition, Event Review Pro captures and allows review of more comprehensive response and patient data, including Basic Life Support (BLS) and Advanced Life Support (ALS) responder observations and interventions. It also manages 12-lead ECG information and monitoring data collected by the Philips HeartStart MRx and XL Monitors/Defibrillators.

What's new in this release

Event Review Pro 4.2 provides the following new features:

- Support for the Philips HeartStart FR3 defibrillator.
- Support for the Philips HeartStart MRx monitor/defibrillator, version F.02/R.02
- Support for the Philips HeartStart XL+ defibrillator
- Upgraded database support: SQL Express 2008 R2 and SQL Server 2008 R2

- Improved waveform display: All ECG waveforms are now rendered with a uniform scale (initially, +/-2mV (1mV/cm). This removes the misperception that the ECG signal amplitude changed, when only the rendering scale recorded by the device changed. Furthermore, Event Review Pro can now display many segments as longer, contiguous segments.
- Support for Windows 7 64-bit as well as 32-bit.
- Automatic saving of all user settings, including ECG waveform order. You can restore defaults by individual UI section.
- Multiple file import/export: When you import data, you can choose multiple files. Each file is imported to a separate new case. You can also choose multiple cases and export them, and the case IDs will be used to keep their names unique.
- User-entered timeline annotations now use real time (as opposed to case time), which is more accurate.
- Customer support for device history data via email from FRx, HS1, FR2, and FR3 defibrillators. You can send the data file to Philips Customer Support for assistance.
- Q-CPR detail reports now can include information from multiple cases.

In addition, the User Account administration feature (including logon and logoff) and the ECG Views feature have been removed.

Any late-breaking news on features is in the Read Me document, available from the Help menu.

Getting assistance

Event Review Pro offers the following sources of help.

Application tooltips

When you hover the cursor over a toolbar button or event on the ECG waveform, text appears in a pop-up window. The text identifies the toolbar button and describes the event.

Wizards

A wizard guides you through common tasks.

The wizard that is available on the Getting Started navigation pane and File menu guides you through creating a case.

The IrDA wizard is available from the Attach ECG window in Cases. This wizard guides you through the process that downloads and attaches ECG information using the infrared transfer method.

Event Review Pro essentials

This document, in PDF format, is designed for new users. It focuses on the simplest and most common workflows: adding and documenting a case, reviewing a case, generating reports, and sending reports to other medical personnel. It includes only basic procedures to accomplish these tasks. To see this document, on the Help menu, click Event Review Pro Essentials.

Online help

To get help with the current feature, press F1.

To read an overview of a topic, on the Help menu, click **Help**. Select a topic from Help Contents. You can use the browse buttons to move from topic to topic.

User guide

The printable user guide, which includes all the information in the online Help system, is in Portable Document Format (PDF). You can view it using a PDF reader such as Adobe Acrobat.

You can view the user guide online, print sections, or print the whole guide.

For customer support, see [Customer support on page 132](#).

Conventions used in this guide

This document uses the following conventions to help identify information.

Conventions	Used for	Example
Bold	Menu options, buttons, field names, and list box names that you need to click.	On the File menu, click Import .
<i>Italic</i>	Manual names.	<i>Instructions for administrator's guide.</i>

Conventions	Used for	Example
Screen message	To distinguish defibrillator voice prompts and HeartStart FR3 screen messages from general text.	The defibrillator voice prompt announces Administration Press the center soft key (soft key #2). The screen message reads Ready for Data Transfer
<i>Italics with leading period</i>	File extension types.	The configuration is exported as an <i>.xml</i> file.
Lower-case <i>italics</i>	Specific file names.	Name the configuration file <i>sample.xml</i> .
<i>Italics</i>	Variables and placeholders	The default folder is located at C:\Document and Settings\ <i>user name</i> \My Documents.



A yellow box with a triangular caution icon and exclamation mark identifies circumstances that can result in data corruption or information loss.



A blue box with a note icon contains information on how features are used.



A green box with a light bulb tip icon contains information to complete a task.

Using online help

Use online help system to search for information using the full-text search feature, glossary, or an index.

You can access online help in these ways::

- Press F1 to open online help for the selected feature.
- On the Help menu, click **Event Review Pro Help**.

Starting the application

Event Review Pro does not place an icon on the desktop. To start Event Review Pro, use the shortcut on the Start menu or the All Programs menu.

To start Event Review Pro

1. Click the Windows **Start** button.
2. Click **All Programs**.
3. **Click Philips HeartStart Event Review Pro 4.2.**

You can also start Event Review Pro directly from a *.cod* or *.mic* file. Navigate to the folder containing the *.cod* or *.mic* file, and double-click the file.

To use Event Review Pro, click a navigation button on the left pane of the application window.

Understanding the application window

When you start Event Review Pro, the Getting Started pane is above the navigation buttons on the left. A welcome page is on the right. The Getting Started pane helps you to start working in the application.

Menu options and toolbar buttons change based on the navigation button and functions that you select.


- The menus list options to access the application features.
- The toolbar buttons provide quick access to frequently used menu options.

Beneath the menus and toolbar buttons, the navigation pane is on the left and the workspace is on the right.

- The navigation pane has navigation buttons and a pane above the buttons.
- A navigation button groups major application tasks by their function. You can hide or display navigation buttons.
- A navigation pane is above the navigation buttons. The pane lists navigation tasks.

The workspace on the right pane changes based on the navigation button and the task that you select. For example, if you click **Cases**, and then click the **All Records** navigation tool, all the cases stored in the database appear in the workspace. When you double-click a row on the workspace, the details of the case appear for viewing or editing.

The following illustration shows a typical Cases tab.

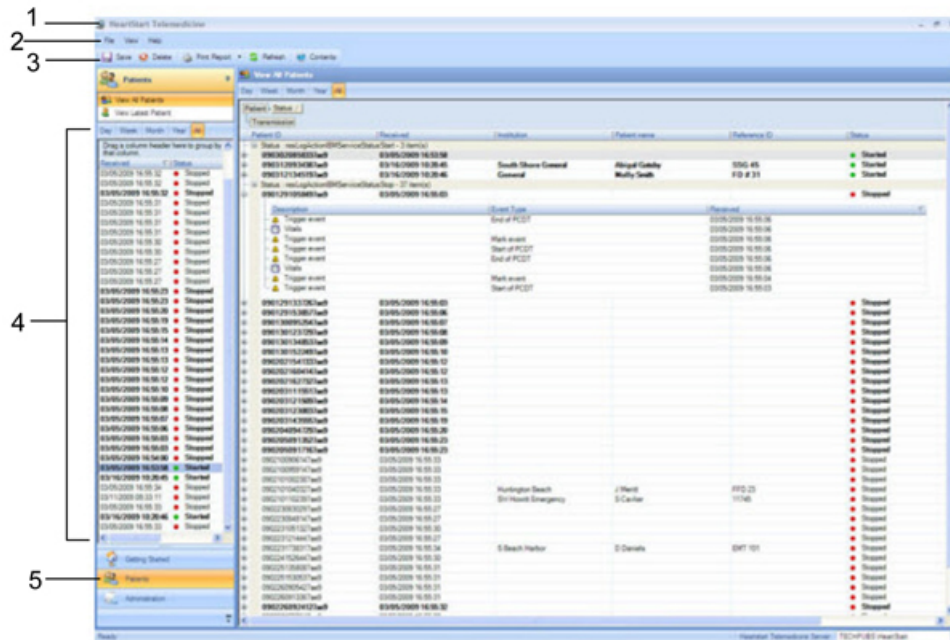

Key to illustration numbers

1	The title bar displays the product name and software version.
2	The toolbar provides quick access to frequently used actions that you can perform for the selected workspace, as represented by icons.
3	The menu bar displays the following menus: File, Edit, View, Tools, and Help. Use the menu bar to select the actions that you want to perform.
4	The navigation pane lists tasks that you can perform in the workspace. The pane content changes depending on which button is clicked and the workspace that is active.
5	Use the navigation buttons in the navigation pane to move among the workspaces.

Completing fields

Use a pointing device, such as a mouse or stylus, to select field values from a window and to change the way information on tables and logs appears.

In Cases, when you use the mouse to click or point to a field, a menu with options for that field might appear.

- If check boxes appear to the left of the options, you can select *multiple* options for the field.
- If check boxes do not appear to the left of the options, you can select *only one* option for the field.
- If an option is a text box, type or write the information using your keyboard or stylus.
- If you want to remove the selected option, click **Clear**.
- If a menu does not appear, type or write the information using the keyboard or stylus.

Using the application tables and logs

In Cases, Administration, and Defibrillators, information appears in tables and logs to help you manage your data.

Each row represents one record in the database. A record can be a case, an event, an attachment, or a system activity. Each column has values for a specific field.

You can modify your view of the information by grouping and sorting the rows. You can also print a table or log from the File menu or toolbar; click Print, then Table, and then the table or log.

All Records table

The All Records table (in the Cases workspace) lists each case in the Event Review Pro database. The cases show the reference ID, case ID, patient ID, case date and time, and patient name.

You can use the All Records table to do the following tasks:

- Sort and group cases. See [Sorting and grouping cases on page 43](#).

- Search for cases based on a variety of values. For more information, see [Filtering entries on page 122](#).
- Hide or display columns. See [Working with columns on page 121](#).
- Open and delete cases. For more information, see [Displaying case details on page 43](#) and [Deleting cases on page 52](#).
- View and print reports. For more information, see [Working with reports on page 106](#).

Once you open a case, you see a set of screen tabs with case details. Select a tab to view, modify, and delete the information. You can use the case tabs in any order based on your tasks within your organization.

Event Log

The Event Log is one of the tabs for a selected case. It lists all defibrillator events and notes written for that case.

Once you open a case and select the Event Log tab, you can see details for an event. You can sort, group, and filter the defibrillator events and notes that a responder or reviewer added to the case. You can also print the events for the open case.

For more information, see [Reviewing case events on page 100](#).

Attachments table

The Attachments table, in the Cases workspace, lists each file that is attached to the case. It identifies each attachment by the file name, file type, and date that the file was attached. Use the Attachments tab to manage the files that you attach to the open case.

Once you open a case and click the Attachments tab, you can add and delete an attached file. You can also sort, group, and filter the list of attached files.

For more information, see [Attaching \(adding\) files to a case on page 102](#).

System Log

The System Log table identifies all system activities based on the type, action, description, date and time, and the user name. The System Log table is on the Administration workspace.

For more information, see [Working with the System Log on page 120](#).

Resizing panes and workspaces

You can change the default layout of the application window. The revised layout persists the next time you use the navigation pane or workspace.

To resize a pane or workspace

1. Use the mouse to click and hold a border on the pane that you want to resize.
2. Drag the border in the appropriate direction to a new location.

The size of the adjacent panes adjusts.

Managing mouse behavior

If you work on a personal computer, you can change the mouse to behave like a stylus. This is useful when your mouse tends to change its location after you remove your hand from the mouse.

To change the mouse behavior, on the Tools menu, click **Options**, and then click **Click Mode**. (This used to be called Tablet mode.)

In Click mode, you must first point and click the field, and then click the value. In the standard mode, you can point or click the field plus its value at the same time.

Using a Tablet PC


You can also work in Event Review Pro using a Tablet PC. Use the stylus in the same way as you would use a mouse or pointing device. You can also use the stylus with the Tablet PC Input Panel to enter text.

To complete a text entry field in tablet mode

1. Use the stylus to tap a field.
The tablet displays a floating icon.
2. Tap the icon to display the Tablet PC Input Panel.
3. Click one of the following three modes to enter the information:
 - Writing pad
 - Character pad
 - On-screen keyboard

4. Type the information. For more information, click the **Help** icon on the input panel.
5. Click **Insert** or tap the **Enter** key on the input panel.

Saving your work

Philips recommends that you save your modified patient details periodically. To save your work, on the **File** menu or toolbar, click **Save** .

Restoring confirmation messages

A confirmation message does one of two things:

- Informs you that an action might result in a change to a configuration
- Warns you that an action may cause a permanent change.

By default, confirmation messages are enabled.

If you have disabled confirmation messages (by checking the Do Not Show This Message Again box when such a message appears), you can re-enable them at any point.

To restore confirmation messages, on the **Tools** menu, click **Options**, and then **Restore Confirmations**.

Using the Navigation pane

The navigation pane on the left of the Event Review Pro display provides access to the application tasks.

To display a set of tasks, click a navigation button.


Following is a description of the navigation buttons.

- **Cases** has tasks to work with defibrillator patient data. It includes tools to download patient data from a supported defibrillator. It can also import it from a file that was previously downloaded by another HeartStart product such as Data Messenger. Use Cases to record information such as patient demographics, initial treatment, and outcome. You can also add comments to clarify a defibrillator event or treatment. Responders and the medical staff can add notes on the case.

- **Reports** includes a list of all available reports, plus a toolbar for printing, navigation, view changing, and export.
- **Administration** helps the system administrator monitor and troubleshoot system activity.
- **Getting Started** helps you create a case quickly, import an existing patient case, and open a case already in your database.



You can hide the left navigation pane to increase the size of the workspace on the right pane. You can always return to the navigation pane while you work on the set of features.

To hide the navigation pane, at the top of the navigation, click the **Auto Hide**  icon.

The navigation pane disappears, and the Navigation Panel tab along the left edge of the pane appears.

To display the navigation pane temporarily

1. Point to the tab along the left edge of the pane.
The navigation pane appears temporarily.
2. Point to or tap the right pane.
The navigation pane disappears.
3. Point to the tab again to redisplay the navigation pane.

To restore the navigation pane, click the tab, and then click the **Stay Open**  icon.

Working with cases

The word "case" refers to all the information relating to a specific cardiac or monitoring event. This includes information recorded by supported defibrillators and information entered by a user. For information about supported defibrillators, see [Supported defibrillators on page 137](#).

The first step in working with ECG information (ECG) is to create a case. You can do this in the following ways:

- Use the case wizard to download data from a defibrillator or its data card
- Import a case with an ECG from a *.mic*, *.hic*, or *.cod* file created by another application, such as Data Messenger

.Mic files are created by the export function of products such as HeartStart Data Messenger and Event Review Pro.

.Hic files are created by third-party ePCR software using the Philips Software Development Kit (SDK).

.Cod files are created by the export function of older products such as Review Express Connect and Event Review Pro 3.5.

- Create a case that does not have ECG information in it yet and add it later

Starting the cases feature

To create, review, and manage cases, click the **Cases** navigation button.

On the Navigation Pane, the database records appear, and on workspace, the All Records table appears.

For more information, see [Displaying case details on page 43](#).

Understanding the Cases navigation pane

Use the Cases navigation pane to work with ECG information. You can complete the following tasks in Cases:

- Download information from a supported defibrillator or its data card

- Import information from a patient data file created by another HeartStart product, such as Data Messenger
- Add information such as patient demographics, initial treatment, and notes to clarify or comment on a defibrillator event

If you select information from a drop-down list, your choice persists the next time you use the field. At that time, you can select another value.

Use the Database Records pane to access cases in the database.

Creating a case (including ECG) with the Case Wizard

The wizard that is available on the Getting Started pane and the File menu guides you through the steps required to add/download patient data from a defibrillator or its data card, create a case, and complete some basic tasks.

After briefly displaying a message listing the completed tasks, the wizard displays the page for the next task.

To start the wizard

1. On the navigation pane, click the **Getting Started** navigation button.
2. On the **Getting Started** navigation pane, click **Case Wizard**.
3. In the **Welcome Page** window, click **Next** to continue.

The Attach ECG page appears.

Attaching (adding) an ECG

To attach an ECG from a Philips defibrillator or its data card

1. Use the table in [Adding ECGs table on page 35](#) to identify the method used to add/download an ECG from your defibrillator.
2. Follow the instructions in [Adding ECGs on page 34](#) to attach the ECG.
3. Click **Next**.

The Case Summary page appears.

Entering case summary information

On the Case Summary page, you can handle case information three ways:

- Add new information
- Accept the information that the defibrillator provides
- Override the information from the defibrillator

All information is saved to the database.



Different defibrillators may capture different patient information as part of recorded ECG data. If you use data from the HeartStart MRx, the ECG information might include patient information such as name, age, and gender.

To complete the Case Summary fields

1. Event Review Pro completes the Case ID field. You can change the case ID.
For defibrillators that generate an incident ID, such as HeartStart MRx, the case ID corresponds to this incident ID. Each device has its own system for generating an ID.
2. Fill in the **Patient ID**, **Last Name**, and **First Name** fields.
3. In **Reference ID**, type a reference ID for the case that you can use to identify it later. For example, "ER Room 3" or "Engine #3, Station 27". This information may already appear if the case was imported from an application that provided it.
4. In **Institution**, type an institution name for the case. For example, "Valley Hospital" or "District 5 Fire Department".
5. Click **Next**.
The Export Case page appears.

Exporting a file

You can export a case to a *.mic* file at any point. The file is then available to a colleague who uses or Event Review Pro. To create a file containing the ECG and case information, check the **Export to a File** box and verify the file name and folder where you want to save it.

If you want to password-protect the exported file, in the **Password** box, type the password. Passwords can have no more than 16 characters. They are case-sensitive.

Remember to provide this password to anyone who needs to open the exported file.

If you need to remove patient identity, check the **Remove Patient Identity** box.

Click **Next**.

Printing a report

You can print a case report at any point. If you want to print it as part of the Case Wizard, check the **Print a Report** box, click the buttons for any available reports, and click **Next**. The Email page opens.

Emailing case data

You can email a case at any point; the case is added to an email as a *.mic* file. If you want to email it as part of the Case Wizard, check the **Email Case** box, type or click the appropriate address or addresses, change the subject if you want, and type a note to accompany the case if you want. You can also password-protect the file or remove patient identify information.

Click **Next**. The Execute page opens.

Opening the case

At this point, you can specify that the case will open after all tasks are completed, so that you can make further changes and review them. Check the **Open Case** box and click **Finish**. The status of each task within the Case Wizard displays in the box.

If your email system needs verification before it can send the email, a confirmation message appears. Click **Allow**, and the process continues.,

When the process is complete, click **Exit**. The case opens.

Creating a case

You can create a case manually or automatically, with or without ECG information.

To create a case that has ECG information manually

1. Use one of the following methods:
 - Click the **Getting Started** navigation button, and then click **Case Wizard**.
 - On the **File** menu, click **Case Wizard**.

The Case Wizard window opens. For more information, see [Creating a case \(including ECG\) with the Case Wizard on page 31](#).

2. Follow the on-screen instructions.

To create a case that does not initially contain ECG information manually

1. On the **File** menu or toolbar, click **New**.

The tabs for the case in the workspace appear.

2. Complete fields to briefly describe or identify the case. For example, add the case date and time, or select the site of collapse.
3. Click **Save**.

At a later time, you can open the case, attach ECG information, and complete the tabs.

To create a case automatically, see [Importing case files automatically on page 46](#).

Adding ECGs

"Attaching" an ECG means reading ECG information from a defibrillator (downloading it), a defibrillator data card, or a file and adding it to a case.

You can attach an ECG to a case you are creating, or you can open an existing case and attach an ECG to it. You can add extra ECGs from other defibrillators as well. Depending on the defibrillator, you can attach an ECG information and 12-lead information from a data card or a file, or through a *Bluetooth* or an infrared transmission.

If you need help identifying the defibrillator or inserting a data card, see [Working with defibrillators on page 137](#).

The defibrillator that you use determines the method that you can use to attach an ECG. The Attach ECG window lists the types of transfer methods that you can use to attach an ECG or 12-lead report on the left pane. The window lists data files in the right pane. The data files include those that were transmitted or are stored on a removable device.

To display the Attach ECG window for an existing case

1. Click the **Cases** navigation button.
2. On the **All Records** table, double-click a case.
3. On the toolbar, click **Attach ECG File**.






Any time you have a case open for reviewing or editing, you can attach an ECG. On the ECG menu, click **Attach**, or on the toolbar, click **Attach ECG File**.

Choosing a transfer method

Event Review Pro activates the following methods automatically. (If you are transferring data from the MRx, you might have to click the *Bluetooth* icon the first time that you use the feature.)

You can also click the following icons to attach ECG data.


-  *Bluetooth*—begins the process of gaining access to ECG data from MRx and FR3 devices. When you click this icon, a wizard appears with instructions for establishing both MRx and FR3 connections. If you click MRx, you then navigate to the MRx Bluetooth Exchange Folder that contains (or will contain) the transmitted *Bluetooth* files. If you click FR3, a wizard appears with instructions for how to put the FR3 into Administration mode so that you can see and download all available ECG files. (You might have to click the *Bluetooth* icon the first time that you use the feature.)
-  IrDA—opens help on how to transmit data from HS1 and FRx defibrillators.
-  File System—displays the Open window that you can use to browse and attach ECG information from a device's data card or your computer.

Adding ECGs table

Use the following table to identify the method and instructions used to read and/or download an ECG from a defibrillator.

Defibrillator model	Method used to read data	Instructions
FR2 Series	Card reader	See Adding an ECG from a data card on page 36 .
FR3 Series	Card reader <i>Bluetooth</i> connection	See Adding an ECG from a data card on page 36 . See Downloading an ECG using an FR3 Bluetooth transmission on page 37 .
XL	Card reader	See Adding an ECG from a data card on page 36 .
XL+	USB data drive	See Transferring XL+ data on page 146 ;
HeartStart MRx	Card reader	See Adding an ECG from a data card on page 36 .

Defibrillator model	Method used to read data	Instructions
HeartStart MRx	HeartStart MRx, version 9.0 or later <i>Bluetooth</i> connection	See Downloading an ECG using an MRx Bluetooth transmission on page 38 .
HeartStart FRx	Infrared connection	See Downloading an ECG from an infrared connection on page 39 .
HS1 family of defibrillators, including HeartStart Home and HeartStart OnSite Defibrillators	Infrared connection	See Adding an ECG from a file on page 40 .
Imported case file	File system	See Adding an ECG from a file on page 40 .

 Event Review Pro documentation uses the HS1 or HS1 family of defibrillators to refer to HeartStart Home, HeartStart OnSite, and HeartStart HS1 Defibrillators.

If you need help identifying the defibrillator or using a transfer method, see [Working with defibrillators on page 137](#).

Adding an ECG from a data card

The FR2, FR3, XL, and HeartStart MRx defibrillators can use a data card to transfer the recorded ECG and other data from the defibrillator to Event Review Pro.

To attach an ECG from a data card


1. Insert the data card into the card reader.
2. Use one of the following methods:
 - On the ECG menu, click **Attach**.
 - On the toolbar, click **Attach ECG File**.

The Attach ECG window opens, listing ECG files grouped by the source of the file.

3. On the right pane of the Attach ECG window, find the file that is associated with the ECG information or 12-lead report. The medium for the data source appears under the Source column.

4. (Optional for FR2 Series defibrillators) If the pediatric pads were used, choose the FR2 Pediatric option. This allows Event Review Pro to properly calculate the pads impedance and the shock energy desired. If you have any doubt about which pads the responder used, contact the responder.

Click the ECG you want to attach. In the Device column, right-click the **FR2**

icon and click **FR2 Pediatric**. The pink teddy bear  icon for FR2 Pediatric appears.

5. Use one of the following methods:
 - Click the file for the ECG you want to attach, and then click **Open**.
 - Double-click the file for the ECG you want to attach.

Event Review Pro downloads the file and attaches the ECG to the case. It appears on the ECG tab.

Adding an ECG from a USB port

The HeartStart XL+ defibrillator uses a USB data drive to transfer the recorded ECG from the defibrillator to Event Review Pro.

To import data from the XL+ USB data drive to Event Review Pro

1. Insert the XL+ USB data drive into the USB data port on the computer running Event Review Pro.
2. Ensure that Windows recognizes the XL+ USB data drive as a connected device, and that Event Review Pro is running.

Windows detects the USB data drive and begins transferring data. All XL+ data is transferred to Event Review Pro.

3. When data transfer is completed, follow Windows procedures to safely remove the XL+ USB data drive from your computer.

Downloading an ECG using an FR3 *Bluetooth* transmission

You can transfer files from the FR3 to the computer using *Bluetooth* transmission, if the FR3 is equipped with *Bluetooth* capability.

To attach an ECG using a *Bluetooth* transmission

1. Set up the *Bluetooth* software and adapter on the computer. For information, see your defibrillator documentation and [Setting up Bluetooth transmissions for the](#)

[MRx and FR3 on page 141.](#)

2. Prepare to send the *Bluetooth* transmission by placing the FR3 (turned on, after the voice acknowledgment) into Administration mode, pressing the option button, and selecting **Wireless Data Transfer**.
3. In Event Review Pro, use one of the following methods:
 - On the ECG menu, click **Attach**.
 - On the toolbar, click **Attach ECG File**.
4. On the right pane of the Attach ECG window, find the ECG recording that you want to attach.

The medium for the data source (in this case, *Bluetooth*) appears under the Source column. Keep in mind that if your ECG recording is very long or contains audio, it may take as long as 5 to 10 minutes (or even longer) before it appears in the window.

5. Use one of the following methods:
 - Click the file for the ECG you want to attach, and then click **Open**.
 - Double-click the file for the ECG you want to attach.

Event Review Pro attaches the ECG to the case, and the ECG tab appears with the attached ECG waveform(s).

Downloading an ECG using an MRx *Bluetooth* transmission

Depending on the version of your MRx, and its options, you can transfer files from the MRx to the computer using *Bluetooth* transmission.

To attach an ECG or 12-lead report using a *Bluetooth* transmission

1. Set up the *Bluetooth* software and adapter on the computer. For information, see your defibrillator documentation and [Setting up Bluetooth transmissions for the MRx and FR3 on page 141.](#)
2. Send the *Bluetooth* transmission. For information, see [Sending the MRx Bluetooth transmission on page 144.](#)
3. Use one of the following methods:
 - On the ECG menu, click **Attach**.
 - On the toolbar, click **Attach ECG File**.
4. On the right pane of the Attach ECG window, find the file that is associated with the ECG or 12-lead report.

The medium for the data source appears under the Source column.

5. Use one of the following methods:

- Click the file for the ECG you want to attach, and then click **Open**.
- Double-click the file for the ECG you want to attach.

Event Review Pro attaches the ECG to the case.

6. If the file does not appear in the right pane, click **Bluetooth** on the left pane of the Attach ECG window and use the Bluetooth Wizard to navigate to the Windows *Bluetooth* Exchange Folder.
 - On the **Browse for Folder** window, navigate to the **Windows Bluetooth Exchange Folder** for the user. The default folder for Windows XP is located at C:\Document and Settings*user name*\My Documents\Bluetooth Exchange Folder. The default folder for Windows 7 is located at C:\Users*user name*\My Documents.
 - If you specified a different folder when setting up your *Bluetooth* stack or Neighborhood, navigate to that folder.
 - Event Review Pro scans the location you specify to find the ECG. It enables the *Bluetooth* icon to the right of the status bar while it reads the *Bluetooth* files.
7. One or more ECG tabs appear with the attached ECG waveform(s).



If the ECG tab does not display an ECG waveform, contact Customer Support. See [Customer support on page 132](#).

Downloading an ECG from an infrared connection

The HeartStart HS1 family of defibrillators and the FRx defibrillator use an IrDA connection to attach an ECG to Event Review Pro.

To attach an ECG from an infrared connection


To retain the actual event times, do not remove the battery from HS1 or FRx defibrillators. If you remove the battery before you download the case information, the files created before you removed the battery might appear in



the list without a date and time (though they will show elapsed time for events and ECG data when you open them). For information on HSI and FRx defibrillators' case date and time, see [Determining the HS1 and FRx case date and time on page 147](#).

1. Set up the infrared connection between the defibrillator and the computer. For more information, see [Using infrared connections for the HS1 and FRx on page 148](#).
2. Use one of the following methods:
 - On the ECG menu, click **Attach**.

- On the toolbar, click **Attach ECG File**.
3. Put the defibrillator in Administration mode. For more information, see [Setting up the infrared connection on page 148](#).
 4. (Optional) If you need help with setting up the IrDA handshake for HS1 or FRx defibrillators, use the IrDA Wizard.

On the left pane of the Attach ECG window, click IrDA . On the IrDA Wizard window, follow the on-screen instructions. The IrDA icon to the right of the status bar appears enabled while it transfers the data.


5. When the defibrillator connects, the voice message says Sending (for the HS1) or Transferring data (for the FRx).
6. On the right pane of the Attach ECG window, find the file that contains the ECG. The data source appears in the **Source** column.
7. Use one of the following methods:
 - Click the file for the ECG you want to attach, and then click **Open**.
 - Double-click the file for the ECG you want to attach.
 Event Review Pro downloads the file and attaches the ECG to the case.
8. A message appears requesting the case start time when the HS1 or FRx case does not have a start time. Event Review Pro requests the start time.

Adding an ECG from a file

When you attach an ECG or 12-lead information from a file, the file can be on the hard drive, external media, or the network.

This might occur when a remote user downloads a case via another HeartStart application, like Data Messenger, and saves it to a file for your subsequent use.

To attach an ECG from a file

1. Use one of the following methods:
 - On the ECG menu, click **Attach ECG**.
 - On the toolbar, click **Attach ECG File**.
2. On the left pane of the Attach ECG window, click **File System** . The Open window lists file extensions in the Files of Type field.
3. In the **Files of Type** field, click the file type.
4. Click **Device Data** for files with the *.s01*, *.index*, *.cod*, *.xml*, or *.mic* extensions—that is, all files except those that include MRx 12-lead reports.

- Click **FR2 Series Peds** for files that were recorded using FR2 defibrillator pediatric pads.
 - Click **MRx 12- lead (XML)** for 12-lead report files that were recorded with the MRx.
5. Navigate to the file you want and double-click it.

Event Review Pro downloads the file and attaches the ECG to the case.



If the ECG tab does not display an ECG waveform, contact Customer Support. See [Customer support on page 132](#).

After you add an ECG

After you attach the ECG information to a case, use the tabs on the Cases workspace to add details to the case.

The following table lists common tasks and related instructions.

Task	Topic
Add details to a case	Adding case details on page 54
Save a case	Saving cases on page 46
Print reports	Printing case information on page 49
Export a case	Exporting cases on page 49
Email a case	Emailing a case on page 51

Adding multiple ECGs

If the emergency response for a single patient involved multiple ECGs on supported defibrillators, you can attach all ECG files to the same case. On the ECG tab, a tab appears for each attached ECG. The tab text identifies the defibrillator, date, and time transmitted by the defibrillator. For example: HeartStart MRx 3/10/2007 3:10:44 PM.

To attach the first ECG, follow the instructions in [Adding ECGs on page 34](#).

To attach additional ECGs

1. Make sure that the case is open. For more information, see [Displaying case details on page 43](#).


2. Use one of the following methods:
 - If you are attaching an ECG from an HS1 or FRx defibrillator, make sure that the defibrillator is set up to communicate using an infrared connection.
 - If you are attaching an ECG from a data card, insert the card into the card reader.
 - If you are attaching an ECG using a *Bluetooth* connection, make sure that the defibrillator is set up for transmission.
3. Use one of the following methods:
 - On the ECG menu, click **Attach**.
 - On the toolbar, click **Attach ECG File**.
4. On the right pane of the Attach ECG window, find the file that generated the ECG.

The medium for the data source appears in the Source column.
5. Use one of the following methods:
 - Click the file for the ECG you want to attach, and then click **Open**.
 - Double-click the file for the ECG you want to attach.

The ECG appears on a separate tab on the ECG tab.

Removing ECGs from a case

Use the Detach option, on the ECG menu and toolbar, to remove an ECG from the case.

 The ECG menu is available when you open a case and work on the ECG table.

To detach an ECG


1. Open the case in the Cases workspace.

For more information, see [Displaying case details on page 43](#).
2. Click the **ECG** tab.
3. On the ECG tab, click the **ECG ID** tab for the ECG you want to detach.
4. Use one of the following methods:
 - On the ECG menu, click **Detach**.
 - On the toolbar, click **Detach ECG**.

A confirmation message might appear. For more information, see [Restoring confirmation messages on page 28](#).
5. Click **Yes** to remove the ECG from the case.

Displaying case details

Use the table on the Cases workspace to display patient identifiers and case details. The table lists each case in the database. The following appear for each case: the reference ID, case ID, patient ID, date and time, patient name, and institution.

 You can customize how the information appears. For more information, see [Working with columns on page 121](#) and [Grouping and sorting entries on page 121](#).

Once you open a case, you select a tab to view and modify the detailed information. You can also view and print reports for the open case and export the open case to a different location as a file.

To open a saved case

1. Click the **Cases** navigation button.
2. Use one of the following methods:
 - On the **All Records** table, double-click the case you want to open.
 - On the **File** menu or toolbar, click **Open**.

The case appears in the Cases workspace.

Sorting and grouping cases

You can change the way information displays on the All Records table. You can sort cases and group them based on a shared parameter value.

You can filter or select cases that meet the criteria you select. You can also hide or display groups of cases.

For information on how to add columns to and delete columns from the display, see [Working with columns on page 121](#).

For information on how to group and sort information in columns, see [Grouping and sorting entries on page 121](#).

For information on how to filter the information that appears in the list of cases, see [Filtering entries on page 122](#).

Hiding and displaying cases

You can hide or display the list of cases in the groups. You can also view and refresh the list of cases from the View menu or shortcut menu on the All Records table. The shortcut menu lists the Refresh, Expand All, and Collapse All options.

To use the shortcut menu

1. Open **Patient Cases** and then **All Records**. For more information, see [Displaying case details on page 43](#).
2. Right-click the **All Records** table.
3. Click an option.

To refresh, hide, or display cases

- To refresh the list of cases, on the **View** menu or shortcut menu, click **Refresh**.
- To hide cases in a group, click **Collapse All** to the left of the grouping.
- To hide cases in each group, on the **View** menu or shortcut menu, click **Collapse All**.
- To display cases in a group, click **Expand All** to the left of the group.
- To display all cases in the table, on the **View** menu or shortcut menu, click **Expand All**.

Importing case files manually

You can import one or more case files that were exported by another HeartStart software application such as Data Messenger.

For each file you import, Event Review Pro creates a case in the Event Review Pro database, adds the case to the All Records table, and adds the information from the file to the case. You can use the file in the same way as any other case you created. You can save, modify, print, export, and email an imported file.

You can import using the Import option (available on the File menu and toolbar) or using drag-and-drop methods.



When you import a case file, Event Review Pro imports all the details that were in the case, such as patient demographics, responder and reviewer notes, and physician notes. When you use the Attach ECG feature, Event Review Pro attaches only the ECG to the case.

To import files using drag and drop

1. On the navigation pane, click the **Patient Cases** navigation button.
2. Open the folder where your *.mic.*, *.hic*, or *.cod* files are stored, and drag one or more into the All Records table.
3. Event Review Pro imports the cases, saves them to the database, and lists them on the **All Records** table.

If any of the cases failed to import successfully (because they required a password, for example), the outcome for each appears in the System Log. See [Working with the System Log on page 120](#).

To import files using the Import option

1. On the navigation pane, click the **Patient Cases** navigation button.
2. On the **File** menu or toolbar, click **Import**.
The Import Case window opens.
3. Navigate to the folder that contains the files you want to import.

The list of files includes only those with a *.cod*, *.hic*, and *.mic* extension.



A *.hic* file is a defibrillator device file that is forwarded by an electronic patient care record (ePCR) application or HeartStart Data Messenger.

4. If you want to import only one file, click it. If you want to import multiple files, press CTRL and click each one. Click **Open**.
5. If any of the files were password-protected, type the file password in the Password field before you attempt to open the file. Passwords can have up to 16 characters and are case sensitive. If some of the files are associated with a different password, they will not import, and you will need to import them separately with their own password.
6. Event Review Pro imports the cases, saves them to the database, and lists them on the **All Records** table.

If you were working with an open case when you imported more than one case, the last case in the list that you imported opens. If you were working in the All Records table, the same table appears after the import.

If any of the cases failed to import successfully (because they had a different password from the others, for example), the outcome for each appears in the System Log. See [Working with the System Log on page 120](#).

Importing case files automatically

Use the Import Service on the Administration navigation pane to set up inboxes to import cases automatically. When the Import Service is running, Event Review Pro monitors the inboxes for case files with *.hic* , *.mic* , and *.cod* extensions.

If an imported case contains a duplicate ECG, the case appears on the **Duplicated ECGs** table. See [Reviewing cases with duplicate ECGs on page 52](#).

For more information, see [Using the Import Service on page 127](#).

Saving cases

When you save a case, Event Review Pro adds it to the database. All the information is available for reports. Use the All Records table to manage the case.

To save an open case, on the **File** menu or toolbar, click **Save**. The case is now in the Event Review Pro database.

Adjusting the date and time of the defibrillator data

It is possible for the internal clock on your defibrillator to drift slightly over time. In addition, if someone removed the battery from an HS1 or FRx defibrillator before downloading the patient data, case and event real time could be lost, leaving only event elapsed time. In these cases, you may want to adjust the times of the downloaded case and its individual events.

Before you transfer information from a defibrillator, confirm that your computer clock is accurate. For example, if the computer does not connect to a time server, you can compare the date and time with an atomic clock on the Internet. Also, confirm that the battery was not pulled from the HS1 or FRx defibrillator since the cardiac emergency. For information on HS1 and FRx defibrillators case date and time, see [Determining the HS1 and FRx case date and time on page 147](#).

To adjust the defibrillator date and time based on the computer date and time

1. Open the case in the Cases workspace. For more information, see [Displaying case details on page 43](#).
2. On the **Overview** tab, navigate to the Defibrillators area.

Confirm that the date and time for the defibrillator are correct in the Defibrillators area.



If you transfer the case from a defibrillator through a *Bluetooth* transmission, Event Review Pro automatically adjusts the case event date on the received data based on the computer clock.

3. After you attach an ECG, you can change the date and time for the appropriate case in the **Adjusted “On” Time** column.

The date and time automatically adjusts for each event in that case. The **ECG** tab displays the adjustment that you enter.

4. On the **File** menu or toolbar, click **Save**.

Erasing the data source

After you attach an ECG and save the case to the database, consider erasing the data from the source. This ensures that the device has full capacity for the next use.

You can erase ECG information from the data cards used by the MRx, the FR2 and FR3 series, and XL defibrillators. You can also erase the ECG information from an HS1, XL+, or FRx defibrillator.

You can also use *Bluetooth* to erase files that were transmitted via *Bluetooth* from an MRx.



Make sure that you transfer all the information to the database or an archive before erasing the information. Once the data is erased, it is lost and cannot be recovered.

To erase a data card

1. Insert the data card into the card reader.
2. On the **Tools** menu, click **Erase Data on Card or Device**.

The Erase Data on Card or Device window opens.

3. Click the card or device name.
4. Click **Erase**.

A confirmation message might appear. For more information, see [Restoring confirmation messages on page 28](#).

5. Click **Yes** in the message box.

To erase data from an HS1 or FRx Defibrillator

1. Set up the defibrillator to communicate with Event Review Pro. For more information, see [Using infrared connections for the HS1 and FRx on page 148](#).
2. Put the defibrillator in the Administration mode.
3. On the **Tools** menu, click **Erase Data on Card or Device**.
The Erase Data on Card or Device window opens.
4. Click the device name.
5. Click **Erase**.
A confirmation message might appear.
6. Click **Yes** in the message box.

To erase data from the Bluetooth Exchange folder (MRx only)

1. On the **Tools** menu, click **Erase Data on Card or Device**.
The Erase Data on Card or Device window opens.
2. **Click Bluetooth Exchange Folder**.
3. Click **Erase**.
A confirmation message might appear.
4. Click **Yes**.

To erase data from an FR3 defibrillator using Bluetooth

1. Once the unit is on and you have heard the voice prompt, place the FR3 in Administration mode, press the option button, and select **Wireless Data Transfer**.
2. In Event Review Pro, on the **Tools** menu, click **Erase Data on Card or Device**.
The Erase Data on Card or Device window opens.
3. In the *Bluetooth* section, click the FR3 icon whose serial number matches the device that you want to erase.
4. Click **Erase**.
A confirmation message might appear.
5. Click **Yes**.
The clock on the defibrillator is now synched with the clock on the computer. This synchronization occurs even if you have no data to erase.

To erase data from a USB data drive

1. Use the XL+ defibrillator to erase data from the USB data drive. For more information, see the HeartStart XL+ instructions for use.

Printing case information

You can print case information about an open case and about a case that is stored in the database. For more information, see [Working with reports on page 106](#).

You can also print the tables that appear on the Case Records table, Event Log tab, and Attachments tab. The information that prints is about the open case.

To print information listed in a table

1. Open the case and display the Case Records table. For more information, see [Displaying case details on page 43](#).
2. (Optional) You can sort, group, or filter the list of cases. For more information, see [Sorting and grouping cases on page 43](#).
3. On the **File** menu or toolbar, click **Print**.
4. On the **Print** menu, click **Tables**.
5. Click the table you want to print.
6. On the Print window, click the printer, page range, and number of copies.
7. Click **Print**.

Exporting cases

Use Export to create files outside the application (that is, not in the application database) so that you can share information with other Event Review or Event Review Pro users. You can also use Export to back up information.

You can export a single case or multiple cases at a time. The Export option is available on the File menu whenever the Cases list appears or a specific case is open. In addition, if you want to export just one case, you can do this (from the File menu or the toolbar) when that case is open.

You can export files with or without a password. If you add a password, anyone attempting to open the file must have the password to open it.



Philips recommends that you record the password and save it in a secure location. If you forget the password, Customer Support cannot “unlock” the file.

To export one or more case files from the Cases list

1. On the **Navigation** pane, click the **Patient Cases** navigation button.
All cases appear.
2. Choose the case or cases to export:
 - a. Click a single case.
 - b. CTRL-click multiple cases.
3. On the **File** menu, click **Export**.
The Export Selected Cases dialog box appears, showing the last folder you used to export a file. If you want a different location, click the **Save in** field at the top of the dialog box and navigate to that location. The list there shows *.mic* files, which is the only supported file format.
4. Decide how to name the exported file or files. The default name includes a placeholder for the case ID, expressed as {caseID}. You can accept the case ID as the filename, or you can use this placeholder to create more descriptive names; for example, if you type “my data {caseID} for Joe.mic” for the filename and export two cases, one with a case ID of “abc” and one with “xyz,” then the two *.mic* files will have the names “my data abc for Joe.mic” and “my data xyz for Joe.mic.” If you type something entirely different for the name that doesn’t include a case ID, the case ID will be appended to the filename by the software.
5. If you want to remove confidential patient identifiers such as name, patient ID, and age (if over 90) from the exported file (called “redaction”), check the **Remove Patient Identity** box.
6. Decide whether you want to protect the file or files with a password. If so, type it in the **Password** field (for security, it appears as *****). Remember to provide the password for the recipient.
7. Click **Save**.
A progress bar appears as the export operation proceeds, and a confirmation message tells you that the export operation succeeded.

To export a single case file from the case display

1. On the navigation pane, click the **Patient Cases** navigation button.
2. In the list of cases, click a case to open it.
3. On the **File** menu or toolbar, click **Export**.
4. Follow steps 3 through 7 in the procedure above.

Emailing a case

You can use the Email option or your email application to send an exported ECG or an entire case to another Event Review Pro user.

Event Review Pro attaches the exported case file to the email when it displays the email window.

Before you use the Email option, complete the following prerequisites:

- Correctly configure Microsoft Outlook or a MAPI-compliant email client.
- Set up the email profile in Windows, if needed. Click Help on the email profile window for information on setting up your profile. After you complete the profile, the profile does not appear again.

To use the Email option

1. On the **Cases** workspace, display the case details. For more information, see [Displaying case details on page 43](#).

2. On the **File** menu or toolbar, click **Email**.

The Email window opens and attaches the file to the email form.

3. Complete the **To** field. The contact list for your email application is not available for the Email tool.

Email saves the address. The next time you send an email to that recipient, you can select the same address by clicking the arrow at the end of the To field.

4. If you want, add a message.
5. Click **Send**.


To use your email application

1. Export the file. For more information, see [Exporting cases on page 49](#).
2. Note the location of the exported file. The file has the supported extension.
3. Open your email application and create an email to the recipient of your choice.
4. Attach the file created by Export in Step 1 to the email form.
5. Send the email.

Deleting cases

You can delete a case or multiple cases from the database. The case or cases then disappear from the **All Records** table.


To delete a case or cases from the list of cases

1. Click the **Patient Cases** navigation button to start the Cases feature.
2. On the **All Records** table, click the case you want to delete, or press CTRL or SHIFT while clicking to select multiple cases..
3. On the **File** menu or toolbar, click **Delete** .

A confirmation message might appear. For more information, see [Restoring confirmation messages on page 28](#).

4. Click **Yes**.

To delete a case that is open

1. On the **File** menu or toolbar, click **Delete** .

A confirmation message might appear. For more information, see [Restoring confirmation messages on page 28](#).

2. Click **Yes**.

The All Records table appears.

Reviewing cases with duplicate ECGs

The Duplicate ECGs table (in the Cases workspace) lists each case in the Event Review Pro database for which the ECG in the case is also the ECG in another case. The cases are listed based on the reference ID, case ID, patient ID, case date and time, patient name, and institution.

Cases can be inadvertently duplicated when you download them automatically. This table provides an easy way to identify these cases. For more information, see [Using the Import Service on page 127](#),

Use the Duplicate ECGs table to complete the following tasks:

- Delete duplicate cases. For more information, see [Deleting cases on page 52](#).
- Sort cases. See [Sorting and grouping cases on page 43](#).

- Search for cases based on a variety of values. For more information, see [Filtering entries on page 122](#).
- Hide or display columns, See [Working with columns on page 121](#)
- Open cases. For more information, see [Displaying case details on page 43](#). You can click the Cases tab to view, modify, and delete the information. You can use the case tabs in any order based on your tasks within your organization.
- View and print reports. For more information, see [Generating reports on page 107](#).

Adding case details


Use the Patient Cases navigation button to access a case and enter case information that the defibrillator did not record. In Cases, you can do the following tasks:

- Document events
- Document patient outcomes
- Add responder actions and observations
- Add responder notes
- Add reviewer notes
- Document patient history
- Review ECG information
- Review events
- Attach files to a case

For more information, see [Working with ECGs on page 65](#) and [Reviewing case details on page 100](#). You can use the case tabs in any order based on your tasks within your organization.

To add and review case details

1. Click the **Patient Cases** navigation button.
2. On the Cases workspace, double-click a case on the **All Records** table.
3. Click a tab.
4. Add or change information on the tab.
5. Repeat steps 3 and 4 as necessary.
6. On the **File** menu or toolbar, click **Save**.

 Event Review Pro displays case detail tabs based on your installation.

Identifying the case

Use the Overview tab to enter basic information about the patient and the response location. The only required information is the case ID, and the case date and time.

If you use the Case Wizard to create a case, Event Review Pro assigns a case ID, which you can change. It also generates the case ID from the current date and time if Case Wizard or the Import feature did not create the case.



Different defibrillators may capture different patient information as part of recorded ECG data. If you use data from the HeartStart MRx, the ECG information might include patient information (such as name, age, and gender) that was entered by the responder. When you import a case using Case Wizard, the Case Summary page displays the information provided by the MRx, and the Case ID field displays the MRx incident ID.

If you export a redacted case (that is, a case with patient identification removed), Event Review Pro replaces the name and patient ID with a series of asterisks. The Date of Birth field is blank. If the redacted age of the patient is greater than 90, the Age field displays 90. If you then import the case, the redacted information (the information that was removed) does not appear.

After you attach the ECG, you can change any of the information except the information in the following fields, which the defibrillator may supply:

- Recorded “On” Time
- Defibrillator Type
- Serial Number
- Shocks

To complete the fields

1. Click the field.
2. Click an option on the menu or type the information in the field.




If Other appears in the option list, you can type a value in the Other text box. If your entry is the same as an option in the list, your entry appears as shown in the option list. Click **Clear** to remove the information in the field or text box.

Overview section fields table

Group or field	Description
Case ID (Required)	When you attach an ECG, incident ID becomes the case ID. You can change the case ID.
Time zone	Initially, displays the time zone for the computer clock. The time zone displays in the Coordinated Universal Time (UCT) or Greenwich Mean Time (GMT) format in a similar format to the Time Zone tab on the Windows Date and Time Properties window. For example: (+1:00) Amsterdam, Berlin, Bern, Rome, Stockholm, Vienna. Select the time zone of the location where the incident occurred.
Case date and time (Required)	.If you attach an ECG before you save the case, this field is set to the date and time when the device was turned on. To change the date and time, click the value you want to change, and then type a new value.
Site of collapse	The list includes typical kinds of locations such as Ambulatory or outpatient area, Hotel, Workplace, Adult intensive care unit (ICU), and Cardiac catheterization lab (cath lab).
Location Detail	You can use this field to be more specific about the site of collapse. For example, consider using this field to record map coordinates or the name of locations, such as the room number.
Reference ID	You can use this field to help identify the case. For example, Engine 3.
Institution	You can use this field to identify the organization where the case originated. For example, Valley Fire Department.
Code called	The list is: Yes, No, and Unknown
Found by	The list is: Nurse, Physician, Therapist, Med technician, and Other

Defibrillator section fields table

Group or field	Description
Adjusted "On" Time	You can use this field to synchronize the defibrillator-created date and time with the computer date and time. Event Review Pro automatically adjusts the date and time of each event on the ECG tab for the adjustment you enter.

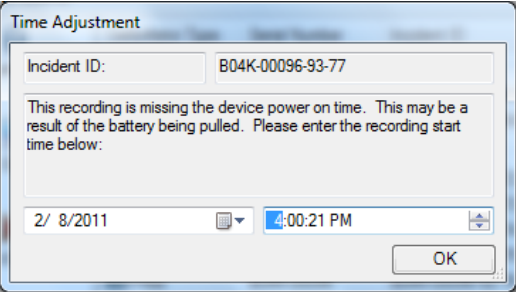Group or field	Description
	<p>Note: Event Review Pro automatically adjusts the date and time based on the clock of the personal computer when it receives a <i>Bluetooth</i> transmission from a HeartStart MRx version 9.0 or later, or from an FR3. This adjustment records the Recorded "On" Time. You can still use this field to override the Recorded "On" Time.</p>
Recorded "On" Time	<p>When you attach an ECG, the field displays the date and time that the defibrillator was turned on.</p> <p>HS1 and FRx only: If you removed the battery before the ECG was attached (downloaded), Event Review Pro displays a message requesting the case date and time. When you click OK, you see this screen:</p>  <p>Set the case date and time to the Recorded "On" Time of the ECG.</p> <p>You can adjust the time in the Adjusted "On" Time field.</p>
Debrillator Type	<p>When you attach an ECG, the field displays the defibrillator model name.</p> <p>You cannot change it.</p>
Serial Number	<p>When you attach an ECG, the field displays the defibrillator serial number.</p> <p>You cannot change it.</p>
Shocks	<p>When you attach an ECG, the field displays the total number of shocks delivered to the patient by this defibrillator. You cannot change it.</p>

Patient section fields table

Group or field	Description
Patient ID	<p>By default, the value in the Patient ID field is added based on the patient ID from the defibrillator (if the device supports entering such data). Otherwise, it is blank, and you can add an ID.</p> <p>The patient ID prints in report headers.</p> <p>A redacted case displays a series of asterisks for the patient name and patient ID.</p>
First name	<p>The first name prints in report headers.</p> <p>If the ECG is from HeartStart MRx, this field can contain information entered on the MRx. You can change it.</p> <p>A redacted case displays a series of asterisks for the patient name and patient ID.</p>
Middle name	<p>The middle name prints in report headers.</p> <p>A redacted case displays a series of asterisks for the patient Event Review name and patient ID.</p>
Last name	<p>The last name prints in report headers.</p> <p>If the ECG is from HeartStart MRx, this field can contain information entered on the MRx. You can change it.</p> <p>A redacted case displays a series of asterisks for the patient name and patient ID.</p>
Date of birth	<p>Type the patient date of birth in the field.</p> <p>A redacted case displays a blank field.</p>

Group or field	Description
Age	<p>This is the age of the patient at the time of the response. Select the format used and type the age in the field. The units of measurements include years, months, or days You can change the format at any time. Event Review Pro calculates the new value for you.</p> <p>If the ECG is attached from HeartStart MRx, this field might contain information entered on the HeartStart MRx. You can change it.</p> <p>If you import a redacted case and the age of the redacted patient is greater than 90, the Age field displays 90.</p> <p>If the date of birth is entered, age is computed automatically, and you may not enter it.</p>
Height	<p>This is the height of the patient.</p> <p>Select the format used by the hospital and type the height in the field. The units of measurement include feet (ft), inches (in), centimeters (cm), or meters (m).</p> <p>You can change the format at any time.</p>
Weight	<p>This is the weight of the patient at the time of the response.</p> <p>Select the format used by your organization and type the weight in the field. The units of measurement include pounds (lb.), ounces (oz.), kilograms (kg), or grams (g).</p> <p>You can change the format at any time.</p>
Gender	<p>The list is:</p> <ul style="list-style-type: none"> Female Male Unknown <p>If the ECG is from HeartStart MRx, this field might contain information entered on the defibrillator.</p>
Race	<p>The list is:</p> <ul style="list-style-type: none"> American Indian or Alaska Native Asian or Pacific Islander Black, Hispanic Black, Non-Hispanic White, Hispanic

Group or field	Description
	White, Non-Hispanic Unknown Other

Describing the conditions at the scene

Use the Scene tab to enter the responder's initial impression of the patient and the scene. The information on this tab appears in the Case Report.

To complete the fields

1. Click the field.
2. Click an option on the menu or type the information.



If Other appears in the option list, you can type a value in the Other text box. If your entry is the same as an option in the list, your entry appears as shown in the option list. Click **Clear** to remove the information in the field or text box.

The following fields can accept more than one input:

- Existing medical treatments
- Past medical history
- Type of first responder
- Defibrillator types used.

Documenting events

Use the Timeline tab to document the events that occurred during the case, the treatments that were administered, and the notes that were collected during the case. The Timeline tab also lists the notes that you add to the ECG on the ECG tab.

Event Review Pro sorts events for a new case in chronological order. You can sort events in ascending (1 to 9, a to z) or descending (9 to 1, z to a) order.

Events and notes appear in the Event Report, on the ECG tab, and on the Event Log tab.



If **Other** appears in the option list, you can type a value in the **Other** text box. If your entry is the same as an option in the list, your entry appears as shown in the option list. Click **Clear** to remove the information in the field or text box.

Events that you add on the ECG tab appear on the Timeline tab as read-only. You can edit the Event Name parameters and the Comment fields, but not the Date and Time field.

Adding, describing, and removing events

You can add, edit, and remove an event in an open case.

To add an event

1. On the **Timeline** tab, click the **Event Name** field.

A menu showing types of events appears. You can click an event that is specific to the phase of events, an administered treatment, and a note added by a responder or reviewer.

2. Click an event from one of the menu columns.

Event Review Pro adds the event to the table in chronological order.

- If an event does not appear in the list of options in the Phase or Treatment columns, click the Custom option in the appropriate column to create one.
- Click **Responder Note** to document notes from a responder.
- Click **Reviewer Note** to document notes from a medical director, code team leader, or operations manager.
- If the event displays additional fields, complete the appropriate information for each field.

Event Review Pro displays a button to the left of the event. Click the button to collapse or expand the additional fields that document the event.

3. (Optional) Type text in the **Comment** field.
4. Repeat steps 1, 2, and 3 for each event, as necessary.

To edit an event

1. On the **Timeline** tab, click the appropriate event in the **Event Name** field.
2. Use one of the following methods:
 - Select the appropriate option on the menu.
 - Type the correction.





To remove one event

1. On the **Timeline** tab, click the **Event Name** field.
2. On the menu, click **Remove this event**.

The event disappears from the event table.

Hiding and displaying event details


You can hide or display the events details for the open case. For more information, see [Displaying case details on page 43](#).

- To refresh the list of cases, on the **View** menu or shortcut menu, click **Refresh**.
- To hide one event row, on the **Timeline** tab, click **Collapse**  to the left of the Event Name field.
- To hide all event rows, on the **View** menu or shortcut menu, click **Collapse All** .
- To display one event row, on the **Timeline** tab, click **Expand**  to the left of the Event Name field.
- To display all rows, on the **View** menu or shortcut menu, click **Expand All** .

Sorting events

You can sort the information based on the column. Events in a new case are first sorted in chronological order.

Click the header to sort the list of values in ascending (1 to 9, or a to z) or descending (9 to 1, or z to a) order. A triangular symbol (arrow) appears to indicate the sort order.

 Event Review Pro adds new events based on the sort order. The location of the event field for the new event changes depending on the sort order. For example, if you sort events in alphabetical order, Event Review Pro inserts the event alphabetically. The blank event field is at the top of the table. If you sort the event by date and time in descending order, Event Review Pro inserts the event at the top of the table. The blank event field is at the bottom of the table.

Changing the date and time of the event

You can change the date and time of an event that you entered on the Timeline tab.

To change the event date and time

1. Click the event.
2. Click **Date and Time**.
3. Click the value you want to change and type the correct date or time.

EMS Edition only: The following fields are used and incorporated in Response Times reports:

- Call Receipt
- Bystander CPR
- Dispatcher Began AED or CPR Instruction
- BLS Unit Notified
- BLS Vehicle Mobile
- BLS Arrival at Scene
- BLS EMS CPR
- ALS Unit Notified
- ALS Vehicle Mobile
- ALS Arrival at Scene
- ALS arrival at Patient
- ALS EMS CPR

EMS Edition only: The following fields are used and incorporated in Utstein reports:

- Admitted ICU or Ward
- Patient Death

Documenting the patient's outcome

Use the Follow-up tab and the Outcome tab to enter information about the patient during and after hospitalization.

To complete the fields

1. Click the field.
2. Click an option on the menu or type the information in the field.



If Other appears in the option list, you can type a value in the Other text box. If your entry is the same as an option in the list, displays your entry as shown in the option list. Click **Clear** to remove the information in the field or text box.

Several fields use scores for Cerebral Performance Categories (CPC) and Overall Performance Categories (OPC). For more information, see [CPC and OPC on page 160](#).

EMS Edition only: Utstein reports use the Any ROSC field:

Working with ECGs

Use the ECG tab in Cases to add, review, and print ECG waveforms. You can also use the ECG tab to enter information and review events.

You can complete the following tasks on the ECG tab:

- View an attached ECG
- Add notes to an ECG
- Modify the display properties of an ECG

The ECG waveform includes recorded defibrillator events as well as reviewer and responder notes. For information on attaching an ECG, see [Adding ECGs on page 34](#).

For details on the displays and functions of the ECG tab, see [Working in the ECG tab on page 65](#).

If the case includes 12-lead information, click the 12-lead tab within the ECG view to see a separate waveform and accompanying information. For details, see [Working with 12-lead ECGs on page 90](#).

Working in the ECG tab

To display an attached ECG, open a case and click the ECG tab. For more information, see [Displaying case details on page 43](#).

Events pane

All defibrillator events and notes appear in chronological order on the Events pane. You can expand events or notes in the Events pane to display additional details. For more information, see [Using the Events pane on page 68](#).

Use the Events toolbar to search for specific events or notes. You can also add a note for an ECG. Use the event filter to control the display of events and notes on the waveform. The event filter displays all possible event types and categories. Check or uncheck them to control whether they appear.

Overview pane

The Overview pane displays the ECG waveform with flagged events and notes for each channel.

You can use the Overview toolbar to:

- Select the channel displayed on the Overview pane
- Select portions of the ECG for detailed view on the Channels pane
- Change the time scale for the ECG

For more information, see [Working with waveforms on page 72](#).

Waveforms (Channels) pane

The Waveform pane provides a more detailed view of waveforms. Event Review Pro can display multiple waveforms in separate channels on this pane.

Use the Waveforms toolbar to

- Zoom in on a channel waveform
- Zoom out on a channel waveform
- Change the time scale for the channel waveforms

You can also use the Waveforms pane to do the following:

- Add, edit, or delete notes
- Control the display of waveforms in the individual channels
- Analyze the quality of the CPR administered

Transport pane


The Transport pane at the bottom of the ECG view displays the ECG lead information, the sweep bar time, and the cursor time relative to the sweep bar

location on the ECG. Use the Play button  to play the entire ECG. For more information, see [Reviewing channel information on page 79](#).

If an audio recording is associated with the ECG, you can use the audio controls to play the audio while viewing correlated events on the waveform on the Channels pane. For more information, see [Reviewing channel information on page 79](#).

You can also use the controls on the Transport pane to find and select events on the waveform for review. For more information, see [Reviewing channel information on page 79](#).

The following illustration identifies the components of the ECG view

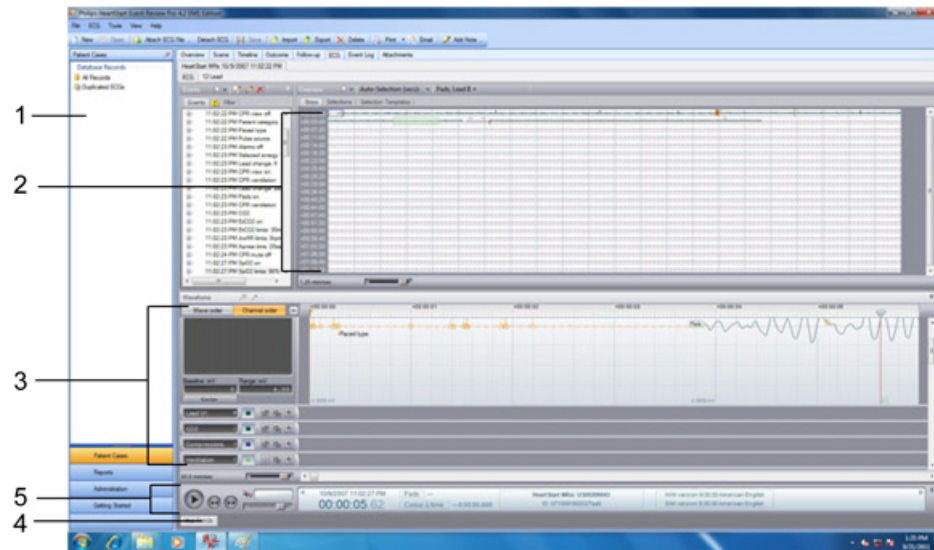

Key to illustration numbers

1	Events
2	Overview
3	Channels
4	Vital Trends
5	Transport

Using the Events pane

A list of defibrillator events and notes associated with an ECG appears in an expandable tree structure on the Events pane. Defibrillator events can include items such as the following:








- Shocks
- Alarms
- Monitor settings
- Equipment status
- Q-CPR events

In the Events pane, you can identify types of events, view them, search for them, and display/hide them. You can also manage the notes that are attached to each case.

Identifying types of events


Event Review Pro displays the notes that were entered when the case was added to the Event pane. For more information, see [Documenting events on page 60](#).

The following icons identify event entries:

Symbol	Meaning
	Note: A user-input observation about the patient's status or the responder's intervention made on the Timeline tab.
	Note: A user-input observation about the patient's status or the responder's intervention made on the ECG tab.
	Shock: The delivery of therapy by the defibrillator.
	Device On: The time the defibrillator was turned on or returned to use.
	Device Off: The time the defibrillator was turned off.
	Defibrillator alarm: An alert provided by the defibrillator. The alarm calls attention to a patient vital sign that is beyond a preset threshold.
	Error condition: An error related to the ability of the defibrillator or part of the defibrillator to perform its intended monitoring function.

Viewing events

When you click an event in the Events pane, the sweep bar automatically moves to the location of the event on the ECG strip on the Overview and Channels panes. You can expand the event in the Events pane to view event details.

To view event details, in the Events pane, click Expand  next to the event you want to view. Event details appear. You can also double-click the event label to view details.

To expand an Events pane node

1. Click an event or event node.
2. Click the PLUS SIGN (+).

To collapse an Events pane node

1. Click an event or event node.
2. Click the MINUS SIGN (-).


To find and view an event using the scroll bar

1. Click and drag the scroll bar down the Events pane to find an event.
2. Click the event .
The sweep bar moves to the location of the event on the waveform.
3. Double-click the event .
The event expands to display event details.

Searching for events

Use the Search  button on the **Events** toolbar to search for events or notes.

To search for an event

1. Type or select the event name in the text box on the **Events** toolbar.
The search text can be any portion (substring) of the event name as it appears in the event label. For example, to search for "patient age," you could type "patient," "patient age," "age," or even just "patien."
2. Click **Search**  .


If the search finds an event, it appears selected in the event tree, and the sweep bar moves to the location of the event on the waveform.

Hiding or displaying events


Use the Filter tab to display or hide events on the waveform on the Overview, Channels, and Events panes according to criteria that you choose.

You can display or hide events either by category or by type. If you filter by category (for example, if you filter by Q-CPR), then all event types in that category (CPR Pads Off, CPR Intubated, etc.) are hidden as well.


By default, all events appear except CPR compression events.

The caution  icon appears on the Filter tab when there is at least one event type or category that was filtered. See the procedure below to display some or all of these events.

To use the event filter

1. On the Events pane, click the **Filter** tab.
A list of event categories appears.
2. Click **Expand**  next to an event category.
The event types in the selected category appear.
3. Uncheck the box for any event type or event category that you want to remove from the display.

The following results occur:

- The event or events disappears from the **Overview** and **Channels** panes.
- The event or events disappears from the Events pane.
- The caution  icon on the **Filter** tab appears.

Managing notes

You can add, modify, and delete notes in the Events pane. You can also add, edit, and delete notes directly on the ECG on the Overview or Channels pane. For more information, see [Managing notes on the waveform on page 86](#).

To add a note in the Events pane

1. Click the part of the Events pane where you want the note to appear.
2. On the toolbar, click **Add Note** .

The Note window appears, with a table similar to the table on the Timeline tab.

3. Click in the **Annotation Type** field.
4. From the list of events in the pop-up window, click and complete an event the same as when you add events on the **Timeline** tab.

For more information, see [Documenting events on page 60](#).

5. Complete the fields, as appropriate.


For example, for a Reviewer Note, type *Reviewer name*, *Reviewer title*, and *Reviewer department*.

If you click **Custom**, type the event name and unit, if appropriate.


6. Click **OK**.

The note now appears on the event tree, the ECG waveform, and the Timeline tab.

To modify a note in the Events pane

1. Click the note you want to modify.
2. On the toolbar, click **Edit Note** .
3. In the Note window, modify the note.
4. Click **OK**.

To delete a note in the Events pane

1. Click the note you want to delete.
2. On the toolbar, click **Delete Note** . The note disappears from the Event pane, the ECG waveform, and the Timeline tab.

Viewing waveforms

You can see the entire ECG as a continuous strip on the Overview pane. It includes markers and labels for recorded events and notes. You can point to event labels on the ECG on the Overview pane to view event details. In the Channel view, you can zoom in to see an ECG segment of interest.


On the Overview and Waveforms (Channels) panes, a vertical sweep bar indicates the current ECG time. The sweep bar provides a visual cue to the current position on the waveform. Click the ECG on the Overview pane to move the sweep bar to the pointer location on the Waveforms pane. For more information, see [Working with waveforms on page 72](#).

You can see a closeup view of a segment of the ECG waveform, as well as a closeup of compression waveforms, on the Waveforms pane. The channel displays the waveform in a grid, indicates the selected scale, and displays events and notes with markers.

Working with waveforms

You can manage notes on the Overview and Waveforms panes.

On the Overview pane, a dashed box (arrow, below) identifies the waveform that appears on the Channels pane.

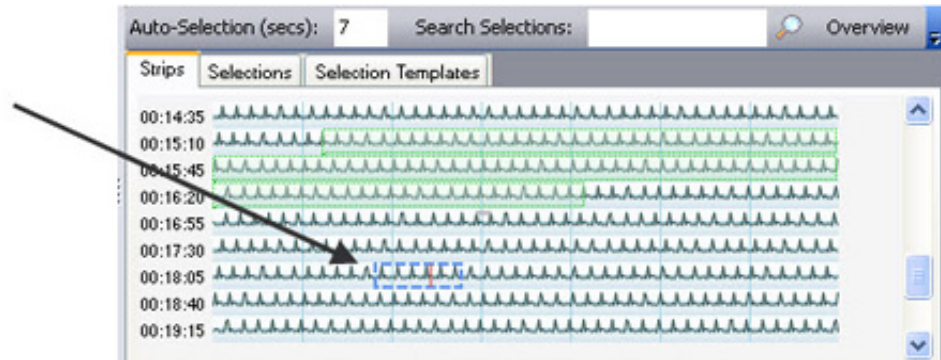

To add and manage notes on the waveform, see [Managing notes on the waveform on page 86](#).

This ECG is intended only for basic rhythm identification. It is not intended for diagnostic and ST segment interpretation,

By default, the primary defibrillator channel appears on the Overview pane.

To select a channel

1. On the **Overview** toolbar, click the **Channel** drop-down list.
2. From the list of all recorded channels, click the channel that you want to view on the Overview pane.

To view event details on the waveform, hover the mouse over the event marker on the ECG waveform. Event details appear in a pop-up window.

To select an event on the waveform, click on or near the event marker on the waveform on the Overview pane. The sweep bar moves to the location of the event on the ECG waveform on the Waveforms pane.

Magnifying waveforms

Use the magnifier to get a detailed view of a section on the waveform.



An asterisk appears for each incomplete release on the compression waveform, similar to the view on HeartStart MRx.

To use the magnifier tool

1. Hold down **SHIFT** and press the mouse button as you point to the area on the waveform that you want to magnify.


A magnified view of the waveform appears, with crosshairs at the location of the pointer.


2. While continuing to press the mouse button, move the pointer to another location on the waveform for closer inspection.

The elapsed time and waveform values appear for the magnified section, as well as the endpoints of each axis.

Displaying and hiding channels

If the ECG has multiple channels, you can select which channels to display or hide on the Channels pane. You can also increase the display area of the Channels pane to accommodate multiple channels.

To display additional channels on the Channels pane, click **Expand**  on the Channel Properties tool for the channel you want to show. The selected channel then appears. You can use the vertical scroll bar to see additional channels. You can also drag the border to increase the size of the Channels pane.

To hide a channel on the Channels pane, click **Collapse**  on the Channel Properties tool for the channel you want to hide.


For information on docking/undocking and pinning/unpinning areas of the ECG display, see [Arranging sections of the ECG display on page 73](#).

Arranging sections of the ECG display


You can use the pinning function to customize the display of sections in the ECG display.

 If you ever need to put the sections back to their default locations, on the **Configuration** menu, click **Reset Panel Layout**.

"Pinning" and "unpinning" refer to the size of an ECG section. A pinned section might look like the illustration below:

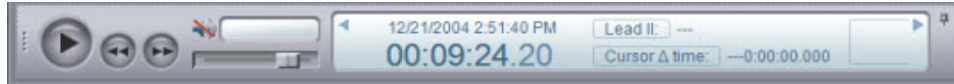

When you click the pin for the section, it is minimized; all that you see is a marker for it at the edge of the display, as in the figure below.

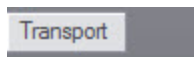


Click the marker to pin it again.

Changing the display time

At the top of the Channels pane, a ruler indicates the time scale of the waveform displayed in elapsed time or real time. Use the channels ruler to change the time scale between elapsed time or real time.

- Real time is the actual clock time based on the Adjusted Time.
- Elapsed time is the time since the defibrillator was turned on.

The time scale appears with divisions for seconds and fractions of seconds.

 If an audio clip is associated with the waveform, the channels ruler displays a blue overlay where the audio is present, in relation to the waveforms.

To change the display time


1. Right-click the channels ruler.
2. On the shortcut menu that appears, click one of the following time formats:
 - Elapsed Time
 - Real Time

The display time changes according to your selection.

Changing the time scale

Use the time scale slider on the bottom left side of the Waveforms (Channels) pane to select the time scale in mm/sec. The default time scale is 25.0 mm/sec. You can also use the mouse wheel to change the time scale if the Channels pane is currently in use.

The time-scale grid adjusts according to the time scale selected. At the default time-scale setting of 25.0 mm/sec, the grid marks indicate 1-second intervals.

 If you decrease the time scale, the grid marks indicate larger intervals. For example, if you decrease the time scale to 12.5 mm/sec, the grid marks indicate 5-second intervals.

To change the time scale using the time scale slider, click the time-scale slider and complete one of the following procedures:

- To increase the scale, drag the slider to the right.
- To decrease the scale, drag the slider to the left.

The adjusted scale appears in mm/sec., and the waveform display adjusts accordingly.

To change the time scale using the mouse wheel, use one of the following methods. Point to the waveform and click the mouse wheel to complete one of the following actions:

- To increase the scale, roll the mouse wheel forward.
- To decrease the scale, roll the mouse wheel backward.

The adjusted scale appears in mm/sec., and the waveform display adjusts accordingly.

To reset the time scale to default, hold down CTRL and click the time-scale slider. The time scale resets to the default scale (25.0 mm/sec.).

Scrolling through the channel

Use the channel scroll bar at the bottom of the Channels pane to scroll through the waveform strip. You can also scroll through the strip using key command shortcuts. For more information, see [Using key command and mouse shortcuts on page 96](#).

To scroll through the waveform, drag the channel scroll bar to the desired location on the waveform.


Customizing channels

Event Review Pro provides a channel properties tool for each channel. This tool displays information about the channel waveform and provides features to modify the waveform display. You can select to display or hide the tool. By default, the tool is closed; click the arrow key on the left (see below) to open the tool.

When the properties area is closed, it looks like this:



When the area is open, it might look like this:

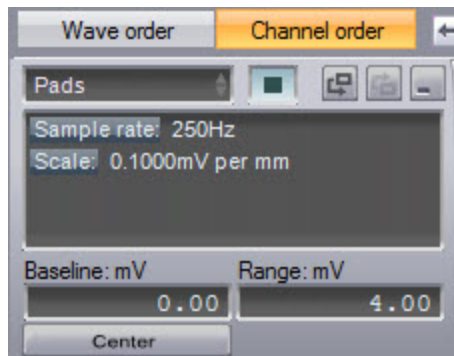





The channel properties tool displays the following information about the current waveform:


- The waveform type or source
- The sample rate in Hz
- The current vertical offset and scale
- The scale in millivolts per millimeter

If multiple waveforms are associated with the channel, you can select the waveform to display in the channel. You can control the vertical offset and can scale the waveform displayed in a channel. You can also control the color of the waveforms displayed on the channels.

If the ECG has multiple channels, you can change the display order for the available channels.


- To display the channel properties tool, click Expand . To hide the tool, click Collapse .
- To move a channel up on the Channels pane, click Move Up .
- To move a channel down on the Channels pane, click Move Down .

To change the color of the waveform

1. Click the color selector  to the right of the waveform ID strip.
2. In the Color window, click the desired color and click **OK**.

The waveform appears using the selected color.

Changing the waveform offset


The Baseline field, available in the ECG channel when you expand the Properties area with the  arrow, controls the vertical offset for the waveform. The Baseline field displays the current vertical offset in the relevant waveform units, for example, *mV* for an ECG. Use the Center button to center the waveform segment displayed in a channel. The waveform appears with the offset determined by the average of the waveform values in the segment currently displayed in the channel.

The image below shows a display before centering:



The image below shows a display after centering:



 You can adjust multidigit scale values based on the digit selected. For example, if the value in the Range field is 100 mV and the third digit is increased by 1, the scale increases to 101; if the second digit is increased by 1, the scale increases to 110.


To change the waveform vertical offset gradually, single-click the specific digit you want to change in the **Baseline** field (for instance, if you want to change 125.63 to 125.69), click the 3 so that it is highlighted, and complete one of the following procedures:


- To increase the vertical offset, turn the mouse scroll wheel forward, or move the mouse so that the cursor moves upward on the screen.
- To decrease the vertical offset, turn the mouse scroll wheel backward, or move the mouse so that the cursor moves downward on the screen.

The adjusted vertical offset value appears in the Baseline field, and modifies the vertical offset of the waveform displayed in the Channels pane, accordingly.

To specify the waveform vertical offset, double click the **Baseline** field and type the offset value. The waveform in the Channels pane appears with the specified offset.

Changing the waveform scale

The Range field, available in the ECG channel when you expand the Properties area with the  arrow, controls the scale of the waveform.

 You can adjust multidigit scale values based on the digit selected. For example, if the value in the Range field is 100 mV and the third digit is increased by 1, the scale increases to 101; if the second digit is increased by 1, the scale increases to 110.

To change the waveform scale gradually, single-click the specific digit that you want to change in the **Range** field (for instance, if you want to change 125.63 to 125.69), click the 3 so that it is highlighted, and complete one of the following procedures:

- To increase the waveform scale, turn the mouse scroll wheel forward, or move the mouse so that the cursor moves upward on the screen.
- To decrease the waveform scale, turn the mouse scroll wheel backward, or move the mouse so that the cursor moves downward on the screen.

The adjusted scale value appears in the Range field, and modifies the scale of the waveform displayed in the Channels pane, accordingly.

To specify the waveform scale, double-click the **Range** field and type the scale value in *mV*. The waveform appears in the Channels pane according to the specified scale.

To restore the default waveform scale, press CTRL , and click in the **Range** field.




Reviewing channel information

Use the transport control (shown below) to navigate through events on an ECG waveform, play the waveform, and play an associated audio clip while reviewing the waveform.





You can also review information about the current sweep bar and pointer locations in an information pane at the bottom of the ECG view.

Navigating ECG events on the waveform

- If you click Next , the sweep bar moves to the location of the next event on the ECG waveform and displays the event time on the information pane.
- If you click Previous , the sweep bar moves to the location of the previous event on the ECG waveform and displays the event time on the information pane.
- If you click Play , the waveform moves at the speed of real time.

Playing the ECG audio

If audio clip was recorded with the ECG, you can play the audio and move along (traverse) the waveform on the Channels pane. The sweep bar indicates the location on the waveform corresponding to the current audio. You can also play and pause the audio by pressing the spacebar if the Channels pane is currently in use. For more information, see [Using key command and mouse shortcuts on page 96](#).

- To play the ECG audio, click **Play** . The audio clip, and the sweep bar moves through the waveform on the Channels pane and indicates the real time on the information pane. If a section of the ECG has no audio, you will see a flashing indicator that displays "No audio."
- To pause the ECG audio, click **Pause** . The audio pauses, and the sweep bar stops at the current location on the ECG waveform.
- To change the volume of the ECG audio, drag the **Volume** slider to the level you want. This adjustment is relative to the volume that you set in the Windows audio control. If your Windows audio is muted or the volume is very low, the Event Review Pro volume setting cannot override it.
- To reset the ECG audio volume to the default level, hold down CTRL and click the **Volume** slider.

Viewing ECG information

The date of the ECG, and the real time and elapsed time for the sweep bar location, appear in an information pane at the bottom of the ECG view.

The waveform type and value for the current pointer location also appears, as well as the cursor delta time. The cursor delta time is the difference in milliseconds between the current pointer location on the waveform and the sweep bar location.

To measure the cursor delta time, click a location on the waveform. The difference between the pointer location and the sweep bar appears in the Cursor time field on the information pane.

To measure the waveform value, click a location on the waveform. Event Review Pro dynamically assigns a name to the waveform type field, according to the waveform type or source, and displays the waveform value in the waveform value field.

The waveform value is displayed in units appropriate to the waveform type, as in the following examples:

- *mV* for leads

- *mmHg* for blood pressure and CO2
- *mm* for compression
- *mOhms* for ventilation

Defining and viewing ECG selections

ECG selections are portions of the ECG that you need to reference for future examination—like bookmarks. You can mark any part of an ECG as a selection and give it a name; you can also assign it a distinctive color to set it off. With selection templates, you can set up a process to create selections automatically when a specific event in the ECG occurs. Once you have more than one selection, you can sort them for quick reference.

You can also print a strip report of the ECG from the File menu or the application toolbar. For more information, see [Generating reports on page 107](#) and [Working with ECG reports on page 110](#)

Shortcut menu commands on the Channels pane assist in these tasks.

Creating ECG selections

If you want to isolate a section of ECG data for separate attention (such as focused analysis during debriefing, or printing it by itself and inserting it in slide presentations, documents, or emails), you can create a "selection." You can do this three different ways:

- To create an ECG selection manually, click where you want the clip to begin on the waveform, and press and hold the mouse button as you move the pointer to where you would like the clip to end. A green overlay (the default color; you can change this) highlights the selection.
- If you want to create a single selection of a predetermined length in the part of the ECG that is on the screen at the moment, use the Auto-selection feature. You can set the length in seconds before you create the selection.
- If you want to set a template to create a selection automatically every time a certain event occurs (like an alarm, a defibrillation, or a ventilation), see [Using selection templates on page 84](#).

To automatically create ECG selections

1. In the **Auto-selection (secs)** text box, type the ECG automatic selection length in seconds.

2. Hold down CTRL and click the ECG where you want to create an automatic selection.

A selection appears that is centered on the pointer location. A colored overlay highlights the selection.

Managing ECG selections

Use the ECG menu or shortcut menu on the Overview or Channels pane to manage ECG selections. To display the shortcut menu, right-click the ECG.

In addition to printing and copying ECG selections, you can do the following tasks:

To create an ECG selection using the shortcut menu

1. Move the sweep bar to the point in the waveform where you want the selection to begin.
2. On the ECG or shortcut menu, click **Selections**, and then click **Set Selection Start**.
3. Move the sweep bar to the point in the waveform where you want the selection to end.
4. On the ECG or shortcut menu, click **Selections**, and then click **Set Selection End**.

To change the color of the ECG selection

1. Click the selection on the waveform.
2. On the ECG or shortcut menu, click **Selections** and then click **Selection Color**.
3. In the Color window, click a color and click **OK**. You can also click **Define Custom Colors** to create your own colors; type values in the boxes for Hue, Saturation, Luminosity, and RGB, preview the results until you have the color you want, click **Add to Custom Colors**, and click **OK**.
4. On the **File** menu or toolbar, click **Save**.

You can now export and import the case with the current ECG view. Event Review Pro also saves the color for future selections.

To convert an ECG selection to a Q-CPR exclusion

1. Click the selection on the waveform.
2. On the ECG or shortcut menu, click **Selections**, and then click **Convert to Q-CPR Exclusion**.

The selection is now a Q-CPR exclusion and is available on the shortcut menu. For more information on Q-CPR reports, see [Working with Q-CPR reports on page 111](#). For more information on Q-CPR exclusions from reports, see [Removing Q-CPR exclusions on page 95](#).

To clear an ECG selection

1. Click the selection on the waveform that you want to remove.
2. Use one of the following methods:
 - On the ECG menu, click **Selections**.
 - Right-click the selection to display the shortcut menu and click **Selections**.
3. From the list of selection commands, click **Remove Selection**.

The selection disappears.

To clear all ECG selections

1. Use one of the following methods:
 - On the ECG menu, click **Selections**.
 - Right-click the selection to display the shortcut menu and click **Selections**.
2. From the list of selection commands, click **Remove All Selections**.

All selections disappear.

Naming, finding, and sorting ECG selections

All ECG selections appear in tabular form on the Selections tab. Use the Selections tab to sort, assign names to, and search for selections.




You can use selection names as bookmarks. The selection name appears at the bottom of a named selection in the Channels pane.

To name ECG selections

1. Click the **Selections** tab.

All selections appear, sorted according to time of creation.
2. Click the name field for a selection to highlight the number assigned at creation.
3. Type a name or unique ID for the selection.
4. Repeat steps 2 and 3 for each selection.

To search for a selection

1. On the Overview toolbar, type search criteria in the **Search selections** field.
Event Review Pro searches for selections based on the name field. The search text must match the selection name in the same order as it appears on the Selections tab and cannot include a "wildcard" (that is, a placeholder for an unknown part of the search string). For example, to search for Shock 1 delivered, you can type "shock" or "Shock 1". Case does not matter. However, if you type only "Shock delivered," you will not get results. If you want to search for Shock 2 delivered, you must specify "Shock 2" to find the subsequent selection.
2. Click **Search** .
The waveform selection matching the search text is now highlighted in the Channels pane.

To change the sort order

1. Click the **Selections** tab.
All selections appear, sorted according to time of creation.
2. Click the heading of the field you want to sort by.
The selections sort themselves by the field that you chose.
3. Click the field heading again to toggle between ascending or descending order.

To move quickly from one waveform selection to the next, use the UP ARROW and DOWN ARROW keys.

Using selection templates

Use selection templates to create event-based ECG selections automatically. For instance, you can create a template that will mark all shocks delivered, showing you 2 seconds before and after the event, and make the selection blue.

To create a selection template

1. Click the **Selection Templates** tab.
The headings for a selection templates table appear.
2. Click **New**.
A new row appears in the selection templates table.
3. From the **Event Category** drop-down list, click an event category.

4. Click in the **Event Type** field for the new selection template.
5. From the **Event Type** drop-down list, click an event type.
6. In the **Pre-seconds** field, type the number of seconds you want to precede the event in the selection.
7. In the **Post-seconds** field, type the number of seconds you want to follow the event in the selection.
8. In the **Color** field, select the color you want to highlight the selection.
9. Click **Create**.

ECG selections appear according to your criteria.

To delete a selection template

1. Click the **Selection Templates** tab.
The headings for a selection templates table appear.
2. Click a template and click **Delete**.

Copying ECGs to the Clipboard

You can copy ECG strips and ECG selections to the Clipboard. You can then paste the strips and selections into documents created by a word processing or graphic application. Copy the ECG strips and selections from the ECG view.

To copy an ECG strip to the Clipboard


1. Right-click the strip in the **Overview** or **Channels** pane.
2. On the menu that appears, click **Copy to Clipboard**.
3. Paste the contents of the clipboard into a document.

To copy an ECG selection to the Clipboard

1. Right-click the selection you want to copy.
2. On the menu that appears, click **Selections**.
3. Click **Copy Selection to Clipboard**.
4. Paste the contents of the clipboard into a document.

Exporting waveform data

To export a case with the *.wfd* format, which is used for data analysis, use the Export option on the File menu or toolbar. For more information on this format, see <http://www.physionet.org>.

 Event Review Pro cannot *import* *.wfd* files. For more information, see [Importing case files automatically on page 46](#).

To export only ECG data

1. Open a case. For more information, see [Displaying case details on page 43](#).
2. Click the **ECG** tab.
3. Click the **ECG view** tab.
4. On the **ECG** menu, click **Export to WFDB**.


The Save As window opens.

5. In the **Save in** box, navigate to the location where you want to save the exported ECG data.
6. In the **File name** box, type a name for the data files.
7. Click **Save**.

Event Review Pro exports the files to the folder that you name.


Managing notes on the waveform

You can add, edit, and delete notes on the waveform. You can also use the application toolbar to add notes.

 If you want to add a note that is "fixed" on the ECG waveform—that is, its time and date do not change—use the Timeline tab. If you want to add a note whose time and date move according to the Adjusted "On" time, set on the Defibrillators section of the Overview tab, use the ECG menu.

To add a note

1. Click the waveform where you want the note to appear.
The sweep bar moves to the pointer location.
2. Use one of the following methods:

- On the ECG menu, click **Add Note**.
- On the Events pane toolbar, click **Add Note** .
- Right-click the waveform to display the shortcut menu, and click **Add Note**.

The Note window opens, with a table similar to the table on the Timeline tab.


3. Click the **Event Name** field.
4. On the list of events, click and complete an event from the list the same as when you add events on the **Timeline** tab.

For more information, see [Documenting events on page 60](#).


5. Complete the fields, as appropriate.
6. Click **OK**.

The note now appears on the waveform, the event tree, and the Timeline tab.

To modify a note

1. On the event tree, click the note that you want to modify.
2. Use one of the following methods:
 - On the Events pane toolbar, click **Edit Note** .
 - Right-click the waveform to display the shortcut menu, and click **Edit Note**.
3. In the Note window, make your modifications to the note.
4. Click **OK**.

To delete a note

1. Click the note that you want to delete.
2. Use one of the following methods:
 - On the Events pane toolbar, click **Delete Note** .
 - Right-click the waveform to display the shortcut menu, and click **Delete Note**.



The note disappears from the waveform.

Zooming in and out of the waveform

Use the ECG menu or shortcut menu to zoom in and out of the waveform.

To zoom in or zoom out

1. Click the waveform.

2. Use one of the following methods:
 - On the ECG menu, click **Zoom In**  or **Zoom Out** .
 - Right-click the waveform to display the shortcut menu, and click **Zoom In** or **Zoom Out**.
 - On the numeric keypad, press PLUS(+) to zoom in or MINUS(-) to zoom out.

To zoom in or zoom out fully

1. Click the waveform.
2. Press SHIFT and, on the numeric keypad, press + to zoom in fully or - to zoom out fully.

Customizing the ECG display

You can configure the display of the waveforms in the Waveforms pane according to your preferences. You can control the display of device events, notes, and labels. You can also configure the display of the grid and the gradient background in the channels.

To remove ECG information from display

1. Click the **Waveforms** pane.
2. Use one of the following methods:
 - Click the **ECG** menu.
 - Right-click the ECG.
3. Click **Display Items**.

The display items menu opens.
4. Click the item that you want removed from display.

The item disappears from the ECG. For additional information, see [Reviewing case events on page 100](#).

To restore ECG information to the display

1. Click the **Waveforms** pane.
2. Use one of the following methods:
 - Click the **ECG** menu.
 - Right-click the ECG.
3. Click **Display Items**.

The display items menu opens.

- Click the item that you want restored.

The item reappears in the ECG.

Reviewing vital trends data

On defibrillator/monitors like the MRx, vital trends data is associated with the ECG. You can view this data in Event Review Pro in a table or on a chart.

In Table view, you can select the time interval for the vital trends information.

	08:40	08:45	08:50	08:55	09:00	09:05
HR (bpm)	0 [^]	73 [^]	73 [^]	73 [^]	73 [^]	--
Temp (°F)	78.9 [^]	79.1 [^]	79.1 [^]	79.1 [^]	79 [^]	--
AwRR (rpm)	0 [^]	0 [^]	0 [^]	0 [^]	20 [^]	--
ABP mean (mmHg)	57 [^]	96 [^]	96 [^]	96 [^]	96 [^]	--
ABP systolic (mmHg)	120	120 [^]	120 [^]	120 [^]	120 [^]	--

On the Vital Trends toolbar, click the **Interval:** drop-down menu and click 5, 10, 15, 30, or 60 minutes from the list that appears.

In Chart view, you can select the duration of the vital trends data.




On the Vital Trends toolbar, click the **Length:** drop-down menu and click 10, 30, or 60 from the list that appears.

You can also choose the color and point style for each trend displayed.

To change the point style in Chart view

- If you are not already in Chart view, click the **Chart** tab.

- On the configuration pane to the left of the vital trends chart, click the style button to the right of the trend you want changed.

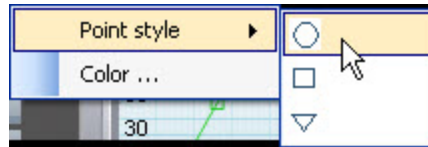

- On the shortcut menu that appears, click **Point style**.
- Click a point style from the list.


To change the trend color in Chart view

- If you are not already in Chart view, click the **Chart** tab.
- On the configuration pane to the left of the vital trends chart, click the style button to the right of the trend you want changed.
- On the shortcut menu that appears, click **Color**.
- In the Color dialog box, click a color from the **Basic Colors** list.

You can also click **Define Custom Colors** to create your own colors; type values in the boxes for Hue, Saturation, Luminosity, and RGB, preview the results until you have the color you want, click **Add to Custom Colors**, and click **OK**.

- Click **OK**.

To remove trends from display in **Chart** view, on the vital trends configuration pane, click the box to the left of the trends (arrow, below).

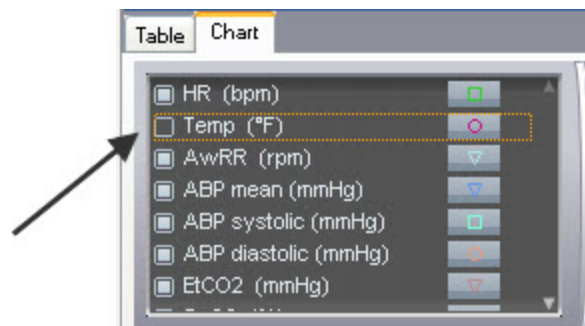

Working with 12-lead ECGs

12-lead data acquired during use of a defibrillator is automatically attached to the case along with continuous ECG. You can also attach 12-lead data that was transmitted through the *Bluetooth* transmission and XML methods.

The 12-lead ECG and associated data appear on the 12-Lead tab. You can also double-click a 12-lead acquisition event in the event pane to display the data on the 12-Lead tab.

This ECG is intended only for basic rhythm identification. It is not intended for diagnostic and ST segment interpretation



Event Review Pro supports HeartStart MRx, version 9.0 or later with the event sum option, for *Bluetooth* transmission. The 12-lead ECG, along with associated summary data, appears in a 12-lead report. When a 12-lead report contains Acute Cardiac Ischemia Time-Insensitive Predictive Instrument (ACI-TIPI) and Thrombolytic Predictive Instrument (TPI) data, the data appears in the summary text box above the waveform.

The 12-lead ECG appears in the middle section of the 12-Lead view. A summary section at the top of the 12-Lead view shows the following 12-lead data:

- Patient demographic information
- Global measurements (Rate, PR, QRSd, QT, QTc)
- Axes (P, QRS, T)
- ECG summary interpretations
- You can use the configuration section at the bottom of the 12-Lead view to control the display of the 12-lead ECG and display additional 12-lead ECGs attached to the case.

You can print 12-lead reports from the File menu or application toolbar. For more information, see [Generating reports on page 107](#).

You can do the following tasks on the 12-Lead tab:

- [Configuring the 12-Lead view on page 91](#)
- [Reviewing the 12-lead ECG on page 92](#)

Configuring the 12-Lead view

Use the configuration section of the 12-Lead view to control the display of the 12-lead ECG.

To select another 12-lead ECG, click the 12-lead drop-down list and click a different ECG. The selected 12-lead ECG and its associate summary information appear.

To change the time format

1. Click the **Timing** drop-down list to change the time format.

The following time formats appear:

- Simultaneous (European standard)
- Time Sequential (United States standard)

2. Click the time segment that you want.

The 12-lead ECG appears in the selected time format.

To scale the report display

You can scale the 12-lead report from 1x, with the Scale control in the far left position, to 3x, with the Scale control in the far right position. You can also use the magnifying tool in the 12-Lead view, when the Scale control is in the 1x position. For more information, see [Magnifying 12-lead ECGs on page 93](#).

- Complete one of the following procedures:
 - To *increase* the scale of the 12-lead ECG displayed, drag the Scale control to the right.
 - To *decrease* the scale of the 12-lead ECG displayed, drag the Scale control to the left.

Reviewing the 12-lead ECG

Use the scroll bar to move through the 12-lead ECG.

If the scale has been increased so that sections of the 12-lead ECG are outside of the current view area, drag the ECG to view the obscured sections. A location key indicates when you have reached the limit of the 12-lead ECG in any direction.

To scroll through the 12-lead ECG, drag the scroll bar to the desired location on the 12-lead ECG.

To use the location key

1. Click a position on the 12-lead display and press the mouse button.
A four-arrow location key appears.
2. Drag the 12-lead ECG display up, down, or to the right or left to display any obscured sections of the 12-lead ECG.

The location key indicates when you have reached the limit of the 12-lead ECG in any direction by dimming the corresponding arrow.

Magnifying 12-lead ECGs

You can use the magnifier tool for a detailed view of the 12-lead ECG waveform.

To use the magnifier tool

1. Hold down **SHIFT** and click the area on the waveform that you want to magnify.
A magnified view of the waveform appears.
2. While continuing to hold the mouse button, move the pointer to another location on the waveform for closer inspection.

Customizing Q-CPR episodes

You can customize the Q-CPR reports that you create. For example, you can use commands on the shortcut menu or the Q-CPR menu to define an episode for Q-CPR reporting and to mark sections of the waveform to exclude from Q-CPR reports.

Q-CPR events appear on the Compression channel.

For information on excluding sections of the waveform from Q-CPR reports, see [Creating Q-CPR exclusions on page 94](#).

For more information on generating and printing Q-CPR Details and !-CPR Report Card reports, see [Working with Q-CPR reports on page 111](#).

To customize a Q-CPR episode

1. If necessary, click the **ECG** tab and display the ECG view.
2. Click the waveform at the point where you want the episode to begin.
3. Use one of the following methods:
 - On the **Q-CPR** menu, click **Set Episode Start**.
 - Right-click the waveform to display the shortcut menu, click **Q-CPR**, and click **Set Episode Start**.

A vertical bar marks the start of the episode.

4. Click the waveform at the point where you want the episode to end.
5. Use one of the following methods:

- On the **Q-CPR** menu, click **Set Episode End**.
- Right-click the waveform to display the shortcut menu, click **Q-CPR**, and click **Set Episode End**.

A vertical bar marks the end of the episode.

Creating Q-CPR exclusions

A Q-CPR exclusion identifies the areas of the waveform that you do *not* want to analyze. Any data *inside* the exclusion area is excluded from the analysis. Use the Q-CPR menu, the shortcut menu on the ECG tab, or shortcut keys to create Q-CPR exclusions.

You can also convert an ECG selection to a Q-CPR exclusion.

To create a Q-CPR exclusion from the Q-CPR menu or shortcut menu

1. If necessary, click the **ECG** tab and display the ECG view.
2. Click the waveform at the point where you want the exclusion to begin.
3. Use one of the following methods:
 - On the **Q-CPR** menu, click **Set Exclusion Start**.
 - Right-click the waveform to display the shortcut menu, click **Q-CPR**, and click **Set Exclusion Start**.
4. Click the waveform at the point where you want the exclusion to end.
5. Use one of the following methods:
 - On the **Q-CPR** menu, click **Set Exclusion End**.
 - Right-click the waveform to display the shortcut menu, click **Q-CPR**, and click **Set Exclusion End**.

A red overlay highlights the exclusion.

To create a Q-CPR exclusion using shortcut keys

1. If necessary, click the **ECG** tab and display the ECG view.
2. Click the waveform where you want the exclusion to begin.
3. Hold down ALT and press the mouse button.
4. Drag the mouse to where you want the exclusion to end.
5. Release the key and mouse.

To create a Q-CPR exclusion from an ECG selection, see [Managing ECG selections on page 82](#).

Removing Q-CPR exclusions

You can remove a single exclusion or all exclusions.

To remove a Q-CPR exclusion

1. If necessary, click the **ECG** tab and display the ECG view.
2. Click the waveform where you want the episode to begin.
3. Use one of the following methods:
 - On the Q-CPR menu, click **Remove Exclusion**.
 - Right-click the waveform to display the shortcut menu, click **Q-CPR**, and click **Remove Exclusion**.

The highlighted exclusion disappears.

To remove all Q-CPR exclusions

1. If necessary, click the **ECG** tab and display the ECG view.
2. Click the waveform.
3. Use one of the following methods:
 - On the Q-CPR menu, click **Remove All Exclusions**.
 - Right-click the waveform to display the shortcut menu, click **Q-CPR**, and click **Remove Exclusions**.

All Q-CPR exclusions disappear.

Adding notes to the Ventilation channel

The notes in the calculations appear in the Q-CPR Report Card report.

To add notes to a ventilation channel

1. If necessary, click the **ECG** tab and display the ECG view.
2. On the **Waveforms** pane, expand the ventilation channel.
3. Click the ventilation waveform at the point that you want to add the note.
4. Right-click the mouse to display the shortcut menu.
5. Click **Q-CPR**, and then **Annotations**.

6. Click one of the following options to add a note, both of which appear at the cursor:
 - **Ventilation**
 - **PR (Perfusing Rhythm)**.

To add notes using keyboard shortcuts

1. If necessary, click the **ECG** tab and display the ECG view.
2. On the **Waveforms** pane, expand the ventilation channel.
3. Click the ventilation waveform at the point that you want to add the note.
4. Press the following key combination:
 - Press **ALT+V** to add a Ventilation note.
 - Press **ALT+R** to add a PR (perfusing rhythm) note at the cursor location.

To move a note to another location on the waveform, click the note on the waveform and drag it to another location using your mouse.

Using key command and mouse shortcuts

You can use the shortcuts in the table below to complete many of the operations on the ECG tab:

Key command shortcuts

Audio

Play/Pause audio

Click the ECG waveform, and then press **SPACEBAR**.

Reset audio volume to default level

Hold down **CTRL** and click the Volume slider

Event tree

Expand selected event tree nodes

SHIFT + PLUS SIGN

Collapse selected event tree nodes

SHIFT + MINUS SIGN

Q-CPR**Annotate compression waveform with ventilation event**

ALT + V

Annotate compression waveform with Perfusing Rhythm.

ALT + R

Selections**Create an automatic selection centered on the pointer location**

CTRL and click waveform

Add a note on a ventilation channel

ALT + V to add a Ventilation note

ALT + R to add a PR (perfusing rhythm) note

Sweep bar**Move sweep bar one minute backward**

PAGE UP

Move sweep bar one minute forward

PAGE DOWN

Move sweep bar one strip backward in the WAVEFORMS pane

SHIFT + PAGE UP

Move sweep bar one strip forward in the WAVEFORMS pane

SHIFT+PAGE DOWN

Move sweep bar to beginning of recording

HOME

Move sweep bar to end of recording

END

Move sweep bar one second backward

LEFT ARROW

Move sweep bar one second forward

RIGHT ARROW

Move sweep bar to previous selection

UP ARROW

Move sweep bar to next selection

DOWN ARROW

Time scale**Time scale expand**

CTRL + PLUS SIGN(+)

Time scale compress

CTRL + MINUS SIGN(-)

Zoom**Zoom in**

PLUS SIGN (+)

Zoom out

MINUS SIGN(-)

Zoom in fully

SHIFT + PLUS SIGN(+)

Zoom out fully

SHIFT + MINUS SIGN(-)

Mouse shortcuts

Magnification**Display magnified view of waveform**

Click the waveform; then hold down SHIFT and press the mouse button.

Q-CPR exclusion**Create Q-CPR Exclusion**

Hold down ALT and drag mouse over waveform .

Time scale**Expand time scale**

Press the mouse wheel and rotate it backward.

Compress time scale

Press the mouse wheel and rotate it forward.

Vertical offset**Increase vertical offset**

Click the BASELINE field and rotate the mouse wheel forward.

Decrease vertical offset

Click the BASELINE field and rotate the mouse wheel backward.

Waveform scale**Increase vertical waveform scale**

Click the RANGE field and rotate the mouse wheel forward.

Decrease vertical waveform scale

Click the RANGE field and rotate the mouse wheel backward.

Reviewing case details

Use the Cases navigation button to add or review case events, plus detailed information about them, and to attach additional files to the case. For example, you can add responder actions and observations, and followup information about the patient.

The events and attached files appear in a table. You can change how the entries appear.

Reviewing case events

Use the Event Log tab to view all the defibrillator events and user notes that are associated with the ECG. The events that appear on the Timeline and ECG tabs are in the Event Log.

In the Event Log, you can do the following:

- View details for an event
- Sort and group the defibrillator events and notes that a responder or reviewer added to the case
- Filter events and notes to display only those that fit your selection criteria
- Print the events for the open case

To display the event log

1. On the **Cases** workspace, click a case to open. For more information, see [Displaying case details on page 43](#).
2. Click the **Event Log** tab.

The table that appears is a log of events for the case in the Event Review Pro database.

Viewing event details

To view event details, double-click an event. The details appear on the ECG tab.

The sweep bar identifies the event. For more information, see [Working with waveforms on page 72](#).



Hiding and displaying events

You can hide or display the list of events for the open case from the shortcut menu.

To use the shortcut menu

1. Open **Cases**. For more information, see [Displaying case details on page 43](#).
2. To display a list of options, use the **View** menu or right-click the **Event Log** table.
3. Click **Refresh**, **Expand All**, or **Collapse All**.

To refresh, hide, and display events

- To refresh the list of events, on the View menu or shortcut menu, click **Refresh**.
- To hide event details, click **Collapse**  to the left of the event.
- To hide all event details, on the **View** menu or shortcut menu, click **Collapse All**.
- To display event details, click **Expand**  to the left of the event.
- To display all event details, on the **View** menu or shortcut menu, click **Expand All**.

Sorting, grouping, and filtering events

You can change the way information displays in the list of events. You can sort and group the list of events, and you can filter them so that only the events that meet the criteria you select appear.

For more information, see [Working with columns on page 121](#), [Grouping and sorting entries on page 121](#), [Filtering entries on page 122](#), and [Removing filters on page 125](#).

Printing case events

You can print the events for the open case or for any case saved in the database.

To print a list of events for the open case

1. On the Event Log tab, decide how you want to format the list before you print it.

For more information, see the procedures above.

2. On the **File** menu or toolbar, click **Print**.
3. On the **Print** menu, click **Tables**.
4. Click **Event Log**.
5. In the **Print** window, click the printer, page range, and number of copies.
6. Click **OK**.

To print a list of events for a case saved in the database

1. Click the **Reports** navigation button.
2. On the **Reports** navigation pane, click a report from the list of reports in the Cases or ECG groups. For more information, see [Working with case reports on page 109](#) and [Working with ECG reports on page 110](#).

Attaching (adding) files to a case

Use the Attachments tab to manage the files that you attach to the open case. For example, you can manage a text document, a *.pdf* file, or a graphics file.

Use the Attachments menu to attach files, open the attachments, and remove attached files from the case.

No attachment file can be larger than 10 Mb.



To view an attachment, you must have the associated application on your computer. For example, to open a picture, you must have an application that can open graphic files.

Displaying the Attachments tab

Use the following procedure to display the Attachments tab.

To display the Attachments tab

1. On the **Cases** workspace, open a case . For more information, see [Displaying case details on page 43](#).
2. Click the **Attachments** tab.

A table of the files that are attached to the case appears. The table is blank when there are no attachments in the Event Review Pro database that are associated with the case.

Adding and removing attachments

Use the following procedures to manage the attachments that are associated with the open case.

To attach a file

1. On the **Attachments** menu, click **Attach**.
2. In the **Attach File** window, navigate to the file location.
3. Click **Open**.

Depending on the size of the file, it can take a few seconds for the file to show in the Attachments list.

To remove an attached file

1. On the **Attachments** tab, click the file description in the **Attachments** table.
2. On the **Attachments** menu, click **Detach**.

Sorting, grouping, and filtering attachments

You can change the way information displays in the list of attached files. You can sort and group the list of files based on the date when the file was attached to the case, the file type, and the file name. For more information, see [Grouping and sorting entries on page 121](#).

You can also decide which columns to display; see [Working with columns on page 121](#).

If you want to limit the number of entries that appear, see [Filtering entries on page 122](#).

A filter shows only those rows in the table that meet the criteria you select.


Hiding and displaying attachment groups

You can also hide or display the group of attachments. The View menu and shortcut menu lists the Refresh, Expand All, and Collapse All options.


To use the shortcut menu

1. Open **Cases**. For more information, see [Displaying case details on page 43](#).
2. Click the **Attachments** tab.
3. Right-click the table.
4. Click **Refresh**, **Expand All**, or **Collapse All**.

To refresh the list of cases, on the **View** menu or short-cut menu, click **Refresh**.

To hide the attached files in a group, click **Collapse**  next to the group.

To hide all event details, on the **View** menu or shortcut menu, click **Collapse All**.

To display the attached files in a group, click **Expand**  next to the group.

To display all attachments, on the **View** menu or shortcut menu, click **Expand All**.

Viewing file contents


To display the contents of an attached file

Use one of the following methods:

- On the **Attachments** menu or toolbar, click **Open**.
- Double-click the row in the **Attachments** table.
The attachment opens in the associated application.

Printing the attachments in an open case

You can print the contents of the files attached to an open case when you have the associated application installed on your computer. Open the attachment in this associated application and print using its features.

 You might not be able to print all attachments, such as a video clip.

You can also print the list of attachments as they appear on the Attachments tab.

To print the list of attached files

1. On the **File** menu or toolbar, click **Print**.
2. On the **Print** menu, click **Tables**.

3. In the **Print** window, click the printer, page range, and number of copies.
4. Click **Print** or **OK**.

Working with reports

Use reports to evaluate the timeliness of, effectiveness of, and trends in your emergency system response. You can generate a report for an open case or directly from the database.

The first step in working with reports is to create a case. Event Review Pro generates the report from the case information.

In Reports, you can also do the following tasks:

- Preview the report as you specify report details
- Print the report
- Export the report

To set up and manage reports, click the **Reports** navigation button. On the Reports navigation pane and workspace, the available reports appear. They are based on your installation. They can appear in the following categories:

Category	Based on
Case	Information on the tabs in Cases
ECG	ECG strip on the ECG tab in Cases
Q-CPR	Quality of CPR
Vital Trends	Vital trends information that appears on the ECG tab in Cases
(EMS Edition only) Response times	Emergency response times
(EMS Edition only) Utstein	Utstein guidelines for reporting on aggregated out-of-hospital cardiac arrests

The workspace provides features to set up, generate, and distribute reports. For more information, see [Generating reports on page 107](#).

Generating reports

You can generate a report for an open case or for a case that is stored in the database. In either case, you use the same steps to generate the reports.

Above the preview of the report, the appropriate parameter settings appear. The default parameter settings appear as selected. You can select or clear the field.

If you select information from a drop-down list, your choice persists the next time you use the field. At that time, you can select another value.



Remember that if you import a redacted case, in any report that includes that case, Event Review Pro replaces the name and patient ID with a series of asterisks. If the redacted age of the patient is greater than 90, the Date of Birth field is blank and the Age field displays 90.

Generating a report for an open case

You can generate a report for an open case on the File menu and toolbar.

To generate a report for an open case

1. Open the case in the **Cases** workspace. For more information, see [Displaying case details on page 43](#).
2. On the **File** menu or toolbar, click **Print**, and then click **Reports**.
3. Click the type of report that you want to generate.
The Report Preview window opens.
4. While you preview the report, you can use the **Report Preview** fields and toolbar to customize the report. For more information, see [Using the Reports toolbar on page 108](#).
5. If you want to add extra reports to this report, check their boxes. These reports may include extra customizable parameters that you can change.
6. On the **Report Preview** toolbar, use a toolbar button to print the report, locate information, export the report, or email the report. For more information, see [Exporting reports on page 118](#).

Generating a report from the database




Use Reports to generate a report for a case saved in the database.






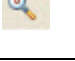






To generate a report from the database

1. Click the **Reports** navigation button.
The Reports navigation pane opens.
2. Double-click the type of report that you want to generate.
On the workspace, above the report preview, you see parameter settings that you can use to specify the report format and details.
3. Depending on the report, do one of the following on the **Filter** tab:
 - a. Click the case from the table and click **Select**.
 - b. Set the filter criteria, and then click **Select**.
4. Depending on the report, click a customization tab.
A preview of the report as you set and change the parameter settings.
5. Complete one or more parameter settings to specify the information that the report displays.
6. To print the report, use one of the following **Report Preview** toolbar buttons:
 - Click **Print** to print the report on the printer of your choice.
 - Click **Quick Print** to print the report on your default printer.
 For more information, see [Printing reports on page 117](#).
7. To generate an export file, click **Export To** on the Report Preview toolbar. For more information, see [Exporting reports on page 118](#).
8. To email the report, click **E-mail As** on the Report Preview toolbar. For more information, see the previous step and [Emailing reports on page 119](#).

Using the Reports toolbar

The following table provides a description of the toolbar buttons on the Reports workspace and Reports Preview window.

Button	Description
	Click Print or press CTRL+P to select a printer, number of copies, or other printing options before printing.
	Click Quick Print to send the document directly to the default printer without making changes.
	Click Find to find text in the document based on your entry in the Find what field and Search options.

Button	Description
	Click First Page or press CTRL+HOME to navigate to the first page of the document.
	Click Previous Page or press PAGEUP to navigate to the previous page of the document.
	Click Next Page or press PAGEDOWN to navigate to the next page of the document.
	Click Last Page or press CTRL+END to navigate to the last page of the document.
	Click Mouse Pointer to show the mouse pointer.
	Click Magnifier to display the full page. Click Magnifier again to zoom to a 100% view of the report page.
	Click Many Pages to select the page layout and to arrange the document pages in preview.
	Click Zoom Out to see more of the page at a reduced size.
	Click Zoom to change the zoom level of the document preview.
	Click Zoom In to get a close-up view of the document.
	Click Export To to export the document in one of the available formats, and save it to a file on a disk. A window to specify parameter settings for the export file appears.
	Click Email As to export the document in one of the available formats, and attach it to the email.

Working with case reports

Case reports include detailed case information based on the information collected on the tabs in Cases. The type of report determines the parameter settings you can specify.

You can generate the reports in the Cases category when a case is open. For information on how to generate a report, see [Generating reports on page 107](#).

After you have generated the report preview, you can choose to append all the available case reports in one report, rather than printing each one separately. On the Case Customization tab, check all the reports that you want to include. Then, on the tab for each of these reports, select the parameter settings you want to include in the report.

For more information, see the following topics:

- [Printing reports on page 117](#)
- [Exporting reports on page 118](#)
- [Emailing reports on page 119](#)

Working with ECG reports

ECG reports are based on the ECG strip on the ECG tab in Cases. See [Working with reports on page 106](#).

ECG Full Disclosure

The ECG Full Disclosure report includes the entire ECG.

The report uses the following parameter settings:

- Case (required)
- Change the scale of the waveform
- Grid
- Events
- Notes
- Beat labels
- Selected channels

Dashes in the report indicate periods during the report when there is no waveform data.

Pre- and Post-shock

The Pre- and Post-shock report documents the ECG segments surrounding each shock event.



The ECG recording begins the presenting rhythm (Presenting seconds) at the start of the ECG, rather than at the time that the defibrillator was turned on.

The report uses the following parameter settings:

- Case (required)
- Change the scale of the waveform
- Grid
- Events
- Notes
- Beat labels
- Selected channels
- Presenting seconds
- Closing seconds
- Preshock seconds
- Postshock seconds

Working with Q-CPR reports

The Q-CPR report card is based on the quality of the CPR that a patient received during the episode. The report uses the information that appears on the ECG tab. You can create a report card based on the CPR results from multiple cases; each report appears in a separate file.

You can choose whether to display the results in millimeters or inches.

You do not have to select parameter settings, and you can choose to exclude ventilation data from the report. See [Working with reports on page 106](#).

The report displays the following information:

- Case ID and patient information
- Episode summary information
- Episode length
- Number of shocks
- Device on and off times
- Ventilation data
- Total number of ventilations

- Total ventilation times before and after intubation
- Average ventilation rate before and after intubation
- No-flow time:
- Flow time
- Average no-flow time before and after shocks
- Use 30 sec. intervals
- Defibrillation data
- Time power was on until the first shock
- Total and average analysis and shock delivery times
- Graphic representation of the quality of CPR for the episode (maximum 30 minutes from the start of the compression waveform))

Working with Q-CPR Details reports

You can see details on all the Q-CPR items that are summarized on the Q-CPR report card and export this report to a spreadsheet. You can also compile details for multiple cases and display them in the same report or in separate reports.

To see what data is included in this report, see [Q-CPR report data on page 112](#).

When you have generated the details report (see [Generating reports on page 107](#)), you can customize the display.

- If you want to show details before and after intubation, in the **Extended Statistics** section, check **Split by Intubation**.
- If you want to change the intervals, click either **No interval**, **Use 30-Sec Interval**, or **Use 60-Sec Interval**.

Q-CPR report data

This section provides a description of the Q-CPR data that appear in the Q-CPR Details report.

You can create a report that includes data from multiple cases; the report can be in one file, with the results separated by case, or in separate reports.

Q-CPR statistic section data label table

Statistic label	Description
Type	The time period for the data: Episode, Period, or Interval
Start	The start time of the statistics coverage, in milliseconds from one second before the first compression event.
Length	The length of the statistics covered time (up to one second after the last compression event), in milliseconds.
NFT	The no-flow time (NFT), in milliseconds.
NFT %	The no-flow time ratio. NFT divided by length, expressed in percent.
FT	The flow time (FT), in milliseconds. It is defined as the length of the statistic minus the NFT, in milliseconds.
FT %	The FT ratio. FT divided by the length, expressed in percent.

Compression section data label table

Statistic label	Description
Comp	The total number of compressions.
Comp Adequate	The total number of compressions that are performed with adequate depth, or are not too shallow
Correct	The total number of compressions that are not too shallow or are not leaning.
Deep	The number of compressions that are too deep.
Shallow	The number of compressions that are too shallow.
Depth	The average depth of compression, in millimeters or inches
Leaning	The number of compressions with incomplete release.
Rate	The average compression rate, per minute. The rate is computed over the active compression time. The time without compression activities is excluded.
Duty Cycle	The compression duty cycle.
Adeq %	The ratio of adequate count to total count.

Ventilation section data label table

Statistic label	Description
Vent	The total number of ventilations.
Vent Rate	The average ventilation rate, per minute. If the statistics are divided into pre- and post-intubation periods, the ventilation rate is for the specific preintubation and postintubation period. If not, the ventilation rate is for the entire episode.

No Sign of Circulation (NSC) section data label table

Statistic label	Description
NSC	The total NSC time, in milliseconds.
Comp	The average compression rate during the NSC time.
NFT/NSC	The ratio of No-Flow Time (NFT) to NSC.

Working with Vital Trends reports

The Vital Trends report provides case information based on the vital trends information that appears on the ECG tab. The information appears in a table format.

The report displays the following information:

- Case ID, device, and patient information
- Vital trends in a table format

The report uses the interval parameter setting. An interval can be 1, 5, 10, 15, 30, or 60 minutes in duration.

Working with Response Times reports


The EMS edition provides Response Times reports.

These reports provide uniform records of your information from out-of-hospital cardiac arrests. They follow the recommendations of the Utstein Consensus Conferences for collecting, summarizing, and reporting response-time information.

You can generate Response Times reports to assess your system's overall response performance and compare performance with service-level commitments.

The type of report determines the parameter settings you can specify. For information on how to generate a report, see [Generating reports on page 107](#).

Event Review Pro saves the parameter settings and displays them the next time you select the report. You can select or clear any of these fields.

 If you import a redacted case, Event Review Pro replaces the name and patient ID with a series of asterisks. If the redacted age of the patient is greater than 90, the Date of birth field is blank and the Age field shows 90.

Average Response Times - Total System

This report summarizes the average response times for key milestones in your system's response protocols. You can specify the date range and site of collapse to limit the cases you want to summarize in the reports.

The report uses the following parameter settings:

- Start date
- End date
- Site of collapse

Percentile Response Times -Total System

This report summarizes response times achieved in a specific percentage of responses. You can specify the date range and site of collapse to limit the cases you want to summarize in the report. You must specify the percentile.

The report uses the following parameter settings:

- Start date
- End date
- Site of collapse
- Percentile (required, median is 50%)

Average Response Times

This report summarizes, for a specific unit, the average response times for key milestones in your response protocols. You can specify the date range and site of

collapse to limit the cases you want to summarize in the reports.

The report uses the following parameter settings:

- Start date
- End date
- Site of collapse
- Unit name

Percentile Response Times

This report summarizes, for a specific unit, the response times achieved in at least a specific percentage of responses. You can specify the date range and site of collapse to limit the cases you want to summarize in the report. You must specify the percentile.

The report uses the following parameter settings:

- Start date
- End date
- Site of collapse
- Unit
- Percentile (Required. The median is 50%.)

Working with Utstein reports

The EMS edition provides Utstein reports.

These reports follow the Utstein guidelines for reporting information on out-of-hospital cardiac arrests. They use the recommended terms and follow the Utstein template approach.

The type of report determines the parameter settings that you can specify. For information on how to generate a report, see [Generating reports on page 107](#).

Event Review Pro saves the parameter settings and displays them the next time you select the report.

Unwitnessed with Bystander CPR

The Unwitnessed with Bystander CPR report summarizes the unwitnessed cardiac

arrests that are in your system.

The report uses the following parameter settings:

- Start date
- End date
- Location
- Population served by EMS system

Witnessed with Bystander CPR

The Witnessed with Bystander CPR report summarizes witnessed cardiac arrest in which a bystander administered CPR.

The report uses the following parameter settings:

- Start date
- End date
- Location
- Population served by EMS system

Witnessed without Bystander CPR

The Witnessed without Bystander CPR report summarizes witnessed cardiac arrests in which no bystander administered CPR.

The report uses the following parameter settings:

- Start date
- End date
- Location
- Population served by EMS system

Printing reports

You can print a report for an open case or for a case that is stored in the database. Before you print a report you must generate the report.

To print a report for an open case

1. From the open case, click **File**, then **Print**, then **Reports**, and then the report that you want to print.
2. On the Customization tab, complete one or more parameter settings to specify the information that the report displays.
3. Click one of the print options.
 - Click **Print** to select a printer, number of copies, and other printing options. When you have made these specifications, click **Print**.
 - Click **Quick Print** to send the report directly to the default printer without making changes.

Exporting reports

You can export a report for an open case or for a case that is stored in the database. Event Review Pro saves a copy of the report to a file outside of the database. You can then email the export file or add the export file to the case on the Attachments tab.



To export a case with the *.wfd* format, which is used by researchers for data analysis, use the Export option on the File menu or toolbar. For more information, see [Exporting cases on page 49](#), (See <http://www.physionet.org> for information on *.wfd* format.)

If you frequently review a report or want to track your system's performance, generate the report, and then use the Export tool to save the report as a file. Sample reports include month-end, quarterly, and year-end reports.

To export a report

1. Generate and view the report. For more information, see [Generating reports on page 107](#).
2. On the **Report Preview** toolbar, click **Export To**.
3. Click the arrow and click the file format for the exported file.

A list of file formats that are appropriate for the type of report appears. Waveform data must be in PDF format.

4. Complete the parameter settings for the file format and click **OK**.

The **Save As** window opens.

Navigate to a location where you want to save this report, or accept the default.

5. In the **File name** field, type a name for the report.
The **Save as Type** field lists the document type that you selected.
6. Click **Save**.
7. If you can open the file, click one of the following buttons:
 - Click **Yes** to preview the file.
 - Click **No** to continue.

Emailing reports

You can email the export file for a report the same way as you email a case.

You can use the E-mail-As option or your email application to send an exported report to another Event Review Pro user.



Email requires that you configure Microsoft Outlook or a MAPI-compliant email client, and an email profile. For more information, see [Emailing a case on page 51](#).

To email a report

1. Generate and view the report. For more information, see [Generating reports on page 107](#).
2. On the **Report Preview** toolbar, click **E-Mail As**.
3. Click the file format for the report. Waveforms must be in PDF format.
4. Depending on the file format that you selected, complete the export options and click **OK**.
The **Save As** window opens.
Navigate to a location where you want to save this report, or accept the default.
5. In the **File name** field, type a name for the report.
The **Save as Type** field lists the document type that you selected.
6. Click **Save**.
Event Review Pro attaches the exported report file when it displays the email window. For more information, see [Emailing a case on page 51](#).
7. Type the email address in the **To** field; type any addresses for copies in the **Cc** fields.
8. Click **Send**.

Working with the System Log

You can view the System Log to see the status of all Event Review Pro activity. You can sort, group, and filter the list of entries. A filter shows only those entries that meet the criteria you select.

This section provides an overview of how to use the System Log .

You can change the default layout of the System Log. This change persists the next time you use the navigation pane or workspace. For more information, see [Resizing panes and workspaces on page 27](#).

To display the system log

1. On the navigation pane, click the **Administration** navigation button.
2. On the **Administration** navigation pane, click **System Log**.

The system log appears.

The System Log lists Event Review Pro activity. For each action, the following information is included:

- The type of activity, for example, Information or Error
- A description of the action
- The date and time when the action occurred
- The user logon name

If a customer support representative asks you to send the System Log, on the **Help** menu, click **Email Customer Support** and send the email message that appears. The System log is attached to it by default.

Each time you start Event Review Pro, the software automatically retains the last 20,000 System Log entries and clears or deletes older entries.

Working with columns

You can customize the columns that appear on the table. You can resize the width of a column. You can also hide or display the columns on the table and the order in which they appear.

- To change a column width, use the mouse to click and hold the column border, and then drag the border to the size you prefer.
- To size a column to fit the width of its current content, use the mouse to click the column border, and then double-click the mouse.
- To change the order in which columns appear, use the mouse to click and hold the column header, and then drag the column to the left or right to a new location.

To hide or display the columns on the table

1. Right-click a column header name to display a shortcut menu.
2. Click **Columns** to display a list of column header names.
3. Click a check box for each column that you want to hide or display on the table.
Column headers with a check mark display on the table.

Grouping and sorting entries

In any grid, such as the Event Log and the System Log, you can arrange entries for your convenience.

Displaying and hiding entries

Right-click any column to display a shortcut menu. The shortcut menu lists all columns available, with check marks next to those that currently display. Click the check boxes to display or hide any of the columns.

Sorting entries

You can click the column header to sort the list of values in ascending (1 to 9, or a to z) or descending (9 to 1, or z to a) order.

Grouping entries

You can group entries so as to sort them by multiple criteria. For example, you can select a primary (major) sort by Type and a secondary (minor) sort by Description, or a primary sort by Device and a secondary one by Event. Use the gray area above the entries to set up groups.

If you want to create further minor sorts, you can then drag more column headings into the gray area to create subgroups.

To group information

1. Click a column header and drag it to the grouping area labeled **Drag a Column Header Here to Group by That Column**.

The column name moves to the grouping area, and the table displays the groups you create.

2. Repeat step 1 as necessary to create more groups. Each additional column name becomes a minor sort below the major sort that you created in step 1. To rearrange any of these column names in the grouping area, drag them to a new position. The table readjusts to the new order of sorts.
3. Click the column header to sort the order in which the grouped information appears.

To ungroup information

1. Click a column header (for example, Type or Event) in the grouping area and drag it below the table header.

If it is the only column head in the grouping area, the table then reverts to its original organization.

If other column headers are still in the grouping area, the table readjusts to sort according to the column heads that are still there.

Filtering entries

Use a filter to limit the number of entries that appear. A filter specifies selection criteria for the entries that appear. Each additional criterion that you specify for the entries further refines the subset of entries that appear. You can combine filtering with sorting and grouping.

Use the upper area above the entries to filter the list of entries. To set up a filter, select a column header and specify a value to use as a criterion. To further refine the list, set up additional criteria.

For example, you might first filter for entries within a certain date range. Then, you might filter these for an institution.

To filter log entries

You can click the arrow to sort the entries in ascending (1 to 9, or a to z) or descending (9 to 1, or z to a) order.

1. Click **Click Here to Add Filter Criteria**.

A row for the criteria condition appears.

2. Click **Choose Field** and click a column name. For example, to limit the display to entries with one specific action, click **Action**.
3. Click **Equal** and click an operator that qualifies the value for the column. You can specify a specific value or a range of values. For a description of each operator that specifies a value, see Comparison Operators, below.

For example, to limit the display to entries with one specific action, click **Equal**.

4. In the next field, click the down arrow and select a value for the column that you selected.
5. (Optional) Click **Add Criteria** to add another criteria condition.
The list of entries updates and displays another criterion row.
6. Click **And** and select a logical operator to evaluate the new criterion in relation to the previous criterion. For a description of each operator used to evaluate criteria, see Logical Operators, below.
7. Repeat steps 2 through 6 for each additional criterion you want to add to the filter.

Event Review Pro compares the entries that meet the last condition to the previous condition. Only the entries that meet the two conditions appear. The process repeats for each preceding criteria condition.

Comparison operators table

You can select a comparison operator to filter a number, a date, or text. For example, to see all entries created after June 1, 2001, click **Greater Than** and type 06/01/2001. You can enter a date using the MM/DD/YYYY format or click the down arrow to choose a date from the calendar. You can filter for specific dates or for a range of dates. You cannot filter for specific times. For the purpose of filtering, the

log uses a time of 12:00 AM for all dates. To select a range of dates, use the **Between** operator.

Use the comparison operators in the following table with these columns: User Name, Type, Action, Description, and Date and Time:

Comparison operator	Lists the entries whose values for the column ...
Equal	...are the same as the selected value.
Not equal	...are not the same as the selected value.
Greater than	...follow or are larger than the selected value.
Less than	...are preceded or are smaller than the selected value.
Greater than or equal to	...are the same, follow, or are larger than the selected value.
Less than or equal to	...are the same, precede, or are smaller than the selected value.
Between	...are between the first and last selected values. Use to specify a range of dates, numbers, or text.
Not between	...are before the first selected value or after the last selected value.
Is null	..do not exist.
Not is null	...exist.
Is empty	...do not appear.
Not is empty	...appear.
Any of	...include the selected items.
None of	...do not include the selected items.

Logical operators table

Select the logical operator from the list of operators that you want to use to compare two criteria. Event Review Pro outlines the criteria and connects the conditions with a line.

Logical operator	Description
And	Displays the entry when the criteria in the first and the second conditions are both met.
And not	Displays the entry when the criteria in the first condition is met

Logical operator	Description
	and the criteria in the second condition is not met.
Or	Displays the entry when either the criteria in the first condition or the criteria in the second condition are met.
Or not	Displays the entry when the criteria in the first condition is met or when the criteria in the second condition is not met.


Removing filters

If you filter the entries and then exit Event Review Pro before removing the filter, the filter settings persist; they appear the next time you use the feature.

You can return a single column to its unfiltered state or remove all the filters.



If you remove a criterion, Event Review Pro removes it and each criterion that appears below it.

To return a single column to its unfiltered state, click **Remove Criteria**  for that criteria condition.


To remove all filters, on the first criterion row, click **Remove Criteria** .

Printing System Log entries

You can print the System Log table in portrait or landscape format.

To print the System Log entries

1. On the **File** menu or toolbar, click **Print**.
2. Click **Table**.
3. In the **Print** window, in the **Select Printer** list, click the printer that you want to use.
4. (Optional) Click **Preferences**.
5. (Optional) In the Print Preferences window, in the **Orientation** area, click an orientation.

 To print the entire width of the system log table, click **Landscape** .

6. Click **OK** to close the window.
7. (Optional) Complete the fields in the **Page** range area.
8. Click **Print** or **OK**.

Using the Import Service

Use the Import Service on the Administration navigation pane to import cases into the Event Review Pro database automatically. When you do this, you can set up inboxes and can use HTTP. See [Configuring the Import Service on page 128](#) for information on your configuration options.

When the Import Service detects that there is a case file with the correct extension, it automatically creates a case using the patient demographic information and the ECG from the case file.

Event Review Pro saves the case to the database and lists it on the Case Records table. The Import Service moves the file to the archive folder in the inbox folder.

If the imported case contains a duplicate ECG, the case appears on the Duplicated ECGs table. For more information, see [Reviewing cases with duplicate ECGs on page 52](#).

For information on how to monitor system usage, see [Working with the System Log on page 120](#).



Case files that are password-protected cannot be imported with the Import Service.

Managing Import Service inboxes

Use the Import Service workspace to view, add, and remove inboxes.

If the Import Service is running when you add or delete inboxes, it detects the change and updates the folders that are monitored. You do not need to restart the Event Review Pro Import Service.

If Import Service is not running when you add or delete inboxes, the changes become effective the next time that Import Service runs.

To view the Import Service inboxes

1. On the navigation pane, click the **Administration** navigation button.
2. On the **Administration** navigation pane, click **Import Service**.

To add an inbox

1. Display the Import Service workspace.
2. In the Inbox folder area, click **Add**.
3. In the Browse for Folder window, use one of the following methods:
 - Navigate to the inbox folder.
 - Navigate to the location for the new folder, click **Make New Folder**, and type a name for the folder.Make sure that the Import Service has read/write permission for this folder.
4. Click **OK**.

The path to the folder in the Inbox folder area now appears.

To remove an inbox

1. Display the Import Service workspace.
2. In the **Inbox** folders area, click the inbox folder location for the inbox that you want to remove.
3. Click **Remove**.
4. Repeat steps 2 and 3 to remove additional inboxes.

Configuring the Import Service

The default installation does not set up the Import Service to start automatically. If you need to import cases automatically, you must set up the Import Service from Windows.

To set up the Import Service from Windows

1. Complete the typical Event Review Pro installation.
2. On the Windows **Start** menu, click **Control Panel**.
3. Double-click **Administrative Tools**.
4. Double-click **Services**.
5. Right-click **Event Review Import Service**.

6. Click **Properties**.
7. In the **Startup type** field on the General tab, click **Automatic** and then **Apply**.
8. Using the information that your IT technician provides, complete the following setup:
 - a. In the **Log on as** field on the Log On tab, click **This account**. The account that you are currently using must have access permission to the Event Review Pro database.
 - b. Complete the account and password information.
 - c. Unless you have a specific reason, do not change the entries on the **Recovery** or **Dependencies** tabs.
 - d. Click **Apply**.
9. Click the **Overview** tab and click **Start**.
10. Click **OK** to close the Properties window.
11. Close the **Service** window.

To start the Import Service manually

1. On the Windows **Start** menu, click **Control Panel**.
2. Double-click **Administrative Tools**.
3. Double-click **Services**.
4. On the **Service** window, click the **Standard** tab.
5. Right-click **Event Review Import Service**.
6. Click **Start**.
7. Close the **Services** window.




The System Log lists any Event Review Import Service activity. For more information, see [Working with the System Log on page 120](#).

Managing Import Service archives

After a file is imported successfully to the database, Import Service moves the file to the archive folder.


Import Service creates the archive folder as a subfolder of the inbox folder. Each inbox folder has its own archive folder. Under the archive folder, Import Service creates separate subfolders to keep the number of files manageable. The name of the daily archive folder is *yyyymmdd*, where *yyyy* is the year, *mm* is the month, and *dd* is the day of the month.

 Import Service does not delete the files or the archive directories. It is the responsibility of the user to periodically delete old archived files and folders.

Using the HTTP Import Service

The Import Service supports the import of case data using HTTP. After you set up HTTP, you can transfer case data from HeartStart Data Messenger. For more information, see the *HeartStart Data Management implementation guide*.

Use the Import Service workspace to enable, change, or disable the working folder for the HTTP import service.

 You must run the Import Service under a named administrator (not the default local service) for HTTP import to work properly. See [Configuring the Import Service on page 128](#).

To enable the HTTP import service

1. Display the Import Service workspace. For more information, see [Managing Import Service inboxes on page 127](#).
2. In the HTTP Import Service area, click **Set Working Folder**.
The **Browse for Folder** window opens.
3. Use one of the following methods to create a working folder:
 - Navigate to the working folder.
 - Navigate to the location for the new folder, click **Make New Folder**, and type a name for the folder.
Make sure that the Import Service has read/write permission for this folder.
4. Click **OK**.
The path to the folder in the HTTP Import Service area appears.

To change the working folder

1. Display the Import Service workspace. For more information, see [Managing Import Service inboxes on page 127](#).
2. In the HTTP Import Service area, click **Set Working Folder**.
3. Use one of the following methods to change the working folder:
 - Navigate to a different working folder.
 - Navigate to the location for the new folder, click **Make New Folder**, and type a name for the folder.

4. Click **OK**.

The path to the folder in the HTTP Import Service area appears.

To disable the HTTP Import Service

1. Display the Import Service workspace. For more information, see [Managing Import Service inboxes on page 127](#).
2. In the HTTP Import Service area, click **Clear Working Folder**.

The path to the folder in the HTTP Import Service area disappears.

Customer support

Philips strives to provide you with excellent customer service and technical support. Software updates for the integrated applications are available from the application Help menu. From the **Help** menu, click **Check for Updates**.

Customer support is available through email, Internet, and telephone.

Email product support is available in English only at: eventreview.support@philips.com. In addition, if a customer support representative asks you to send the System Log, on the **Help** menu, click **Email with System Log** and send the email message that appears. The System log is attached to it by default.

Internet product support is available at the following address:

<http://www.philips.com/DataManagementSupport>

For telephone assistance outside the United States, please call your sales representative or local response center. See the following tables for contact information. You can also navigate to technical support telephone numbers for data management products at this address:

http://www.healthcare.philips.com/main/services/response_center

Region	Address	Telephone number
United States	Philips Healthcare 3000 Minuteman Road Andover, Massachusetts 01810-1099	+1(800) 263-3342
Canada	Philips Healthcare, a Division of Philips Electronics Ltd. 281 Hillmount Road Markham, Ontario, Canada L6C 2S3	+1(800) 291-6743
Authorized EU Representative Europe, Middle East, and Africa	Philips Medizin Systeme Boeblingen GmbH Cardiac and Monitoring Systems Hewlett-Packard Strasse 2	(+49) 7031 463-2254

Region	Address	Telephone number
	71034 Boeblingen, Germany	
Latin America	Philips Medical Systems Ltda. Rua Verbo Divino - 1400 - 7o Andar0 04719-002 - Sao Paulo, Brazil	0800 7017789
Asia Pacific	Philips Electronics Hong Kong Ltd. 6/F, Core Building 1 1 Science Park East Avenue Hong Kong Science Park Shatin. New Territories, Hong Kong	+852 2821 5888

You can also navigate to technical support telephone numbers for data management products at the above addresses.

For telephone support in English only, you can call the following numbers between 9:00 AM and 5:00 PM, Pacific Time:

- (800) 263-3342, inside the United States
- +1.978.659-3332, outside the United States

Philips Response Center telephone numbers

The following table contains telephone numbers for Philips Response Centers around the world.

Area	Country	Telephone number
North America	Canada	800-323-2280
	United States of America	800-722-9377

Area	Country	Telephone number	
Europe	European International Sales	41 22 354 6464	
	Austria	01 60 101 820	
	Belgium	French	02 525 68 80
		Dutch	02 525 68 81
	Finland	09- 615 80 400	
	France	0810 835 624	
	Germany	0180 3333 544	
	Italy	800 232100	
	Luxembourg	+32 2 525 68 80	
	Netherlands	040 27 85600	
	Portugal	800 201766	
	Spain	900 180612	
		902 304050	
	Sweden	08-59 85 2530	
	Switzerland	German	0800 80 3000
French		0800 80 3001	
United Kingdom	0870 532 9741		

Area	Country	Telephone number	
Asia and Pacific	Australia	1800 251 400	
	China	Beijing	800 810 0038
		Hong Kong	852 2821 5888
		Macau	0800 923
	India	18004256788	
	Indonesia	021 794 7542	
	Japan	0120 381 557	
	Korea		080 372 7777
		Seoul	02 3445 9010
	Malaysia	1800 866 188	
	New Zealand	0800 251 400	
	Philippines	02 845 7875	
	Singapore	1800 PHILIPS (1800-744-5477)	
	South Africa	011 471 6000	
	Thailand	02 614 3559	
Taiwan	0800 005 616		

Comments or suggestions?

Please send your feedback and suggestions to: eventreview.support@philips.com

Supported help

Customer support technicians provide help for the following:

- Explaining the proper use of application features and answering your questions about how the application works
- Explaining the proper installation and maintenance of the application
- Assisting you in selecting and configuring card readers

Unsupported help

Customer support technicians do not provide help for the following:

- Interpreting ECG or medical data. Please call your medical director or clinical specialist.
- Repairing hardware. The support technicians can help you determine if you have a hardware problem, but they cannot help you fix problems that are not related to the Philips HeartStart application software.
- Troubleshooting defibrillators. Instead, call Philips Customer Support and ask for defibrillator support.
- Troubleshooting non-Philips products.

Helping us help you

You can help our technicians give you good support by following these steps:

1. Call from a phone near your computer.
2. Start Event Review Pro.
3. Have the following information:
 - Windows version.
 - The application version number. This is available from the Help menu. Click the About option.
 - A written copy of the error message text.
 - The activity and task you did when the error occurred.

Working with defibrillators

Most of the defibrillators that Event Review Pro supports store data on data cards (also called memory cards). Depending on the defibrillator, you need a CompactFlash or SD data card reader to transfer the data.

The HeartStart FR3 defibrillator can transfer information using a data card or a *Bluetooth* connection.

The HeartStart HS1 family of defibrillators and FRx defibrillator use an infrared connection to transfer the information.

The HeartStart MRx Monitor/Defibrillator can transfer information using the following methods:

- A data card to transfer the information
- Defibrillators with version 9.0 or later with the *Bluetooth* option to transfer the information through a *Bluetooth* connection



If you need help choosing a card reader, see [Selecting accessories for data transfer on page 138](#). The quick reference has illustrations of the various data cards and the associated defibrillators. It can help you recognize the data cards and select the appropriate card reader.

Supported defibrillators

You can download product owners' manuals for supported defibrillators from the Philips support Web site;

To download products and services materials

1. Navigate to this site: <http://www.philips.com/productdocs>
2. Click **Resuscitation/Defibrillators products**.

- Click the defibrillator to display a page that lists the product user and service materials.

The following table identifies the method used to read defibrillator data. To understand how to use the method, click the appropriate link in the **Instructions** column.

Data source	Method used to read data	Instructions
HeartStart FR3	Secure digital (SD) card or <i>Bluetooth</i> connection	Reading FR3 series cards on page 140 Downloading an ECG using an FR3 Bluetooth transmission on page 37
HeartStart MRx	Card reader or <i>Bluetooth</i> connection	Reading HeartStart MRx cards on page 144 Downloading an ECG using an MRx Bluetooth transmission on page 38
HeartStart FRx	Infrared (IrDA) connection	Downloading an ECG from an infrared connection on page 39
HeartStart HS1	Infrared (IrDA) connection	Downloading an ECG from an infrared connection on page 39

Selecting accessories for data transfer

Before attempting data transfer from defibrillators and devices to Event Review Pro, be sure that you have all the accessories you need.

Defibrillators and devices vary as to whether they use data cards and card readers, infrared connections, and/or *Bluetooth* connections to transfer data.

Choosing data transfer adapters and readers


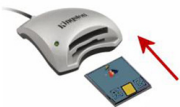





If your computer did not come installed with the appropriate software or adapter for the transfer method, you can use an adapter or reader. Use the following table to identify the requirements for the adapter or drive that you need.

Defibrillator	Data transfer method
HeartStart FR2	CompactFlash data card

Defibrillator	Data transfer method
HeartStart FR3	Custom Philips version of the SD data card <i>Bluetooth</i> data transfer
HeartStart FRx	Infrared data transfer
HeartStart HS1	Infrared data transfer
HeartStart MRx	CompactFlash data card <i>Bluetooth</i> data transfer batch LAN data transfer
HeartStart XL	ATA flash data card

Quick reference to data cards

The following table shows the defibrillators with their associated data cards.

HeartStart FR2 		FR2 defibrillators use the CompactFlash cards with the CompactFlash card reader, or an ATA flash-card reader with an adapter.
HeartStart FR3 		FR3 defibrillators use Philips secure digital SD data cards, with the SD data-card reader, or an ATA flash-card reader with an adapter.
HeartStart MRx 		MRx defibrillators use CompactFlash cards with the CompactFlash card reader, or an ATA flash-card reader with an adapter.
HeartStart XL 		XL defibrillators use ATA flash cards with the ATA flash-card reader.

Choosing *Bluetooth* accessories

HeartStart MRx and FR3 defibrillators support data transfer using *Bluetooth* wireless technology.

Event Review Pro was tested with the following Windows systems and *Bluetooth* stacks:

- Windows XP or Windows 7: Microsoft *Bluetooth* stack (for MRx)
- Windows XP or Windows 7: Widcomm *Bluetooth* stack (for FR3 or MRx)

If your computer did not come with the installed *Bluetooth* option, you must use a *Bluetooth* adapter that supports these stacks.

Reading FR3 series cards

The FR3 series defibrillators store information on a SD card. SD card readers are available through your Philips sales representative or local distributor.

To read a data card

1. One end of the data card has a series of contacts. Plug this end firmly into the card reader.
2. Check in Windows Explorer to see that the card was recognized. When Windows recognizes the card, it assigns a drive letter to the card.
At this point, Event Review Pro can read the data on the card.
3. If the card does not appear in Windows Explorer, remove the card from the reader, and then insert it into the card reader again.

Retrieving FR3 data using *Bluetooth* transmission

You can download data from the FR3 to the computer through a *Bluetooth* transmission. Any data currently contained in the FR3 is transferred.

The *Bluetooth* wizard guides you through the steps necessary to accomplish the transfer. From an open case, click **Attach ECG**, click *Bluetooth*, click **HeartStart FR3**, and then click **Next**.

Follow the instructions from this point.

For further information, see [Downloading an ECG using an FR3 Bluetooth transmission on page 37](#).

Setting up Bluetooth transmissions for the MRx and FR3

You can transmit HeartStart MRx and HeartStart FR3 data to Event Review Pro through *Bluetooth*® wireless transmission.



You must have a supported adapter for Bluetooth transmission. Most Microsoft®-compatible Bluetooth adapters, such as IO Gear GBU421 (or in Japan, the Logitech Bluetooth 3.0 adapter (LBT-13R35861A) are suitable. You can order an adapter from your Philips representative.

To transfer the information from the defibrillator to Event Review Pro, you must set up the *Bluetooth* transmission on the defibrillator and the receiving computer that has Event Review Pro installed.

Consider the following:

- If your computer is not equipped with *Bluetooth*, it requires a *Bluetooth* dongle with a USB connector. The dongle connects the computer and the defibrillator.
- Even if your computer is equipped with *Bluetooth*, you may need other adjustments. For example, Panasonic Toughbooks have a Toshiba *Bluetooth* driver, not a Microsoft driver. This works for the MRx; for the FR3, add a *Bluetooth* dongle that supports the Microsoft or Widcomm driver.

Event summaries and (for MRx) 12-lead reports might contain information that could be considered Patient Healthcare Information (PHI) or patient-identifiable data. Handle the information in accordance with HIPAA or your local patient privacy requirements.

Bluetooth option prerequisites for the MRx and FR3

The prerequisites for *Bluetooth* use are slightly different for the MRx and FR3 defibrillators.

MRx

The *Bluetooth* card on the HeartStart MRx is inside of the defibrillator. Confirm that you have a *Bluetooth* card installed and that the Event Summary *Bluetooth* transmission feature is enabled.

For detailed information on setting up the MRx for *Bluetooth* transmission, see the most current release of the following documentation:

- HeartStart MRx Instructions for Use Addendum
- HeartStart MRx M3535A, M3536A Instructions for Use
- HeartStart MRx Data Transmission Implementation Guide

FR3

Confirm that you have a Bluetooth card installed in the defibrillator.

For detailed information on setting up the FR3 for *Bluetooth* transmission, see the most current release of the following documentation:

- HeartStart FR3 Instructions for Administrators

Pairing and testing the *Bluetooth* option with the computer: MRx and FR3

MRx

Use this information to pair and test the *Bluetooth* option for HeartStart MRx.



Make sure that the application computer is turned on and that the *Bluetooth* software and adapter are installed and visible or discoverable.

To pair the HeartStart MRx *Bluetooth* option with the computer

1. On the HeartStart MRx main menu, click **Bluetooth Devices**.
2. Click **Add Device**.

The defibrillator searches for your computer.

If your computer is not listed after the search, your computer *Bluetooth* is not enabled or set up correctly. Refer to the *Bluetooth* documentation to troubleshoot the setup.

3. Click the name of your computer and then press **Enter**.
4. Type a passkey for the *Bluetooth* option, using the number menu.
The passkey is a user-defined character sequence, such as 000, or 1234. For information on this sequence, see the documents listed in the [Bluetooth option prerequisites for the MRx and FR3 on page 141](#).
5. At the top of the number menu, click **Done**.
6. On your computer, immediately watch for a pop-up message.
7. Click the pop-up message.
8. Type the same *Bluetooth* passkey that you typed on the defibrillator in step 4. The defibrillator and the computer now have a *Bluetooth* connection, and you are ready to test it.
9. Click **Next** and then **Finish**.
10. On the main menu of the MRx, click **File Transfer** and then **Enter**.
HeartStart MRx displays the "transmission test passed" message.
11. Click **Acknowledge**.
12. Scroll to **Exit** and click **Enter**.
13. On the main menu of the MRx, click **Other, Data Management**, and then **Acknowledge**.
14. After changing into Data Management mode, the MRx displays a list of summary case data with the most recent incident at the top. Select the incident that you want to transmit.
15. Click **Send** and then **All Event Data**.
16. On your computer, in the **Attach ECG** dialog box, watch for a display that includes the MRx folder, and open it to see the data that has been transferred.

FR3

Pairing the Bluetooth option for the FR3 occurs automatically. Use the following procedure to verify that this has occurred.

To verify that the FR3 Bluetooth option has paired with the computer:

1. When the FR3 is on and you have heard the voice acknowledgment, put the unit into Administration mode.
2. Press the option button and select wireless data transfer.

3. On your computer, in the **Attach ECG** dialog box, Event Review Pro discovers the PIN number and displays a list of the patient cases currently on the FR3.

Reading HeartStart MRx cards

The HeartStart MRx defibrillator stores information in internal memory.

To open the information in Event Review Pro, transfer the information from the internal memory to an external CompactFlash card and then download it to Event Review Pro.

Follow the instructions in the HeartStart MRx Monitor/Defibrillator *Instructions for use* guide. After you transfer the information to a CompactFlash card, you can read the card in a card reader.



CompactFlash cards and card readers are available through your Philips sales representative or local distributor.

To read a data card

1. Remove the data card from the HeartStart MRx data card tray.
2. One end of the data card has a series of perforations. Insert this end firmly into the card reader.
3. Check in Windows Explorer to see that the card was recognized. When Windows recognizes the card, it assigns a drive letter to the card.
Event Review Pro can read the data on the card.
4. If the card does not appear in Windows Explorer, remove the card from the card reader, and then insert it into the card reader again.

Sending the MRx Bluetooth transmission

You can download a case from HeartStart MRx to the computer through a *Bluetooth* transmission.

Before downloading, make sure that the MRx has been paired with your computer. See [Pairing and testing the Bluetooth option with the computer: MRx and FR3 on page 142](#) for information on pairing and testing the connection.

To make a Bluetooth transmission

1. On the HeartStart MRx, place the defibrillator in the Data Management mode. If necessary, use the MRx documents that are listed in [Bluetooth option prerequisites for the MRx and FR3 on page 141](#).
2. Click **Next Item** or **Prev Item** to navigate to the incident of your choice.
3. Scroll to **Send**.
4. Press **Send**.
5. Select the type of information that you want to transmit to Event Review Pro and click **Enter**.
6. The first time you send a *Bluetooth* transmission, a popup message appears, including an option to **Always Allow Authorization**. If you want to do this, check the box; if you do not, the checkbox will appear every time you send a transmission..
7. After the HeartStart MRx transmits the incident to the computer, you can attach the ECG from the HeartStart MRx in Event Review Pro.

For more information, see [Downloading an ECG using an MRx Bluetooth transmission on page 38](#).

Reading FR2 series cards

The FR2 series defibrillators store information on a CompactFlash card. Card readers are available through your Philips sales representative or local distributor.



The FR2 series of defibrillators can use adult defibrillation pads or infant/child reduced-energy defibrillation pads. When a user connects infant/child reduced-energy defibrillation pads to the FR2, the FR2 delivers a reduced-energy defibrillation shock of 50 joules (nominal).

On the Attach ECG window, you must indicate to the application whether adult or infant/child pads were used on the patient. To do this, click the FR2 Adult or FR2 Pediatric option from a shortcut menu. For more information, see [Adding ECGs on page 34](#).

If you have any doubt about the type of pads that the responder used, contact the responder.

To read a data card

1. One end of the data card has a series of perforations. Plug this end firmly into

the card reader/

2. Check in Windows Explorer to see that the card was recognized. When Windows recognizes the card, it assigns a drive letter to the card.
Event Review Pro can read the data on the card.
3. If the card does not appear in Windows Explorer, remove the card from the reader, and then insert it into the card reader again.

Reading XL cards

The XL defibrillator stores information on an ATA flash card. Read the documentation for the defibrillator to learn how to remove the card from the defibrillator.

You can read these cards in an ATA flash card reader or a PC Card (PCMCIA) reader. If the computer is not equipped with an adapter, you can add one as an external or internal card reader. The external card reader plugs into the desktop USB port.

To read a data card

1. One end of the card has a series of perforations. Insert this end firmly into the card reader.
2. Check in Windows Explorer to see that the card was recognized. When Windows recognizes the card, it assigns a drive letter to the card.

The Philips Data Management Solutions application and Event Review Pro can read the data on the card.

3. If the card does not appear in Windows Explorer, remove the card from the card reader, and then insert it into the card reader again.

Transferring XL+ data

The HeartStart XL+ defibrillator uses a USB data drive to transfer event summary data. The USB data drive connects to a USB port on the computer running Event Review Pro.

You can export data from the XL+ defibrillator to the XL+ USB data drive using a menu selection in the Data Management mode. Once the data is exported, you can transfer it to the computer running Event Review Pro. See [Adding an ECG from a USB port on page 37](#).


To export event-summary data from the XL+ defibrillator to the USB data drive

1. Locate the USB data port on top of the XL+ defibrillator and insert the XL+ USB data drive into it. For more information, see the *HeartStart XL+ instructions for use*.
2. Confirm that you are in **Data Management** mode.
3. Press the **Menu Select** button.
The Internal Memory menu appears.
4. Use the navigation buttons to select the **Export All** option.
5. Press the **Menu Select** button to start the task.
The data is exported to the USB data drive.

Determining the HS1 and FRx case date and time

The HS1 family of defibrillators and FRx defibrillators do not have a real-time clock that keeps the defibrillator date and time. Instead, the defibrillator has a continuously running timer that starts when you insert the battery. Event Review Pro uses this timer and the computer date and time to calculate case and event time. The timer resets when you remove the battery. Therefore, do not remove the battery until after you have downloaded the patient data to the application.

Each time that Event Review Pro communicates with an HS1 and FRx defibrillator through IrDA, it uses the date/time of the computer clock and the timer of the HS1/FRx defibrillator to synchronize the date/time of all events and data since the battery was last inserted. This is called the "time sync point." If Event Review Pro does not find a time sync point, you need to set the case start time.

 The time sync point might not match the date and time of the computer by 15 to 30 seconds.

Using infrared connections for the HS1 and FRx

The HS1 family of defibrillators and the FRx defibrillator use an infrared (IrDA) communications port to transfer information. For information on how to determine the case date and time, see [Determining the HS1 and FRx case date and time on page 147](#).

The infrared port is located on the lower right side of the defibrillator. To transfer information, the receiving computer must be set up for infrared communication.

Some computers have infrared ports. If the computer does not have an infrared port or the port does not work with the defibrillator, you can add an infrared adapter.

An infrared adapter that connects to a USB port, such as the ActiSYS IR4002US (ACT-IR4002US if you purchase it from Philips) is the recommended adapter.

Setting up an infrared adapter

On Windows XP installations, Event Review Pro was tested with ACTiSYS adapters. For best results, use the IR4002US (ACT-IR4002US if you purchase it from Philips) and follow the manufacturer's instructions. More information is available at the following Web site:

<http://www.actisys.com>

To set up the infrared adapter

1. Read the instructions that came with the infrared adapter.
2. Windows XP only: Use the CD that comes with the infrared adapter to run the driver setup program on the Windows operating system.
3. The New Hardware Wizard guides you through setting up the adapter.

Setting up the infrared connection

Downloading an ECG and device history information from an HS1 or FRx defibrillator requires an infrared connection between Event Review Pro and the defibrillator.

To transfer information, you must put the HS1 or FRx defibrillators into the Administration mode.

The following discussion explains how to set up each defibrillator to transfer information. The process is similar for the HS1 family of defibrillators and the FRx defibrillator. However, there are enough differences to warrant separate sections.

To set up the IrDA connection between an FRx defibrillator and Event Review Pro

If the FRx does not receive a transmission within 3 minutes of starting Administration mode, it cancels the mode. At that point, you must start over.

1. In Event Review Pro, complete one of the following procedures:
 - a. Use the Case Wizard to create a case, up to the point where you see the Attach ECG dialog box.
 - b. Create a case for the infrared transmission, and then (on the toolbar) click **Attach ECG file**.
 - c. Open an existing case for the transmission, and then (on the toolbar) click **Attach ECG file**.
2. Remove the pads connector.
3. If you haven't done so already, insert the battery into the defibrillator.

You will hear voice prompts to plug in pads connectors. Disregard these prompts in this instance. These messages end once you place the defibrillator in Administration mode. .
4. On the defibrillator, press and hold the blue **Information** button, and wait for three tones.
5. Release the blue Information button.

The voice message announces Administration. The defibrillator enters data transfer mode.
6. Position the defibrillator to communicate with the computer.

The defibrillator and computer should be between 4 and 24 inches apart and aligned with each other, with an unobstructed path between them. For help aligning the defibrillator and the computer, click the IrDA icon on the Attach ECG window and follow screen instructions.

You can now send and receive information between the defibrillator and the computer. The device announces transferring data.
7. When sending and receiving information is complete, plug the pads connector in again. You can now turn off the defibrillator.

For more information, refer to the defibrillator documentation.

To set up the IrDA connection between an HS1 defibrillator and Event Review Pro

If the HS1 does not receive a transmission within 3 minutes of starting Administration mode, it cancels the mode. At that point, you must start over

1. In Event Review Pro, complete one of the following procedures:
 - a. Use the Case Wizard to create a case, up to the point where you see the Attach ECG dialog box.
 - b. Create a case for the infrared transmission, and then (on the toolbar) click **Attach ECG File**.
 - c. Open an existing case for the infrared transmission, and then (on the toolbar) click **Attach ECG File**.
2. Remove the pads cartridge. Locate the latch at the top, and then slide the latch to the right to release the pads cartridge.
3. If you haven't done so already, insert the battery into the defibrillator.
4. You will hear voice prompts to plug in the pads cartridge. Disregard these prompts in this instance. These messages will be silenced once you place your defibrillator in Administration mode. After the battery is inserted, the defibrillator automatically turns on.
5. On the defibrillator, press and hold the blue **Information** button, and then wait for three tones.
6. Release the blue button.

A voice message announces Administration.
7. Briefly press the blue **Information** button again.

The voice message says Mode 1.
8. Position the defibrillator to communicate with the computer.

The defibrillator and computer should be at least 12 inches apart and aligned with each other. For help aligning the defibrillator and the computer, click the IrDA icon on the Attach ECG window and follow screen instructions.

You can now send and receive information between the defibrillator and the computer. The device announces sending data.
9. When sending and receiving information is complete, reinstall the pads cartridge. You can now turn the defibrillator off.

For more information, refer to the defibrillator documentation.

Understanding voice and system messages

Defibrillators and Windows use messages to announce the current state.

Defibrillator voice messages

Philips designed the defibrillators to ensure that they are always ready for use in an emergency.



Remember to insert the HS1 pads cartridge or plug in the FRx pads connector when you finish transferring information. This will ensure that the defibrillator is ready for use on the next patient.

At several points in the Event Review Pro instruction, you are instructed to ignore the defibrillator voice message that says "Insert the pads cartridge (for HS1)" or "No connector installed. Plug in pads connector (for FRx)." The message is there to ensure that you are aware of the defibrillator's current state. If you remove the HS1 pads cartridge or the FRx pads connector, the defibrillator is not available for use on the patient. When the defibrillator enters the Administration mode, the message stops.



You can put the FR3 into Admin mode without removing pads.

System messages

If Windows detects an active wireless device, it starts the Windows Wireless Link. Event Review Pro does not use this application and cannot disable it. As a result, the Windows Wireless Link might display messages. These messages do not apply to your current task. You can ignore them.

Emailing device history data

If you ask for assistance from Philips Customer Support, you may be asked to retrieve the device history data from a specific defibrillator. Use the following steps (or as directed by Customer Support) to retrieve and email this data.

To email device history data to Customer Support

1. From the **Help** menu, click **Email Device History**.

The Email Device History dialog box appears.

2. Set up the device that you are sending data from:
 - a. If you are sending data from an FRx or HS1 device, put the device into Administration mode. For details on setting up the connection, see [Setting up the infrared connection on page 148](#).
 - b. If you are sending data from a data card from an FR2 or FR3 device, follow the instructions from Customer Support to put the device history data file on the card, and then insert the card into a card reader attached to your computer.
 - c. If you are sending data from a *Bluetooth* connection on an FR3, follow the instructions from Customer Support to put the device history data file on the card; then turn the defibrillator on, wait for the voice prompts, put the device into Administration mode, press the option button, and select Wireless Data Transfer.
3. From the list of devices, click the one that you want to send data for and click **Finish**.

An email message box appears, with the device history data appearing as an attachment.

4. Type any comments you want to add to the message.
5. Click **Send**.

Managing the database

Event Review Pro stores the application database in Microsoft SQL Server 2008 R2 Express Edition. This is a free version that Philips installs by default. It limits the database size to 10GB. The number of cases that can be stored in a database of this size varies depending on the length of captured ECGs and whether or not the file includes audio.

The versions sold by Microsoft are limited only by the size of the hard disk.

The Event Review Pro database is low-maintenance, compared to many other databases. Nonetheless, attention to routine maintenance tasks will help to ensure database integrity.



It is the database administrator's responsibility to set up the appropriate database management features and to check that they execute successfully.

Back up the database on a regular basis (preferably every day) to tape or some other medium. Disaster recovery experts recommend that you store the backup tapes somewhere safe, so that disasters such as fire or theft cannot harm them. Without a recent backup, you have no chance of recovery after a catastrophe, such as a disk failure or fire. Use your organization's backup and recovery policy to set up this procedure.

For information on migrating data from an older version of the application, see [Migrating cases from previous versions on page 157](#).

Using the Event Review Pro database on a remote database server

Event Review Pro supports use of a remote database that is installed on another machine. Multiple instances of Event Review Pro can share the remote database.

Philips supports implementing this feature with the following software:

- SQL Server 2008 R2
- SQL Server 2008

- SQL Server 2005

These can run on the following:

- Windows Server 2008 R2.
- Windows Server 2008
- Windows Server 2003 R2
- Windows Server 2003



Not all versions of SQL Server run on all versions of Windows servers. Check the SQL Server specifications before proceeding.

You must install the same version of Event Review Pro on each client machine connected to a shared database.



Philips recommends that Information Technology (IT) personnel install and configure the remote database server.

Use this feature with care. Before implementing the application on a shared database server, keep the following in mind:

- The application saves only the changes made by the most recent user to save the case. For example, when more than one user opens the same case, makes changes, and saves the case, the application saves the changes made by the *last user* to save the case.
- Event Review Pro installed on a client machine cannot run when the remote database is unavailable, such as when the server or network is shut down or disconnected.

Changing the database server

You can change the Event Review Pro database server in the following ways:

To change the database server from the Tools menu

1. On the **Tools** menu, click **Database Server**.

The Configure SQL Server for Event Review Pro 4.2 window opens, showing the current database name.

2. In the **SQL Server name** box, use the down arrow to display a list of all available SQL Servers on the network.
3. Click the SQL server name that you want.

4. To verify the connection, click **Test**.

A message from the SQL server software appears, noting success or failure. If Event Review Pro displays a failure message, refer to your SQL Server documentation or your local database administrator for instructions.

5. If you see a success message, click **Save**.

A message appears stating that you must restart Event Review Pro to use the new database server.

To verify the database server on the About Event Review Pro window

1. Start Event Review Pro.
2. On the **Help** menu, click **About Event Review Pro**.
3. On the **Server** line, confirm that Event Review Pro is connected to the correct server database.

Using Microsoft SQL

You can purchase the SQL Server software to set up a remote database environment for Event Review Pro.



Philips recommends that IT administrators who use domain controllers create a user group for the Event Review Pro users.

Setting up SQL Server as a remote database server involves the following tasks:

1. Install the SQL Server Software.
2. Set up Event Review Pro on the client machines.
3. Create the Event Review Pro database.
4. Set up access control on the database.
5. Troubleshoot the client machine connections to the database.

1 - Install the SQL Server software

Install Microsoft SQL Server on the server machine operating system.

2 - Set up Event Review Pro on the client machines

1. Navigate to the installation folder for Event Review Pro on the client machine.
2. Double-click **SqlInstaller.exe**.

The Install SQL database for Event Review Pro 4.2 window opens.

3. In the **SQL Server name** box, use the down arrow to display a list of all available SQL Servers on the network.
4. Click the SQL server name that you want.
5. To verify the database connection, click **Test**.
A message from the SQL server software appears, noting success or failure. If Event Review Pro displays a failure message, refer to your SQL Server documentation or your local database administrator for instructions .
6. If you see a success message, in the **Database Name** box, type the database name that you want.
7. Click **Install**.
Event Review Pro is now installed in that database on the server that you selected.
8. Click **Save**.
A message appears stating that you must restart Event Review Pro to use the new database.
9. Start Event Review Pro.
The software connects automatically to the remote database.
10. On the **Help** menu, click **About** to verify the connection.
Event Review Pro lists the database in use with the remote machine name.
11. If the test fails, continue with step 3 - Troubleshoot the client machine connections to the database.

3 - Create the Event Review Pro database

The following must exist before you create the Event Review Pro database on the server database:

- A domain account.
- Log on to the domain using the workstation from which you will run the installation software.
- Administrative privileges on the Event Review Pro machine used to run the installation software. You might need to add your domain account to the local Administrators group on the workstation.

4 - Set up access control on the database

Philips recommends that you create a user group for Event Review Pro users to access the database. This is easy to administer: for example, create a group called "Phillips Data Management" and add users to that group. However you should set up

security and access according to the security policies of your organization.



It is the SQL Server database administrator's responsibility to set up the network security to recognize the user name and password on the domain and workgroup.

Your SQL Server database administrator must complete the following:

- Add the user group to the logons for the database server
- Grant read/write permission to the database created earlier in these procedures to the user group

5 - Troubleshoot the client machine connections to the database

Consider the following:

- If Event Review Pro displays a message, refer to your SQL Server documentation for a description of the message.
- Do you have access privileges on the remote machine?

Event Review Pro uses Windows Authentication for database access.

Migrating cases from previous versions

You can migrate cases from Event Review Pro 3.5, 4.0, or 4.1 to the Event Review Pro 4.2 database. If the source database is on another machine, see your IT professional for assistance.

Migrating the data involves the following steps:

1. You identify the database from which the cases are to be migrated—the *source* database.
2. Event Review Pro creates a list of all the cases from the source database and compares them (using the case ID) to the cases in the current Event Review Pro 4.2 database (the *destination* database). Cases that are not already in the destination database are marked for migration.
3. You migrate the marked cases.

Use the following procedures to complete the migration.

To select and identify the source database

1. In the main workspace, click the **Administration** button.
2. In the **Administration** navigation pane, click one of the following options:
 - Migrate Records from Event Review Pro 3.5
 - Migrate Records from Event Review Pro 4.0
 - Migrate Records from Event Review Pro 4.1

A database selection dialog box appears, displaying only those computers with the expected, well-known instance and database names, as shown in this table:

Event Review Pro version	SQL Server instance name	Database name
3.5	HeartStart	ERPro35
4.0	Philips	ER40
4.1	Philips	ER41

3. Click the appropriate source database.
4. (Recommended) To test the accessibility of the source database, click **Test**.
A message appears reporting the success or failure of the test.
5. Click **OK** to continue the migration process.

To view and confirm the list of cases to be migrated

1. Once you have selected the source database, Event Review Pro connects to this database and displays a list of all the cases. The ones to be migrated have the + sign.
2. If a case with the same case ID already exists in the destination database, the case will not be migrated. If you still want to migrate such a case, cancel the migration, rename or delete the case in the destination database, and then restart the migration process.

To migrate the cases

1. Click **OK**.
The data migration begins. This takes some time, particularly if you have many cases.
2. If you need to cancel migration, click **Cancel**. You can restart the migration process at another time. Because any already migrated cases will have the same case ID as those in the source database, those cases will not be migrated again.

The migration status of each case is logged in the System Log, including any errors.

The successfully migrated cases now appear in the All Records table.

CPC and OPC

The Cerebral Performance Categories (CPC) and Overall Performance Categories (OPC) are two scores that are used when reporting cardiac arrest. The Follow-up tab uses these categories.

EMS Edition only: The two scores are the recommended guidelines for uniform reporting of information from out-of-hospital cardiac arrest in the Utstein-style templates. (Cummins RO, Chamberlain DA, Abramson NS, et al. AHA Medical/Scientific Statement - Recommended guidelines for uniform reporting of information from out-of-hospital cardiac arrest: the Utstein style. *Circulation* 1991;84:960-975.)

This section explains the scores for each category.

Cerebral performance categories table

Score	Description
1	Good cerebral performance: Conscious. Alert, able to work and lead a normal life. May have minor psychological or neurological deficits (mild dysphasia, nonincapacitating hemiparesis, or minor cranial nerve abnormalities).
2	Moderate cerebral disability: Conscious. Sufficient cerebral function for part-time work in sheltered environment or independent activities of daily life (dressing, traveling by public transportation, and preparing food). May have hemiplegia, seizures, ataxia, dysarthria, dysphasia, or permanent memory or mental changes.
3	Severe cerebral disability: Conscious. Dependent on others for daily support because of impaired brain function (in an institute or at home with exceptional family effort). At least limited cognition; includes a wide range of cerebral abnormalities from ambulatory with severe memory disturbance or dementia precluding independent existence to paralytic and able to communicate only with eyes, as in locked-in syndrome.
4	Coma, vegetative state: Not conscious. Unaware of surroundings, no cognition. No verbal or psychological interactions with environment.
5	Death: Certified brain dead by traditional criteria.

Overall performance categories table

Score	Description
1	Good overall performance: Healthy, alert, capable of normal life. Good cerebral performance (CPC1) plus no or only mild functional disability from noncerebral organ system abnormalities.
2	Moderate overall disability: Conscious. Moderate cerebral disability alone (CPC2) or moderate disability from noncerebral system dysfunction alone or both. Performs independent activities of daily life (dressing, traveling, and food preparation). May be able to work part-time in sheltered environment but disabled for competitive work.
3	Severe overall disability: Conscious. Severe cerebral disability alone (CPC3) or severe disability from noncerebral organ system dysfunction alone or both. Dependent on others for daily support.
4	Coma, vegetative state: Not conscious. Unaware of surroundings, no cognition. No verbal or psychological interactions with environment.
5	Death: Certified brain dead by traditional criteria.

Glasgow Coma score

The Glasgow Coma score is a score given to head trauma victims based on the Glasgow Coma Scale (GCS). The total score is based on examination of eye opening, verbal output, and motor (movement) responses to different modalities and intensity of stimulation. You give each of these three attributes a score, and then Event Review Pro totals the three scores.

The choices in the GCS Total field are:

0 (3) (total up to 3)

1 (4-5) (total from 4 to 5)

2 (6-8) (total from 6 to 8)

3 (9-12) (total from 9 to 12)

4 (13-15) (total from 13 to 15)

Unknown



A score of “unknown” indicates that one or more response scores do not have an entry or have an entry of “unknown.”

Eye opening		Verbal output		Motor (movement) responses	
Score	Finding	Score	Finding	Score	Finding
				6	Obeys commands or Spontaneous or purposeful Patient follows commands
		5	Oriented or Coos or babbles Normal, fluent, appropriate speech	5	Purposeful or Withdrawal from touch Patient locates pain on his body

Eye opening		Verbal output		Motor (movement) responses	
4	Spontaneous Eyes open without stimulation	4	Confused or Irritable and continually crying Confused but fluent speech	4	Withdraws from pain Patient pulls away from pain
3	To voice Eye opening to loud noise	3	Inappropriate or Cries to pain Mumbling, occasional word recognizable	3	Flexion to pain or Abnormal flexion to pain Patient flexor postures
2	To pain Eye opening to pain only	2	Incomprehensible or Moans to pain Vocalizations but not verbalizations (no words)	2	Extension to pain Patient extensor postures
1	None No eye opening	1	None No vocalization or verbalization	1	None No motor response

Glossary

ACI-TIPI

For acute cardiac ischemia time insensitive predictive instrument. A Philips software tool in that can provide a second opinion and decrease the time between the onset of a patient's acute cardiac ischemia (ACI) symptoms and the treatment of interventional cardiology. The ACI-TIPI feature computes a patient's likelihood of having an ACI based on the patient's age, gender, chest pain status, and acquired 12-lead ECG. When the defibrillator runs the ACI-TIPI analysis, the interpretative block on the 12-lead report includes ACI-TIPI data.

Administration

Administration mode allows you to apply or change configuration settings in a HeartStart supported defibrillator, and also allows an administrator to view HeartStart Configure activities.

Administration mode

A setting on HeartStart defibrillators that allows the defibrillator to send and receive information.

Administration pane

The Administration pane allows you to view system log entries for the computer running HeartStart Configure.

AED

Stands for Automated External Defibrillator. A defibrillator that automatically performs rhythm analysis of the patient's surface electrocardiogram.

All Configurations

A feature in the Saved Configurations pane that allows you to view all configurations that are available for viewing and editing in HeartStart Configure.

Archive

To store or save information.

Attend to Patient Period

A period following a no shock Advised (NSA) decision during which the responder may perform CPR, if needed, or otherwise attend to the patient.

BLS

Basic Life Support

Bluetooth wireless transmission

Use of an optional FR3 Bluetooth wireless technology transceiver module to transfer a configuration from a HeartStart FR3 defibrillator to a Bluetooth-enabled computer.

Bluetooth®

A short-range wireless technology that uses radio links between devices such as defibrillators and computers, mobile computers, mobile phones, and other portable devices.

Bystander CPR

The attempt to perform basic CPR by someone who is not part of an organized emergency response system. Typically, this person witnessed the arrest. In certain situations, physicians, nurses, and paramedics may perform bystander CPR.

Card reader

Hardware that reads information from a data card.

CFG

Stands for configuration file. This binary file extension type is used for a HeartStart FR3 configuration file.

Configuration

A file that sets the behavior of a defibrillator.

Configuration Tasks

A task pane in HeartStart Configure that allows you to create defibrillator configurations per American Heart Association (AHA) guidelines, and set time, date, and language options for the HeartStart FR3 defibrillator.

Configurations (feature)

A feature set that allows you to configure, manage, and review information about a defibrillator configuration. HeartStart Configure identifies each configuration based on the type, name and the date and time of the configuration. This allows you to list configurations based on a variety of values.

Configurations pane

The Configurations pane allows you to view and manage supported defibrillator configurations.

CPR

Stands for Cardiopulmonary Resuscitation. A technique for providing artificial respiration and heart compressions to maintain life in a victim of Sudden Cardiac Arrest (SCA).

DAT

Stands for data file. This binary file format is used by Event Review 3.5 and Event Review Pro 3.5 for HeartStart FRx and HS1 defibrillator configurations.

Data Card

A computer storage device used for recording and storing information. Some HeartStart defibrillators use data cards to record configuration, ECG, and audio information.

Defibrillator event

Information received from a defibrillator. Examples are alarms, shocks, measurements, and error conditions.

Device

A generic term used for defibrillators and data cards.

ECG

Stands for electrocardiogram. The electrical rhythm of the heart as detected through defibrillator pads.

EMS

Stands for Emergency Medical Services.

Error

Describes a system error or an erroneous input by the user.

Export

Saving a configuration as a file for use on another computer running HeartStart Configure. This is typically done to share a configuration with other users or to back up a configuration.

Feature

A link or button within the HeartStart Configure task pane that allows you to complete a task.

Fibrillation

A disturbance of the normal heart rhythm that results in chaotic, disorganized activity that cannot effectively pump blood. Ventricular fibrillation (fibrillation in the lower chambers of the heart) is the most common cause of sudden cardiac arrest.

FR3

A HeartStart compact, battery-powered automated external defibrillator (AED) designed for use by trained responders to treat ventricular fibrillation (VF), the most common cause of sudden cardiac arrest (SCA).

FRx

A HeartStart compact, battery-powered automated external defibrillator (AED) designed for use by trained responders to treat ventricular fibrillation (VF), the most common cause of sudden cardiac arrest (SCA).

Getting Started

A task pane in HeartStart Configure that allows you to navigate to online help, selecting accessories, and check for software updates.

HS1

A HeartStart compact, battery-powered automated external defibrillator (AED) that is designed for simple and reliable operation by minimally trained users to treat ventricular fibrillation (VF), the most common cause of sudden cardiac arrest (SCA).

HTTP

Stands for Hypertext Transfer Protocol. The set of rules for exchanging files (text, graphic images, sound, video, and other multimedia files) on the World Wide Web.

Import

Opening a defibrillator configuration file that is not stored in HeartStart Configure, and was created by another user on another installation of HeartStart Configure.

Incident

The series of events involved in treating a patient with a defibrillator.

Information

Details action and events of an end-user using HeartStart Configure.

Language Card

A secure digital (SD) card that can be purchased separately for the HeartStart FR3 defibrillator model. It contains languages used by the local culture. You can either change the primary language, or set a second language for bilingual operation.

License Manager

The License Manager allows you to type in your registration key to activate HeartStart Configure software.

Navigation pane

Use the navigation pane to navigate to workspaces and to select a defibrillator and perform defibrillator configuration tasks.

NSA

Stands for no shock advised. A decision that is made by the defibrillator based on analysis of the patient's heart rhythm. The defibrillator will not deliver a shock in this mode.

PEA

Pulseless Electrical Activity.

PR

Perfusing Rhythm.

PST

For Periodic Self Test. A self-administered test mode that a HeartStart defibrillator runs to ensure that all defibrillator components are functional and work properly. Test results are Pass, Warning, or Call customer support.

Recorded "On" Time

The date and time the responder turned on the defibrillator.

Save

Adding a new configuration the information and its associated details to the HeartStart Configure dataset for reuse later on other defibrillators in order to standardize defibrillator behavior.

Saved configurations

The area of the Navigation pane where configurations can be viewed and filtered based on supported defibrillator type.

Shock series

Also called shock stack. A sequence of one or more shocks, each separated by no more than a preset interval. After completion of a shock series, the defibrillator automatically provides a CPR protocol.

SMART analysis

The proprietary algorithm used by the FR3 to analyze the patient's heart rhythm and determine whether a shock is advised.

Smart CPR

A configurable feature that, either automatically or by manual selection, enables a CPR interval before rhythm analysis and shock decision for patients with a shockable rhythm.

System log

The system log lists all monitored application activity. Use the system log to review application usage.

Task wizard

A software feature that guides the user through a task.

URL

Stands for Uniform Resource Locator. An address for a resource on the Internet. URLs are used by Web browsers to locate Internet resources. A URL specifies the protocol to be used in accessing the resource (such as http: for a World Wide Web page or ftp: for an FTP site), the name of the server on which the resource resides (such as //www.whitehouse.gov), and, optionally, the path to a resource (such as an HTML document or a file on that server).

Workspace

The pane to the right of the HeartStart Configure navigation pane. It displays the active window for entering and viewing information.

XML

Stands for eXtensible Markup Language. This file format has a set of rules for encoding documents in machine-readable form. It is used for HeartStart Configure 3.0 configuration files. A HeartStart Configure 3.0 file is exported in the XML file format.

	Average Response Times report	116
.		
.cod file		30
.hic file		30, 45
.mic file		30
.wfdb file format		86
1		
12-lead information		41, 90-91
attaching		34
magnifying		93
reviewing		92
A		
accessories, selecting		138
activating the product		14
email		15
Internet		15
within 60 days		17
All Records table		25
ATA flash card reader		146
attachments		26, 35
Bluetooth		35
ECGs		31, 34-41, 146
multiple ECGs		41
printing		104
sorting and filtering		103
viewing		103
Attachments table		26
audio, playing on ECGs		80
Average Response Times–Total System report		115
B		
Bluetooth		90, 137, 140-141, 144
accessories		140
attaching ECGs		35, 37-38
automatic synchronization		47
erasing data source		47
exchange folders		39
pairing and testing		142
prerequisites		141
C		
card readers		139
case		30
Case Capture 1.0		20
Case Details report		110
case ID		50
Case Summary		32
Case Wizard		31, 33
cases		30-32, 43
adding details		54
adding patient information		63
attaching files		102
creating		30-31, 33
defined		19, 30
deleting		52
duplicate ECGs		52
emailing		51
exporting		49
hiding and displaying		44
identifying		55
importing automatically		46
importing manually		44

migration from earlier versions	157	database	153
opening	33	remote server	153
printing	49	supported	19
reviewing details	100	verifying location	154
reviewing events	100	date and time	
saving	46	adjusting defibrillator	46
sorting and grouping	43	events, changing	62
Cases navigation pane	30	defibrillators	137
CDR file, sending to Customer Support	151	adjusting date and time	46
Cerebral Performance Categories (CPC)	160	device history data	20
channels		self-test	20
changing display time	74	supported	137
customizing	76	voice and system messages	151
displaying and hiding	73	detaching ECGs	42
magnifying waveforms	73	device history data	20
selecting	72	device tests	151
viewing ECG information	80	Do Not Show This Message Again	28
Channels pane	66, 75	document conventions	21
Click mode	27	downloading ECGs	34
color		Duplicate ECGs table	127
ECG selections	82		
vital trends display	90	E	
CompactFlash card reader	144	ECG Pre- and Post-shock report	110
compatibility	18	ECG selections	81
confirmation messages, restoring	28	color	82
cumulative device records, sending to Customer Support	151	creating	81
cursor delta time	80	managing	82
customer support	132	naming and sorting	83
emailing System Log to	120	Q-CPR	82
		templates	84
D		ECGs	65
data card for attaching ECGs	36, 42, 137, 139	12-Lead view	91-92
data migration from older versions	18	adding	34-41
data source, erasing	47	attaching	31, 34-41, 146
data transfer adapters and readers	138	changing time scale	75
data, migration from earlier versions	157	Channels pane	66
		copying to clipboard	85
		customizing display	88
		detaching	42
		duplicates	52, 127

Events pane	65		
multiple attachments	41		
navigating on waveform	79		
Overview pane	66		
playing audio	80		
removing	42		
reports	110		
shortcuts	96		
tab	65		
Transport pane	66		
view	65		
viewing information	80		
views (deleted)	20		
vital trends	89		
email			
cases	51		
reports	33		
erasing data source	47		
Event Log	26, 100		
events	26		
adding	61		
changing date and time	62		
documenting	60		
editing	61		
expanding and collapsing nodes	69		
filtering	101		
finding and viewing	69		
grouping	101, 122		
hiding and displaying	70, 101		
icons	68		
managing notes	70		
printing	101		
removing	61		
searching for	69		
sorting	62, 122		
tree structure	68		
viewing on waveform	72		
Events pane	65		
exclusions for Q-CPR			
creating	94		
removing	95		
exporting files	32		
multiple	20		
password protection	32		
		F	
		F1, using for Help	21
		fields, completing	25
		files, adding to case	102
		files, attaching ECGs from	40, 102
		filtering events	101
		FR2 defibrillators	37
		attaching ECGs	35
		data cards	145
		FR3 defibrillators	35, 37
		Bluetooth	140, 142
		data cards	140
		synchronizing clock via data erasure	48
		FRx defibrillators	36, 39, 42
		Admin mode	149
		case date and time	147
		G	
		Glasgow Coma Score (GCS)	162
		grouping events	101
		H	
		HeartStart Home defibrillators	36
		HeartStart OnSite defibrillators	36
		help	20, 22, 132
		hiding and displaying cases	44
		HS1 defibrillators	36, 39, 42
		Admin mode	150
		case date and time	147
		HTTP, using for Import Service	130

I

icons for events	68
Import Service	46, 127
archives	129
configuring	128
managing inboxes	127
using HTTP	130
importing files	20, 127
automatically	46
manually	44
infrared connections	35, 39, 42, 137, 147-148
installation	8
introduction	21
IrDA connections	35, 39, 42, 137, 147

K

key command shortcuts	96
-----------------------	----

M

main window	23
migration of older data	8, 157
mouse	27, 96
MRx defibrillators	19, 32, 35, 38, 55, 90
Bluetooth	142, 144
data cards	144
multiple files	
exporting	20
importing	20

N

navigation pane	23, 28
minimizing	29

new features	19
notes, managing	70

O

offset, waveform	77
online help	21
Overall Performance Categories (OPC)	160
Overview pane	66

P

password protection	32, 45, 49, 127
patient identity, removing	32
patient information, adding	63
PC card reader	146
PCMCIA card reader	146
pediatric option for FR2	37, 41
Percentile Response Times–Total System report	115
Percentile Response Times report	116
pinning and unpinning	73
point style, vital trends view	90
printing	
case information	49
reports	33, 117
product support	132

Q

Q-CPR	
adding notes to ventilation channel	95
creating exclusions	94
customizing episodes	93
removing exclusions	95
report data	112
Q-CPR Report Card report	111

R

redaction	50, 55, 107, 115
registration	17
removing ECGs	42
removing the application	17
reports	106
case details	110
ECGs	110
emailing	33, 119
exporting	118
generating	107
generating for open case	107
generating from database	107
printing	33, 117
Q-CPR	111
Response Times	114
toolbar icons	108
Vital Trends	114
resetting ECG display layout	74
resizing panes and workspaces	27
Response Times report	114

S

saving cases	28, 46
scale, waveform	78
selections templates, ECGs	84
self-test of defibrillators	20
shortcuts	96
SQL Server 2005, 2008	153
starting the product	23
synchronization, defibrillator and computer date and time	48
System Log	26, 120
emailing to Customer Support	120
filtering entries	122

grouping and sorting entries	122
printing entries	125
removing filters	125
system requirements	9

T

tables	
duplicate ECGs	52
printing	49
Tablet Mode	27
Tablet PC	27
telephone support	133
time	
changing on waveform	74
changing scale	75
Timeline tab	60
tooltips	20
Transport pane	66

U

uninstallation	17
version 3.5	18
Unwitnessed with Bystander CPR report	117
updates, checking for	12
upgrade	8
USB data drive, attaching ECGs from	37, 146
user accounts	20
Utstein guidelines	114, 116, 160
Utstein reports	116

V

ventilation channel, adding notes to	95
vertical offset, changing	77

vital trends	89
color	90
point style	90
Vital Trends report	114

W

waveform display	20
waveforms	72
audio	74
changing color	77
changing display time	74
changing scale	78
displaying in channel	77
ECG information	80
exporting data	86
incomplete releases	73
magnifying	73
measuring value	80
navigating ECG events	79
notes	86
offsets	77
viewing	71
zooming in and out	87
Waveforms pane	75
Windows 7, 64-bit	20
Windows Server	154
Witnessed with Bystander CPR report	117
Witnessed without Bystander CPR report	117
wizards	20, 40

X

XL defibrillators	
attaching ECGs	35
data cards	146
XL+ defibrillators	37, 49, 146
attaching ECGs	35