

Register your product and get support at
www.philips.com/welcome

SC2001

PHILIPS
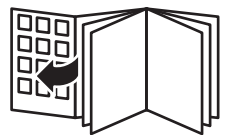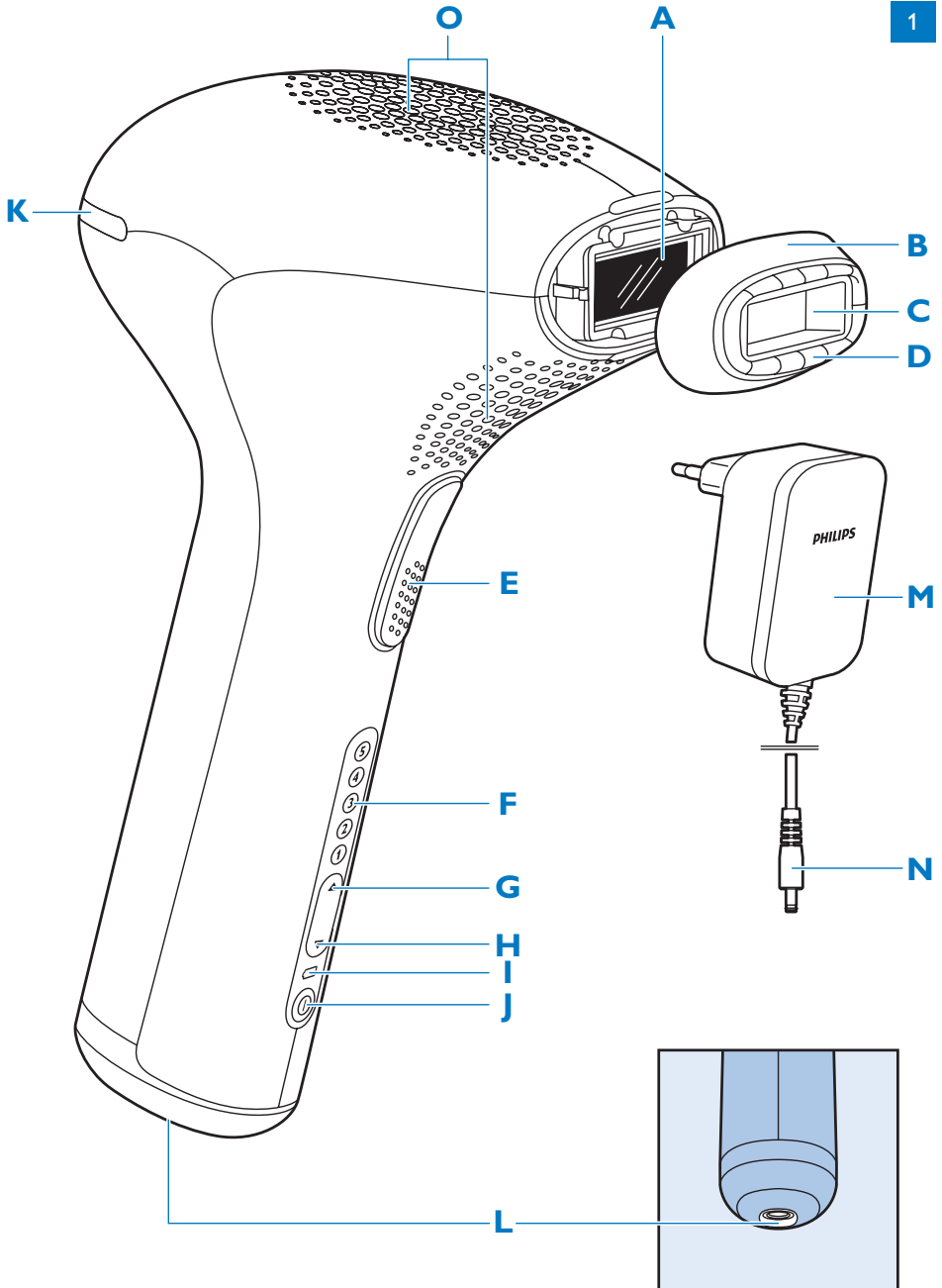
| | |
|-----------------|------------|
| ENGLISH | 6 |
| DANSK | 28 |
| ΕΛΛΗΝΙΚΑ | 50 |
| SUOMI | 78 |
| NORSK | 100 |
| SVENSKA | 121 |
| TÜRKÇE | 143 |

Table of contents

| | |
|---|----|
| Introduction | 6 |
| Benefits | 7 |
| Intense Pulsed Light (IPL) technology for home use | 7 |
| Effective hair re-growth prevention for everyday smoothness | 7 |
| Gentle treatment, even on sensitive body areas | 7 |
| Cordless operation for maximum freedom and flexibility | 7 |
| Lamp does not require replacement | 7 |
| General description | 7 |
| Important | 8 |
| Danger | 8 |
| Warning | 8 |
| Caution | 9 |
| For whom is Philips Lumea not suitable? | 9 |
| Contraindications | 10 |
| How Philips' Intense Pulsed Light technology works | 11 |
| Hair growth | 11 |
| Working principle | 11 |
| Charging | 12 |
| Preparing for use | 12 |
| Preparing treatment areas | 12 |
| Setting the light intensity | 13 |
| Recommended light intensities (1-5) | 14 |
| Tanning advice | 14 |
| Using the appliance | 15 |
| First use | 15 |
| Subsequent use | 16 |
| Guidelines for treatment | 17 |
| Average treatment time | 17 |
| Using the appliance on the legs | 18 |
| How to achieve optimal results | 18 |
| Possible side effects and skin reactions | 19 |
| After use | 19 |
| Cleaning and maintenance | 19 |
| Storage | 20 |
| Environment | 21 |
| Removing the rechargeable batteries | 21 |
| Guarantee and service | 22 |
| Replacement parts | 22 |
| Technical specifications | 22 |
| Troubleshooting | 23 |

Introduction

Congratulations on your purchase and welcome to Philips! To fully benefit from the support that Philips offers, register your product at www.philips.com/lumea.

Hair removal with light is one of the most effective methods to get rid of unwanted body hair with longer-lasting results. It is completely different from current home hair removal methods and Philips has been investigating this method together with leading dermatologists for almost 10 years. You can now use this innovative technology comfortably in the privacy of your own home, whenever and wherever you want.

This appliance has been designed and developed for use by women only and is intended for removing body hair on areas mentioned below:

- underarms
- bikini area
- legs

Do not use the Philips Lumea IPL hair removal system on the face. Philips Lumea is designed for use by one person only. Philips Lumea is not designed for use by men.

For more information about Philips Lumea, see www.philips.com/lumea. Check our homepage regularly for updates.

Benefits

This appliance provides the following benefits:

Intense Pulsed Light (IPL) technology for home use

Philips Lumea uses a light-based technology called 'Intense Pulsed Light' (IPL). IPL is also used in the professional beauty market for hair removal and has been in use for the last fifteen years. Philips Lumea now brings this innovative technology to the comfort of your home.

Effective hair re-growth prevention for everyday smoothness

This method provides you with long-lasting hair regrowth prevention and silky-smooth skin. Our studies have shown a significant hair reduction already after two treatments. Optimal results were generally achieved after four to five treatments. A reduction in hair density is visible from treatment to treatment. How fast the effect becomes visible and how long the results last varies from person to person. To treat all hairs successfully and prevent the hair root from becoming active again, the treatments have to be repeated once every two weeks for the first two months. Your skin should be smooth and hair-free by then. To ensure that you keep these results, we advise you to repeat treatments on the underarms and bikini line every two weeks. On the legs, regrowth is usually slower. Therefore, a treatment once every four weeks is probably sufficient to keep your legs smooth and hair-free. If you see too many hairs grow back within these four weeks, simply shorten the interval between treatments again, but do not treat more often than once every two weeks.

Gentle treatment, even on sensitive body areas

Philips Lumea has been developed in close cooperation with leading dermatologists for an effective and gentle treatment, even on the most sensitive body areas.

Cordless operation for maximum freedom and flexibility

Because the appliance runs on rechargeable batteries, it is portable and can be used anywhere.

Lamp does not require replacement

Philips Lumea has a lamp with a capacity of 40,000 flashes, which is enough to last the product's lifetime if used according to the instructions.

General description (Fig. 1)

- A** Light exit window with filter glass
- B** Removable window
- C** Metallic surface inside the removable window
- D** Safety system (safety ring with contact switches)

- E** Flash button
- F** Intensity lights (1- 5)
- G** ▲ Intensity increase button
- H** ▼ Intensity decrease button
- I** Charging light and battery low indication
- J** On/off button
- K** 'Ready to flash' light
- L** Appliance socket
- M** Adapter
- N** Small plug
- O** Air vents

Also included: pouch and cleaning cloth

Important

Read this user manual carefully before you use the appliance and save it for future reference.

Danger

- Water and electricity are a dangerous combination. Do not use this appliance in wet surroundings (e.g. near a bath, a shower or a swimming pool).
- Keep the appliance and the adapter dry.

Warning

- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- The appliance is not intended for children younger than 14 years. Teenagers aged between 15 and 18 years can use the appliance with the consent and/or assistance of their parents or the persons who have parental authority over them. Adults from 19 years and older can use the appliance freely.
- Only charge the appliance with the adapter supplied.
- Do not use the appliance or the adapter if it is damaged.
- Do not use the appliance if the filter glass is broken.
- If the adapter is damaged, always have it replaced with one of the original type in order to avoid a hazard.
- The adapter contains a transformer: Do not cut off the adapter to replace it with another plug, as this causes a hazardous situation.
- To prevent damage:
- Never cover the ventilation slots during use.

- Make sure nothing can obstruct the air flow through the ventilation slots.
- Never subject the appliance to heavy shocks and do not shake or drop it.
- If you take the appliance from a very cold environment to a very warm environment, wait approximately 3 hours before you use it.
- Never use or store the appliance in a dusty environment.
- Never leave the appliance unattended when it is switched on. Always switch off the appliance after use.
- Do not expose the appliance to direct sunlight or UV light for a longer time.
- Do not flash against white, smooth or reflective surfaces like paper; mirrors, floors, walls etc. This can cause severe damage to the removable window. Only flash when the appliance is in contact with skin.

Caution

- Please use this appliance only at settings suitable for your skin type. Use at higher settings than those recommended can increase the risk of skin reactions and side effects.
- This appliance is only intended for removing women's body hair from areas below the neck. Do not use it for any other purpose. This appliance is not designed for use by men.
- **Keep the following things in mind to avoid compromising the lifetime of Philips Lumea: Use by more than one person and more frequent use than indicated in this user manual can severely hamper the efficacy of the appliance and can cause certain components to stop functioning.**
- Do not expose the appliance to temperatures lower than 15°C or higher than 35°C during use.
- The appliance is equipped with a built-in overheat protection. If the appliance overheats, it does not flash. Do not switch off the appliance, but let it cool down for approx. 15 minutes before you continue with your treatment.
- Always return the appliance to a service centre authorised by Philips for examination or repair. Repair by unqualified people could cause an extremely hazardous situation for the user.
- Never use the appliance in any of the cases mentioned under Contraindications.

For whom is Philips Lumea not suitable?

Philips Lumea is not designed for everyone. If any of the following is true for you, then this appliance is **not** suitable for you to use!

- If your natural hair colour in the treatment areas is light blond, white, grey or red.
- If you are of skin type VI (you rarely or never get sunburnt and you tan very darkly). In this case you run a high risk of developing skin reactions such as strong redness, burns or discolourations when you use Philips Lumea.
- There are medications and health conditions that may lead to adverse reactions when you use Philips Lumea. Please read section 'Contraindications' to determine whether any of these apply to you.

- All types of tanning can influence the safety of the treatment. This applies to exposure to natural sunlight as well as artificial tanning methods such as tanning lotions, tanning appliances etc. If you are tanned or are planning to tan, please see chapter 'Preparing for use', section 'Tanning advice'.

Contraindications

Your skin may become a little red and irritated after you have used the appliance. This phenomenon is absolutely normal and quickly disappears. If the irritation does not disappear within three days, we advise you to consult a doctor.

Never use the appliance:

- If you have an unsuitable skin and/or hair type (consult the skin/hair colour table to determine whether this appliance is suitable for you).
- If you have a skin disease such as active skin cancer; if you have a history of skin cancer or any other localised cancer in the areas to be treated.
- If you have pre-cancerous lesions or multiple atypical moles in the areas to be treated.
- If you have a history of collagen disorder; including a history of keloid scar formation or a history of poor wound healing.
- If you have a history of vascular disorder; such as the presence of varicose veins or vascular ectasia in the areas to be treated.
- If your skin is sensitive to light and easily develops a rash or an allergic reaction.
- If you have infections, eczema, burns, inflammation of hair follicles, open lacerations, abrasions, herpes simplex, wounds or lesions and haematomas in the areas to be treated.
- If you have had surgery in the areas to be treated.
- If you have epilepsy with flashlight sensitivity.
- If you have diabetes, lupus erythematoses, porphyria or congestive heart disease.
- If you have any bleeding disorder.
- If you have a history of immunosuppressive disease (including HIV infection or AIDS).
- On areas of your skin where you are currently being treated with or have recently been treated with Alpha-Hydroxy Acids (AHAs), Beta-Hydroxy Acids (BHAs), topical isotretinoin and azelaic acid.
- If you have taken oral isotretinoin Accutane® or Roaccutane® in the last six months. This treatment can make skin more susceptible to tears, wounds and irritations.
- When you are on painkillers, which reduce skin's sensitivity to heat.
- If you are taking photosensitising agents or medications, check the package insert of your medicine and never use the appliance if it is stated that it can cause photo-allergic reactions, photo-toxic reactions or if you have to avoid sun when taking this medicine.
- If you take anticoagulation medications, including heavy use of aspirin, in a manner which does not allow for a minimum 1-week washout period prior to each treatment.
- If you take immunosuppressive medications.
- Over or near anything artificial like silicon implants, Implanon contraceptive implants, pacemakers, subcutaneous injection ports (insulin dispenser) or piercings.

- On moles, freckles, large veins, darker pigmented areas, scars, skin anomalies without consulting your doctor. This can result in a burn and change of colour, which makes it potentially harder to identify skin-related diseases (e.g. skin cancer).
- On nipples, areolas and mucous areas (inner vagina and anus).
- On warts, tattoos or permanent make-up. This can result in a burn and change of colour.
- After intensive sunbathing and/or on sunburnt skin or if you use tanning accelerators. See chapter 'Preparing for use', section 'Tanning advice' for more information.
- If you use long-lasting deodorants. This can result in skin reactions. See chapter 'Using the appliance', section 'Possible side effects and skin reactions'.
- If you are pregnant or breastfeeding.

How Philips' Intense Pulsed Light technology works

Hair growth

Depending on age, metabolism and other factors, everyone's hair grows differently, but hair normally grows in three phases:

1 Anagen phase (growing phase)

The hair grows actively from the root. In this phase the concentration of melanin is highest. Melanin is responsible for the pigmentation of the hair. The more melanin there is in the root, the darker the hair is. Melanin concentration is also the determining factor for the efficacy of this method. Only hairs in the anagen phase are susceptible to treatment with light.

2 Catagen phase (degradation phase)

The hair growth stops and the root shrinks before the hair sheds.

3 Telogen phase (resting phase)

The old hair separates from the root and sheds. The root remains at rest until the biological clock tells it to become active again and grow a new hair. Hairs on an area of skin are usually in different stages of the cycle of hair growth.

Working principle


The appliance works by heating up the hair and root beneath the skin.

The melanin in the hair and hair root absorbs the released light pulses. The darker the colour of the hair, the more heat can be absorbed. This process stimulates the hair to go into the resting phase.

The hair now sheds naturally and the hair growth is inhibited. Hairs can take one or two weeks to shed. Optimal hair removal cannot be obtained with a single session, since only hairs in the anagen phase are susceptible to treatment with light. During a treatment, hairs on the treated area are usually in different stages of the cycle of hair growth.

To treat all hairs successfully and prevent the hair root from becoming active again, the treatments have to be repeated once every two weeks for the first two months. Your skin should be smooth and hair free by then.

To ensure that you keep these results, it is advised that you keep repeating treatments every two weeks on the underarms and bikini line. On the legs, however, regrowth of hair is usually slower. Therefore, a treatment once every four weeks should also be sufficient to keep your legs smooth and hair free. If too many hairs grow back in these four weeks, then simply shorten the interval between treatments again, but do not treat more often than once every two weeks.

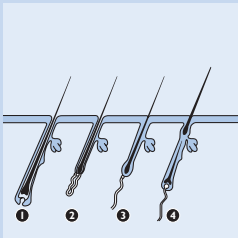An integrated optical filter makes sure no UV light reaches the skin and therefore leaves it unharmed.

Electromagnetic fields (EMF)

This Philips appliance complies with all standards regarding electromagnetic fields (EMF). If handled properly and according to the instructions in this user manual, the appliance is safe to use based on scientific evidence available today.

Charging

Fully charge the batteries before you use the appliance for the first time and when the batteries are empty. Fully charging the batteries takes up to 1.5 hours.

Charge the appliance when the charging light lights up orange during use to indicate that the battery is low. When this happens, the remaining capacity is at least 30 flashes.

Fully charged batteries provide at least 150 flashes at light intensity 5, which amounts to approx. 15 minutes of cordless operation time.

Fully charge the appliance every 3 to 4 months, even if you do not use the appliance for a longer time.

Charge the appliance in the following way:

- 1** Switch off the appliance.
- 2** Insert the small plug into the appliance and put the adapter in the wall socket.
 - ▶ The charging light flashes green to indicate that the appliance is charging.
 - ▶ When the batteries are fully charged, the charging light lights up green continuously.

Notes:

- The adapter feels warm during charging. This is normal.
- You cannot use the appliance during charging.
- This appliance is equipped with battery overheat protection and does not charge if the room temperature exceeds 40°C.

Never cover the appliance and adapter during charging.

- 3** After charging, remove the adapter from the wall socket and pull the small plug out of the appliance.

Preparing for use

For optimal results and highest efficacy, prepare the areas on which you intend to use the appliance in the manner described below.

Preparing treatment areas

- 1** Shave the areas you intend to treat as long as you still experience hair regrowth.

Note: If your last depilation method involved removing hairs by the roots (e.g. epilating, waxing etc.), wait until you notice significant hair regrowth before you use Philips Lumea. Make sure you shave the area before you use Philips Lumea. Shaving is necessary to ensure that the light energy is channelled to the hair root in the most efficient way.

Please note that treatment on unshaved or improperly shaved areas has a number of undesired effects:

- Contamination (burns/black spots) of the removable window with stray hairs and dirt. A contaminated removable window can reduce the lifetime and efficacy of the appliance. In addition, burns/black spots on the removable window can cause the treatment to become painful or cause skin reactions like redness and discolouration.
- Hairs on your skin can burn and this can result in an unpleasant smell.

2 Clean your skin and make sure it is hair-free, entirely dry and free from oily substances (e.g. present in cosmetic products like deodorants, lotions, perfumes, tanning creams and sunscreens).

3 Check the appliance, especially the removable window and filter glass for trapped hair, dust and fibres. If you see any contamination, clean the appliance according to the instructions in chapter 'Cleaning and maintenance'.

Notes:

- If shaving causes wounds or lesions, your skin may be more sensitive and the light treatment might feel uncomfortable sometimes.
- Once the hairs stop growing back - which usually happens after 4-5 treatments - you no longer have to shave before you use the appliance.
- If shaving causes skin irritation, you can also use a trimmer to trim the hair as short as possible. Light treatment after trimming might make it a little uncomfortable and you may get the smell of burning hair. Make sure you always use a setting that is comfortable for you. Using the appliance on trimmed hair might increase the chances of hairs getting trapped in the recessed window. Make sure you clean the appliance carefully after each treatment.

Setting the light intensity

Philips Lumea is an appliance developed for consumer use at home and therefore NOT COMPARABLE with professional light-based appliances with respect to pain sensation. Pulses of light delivered by Philips Lumea are gentler than those from professional appliances and therefore you can expect a convenient and effective treatment.

You have to adjust the light intensity to your skin and body hair colour and to a level you find comfortable.

The table with recommended light intensities below helps you to establish the right light intensity.

1 Consult the table below to determine which light intensities are most suitable for your skin and body hair colour and to check if this method is suitable for you (if it is not, this is indicated with x in the table).

Persons with darker skin run a higher risk of skin irritation because their skin absorbs more light. Therefore lower intensities are recommended for persons with darker skin.

2 For each body area you want to treat, select a suitable intensity in the recommended range that is as high as possible but does not cause any discomfort on the skin. Each body area and possibly even parts of a body area may require a separate appropriate light

intensity. Therefore you should lower the light intensity if you experience any discomfort anywhere.

The light intensity you use can feel warm or hot on the skin, but should never feel uncomfortable. To prevent skin reactions, see chapter 'Using the appliance', section 'Possible side effects and skin reactions'.

Recommended light intensities (1-5)

In the table below, 'x' means that the appliance is not suitable for you.

| Skin colour & characteristics: | White (always sunburn, no tanning) | Beige (easy sunburn, minimal tanning) | Light brown (sometimes sunburn, good tanning) | Mid brown (rarely sunburn, fast and good tanning) | Dark brown (rarely sunburn, very good tanning) | Black (rarely to never sunburn, very dark tanning) |
|--|------------------------------------|---------------------------------------|---|---|--|--|
| hair colour: dark blond/light brown/ brown/dark brown/ black | 4/5 | 4/5 | 3/4/5 | 3/4 | 1/2/3 | x |
| hair colour: white/ grey/red/light blond | x | x | x | x | x | x |

Important:

- To check whether Philips Lumea is suitable for you, see section 'For whom is Philips Lumea not suitable?'.
- At higher light intensities:
 - the method is more effective.
 - you have a higher risk of skin reactions. See chapter 'Using the appliance', section 'Side effects and skin reactions'.
 - the pain sensation can be stronger.
- If your skin has been exposed to natural sunlight or an artificial tanning method, your skin may be darker than during the previous treatment. Perform a skin test (see chapter 'Using the appliance', section 'First use') on the area to be treated to determine if you have to use a lower light intensity than usual. For more information, see section 'Tanning advice' below.

Tanning advice

Sun exposure before treatment

- Wait at least 48 hours after sunbathing before you use the appliance. Sunbathing can cause skin to darken (tan) and/or can cause sunburn (reddening of the skin).
- After these 48 hours, check the skin for any remaining sunburn. Do not use Philips Lumea as long as sunburn persists.
- When the sunburn has disappeared, perform a skin test (see chapter 'Using the appliance', section 'First use') on the area to be treated to determine the appropriate light intensity.
- If your skin is tanned, you run a very high risk of developing skin reactions after treatment with Philips Lumea. Perform a skin test (see chapter 'Using the appliance', section 'First use') on the area to be treated to determine the appropriate light intensity.

Sun exposure after treatment

- Wait at least 24 hours after treatment before you expose treated areas to the sun. Even after 24 hours, make sure that the treated skin does not show any redness from the treatment anymore before you expose it to the sun.
- Cover treated areas when you go out in the sun or use a non-perfumed sun block in the two weeks after the treatment.
- Exposing treated areas to the sun immediately after treatment without protection can increase the risk of adverse side effects and skin reactions. See chapter 'Using the appliance', section 'Possible side effects and skin reactions'.

Tanning with artificial light

- Follow the instructions in 'Sun exposure before treatment' and 'Sun exposure after treatment' above.

Tanning with creams

- If you have used artificial tanning lotion, wait until the artificial tan has fully disappeared before you use Philips Lumea.

Note: For more details on unwanted skin reactions, see chapter 'Using the appliance', section 'Possible side effects and skin reactions'.


Using the appliance

Shave the area you intend to treat before you use the appliance (see chapter 'Preparing for use').


If the treatment becomes intolerably painful at any point in time, reduce the light intensity setting until the treatment is comfortable again. If you notice skin reactions stronger than slight redness, stop the treatment immediately and read section 'Possible side effects and skin reactions' in this chapter.

First use

When you use the appliance for the first time, we advise you to perform a skin test on the area to be treated to check your skin's reaction to the treatment and to get accustomed to light-based treatment.

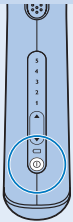
- 1 Choose an area close to the intended treatment area.
- 2 Press the on/off button to switch on the appliance.

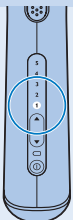- 3 Intensity light 1 goes on automatically to indicate that the appliance is switched on at default intensity 1.
- 3 Apply one flash at the lowest recommended intensity for your skin type.
- 4 As long as it feels comfortable (not painful), increase the setting by one level within the recommended range for your skin type and apply one flash for each setting. Do not apply more than one flash on the same spot.


- 5** After the skin test, wait for 24 hours and check your skin for any reaction. If the skin shows reactions, choose the highest setting that did not result in any skin reaction for subsequent use.

Subsequent use

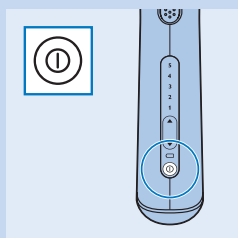- 1** Before every treatment, clean the removable window, filter glass and metallic surface of Philips Lumea.

- 2** Press the on/off button to switch on the appliance.

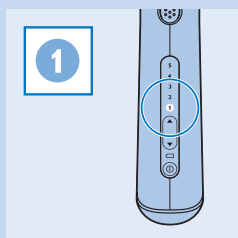
- 1** Intensity light 1 goes on automatically to indicate that the appliance is switched on at default intensity 1.

Note: It is normal that the appliance becomes warm during use.


- 3** Start each treatment by selecting the light intensity very carefully, based on your skin's sensitivity and comfort at that moment.

Your skin may react differently on different days/occasions for a number of reasons. See section 'Possible side effects and skin reactions' in this chapter for more information.

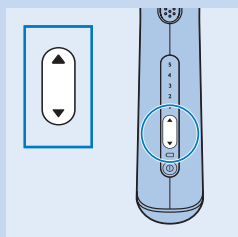
- 4** Press the ▲ button to increase the intensity. To lower the intensity, press the ▼ button.

- 1** Every time you press the button, the corresponding intensity light starts to flash. This takes a few seconds.

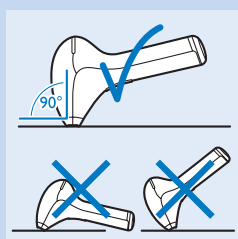
- 5** Place the appliance perpendicularly on the skin so that the removable window and the safety ring are in contact with the skin.

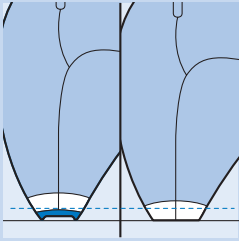- 6** Press the safety ring fully onto the skin with slight pressure.


Note: The safety ring has contact switches that together are the safety system of the appliance. This safety ring prevents unintentional flashing without skin contact.

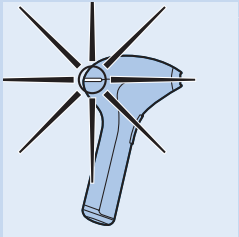- When all contact switches of the safety ring have been pushed into the appliance, the 'ready to flash' light on the back of the appliance goes on. This may take up to 5 seconds. When the 'ready to flash' light is on, the appliance is ready to release a flash.

Tip: For easier use, stretch the skin you want to treat. On areas with softer skin, you may need to exert more pressure on the appliance to push in all contact switches.

Note: Bony areas like shinbone and ankles are more sensitive than other areas of the body. When you press the window hard onto these areas, the skin can be more susceptible to skin reaction like irritation and redness. See section 'Possible side effects and skin reactions' in this chapter. Press the window onto the skin just hard enough to make the 'ready to flash' light go on, but no harder.

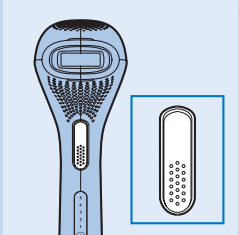- 7** Press the flash button to release a flash. After the appliance has flashed once, release the flash button.

Note: The visible light produced by the appliance is the reflection of the flash on the skin and is harmless to your eyes. It is not necessary to wear goggles during use.

Note: Avoid looking straight into the light that comes out of the window when the appliance is switched on.


- 8** For the next flash, repeat steps 5, 6 and 7. After each flash it takes 3 to 5 seconds until the appliance is ready to flash again.

Never flash the same area twice without lifting the appliance off the skin. This prevents excessive heating and possible burning of the skin.

- 9** To avoid untreated areas, always make sure there is some overlap with the previously treated area when you place the appliance on the skin. The effective light only comes out of the light exit window. Make sure the flashes are produced close to each other.

- 10** Switch off the appliance after you have finished the treatment.

- 11** Check the removable window and filter glass for hairs and dirt and clean the appliance after use. See chapter 'Cleaning and maintenance'.

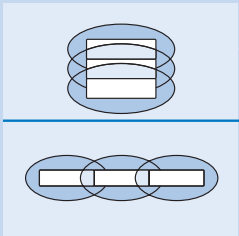Guidelines for treatment

Average treatment time

| Area | Approx. treatment time |
|--------------------|------------------------|
| Underarms | 2 minutes per underarm |
| Bikini line | 2 minutes per side |
| Entire bikini area | 8 minutes |
| One lower leg | 15 minutes |
| One full leg | 30 minutes |

Note: The average treatment times above have been observed during research and tests. Individual treatment times may vary from person to person.

Note: Fully charged batteries offer at least 150 flashes at intensity 5, which amounts to an operating time of approx. 15 minutes.

Using the appliance on the legs

If you use the appliance to treat your legs, you may have to recharge the appliance during the treatment. This is normal. How long you can use the appliance depends on the size of the treatment area and the intensity setting used.

Tip: Many users combine leg treatments with other activities around the house such as watching TV or listening to music. You can also try this to make the experience more enjoyable and relaxing.

How to achieve optimal results

- Optimal results cannot be obtained with one single session. To remove all hairs successfully and prevent the hair root from becoming active again, the treatments have to be repeated once every two weeks for the first two months. Your skin should be smooth and hair-free by then.
- To ensure that you keep these results, we advise you to repeat the treatments on the underarms and bikini line every two weeks. On the legs regrowth is usually slower. Therefore, treatment of the legs once every four weeks is probably sufficient to keep your legs smooth and hair-free. If too many hairs grow back within these four weeks, simply shorten the interval between treatments again, but do not treat more often than once every two weeks. It is recommended to adhere to this frequency for optimum results. When you use the appliance more frequently than recommended, you also increase the risk of skin reactions. See section 'Possible side effects and skin reactions' in this chapter.
- Avoid treating exactly the same area multiple times during one session as this does not improve the effectiveness of the appliance. However, this does increase the risk of skin reactions.
- Never flash the same area twice without lifting the appliance off the skin. This prevents excessive heating and possible burning of the skin.
- It is normal that hairs are still visible after the first few treatments. See chapter 'Important', section 'How photo-epilation works' for more information.

- Some hairs grow back even though you use the appliance frequently. This happens because hairs are in different stages of the cycle of hair growth during treatment. The hairs that grow back tend to be softer and thinner.

Possible side effects and skin reactions

- Your skin shows slight redness and/or prickles, tingles or feels warm. This phenomenon is absolutely harmless and disappears quickly.
- A skin reaction similar to sunburn may occur on the treated area for a short while. If this does not disappear within 3 days, we advise you to consult a doctor.
- Dry skin and itching may occur because of the combination of shaving and light treatment. This is harmless and disappears within a few days. You can cool the area with an ice pack or a wet facecloth. If dryness persists, you can apply a non-scented moisturiser on the treated area 24 hours after the treatment.
- Burns, excessive redness and swelling: these reactions occur very rarely. They are the result of using a light intensity that is too high for your skin colour. If these reactions do not disappear within 3 days, we advise you to consult a doctor. Wait with the next treatment until the skin has healed completely and make sure you use a lower light intensity.
- Skin discolouration: this occurs very rarely. Skin discolouration manifests itself as either a darker or lighter patch than the surrounding area. This is the result of using a light intensity that is too high for your skin colour. If the discolouration does not disappear within 2 weeks, we advise you to consult a doctor. Do not treat discoloured areas until the discolouration has disappeared and your skin has regained its normal skin colour.
- Skin infection and inflammation: This occurs very rarely and may be caused by using the appliance on wounds or cuts that resulted from shaving, on already existing wounds or on ingrown hairs.
- Excessive pain: this can occur during or after treatment if you have used the appliance on unshaved skin, if you use the appliance at a light intensity that is too high for your skin colour, if you flash the same area more than once and if you use the appliance on open wounds, inflammations, infections, tattoos, burns, etc. See chapter 'Important', section 'Contraindications' for more information.

After use

- Do not apply any cosmetic substances to the treated areas immediately after treatment.
- Do not use deodorant after the treatment until any redness has disappeared completely.
- See chapter 'Preparing for use', section 'Tanning advice' for information about sun exposure and artificial tanning.

Cleaning and maintenance

To ensure optimal results and a long life for the appliance, clean the appliance before and after each use and, if necessary, also during the treatment. The appliance loses its effectiveness if you do not clean it properly.

Never clean the appliance or any of its parts under the tap or in the dishwasher.

Never use scouring pads, abrasive cleaning agents or aggressive liquids such as petrol or acetone to clean the appliance.

Never scratch the filter glass or the metallic surface inside the removable window.

Note: Make sure that the contact switches of the safety ring do not get jammed with dirt.

Note: Stop using the appliance when it is no longer possible to clean the filter glass or the removable window.

- 1 Switch off the appliance, unplug it and let it cool down.

Note: The filter glass becomes hot during use. Make sure it has cooled down before you clean it.

- 2 To remove the removable window, put your fingers in the recesses at the top and bottom and pull gently. The window should come off easily.



- 3 Moisten the soft cloth supplied with the appliance with a few drops of water and use it to clean the following parts:

- the surface of the filter glass
- the outside surface of the removable window
- the metallic frame inside the removable window

Note: If water does not clean effectively, use a few drops of 70% isopropyl alcohol to clean the parts mentioned above.

Note: A contaminated recessed window can look like this picture. Make sure you clean the appliance before and after each use and, if necessary, also during the treatment according to the instructions in this section.

- 4 If necessary, clean the outside of the appliance with the dry soft cloth supplied with the appliance.

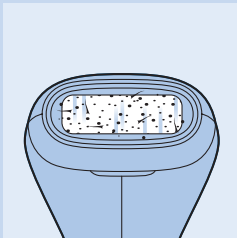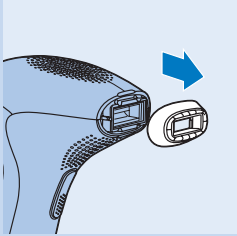Storage

- 1 Switch off the appliance, unplug it and let it cool down.
- 2 Clean the appliance before you store it.
- 3 Put the appliance in the pouch supplied.
- 4 Store the appliance in a dry place at a temperature between 0°C and 60°C.

Environment


- Do not throw away the appliance with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment.


- The built-in rechargeable batteries contain substances that may pollute the environment. Always remove the batteries before you discard and hand in the appliance at an official collection point. Dispose of the batteries at an official collection point for batteries. If you have trouble removing the batteries, you can also take the appliance to a Philips service centre. The staff of this centre will remove the batteries for you and will dispose of them in an environmentally safe way.

Removing the rechargeable batteries


Only remove the rechargeable batteries if they are completely empty.

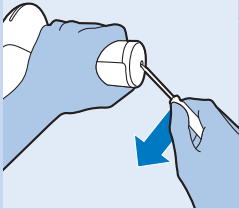
- 1 Insert an awl or another pointed tool into the bottom cap of the appliance and pry off the bottom cap.

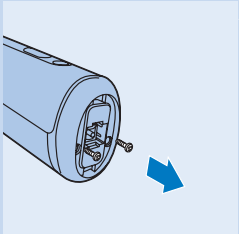
- 2 Undo the two screws.

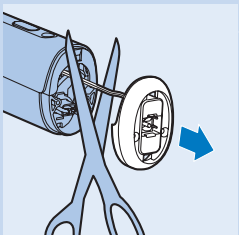
- 3 Remove the bottom part of the handle and cut the two wires that connect the bottom part to the appliance.

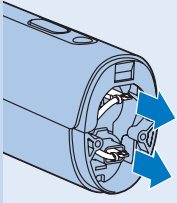4 Pull the battery connectors off the battery terminals.

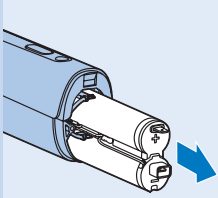5 Shake the rechargeable batteries out of the appliance.

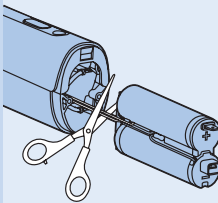6 Cut the two wires one at a time to ensure no mishaps due to possible residual charge.

Guarantee and service

If you need service or information or if you have a problem, please visit the Philips website at www.philips.com or contact the Philips Consumer Care Centre in your country (please look for the phone number in the worldwide guarantee leaflet). If there is no Consumer Care Centre in your country, go to your local Philips dealer.

Replacement parts

The following replacement parts are available:

- Removable window
- Adapter
- Pouch

Technical specifications

Model SC2001

| | |
|-----------------------------------|--|
| Rated voltage | 100V-240V |
| Rated frequency | 50Hz-60Hz |
| Rated input | 7.5W |
| Protection against electric shock | Class II  |
| Protection rating | IP 30 (EN 60529) |
| Operating conditions | Temperature: +15° to +35°C Relative humidity: 25% to 75% |
| Storage conditions | Temperature: 0° to +60°C Relative humidity: 5% to 95% |
| Intense pulsed light lamp | Spectrum: >570nm |
| Lithium-ion battery | 2 x 3.7 Volt 1100 mAh |

Troubleshooting

This chapter summarises the most common problems you could encounter with the appliance. If you are unable to solve the problem with the information below, contact the Consumer Care Centre in your country (please look for the phone number in the worldwide guarantee leaflet).

| Problem | Possible cause | Solution |
|--|--|--|
| The appliance does not work. | The rechargeable batteries are empty. | Charge the appliance (see chapter 'Charging'). |
| | The appliance is broken. | Contact the Consumer Care Centre in your country, your Philips dealer or a Philips service centre. |
| The appliance suddenly switches off. | The rechargeable batteries are empty. | Charge the appliance (see chapter 'Charging'). |
| The adapter becomes warm during charging. | This is normal. | No action required. |
| The charging light does not go on when I put the small plug in the appliance socket. | You have not put the adapter in the wall socket. | Insert the small plug into the appliance socket and put the adapter in the wall socket. |

| Problem | Possible cause | Solution |
|--|---|---|
| | The wall socket is not live. | Connect another appliance to the wall socket to check if the wall socket is live. If the socket is live but the appliance still does not charge, contact the Consumer Care Centre in your country, your Philips dealer or a Philips service centre. |
| | The adapter is not inserted properly into the wall socket and/or the small plug is not inserted properly into the appliance socket. | Make sure the adapter is inserted properly into the wall socket and/or the small plug is inserted properly into the appliance socket. |
| | The appliance is broken. | Contact the Consumer Care Centre in your country, your Philips dealer or a Philips service centre. |
| I have switched on the appliance but I cannot increase or decrease the light intensity. | The appliance has to be reset. | To reset the appliance, switch off the appliance and then switch it on again. If you still cannot adjust the light intensity, contact the Consumer Care Centre in your country, your Philips dealer or a Philips service centre. |
| The 'ready to flash' light does not go on when I place the appliance on the skin. | You have not put the removable window with safety ring properly on the skin. | Place the appliance perpendicularly onto the skin in such a way that all contact switches touch the skin. Then check if the 'ready to flash' light is on and press the flash button. |
| | | If this does not work, place the appliance on a part of your body where it is easy to establish full skin contact, for example on your lower arm. Then check if the 'ready to flash' light works. If the 'ready to flash' light still does not go on when you place the appliance on the skin, contact the Consumer Care Centre in your country, your Philips dealer or a Philips service centre. |
| | The safety ring is dirty. | Clean the safety ring carefully. If you are unable to clean the safety ring properly, contact the Consumer Care Centre in your country, your Philips dealer or a Philips service centre to replace the flash screen. |
| The 'ready to flash' light does not go on when I place the removable window on my skin, but the cooling fan works. | The overheat protection has been activated. | When the overheat protection has been activated, the fan still works. Do not switch off the appliance and let it cool down for approx. 15 minutes before you continue to use it. If the 'ready to flash' light still does not go on when you place the appliance on the skin, contact the Consumer Care Centre in your country, your Philips dealer or a Philips service centre. |

| Problem | Possible cause | Solution |
|---|--|--|
| The 'ready to flash' light is green but the appliance does not produce a flash when I press the flash button. | The appliance needs to be reset. | To reset the appliance, switch off the appliance and switch it on again. If the 'ready to flash' light is green but the appliance still does not produce a flash when you press the flash button, contact the Consumer Care Centre in your country, your Philips dealer or a Philips service centre. |
| The appliance becomes warm during use. | This is normal. | No action required. |
| The appliance does not show light intensity 1 when I switch it on. | The appliance needs to be reset. | To reset the appliance, switch it off and on again. If the appliance still does not show light intensity 1, contact your Customer Care Centre, your Philips dealer or a Philips Service Centre. |
| The appliance produces a strange smell. | The removable window is dirty. | Clean the removable window carefully. If you are unable to clean the removable window properly, contact the Consumer Care Centre in your country, your Philips dealer or a Philips service centre to replace the removable window. |
| | You have not shaved the area to be treated properly. | If there are hairs on the area to be treated, these hairs may get burnt when you use the appliance. As a result, you notice a strange smell. Shave the area to be treated properly before you use the appliance. If shaving causes skin irritation, trim hairs as short as possible and use a comfortable setting. |
| | You stored the appliance in a dusty environment. | The strange smell disappears after a few flashes. |
| The skin feels more sensitive than normal during treatment. | The light intensity you use is too high. | Check if you have selected the right light intensity. If necessary, select a lower light intensity. |
| | You did not shave the areas to be treated. | Shave the areas to be treated before you use the appliance. If shaving causes skin irritation, trim hairs as short as possible and use a comfortable setting. |
| | The filter glass in the light exit window is broken. | If the filter glass in the light exit window is broken, do not use the appliance anymore. Contact the Consumer Care Centre in your country, your Philips dealer or a Philips service centre. |
| | The appliance is defective. | Contact the Consumer Care Centre in your country, your Philips dealer or a Philips service centre. |
| I experience an unacceptable pain sensation when I use the appliance. | You did not shave the areas to be treated | Shave the areas to be treated before you use the appliance. If shaving causes skin irritation, trim hairs as short as possible and use a comfortable setting. |

| Problem | Possible cause | Solution |
|---|---|---|
| | You treated an area for which the appliance is not intended. | Only use the appliance on the legs, underarms and bikini area. |
| | You have used a light intensity which is too high for you. | Reduce the light intensity to a level that is comfortable for you. See chapter 'Preparing for use', section 'Setting the light intensity'. |
| | The filter glass of the light exit window is broken. | If the filter glass of the light exit window is broken, do not use the appliance any more. Contact the Consumer Care Centre in your country, your Philips dealer or a Philips service centre. |
| | The appliance is not suitable for your skin colour. | Do not use the appliance if you have a very dark skin. Do not use it either if your body hair is light blond, red, grey or white. |
| The treated areas become red after the treatment. | Slight redness is harmless and normal and disappears quickly. | No action required. |
| The skin reaction after the treatment lasts longer. | You have used a light intensity which is too high for you. | Select a lower light intensity the next time. See chapter 'Preparing for use', section 'Setting the light intensity'. |
| | | If the skin reaction lasts longer than 3 days, contact your doctor. |
| The hair reduction results are not satisfactory. | You have used a light intensity which is too low for you. | Select a higher light intensity the next time. |
| | You did not make enough overlaps when you used the appliance. | For proper hair-removal results, you have to make overlaps when you use the appliance. See chapter 'Using the appliance', step 7. |
| | You do not use the appliance as often as recommended. | To remove all hairs successfully and to prevent the hair root from becoming active again, the treatments have to be repeated once every two weeks for the first two months. Your skin should be smooth and hair-free by then. To ensure that you keep these results, we advise you to repeat the treatments on the underarms and bikini line every two weeks. On the legs, regrowth is usually slower. Therefore, a treatment once every four weeks is probably sufficient to keep your legs smooth and hair-free. If you see too many hairs grow back within these four weeks, then simply shorten the interval between treatments again, but do not treat more often than once every two weeks. |

| Problem | Possible cause | Solution |
|---|--|--|
| | The appliance is not suitable for your hair or skin colour. | Do not use the appliance if your body hair is light blond, red, grey or white. Do not use the appliance either if you have very dark skin. |
| Hair starts to grow again on some spots of the treated areas. | You did not make enough overlaps when you used the appliance. | For proper hair-removal results, you have to make overlaps when you use the appliance. See chapter 'Using the appliance', step 7. |
| | Hair starts to grow again as part of the natural cycle of hair growth. | This is normal. It is part of the natural cycle of hair growth that hair starts growing again. The hairs that come back, however, tend to be softer and thinner. |

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