Register your product and get support at


[www.philips.com/SatinLux](http://www.philips.com/SatinLux)


or call 00800-S-A-T-I-N-L-U-X (7-2-8-4-6-5-8-9)

SC2000

**PHILIPS**
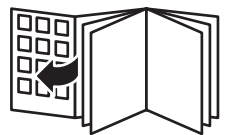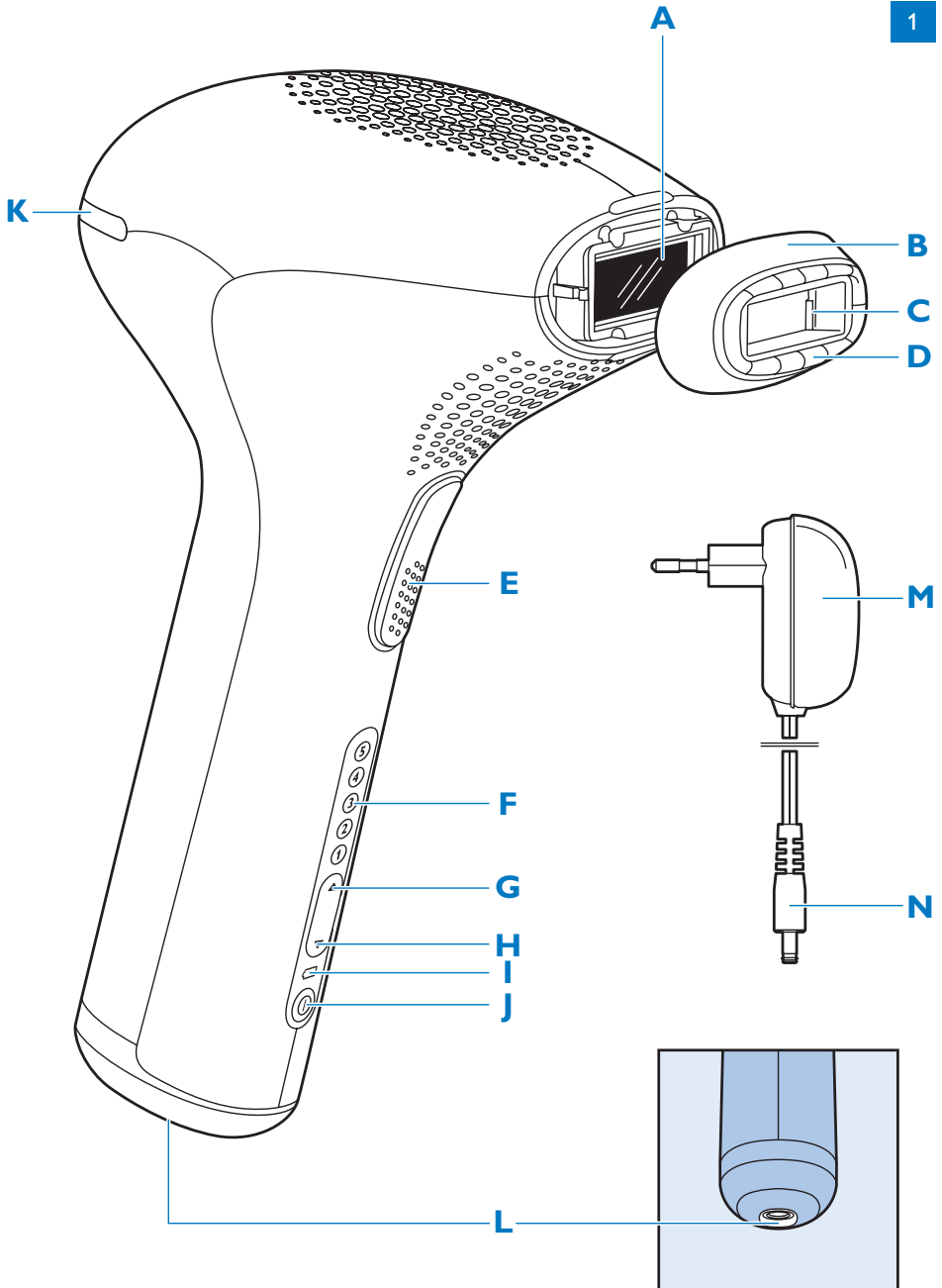
<b>ENGLISH</b>	<b>6</b>
<b>DEUTSCH</b>	<b>22</b>
<b>ESPAÑOL</b>	<b>40</b>
<b>FRANÇAIS</b>	<b>57</b>
<b>ITALIANO</b>	<b>74</b>
<b>NEDERLANDS</b>	<b>92</b>
<b>PORTUGUÊS</b>	<b>110</b>

## Introduction

Congratulations on your purchase and welcome to Philips! To fully benefit from the support that Philips offers, register your product at [www.philips.com/satinlux](http://www.philips.com/satinlux).

Hair removal with light is one of the most effective methods to get rid of unwanted body hair with longer lasting results. It is completely different from current home hair removal methods and Philips has been investigating this method together with leading dermatologists for almost 10 years. You can now use this innovative technology comfortably in the privacy of your own home, whenever you want.

This appliance is intended for removing women's body hair on areas below the neck:

- underarms
- bikini area
- legs

The appliance is not intended for:

- persons with a too light body hair colour (light blond, red, grey or white body hair) or a very dark skin. See chapter 'Preparing for use', section 'Setting the light intensity'.

For more information about the SatinLux, see [www.philips.com/satinlux](http://www.philips.com/satinlux).

## How photo-epilation works

### Hair growth

Depending on age, metabolism and other factors, everyone's hair grows differently, but hair normally grows in three phases:

1 Anagen phase (growing phase)

The hair grows actively from the root. In this phase the concentration of melanin is highest. Melanin is responsible for the pigmentation of the hair: The more melanin there is in the root, the darker the hair is. Melanin concentration is also the determining factor for the efficacy of this method. Only hairs in the anagen phase are susceptible to treatment with light.

2 Catagen phase (degradation phase)

The hair growth stops and the root shrinks before the hair sheds.

3 Telogen phase (resting phase)

The old hair separates from the root and the new hair slowly pushes out the old hair. The root remains at rest until the biological clock tells it to become active again and grow a new hair (4). Hairs on an area of skin are usually in different stages of the cycle of hair growth.


### Working principle

The appliance removes hair by heating up the hair and root beneath the skin. The melanin in the hair and hair root absorbs the released light pulses. The darker the colour of the hair, the more heat can be produced. This process stimulates the hair to go into the resting phase.

The hair now sheds naturally and the hair growth is inhibited.

Optimal hair removal cannot be obtained with a single session, since only hairs in the anagen phase are susceptible to treatment with light. During a treatment, hairs on the treated area are usually in different stages of the cycle of hair growth.

To treat all hairs successfully and to prevent the hair root from becoming active again, the treatments have to be repeated **once every two weeks**.

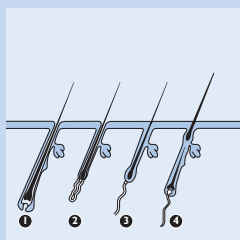An integrated optical filter makes sure no UV light reaches the skin and therefore leaves it unharmed.

### Benefits

This appliance has the following advantages:

#### Professional IPL technology for home use

SatinLux uses a light-based technology called 'Intense Pulsed Light' (IPL). IPL is also used in the professional beauty market for hair removal. SatinLux now brings this innovative technology to the comfort of your home.

#### Effective hair regrowth prevention for everlasting smoothness

This method prevents hair from growing back for a long time and provides a silky-smooth skin. Our studies have shown a significant hair reduction already after 2 treatments. Optimal results were generally achieved after 4 to 5 treatments. A reduction in hair density is visible from treatment to treatment. How fast the effect becomes visible and how long the results last varies from person to person. Repeat this treatment every two weeks to achieve optimal results and to keep your skin smooth for ever.

#### Gentle treatment, even on sensitive body areas

SatinLux has been developed in close cooperation with leading dermatologists for an effective and gentle treatment, even on the most sensitive body areas.

#### Cordless operation for maximum freedom and flexibility

Because the appliance runs on rechargeable batteries, it is portable and can be used anywhere.

### General description (Fig. 1)

- A** Light exit window with filter glass
- B** Removable window
- C** Metallic surface inside the removable window
- D** Safety system (safety ring with contact switches)
- E** Flash button
- F** Intensity lights (1- 5)
- G** ▲ Intensity increase button
- H** ▼ Intensity decrease button
- I** Charging light and battery low indication
- J** On/off button
- K** 'Ready to flash' light
- L** Appliance socket
- M** Adapter
- N** Small plug

Also included: pouch and cleaning cloth

### Important

Read this user manual carefully before you use the appliance and save it for future reference.

### Danger

- Water and electricity are a dangerous combination. Do not use this appliance in wet surroundings (e.g. near a bath, a shower or a swimming pool).
- Keep the appliance and the adapter dry.

### Warning

- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- Only use the appliance in combination with the adapter supplied.
- Do not use the appliance if the adapter or the appliance itself is damaged.
- Do not use the appliance if the filter glass is broken.
- If the adapter is damaged, always have it replaced with one of the original type in order to avoid a hazard.
- The adapter contains a transformer. Do not cut off the adapter to replace it with another plug, as this causes a hazardous situation.
- To prevent damage:
  - 1 never cover the ventilation slots during use.
  - 2 make sure nothing can obstruct the air flow through the ventilation slots
  - 3 never subject the appliance to heavy shocks and do not shake or drop it.
- If you take the appliance from a very cold environment to a very warm environment, wait approximately 3 hours before you use it.
- Never use or store the appliance in a dusty environment.
- Never leave the appliance unattended when it is switched on. Always switch off the appliance after use.
- Do not expose the appliance to direct sunlight or UV light for a longer time.

### Caution

- This appliance is only intended for removing women's body hair. Do not use it for any other purpose.
- The appliance is only intended for removing human hair from areas below the neck. Do not use it for any other body part.
- This appliance is not designed for commercial or professional use.
- If the appliance is used for commercial or professional purposes, if it is used on other body parts than for which it is explicitly intended or if it is used by persons explicitly excluded from use, the guarantee becomes invalid and Philips refuses any liability for damage caused.

- Your skin may become a little red and irritated after using the appliance. This phenomenon is absolutely normal and quickly disappears. If the irritation does not disappear within three days, we advise you to consult a doctor.
- Do not expose the appliance to temperatures lower than 15°C or higher than 35°C during use.
- The appliance is equipped with a built-in overheat protection. If the appliance overheats, it does not flash. Do not switch off the appliance, but let it cool down for approx. 15 minutes before you continue with your treatment.
- Always return the appliance to a service centre authorised by Philips for examination or repair; Repair by unqualified people could cause an extremely hazardous situation for the user.
- Never use the appliance in any of the cases mentioned under Contraindications.

### **Contraindications**

Never use the appliance:

- If you have a skin disease such as active skin cancer; if you have a history of skin cancer or any other localised cancer in the areas to be treated, if you have pre-cancerous lesions or multiple atypical moles in the areas to be treated.
- If you have epilepsy with flashlight sensitivity.
- If you have a history of collagen disorder; including a history of keloid scar formation or a history of poor wound healing.
- If you have a history of vascular disorder; such as the presence of varicose veins or vascular ectasia in the areas to be treated.
- If your skin is sensitive to light and causes a rash or an allergic reaction. If you are taking photosensitising agents or medications, check the package insert of your medicine and never use the appliance if it is stated that it can cause photo-allergic reactions, photo-toxic reactions or if you have to avoid sun when taking this medicine.
- If you have diabetes, lupus erythematoses, porphyria or congestive heart disease.
- On areas of your skin where you are currently being treated with or have recently been treated with Alpha-Hydroxy Acids (AHAs), Beta-Hydroxy Acids (BHAs), topical isotretinoin and azelaic acid.
- If you have taken oral isotretinoin Accutane® or Roaccutane® in the last six months. This treatment can make skin more susceptible to tears, wounds and irritations.
- If you have any bleeding disorder or take anticoagulation medications, including heavy use of aspirin, in a manner which does not allow for a minimum 1-week washout period prior to each treatment.
- If you have infections, eczema, burns, inflammation follicles, open lacerations, abrasions, surgeries, herpes simplex, wounds or lesions and haematomas in the areas to be treated.
- If you have a history of immunosuppressive disease (including HIV infection or AIDS) or when you take immunosuppressive medications.
- When you are on painkillers, which reduce the sensitivity to heat.
- If you are pregnant or breastfeeding.
- On moles, freckles, large veins, darker pigmented areas, scars, skin anomalies without consulting your doctor. This can result in a burn and

change of colour, which makes it potentially harder to identify skin-related diseases (e.g. skin cancer).

- On nipples and mucous areas.
- After intensive sunbathing and/or on sunburnt skin or if you use suntan accelerators.
- On warts, tattoo or permanent make-up. This can result in a burn and change of colour.
- If you use long-lasting deodorants. This can result in skin reactions.
- Over or near anything artificial like silicon implants, Implanon contraceptive implants, pacemakers, subcutaneous injection ports (insulin dispenser) or piercings.

### Electromagnetic fields (EMF)

This Philips appliance complies with all standards regarding electromagnetic fields (EMF). If handled properly and according to the instructions in this user manual, the appliance is safe to use based on scientific evidence available today.

### Charging

Fully charge the batteries before you use the appliance for the first time and when the batteries are empty. Fully charging the batteries takes up to 4 hours.

Charge the appliance when the charging light lights up orange during use. When the charging light goes on, the remaining capacity is at least 30 flashes.

Fully charged batteries provide at least 160 flashes at light intensity 5, which amounts to approx. 15 minutes of cordless operation time. The actual number of flashes varies, depending on the selected light intensity (see 'Efficacy at different light intensities' in chapter 'Preparing for use').

Fully charge the appliance every 3 to 4 months, even if you do not use the appliance for a longer time.

- 1** Switch off the appliance.
- 2** Insert the small plug into the appliance and put the adapter in the wall socket.
  - ▶ The charging light flashes green to indicate that the appliance is charging.
  - ▶ When the batteries are fully charged, the charging light lights up green continuously.

**Caution:** Never cover the appliance and adapter during charging.

*Note:* During charging, the appliance cannot be used.

*Note:* The adapter feels warm during charging. This is normal.

- 3** After charging, remove the adapter from the wall socket and pull the small plug out of the appliance.

### Preparing for use

For optimal results and highest efficacy, prepare the areas on which you intend to use the appliance in the manner described below.

- 1 Shave the areas you intend to treat as long as you still experience hair regrowth.

Never remove the hairs by means of wax, epilation, tweezers or depilation cream. These methods counteract photo-epilation because they remove hairs by the root. In these cases, the light cannot be successfully absorbed by the root and the hairs will grow back again.

*Note: Shaving is necessary to ensure that the light energy is channelled to the hair root in the most efficient way. Furthermore, visible hairs in the treatment area can be burnt by the light and can thus cause an unpleasant smell. In general this is harmless, but it may be a bit painful when you use the appliance.*

- 2 Clean your skin and make sure it is entirely dry and free from grease and oil (e.g. present in cosmetic products like deodorants, lotions, perfumes, tanning creams and sunscreens).

*Note: If shaving causes wounds or lesions, the skin may be more sensitive as a result. In that case, the light treatment may be more perceptible and sometimes even uncomfortable.*

*Note: Once the hairs stop growing back - which usually happens after a few treatments - you no longer have to shave before you use the appliance.*

*Tip: If shaving causes skin irritation, you can also use a trimmer to trim the hair as short as possible. In that case, the light treatment may be more perceptible and sometimes even slightly uncomfortable.*

### **Setting the light intensity**

You have to adjust the light intensity to your skin and body hair colour and to a level you find comfortable.

The table with recommended light intensities below helps you to establish the right light intensity.

- 1 Consult the table below to determine which light intensities are most suitable for your skin and body hair colour and to check if this method is suitable for you at all (if it is not, this is indicated with x in the table).

*Note: Persons with darker skin run a higher risk of skin irritation because their skin absorbs more light. Therefore lower intensities are recommended for persons with darker skin.*

- 2 For each body area you want to treat, select a suitable intensity in the recommended range that is as high as possible but does not cause any discomfort on the skin.

The light intensity you use can feel warm or hot on the skin, but should never feel uncomfortable.

*Note: Before every treatment, check if the appliance is set to the right intensity.*

### Recommended light intensities (1-5)

Skin colour & characteristics :	White (always sunburn, no tanning)	Beige (easy sunburn, minimal tanning)	Light brown (sometimes sunburn, good tanning)	Mid brown (rarely sunburn, fast and good tanning)	Dark brown (rarely sunburn, very good tanning)	Black (rarely to never sunburn, very dark tanning)
Body hair colour: white/ grey	x	x	x	x	x	x
red	x	x	x	x	x	x
light blond	x	x	x	x	x	x
blond	4/5	4/5	4/5	3/4/5	1/2/3/4	x
dark blond/ light brown	4/5	4/5	4/5	3/4/5	1/2/3/4	x
brown	4/5	4/5	4/5	3/4/5	1/2/3/4	x
dark brown	4/5	4/5	4/5	3/4/5	1/2/3/4	x
black	4/5	4/5	4/5	3/4/5	1/2/3/4	x

x means that the appliance is not suitable for use.

*Note: The appliance is not suitable for persons with very dark to black skin because of the higher risk of skin irritations.*

*Note: The appliance is not suitable for persons with light blond, red, grey and white body hair because the low melanin concentration in the hair results in very poor hair removal efficacy.*

*Note: At higher light intensities the method is not only more efficacious, but it is also more perceptible.*

#### Efficacy at different light intensities

Light intensity	Efficacy	Number of flashes
5	*****	160
4	****	190
3	***	250
2	**	350
1	*	520

\* indicates the level of efficacy

*Note: This table shows the number of flashes the appliance can give at each light intensity and the corresponding expected efficacy. Please use the intensity appropriate for your skin colour and body hair colour.*

## Using the appliance

The light intensity you use can feel warm or hot on the skin, but should never feel uncomfortable.

This appliance is only intended for removing women's body hair on areas below the neck: underarms, bikini area and legs.

*Note: Shave the area you intend to treat before you use the appliance (see chapter 'Preparing for use').*

*Tip: The best time for treatment is before you go to bed at night. This diminishes any skin reaction that may occur overnight.*

**1** Press the on/off button to switch on the appliance.

Intensity light 1 goes on automatically to indicate that the appliance is switched on at default intensity 1.

*Note: You can switch off the appliance any time by pressing the on/off button again.*




*Note: It is normal that the appliance becomes warm during use.*


**2** Press the ▲ button to increase the intensity. To lower the intensity, press the ▼ button.

Every time you press the button, the corresponding intensity light starts to flash. This takes a few seconds.

**3** Place the appliance perpendicularly on the skin so that the removable window and the safety ring are in contact with the skin.


*Tip: When you use the appliance for the first time, try it out, for instance on your hip, to get accustomed to the light-based method. Start with the lowest recommended intensity and increase it until you have reached the highest recommended intensity. This test tells you how sensitive your skin is to light. Check if your skin shows any reactions 20 minutes after the test. If it does, use a lower recommended intensity for the areas to be treated.*

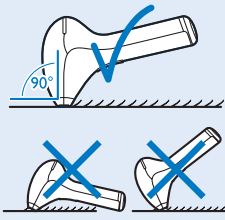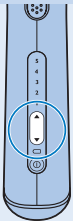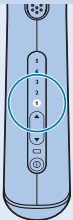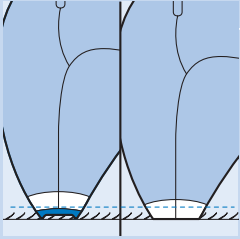**4** Press the safety ring fully onto the skin with slight pressure.

*Note:* The safety ring has contact switches that together are the safety system of the appliance. This safety ring prevents unintentional flashing without skin contact.

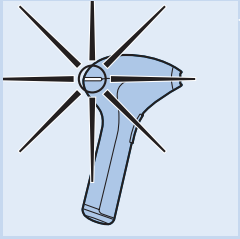**D** When all contact switches of the safety ring have been pushed into the appliance, the 'ready to flash' light on the back of the appliance goes on. This may take a few seconds. When the 'ready to flash' light is on, the appliance is ready to release a flash.


*Tip:* For easier use, stretch the skin you want to treat. On areas with softer skin, you may need to exert more pressure on the appliance to push in all contact switches.

**5** Press the flash button to release a flash. After the appliance has flashed once, release the flash button.

**6** For the next flash, repeat steps 3, 4 and 5. This can take a few seconds.

Never flash the same area twice without lifting the appliance off the skin. This prevents excessive heating and possible burning of the skin.

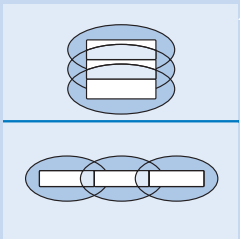**7** To avoid untreated areas, always make sure there is some overlap with the previously treated area when you place the appliance on the skin. The effective light only comes out of the light exit window. Make sure the flashes are produced close to each other.

**8** Switch off the appliance after you have finished the treatment.

*Note:* Avoid looking straight into the light coming out of the window when the appliance is switched on.

*Note:* The visible light produced by the appliance is the reflection of the flash on the skin and is harmless to your eyes. It is not necessary to wear goggles during use.

**Guidelines for treatment**

**Recommended number of flashes and treatment time**

Area	Approx. number of flashes	Approx. treatment time
Underarms	25 per underarm	3 minutes per underarm
Bikini line	25 per side	3 minutes per side
Entire bikini area	90 for entire area	10 minutes
One lower leg	160-190	15 minutes

Area	Approx. number of flashes	Approx. treatment time
One full leg	320-380	30 minutes

*Note: Fully charged batteries offer at least 160 flashes at intensity 5, which amounts to a cordless operating time of approx. 15 minutes.*

### How to achieve optimal results

Optimal results cannot be obtained with a single session, since only hairs in the anagen phase are susceptible to treatment with light. During a treatment, hairs on a treated area are usually in different stages of the cycle of hair growth. To treat all hairs successfully and to prevent that the hair root becomes active again, you have to repeat the treatment every 2 weeks. It is recommended to adhere to this frequency. Better results are not expected if you increase or decrease the treatment frequency. Clinical trials have shown that this is the optimal frequency of use.

### After treatment

#### After treatment it may happen that:

- Your skin shows slight redness and/or prickles, tingles or feels warm. This reaction is harmless and disappears quickly.
- A skin reaction similar to sunburn may occur on the treated area for a short while.
- Your skin may appear to be dry, as the combination of shaving and light treatment may cause dry skin.
- Some hairs grow back although you use the appliance frequently. This happens because hairs are in different stages of the cycle of hair growth during treatment. The hairs that grow back tend to be softer and thinner.

#### It is recommended to:

- Consult your doctor if skin irritation does not disappear within three days (see chapter 'Important', section 'Caution').
- Wait with the next treatment until skin irritation has disappeared completely. Perhaps you have used an intensity that was too high for you. Use a lower intensity the next time.
- Do not take an intensive sunbath or suntan without proper protection immediately after a treatment. If you want to go sunbathing after a treatment, apply a sun blocker to the treated areas.
- Avoid swimming in chlorinated water for 48 hours after treatment. The chlorine in the water can cause itchiness.
- Avoid applying any cosmetic substances to the treated areas immediately after treatment.


## Cleaning and maintenance

Never clean the appliance or any of its parts under the tap or in the dishwasher.

Never use scouring pads, abrasive cleaning agents or aggressive liquids such as petrol or acetone to clean the appliance.

Never scratch the filter glass or the metallic surface inside the removable window.

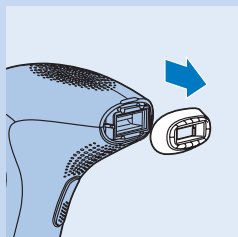*Note: Make sure that the contact switches of the safety ring do not get jammed with dirt.*

*Note: Stop using the appliance when it is no longer possible to clean the filter glass or the removable window.*

To ensure optimal results and a long life for the appliance, clean the appliance after each treatment.

- 1 Switch off the appliance, unplug it and let it cool down.

*Note: The filter glass becomes hot during use. Make sure it has cooled down before you clean it.*

- 2 To remove the removable window, put your fingers in the recesses on the housing of the appliance and pull the removable window off the appliance.

- 3 Use the dry, soft cloth supplied with the appliance to clean the following parts:


- the surface of the filter glass
- the outside surface of the removable window
- the metallic frame inside the removable window
- the housing of the appliance (if necessary)

*Note: If necessary, wet the cloth with some drops of water.*

### Storage

- 1 Switch off the appliance, unplug it and let it cool down.
- 2 Clean the appliance before you store it.
- 3 Put the appliance in the pouch supplied.
- 4 Store the appliance in a dry place at a temperature between 0°C and 60°C.

### Environment

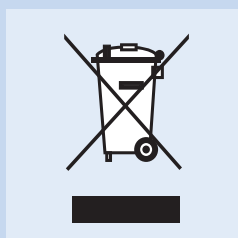
- Do not throw away the appliance with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment.
- The built-in rechargeable batteries contain substances that may pollute the environment. Always remove the batteries before you discard and hand in the appliance at an official collection point. Dispose of the batteries at an official collection point for batteries. If you have trouble removing the batteries, you can also take the appliance to a Philips service centre. The staff of this centre will remove the batteries for you and will dispose of them in an environmentally safe way.


### Removing the rechargeable batteries

Only remove the rechargeable batteries if they are completely empty.


- 1 Insert an awl or another pointed tool into the appliance tool and pry the bottom cap off the handle.

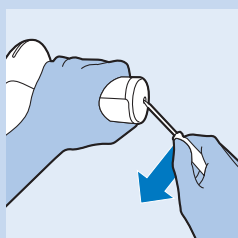
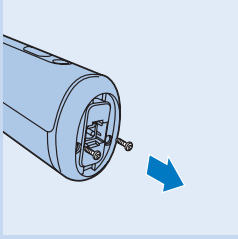
**2** Undo the two screws.

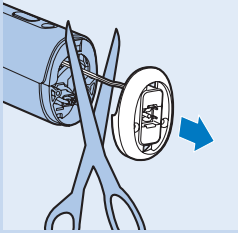
**3** Remove the bottom part of the handle and cut the two wires that connect the bottom part to the appliance.

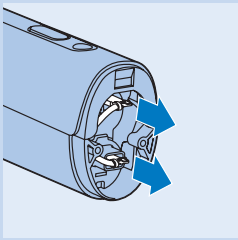**4** Pull the battery connectors off the battery terminals.

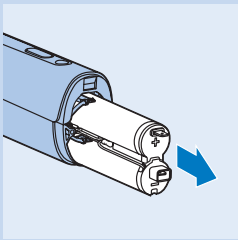**5** Shake the rechargeable batteries out of the appliance.

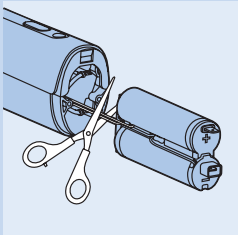**6** Cut the two wires that connect the rechargeable batteries to the appliance.

### Guarantee and service

If you need service or information or if you have a problem, please visit the Philips website at [www.philips.com/satinlux](http://www.philips.com/satinlux) or call 00800-SATINLUX (72846589).

### Replacement parts

The following replacement parts are available:

- Removable window
- Adapter
- Pouch

### Technical specifications

#### Model SC2000

Rated voltage	100V-240V
Rated frequency	50Hz-60Hz
Rated input	7.5W
Protection against electric shock	Class II 
Protection rating	IP 30 (EN 60529)
Operating conditions	Temperature: +15° to +35°C
	Relative humidity: 25% to 75%
Storage conditions	Temperature: 0° to +60°C
	Relative humidity: 5% to 95%
Intense pulsed light lamp	Spectrum: >570nm
Lithium ion battery	2 x 3.7 Volt
	1100 mAh

### Troubleshooting

This chapter summarises the most common problems you could encounter with the appliance. If you are unable to solve the problem with the information below, contact the Consumer Care Centre in your country. See the telephone number on the back of the user manual or visit [www.philips.com/satinlux](http://www.philips.com/satinlux) for further information.

Problem	Possible cause	Solution
The appliance does not work.	The rechargeable batteries are empty.	Charge the appliance (see chapter 'Charging').
	The appliance is broken.	Contact the Consumer Care Centre in your country, your Philips dealer or a Philips service centre.
The appliance suddenly switches off.	The rechargeable batteries are empty.	Charge the appliance (see chapter 'Charging').

Problem	Possible cause	Solution
The adapter becomes warm during charging.	This is normal.	No action required.
The charging light does not go on when I put the small plug in the appliance socket.	You have not put the adapter in the wall socket.	Insert the small plug into the appliance socket and put the adapter in the wall socket.
		Connect another appliance to the wall socket to check if the wall socket is live.
		If the appliance still does not charge, contact the Consumer Care Centre in your country, your Philips dealer or a Philips service centre.
	The adapter is not inserted properly into the wall socket and/or the small plug is not inserted properly into the appliance socket.	Make sure the adapter is inserted properly into the wall socket and/or the small plug is inserted properly into the appliance socket.
	The appliance is broken.	Contact the Consumer Care Centre in your country, your Philips dealer or a Philips service centre.
I have switched on the appliance but I cannot increase or decrease the light intensity.	The appliance has to be reset.	To reset the appliance, switch off the appliance and then switch it on again. If you still cannot adjust the light intensity, contact the Consumer Care Centre in your country, your Philips dealer or a Philips service centre.
The 'ready to flash' light does not go on when I place the appliance on the skin.	You have not put the removable window with safety ring properly on the skin.	Place the appliance perpendicularly onto the skin in such a way that all contact switches touch the skin. Then check if the 'ready to flash' light is on and press the flash button.
		If this does not work, place the appliance on a part of your body where it is easy to establish full skin contact, for example on your lower arm. Then check if the 'ready to flash' light works. If the 'ready to flash' light still does not go on when you place the appliance on the skin, contact the Consumer Care Centre in your country, your Philips dealer or a Philips service centre.
	The safety ring is dirty.	Clean the safety ring carefully. If you are unable to clean the safety ring properly, contact the Consumer Care Centre in your country, your Philips dealer or a Philips service centre to replace the flash screen.

Problem	Possible cause	Solution
The 'ready to flash' light does not go on when I place the removable window on my skin, but the cooling fan works.	The overheat protection has been activated.	When the overheat protection has been activated, the fan still works. Do not switch off the appliance and let it cool down for approx. 15 minutes before you continue to use it. If the 'ready to flash' light still does not go on when you place the appliance on the skin, contact the Consumer Care Centre in your country, your Philips dealer or a Philips service centre.
The 'ready to flash' light is green but the appliance does not produce a flash when I press the flash button.	The appliance needs to be reset.	To reset the appliance, switch off the appliance and switch it on again. If the 'ready to flash' light is green but the appliance still does not produce a flash when you press the flash button, contact the Consumer Care Centre in your country, your Philips dealer or a Philips service centre.
The appliance becomes warm during use.	This is normal.	No action required.
The appliance does not show light intensity 1 when I switch it on.	The appliance needs to be reset.	To reset the appliance, switch it off and on again. If the appliance still does not show light intensity 1, contact your Customer Care Centre, your Philips dealer or a Philips Service Centre.
The appliance produces a strange smell.	The removable window is dirty.	Clean the removable window carefully. If you are unable to clean the removable window properly, contact the Consumer Care Centre in your country, your Philips dealer or a Philips service centre to replace the removable window.
	You have not shaved the area to be treated properly.	If there are hairs on the area to be treated, these hairs may get burnt when you use the appliance. As a result, you notice a strange smell. Shave the area to be treated properly before you use the appliance.
The treatment is more perceptible than usual.	You used a too high light intensity.	Check if you have selected the right light intensity. If necessary, select a lower light intensity.
	You did not shave the areas to be treated.	Always shave the areas to be treated before you use the appliance.
	The appliance is defective.	Contact the Consumer Care Centre in your country, your Philips dealer or a Philips service centre.
	You stored the appliance in a dusty environment.	The strange smell disappears after a few flashes.
	The filter glass in the light exit window is broken.	If the filter glass in the light exit window is broken, do not use the appliance anymore. Contact the Consumer Care Centre in your country, your Philips dealer or a Philips service centre.

Problem	Possible cause	Solution
I experience an unacceptable pain level when I use the appliance.	You did not shave the areas to be treated	Always shave the area to be treated before you use the appliance.
	You treated an area for which the appliance is not intended.	Only use the appliance on the legs, underarms and bikini area.
	You used a too high light intensity.	Reduce the light intensity to a level that is comfortable for you. See chapter 'Preparing for use', section 'Setting the light intensity'.
	The filter glass of the light exit window is broken.	If the filter glass of the light exit window is broken, do not use the appliance any more. Contact the Consumer Care Centre in your country, your Philips dealer or a Philips service centre.
	The appliance is not suitable for your hair or skin colour.	Do not use the appliance if your body hair is light blond, red, grey or white. Do not use the appliance either if you have very dark skin.
The treated areas become red after the treatment.	Slight redness is harmless and normal and disappears quickly.	No action required.
The skin reaction after the treatment lasts longer.	You used a too high light intensity.	Select a lower light intensity the next time. See chapter 'Preparing for use', section 'Setting the light intensity'.
		If the skin reaction lasts longer than 3 days, contact your doctor.
The hair reduction results are not satisfactory.	You used a too low light intensity.	Select a higher light intensity the next time.
	You did not let the appliance give enough flashes.	Make sure you let the appliance give enough flashes. See chapter 'Preparing for use', section 'Setting the light intensity'.
	You did not make enough overlaps when you used the appliance.	To obtain proper hair-removal results, you have to make overlaps when you use the appliance. See chapter 'Using the appliance', step 7.
	You do not use the appliance often enough.	Use the appliance once every 2 weeks.
	The appliance is not suitable for your hair or skin colour.	Do not use the appliance if your body hair is light blond, red, grey or white. Do not use the appliance either if you have very dark skin.
Hair starts to grow again on some spots of the treated areas.	You did not make enough overlaps when you used the appliance.	For proper hair removal results, you have to make overlaps when you use the appliance. See chapter 'Using the appliance'.
	Hair starts to grow again as part of the natural cycle of hair growth.	It is part of the natural cycle of hair growth that hair starts growing again. The hairs that come back, however, tend to be softer and thinner.

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