Always here to help you

Register your product and get support at
www.philips.com/welcome


SC198X

User manual

PHILIPS

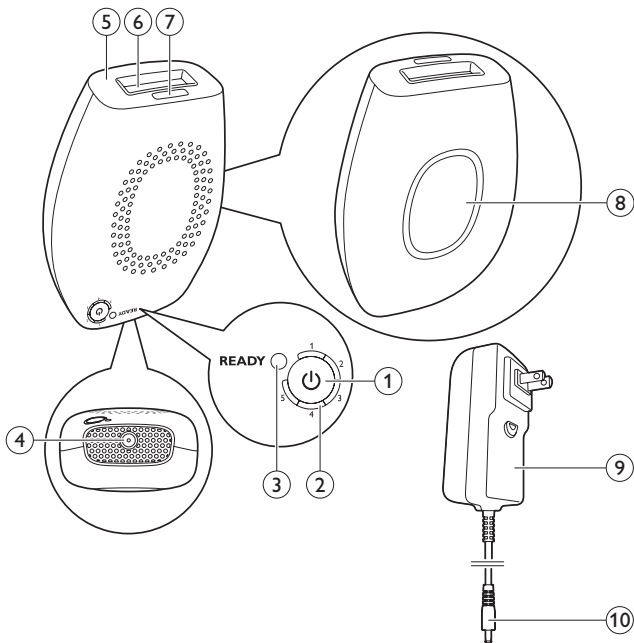SC198X

ENGLISH 8

ESPAÑOL 37

Table of contents

Introduction	8
What is Philips Lumea Comfort?	9
Intended use of Philips Lumea Comfort	9
Safety with Philips Lumea Comfort	10
Contraindications	10
When to avoid using Philips Lumea Comfort	12
Precautions - How to use Philips Lumea Comfort safely	14
Reducing the risk of injury	16
Risks of using Philips Lumea Comfort	18
Treatment plan and effectiveness	21
How to use Philips Lumea Comfort	23
Cleaning and Maintenance of Philips Lumea Comfort	28
Disposal	29
Troubleshooting	29
Assistance	30
Full Two-Year Warranty	31
45-Day Money-Back Guarantee	32
Specifications	33
Labeling	35

Introduction

Congratulations on your purchase and welcome to Philips! To benefit fully from the support that Philips offers, register your product at www.philips.com/welcome.

Before you use Philips Lumea Comfort for the first time, please read the whole user manual. Pay special attention to all warnings and safety precautions and follow them strictly during use.

What is Philips Lumea Comfort?

Philips Lumea Comfort is a light-based device for long-term hair removal, designed for home use.

Device description (Fig. 1)

- 1 Intensity button
- 2 Intensity lights
- 3 READY/skin color warning light
- 4 Device socket
- 5 Skin contact surface
- 6 Light exit window
- 7 Skin color sensor
- 8 Flash button
- 9 Adapter
- 10 Small plug of AC adapter cord

Box contents

When you open the box, you find it contains the following parts:

- 1 Philips Lumea Comfort device
- 2 AC adapter cord
- 3 This user manual and a quick start guide

Intended use of Philips Lumea Comfort

Philips Lumea Comfort is intended for the removal of unwanted hair. Philips Lumea Comfort is also intended for permanent reduction in hair regrowth, defined as a long-term stable reduction in the number of hairs regrowing when measured at 6, 9, and 12 months after a treatment regime.

Philips Lumea Comfort may be used to remove unwanted body hair. Ideal body areas for use of Philips Lumea Comfort include the underarms, bikini line, arms and legs.

Safety with Philips Lumea Comfort

With Philips Lumea Comfort safety comes first.

IPL technology offers superior safety with lower energy level

Intense Pulsed Light technology is able to achieve long-term hair removal results at a fraction of the energy level used in other light-based hair removal equipment. The low energy used in Philips Lumea Comfort reduces its potential to cause harm or complications, and contributes to your overall safety.

Philips Lumea Comfort protects your eyes

The Philips Lumea Comfort has a **built-in safety feature** for eye protection. It has been designed so that a flash cannot be emitted when the treatment surface is facing open air. The flash button only works when the treatment surface is in full contact with the skin.

Philips Lumea Comfort protects your skin

Philips Lumea Comfort comes with a skin color sensor that only enables use on lighter skin colors to ensure skin safety.

The 0.35 × 1.18in (0.9 × 3cm) light exit window through which flashes are delivered is recessed in the treatment surface. This enables Philips Lumea Comfort to protect your skin by avoiding direct contact between the light exit window and the skin.

Contraindications

IMPORTANT SAFETY INFORMATION - Read Before Use!

Philips Lumea Comfort is not designed for everyone. Please read and consider the information in the following section before use. For further information and personalized advice, you may also visit www.philips.com/lumea.

If you are unsure Philips Lumea Comfort is safe for you to use, please consult your physician or dermatologist.

DO NOT USE Philips Lumea Comfort if you have a naturally dark skin color!

Philips Lumea Comfort removes unwanted hair by selectively addressing hair pigment. Varied quantities of pigment also exist in the surrounding tissue of skin. The quantity of pigment in a particular person's skin, which is manifested by their skin color, determines the degree of risk they are exposed to when they use this device.

Treating dark skin with Philips Lumea Comfort can result in adverse effects such as burns, blisters, and skin color changes (hyperpigmentation or hypopigmentation).

- **DO NOT USE Philips Lumea Comfort on naturally dark skin!**

A unique skin color sensor is incorporated in Philips Lumea Comfort. It measures the color of the skin to be treated at the beginning of each session and occasionally during the session. This skin color sensor ensures that a flash will be emitted only on suitable skin colors.

DO NOT USE Philips Lumea Comfort on tanned skin or after recent sun exposure!

Tanned skin, particularly after sun exposure, contains large quantities of the pigment melanin. This applies to all skin types and colors, including those that do not seem to tan quickly. The presence of large quantities of melanin exposes the skin to higher risk when using Philips Lumea Comfort.

Using Philips Lumea Comfort on skin that has been recently exposed to sunlight can result in adverse effects such as burns, blisters, and skin color changes (hyperpigmentation or hypopigmentation).

- **DO NOT USE Philips Lumea Comfort on tanned skin or after sun exposure! Such use can cause burns or skin injury.**
- **Avoid exposure to the sun for 4 weeks before your treatment with Philips Lumea Comfort!**

Never use Philips Lumea Comfort around or near the eyes - Philips Lumea Comfort is not recommended for facial treatment.

Philips Lumea Comfort is not effective on naturally white, gray, blond and red body hair. If your body hair is of these colors, the device does not work on you. Hair removal by lasers or intense pulsed light can cause increased hair growth in some individuals. Based upon currently available data, the highest risk for this response occurs in females of Mediterranean, Middle-Eastern, and South-Asian descent, when they use laser or intense pulsed light on their face or neck.

When to avoid using Philips Lumea Comfort

Certain conditions may prevent the use of Philips Lumea Comfort temporarily. **DO NOT USE THE DEVICE** if any of the following currently applies to you:

- If you are pregnant or nursing (lactating).
- If you were exposed to strong sunlight or an artificial tanning machine during the past 4 weeks.
- If you have a tattoo or permanent makeup on the area to be treated.
- If you have dark spots, such as large freckles, birthmarks, moles or warts in the area to be treated.
- If you have eczema, psoriasis, lesions, open wounds or active infections, such as cold sores, in the area to be treated. Wait for the area to heal before you use the device.

- If you have a history of keloidal scar formation, a known sensitivity to light (photosensitivity) or are taking medication that makes the skin more sensitive to light, including non-steroidal anti-inflammatory agents, (e.g., aspirins, ibuprofens, acetaminophen), tetracyclines, phenothiazines, thiazide, diuretics, sulfonluraes, sulfonamides, DTIC, fluorouracil, vinblastine, griseofulvin, Alpha-Hydroxi Acids (AHAs), Beta-Hydroxi Acids (BHAs), Retin-A[®], Accutane[®] and/or topical retinoids.
- If you have abnormal skin conditions due to diabetes, for example. or another systemic or metabolic disease.
- If you are currently or have recently been treated with Alpha-Hydroxi Acids (AHAs), Beta-Hydroxi Acids (BHAs), Retin-A[®], Accutane[®] (isotretinoin), topical retinoids, or azelaic acid.
- If you have been treated with Accutane[®] (isotretinoin) within the past 6 months.
- If you have been on a steroid regimen within the past 3 months.
- If you have a history of herpes outbreaks in the treatment area, unless you have consulted your physician and received preventive treatment before you start using the device.
- If you suffer from epilepsy.
- If you have an active implant, such as a pacemaker, incontinence device, insulin pump, etc.
- If you have a disease related to photosensitivity, such as porphyria, polymorphic light eruption, solar urticaria, lupus, etc.
- If you have a history of skin cancer or potential skin malignancies.
- If you have received radiation therapy or chemotherapy treatments within the past 3 months.

- If you have any other condition which in your physician's opinion would make it unsafe for you to be treated.

If you are unsure Philips Lumea Comfort is safe for you to use, please consult your physician or dermatologist.

Precautions - How to use Philips Lumea Comfort safely

Choose the intensity level carefully

Intensity level refers to the intensity of the flash emitted onto your skin during use, from the lowest level to the highest level. The intensity lights around the intensity button indicate the intensity level at which the device is set. As the intensity level increases, so do the results as well as the risk of side effects (see chapter 'Risks of using Philips Lumea Comfort').

Always start using Philips Lumea Comfort at the lowest intensity level!

Only if you experience little or no discomfort during and after use of Philips Lumea Comfort at the lowest intensity level, raise the intensity by one level the next time you use the device and so on for each subsequent session.

For detailed instructions on setting the intensity level, see section 'Treating with Philips Lumea Comfort' in chapter 'How to use Philips Lumea Comfort'.

Avoid adverse effects!

- Do not treat the same area of skin more than once during a session!
- Avoid overlapping flashes!
- If your skin blisters or burns, **STOP USING THE DEVICE IMMEDIATELY!**

Avoid complications after use!

- Do not expose treated areas of skin to the sun. Sun exposure includes constant unprotected exposure to direct sunlight for over 15 minutes, or constant unprotected exposure to diffuse sunlight for over 1 hour. To protect recently treated skin when exposed to sunlight, be sure to apply a sunscreen with SPF 30 or higher thoroughly for 2 weeks after each hair removal session.
- Always shave the area to be treated and make sure that the skin is clean and dry before you use the device.
- Cover birthmarks and tattoos before treatment.
- Cover dark brown or black spots, such as large freckles, birthmarks, moles or warts before treatment.
- Never look directly at light coming out of the light exit window of the device.
- Do not use the device on nipples and genitals (male or female).
- Do not use the device on any body site where you might later want hair.
- Do not use the device for any other purpose than hair removal.
- Never emit a flash into open space. Always make sure the skin contact surface is in full contact with the skin during treatment.
- Remove the skin contact surface from the skin if either the skin or the surface is too hot.
- Never use flammable liquids such as alcohol (including perfumes, sanitizers, or other substances containing alcohol) or acetone to clean the skin before you use the device.

- Use of the device may cause temporary pigmentation changes (see chapter 'Risks of using Philips Lumea Comfort').
- Keep this device out of the reach of children. Do not use the device on children or allow children to use it.

Reducing the risk of injury

As with any electrical device, certain precautions must be taken to ensure your safety when using Philips Lumea Comfort.

- **No modification of the device is allowed!**

- **Keep Philips Lumea Comfort away from water!**

- Keep Philips Lumea Comfort away from water. This is an electrical device. As such, it should always be kept away from water. Do not place or store Philips Lumea Comfort where it can fall or be pushed into a tub, sink or any other vessel containing water. Do not place in, or drop into water or any other liquid. This may cause electrocution!
- Do not use Philips Lumea Comfort while bathing.
- Do not use Philips Lumea Comfort if it becomes damp or wet.
- Do not reach for Philips Lumea Comfort if it has fallen into water. Unplug Philips Lumea Comfort immediately if it has fallen into water.
- **Never open Philips Lumea Comfort!**
- Do not attempt to open or repair Philips Lumea Comfort. Opening the device may expose you to dangerous electrical components and to intense light, either of which may cause bodily damage and/or eye injury.

- **Do not attempt to open or repair your Philips Lumea Comfort. For assistance call 1-800-243-7884.**
- If you try to open Philips Lumea Comfort, you may also damage the device. Opening the device will void your warranty. Please visit www.philips.com/support or call 1-800-234-7884 for assistance if you have a broken or damaged device in need of repair.
- Use Philips Lumea Comfort only for its intended purpose and as described in this user manual.
- Never leave the device unattended when plugged into an electrical outlet.
- Do not operate the device if it has a damaged cord or adapter. Keep the cord away from heated surfaces.
- Do not use the device if you see or smell smoke during use.
- Do not use the device if it does not work properly or if it appears damaged.
- Do not use the device if the fan vent in the skin contact surface is cracked, detached or missing.
- Do not use the device if the skin color sensor is cracked or broken.
- Do not use the device if the outer shell is cracked or coming apart.
- Always unplug the device from the electrical outlet immediately after use.
- Always unplug the device before cleaning.
- Do not use the device with any attachments or accessories not recommended by Philips.
- Philips Lumea Comfort must be operated only with its designated adapter. For details on the right adapter, see chapter 'Specifications'.
- This device is not suitable for use in the presence of flammable anaesthetic mixtures with air, oxygen or nitrous oxide.

Risks of using Philips Lumea Comfort

Side effects and complications associated with the use of Philips Lumea Comfort are uncommon. However, every cosmetic procedure, including those designed for home use, involves some degree of risk. Therefore, it is important that you understand and accept the risks and complications that can occur with intense pulsed light hair removal systems designed for home use.

Adverse event	Likelihood of adverse event	Adverse effect	Likelihood of the adverse effect
Stacking or overlapping of multiple flashes	Minor	Minor skin discomfort	minor
		Skin redness	minor
		Increased skin sensitivity	minor
		Skin wounds and burns	rare
		Scarring	rare
		Pigment changes	rare
		Excessive redness and swelling	rare
		Infection	negligible
		Bruising	negligible

Minor skin discomfort

Although intense pulsed light hair removal is generally very well-tolerated, most users feel some mild discomfort during use. This minor discomfort is usually described as being a mild stinging sensation on the treated skin areas. The stinging sensation usually lasts during release of the flash or for a few minutes after release of the flash. Anything beyond this minor discomfort is abnormal and means that either you should not continue to use the device because you are unable to tolerate its hair removal method, or that the intensity level is set too high.

Skin redness

Your skin may become red right after using Philips Lumea Comfort or within 24 hours of using Philips Lumea Comfort. Redness generally clears up within 24 hours. Consult your physician if redness does not go away within 2 to 3 days.

Increased skin sensitivity

The skin of the treated area is more sensitive so you may encounter dryness or flaking of the skin.

Skin wounds and burns

Very rarely, burns or wounds to the skin can occur following use of the device. The burn or wound can require a few weeks to heal and, extremely rarely, may leave a noticeable scar.

Scarring

Although very rare, scarring may occur. Usually when scarring occurs, it is in the form of a flat and white lesion on the skin (hypotrophic). However, it can be large and red (hypertrophic) or large and extend beyond the margins of the injury itself (keloid). Subsequent aesthetic treatments may be required to improve the appearance of the scar.

Pigment changes

Philips Lumea Comfort targets the hair shaft, in particular the pigmented cells in the hair follicle and the hair follicle itself. Nevertheless, there is risk of temporary hyperpigmentation (increased pigment or brown discoloration) or hypopigmentation (whitening) to the surrounding skin. This risk of changes in skin pigmentation is higher for people with darker skin color.

Usually discoloration or changes to skin pigment are temporary and hyperpigmentation or hypopigmentation rarely occur.

Excessive redness and swelling

In rare cases, treated skin may become very red and swollen. This is more common in sensitive areas of the body. The redness and swelling should subside within 2 to 7 days. Application of ice may help to reduce the swelling. Gentle cleansing of the skin is allowed, but exposure to the sun should be avoided.

Infection

Infection of the skin is exceedingly rare. However, it is a possible risk following a skin burn or wound caused by treatment with Philips Lumea Comfort.

Bruising

Very rarely, use of Philips Lumea Comfort may cause blue-purple bruising that may last 5 to 10 days. As the bruising fades, there may be a rust-brown discoloration of the skin (hyperpigmentation).

Treatment plan and effectiveness

A typical full hair growth cycle may take 18-24 months. During this time, multiple treatments with Philips Lumea Comfort may be required in order to achieve permanent hair reduction.

The effectiveness of hair removal varies from person to person according to body area, hair color and how the Philips Lumea Comfort is used.

Treatment plan

The typical treatment plan during a full hair growth cycle:

- Treatments 1 - 4 : Plan 2 weeks apart
- Treatments 5 - 7: Plan 4 weeks apart
- Treatments 8 and after: Treat as needed until decided results are achieved.

Expected treatment outcome for hair removal

Skin type	Hair color	Light fluence [J/cm ²]	Anatomic area of the body	Average number of treatments	Effectiveness (% of hair removal/regrowth)
I-IV	Light brown to brown	max. 4	Lower or upper legs	Session 10-12, approx. 4 weeks apart	60% hair reduction
			Arms		50% hair reduction
			Bikini line		60% hair reduction

Skin type	Hair color	Light fluence [J/cm ²]	Anatomic area of the body	Average number of treatments	Effectiveness (% of hair removal/regrowth)
			Underarms		55% hair reduction
I-IV	Brown to black	max. 4	Lower or upper legs	Session 8-10, approx. 4 weeks apart	70% hair reduction
			Arms		60% hair reduction
			Bikini line		70% hair reduction
			Underarms		65% hair reduction

Notes

When looking at the tables above, please bear in mind the following:

- Individual response depends on hair type as well as biological factors that may affect hair growth patterns. Some users may respond quicker or slower than the average number of treatments.
- You cannot expect permanent hair reduction in one or two treatment sessions. The duration of the resting period for hair follicles also depends on the body area.

How to use Philips Lumea Comfort

Philips Lumea Comfort setup

Remove the Philips Lumea Comfort device and other parts from the box.

Preparing the adapter (specific types only)





- 1 Pull down the open switch and hold it.
- 2 Attach the plug in the correct direction.

- 3 Release the open switch. The plug is attached now.

Note: To release the plug, pull down the switch and hold it. You can now release the plug.


- 4 Plug the AC cord into the device socket of the Philips Lumea Comfort device.
- 5 Plug the adapter into an electrical outlet.

Your Philips Lumea Comfort is now ready to start.

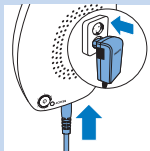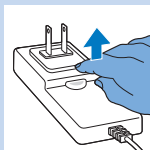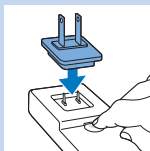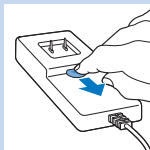
Treating with Philips Lumea Comfort for the first time

Use Philips Lumea Comfort only for its intended purpose and as described in this user manual. Make sure that your skin is shaved, clean, dry and free of any powders, creams, antiperspirants or deodorants.

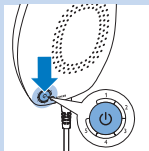
- 1 Press the intensity button briefly. The green intensity level 1 light turns on.

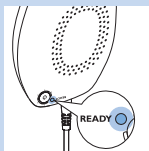- 2 About 1 second later, the yellow READY light turns on. The device is now ready to deliver a first flash at the lowest intensity level.

- 3 If you wish to change the intensity level, press the intensity button quickly and repeatedly to increase the intensity level from 1 to 5 and from 5 directly back to 1, until you reach the desired level. The green level lights indicate the level selected (see section 'Intensity level' below).

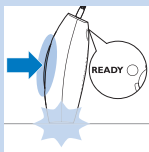- 4 Place the skin contact surface on the skin. Make sure that the skin is spread evenly and smoothly and the skin contact surface is in full contact with the skin. Make sure that the yellow READY light starts blinking slowly.

- 5 Press the flash button. The device first determines the color of your skin. If the color of the skin is light enough for safe treatment, the device emits a flash of light and produces a soft popping sound at the same time. This popping sound is normal for this device. You may also feel a mild warm and tingling sensation.

If you see no flash and the yellow READY light blinks rapidly, the skin color warning light warns you that your skin color is too dark for safe treatment. Try using the device on a different body part. For assistance, please visit www.philips.com/support or call 1-800-243-7884.

- 6** After emitting a flash, Philips Lumea Comfort recharges immediately for the next flash. Recharging takes a few seconds, depending on the intensity level set. When recharging is finished, the READY light turns on again.
- 7** Move the skin contact surface to another area of skin. When there is full skin contact, the READY light starts blinking slowly again.

Be careful not to treat the same area of skin more than once per treatment session! Treating the same area of skin more than once per session increases the likelihood of adverse effects.

- 8** Repeat the process, starting at step 5.

Intensity level

The intensity level setting determines the intensity of the flash delivered to your skin. There are 5 levels, from the lowest (1) to the highest (5).

Corresponding intensity lights show the intensity selected. As the intensity level increases, so do the hair removal results and the risk of possible side effects and complications.

Always start using Philips Lumea Comfort at the lowest intensity level.

If you experience little or no discomfort during and after your first hair removal session at the lowest intensity level, raise the intensity level by 1 intensity level at the next session. Do the same for each subsequent session until you reach the highest intensity that feels comfortable.

- **DO NOT raise the intensity level if you experience abnormal discomfort during or after treatment with Philips Lumea Comfort (see chapter 'Risks of using Philips Lumea Comfort').**
- **DO NOT raise the intensity level during a hair removal session even if you experience no discomfort. Discomfort may also appear some time after the session.**

When you turn on Philips Lumea Comfort, its intensity level is automatically reset to level 1. At the start of each session, only one intensity light is on: the intensity light for level 1. This means that, at the start of each session, you have to increase the intensity to the highest level that still felt comfortable in the previous session.

Slide & Flash or Step & Flash

Slide & Flash

Keep the flash button pressed. The appliance releases one flash after another as long as the skin contact surface is in full contact with the skin. After each flash, slide the appliance to the next area. Continue to press the skin contact surface onto the skin. This Slide & Flash mode is especially convenient for treating larger areas like the legs.

Step & Flash

Let go of the flash button after each flash. For the next flash, place the appliance on the area right next to the area just treated. Make sure that the skin contact surface is fully pressed onto the skin. This Step & Flash mode guarantees a precise treatment of e.g. the knees or ankles.

- **Try to avoid overlapping flashes! Do not treat the same area of skin more than once during a hair removal session!**

- **If your skin blisters or burns, STOP USING THE DEVICE IMMEDIATELY!**

What to expect when treating with Philips Lumea Comfort

For many people, using Philips Lumea Comfort may be their first experience with a light-based hair removal device designed for home use. Philips Lumea Comfort is simple to use, and hair removal sessions proceed quickly. During a treatment session with Philips Lumea Comfort, it is normal to experience and feel:

- A **flash of light** - The bright light emitted by Philips Lumea Comfort does not harm the eyes when directed at non-facial areas. Special eye protection is not necessary when using Philips Lumea Comfort.
- A **fan noise** - The cooling fan inside the device sounds similar to a hairdryer. This is normal.
- A **popping sound** at each flash - When a flash is activated, it is normal to hear a soft popping sound and see a flash at the same time when a flash is being emitted.

- **Moderate pressure of the treatment area** - This helps to ensure the adjacent administration of light flashes. It is part of the unique safety features of Philips Lumea Comfort.
- **A warm and tingling sensation** - When a flash is administered, it is normal to feel a mild warm and tingling sensation. Remember that it is important to always use low settings for initial hair removal sessions. You may continue to feel some warmth up to an hour after the treatment.
- **Mild red or pink color** - During or just after your treatment session, it is not uncommon to see very mild pink coloration of the skin. This pink color is usually most noticeable around the hairs themselves. If your skin is completely red or shows blisters or burns, stop using the device immediately.

After treatment with Philips Lumea Comfort

- 1 When you have completed your treatment session, press and hold the intensity button to turn off the device. Be sure to remember the last intensity level setting you used, as the device is automatically reset to the lowest setting when you turn it off, and it does not remember the previously selected setting.
- 2 Unplug the adapter from the electrical outlet and the AC cord from the device socket.
- 3 After each treatment session, you are recommended to clean your Philips Lumea Comfort, especially its skin contact surface (see chapter 'Maintenance').
- 4 After cleaning, store your Philips Lumea Comfort in its original box and keep it away from water.

Skin care after hair removal session

- Do not expose the treated areas of the skin to the sun. To protect recently treated skin when exposed to sunlight, be sure to apply a sunscreen with SPF 30 or higher thoroughly for 2 weeks after each hair removal session.

Side effects and complications

Some users may experience pigmentation changes as a result of treatment with Philips Lumea Comfort. These effects, if they occur, are generally mild and temporary.

- **If you experience any complication, please contact your physician immediately!**
See chapter 'Risks of using Philips Lumea Comfort' for risks that may occur.

Cleaning and Maintenance of Philips Lumea Comfort

It is recommended that you clean the device and especially its skin contact surface after each hair removal session.

- 1** Turn off the device and unplug it from the electrical outlet.
- 2** Use a dry, clean cloth and a specially formulated cleaner for electronic equipment to gently wipe the outer shell and the skin contact surface of the device.

Never immerse the device or any of its parts in water.

Disposal

- Contact your local town or city officials for recycling information.

Troubleshooting

My Philips Lumea Comfort does not start.

- Make sure that the power cord is plugged into the device socket properly.
- Make sure that the adapter is plugged into the electrical outlet properly.

No flash is emitted when I press the flash button

- Make sure that the yellow READY light is blinking.
- Make sure that you have good skin contact and that the skin contact surface is evenly and firmly pressed onto the skin. For your safety, the flash button will only activate a flash if the skin contact surface is in full contact with the skin.
- Check if the skin color warning light blinks orange. If it does, the skin color sensor has stopped the emission of a flash. Your skin may be too dark to use Philips Lumea Comfort. Try using the device on a different part of your body.
- If the READY light is on constantly together with any number of other lights, there is a system error. Turn off the device and unplug it. Then plug in the device and try to turn it on again. If the problem persists, please visit www.philips.com/support or call 1-800-234-7884 for assistance.

- Check if the system warning light is on and all other lights are blinking. If so, unplug the device. The device has reached the end of its life. Please refer to chapter 'Disposal' for instructions on how to dispose of the device.
- **Do not attempt to open or repair your Philips Lumea Comfort. For assistance call 1-800-243-7884.**
- Do not attempt to open or repair Philips Lumea Comfort. Opening the device may expose you to dangerous electrical components and to intense light, either of which may cause bodily damage and/or eye injury.
- If you try to open Philips Lumea Comfort, you may also damage the device. Opening the device will void your warranty. Please visit www.philips.com/support or call 1-800-243-7884 for assistance if you have a broken or damaged device in need of repair.

Assistance

For assistance, visit our website: **www.philips.com/support** or call toll-free 1-800-243-7884. Online information is available 24 hours a day, 7 days a week.

Full Two-Year Warranty

Full Two-Year Warranty

Philips Electronics North America Corporation warrants each new Philips product, Model SC1982, SC1981, against defects in materials or workmanship for a period of two years from the date of purchase, and agree to repair or replace any defective product without charge.

IMPORTANT: This warranty does not cover damage resulting from accident, misuse or abuse, lack of reasonable care, the affixing of any attachment not provided with the product or loss of parts or subjecting the product to any but the specified voltage.* Use of unauthorized replacement parts will void this warranty.

NO RESPONSIBILITY IS ASSUMED FOR ANY SPECIAL, INCIDENTAL OR CONSEQUENTIAL DAMAGES.

In order to obtain warranty service, simply go to www.philips.com or call 1-800-243-7884 for assistance. It is suggested that for your protection you return shipments of product by insured mail, insurance prepaid. Damage occurring during shipment is not covered by this warranty.

NOTE: No other warranty, written or oral, is authorized by Philips Electronics North America Corporation. This warranty gives you specific legal rights, and you may also have other rights which vary from state to state or province to province. Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above exclusion and limitations may not apply to you.

*Read enclosed instructions carefully.

Manufactured for:

Philips Consumer Lifestyle, A Division of Philips Electronics North America Corporation,
P.O. Box 10313, Stamford, CT 06904

PHILIPS and Philips Shield are Registered Trademarks of Koninklijke Philips N.V.

© 2014 Philips Electronics North America Corporation. All Rights Reserved.

45-Day Money-Back Guarantee

If you are not fully satisfied with your Philips Lumea send the product back and we'll refund you the full purchase price.

The Lumea must be shipped prepaid by insured mail, insurance prepaid, have the original sales receipt, indicating purchase price and date of purchase, and the money-back guarantee return authorization form enclosed. We cannot be responsible for lost mail. The Lumea must be postmarked no later than 45 days after the date of purchase. Philips reserves the right to verify the purchase price of the product and limit refunds not to exceed suggested retail price. To obtain money-back guarantee return authorization form, call 1-800-243-7884 for assistance.

Delivery of refund check will occur 6–8 weeks AFTER receipt of returned product.

Specifications

Model	SC1982, SC1981	SC1986, SC1985
Spot size	0.35 × 1.18in [0.42 sq. in.] (0.9cm × 3cm [2.7cm ²])	0.9 cm × 3cm [2.7cm ²]
Technology used	Intense Pulsed Light	Intense Pulsed Light
Flash duration	500 - 800 µsec	500 - 800 µsec
Wavelength	475 - 1200 nm	475 - 1200 nm
Power output	max. 4 J/cm ²	max. 4 J/cm ²
Charging time/ Power source	Continuous operation	Continuous operation

Model	SC1982, SC1981	SC1986, SC1985
Time needed to treat one lower leg	11 minutes	6,5 minutes
Safety during operation	Skin color sensor ensures that a flash is only delivered to appropriate skin colors.	Skin color sensor ensures that a flash is only delivered to appropriate skin colors.
Classified according to standard	IEC 60601-2-57:2011	IEC 60601-2-57:2011
- Risk group	exempt	exempt
Package size	7.1 × 7.1 × 3.9in (180 × 180 × 100mm) HxWxD	8.7 × 8.8 × 3.5in (220 × 223 × 90mm) HxWxD
Product weight	7.1oz. (200g)	7.1oz. (200g)
Adapter	Model: KSA24A1200150HU	Model: KSAPV0361200300D5
- Rated	100-240V; 50-60Hz; 0.5A	100-240V; 50-60Hz; 0.5A


Operating conditions

Temperature	+50°F to +95°F (+10°C to +35°C)
Relative humidity	30% to 75% rH
Atmospheric pressure	700 to 1060hPa

Storage and transport conditions

Temperature	-40°F to +158°F (-40°C to +70°C)
Relative humidity	10 to 90% rH
Atmospheric pressure	500 to 1060hPa

Labeling

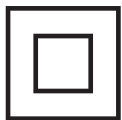- This symbol indicates that this device is double insulated (Class II).

- This symbol means: Contact your local town or city officials for recycling information.

- This symbol means: Conforms to EC Directives. CE stands for 'Conformité Européenne'.

- Degree of protection against ingress of water: ordinary.

- This symbol means: Read the user manual before you start using the device.

- This symbol means: Type BF (Body Floating) applied part (i.e. electrode) complying with IEC 60601-1.

- This is the CSA mark for the USA and Canada.

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© 2014 Philips Electronics North America Corporation. Reservados todos los derechos.