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EN Always perform a skin test before use. After the test, wait 24 hours. If you see skin reactions, use the highest intensity level that did not give skin reactions.

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
EN If you recently sunbathed or used artificial tanning equipment, wait at least 2 weeks before you use the appliance.

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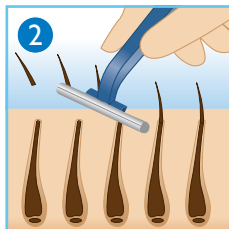
EN Shave the skin before use. When the hairs stop growing back, you do not have to shave anymore.

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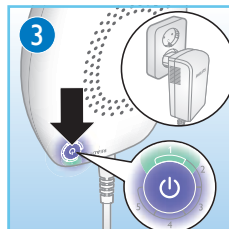
EN Plug in the appliance and switch it on. The light of intensity 1 goes on and the fan inside starts.

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SC1982, SC1981

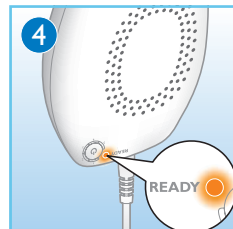
EN After 1 second, the READY/skin colour warning light lights up. The appliance is ready to flash.

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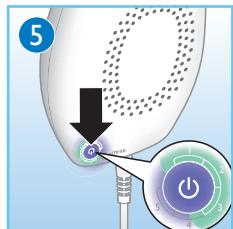
EN Press the on/off/intensity button to increase the intensity. If the treatment becomes painful, reduce the intensity.

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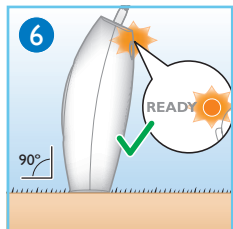
EN Place the skin contact surface on your skin. When there is full skin contact, the READY/skin colour warning light starts to blink slowly.

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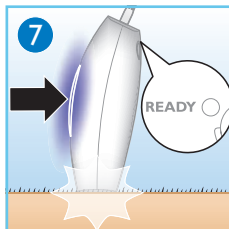
EN Press the flash button. If your skin colour is light enough, the READY/skin colour warning light goes off and the appliance releases a flash.

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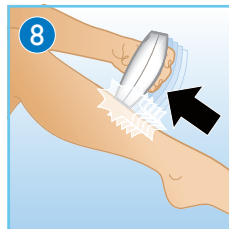
EN In the Slide & Flash mode, you keep the flash button pressed. The appliance continues to release flashes as you slide it across your skin.

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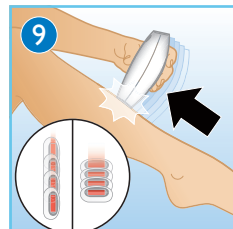EN In the Step & Flash mode, release the flash button after each flash. Then move the appliance to an area right next to the area you just treated and release another flash.

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EN - ⊗ means that this method is not suitable for you.

- 1 to 5 are the light intensity settings. The table states the suitable light intensity for each skin colour.

Note: If the treatment becomes painful at any point in time, reduce the light intensity setting until the treatment is comfortable again.

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