Always here to help you

Register your product and get support at
www.philips.com/welcome


SC1982
SC1981

User manual

PHILIPS

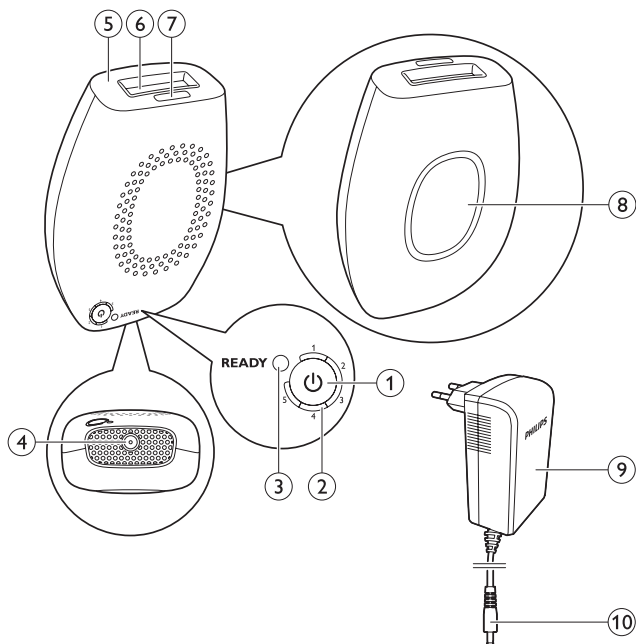ENGLISH 6

ESPAÑOL 35

Table of contents

Introduction	6
What is Philips Lumea Comfort?	7
Intended use of Philips Lumea Comfort	8
Safety with Philips Lumea Comfort	8
Contraindications	9
When to avoid using Philips Lumea Comfort	11
Precautions - How to use Philips Lumea Comfort safely	12
Reducing the risk of injury	14
Risks of using Philips Lumea Comfort	16
Treatment plan and effectiveness	19
Treatment plan	19
Expected treatment outcome	19
Expected treatment outcome for hair removal	22
How to use Philips Lumea Comfort	23
Philips Lumea Comfort setup	23
Treating with Philips Lumea Comfort for the first time	23
Intensity level	25
Covering skin areas	26
What to expect when treating with Philips Lumea Comfort	27
After treatment with Philips Lumea Comfort	28
Maintenance of Philips Lumea Comfort	29
Disposal	30
Troubleshooting	30
Customer support	31
Warranty	31
Specifications	32
Labeling	34


Introduction

Congratulations on your purchase and welcome to Philips! To benefit fully from the support that Philips offers, register your product at www.philips.com/lumea.

Before you use Philips Lumea Comfort for the first time, please read the whole user manual. Pay special attention to all warnings and safety precautions and follow them strictly during use.

What is Philips Lumea Comfort?

Philips Lumea Comfort is a light-based device for long-term hair removal, designed for home use.

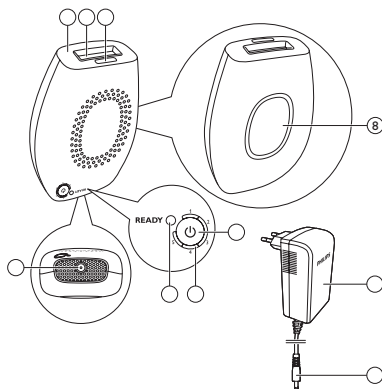Device description (Fig. 1)

- 1 Intensity button
- 2 Intensity lights
- 3 READY/skin color warning light
- 4 Device socket
- 5 Skin contact surface
- 6 Light exit window
- 7 Skin color sensor
- 8 Flash button
- 9 Adapter
- 10 Small plug of AC adapter cord

Box contents

When you open the box, you find it contains the following parts:

- 1 Philips Lumea Comfort device
- 2 AC cord adapter
- 3 This user manual

Intended use of Philips Lumea Comfort

Philips Lumea Comfort is intended for the removal of unwanted hair. Philips Lumea Comfort is intended to effect long-term or permanent results. Permanent results are defined as a long-term stable reduction in the number of hairs regrowing after a treatment regime.

Philips Lumea Comfort may be used to remove unwanted body hair. Ideal body areas for use of Philips Lumea Comfort include the underarms, bikini line, arms and legs.

Safety with Philips Lumea Comfort

With Philips Lumea Comfort safety comes first.

HPL technology offers superior safety with lower energy level

Home Pulsed Light technology is able to achieve long-term hair removal results at a fraction of the energy level used in other light-based hair removal equipment. The low energy used in Philips Lumea Comfort reduces its potential to cause harm or complications, and contributes to your overall safety.

Philips Lumea Comfort protects your eyes

The Philips Lumea Comfort has a **built-in safety feature** for eye protection. It has been designed so that a light pulse cannot be emitted when the treatment surface is facing open air. The pulse button only works when the treatment surface is in full contact with the skin.

Philips Lumea Comfort protects your skin

Philips Lumea Comfort comes with a skin color sensor that only enables use on lighter skin colors to ensure skin safety.

The 0.35x1.2in (0.9X3cm) light exit window through which pulses of light are delivered is recessed in the treatment surface. This enables Philips Lumea Comfort to protect your skin by avoiding direct contact between the light exit window and the skin

Contraindications**IMPORTANT SAFETY INFORMATION
READ BEFORE USE!**

Philips Lumea Comfort is not designed for everyone. Please read and consider the information in the following section before use. For further information and personalized advice, you may also visit www.philips.com/lumea.

If you are unsure Philips Lumea Comfort is safe for you to use, please consult your physician or dermatologist.

DO NOT USE Philips Lumea Comfort if you have a naturally dark skin color!

Philips Lumea Comfort removes unwanted hair by selectively addressing hair pigment. Varied quantities of pigment also exist in the surrounding tissue of skin. The quantity of pigment in a particular person's skin, which is manifested by their skin color, determines the degree of risk they are exposed to when they use this device.

Treating dark skin with Philips Lumea Comfort can result in adverse effects such as burns, blisters, and skin color changes (hyper- or hypo-pigmentation).

- **DO NOT USE Philips Lumea Comfort on naturally dark skin!**

A unique skin color sensor is incorporated in Philips Lumea Comfort. It measures the color of the skin to be treated at the beginning of each session and occasionally during the session. This skin color sensor ensures that a pulse will be emitted only on suitable skin colors.

DO NOT USE Philips Lumea Comfort on tanned skin or after recent sun exposure!

Tanned skin, particularly after sun exposure, contains large quantities of the pigment melanin. This applies to all skin types and colors, including those that do not seem to tan quickly. The presence of large quantities of melanin exposes the skin to higher risk when using Philips Lumea Comfort.

Using Philips Lumea Comfort on skin that has been recently exposed to sunlight can result in adverse effects such as burns, blisters, and skin color changes (hyper- or hypo-pigmentation).

- **DO NOT USE Philips Lumea Comfort on tanned skin or after sun exposure! Such use can cause burns or skin injury.**

- **Avoid exposure to the sun for 4 weeks before your treatment with Philips Lumea Comfort!**

Never use Philips Lumea Comfort around or near the eyes - Philips Lumea Comfort is not recommended for facial treatment.

Philips Lumea Comfort is not effective on naturally white, gray, blond and red body hair. If your body hair is of these colors, the device does not work on you.

When to avoid using Philips Lumea Comfort

Certain conditions may prevent the use of Philips Lumea Comfort temporarily. **DO NOT USE THE DEVICE** if any of the following currently applies to you:

- If you are pregnant or nursing.
- If you were exposed to strong sunlight or an artificial tanning machine during the past 4 weeks.
- If you have a tattoo or permanent makeup on the area to be treated.
- If you have dark spots, such as large freckles, birthmarks, moles or warts in the area to be treated.
- If you have eczema, psoriasis, lesions, open wounds or active infections in the area to be treated.
- If you have a history of keloidal scar formation, a known sensitivity to light (photosensitivity) or are taking medication that makes the skin more sensitive to light, including non-steroidal anti-inflammatory agents, (e.g., aspirins, ibuprofens, acetaminophen), tetracyclines, phenothiazines, thiazide, diuretics, sulfonamides, DTIC, fluorouracil, vinblastine, griseofulvin, Alpha-Hydroxi Acids (AHAs), Beta-Hydroxi Acids (BHAs), Retin-A[®], Accutane[®] and/or topical retinoids.
- If you have abnormal skin conditions due to diabetes, for example. or another systemic or metabolic disease.
- If you are currently or have recently been treated with Alpha-Hydroxi Acids (AHAs), Beta-Hydroxi Acids (BHAs), Retin-A[®], Accutane[®] (isotretinoin), topical retinoids, azelaic acid.
- If you have been treated with Accutane[®] (isotretinoin) within the past 6 months.
- If you have been on a steroid regimen within the past 3 months.

- If you have a history of herpes outbreaks in the treatment area, unless you have consulted your physician and received preventive treatment before you start using the device.
- If you suffer from epilepsy.
- If you have an active implant, such as a pacemaker, incontinence device, insulin pump, etc.
- If you have a disease related to photosensitivity, such as porphyria, polymorphic light eruption, solar urticaria, lupus, etc.
- If you have a history of skin cancer or potential skin malignancies.
- If you have received radiation therapy or chemotherapy treatments within the past 3 months.
- If you have any other condition which in your physician's opinion would make it unsafe for you to be treated

If you are unsure Philips Lumea Comfort is safe for you to use, please consult your physician or dermatologist.

Precautions - How to use Philips Lumea Comfort safely

Choose the intensity level carefully

Intensity level refers to the intensity of the pulse emitted onto your skin during use, from the lowest level to the highest level. The intensity lights around the intensity button indicate the intensity level at which the device is set. As the intensity level increases, so do the results as well as the risk of side effects (see chapter 'Risks of using Philips Lumea Comfort').

Always start using Philips Lumea Comfort at the lowest intensity level!

Only if you experience little or no discomfort during and after use of Philips Lumea Comfort at the lowest intensity level, raise the intensity by one level the next time you use the device and so on for each subsequent session.

For detailed instructions on setting the intensity level, see section 'Treating with Philips Lumea Comfort' in chapter 'How to use Philips Lumea Comfort'.

Avoid adverse effects!

- Do not treat the same area of skin more than once during a session!
- Avoid overlapping pulses!
- If your skin blisters or burns, **STOP USING THE DEVICE IMMEDIATELY!**

Avoid complications after use!

- Do not expose treated areas of skin to the sun. Sun exposure includes constant unprotected exposure to direct sunlight for over 15 minutes, or constant unprotected exposure to diffuse sunlight for over 1 hour. To protect recently treated skin when exposed to sunlight, be sure to apply a sunscreen with SPF 30 or higher thoroughly for 2 weeks after each hair removal session.
- Always shave the area to be treated and make sure that the skin is clean and dry before you use the device.
- Cover birthmarks and tattoos before treatment.
- Cover dark brown or black spots, such as large freckles, birthmarks, moles or warts before treatment.

- Never look directly at light coming out of the light exit window of the device.
- Do not use the device on nipples and genitals (male or female).
- Do not use the device on any body site where you might later want hair.
- Do not use the device for any other purpose than hair removal.
- Never emit a pulse into open space. Always make sure the skin contact surface is in full contact with the skin during treatment.
- Remove the skin contact surface from the skin if either the skin or the surface is too hot.
- Never use flammable liquids such as alcohol (including perfumes, sanitizers, or other substances containing alcohol) or acetone to clean the skin before you use the device.
- Use of the device may cause temporary pigmentation changes (see chapter 'Risks of using Philips Lumea Comfort').
- Keep this device out of the reach of children. Do not use the device on children or allow children to use it.

Reducing the risk of injury

As with any electrical device, certain precautions must be taken to ensure your safety when using Philips Lumea Comfort.

- **No modification of the device is allowed!**
- Keep Philips Lumea Comfort away from water. This is an electrical device. As such, it should always be kept away from water. Do not place or store Philips Lumea Comfort where it can

fall or be pushed into a tub, sink or any other vessel containing water. Do not place in, or drop into water or any other liquid. This may cause electrocution!

- Do not use Philips Lumea Comfort while bathing.
- Do not use Philips Lumea Comfort if it becomes damp or wet.
- Do not reach for Philips Lumea Comfort if it has fallen into water. Unplug Philips Lumea Comfort immediately if it has fallen into water.

- **Keep Philips Lumea Comfort away from water!**

- Never open Philips Lumea Comfort! Do not attempt to open or repair Philips Lumea Comfort. Opening the device may expose you to dangerous electrical components and to intense light, either of which may cause bodily damage and/or eye injury.

- **Do not attempt to open or repair your Philips Lumea Comfort. Always have the device repaired by an authorized Philips repair center.**

- If you try to open Philips Lumea Comfort, you may also damage the device. Opening the device will void your warranty. Please visit www.philips.com/support or call 1-800-243-7884 for assistance if you have a broken or damaged device in need of repair.
- Use Philips Lumea Comfort only for its intended purpose and as described in this user manual.
- Do not use the device if you see or smell smoke during use.

- Do not use the device if it does not work properly or if it appears damaged.
- Do not use the device if the fan vent in the skin contact surface is cracked, detached or missing.
- Do not use the device if the skin color sensor is cracked or broken.
- Do not use the device if the outer shell is cracked or coming apart.
- Always unplug the device from the electrical outlet immediately after use.
- Always unplug the device before cleaning.
- Do not use the device with any attachments or accessories not recommended by Philips.
- Philips Lumea Comfort must be operated only with its designated adapter. For details on the right adapter, see chapter 'Specifications'.
- Never leave the device unattended when plugged into an electrical outlet.
- Do not operate the device if it has a damaged cord or adapter. Keep the cord away from heated surfaces.
- This device is not suitable for use in the presence of flammable anaesthetic mixtures with air, oxygen or nitrous oxide.

Risks of using Philips Lumea Comfort

Side effects and complications associated with the use of Philips Lumea Comfort are uncommon. However, every cosmetic procedure, including those designed for home use, involves some degree of risk. Therefore, it is important that you understand and accept the risks and complications that can occur with pulsed light hair removal systems designed for home use.

Adverse event	Likelihood of adverse event	Adverse effect	Likelihood of the adverse effect
Stacking or overlapping of multiple pulses	Minor	Minor skin discomfort	minor
		Skin redness	minor
		Increased skin sensitivity	minor
		Skin wounds and burns	rare
		Scarring	rare
		Pigment changes	rare
		Excessive redness and swelling	rare
		Infection	negligible
		Bruising	negligible

Minor skin discomfort

Although home pulsed light hair removal is generally very well-tolerated, most users feel some mild discomfort during use. This minor discomfort is usually described as being a mild stinging sensation on the treated skin areas. The stinging sensation usually lasts during release of the pulse or for a few minutes after release of the pulse. Anything beyond this minor discomfort is abnormal and means that either you should not continue to use the device because you are unable to tolerate its hair removal method, or that the intensity level is set too high.

Skin redness

Your skin may become red right after using Philips Lumea Comfort or within 24 hours of using Philips Lumea Comfort. Redness generally clears up within 24 hours. Consult your physician if redness does not go away within 2 to 3 days.

Increased skin sensitivity

The skin of the treated area is more sensitive so you may encounter dryness or flaking of the skin.

Skin wounds and burns

Very rarely, burns or wounds to the skin can occur following use of the device. The burn or wound can require a few weeks to heal and, extremely rarely, may leave a noticeable scar.

Scarring

Although very rare, scarring may occur. Usually when scarring occurs, it is in the form of a flat and white lesion on the skin (hypotrophic). However, it can be large and red (hypertrophic) or large and extend beyond the margins of the injury itself (keloid). Subsequent aesthetic treatments may be required to improve the appearance of the scar.

Pigment changes

Philips Lumea Comfort targets the hair shaft, in particular the pigmented cells in the hair follicle and the hair follicle itself. Nevertheless, there is risk of temporary hyperpigmentation (increased pigment or brown discoloration) or hypopigmentation (whitening) to the surrounding skin. This risk of changes in skin pigmentation is higher for people with darker skin color.

Usually discoloration or changes to skin pigment are temporary and hyperpigmentation or hypopigmentation rarely occur.

Excessive redness and swelling

In rare cases, treated skin may become very red and swollen. This is more common in sensitive areas of the body. The redness and swelling should subside within 2 to 7 days. Application of ice may help to reduce the swelling. Gentle cleansing of the skin is allowed, but exposure to the sun should be avoided.

Infection

Infection of the skin is exceedingly rare. However, it is a possible risk following a skin burn or wound caused by treatment with Philips Lumea Comfort.

Bruising

Very rarely, use of Philips Lumea Comfort may cause blue-purple bruising that may last 5 to 10 days. As the bruising fades, there may be a rust-brown discoloration of the skin (hyperpigmentation).

Treatment plan and effectiveness

A typical full hair growth cycle may take 18-24 months. During this time, multiple treatments with Philips Lumea Comfort may be required in order to achieve permanent hair removal.

The effectiveness of hair removal varies from person to person according to body area, hair color and how the Philips Lumea Comfort is used.

Treatment plan

The typical treatment plan during a full hair growth cycle:

- Treatments 1 - 4 : Plan 2 weeks apart
- Treatments 5 - 7: Plan 4 weeks apart
- Treatments 8 and after: Treat as needed until decided results are achieved.

Expected treatment outcome

The treatments take place at a light fluence of 3 - 4.5 J/cm².

Expected treatment outcome for light brown to brown hair

Body area	Average no. of treatments, approx 4 weeks apart	Effectiveness (% of hair removal/regrowth)
Lower or upper legs	10 - 12	60%
Arms	10 - 12	50%
Bikini line	10 - 12	60%
Underarms	10 - 12	55%

Expected treatment outcome for brown to black hair

Body area	Average no. of treatments, approx 4 weeks apart	Effectiveness (% of hair removal/regrowth)
Lower or upper legs	8 - 10	70%
Arms	8 - 10	60%
Bikini line	8 - 10	70%
Underarms	8 - 10	65%

Notes

When looking at the tables above, please bear in mind the following:

- Individual response depends on hair type as well as biological factors that may affect hair growth patterns. Some users may respond quicker or slower than the average number of treatments.
- You cannot expect permanent hair removal in one or two treatment sessions. The duration of the resting period for hair follicles also depends on the body area.

Expected treatment outcome for hair removal

Fitzpatrick skin type	Hair color	Light fluence [Jcm ²]	Anatomic area of the body	Average number of treatments	Effectiveness (% of hair removal/regrowth)
I-VI	Light brown to brown	3-4.5	Lower or upper legs	Session 10-12, approx. 4 weeks apart	60% hair reduction
			Arms		50% hair reduction
			Bikini line		60% hair reduction
			Under-arms		55% hair reduction
I-VI	Brown to black	3-4.5	Lower or upper legs	Session 8-10, approx. 4 weeks apart	70% hair reduction
			Arms		60% hair reduction
			Bikini line		70% hair reduction
			Under-arms		65% hair reduction

Notes

When looking at the tables above, please bear in mind the following:

- Individual response depends on hair type as well as biological factors that may affect hair growth patterns. Some users may respond quicker or slower than the average number of treatments.
- You cannot expect permanent hair removal in one or two treatment sessions. The duration of the resting period for hair follicles also depends on the body area.

How to use Philips Lumea Comfort

Philips Lumea Comfort setup

- 1** Remove the Philips Lumea Comfort device and other parts from the box.
- 2** Plug the AC cord into the device socket of the Philips Lumea Comfort device.
- 3** Plug the adapter into an electrical outlet.
- 4** All 6 lights blink on and off, one after the other, in a circle.

Your Philips Lumea Comfort is now ready to start.

Treating with Philips Lumea Comfort for the first time

Make sure that your skin is shaved, clean, dry and free of any powders, creams, antiperspirants or deodorants.

- 1** Press the intensity button briefly. The green intensity level 1 light turns on.

- 2** About 1 second later, the yellow READY light turns on. The device is now ready to deliver a first pulse at the lowest intensity level.
- 3** If you wish to change the intensity level, press the intensity button quickly and repeatedly to increase the intensity level from 1 to 5 and from 5 directly back to 1, cyclically, until you reach the desired level. The green level lights indicate the level selected (see section 'Intensity level' below).
- 4** Place the skin contact surface on the skin. Make sure that the skin is spread evenly and smoothly and the skin contact surface is in full contact with the skin. Make sure that the yellow READY light starts blinking slowly.
- 5** Press the pulse button. The device first determines the color of your skin. If the color of the skin is light enough for safe treatment, the device emits a pulse of light and produces a soft popping sound at the same time. This popping sound is normal for this device. You may also feel a mild warm and tingling sensation.

If you see no light pulse and the yellow READY light blinks rapidly, this light is now the skin color warning light that warns you that your skin color is too dark for safe treatment. Try using the device on a different body part. For assistance, please visit www.philips.com/support or call 1-800-243-7884.

- 6** After emitting a pulse, Philips Lumea Comfort recharges immediately for the next pulse. Recharging takes a few seconds, depending on the intensity level set. When recharging is finished, the READY light turns on again.
- 7** Move the skin contact surface to another area of skin. When there is full skin contact, the READY light starts blinking slowly again. Be careful not to treat the same area of skin more than once per treatment session! Treating the same area of skin more than once per session increases the likelihood of adverse effects.
- 8** Repeat the process, starting at step 5.

Intensity level

The intensity level setting determines the intensity of the light pulse delivered to your skin. There are 5 levels, from the lowest (1) to the highest (5). Corresponding intensity lights show the intensity selected. As the intensity level increases, so do the hair removal results and the risk of possible side effects and complications.

Always start using Philips Lumea Comfort at the lowest intensity level.

If you experience little or no discomfort during and after your first hair removal session at the lowest intensity level, raise the intensity level by 1 intensity light at the next session. Do the same for each subsequent session until you reach the highest intensity that feels comfortable.

- **DO NOT** raise the intensity level if you experience abnormal discomfort during or after treatment with Philips Lumea Comfort (see chapter 'Risks of using Philips Lumea Comfort').
- **DO NOT** raise the intensity level during a hair removal session even if you experience no discomfort. Discomfort may also appear some time after the session.

When you turn on Philips Lumea Comfort, its intensity level is automatically reset to level 1. At the start of each session, only one intensity light is on: the intensity light for level 1. This means that, at the start of each session, you have to increase the intensity to the highest level that still felt comfortable in the previous session.

Covering skin areas

You have to administer the pulses in rows. Start at one end of each row and progress sequentially towards the other end. This technique allows better control over skin coverage and helps avoid treating the same area more than once or overlapping areas treated. When you place Philips Lumea Comfort on the skin, the device creates temporary pressure marks on the area treated. These marks can be used to position the device for administering the subsequent pulse.

- **Try to avoid overlapping pulses! Do not treat the same area of skin more than once during a hair removal session!**

- **If your skin blisters or burns, STOP USING THE DEVICE IMMEDIATELY!**

What to expect when treating with Philips Lumea Comfort

For many people, using Philips Lumea Comfort may be their first experience with a light-based hair removal device designed for home use. Philips Lumea Comfort is simple to use, and hair removal sessions proceed quickly. During a treatment session with Philips Lumea Comfort, it is **normal** to experience and feel:

- A **flash of light** - The bright light emitted by Philips Lumea Comfort does not harm the eyes when directed at non-facial areas. Special eye protection is not necessary when using Philips Lumea Comfort.
- A **fan noise** - The cooling fan inside the device sounds similar to a hairdryer. This is normal.
- A **popping sound** at each pulse - When a pulse is activated, it is normal to hear a soft popping sound and see a flash at the same time when a pulse is being emitted.
- **Moderate pressure of the treatment area** - This helps to ensure the adjacent administration of light pulses. It is part of the unique safety features of Philips Lumea Comfort.

- A **warm and tingling sensation** - When a pulse is administered, it is normal to feel a mild warm and tingling sensation. Remember that it is important to always use low settings for initial hair removal sessions. You may continue to feel some warmth up to an hour after the treatment.
- **Mild red or pink color** - During or just after your treatment session, it is not uncommon to see very mild pink coloration of the skin. This pink color is usually most noticeable around the hairs themselves. If your skin is completely red or shows blisters or burns, stop using the device immediately.

After treatment with Philips Lumea Comfort

- 1** When you have completed your treatment session, press and hold the intensity button to turn off the device. Be sure to remember the last intensity level setting you used, as the device is automatically reset to the lowest setting and does not remember the previously selected setting.
- 2** Unplug the adapter from the electrical outlet and the AC cord from the device socket.
- 3** After each treatment session, you are recommended to clean your Philips Lumea Comfort, especially its skin contact surface (see chapter 'Maintenance').
- 4** After cleaning, store your Philips Lumea Comfort in its original box and keep it away from water.

Skin care after hair removal session

- Do not expose the treated areas of the skin to the sun. To protect recently treated skin when exposed to sunlight, be sure to apply a sunscreen with SPF 30 or higher thoroughly for 2 weeks after each hair removal session.

Side effects and complications

Some users may experience pigmentation changes as a result of treatment with Philips Lumea Comfort. These effects, if they occur, are generally mild and temporary.

- **If you experience any complication, please contact your physician immediately!** See chapter 'Risks of using Philips Lumea Comfort' for risks that may occur.

Maintenance of Philips Lumea Comfort

It is recommended that you clean the device and especially its skin contact surface after each hair removal session.

- 1** Turn off the device and unplug it.
- 2** Use a dry, clean cloth and a specially formulated cleaner for electronic equipment to gently wipe the outer shell and the skin contact surface of the device.

Never immerse the device or any of its parts in water.

Disposal

- Do not throw away the device with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment.

Troubleshooting

My Philips Lumea Comfort does not start.

- Make sure that the AC cord is plugged into the device socket properly.
- Make sure that the adapter is plugged into the electrical outlet properly.

No light pulse is emitted when I press the pulse button

- Make sure that the yellow READY light is blinking.
- Make sure that you have good skin contact and that the skin contact surface is evenly and firmly pressed onto the skin. For your safety, the pulse button will only activate a pulse if the skin contact surface is in full contact with the skin.
- Check if the skin color warning light blinks orange. If it does, the skin color sensor has stopped the emission of a pulse. Your skin may be too dark to use Philips Lumea Comfort. Try using the device on a different part of your body.
- If the READY light is on constantly together with any number of other lights, there is a system error. Turn off the device and unplug it. Then plug in the device and try to turn it on again. If the problem persists, please visit www.philips.com/support or call 1-800-243-7884 for assistance.

- Check if the system warning light is on and all other lights are blinking. If so, unplug the device. The device has reached the end of its life. Please refer to chapter 'Disposal' for instructions on how to dispose of the device.
- **Do not attempt to open or repair your Philips Lumea Comfort. Always have the device repaired by an authorized Philips repair center.**
- Do not attempt to open or repair Philips Lumea Comfort. Opening the device may expose you to dangerous electrical components and to intense light, either of which may cause bodily damage and/or eye injury.
- If you try to open Philips Lumea Comfort, you may also damage the device. Opening the device will void your warranty. Please visit www.philips.com/support or call 1-800-243-7884 for assistance if you have a broken or damaged device in need of repair.

Customer support

If you need information, support, or if you have a problem, please visit **www.philips.com/support** or call 1-800-243-7884 for assistance.

Warranty

Philips warrants that this product shall be free from defects in material, workmanship and assembly, under normal use, in accordance with the specifications and warnings, for a period of one year from the date of purchase. This warranty extends

only to the original purchaser of the product, and is not transferable. To exercise your rights under this warranty, you must provide proof of purchase in the form of an original sales receipt that shows the product name and the date of purchase. For customer support or to obtain warranty service, please visit our website www.philips.com/support. THERE ARE NO OTHER EXPRESS OR IMPLIED WARRANTIES. Philips' liability is limited to repair or, at its sole option, replacement of the product. Incidental, special and consequential damages are disclaimed where permitted by law. This warranty gives you specific legal rights. You may also have other rights that vary from state to state. For support, please visit www.philips.com/support

Manufactured for:
 Philips Consumer Lifestyle
 A Division of Philips Electronics North America Corp.
 P.O. Box 10313, Stamford, CT 06904

Specifications

Model	SC1981
Spot size	0.35 × 0.46in [0.42 sq. in.] (0.9cm × 3cm [2.7cm ²])
Technology used	New Home Pulsed Light
Pulse duration	500 - 800 μsec
Wavelength	475 - 1200 nm
Charging time/Power source	Continuous operation
Time needed to treat lower legs	30 minutes


Model SC1981	
Safety during operation	Skin color sensor ensures that a pulse is only delivered to appropriate skin colors.
Classified according to standard	IEC 60601-2-57:2011
- Risk group	exempt
Ingress of water	IP30
Package size	9 × 9 × 5.2in (180 × 180 × 100mm) HxWxD
Product weight	7.1oz. (200g)
Adapter	Model: KSA24A1200150HU
- Rated	100-240V; 50-60Hz; 0.5A

Operating conditions	
Temperature	+50°F to +95°F (+10°C to +35°C)
Relative humidity	30% to 75% rH
Atmospheric pressure	700 to 1060hPa

Storage and transport conditions	
Temperature	-104°F to +158°F (-40°C to +70°C)
Relative humidity	10 to 90%
Atmospheric pressure	700 to 1060hPa

Level	Energy [J]	Floucnce [J/cm ²]	Repetition [sec]	Tolerance
1	3.5	1.3	1.1	+ / - 20
2	5.0	1.9	1.5	+ / - 20
3	6.5	2.4	2.3	+ / - 20
4	8.0	3.0	3.0	+ / - 20
5	9.5	3.5	3.5	+ / - 20

Labeling

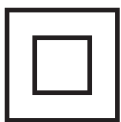- This symbol indicates that this device is double insulated (Class II).

- This symbol means: Do not throw away with the normal household waste. For further instructions, see chapter 'Disposal'.

- This symbol means: Conforms to EC Directives. CE stands for 'Conformité Européenne'.

- This symbol means: Type BF (Body Floating) applied part (i.e. electrode) complying with IEC 60601-1.

- This symbol means: Read the user manual before you start using the device.

- This is the CSA mark for the USA and Canada.

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