PHILIPS


Lumea

BRI861, BRI860

1

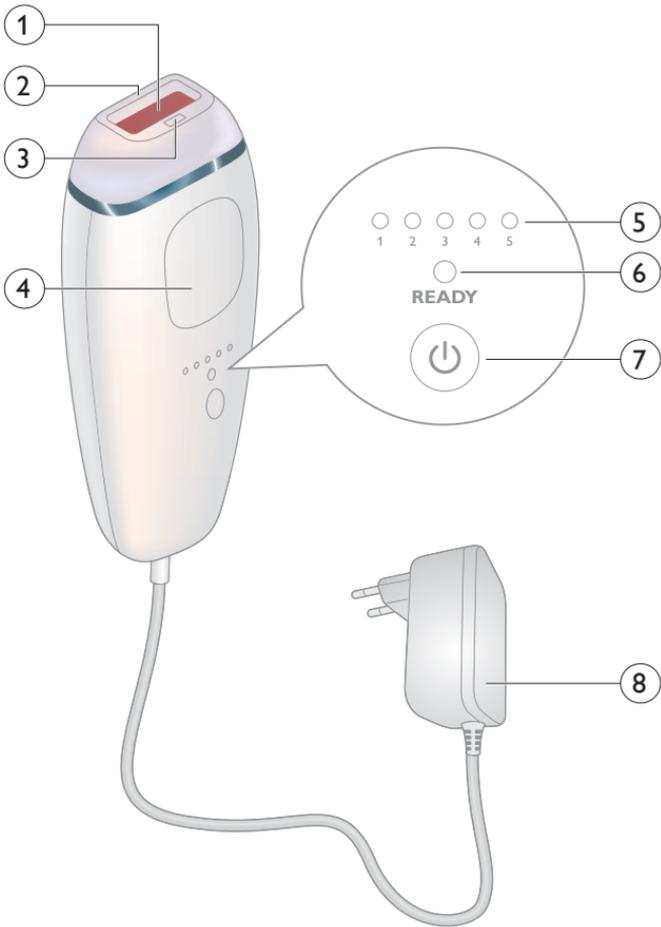English 6
Español 40
Español Latino 78
Português do Brasil 117
Türkçe 155

6 English

Contents

Welcome	7
General description (Fig. 1)	8
Who should not use Lumea? Contraindications	8
Important	15
Danger	15
Warning	16
To prevent damage	19
Caution	20
Electromagnetic fields (EMF)	22
How Intense Pulsed Light- technology works	23
Suitable body hair colors	23
What to expect	23
Immediately after the 1st treatment	23
After 2-3 treatments	24
After 4-5 treatments	24
Tanning advice	24
Tanning with natural or artificial sunlight	24
Tanning with creams	25
Before you use your Lumea	25
Pretreating your skin	25
Skin test	26
Using your Philips Lumea	27
Skin tone sensor	27
Selecting the right light intensity	27
Handling the device	29
Two treatment modes: Stamp & Flash and Slide & Flash	31
Recommended treatment schedule	31
Initial phase	31
Touch-up phase	31
Treatment time per area	32

After use	33
Common skin reactions	33
Rare side effects	33
Aftercare	34
Cleaning & storage	35
Guarantee and support	35
Recycling	35
Technical specifications	36
Troubleshooting	37

Welcome

Welcome to the beauty world of Lumea! You are only a few weeks away from silky-smooth skin. Philips Lumea uses Intense Pulsed Light (IPL) technology, known as one of the most effective methods to continuously prevent hair regrowth. In close collaboration with skin experts we adapted this light-based technology, originally used in professional beauty salons, for easy and effective use in the safety of your home. Philips Lumea is gentle and offers convenient and effective treatment at a light intensity that you find comfortable. Unwanted hairs are finally a thing of the past. Enjoy the feeling of being hair-free and look and feel amazing every day.

To fully benefit from the support that Philips offers, register your product at www.philips.com/welcome. For further information, please go to www.philips.com/lumea to find our experts' advice, tutorial videos and FAQs and make the most of your Lumea.

Note: Keep these instructions with your product at all times

General description (Fig. 1)

- 1 Light exit window with integrated UV filter
- 2 Integrated safety system
- 3 Skin tone sensor
- 4 Flash button
- 5 Intensity lights (1-5)
- 6 'Ready to flash' light
- 7 On/off button
- 8 Adapter

Who should not use Lumea? Contraindications

- Never use the device if you have skin type V and VI (You rarely to never sunburn, very dark tanning). In this case you run a high risk of developing skin reactions, such as hyperpigmentation and hypopigmentation, strong redness or burns.

- Never use the device if you are pregnant or breast feeding as the device was not tested on pregnant or breast feeding women.
 - Never use the device if you have any active implants such as pacemaker, neurostimulator, insulin pump etc.
 - Never use the device if you suffer from any of the diseases listed below:
 - If you have a skin disease such as active skin cancer, you have a history of skin cancer or any other localized cancer in the areas to be treated.
 - If you have a history of collagen disorder, including a history of keloid scar formation or a history of poor wound healing.

- If you have a history of vascular disorder, such as the presence of varicose veins or vascular ectasia in the areas to be treated.
- If your skin is sensitive to light and easily develops a rash or an allergic reaction.
- If you have infections, eczema, burns, inflammation of hair follicles, open lacerations, abrasions, herpes simplex, wounds or lesions and haematomas in the areas to be treated.
- If you have had surgery in the areas to be treated in the last three weeks.
- If you have received radiation therapy or chemotherapy within the past 3 months.
- Never use the device if you take any of the medications listed below:

- If your skin is currently being treated with or has recently been treated in the past week with Alpha-Hydroxy Acids (AHAs), Beta-Hydroxy Acids (BHAs), topical isotretinoin and azelaic acid.
- If you have taken any form of isotretinoin Accutane or Roaccutane in the last six months. This treatment can make skin more susceptible to tears, wounds and irritations.
- If you are on painkillers which reduce the skin's sensitivity to heat.

- If you are taking photosensitising agents or medications, check the package insert of your medicine and never use the device if it is stated that it can cause photo-allergic reactions, photo-toxic reactions or if you have to avoid sun when taking this medicine.
- If you take anticoagulation medications, including heavy use of aspirin, in a manner which does not allow for a minimum 1-week washout period prior to each treatment.
- If you take immunosuppressive medications.

- If you have had surgery in the areas to be treated in the last three weeks.
- Never use the device on the following areas:
 - Around the eyes and near the eyebrows.
 - On lips, nipples, areolas, labia minora, vagina, anus and the inside of the nostrils and ears.
 - Men must not use it on the face and neck including all beard growing areas and whole genital area.
 - Over or near anything artificial like silicone implants, subcutaneous injection ports (insulin dispenser) or piercings.

- On moles, freckles, large veins, darker pigmented areas, scars, skin anomalies without consulting your doctor. This can result in a burn and a change in skin color, which makes it potentially harder to identify skin-related diseases.
- On warts, tattoos or permanent make-up. This can result in a burn and a change in skin color.
- On areas, where you use long-lasting deodorants. This can result in skin reactions.
- Never use the device if your skin is irritated (red or cut), sunburned or exposed to strong sunlight or artificial tanning during the past 2 weeks.

Note: This list is not exhaustive. If you are not sure whether you can

use the device, we advise you to consult a doctor.

Important Danger

- Keep the device and the adapter dry.
- If the device is broken, do not touch any inner part to avoid electric shock.
- Never insert any objects into the device.
- Water and electricity are a dangerous combination. Do not use this appliance in wet surroundings (e.g. near a filled bath, a running shower or a filled swimming pool).

Warning

- This device is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, unless they have been given supervision or instruction concerning use of the device by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the device.
- The device is not intended for children under the age of 15 years. Teenagers aged between 15 and 18 years can use the device with the consent and/or assistance of their parents or the persons who have parental authority over them. Adults of 18 years and older can use the device freely.

- Always check the device before you use it. Do not use the device or adapter if it is damaged. Always replace a damaged part with one of the original type.
- The supply cord cannot be replaced. If the cord is damaged the appliance should be scrapped.
- The adapter contains a transformer. Do not cut off the adapter to replace it with another plug, as this causes a hazardous situation.
- Do not use the device if the UV filter of the light exit window is broken.
- The adapter cord is always attached to device. Don't try to separate the adapter from the device. This might damage the cord.

- Do not use any pencil or pen to mark the areas to be treated. This may cause burns on your skin.
- If you have a dark colored skin, be careful with treating a darker area immediately after a lighter area. The skin tone sensor may not block immediately the treatment on the darker body area.
- Hair removal by lasers or intense pulse light sources can cause increased hair growth in some individuals. Based upon currently available data, the highest risk groups for this response are females of Mediterranean, Middle Eastern and South Asian heritage treated on the face and neck.
- Light exit window can become very hot after usage. Do not touch the inner part of the light exit window.

Note: Perform a skin test if you notice a skin tone change from last treatment and wait 30 minutes before treatment.

To prevent damage

- Make sure that nothing obstructs the air flow through the ventilation slots.
- Never subject the appliance to heavy shocks and do not shake or drop it.
- If you take the appliance from a very cold environment to a very warm environment or vice versa, wait approximately 3 hours before you use it.
- Store the appliance in a dust free and dry place
- Do not expose the appliance to temperatures lower than 5°C or higher than 35°C during use.
- To prevent damage, do not expose the appliance to direct sunlight or UV light for several hours.

Caution

- This device is only intended for removing unwanted body hair from areas below the cheekbones. Do not use it for any other purpose.
- This device is not washable. Never immerse the device in water and do not rinse it under the tap.
- For hygienic reasons, the device should only be used by one person.
- Use the device only at settings suitable for your skin type. Use at higher settings than those recommended can increase the risk of skin reactions and side effects.
- Only use the device for its intended purpose as shown in the user manual.

- Never use compressed air, scouring pads, abrasive cleaning agents or aggressive liquids such as petrol or acetone to clean the device.
- Always return the device to a service centre authorised by Philips for examination or repair. Repair by unqualified people could cause an extremely hazardous situation for the user.
- Never leave the appliance unattended when it is switched on. Always switch off the appliance after use.
- Do not use the appliance if any of the conditions mentioned in chapter 'Contraindications, precautions and warnings', section 'Contraindications' apply to you.

- Tanning with natural or artificial sunlight might influence the sensitivity and color of your skin. Perform a skin test to determine the appropriate light intensity setting.
- Before you use Lumea, you should clean your skin and make sure it is hair-free, entirely dry and free from oily substances.
- Do not treat the same area of skin more than once during a session. This does not improve effectiveness of the treatment, but increases the risk of skin reactions.
- Lumea should never be painful. If you experience discomfort, reduce the light intensity setting.


Electromagnetic fields (EMF)

This Philips device complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

How Intense Pulsed Light- technology works

With IPL technology, gentle pulses of light are applied to the skin and absorbed by the hair root. The darker the hair, the better the pulses of light are absorbed.

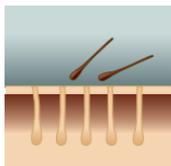The pulses of light stimulate the hair follicle to go into a resting phase. As a consequence, the hair sheds naturally and hair regrowth is prevented, leaving your skin continuously silky-smooth.

The cycle of hair growth consists of different phases. IPL technology is only effective when the hair is in its growing phase. Not all hairs are in the growing phase at the same time. This is why we recommend an 8-week treatment schedule followed by touch-ups to make sure all hairs are effectively treated in the growing phase.

Note: Treatment with Lumea is not effective if you have light blond, grey, red or white hairs as light hairs do not absorb enough light. Below you can see the hair colors for which Lumea is suitable and effective.

Suitable body hair colors

What to expect

Immediately after the 1st treatment

After the first treatment, it can take 1 to 2 weeks for the hairs to fall out. In the first weeks following the

initial treatments, you still see some hairs growing. These are likely to be hairs that were not in their growing phase during the first treatments.

After 2-3 treatments


After 2-3 treatments, you should see a noticeable reduction in hair growth. However, to effectively treat all hairs, it is important to keep on treating according to the recommended treatment schedule.

After 4-5 treatments

After 4-5 treatments, you should see a significant reduction of hair growth in the areas that you treated with Lumea. A reduction of hair density should be visible as well. Keep on treating with frequent touch-ups (every 4 weeks) to maintain the result.

Tanning advice

Tanning with natural or artificial sunlight

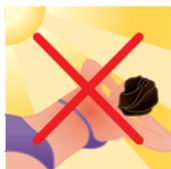Intentionally exposing your skin to natural or artificial sunlight with the aim of developing a tan influences the sensitivity and color of your skin.

Therefore the following is important:

- After each treatment, wait at least 48 hours before tanning. Even after 48 hours, make sure that the treated skin does not show any redness from the treatment anymore.
- In case of exposing your skin to the sun (without tanning intentionally) in the 48 hours after treatment, use a sunblock SPF 50+ on the treated areas. After this period, you can use a sunblock SPF 30+ for two weeks.
- After tanning, wait at least 2 weeks before you use Lumea.

- After recent tanning, perform a skin test to determine the appropriate light intensity setting. For instructions see chapter 'Skin test'.
- Do not use Lumea on sunburned body areas.

Note: Occasional and indirect sun exposure does not qualify as tanning.

Tanning with creams

If you have used an artificial tanning lotion, wait until the artificial tan has disappeared completely before you use the device.

Before you use your Lumea

Pretreating your skin

Before you use Lumea, you should pretreat your skin by removing hairs on the surface of your skin. This allows the light to be absorbed by the hair parts below the skin surface to ensure effective treatment. You can either shave, short-trim, epilate or wax. Do not use depilatory creams, as chemicals may cause skin reactions.

If you choose to wax, please wait 24 hours before using Lumea to let your skin rest. We recommend that you take a shower before the treatment to ensure that all possible residue of wax has been removed from your skin.

- 1 Pretreat the areas you intend to treat with Lumea.
- 2 Clean your skin and make sure it is hair-free, entirely dry and free from oily substances.

Note: if shaving causes skin irritation, we advise you not to use the device until skin irritation is resolved.

Skin test

When you use Lumea for the first time or after recent tanning, perform a skin test on each area to be treated. The skin test is necessary to check your skin's reaction to the treatment and to determine the correct light intensity setting for each body area.

- 1 Choose an area close to the area you intend to treat.

Note: do not try the device on difficult or sensitive areas (ankle and bony area).

- 2 Turn on the device. Make sure you select setting 1.
- 3 Place the device at a 90° angle on the skin so that the integrated safety system is in contact with your skin.

The integrated safety system prevents unintentional flashing without skin contact.

- 4 Press the flash button to release a flash.
- 5 Slide the device over the skin to the next area to be treated.
- 6 Increase the setting by one level, apply a flash and slide the device to the next area. Repeat this for all levels within the recommended range for your skin type.
- 7 After the skin test, wait 24 hours and check your skin for any reaction. If your skin shows reactions, choose the highest setting that did not result in any skin reaction for subsequent use.

Note: Perform a skin test if you notice a skin tone change from last treatment and wait 30 minutes before treatment.

Using your Philips Lumea

Skin tone sensor

For extra safety, the integrated skin tone sensor measures the skin tone at the beginning of each session and occasionally during the session. If it detects a skin tone that is too dark for treatment with Lumea, the 'ready to flash' light starts blinking orange and the device automatically disables to prevent you from developing skin reactions. This means that it does not emit pulses when you press the flash button.

Selecting the right light intensity

Lumea provides 5 different light intensities. Depending on your skin type and on what you find comfortable, you can select the right light intensity setting.

- 1 Consult the table below to select the right setting.
- 2 To adjust the light intensity setting, press the on/off button one or more times until you have reached the required setting. The corresponding intensity light level lights up green. You can change the setting at all times.

Note: The integrated skin tone sensor does not automatically select a light intensity setting for you, but only warns you when your skin is too dark for treatment.

Note: The device automatically disables when your skin tone is too dark, to prevent you from developing skin reactions. The 'ready to flash' light will blink orange to indicate if the skin tone is too dark.

- 3 Lumea should never be painful. If you experience discomfort, reduce the light intensity setting. You can do this by using the on/off button.
- 4 After moving to another body area or recent tanning, perform a skin test to determine the right light intensity setting.


Note: Since body areas might differ in skin color, you need to select the right setting for each body area separately. The Philips Lumea does not automatically select the correct setting.

Skin type	Skin tone	Light intensity setting
I	White: you always sunburn, never tan.	4/5
II	Beige: you easily sunburn, tan minimally.	4/5
III	Light brown: you sunburn easily, tan slowly to light brown.	3/4
IV	Mid brown: you rarely sunburn, tan easily.	1/2/3
V	Dark brown: you rarely sunburn, tan very easily.	You cannot use the device
VI	Brownish black or darker: you rarely or never sunburn, very dark tanning.	You cannot use the device

Note: Your skin may react differently on different days/occasions for a number of reasons.

Handling the device

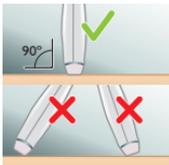- Before use, clean the light exit window.
- Put the adapter in the wall socket.
- Switch on the device and select the right light intensity for your skin tone. See chapter: 'Selecting the right light intensity'.
- Place the device at a 90° angle on the skin so that the integrated safety system is in contact with your skin.
The integrated safety system prevents unintentional flashing without skin contact.
- Press the device firmly onto your skin to ensure proper skin contact. The 'ready to flash' light lights up green to indicate that you can proceed with the treatment.

Note: When the 'ready to flash' light starts blinking orange, your skin is not suitable for treatment. In this case, the device automatically disables. You can try to use the device on other areas with a lighter skin tone.

Note: If the 'ready to flash' light does not light up green, the device is not completely in contact with your skin.

TIP: If you use Lumea on your bikini area and want to leave some hair in place, make sure that the skin tone sensor is not placed on the

remaining hair during the treatment, as this can cause the device to block

- 6 Press the flash button to release a flash. The device makes a soft popping sound. You should feel a warm sensation because of the flash.

Note: The scattered light produced by the device is harmless to your eyes. It is not necessary to wear goggles during use. Use the device in a well-lit room so that the light is less glaring to your eyes.

- 7 Place the device on the next area to be treated. After each flash, it takes up to 2.9 seconds until the device is ready to flash again. You can release a flash when the 'ready to flash' light lights up green.

- 8 To make sure that you have treated all areas, release the flashes close to each other. Effective light only comes out of the light exit window. Part of the device that is in contact with your skin is slightly larger therefore there should be some overlap. However, make sure you flash the same area only once.

Note: Do not use any pencil or pen to mark the areas to be treated, as this may cause side effects to your skin.

Note: Flashing the same area twice does not improve the effectiveness of the treatment, but increases the risk of skin reactions.

- 9 When you have finished the treatment, press and hold the on/off button for 2 seconds to switch off the device. Remove the adapter from the wall socket.

Two treatment modes: Stamp & Flash and Slide & Flash

Your Philips Lumea has two treatment modes for more convenient use on different body areas:

- The Stamp & Flash mode is ideal to treat small or curvy areas like knees and underarms. Simply press and release the flash button to release a single flash.

- The Slide & Flash mode offers convenient use on larger areas like legs. Keep the flash button pressed while you slide the device over your skin to release several flashes in a row.

Note: The device needs about 2.9 seconds in between two flashes.

Recommended treatment schedule

Initial phase

For the first 4 to 5 treatments, we advise you to use Philips Lumea once every two weeks to ensure that all hairs are treated.

Note: Replacing one of the IPL treatments with another hair removal method (waxing, epilating, etc) will not help to reach the desired hair reduction.

Touch-up phase


After the initial treatment phase (4-5 treatments), we recommend touch-ups every 4 weeks, when you see hairs growing back. This is to maintain results and enjoy smooth skin for months. The time between treatments may vary based on your individual hair regrowth and also across different body areas.

Tip: You can write the treatment schedule in your agenda to remind yourself to the treatments so you will not forget.

Note: Using the device more often does not enhance the effectiveness.

Treatment time per area

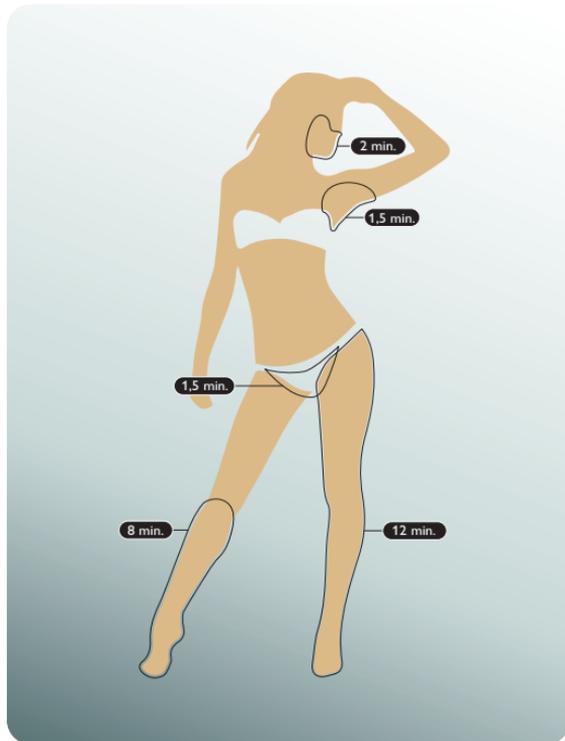After use

Common skin reactions

Your skin may show slight redness and/or may prickle, tingle or feel warm. This reaction is absolutely harmless and disappears quickly.

Dry skin and itching may occur due to shaving or a combination of shaving and light treatment. You can cool the area with an ice pack or a wet face cloth. If dryness persists, you can apply a non-scented moisturizer on the treated area.

Rare side effects

- Burns, excessive redness (e.g. around hair follicles) and swelling: these reactions occur rarely. They are the result of using a light intensity that is too high for your skin tone. If these reactions do not disappear within 3 days, we advise you to consult a doctor. Wait with the next treatment until the skin has healed completely and make sure you use a lower light intensity.
- Skin discoloration: this occurs very rarely. Skin discoloration manifests itself as either a darker patch (hyperpigmentation) or a lighter patch (hypopigmentation) than the surrounding area. This is the result of using a light intensity that is too high for your skin tone. If the discoloration does not disappear within 2 weeks, we advise you to consult a doctor. Do not treat discolored areas until the discoloration has disappeared and your skin has regained its normal skin tone.
- Skin infection is very rare but is a possible risk following a (micro)wound, a skin burn, skin irritation etc.

- Epidermal heating (a sharply defined brownish area which often occurs with darker skin tones and is not accompanied with skin dryness): This reaction occurs very rarely. In case this reaction does not disappear within 1 week, we advise you to consult a doctor. Wait with the next treatment until the skin has healed completely and make sure you use a lower light intensity.
- Blistering (looks like small bubbles on the surface of the skin): this occurs very rarely. In case this reaction does not disappear within 1 month or when the skin gets infected, we advise you to consult a doctor. Wait with the next treatment until the skin has healed completely and make sure you use a lower light intensity.
- Scarring: often the secondary effect of a burn, which can take longer than a month to heal.
- Folliculitis (swelling around hair follicles combined with pustule formation): this reaction occurs very rarely and is the result of bacteria penetrating the damaged skin. In case this reaction occurs, we advise you to consult a doctor as folliculitis may need antibiotic ointment.
- Excessive pain: this can occur during or after treatment if you have used the device on skin that is not hair-free, if you use the device at a light intensity that is too high for your skin tone, if you flash the same area more than once and if you use the device on open wounds, inflammations, infections, tattoos, burns, etc.


Aftercare

After use, you can safely apply lotions, creams, deodorant, moisturizer or cosmetics to the treated areas. If you experience skin irritation or skin redness after treatment, wait until it disappears

before applying any product to your skin. If you experience skin irritation after applying a product to your skin, wash it off with water.

Note: when following all instructions you may still get skin reactions. In this case stop using the device and contact the consumer care center in your country.

Cleaning & storage

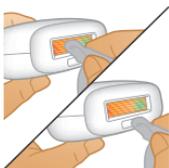- 1 After use, switch off the device, unplug it and let it cool down.
- 2 Moisten a soft cloth with a few drops of water and use it to clean the following parts:
 - the light exit window including the reflectors on the side
 - the skin tone sensor
- 3 Store the device in a dust free and dry place at a temperature between -25°C and 75°C , with a storage humidity up to 75%.

Guarantee and support

If you need information or support, please visit www.philips.com/support or read the separate worldwide guarantee leaflet.

Recycling

- This symbol means that this product shall not be disposed of with normal household waste (2012/19/EU).
- Follow your country's rules for the separate collection of electrical and electronic products. Correct disposal helps prevent negative consequences for the environment and human health.

Technical specifications

BRI861, BRI860

Rated voltage	100V-240V	
Rated frequency	50Hz-60Hz	
Rated input	20W	
Protection against electric shock	Class II	
Protections rating	IP 30 (EN 60529)	
Operating conditions	Temperature: 5 °C to 35 °C Relative humidity: 15% to 75%	
Spectrum	>530nm	
Maximum variation of the light output over the treatment area	<20%	
Operating pressure	700hPa-1060hPa	
Altitude	max. 3000m	
Pulse duration	0.8 +/- 0.2 mSeconds	
Duration of the pulsetrain	Single pulse	
Pulse interval	Min 2.9 Second	
	The minimal time between flashes	Maximum output of radiation
Setting 1	1.3s	5.76 +/- 1.28 J
Setting 2	1.7s	7.68 +/- 1.28 J
Setting 3	2.3s	10.24 +/- 1.28 J
Setting 4	2.6s	11.52 +/- 1.28 J
Setting 5	2.9s	12.8 +/- 1.28 J

Troubleshooting

This chapter summarizes the most common problems you could encounter with the device. If you are unable to solve the problem with the information below, visit www.philips.com/support for a list of frequently asked questions or contact the Consumer Care Center in your country.

Problem	Possible cause	Solution
The device/adapter becomes warm during use.	This is normal.	No action required.
When I place the device on the skin, it does not release a flash. The 'ready to flash' light blinks orange.	Your skin tone in the area to be treated is too dark.	Treat other body areas with lighter skin tones with Lumea.
The ready to flash light blinks orange and all 5 intensity lights blink as well.	The device needs to be reset.	To reset the device take the plug out of the socket, wait for 30 minutes to let the device cool down. The device should function normally again. In case it doesn't work again, contact the Consumer Care Center in your country.
The ready to flash light does not light up green.	The device is not completely in contact with your skin.	Place the device at a 90 degree angle on the skin so that the integrated safety system is in contact with your skin.

Problem	Possible cause	Solution
The device produces a strange smell.	The light exit window or the skin tone sensor is dirty.	Clean the light exit window and the skin tone sensor carefully.
	You have not removed the hairs on the area to be treated properly. These hairs may get burned and can cause the smell.	Pretreat your skin before you use Lumea.
The skin feels more sensitive than usual during treatment. I experience an unacceptable pain sensation when I use the device.	The light intensity setting you use is too high.	Check if you have selected the right light intensity setting. If necessary, select a lower setting.
	You did not remove the hairs on the areas to be treated.	Pretreat your skin before you use Lumea.
	The UV filter of the light exit window is broken.	If the UV filter is broken, do not use the device anymore. Contact the Consumer Care Center in your country, your Philips dealer or a Philips service center.
	You treated an area for which the device is not intended.	Never use the device on the following areas: inner labia, anus, nipples, areolas, lips, moles, freckles, tattoos, piercings, inside the nostrils and ears, around the eyes and near the eyebrows. Men must not use it on the face, neck or on the scrotum.

Problem	Possible cause	Solution
The skin reaction after treatment lasts longer than usual.	You have used a light intensity setting which is too high for you.	Select a lower intensity next time. See chapter 'Using your Philips Lumea', section 'Selecting the right light intensity'.
The results of the treatment are not satisfactory.	You have used a light intensity setting which is too low for you.	Select a higher setting next time.
	You did not flash an area adjacent to an area you treated before.	You have to release the flashes close to each other and there should be some overlap between the flashes.
	The device is not effective on your body hair color.	If you have light blond, grey, red or white hair, the treatment is not effective.
	You do not use the device as often as recommended.	To remove all hairs successfully, we advise you to follow the recommended treatment schedule. You can reduce the time between treatments, but do not treat more often than once every two weeks.
	You respond more slowly to IPL treatment.	Continue using the device for at least 6 months, as hair regrowth can still decrease over the course of this period.

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©2016 Koninklijke Philips N.V. All rights reserved. Specifications are subject to change without notice. Trademarks are the property of Koninklijke Philips N.V. (Royal Philips) or their respective owners.