

**PHILIPS**

# SatinLux

The next generation in hair removal



Introduction	5
Key benefits of SatinLux	6
Who can use SatinLux?	7
How does SatinLux work?	9
Treatment areas	12
Treatment time	13
The product overview	14
Warnings	16
FAQs	18



All current hair removal methods have one common downside: hair needs to grow back before you can remove it. Philips now grants women's ultimate beauty wish: to feel free from unwanted body hair and to enjoy smoothness forever.

**Philips launches SatinLux, the first at home photo-epilator for effective hair re-growth prevention**



### Professional IPL technology for home use

SatinLux uses a light-based technology called 'Intense Pulsed Light' (IPL). IPL is also used in the professional beauty market for hair removal.

### Effective hair re-growth prevention for everlasting smoothness

This method prevents hair from growing back. Philips Clinical Studies show a significant hair reduction already after 2 treatments and optimal results were generally achieved after 4-5 treatments.

### Gentle treatment, even on sensitive body areas

SatinLux has been developed in close cooperation with leading dermatologists for an effective and gentle treatment, even on the most sensitive body areas.

Light based depilation works on almost all skin tones and body hair colours, but is not effective on very light body hair (light blond, red, grey and white). Light based depilation is not suitable for black skin.

	White (always sunburn, no tanning)	Beige (easy sunburn, minimal tanning)	Light brown (sometimes sunburn, good tanning)	Mid brown (rarely sunburn, fast and good tanning)	Dark brown (rarely sunburn, very good tanning)	Black (rarely to never sunburn, very dark tanning)
White/grey	x	x	x	x	x	x
Red	x	x	x	x	x	x
Light blond	x	x	x	x	x	x
Blond	4/5*	4/5*	4/5*	3/4/5*	1/2/3/4*	x
Dark blond	4/5*	4/5*	4/5*	3/4/5*	1/2/3/4*	x
Brown	4/5*	4/5*	4/5*	3/4/5*	1/2/3/4*	x
Dark brown	4/5*	4/5*	4/5*	3/4/5*	1/2/3/4*	x
Black	4/5*	4/5*	4/5*	3/4/5*	1/2/3/4*	x

\*recommended setting to be used

### Preparation



#### Charge the battery

Fully charge the battery before you use the appliance for the first time and when the battery is empty.



#### Prepare your skin

For optimal results and highest efficacy, shave the areas you intend to treat before each treatment. This is only needed for the first few treatments as long as you still experience hair re-growth.

### In-use



#### Step 1

Press the on/off button to switch on the appliance.



#### Step 2

Press the intensity increase / decrease button up or down to change intensity (from 1 to 5).



#### Step 3

Place the appliance perpendicularly on the skin so that the removable window and safety ring are in contact with the skin.



#### Step 4

Press the safety ring fully onto the skin with slight pressure. When correct, the 'ready to flash' light on the back of the appliance goes on.



#### Step 5

Press the ready to flash button. Repeat steps 3, 4 and 5.



#### Step 6

To avoid untreated areas, always make sure there is some overlap.



#### Step 7

Switch off the appliance after you have finished the treatment.

SatinLux is intended for the use on underarms, bikini area and legs



Area	Approx. number of flashes	Approx. treatment time
Underarms	25 per underarm	3 minutes per underarm
Bikini line	25 per side	3 minutes per side
Entire bikini area	90 for entire area	10 minutes
One lower leg	160-190	15 minutes
One full leg	320-380	30 minutes

On average, women have approximately 16,800 hairs on their legs and 600 hairs under their arms. Naturally, it requires more treatment time (and possibly recharging), to cover the leg areas.

A fully charged battery offers at least 160 flashes at intensity 5, which amounts to a cordless operating time of approx. 15 minutes. On a lower setting (though less efficacious!) the battery could last up to 465 flashes.

The product overview



**'ready to flash' light**  
indicates a correct application on the skin

**rechargeable battery**  
for easy cordless handling

retail price € 499

**5 light energy settings**  
adjustable to skin type and personal comfort

**3 cm<sup>2</sup> application window**  
for an effective treatment of large areas

**removable window with UV-filter**  
protects the skin from UV light

**integrated safety system**  
prevents unintentional flashing

storage pouch and cleaning cloth included





Do not use the device under the following circumstances or body areas.



Facial area



On moles,  
freckles and scars



On tattoos  
or piercings



After recent  
sunbath



On light body  
hair colours



On dark skin




On nipples



Children



Pregnant

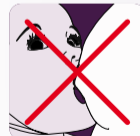

Breastfeeding

**Is it suitable for all skin tones and hair types?**

SatinLux is suitable for most women. All light-based depilation methods work best on dark blond, brown or black body hair and are not effective on light blond, red, grey and white hair. Light-based depilation is not suitable for black skin.

**Why do I have to pre-shave before using SatinLux?**

Shaving is necessary to ensure that the light energy is channeled to the hair root in the most efficient way. Furthermore, visible hairs in the treatment area can be burnt by the light and this can feel and smell uncomfortable.

**Before using SatinLux, can the re-grown hair be epilated or waxed instead of shaved?**

No. The short hair root is necessary to channel the light energy to the follicle. An epilator or wax would remove the hair completely.

**Is it necessary to still use the appliance every 2 weeks, even after having reached an almost full hair free result?**

Yes, it is. It is proven that to treat all hairs successfully and to prevent that the hair root becomes active again, you have to repeat the treatment every 2 weeks. Better results are not expected if you increase or decrease the treatment frequency.

**Does it really work/is it really effective?**

Our clinical studies show a significant hair reduction after 2 treatments and optimal results were generally achieved after 4-5 treatments. At this stage, nearly all women were already extremely satisfied with this result and considered themselves as being smooth.



© 2008 Koninklijke Philips Electronics N.V.