

# PHILIPS

IPL - Hair removal device

Lumea Advanced

For use on body, face, bikini  
15 minutes to treat lower legs  
>250,000 lamp flashes  
Extra long cord

SC1999/00



## Gentle, long-term hair reduction

Philips Lumea Advanced works wonders to prevent the reappearance of hair on body and face. Gentle pulses of light, applied regularly, keep skin silky-smooth every day.

### Long lasting hair removal with IPL technology

- Enjoy smooth skin every day
- Effortlessly effective
- Adapted for safe and effective use in your own home
- Suitable for a wide variety of hair and skin types

### Safe, fast and effective on body, face and bikini

- Safe and effective even on sensitive areas
- Large treatment window for fast application
- Precision attachment for extra safety on the face
- Integrated skin tone sensor
- Gentle treatment
- Attachment for specialised treatment of bikini area hair

### Maximum freedom, minimum hassle

- No replacement parts, no hidden costs
- Slide and Flash mode for a more convenient application


# Highlights

## IPL technology



Philips Lumea uses an innovative light-based technology called IPL (Intense Pulsed Light), derived from the technology used in professional beauty salons. Philips has adapted this technology for safe and effective use in the comfort of your own home. Philips worked closely with leading dermatologists to develop our breakthrough hair removal system. For more than 10 years we carried out extensive consumer research with more than 2000 volunteers.

## Enjoy smooth skin every day

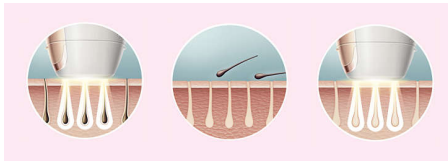

Philips Lumea applies gentle pulses of light to the hair root. As a consequence, the amount of hair your body grows gradually decreases. Repeating the treatment leaves your skin beautifully hair free and touchably smooth.

## Effortlessly effective

Clinical results after treatment: **example 1 (leg)**



Clinical results after treatment: **example 2 (bikini line)**



Our clinical studies have shown significant hair reduction in as little as four fortnightly treatments resulting in smooth hair free skin. To maintain these results, simply repeat the treatment when needed. Time between

treatments may vary based on your individual hair regrowth. - To enlarge this picture, please click on the image in the image gallery on the top of this page

## Suitable hair and skin types




Philips Lumea works effectively on (naturally) dark blonde, brown and black hairs. As with other light-based treatments, Philips Lumea is NOT effective on red, light-blonde or white/grey hair. This product is also NOT suitable for dark skin.

## Application areas



Use with confidence to remove hair on face (upper lip, chin and sideburns) and body parts including legs, underarm, bikini area, belly and arms.

## Fast treatment on legs

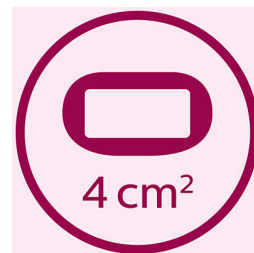

Big treatment window for fast application on large areas like legs.

## Extra safety for the face



With the additional integrated light filter, it can be used for safe treatment of facial hair on the upper lip, chin, cheeks and other sensitive areas.

## Bikini area attachment



Bikini area attachment for effective treatment of bikini hair. Hair in this area tends to be typically stronger and thicker than leg hair.

## Skin Tone Sensor



The integrated skin tone sensor measures the treated skin complexion at the beginning of each session and occasionally during the session. If it detects a skin tone that is too dark for this product it will automatically stop emitting pulses.

# Specifications

## Power

**Voltage:** 100-240 V

## Service

**Warranty:** 2 years global warranty

## Technical specifications

**High performance lamp:** Delivers <gt/> 250,000 flashes

## Application areas

**Body areas:** Arms, Legs, Belly, Bikini, Armpits

**Face areas:** Chin, Sideburns, Upper lip

## Attachments

**Precision attachment (2 cm<sup>2</sup>):** For safe application on face

**Body attachment (4 cm<sup>2</sup>):** For use below the

neckline

**Bikini area attachment (2 cm<sup>2</sup>):** Special treatment bikini hair

## Items included

**Instructions for use:** User manual

**Storage:** Pouch

**Adapter:** 24V / 1500mA

**Accessory:** Cleaning cloth

## Power

**Corded use**

## Safety and adjustable settings

**Integrated UV filter:** Protects skin from UV light

**5 light energy settings:** Adjustable to your skin type

**Integrated safety system:** Prevent unintentional flashing

**Skin Tone sensor:** Detects your skin tone

## Application time

**Bikini:** 4min

**Lower leg:** 8 min

**Armpit:** 1 min

**Face areas:** 2 min

## Application mode

**Slide and Flash:** For convenient application

**Corded/cordless use:** Corded

**Stamp and Flash:** For treatment on curvy areas

