GB Call FREE on **0800 289 064** (GB Only)

Call FREE on **1800 509 021** (IRL Only)

Call on **1 300 364 474** (AUS)

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www.philips.com/AVENT

AVENT, Suffolk, CO10 7QS, England

Philips of Holland. All logos brands and product names are trademarks or registered
trademarks of their respective holders and are hereby recognised as such.

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HR ISIS Ručna izdajalica

CLINICALLY PROVEN
AS EFFECTIVE AS HOSPITAL
ELECTRIC PUMPS*

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







* Randomised Trial Comparing the
Efficacy of a Novel Manual Breast
Pump With a Standard Electric
Breast Pump in Mothers Who
Delivered Preterm Infants.
M S Fewtrell, P Lucas et al.
MRC Childhood Nutrition
Research Centre, London.
Paediatrics June 2001.

2.

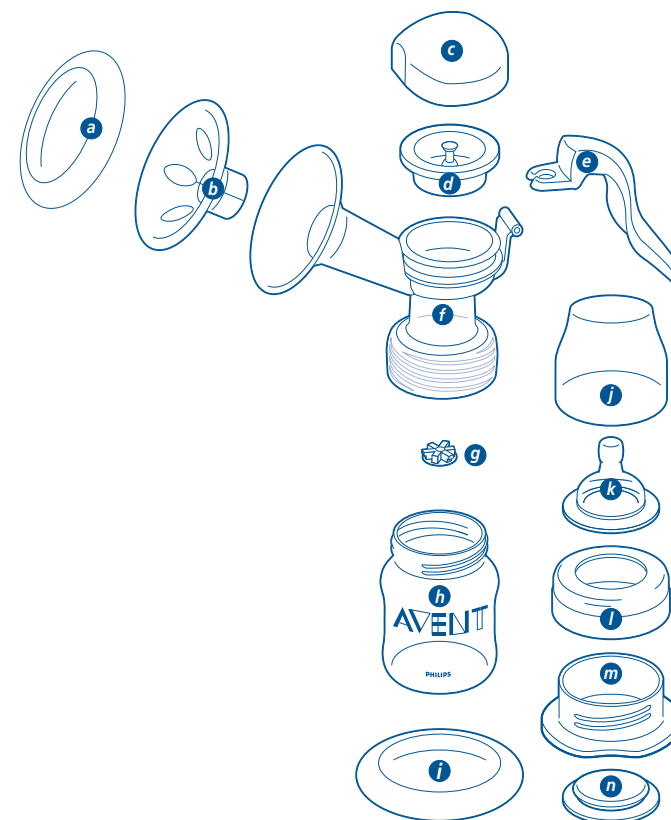6. (1) – (7)

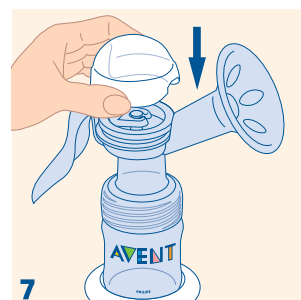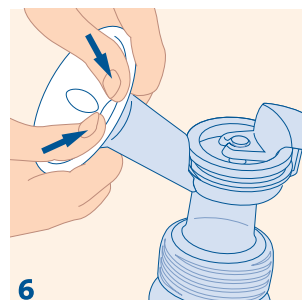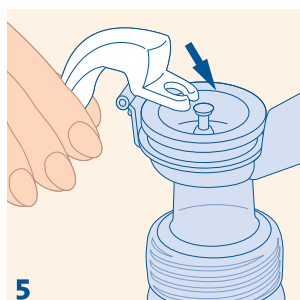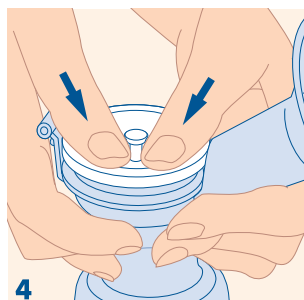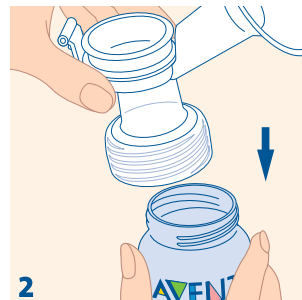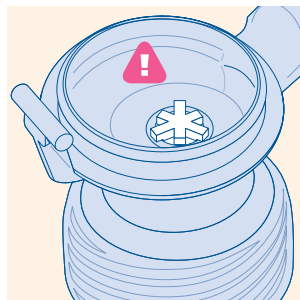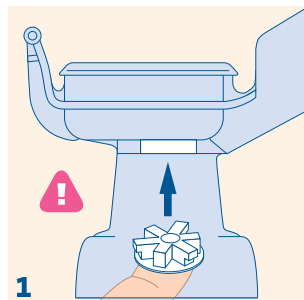IMPORTANT: DO NOT LOSE THE WHITE VALVE. YOUR PUMP WILL NOT WORK WITHOUT IT OR IF IT IS FITTED INCORRECTLY.
If lost, spare valves are available direct from AVENT.

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GB

Please take five minutes to read through the **WHOLE** of this leaflet **BEFORE** you use your Philips AVENT ISIS Manual Breast Pump for the first time.

Contents

1. Introducing the ISIS Manual Breast Pump
2. Detailed parts guide
3. Cleaning and sterilisation
4. Keeping your breast pump sterile for travel
5. Using the ISIS Manual Breast Pump with the AVENT VIA Feeding System
6. How to assemble the ISIS Manual Breast Pump
7. When to express milk
8. Hints to help you succeed
9. How to use the ISIS Manual Breast Pump
10. Storing your breast milk
11. Feeding your baby your expressed breast milk
12. Choosing the right teat for your baby
13. Troubleshooting guide
14. Other helpful products from Philips AVENT

Should you ever lose these instructions, they are available on our website: www.philips.com/AVENT

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1. Introducing the ISIS Manual Breast Pump

CLINICALLY PROVEN AS EFFECTIVE AS HOSPITAL ELECTRIC PUMPS*

The ISIS Manual Breast Pump is a breakthrough in breast pump design. Its silicone diaphragm provides 100% reliable vacuum, whilst the Let-down Massage Cushion naturally imitates your baby's suckling action to provide fast milk flow – silently, comfortably and gently. The pump is easy to assemble, discreet to use and all parts are dishwasher safe.

Doctors recommend that breast milk is the best nutrition for babies during the first year, combined with solid food after the first 6 months. Your milk is specially adapted to your baby's needs and contains antibodies which help protect your baby against infection and allergies.

Your ISIS Manual Breast Pump will help you breastfeed longer. You can express and store your milk so that – whether you are out for a few hours, taking a well deserved rest, or back at work – your baby can still enjoy the benefits of your breast milk, even if you can't be there to give it yourself. As the pump is so silent and discreet you can take it with you anywhere, allowing you to express milk at your own convenience and maintain your milk supply.

Since the AVENT Bottle and Airflex Teat are specially designed to mimic breastfeeding, your baby will switch easily between breast and bottle. This allows mothers to continue providing breast milk to baby even when they return to work.

* Randomised Trial Comparing the Efficacy of a Novel Manual Breast Pump With a Standard Electric Breast Pump in Mothers Who Delivered Preterm Infants. **M S Fewtrell, P Lucas et al. MRC Childhood Nutrition Research Centre, London. Paediatrics June 2001.**

2. Detailed parts guide

ISIS MANUAL BREAST PUMP

- a) Funnel cover (same as i)
- b) Let-down Massage Cushion
- c) Pump cover
- d) Silicone diaphragm and stem*
- e) Handle
- f) Pump body
- g) White valve*

BOTTLE

- h) Airflex Feeding Bottle (125ml/4oz)
- i) Stand (same as a)
- j) Dome cap
- k) Extra Soft Teat, Newborn Flow (0m+)
- l) Screw ring
- m) Sealing base
- n) Sealing disc
- *spares included

3. Cleaning and sterilisation

Before you use your breast pump for the first time and every time before use:

Separate all parts, then wash in the dishwasher or in warm soapy water and rinse. Do not use antibacterial or abrasive cleaners/detergents. Sterilise in a Philips AVENT Steam Steriliser (Electronic, Electric or Microwave), or by boiling for 5 minutes. If using a Philips AVENT Steam Steriliser, the pump will remain sterile in the steriliser for a minimum of 6 hours, provided the lid is not removed.

4. Keeping your breast pump sterile for travel

After sterilising all parts, insert the white valve (g) and fit the diaphragm with stem (d) onto the pump body (f). Clip the cover (a) over the funnel and Let-down Massage Cushion (b) and screw on the sealing base (m) in place of the Airflex Feeding Bottle (h).

5. Using the ISIS Manual Breast Pump with the AVENT VIA Feeding System

Simply replace the Airflex Feeding Bottle with a VIA Cup and Adaptor. To make the pump completely stable insert the VIA Cup into the VIA Lid/Stand.

6. How to assemble the ISIS Manual Breast Pump

Separate all parts and ensure pump has been cleaned and sterilised as in section 3:

- 1) Wash your hands thoroughly. Insert the white valve (g) into the pump body (f) from underneath making sure that the flat side of the valve is facing downwards and the star side is facing up.
- 2) Place the pump body (f) onto the Airflex Feeding Bottle (h).
- 3) Twist gently clockwise until secure. **DO NOT OVERTIGHTEN!**
- 4) Place the silicone diaphragm and stem (d) into the pump body (f). Make sure it fits securely around the rim by pressing down with your fingers to ensure a perfect seal. (This is easier if assembled whilst wet.)
- 5) Place the forked end of the handle (e) under the silicone diaphragm and stem (d) and push down gently on the handle until it clicks into place.
- 6) Carefully insert the Let-down Massage Cushion (b) into the pump funnel (f) making sure it is perfectly sealed all round the rim of the funnel. (This is easier if assembled whilst wet.)
- 7) Click the pump cover (c) onto the pump body (f). To make the pump completely stable insert the bottle into the stand (i).

7. When to express milk

- 1) If possible, wait until your milk supply and breastfeeding schedule are well established, normally at least 2-4 weeks after birth, unless told otherwise by your doctor or breastfeeding advisor.

- 2) The exceptions to the above are:
 - a) If you are expressing milk for your baby to be given in hospital.
 - b) If your breasts are engorged (painful or swollen), you can express a small amount of milk before or between feeds to ease the pain and to help your baby latch on more easily.
 - c) If you have sore or cracked nipples, you may wish to express your milk until they are healed.
 - d) If you are separated from your baby and wish to continue to breastfeed when you are reunited, you should express your milk regularly to stimulate your milk supply.
- 3) You will need to find the right times of the day for you to express your milk, e.g.
 - a) Early in the morning when your breasts are full, just before or after your baby's first feed.
 - b) After a feed, if your baby has not emptied both breasts.
 - c) Between feeds, or in a break if you have returned to work.

8. Hints to help you succeed

Using a breast pump takes practice – you may need to make several attempts before you succeed, but because ISIS is so simple and natural to use, you will soon get used to expressing your milk.

- 1) Choose a time when you're not rushed and won't be interrupted.
- 2) Having your baby or a photograph of your baby nearby can help encourage 'let-down'.
- 3) Placing a warm cloth on the breasts for a few minutes before you start pumping can encourage milk flow and soothe painful breasts.
- 4) Warmth and relaxation can encourage milk flow. Try expressing after a bath or shower.
- 5) Try expressing from one breast while your baby is feeding from the other, or continue expressing just after a feed.

- 6) Repositioning the pump on the breast from time to time can help stimulate your milk ducts.

Practice with your pump to find the best technique for you. However, if the process becomes very painful, stop and consult your breastfeeding advisor.

9. How to use the ISIS Manual Breast Pump

- 1) Wash your hands thoroughly and make sure your breasts are clean. Gently squeeze a little milk from each nipple to ensure milk ducts are not blocked.
- 2) Make sure you have sterilised and carefully assembled your pump exactly as shown in sections 3 and 6.
- 3) Consult the 'Hints to help you succeed' section.
- 4) Relax in a comfortable chair, leaning slightly forward (use cushions to support your back). Make sure you have a glass of water nearby.
- 5) Press the funnel and Let-down Massage Cushion of the pump firmly against your breast, making sure that no air can escape – otherwise there will be no suction.
- 6) As you gently push down on the pump handle, you will feel the suction on your breast. You don't need to depress the handle fully to create a vacuum – only as much as is comfortable. Your milk will soon start flowing, even though you may not use all the suction the pump can generate.
- 7) Begin by pumping 5-6 times rapidly to initiate let-down. Then, hold the handle down for 2-3 seconds, and allow it to return to its resting place. These 2-3 second cycles imitate your baby's natural suckling pattern and allow the milk to flow between strokes.
- 8) Milk should begin to flow after the first few strokes. Don't worry if your milk does not flow immediately. Relax and continue pumping.

Do not continue pumping for more than 5 minutes at a time if NO result is achieved. Try expressing at another time during the day.

Please remember that, should the process become very painful or uncomfortable, you should stop using the pump and consult your health professional or breastfeeding advisor.

- 9) Whilst some women may prefer to use the pump without the Let-down Massage Cushion, extensive trials show that, in most cases, the massaging effect of the petals flexing in and out around the areola stimulates the let-down reflex and allows you to express more milk, faster.
- 10) On average you will need to pump for 10 minutes to express 60-125ml or 2-4oz. However, this is just a guide and varies from woman to woman. If you find that you regularly express more than 125ml/4oz per session, please use a 260ml/9oz AVENT Bottle or 240ml/8oz VIA Cup.


Do not overfill Airflex Feeding Bottles or VIA Cups and keep upright, otherwise milk may leak from below the pump body.

10. Storing your breast milk

- Breast milk can be stored in the fridge or freezer in either sterile breast milk containers (125ml/4oz or 260ml/9oz bottles fitted with sealing discs) or in VIA Cups with Lids.
- Only milk collected with a sterile pump should be stored to be fed to your baby.
- Expressed milk should be refrigerated immediately.
- Expressed breast milk can be stored in the refrigerator for up to 48 hours (not in the door) or in the freezer for up to 3 months.
- If you are storing milk in the fridge to add to during the day, only add milk which has been expressed into a sterile container. Milk can only be kept like this for a maximum of 48 hours (from first expression) after which it must be used immediately or frozen for future use.
- For products other than breast milk make feeds up fresh when needed, following the manufacturer's instructions.

Storing your milk in the fridge

If your expressed breast milk will be fed to your baby within 48 hours, it can be stored in the fridge (not in the door) in an assembled Airflex

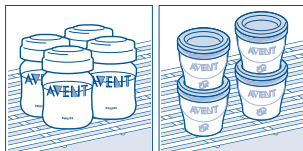Feeding Bottle or VIA Feeding System. Assemble a sterilised teat, screw ring and dome cap and attach onto the body of the Airflex Bottle or the VIA Adaptor and Cup (see below).

Storing your milk in the freezer

To store in the freezer, use a sealing disc in place of the teat and dome cap or use a sterilised VIA Lid to seal the VIA Cup. The storage containers should be clearly labelled with the date and time of expression and older milk used up first.

NEVER refreeze breast milk. DO NOT add fresh breast milk to already frozen milk. ALWAYS discard any milk that is left over at the end of a feed.

11. Feeding your baby your expressed breast milk

If breast milk has been stored in an AVENT Breast Milk Container, replace the sealing disc with a sterilised Airflex Teat.

Alternatively replace the VIA Lid with a sterilised VIA Adaptor, screw ring, dome cap and Airflex Teat.

Frozen breast milk should be defrosted overnight in the refrigerator. However, in an emergency, it can be defrosted in a jug of hot water. Once thawed use within 24 hours.

Breast milk can be heated in the same way as any other bottle feed. Stand the Airflex Bottle or VIA Feeding System in a bowl of hot water, or alternatively, for a safe and fast way to warm feeds, use the Philips AVENT Bottle and Baby Food Warmer. IMPORTANT! Milk or food will heat faster in VIA Cups than in standard feeding bottles/containers. Always check milk or food temperature before feeding.


Never heat milk in a microwave as this may destroy valuable nutrients and antibodies. In addition, the uneven heating may cause hot spots in the liquid which could cause burns to the throat.

Never immerse a frozen milk container or VIA Cup in boiling water as this may cause it to crack/split.

Never put boiling water directly into a VIA Cup, allow to cool for around 20 minutes before filling.

12. Choosing the right teat for your baby

AVENT Airflex Teats can be used on both the Airflex Feeding Bottle and VIA Feeding System. AVENT Teats are clearly numbered on the side to indicate flow rate.

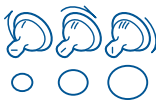0m+		NEWBORN FLOW EXTRA SOFT SILICONE TEAT Ideal for newborns and breastfed babies of all ages.
1m+		SLOW FLOW EXTRA SOFT SILICONE TEAT Ideal for newborns and breastfed babies of all ages.
3m+		MEDIUM FLOW SOFT SILICONE TEAT For bottle fed babies at 3 months of age and up.
6m+		FAST FLOW SOFT SILICONE TEAT Offers a faster flow for older babies.
3m+		VARIABLE FLOW SOFT SILICONE TEAT – SLOT CUT Provides an extra fast flow and accommodates thicker liquids. Flow rate can be varied by turning the bottle to align the I, II or III markings on the teat with baby's nose.

13. Troubleshooting guide

Lack of suction	Check that all the pump components have been assembled correctly and that the white valve is the right way up. Ensure that the diaphragm and Let-down Massage Cushion are fitted securely to the pump body and a perfect seal is created.
No milk being expressed	Ensure suction is created and the pump is correctly assembled. Relax and try again, expression may improve after practice. Please refer to the 'Hints to help you succeed' section.
Pain in the breast area when expressing	You may be pumping too hard. You do not need to use all the suction the pump can generate. Try depressing the handle halfway and remember the 2-3 second pumping rhythm. Consult your health professional/breastfeeding advisor.
Milk is drawn up the Let-down Massage Cushion	Remove and refit the Let-down Massage Cushion ensuring it is firmly attached and pushed far enough down the funnel. Try leaning slightly forward.
Cracking/discolouration of pump parts	Avoid contact with antibacterial or abrasive detergents/cleaners as these can damage the plastic. Combinations of detergents, cleaning products, sterilising solution, softened water and temperature fluctuations may, under certain circumstances, cause the plastic to crack. If this occurs, do not use. Contact Philips AVENT Customer Services for spare parts. The pump is dishwasher safe but food colourings may discolour components.
Lost or broken parts	Spare parts are available through Philips AVENT Customer Services (see back cover for contact details)

If you are still experiencing problems, please contact Philips AVENT Customer Services. (see back cover for contact details). Alternatively contact your health professional/breastfeeding advisor.

14. Other helpful products from Philips AVENT

Philips AVENT ISIS Out & About Manual Breast Pump Set – for expressing and transporting breast milk.

AVENT Breast Milk/Baby Food Containers – for storage of breast milk or baby food in the fridge or freezer.

AVENT VIA Breast Milk Containers – for expressed milk storage and transport.

AVENT ISIS Breast Shell Set (2 pairs) with Ultrasoft Backing Cushion – to protect sore/cracked nipples, help ease engorgement and collect leaking breast milk.

AVENT Nipple Protectors – to protect sore/cracked nipples during breastfeeding (small and standard sizes).

AVENT Ultra Comfort Disposable Breast Pads – unique ULTRA absorbent core guarantees total dryness and comfort.

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