

ENGLISH

Introduction

Congratulations on your purchase and welcome to Philips!
To fully benefit from the support that Philips offers, register your product at www.philips.com/welcome.

IMPORTANT SAFEGUARDS

READ ALL INSTRUCTIONS BEFORE USE

- Dangers**
 - Do not immerse, place or store the charger over or near water: Do not reach for a charger that has fallen into water: Unplug immediately. After cleaning, make sure the charger is completely dry before you connect it to the mains.
- Warnings**
 - The mains cord cannot be replaced. If the mains cord is damaged, discard the charger. Always have the charger replaced with one of the original type in order to avoid a hazard.
 - Discontinue use if this product appears damaged in any way (brush head, charger, toothbrush handle). Cord, plugs and batteries cannot be replaced. This device contains no user-serviceable parts.
 - Do not use the charger outdoors or near heated surfaces.
 - Do not use attachments and/or brush heads other than those recommended by the manufacturer.
 - This device is not intended for use by children or persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless supervised or instructed in use by a person responsible for their safety.
 - Children should be supervised to ensure that they do not play with the device.
- Electromagnetic fields (EMF)**

This device complies with current standards regarding electromagnetic fields (EMF). If used as directed this device is safe.

- The Sonicare (Fig. 1)**
 - 1 Replaceable brush head
 - 2 Power on/off button
 - 3 Handle
 - 4 Hygienic travel cap
 - 5 Charger

Getting started

Attaching the brush head

- 1 Align the front of the brush head with the front of the toothbrush handle.
- 2 Push the brush head onto the metal shaft. (Fig. 2)

Note: When the brush head is attached, it is normal that there is a small gap between the brush head and the handle.

Charging the Sonicare



- 1 Put the plug of the charger in the wall socket.
- 2 Place the handle on the charger.

Note: To keep the batteries fully charged at all times, you may keep your toothbrush on the charger when not in use. It takes approximately 30 hours to fully charge the batteries. When the toothbrush is used in clinical studies, the handle should be fully charged.

Using the Sonicare

- 1 Wet the bristles and apply a small amount of toothpaste (Fig. 3).
- 2 To avoid splattering, place the brush head in your mouth before switching the toothbrush on.
- 3 Press the power on/off button to switch on the toothbrush.
- 4 Brush like you normally would do. (Fig. 4)

Note: To achieve optimal results, thoroughly brush your entire mouth. Brush the outside, the inside and chewing surfaces of your teeth, as well as your gums.

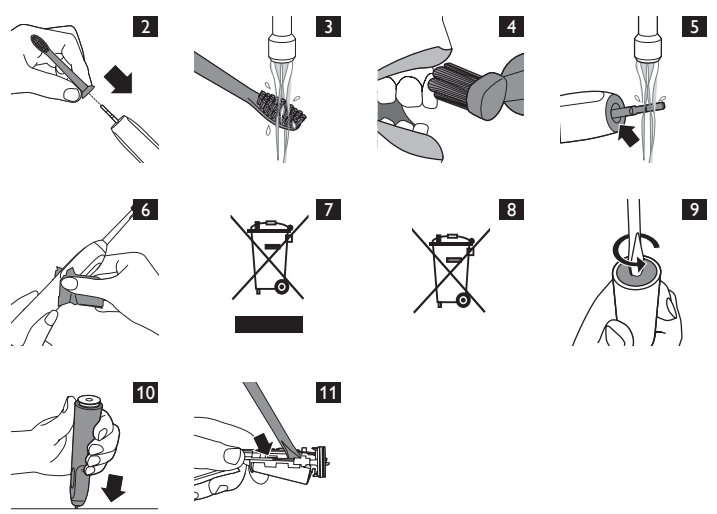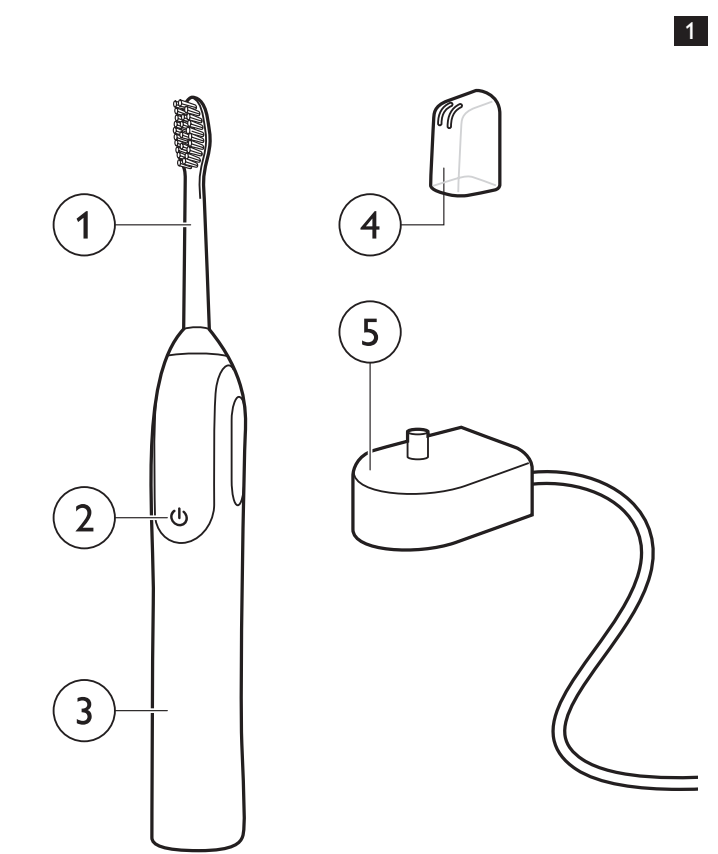

www.philips.com

100% recycled paper
100% papier recyclé
4235.020.5362.3

Storage

If you are not going to use the toothbrush for a long time (a month or longer), unplug the charger, clean the device and store it in a cool and dry place away from direct sunlight.

Replacement

Brush head

- Replace Sonicare brush heads every 3 months to achieve optimal results.
- Do not use a brush head with crushed or bent bristles. Replace the brush head every 3 months or sooner if signs of wear appear.
- Use only Philips Sonicare brush heads.

Note: Use of Sonicare for Kids brush heads are not recommended with this toothbrush.

Disposal

- Do not throw away the appliance with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment (Fig. 7).
- The built-in rechargeable batteries contain substances that may pollute the environment. Always remove the batteries before you discard and hand in the appliance at an official collection point. Dispose of the batteries at an official collection point for batteries. If you have trouble removing the batteries, you can also take the appliance to a Philips service centre. The staff of this centre will remove the batteries for you and will dispose of them in an environmentally safe way (Fig. 8).

Register your product and get support at www.philips.com/welcome



Philips Sonicare PowerUp HX3100 series

PHILIPS

Removing the rechargeable batteries

- 1 To deplete the rechargeable batteries of any charge, remove the handle from the charger, switch on the toothbrush and let it operate until it stops. Repeat this step until you can no longer switch on the toothbrush.
- 2 Insert a flat-head (standard) screwdriver into the slot located at the bottom of the handle. Turn the screwdriver anti-clockwise to release the bottom cap (Fig. 9).
- 3 Hold the handle upside down and push down on the shaft to release the internal components of the handle (Fig. 10).
- 4 Cut the wires to disconnect the batteries from the housing.
- 5 Place a flat-head (standard) screwdriver onto each of the three snaps and tap with a hammer to break the snaps. (Fig. 11)
- 6 Separate the bottom cap from the housing and remove the rechargeable batteries.
- 7 The rechargeable batteries can now be recycled and the rest of the appliance discarded appropriately.

Guarantee and service

Philips warrants its products for two years after the date of purchase. Defects due to faulty materials and workmanship will be repaired or replaced at Philips expense provided that convincing proof of purchase in the qualifying period is provided. Use of unauthorized replacement parts will void this warranty. Contact our Consumer Care Center at 1-800-682-7664 (North America), outside North America contact your local Philips Consumer Care Center. Internet information: www.sonicare.com (North America) or www.philips.com (outside North America).

Guarantee restrictions

- The terms of the guarantee do not cover the following:
- Brush heads
 - Damage caused by use of unauthorized replacement parts.
 - Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
 - Normal wear and tear, including chips, scratches, abrasions, discolouration or fading.