
FlexCare
900 series
Rechargeable
sonic toothbrush

PHILIPS
sonicare
the sonic toothbrush

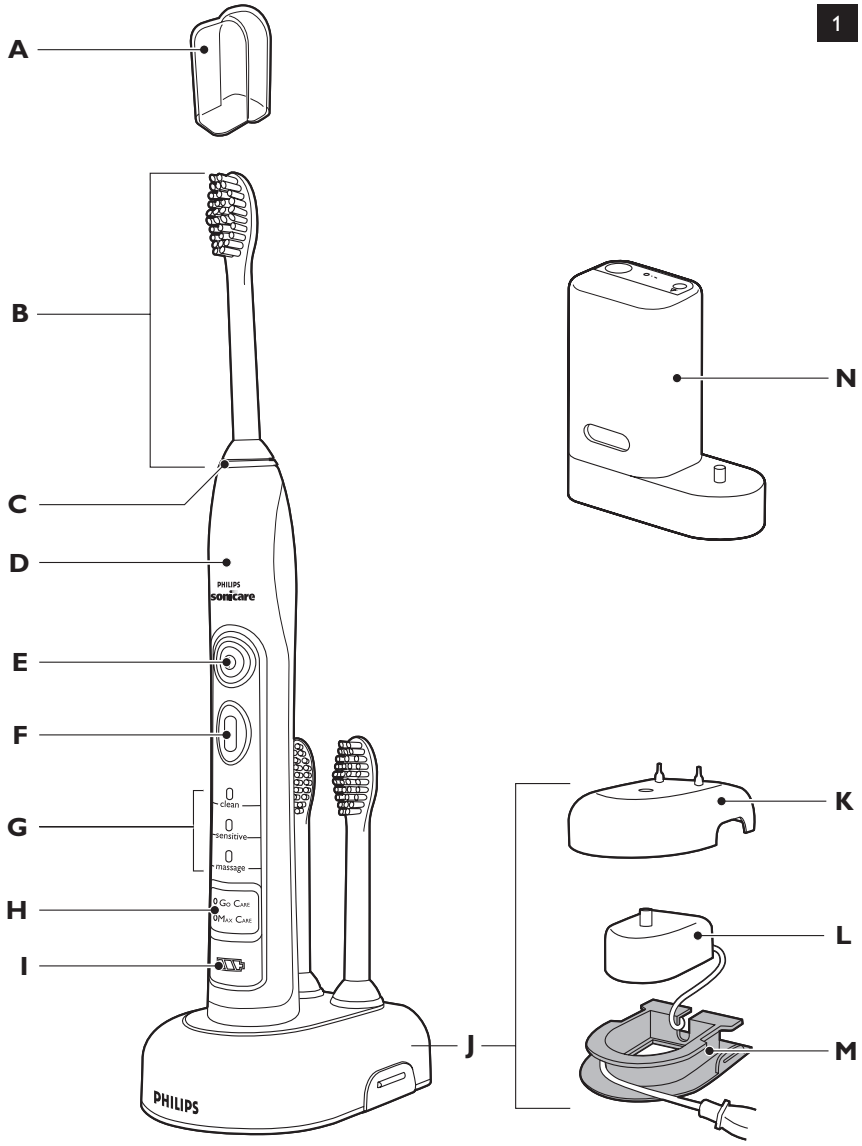HX6985, HX6983, HX6982, HX6942, HX6933, HX6932, HX6902,

ENGLISH	4
DEUTSCH	15
FRANÇAIS	27
ITALIANO	40
NEDERLANDS	51

Important

Read this user manual carefully before you use the appliance and save it for future reference.

Danger

- Keep the charger and/or sanitiser away from water. Do not place or store it over or near water contained in a bathtub, washbasin, sink etc. Do not immerse the charger and/or sanitiser in water or any other liquid. After cleaning, make sure the charger and/or sanitiser is completely dry before you connect it to the mains.

Warning

- Check if the voltage indicated on the bottom of the charger and/or sanitiser corresponds to the local mains voltage before you connect the appliance.
- The mains cord cannot be replaced. If the mains cord is damaged, discard the charger and/or sanitiser.
- Always have the charger and/or sanitiser replaced with one of the original type in order to avoid a hazard.
- Do not use the charger and/or sanitiser outdoors or near heated surfaces.
- If the appliance is damaged in any way (brush head, toothbrush handle, charger and/or sanitiser), stop using it. This appliance contains no serviceable parts. If the appliance is damaged, contact the Customer Care Centre in your country (see chapter 'Guarantee & service').
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- Discontinue use of the sanitiser if the UV light bulb remains on when the door is open. UV light can be harmful to the human eye and skin. Keep the sanitiser out of the reach of children at all times.

Caution

- Do not clean the brush head, the handle, the charger, the charger cover and/or the sanitiser in the dishwasher.
- If you have had oral or gum surgery in the previous 2 months, consult your dentist before you use the toothbrush.
- Consult your dentist if excessive bleeding occurs after using this toothbrush or if bleeding continues to occur after 1 week of use.
- The Sonicare toothbrush complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician or the device manufacturer of the implanted device prior to use.
- If you have medical concerns, consult your doctor before you use the Sonicare.
- This appliance has only been designed for cleaning teeth, gums and tongue. Do not use it for any other purpose. Stop using the appliance and contact your doctor if you experience any discomfort or pain.
- The Sonicare toothbrush is a personal care device and is not intended for use on multiple patients in a dental practice or institution.

- Stop using a brush head with crushed or bent bristles. Replace the brush head every 3 months or sooner if signs of wear appear.
- Do not use other brush heads than the ones recommended by the manufacturer.
- If your toothpaste contains peroxide, baking soda or bicarbonate (common in whitening toothpastes), thoroughly clean the brush head with soap and water after each use. This prevents possible cracking of the plastic.
- The UV light bulb is hot during and immediately after the sanitising cycle. Do not touch the UV light bulb when it is hot.
- Do not operate the sanitiser without the protective screen in place to avoid contact with a hot bulb.

Electromagnetic fields (EMF)

This Philips appliance complies with all standards regarding electromagnetic fields (EMF). If handled properly and according to the instructions in this user manual, the appliance is safe to use based on scientific evidence available today.

The Sonicare (Fig. 1)


- A** Hygienic travel cap
 - B** Brush head
 - C** Interchangeable colour code ring
 - D** Handle with soft grip
 - E** Power on/off button
 - F** Personalised Brushing button
 - G** Brushing modes
 - H** Brushing routines
 - I** Deluxe recharge gauge
 - J** Deluxe charger (specific types only)
 - K** Charger cover with brush head holder (specific types only)
 - L** Travel charger
 - M** Charger base with cord wrap
 - N** UV sanitiser with integrated charger and cord wrap (specific types only)
- Not shown: UV light bulb
 - Not shown: Sanitiser drip tray
 - Not shown: Protective screen for UV light bulb

Getting started


Changing the colour code ring

Sonicare brush heads come with interchangeable colour code rings to identify your brush head.
To change the colour code ring:

6 ENGLISH

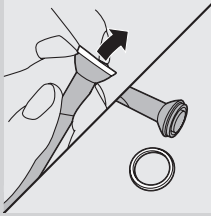
- 1 Pull the colour code ring from the bottom of the brush head.

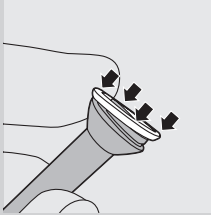- 2 Slip one edge of the new ring over the bottom of the brush head. Then press down on the other side to snap the ring in place.


Attaching the brush head

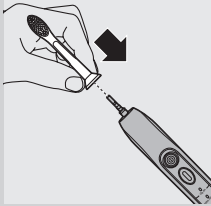- 1 Align the brush head so the bristles face the front of the handle.
- 2 Firmly press the brush head down on the metal shaft until it stops.

Note: There is a small gap between the colour code ring and the handle.

Charging your Sonicare

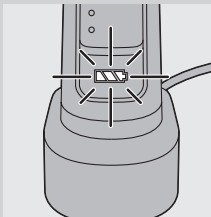- 1 Plug the charger or sanitiser into a wall socket.
- 2 Place the handle on the charger or sanitiser.
 - ▶ The flashing light on the battery gauge indicates that the toothbrush is charging.

Deluxe recharge gauge:

Indicates the amount of charge left in the battery




- 3 green LEDs: 75-100%
- 2 green LEDs: 50-74%
- 1 green LED: 25-49%
- 1 flashing yellow LED: less than 25%

Note: If the battery charge of your Sonicare is low, you hear 3 beeps and 1 LED on the recharge gauge flashes yellow for 30 seconds.

Note: To keep the battery fully charged, we recommend that you keep your Sonicare on the charger or sanitiser when not in use. It takes at least 24 hours to charge the battery fully.

Using the Sonicare

Brushing instructions

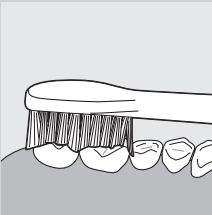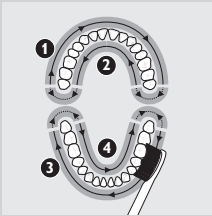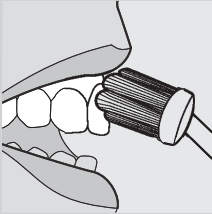- 1** Wet the bristles and apply a small amount of toothpaste.
- 2** Place the toothbrush bristles against the teeth at a slight angle towards the gumline.
- 3** Press the power on/off button to switch on the Sonicare.
- 4** Apply light pressure to maximise the Sonicare's effectiveness and let the Sonicare toothbrush do the brushing for you.
- 5** Gently move the brush head slowly across the teeth in a small back and forth motion so the longer bristles reach between your teeth. Continue this motion throughout your brushing cycle.

Note: To make sure you brush evenly throughout the mouth, divide your mouth into 4 sections using the Quadpacer (see chapter 'Features').

- 6** Begin brushing in section 1 (outside of upper teeth) and brush for 30 seconds before moving to section 2 (inside of upper teeth). Continue brushing in section 3 (outside of lower teeth) and brush for 30 seconds before moving to section 4 (inside of lower teeth).

- 7** After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs. You may also brush your tongue, with the toothbrush switched on or off, as you prefer.

Your Sonicare is safe to use on:

- Braces (brush heads wear out sooner when used on braces)
- Dental restorations (fillings, crowns, veneers)

Personalising your brushing experience

The Sonicare automatically starts in the default Clean mode. To personalise your brushing:

- 1** Before you switch on the Sonicare, press the Personalised Brushing button to toggle between modes and routines.
 - ▶ The green LED indicates the selected mode or routine.

Note: When the toothbrush is switched on, you can switch between modes but not between routines. Routines have to be selected before you switch on the appliance.

Brushing modes

Clean mode

Standard mode for superior teeth cleaning.

Sensitive mode

Gentle, yet thorough cleaning for sensitive gums and teeth.

Massage mode

Gentle gum stimulation.

Brushing routines

Go Care routine

1-minute brushing cycle in the Clean mode, for a quick clean. You hear the Quadpacer signal at 15-second intervals.

Max Care routine

3-minute brushing cycle that combines the Clean and Massage modes in one routine for a thorough mouth clean. There are 30 seconds of Clean mode and 15 seconds of Massage mode for each of the 4 sections of your mouth. You hear the Quadpacer signal at 45 second intervals.

Note: When the Sonicare is used in clinical studies, the default 2-minute Clean mode must be selected. The handle has to be fully charged. Deactivate the Easy-start feature. For areas where excess staining occurs, an additional 30 seconds of brushing time can be spent to help remove stains.

Features

Easy-start

- This Sonicare model comes with the Easy-start feature activated.
- The Easy-start feature gently increases power over the first 14 brushings to ease into the Sonicare experience.

Note: Each of the first 14 brushings must last at least 1 minute to properly move through the Easy-start ramp-up cycle.

Deactivating or activating the Easy-start feature

1 Attach your brush head to the handle.

2 Place the handle on the plugged-in charger or sanitiser.

- To deactivate Easy-start:

Press and hold the on/off button for 5 seconds. You hear 1 beep to indicate that the Easy-start feature has been deactivated.

- To activate Easy-start:


Press and hold the on/off button for 5 seconds. You hear 2 beeps to indicate that the Easy-start feature has been activated.

Note: Using the Easy-start feature beyond the initial ramp-up period is not recommended and reduces the Sonicare's effectiveness in removing plaque.

Smartimer

The Smartimer indicates that your brushing cycle is complete by automatically switching off the toothbrush at the end of the brushing cycle. Dental professionals recommend brushing no less than 2 minutes twice a day.

Quadpacer

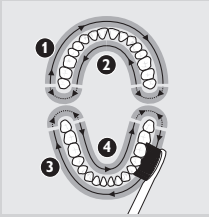- The Quadpacer is an interval timer that has a short beep and pause to remind you to brush the 4 sections of your mouth. Depending on the brushing mode or brushing routine you have selected, the Quadpacer beeps at different intervals during the brushing cycle. See the 'Brushing instructions' and 'Personalising your brushing experience' sections.


The Quadpacer on this model has been activated. To deactivate or reactivate the Quadpacer:

- 1** Place the handle with a brush head attached on the plugged-in charger or sanitiser
- 2** To deactivate the Quadpacer: press and hold the Personalised Brushing button for 5 seconds. You hear 1 beep to indicate that the Quadpacer has been deactivated.
- 3** To reactivate the Quadpacer: press and hold the Personalised Brushing button for 5 seconds until you hear 2 beeps. This indicates that the Quadpacer is active again.


Note: The Quadpacer cannot be deactivated on the pre-programmed Go Care and Max Care brushing routines.

Cord wrap on deluxe charger (specific types only)


If your model includes the deluxe charger cover and base, the travel charger is pre-installed inside the cover. If you want to shorten the mains cord, you can store excess cord in the cord-wrap feature built into the charger base.

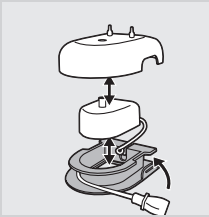- 1** To separate the charger cover from the charger base, press the two grey snap bars on the charger base and pull the white charger cover upwards.

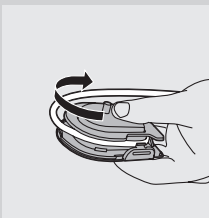
- 2** Wrap excess cord around the grey charger base as shown in the picture. Be sure to wrap the cord on the inside of the two snap levers.

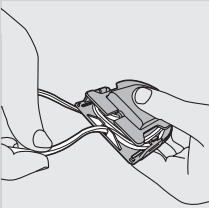
- 3** When you have stored the excess cord, guide the cord through the small groove in the rear of the grey charger base.
- 4** Reattach the charger cover by pressing it down over the charger base until it snaps into place.

Tip: For extra convenience during travel, you can remove the travel charger and use it without the charger cover and charger base.

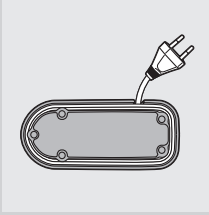Cord wrap on sanitiser (specific types only)


- If your model includes the sanitiser, you can store excess cord with the cord-wrap feature built into the bottom of the sanitiser.

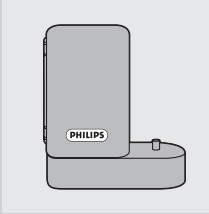Sanitising (specific types only)

- With the sanitiser, you can clean your brush head after every use.

Discontinue use of the sanitiser if the UV light bulb remains on when the door is open or if the Philips logo is broken or missing from the sanitiser. UV light can be harmful to the human eye and skin. Keep the sanitiser out of the reach of children at all times.


Unplug the sanitiser and call Customer Service if the sanitiser gives off smoke or a burning smell while it is operating.

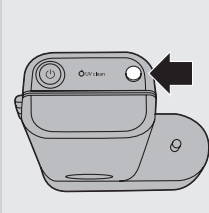- 1 After brushing, rinse the brush head and shake off excess water.

Do not place the travel cap on the brush head during sanitisation.

- 2 To open the sanitiser door, press the door release button.

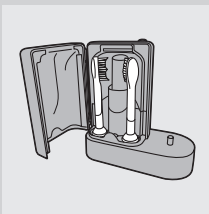- 3 Place the brush head on one of the 2 pegs in the sanitiser.
 - Make sure the bristles of the brush head directly face the light bulb.

Note: Only clean Sonicare ProResults brush heads in the sanitiser.

- 4 Make sure the sanitiser is plugged into wall socket.

- 5 Close the door and press the green power on/off button once to select the UV clean cycle.

Note: You can only switch on the sanitiser when the door is properly closed.

Note: If you open the door during the sanitising cycle, the sanitiser stops.


Note: The sanitiser cycle runs for 10 minutes and then automatically switches off.

- ▶ The sanitiser is in operation when the blue light glows through the Philips logo and the UV clean LED flashes slowly.
- ▶ When the sanitising cycle is complete, the UV clean LED lights up green continuously and the sanitiser automatically switches off.

Cleaning

Do not clean brush heads, the handle, the travel charger, the charger cover and the UV sanitiser in the dishwasher.

Toothbrush handle

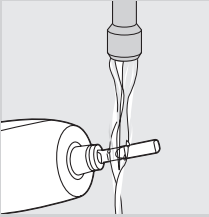
- 1 Remove the brush head and rinse the metal shaft with warm water.

Do not push on the rubber seal round the metal shaft with sharp objects, as this may cause damage.

- 2 Use a damp cloth to wipe the entire surface of the handle.

Brush head

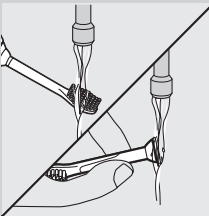
- 1 Rinse the brush head and bristles after each use.
- 2 Remove the brush head from the handle and rinse the brush head connection at least once a week with warm water.

Travel charger and deluxe charger (specific types only)


- 1 Unplug the charger.
- 2 Specific types only: Disassemble the deluxe charger and take out the travel charger.
- 3 Use a damp cloth to wipe the surface of the charger.

Sanitiser (specific types only)

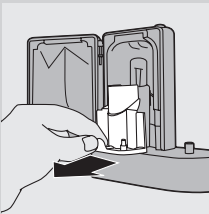Do not clean the sanitiser when the UV light bulb is hot.
For optimal effectiveness, it is recommended to clean the sanitiser weekly.

- 1 Unplug the sanitiser.
- 2 Slightly lift the drip tray and pull it out. Rinse the drip tray and wipe it clean with a damp cloth.
- 3 Clean all reflector surfaces with a damp cloth.

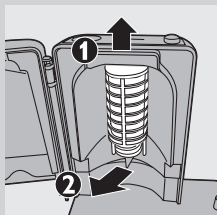**4 Remove the protective screen in front of the UV light bulb.**

To remove the screen, gently lift it up (1) and pull it out (2).

5 Remove the UV light bulb.

To remove the light bulb, grasp it and pull it out of the metal clasp.

6 Clean the protective screen and the UV light bulb with a damp cloth.**7 Reinsert the UV light bulb.**

To reinsert the light bulb, align the bottom of the light bulb with the metal clasp and push the bulb into the clasp.

8 Reinsert the protective screen.

To reinsert the screen, align the pegs on the screen with the slots on the reflective surface near the UV light bulb. Then insert the pegs into the slots and slide down the screen to secure it in the sanitiser.

Storage

- If you are not going to use the appliance for a long time, remove the mains plug from the wall socket. Then clean the appliance and store it in a cool and dry place away from direct sunlight. You can use the cord wrap to store the mains cord neatly.
- You can store brush heads on the pegs at the back of the deluxe charger (specific types only).

Replacement

Brush head

- Replace Sonicare brush heads every 3 months to achieve optimal results.
- Use only Sonicare ProResults replacement brush heads.

UV light bulb


- You can order replacement UV light bulbs from the Customer Care Centre in your country or an authorised Philips service centre.

Disposal

- Do not throw away the appliance with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment.
- The lamp of the sanitiser contains mercury, which is a substance that can be harmful to the environment. When you discard the lamp, do not throw it away with the normal household waste but hand it in at an official collection point.
- The built-in rechargeable battery contains substances that may pollute the environment. Always remove the battery before you discard and hand in the appliance at an official collection point. Dispose of the battery at an official collection point for batteries. If you have trouble removing the battery, you can also take the appliance to a Philips service centre. The staff of this centre will remove the battery for you and will dispose of it in an environmentally safe way.


Removing the rechargeable battery

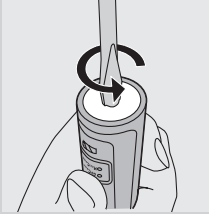
Please note that this process is not reversible.

1 To deplete the rechargeable battery of any charge, remove the handle from the charger or sanitiser, switch on the Sonicare and let it run until it stops.

Repeat this step until you can no longer switch on the Sonicare.

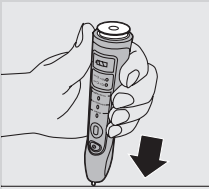2 Insert a flat-head (standard) screwdriver into the slot located at the bottom of the handle. Turn the screwdriver anticlockwise to release the bottom cap.

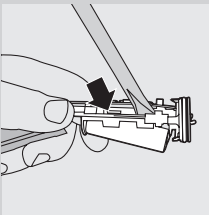3 Hold the handle upside down and push down on the shaft to release the internal components of the handle.

4 Insert the screwdriver under the circuit board, next to the battery connections, and twist to break the connections. Remove the circuit board and pry the battery from the plastic carrier.

Do not connect the appliance to the mains again after you have removed the battery.

Guarantee & service

If you need information or if you have a problem, please visit the Philips website at www.philips.com or contact the Philips Customer Care Centre in your country (you find its phone number in the worldwide guarantee leaflet). If there is no Customer Care Centre in your country, go to your local Philips dealer or contact the Service Department of Philips Domestic Appliances and Personal Care BV.

Guarantee restrictions

The terms of the international guarantee do not cover the following:

- Brush heads
- Damage caused by misuse, abuse, neglect or alterations
- Normal wear and tear, including chips, scratches, abrasions, discolouration or fading
- UV light bulb

Frequently asked questions

This chapter lists the questions most frequently asked about the appliance. If you cannot find the answer to your question, contact the Customer Care Centre in your country.

14 ENGLISH

Question	Answer
Why doesn't the Sonicare toothbrush work?	You may need to recharge the Sonicare. Charge the Sonicare for at least 24 hours. The socket may not be live. The electricity supply to sockets in mirror cabinets over washbasins or shaver sockets in bathrooms may be cut off when the light is switched off. You may need to reattach or replace the brush head.
Why do I feel a tickling sensation when I use the Sonicare?	You may experience a slight tickling/tingling sensation when you use Sonicare for the first time. As you become accustomed to brushing with the Sonicare, this sensation diminishes.
The brush head chatters against my teeth. What can I do about this?	To avoid this, keep the brush head properly oriented, with the bristles placed on the gumline at a slight angle.
Why does the toothbrush seem less powerful?	You may need to charge the Sonicare. The Easy-start feature may be on. Deactivate the Easy-start feature (See chapter 'Features'). You may need to replace the brush head.
Why doesn't the Personalised Brushing button work during brushing?	You may have selected a brushing routine. If this is the case, the Personalised Brushing button is not active (See chapter 'Using the Sonicare').
Why doesn't the sanitiser light go on?	The socket may not be live. The electricity supply to sockets in mirror cabinets over washbasins or shaver sockets in bathrooms may be cut off when the light is switched off. You may need to shut the sanitiser door. You may need to replace the UV light bulb (see chapter 'Replacement').
Why does the toothbrush seem to lose power during brushing?	You may be in the Max Care brushing routine. If this is the case, the routine alternates between the Clean and Massage modes throughout the 3-minute brushing cycle. The Clean and Massage modes have a different sound and brush movement.
Do I need a special charger if I travel with my Sonicare?	Travel tip: the travel charger, the deluxe charger and the sanitiser can be used at common worldwide voltages between 100-240V AC, at either 50 or 60Hz. All you need is a simple plug adapter.

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© 2008 Koninklijke Philips Electronics N.V. (KPEENV). All rights reserved.
PHILIPS and the Philips shield are registered trademarks of KPEENV.
Quadpacer, Sonicare and the Sonicare logo are trademarks of Philips Oral
Healthcare, Inc. and/or KPEENV
Printed in the U.S.A.