

Quick start guide

Review **IMPORTANT SAFEGUARD** section of your user manual prior to the first use.

Contents: 2 Sonicare HealthyWhite handles
2 DiamondClean brush heads (Standard)
1 UV Sanitizer with integrated charger, cord wrap and brush head storage
1 Multi-voltage, compact travel charger
2 Travel cases



Before you start



Charge for 2 hours minimum

1 HealthyWhite brushing modes and brush heads:

Your Sonicare HealthyWhite is set to Clean mode. Push the on/off button to toggle between modes. A green LED indicates the selected mode.

Modes:


Clean mode: Standard mode for complete cleaning.

Sensitive mode: Gentle, yet thorough cleaning for sensitive gums and teeth.

Clean and White mode: 2 minutes of Clean mode with an additional 30 seconds of White mode to focus on your visible front teeth.


Brush heads:

DiamondClean brush head: Diamond shaped bristles provide exceptional cleaning and whiter teeth.

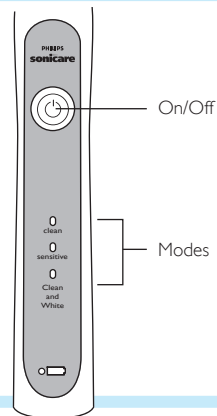



Questions? Call Philips Consumer Care at 1-800-682-7664

2 Brushing instructions:

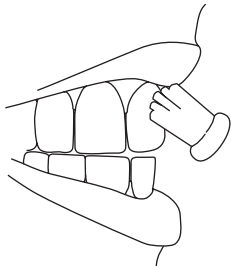

Before turning your Sonicare on, wet the bristles and apply a small amount of toothpaste. Place bristles against your teeth at a slight angle towards the gumline. Press the on/off button. Hold the handle with a light grip and gently move the brush head slowly across the teeth in a small back and forth motion so the longer bristles reach between your teeth. Continue this motion throughout your brushing cycle. You may experience a slight tickling/tingling sensation when using your Sonicare toothbrush for the first time. As you become accustomed to brushing with Sonicare, this sensation will diminish.

3 UV Sanitizer:

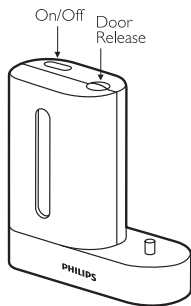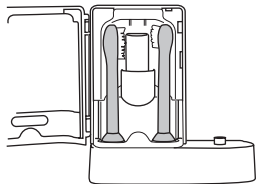

The UV sanitizer uses a specially designed UV light bulb to eliminate germs on your brush head.

After brushing, rinse the brush head and shake off excess water. Place the brush head(s) on the peg(s) located inside the UV sanitizer. Make sure bristles face the UV light bulb. Do not place the travel cap on the brush head during sanitization.

Close the door and press the green power on/off button once to select the UV clean cycle.

The sanitizer is in operation when the blue light glows through the window. When the sanitizing cycle is complete, the sanitizer automatically shuts off.