
UV Sanitiser

HX7990

PHILIPS
sonicare
the sonic toothbrush

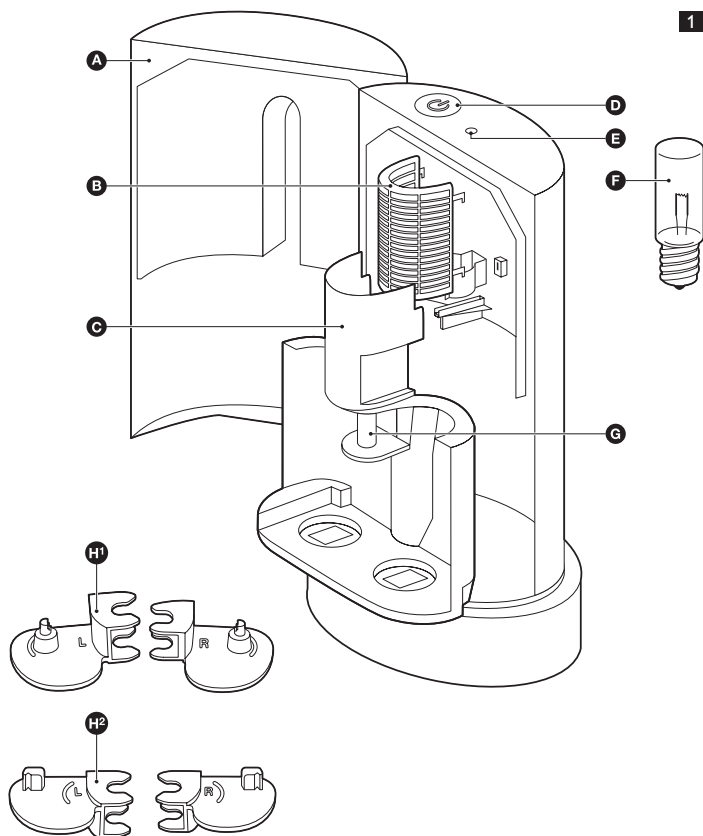HX7990

ENGLISH 4
DEUTSCH 15
FRANÇAIS 26
ITALIANO 37
NEDERLANDS 48

Introduction

Congratulations on your purchase and welcome to Philips! To fully benefit from the support that Philips offers, register your product at www.philips.com/welcome.

The Sonicare UV Sanitiser uses a specially designed UV light bulb that helps kill germs on the brush head of your rechargeable toothbrush in a 10-minute sanitising cycle. When used properly on Sonicare ProResults, Sonicare e-series, Oral-B Flexisoft* and Oral-B FlossAction** brush heads, this appliance is designed and demonstrated to act upon the following germs: Streptococcus mutans, Herpes simplex virus and Escherichia coli.

*Flexisoft is a registered trademark of Braun Aktiengesellschaft Joint Stock Company Fed Rep Germany.

**Oral-B and FlossAction are registered trademarks of Gillette Canada Company, dba Oral-B Laboratories Corporation Canada.

Important

Read this user manual carefully before you use the appliance and save it for future reference.

Danger

- Keep the sanitiser away from water. Do not place or store the sanitiser over or near water contained in a bathtub, washbasin, sink etc. Do not immerse the sanitiser in water or any other liquid. After cleaning, make sure the sanitiser is completely dry before you connect it to the mains.

Warning

- Check if the voltage indicated on the bottom of the sanitiser corresponds to the local mains voltage before you connect the appliance.
- The mains cord of this appliance cannot be replaced. If the mains cord is damaged, the appliance must be discarded.
- If the appliance is damaged in any way, stop using it. This appliance contains no serviceable parts. If the appliance is damaged, contact the Customer Care Centre in your country (see chapter 'Guarantee & service').
- Keep the sanitiser away from heated surfaces. Do not use the sanitiser outdoors or in areas with direct sunlight.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- Discontinue use of the sanitiser if the UV light bulb remains on when the door is open. UV light can be harmful to the human eye and skin. Keep the sanitiser out of the reach of children at all times.

Caution

- Do not use the sanitiser or any of its parts for purposes other than described in this manual.
- Do not clean the sanitiser in the dishwasher.
- The UV light bulb is hot during and immediately after the sanitising cycle. Do not touch the UV light bulb when it is hot.

6 ENGLISH

- Do not operate the sanitiser without the protective screen in place to avoid contact with a hot bulb.
- If you have used the sanitiser for 3 UV clean cycles in a row, switch off the appliance for at least 30 minutes before you start another UV clean cycle.

Electromagnetic fields (EMF)


This Philips appliance complies with all standards regarding electromagnetic fields (EMF). If handled properly and according to the instructions in this user manual, the appliance is safe to use based on scientific evidence available today.

General description (Fig. 1)


- A** Door recess
 - B** Removable protective screen for light bulb
 - C** Removable drip tray
 - D** On/off button
 - E** UV clean LED
 - F** UV light bulb
 - G** Centre post to attach brush head holders
 - H** Brush head holders
 - H1 brush head holder for Sonicare ProResults brush heads
 - H2 brush head holder for Oral-B Flexisoft*, Oral-B FlossAction**
- No brush head holder needed for Sonicare e-series brush heads.

Before first use


Note: Your sanitiser may come with the UV light bulb installed. If this is the case, go on to chapter 'Preparing for use'.

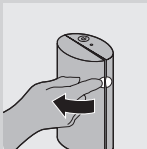
- 1** Place your finger in the door recess on the right-hand side of the sanitiser and pull open the door.

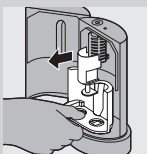- 2** Remove the drip tray.

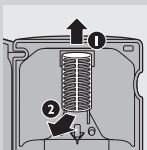- 3** To remove the protective screen, gently lift it up (1) and pull it forward (2).

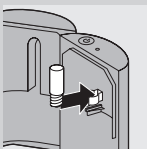- 4** Snap the UV light bulb into the lamp holder.

Do not crush or break the bulb. UV light bulbs contain a small amount of mercury.


Note: Use only the UV light bulb supplied or a replacement UV light bulb ordered from Philips or service centres authorised by Philips.

Note: To remove the UV light bulb, pull it out of the lamp holder.

- 5** Align the pegs on the protective screen with the slots in the reflective surface of the sanitiser on either side of the UV light bulb.


- 6** Insert and gently push down the screen to secure it in the sanitiser.

8 ENGLISH

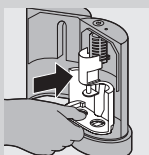
- 7 Reattach the drip tray.

Preparing for use

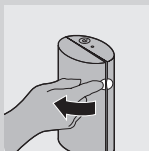Sanitising Sonicare e-series brush heads


- 1 Place your finger in the door recess on the right-hand side of the sanitiser and pull open the door.

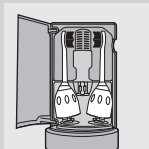
- 2 Place the brush heads directly on the drip tray with the bristles facing the UV light bulb.

Sanitising Sonicare ProResults, Oral-B Flexisoft and Oral-B FlossAction brush heads


Note: This sanitiser is NOT suitable for other brush heads.

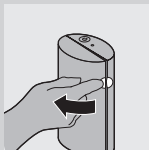- 1 Place your finger in the door recess on the right-hand side of the sanitiser and pull open the door.

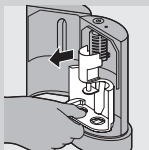
- 2 Remove the drip tray.
- 3 Select the appropriate brush head holder (see chapter 'General description').
- 4 Attach the brush head onto the appropriate brush head holder (see chapter 'General description').

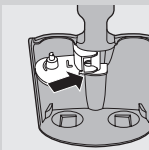
- 5** Attach the left brush head holder (L) to the left side of the centre post until it clicks into place.

Note: The brush head holders do not fit properly when you place the right brush head holder first. Always start with the left brush head holder.

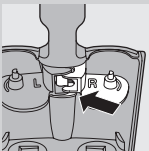- 6** Attach the right brush head holder (R) to the right side of the centre post until it clicks into place.


- 7** Place the drip tray back into the sanitiser.

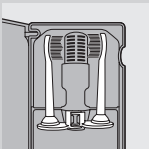
- 8** Place the brush heads onto the brush head holders with the bristles pointing towards the UV light bulb.

Sanitising brush heads of two different types

If you want to sanitise brush heads of two different types, use the appropriate brush head holder for each type.

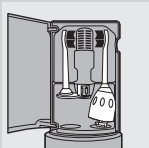
- 1** Select the appropriate brush head holder for each type of brush head (see chapter 'General description'. Sonicare e-series brush heads do not require a brush head holder).

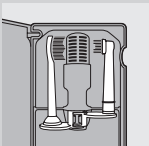- 2** Follow the instructions given under 'Sanitising Sonicare ProResults, Oral-B Flexisoft and Oral-B FlossAction brush heads'.

Using the appliance

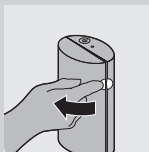
Do not place the travel cap on the brush head during sanitisation.

- 1** After brushing, rinse the brush head and shake off excess water.
- 2** Place your finger in the door recess on the right-hand side of the sanitiser and pull open the door.
- 3** Put the brush heads into the sanitiser (see chapter 'Preparing for use').

Note: Make sure the bristles of the brush head face the light bulb.

- 4** Close the sanitiser door.

Note: The sanitiser does not work when the door is not properly closed.

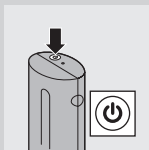- 5** Put the mains plug in a wall socket.
- 6** Press the on/off button to start the UV clean cycle.

Note: The sanitiser is in operation when the blue light glows through the window on the front of the sanitiser and the UV clean LED flashes slowly.

Note: The sanitiser needs to run for the full 10-minute UV clean cycle to get optimal results.

Note: When the sanitising cycle is complete, the UV clean LED turns solid green and the sanitiser automatically switches off.

Note: You can stop the sanitiser during the sanitising cycle by pressing the on/off button or by opening the door.

Cleaning and maintenance

Do not clean any of the sanitiser accessories in a dishwasher.


Do not clean the sanitiser when the UV light bulb is hot.

For optimal effectiveness, it is recommended to clean your sanitiser as follows:

- 1** Unplug the sanitiser.
- 2** Remove the drip tray, protective screen and the light bulb (see chapter 'Before first use').
- 3** Thoroughly rinse the brush head holders and the drip tray with warm soapy water and wipe them clean with a damp cloth.
- 4** Wipe the bulb, the protective screen and all internal and external surfaces with a damp cloth.
- 5** Let all items dry before you reinstall them.
- 6** Reattach the drip tray, protective screen and light bulb (see chapter 'Before first use').

Storage

- If you are not going to use the sanitiser for an extended period of time, unplug it, clean it and store it in a cool and dry place away from direct sunlight.
- You can use the cord-wrap feature to store the mains cord in an organised manner.

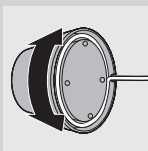Replacement

Replacement parts are available through your Philips dealer or the Philips Customer Care Centre in your country.

The parts available for ordering are:

- UV light bulb
- Brush head holders
- Protective screen

Replacing the UV light bulb

When the UV clean LED flashes quickly, the UV light bulb must be replaced.

See chapter 'Before first use' for replacing the UV light bulb.

Environment

- The UV light bulb contains mercury and may be harmful to the environment. Do not crush, break or attempt to disassemble the bulb. Do not throw away the UV light bulb with the normal household waste when you replace it, but hand it in at an official collection point for recycling.
- Do not throw away the appliance with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment.

Guarantee and service

If you need information or if you have a problem, please visit the Philips website at www.philips.com or contact the Philips Customer Care Centre in your country (you find its phone number in the worldwide guarantee leaflet). If there is no Customer Care Centre in your country, go to your local Philips dealer or contact the Service

Department of Philips Domestic Appliances and Personal Care BV.

Guarantee restrictions

The terms of the international guarantee do not cover the following:

- UV light bulb
- Damage to brush heads
- Damage caused by misuse, abuse, neglect, alterations or unauthorised repair.
- Normal wear and tear, including chips, scratches, abrasions, discolouration or fading.

Troubleshooting

This chapter summarises the most common problems you may encounter with your sanitiser. If you are unable to solve the problem, contact the Customer Care Centre in your country.

Problem	Solution
The UV clean LED flashes quickly.	The UV light bulb must be replaced (see chapter 'Replacement').
The UV light does not go on.	Make sure the door is properly closed.
	Make sure the sanitiser is plugged into a live wall socket.
	The UV light bulb must be replaced (see chapter 'Replacement').
The brush head of my toothbrush does not fit in the sanitiser.	Check the list of brush heads that fit in the sanitiser (see chapter 'Introduction').
	Make sure you use the correct brush head holder (see chapter 'General description').

14 ENGLISH

Problem	Solution
The bristles of my brush head change colour.	Make sure the bristles point towards the UV light bulb. Exposure to UV light may cause bristles to change colour slightly over time.

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© 2007 Philips Oral Healthcare, Inc. All rights reserved.
PHILIPS and the Philips shield are registered trademarks of
Sonicare logo are trademarks of Philips Oral Healthcare, Inc.
Printed in the U.S.A.