Register your product and get support at
www.philips.com/welcome

**Philips Sonicare
AirFloss
100 series**

PHILIPS

HX8181, HX8141, HX8111

ENGLISH 4

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1	Important	5
2	What's in the box	7
3	Before first use	8
4	Preparing for use.....	9
5	Using	10
6	Auto off	11
7	Cleaning	12
8	Storing your AirFloss	13
9	Environment	14
10	Guarantee and service	15

1

IMPORTANT

Read this user manual carefully before you use the appliance and save it for future reference.

DANGER

- Keep the charger away from water. Do not place or store it over or near water contained in a bathtub, washbasin, etc. Do not immerse the charger in water or any other liquid. After cleaning, make sure the charger is completely dry before you connect it to the mains.

WARNING

- The mains cord cannot be replaced. If the mains cord is damaged, discard the charger. Always have the charger replaced with one of the original type in order to avoid a hazard.
- If the appliance is damaged in any way (nozzle, handle and/or charger), stop using it.
- This appliance contains no user-serviceable parts. If the appliance is damaged, contact the Consumer Care Centre in your country (see chapter 'Guarantee and service').
- Do not use the charger outdoors or near heated surfaces.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.


CAUTION

- If you have had oral or gum surgery in the previous 2 months, consult your dentist before you use this appliance.
- Consult your dentist if excessive bleeding occurs after using this appliance or if bleeding continues to occur after 1 week of use.
- If you have medical concerns, consult your doctor before you use this appliance.
- This Philips appliance complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician or the appliance manufacturer of the implanted appliance prior to use.
- Do not use attachments other than the ones recommended by the manufacturer.
- Do not clean the nozzle, the handle, the charger, the water reservoir and its cover in the dishwasher or a microwave.

ELECTROMAGNETIC FIELDS (EMF)

This Philips appliance complies with all standards regarding electromagnetic fields (EMF). If handled properly and according to the instructions in this user manual, the appliance is safe to use based on scientific evidence available today.

2

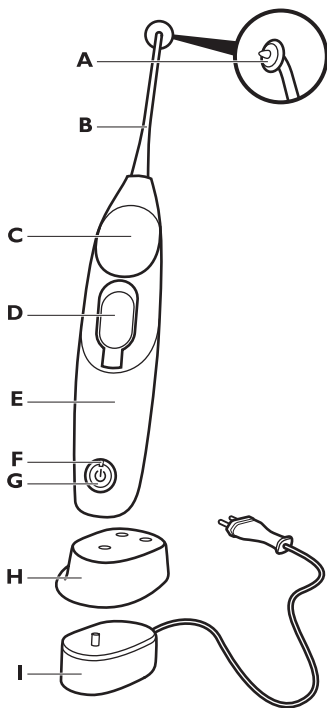A Guidance tip

B Nozzle of AirFloss

C Activation button

D Water reservoir and cover

E Handle

F Charge indicator

G Power on/off button





H Removable charger cover
(specific types only)

I Charger

3

BEFORE FIRST USE







Charging - Charge the AirFloss for 24 hours before first use. The charge indicator flashes green until fully charged. Rapid yellow flashing indicates the battery is low (fewer than 3 uses left). A full charge lasts 2 weeks (or 14 uses).

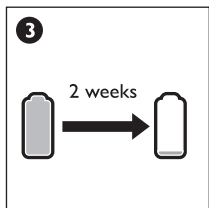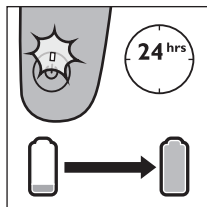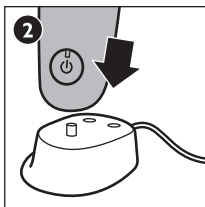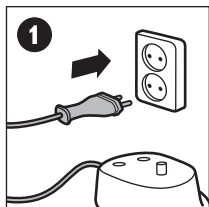4

PREPARING FOR USE






Preparing for use - Attach the nozzle, fill the reservoir with water or mouthwash. After filling the reservoir, press the power on/off button to switch on. Then press the activation button approximately 6 times or until spray is produced. The reservoir holds enough liquid for approximately 2 uses. Replace the nozzle every 6 months for optimal results.

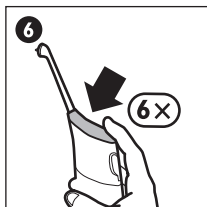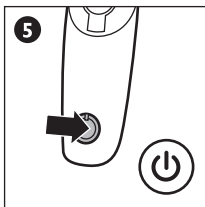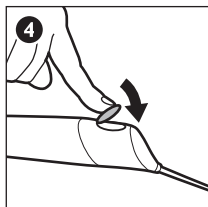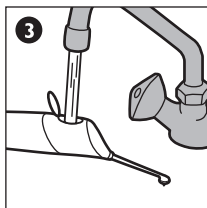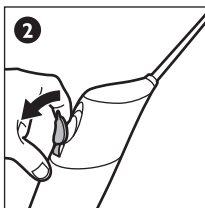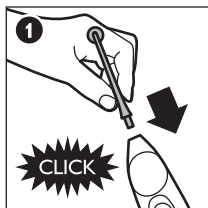
5

USING

Using the AirFloss - The AirFloss is intended to be used on the outside of the teeth. Place the tip between teeth and apply gentle pressure. Press the activation button to deliver a burst of air and water. Glide the guidance tip along the gum line until you feel it settle between the next two teeth. Press the activation button again and repeat for all spaces, including behind your back teeth. Press the power on/off button to switch off after use.

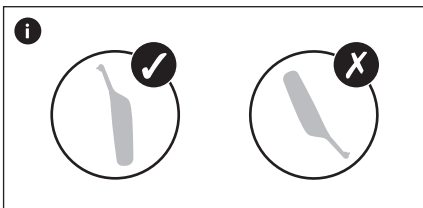
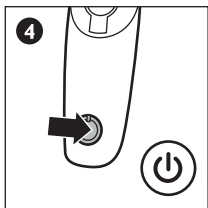
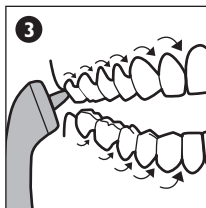
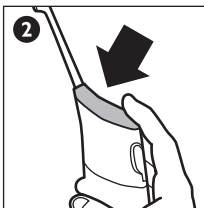
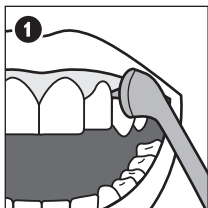
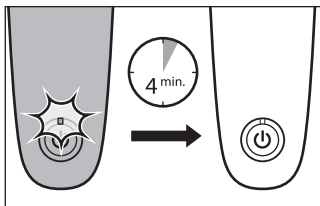
6**AUTO OFF**

7

CLEANING





Cleaning - Remove the nozzle and rinse with water. Use a cotton bud to wipe the reservoir, if desired. Make sure no foreign objects get inside. Do not use cleaning agents to clean the appliance. Do not clean in the dishwasher or any other sanitising device.

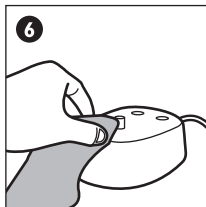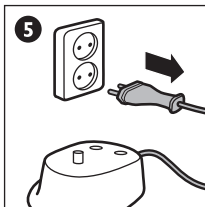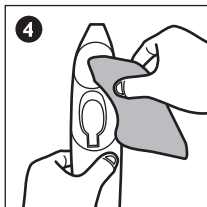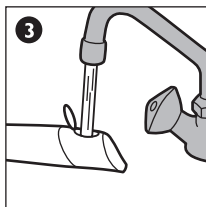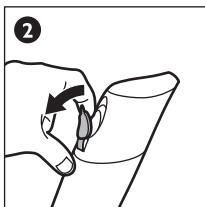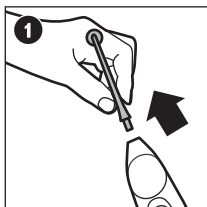8

STORING YOUR AIRFLOSS

Storing - If the AirFloss is not used for a longer period of time, empty the reservoir, press the activation button until no spray is produced, unplug the charger, clean the appliance and store it in a cool, dry place away from direct sunlight.




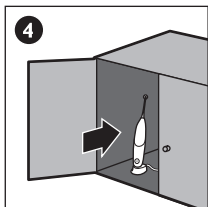
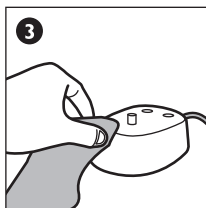
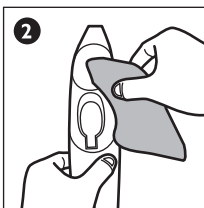
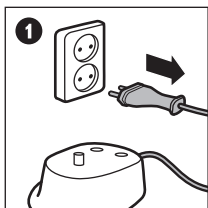
9

ENVIRONMENT

Do not throw away the appliance with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment.

The built-in rechargeable battery contains substances that may pollute the environment. Always remove the battery before you discard and hand in the appliance at an official collection point. Dispose of the battery at an official collection point for batteries. If you have trouble removing the battery, you can also take the appliance to a Philips service centre. The staff of this centre will remove the battery for you and will dispose of it in an environmentally safe way.

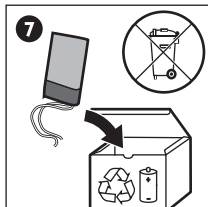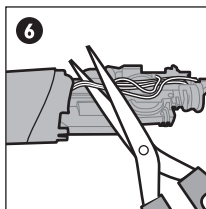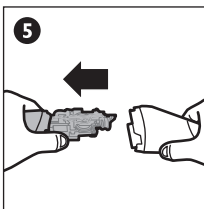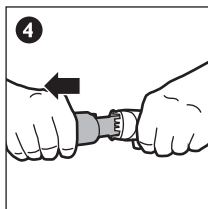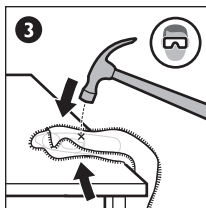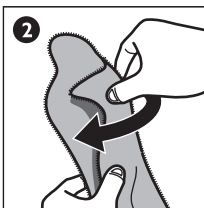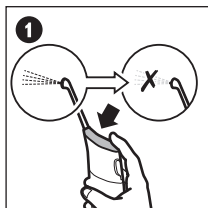REMOVING THE RECHARGEABLE BATTERY

Please note that this process is not reversible. Press the activation button repeatedly until the AirFloss no longer produces any bursts of air. Use a screwdriver to pry open the handle and remove the battery. Observe basic safety precautions.

If you need service or information or if you have a problem, please visit the Philips website at www.philips.com/sonicare or contact the Philips Consumer Care Centre in your country. You find its phone number in the worldwide guarantee leaflet. If there is no Consumer Care Centre in your country, go to your local Philips dealer.

GUARANTEE RESTRICTIONS

The terms of the international guarantee do not cover the following:

- Nozzle of AirFloss
- Damage caused by misuse, abuse, neglect, alterations or unauthorised repair
- Normal wear and tear, including chips, scratches, abrasions, discolouration or fading

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©2012 Koninklijke Philips Electronics NV (KPENV). All rights reserved.
PHILIPS and the Philips shield are registered trademarks of KPENV.
AirFloss, Sonicare and the Sonicare logo are trademarks of Philips Oral
Healthcare, Inc. and/or KPENV.