PHILIPS

sensiflex

ENGLISH 4

FRANÇAIS 12

DEUTSCH 21


NEDERLANDS 29

ITALIANO 37

ESPAÑOL 45

PORTUGUÈS 54


First use

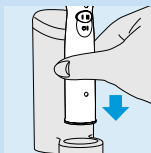1 Put the plug in the wall socket.

2 Place the toothbrush in the charging unit. The brush will now be charged automatically.

Let the appliance charge for at least 16 hours before you use it for the first time.

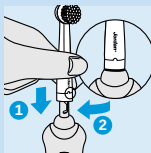3 Mount the brush head on the toothbrush by inserting it into the bayonet catch and turning it.


► Put some toothpaste on the brush.

4 Place the brush against your teeth and press the on/off button.

5 Brush one tooth at a time, until you have cleaned all teeth completely (minimum 2 minutes, see also the brushing instructions).

► As a result of the Gum Protection System, you will hear and feel the brush flexing backwards when too much pressure is applied.

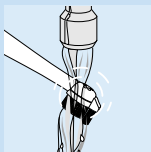6 After use, always clean the inside and the outside of the brush head as well as the toothbrush itself.

► You can put the toothbrush back in the charging unit after use.



For more information on how to use the appliance, read the section 'instructions for use'.

Important

- ▶ Be sure to check if the wall socket to which you connect the charging unit is live. In some bathroom cabinets, the electricity supply to the built-in socket may be interrupted as soon as the light is switched off.
- ▶ Charge the appliance for at least 16 hours before you use it for the first time.
- ▶ The mains cord of the charging unit cannot be replaced. If the cord is damaged, the unit should be discarded. Please apply to your Philips dealer or to the Philips organisation in your country (refer to the World Wide Guarantee leaflet).
- ▶ The built-in battery contains substances that can be harmful to the environment. If you discard the appliance in due course, remove the built-in battery and hand it in at an official collection point or dispose of it in accordance with the regulations. See section 'Environmental aspects'.
- ▶ *If your appliance is marked with  the following text applies:*
The built-in rechargeable (NiCd) battery contains substances which may pollute the environment.
When you discard your appliance in due course, please ensure that the battery is kept separate from normal household refuse and disposed of at an officially assigned collection point. Please refer to section 'Environmental aspects'.


Instructions for use

Mounting the appliance on the wall

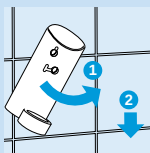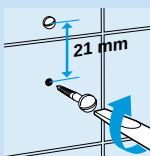- 1 Mount the charging unit to the wall.
 - Do not allow the charging unit to come into contact with water. The toothbrush, however, is watertight and can therefore be used safely anywhere.
- 2 To mount the charger to the wall, you will need 2 wood screws (3,5 x 25 mm) and two wall plugs (S5). They are available in any home centre.

Charging

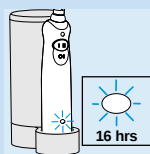
- 1 Leave the appliance to charge for at least 16 hours before you use it for the first time.
 - The toothbrush can be put back in the charging unit after use.
 - Only type HX2565: If the battery in the toothbrush is being charged, the pilot light is green and blinks. If the battery has been fully charged, the pilot light burns continuously.
 - If the battery is almost empty, the pilot light is red and blinks for 10 seconds after the toothbrush has been switched off.


Tip: The toothbrush can be left to charge continuously. When you start using the appliance, we recommend you to use it until the battery is empty and do so three times in a row. You can do this by using the appliance in the normal way without recharging it between brushing session until the motor almost stops running. To extend

the life time of the battery, repeat this procedure every six months.


How to use

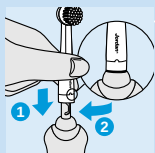
- 1 Mount the brush head onto the toothbrush by inserting it into the bayonet catch and turning it. This will make it impossible for the brush head to become detached due to vibration.

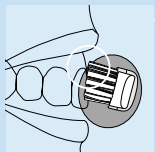- 2 Place the brush in your mouth.

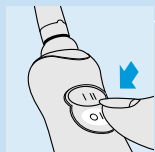- 3 Press the on/off button to switch the toothbrush on.

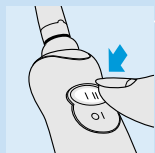
- 4 The toothbrush will start rotating at the highest speed. To switch to a lower speed, press the 2-speeds button.

Brushing



- 1 Place the brush against the tooth at an angle, in such a way that it covers the surface of the tooth as well as the gum line. Hold the brush in that position for a few seconds and then move on to the next tooth.

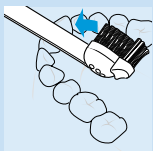- ▶ The round brush head removes plaque from the teeth and beneath the gum margin.
- ▶ The individual moving Active Tip cleans between the teeth.

- 2 The brush will do the work for you. All you have to do is guide the brush carefully along the teeth and molars.
- ▶ Make sure that you reach all surfaces of the teeth (inside, outside and chewing surface) as well as the gums properly. A fixed brushing pattern will prevent you from skipping certain areas.

Safe brushing

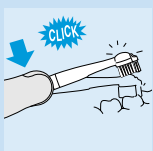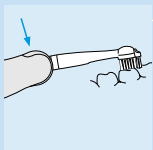
- 1 The Gum Protection System will keep a check on the pressure you exert on the brush.
- 2 If you press the brush too hard against your gums, you will feel and hear a click when the brush head flexes backwards. As soon as you reduce the pressure, the brush will return to its normal position.
- ▶ This will probably happen regularly when you first start using the brush, but after some time you will have discovered the optimal brushing pressure, and the brush will flex backwards less often.

Tip: If plaque is removed properly, the gums will have a pink colour and firmly surround the teeth. Healthy gums do not bleed when you eat or brush your teeth (see also the section 'taking optimal care of your teeth').


If, nevertheless, your gums bleed slightly the first few times you use the brush, this bleeding will generally stop after a few days. If it does not, consult your dentist or dental hygienist.

Minimum brushing time: 2 minutes

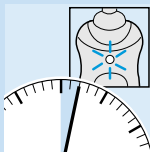
- ▶ In general, 2 minutes will be enough to clean your teeth properly. Therefore the brush is equipped with a red pilot light that will blink for 40 seconds when you have been brushing for 2 minutes.
- ▶ Naturally, you can also brush longer than 2 minutes.


After brushing

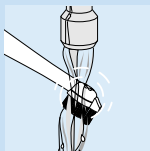
- 1 Let the toothbrush run while you rinse the brush head under the tap. Use lukewarm water (max. 60°C).

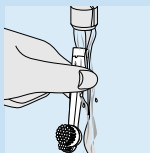- 2 Also rinse the inside of the brush head.

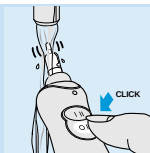- 3 Switch the toothbrush on again to rinse toothpaste residues from the drive shaft.

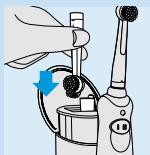
- 4 Store the brush head upside down in the storage compartment, with the coloured ring visible on top.

- 5 You can clean the charging unit with a moist cloth.

► First unplug the charging unit!




Brush head replacement

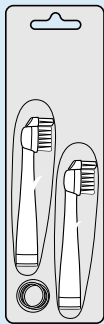- 1 Replace the brush head every three months for the best brushing results.
- The brush head should also be replaced when the bristles start to bend outwards.
 - New brush heads are available from the shop in which you have purchased your Philips Jordan toothbrush.

Tip: Worn bristles will move plaque around rather than remove it. They may also damage the gums. It is therefore important that you replace the brush head every three months.

Environmental aspects

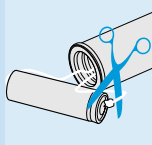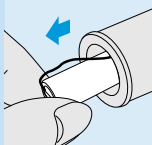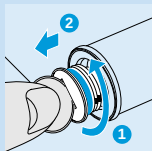
If you discard the appliance in due course, remove the battery and dispose of it in an environmentally safe way.

- ▶ NB: You can also hand the appliance in at a Philips Service Centre, which will dispose of the battery for you.
- ▶ Before you remove the battery, it must be completely empty. Let the appliance run until it stops before removing the battery.
- ▶ NB: Only unscrew the bottom of the appliance if you want to remove the battery.

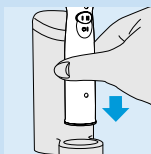After the bottom has been unscrewed, the appliance may not be watertight anymore. Therefore the guarantee will become invalid if the bottom of the appliance has been unscrewed. If undamaged, the sticker on the bottom of the appliance will prove that the appliance has not been opened.

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