
Always here to help you / 始终如一地为 您提供帮助  
Register your product and get support at / 请登录以下网址 并得到相应的帮助  
[www.philips.com/welcome](http://www.philips.com/welcome)

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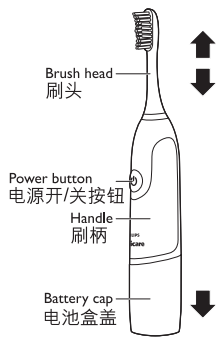READ ALL INSTRUCTIONS BEFORE USE

#### Warnings

- Use only as intended and with manufacturer attachment.
- Use only AA Alkaline batteries. Always replace with two new batteries; never mix old and new batteries.
- Supervise use by children and persons with reduced physical, sensory or mental capabilities.
- Discontinue use if any part of this product appears damaged in any way (brush head and handle).
- Use this product only for its intended use. This is not a toy.
- Consult dentist or physician if you experience

significant bleeding, discomfort, pain or medical concerns.

### SAVE THESE INSTRUCTIONS

#### Electromagnetic fields (EMF)

This appliance complies with current standards regarding electromagnetic fields (EMF). If used as directed, this appliance is safe.

#### OPERATING INSTRUCTIONS

*Note: Not intended for use on multiple patients in a dental practice or institution.*

- 1 Pull battery cap down. Install 2 size AA Alkaline batteries. Proper battery direction is indicated under cap. Do not use the product without battery cap installed.
- 2 Install brush head by placing on metal shaft and pushing down. Remove brush head by pulling up.

*Note: A gap between brush head and handle is normal.*

- 3 Avoid splattering: place brush in mouth before turning on.
- 4 Press power on/off button to turn on.
- 5 Brush normally with light pressure for 2 minutes.
- 6 Automatically turns off after two minutes.
- 7 Replace brush heads every 3 months or sooner if needed.

The toothbrush is safe to use on:

- Braces (brush heads wear out sooner when used on braces)
- Dental restorations (fillings, crowns, veneers)

#### CLEANING AND STORAGE

Clean head with water and wipe handle regularly with a damp cloth only.

Store in a cool dry location when not used.

#### DISPOSAL

At product end of life - do not throw away in trash can; recycle products and batteries appropriately.

#### WARRANTY

This product is warranted for 1 years from the date of purchase. What is not covered under the warranty:

- Brush heads
- Damage from misuse, abuse, neglect alterations, unauthorized repair/parts
- Normal wear and tear including discoloration or fading

For complete warranty details, visit [www.philips.com/support](http://www.philips.com/support) and download the operator's manual with full warranty details.

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