
  
**sonicare**<sup>®</sup>  
the **sonic** toothbrush

elite<sup>®</sup>  
e7000  
series

**PHILIPS**

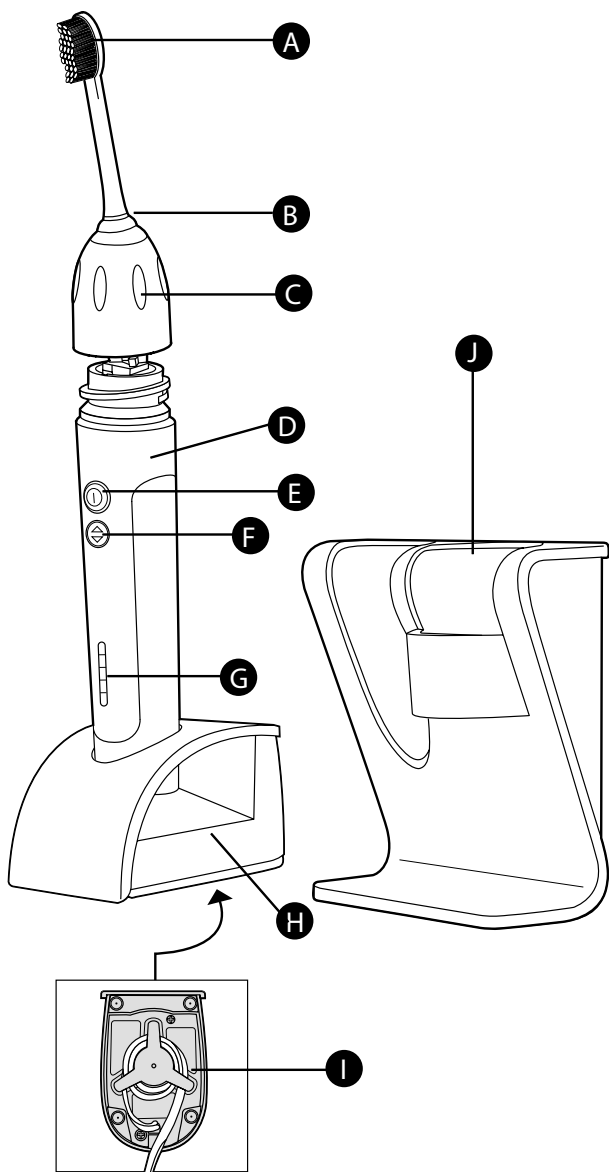**SONICARE ELITE 7000**

**ENGLISH 4**

**NORSK 16**

**SVENSKA 28**

**SUOMI 40**

**DANSK 52**

## General description

- A** Brush head
- B** Colour code
- C** Removable nut
- D** Sonic power handle
- E** Power on/off button
- F** Dual Speed Control button (specific models only)
- G** Charge indicator LEDs (specific models only)
- H** Charger
- I** Cord wrap facility
- J** Luxury holder (specific models only)

## Important

Read these instructions for use carefully before using the appliance and save them for future reference.

If you have any questions after reading this manual, please contact the Philips Customer Care Centre in your country or visit our website ([www.philips.com](http://www.philips.com)).

Sonicare® has been carefully engineered for both effectiveness and safety. However, basic safety precautions should always be followed when using any electrical product, especially when children are present.

## Danger

To reduce the risk of electrocution:

- 1) Do not place or store the charger where it can fall or be pulled into a bathtub or washbasin.
- 2) Do not place or drop the charger into water or other liquid.
- 3) Do not reach for a charger that has fallen into water. Unplug immediately.
- 4) Do not use while bathing.
- 5) Do not plug the charger into a socket that supplies a voltage for which the charger is unsuitable. It has been designed for use with a specific electrical voltage.
- 6) Refer to the bottom of the charger for voltage specification. Voltage converters or plug adapters DO NOT guarantee voltage compatibility. Never force the plug into a socket.

## Warning

To reduce the risk of burns, electrocution, fire or injury:

- 1) Close supervision is necessary when this product is used by or on children or the disabled.
- 2) Use this product only for its intended purpose as described in this booklet. Do not use attachments not recommended by the manufacturer: See voltage information above.


- 3) Never operate the charger if it has a damaged cord or plug, if it is not working properly, if it has been dropped or if it has been dropped into water. For further instructions, contact your local Philips Customer Care Centre.
- 4) Keep the cord away from heated surfaces.

### Medical warnings


- 1) Avoid using a brush head with crushed or bent bristles. Damaged bristles may break off while brushing. In order to prevent this, replace the brush head with a new one every 6 months or sooner if signs of wear appear.
- 2) Similar to starting with any new manual or electric oral hygiene product, consult your dental professional before using this appliance if you have had oral or gum surgery in the previous 2 months.
- 3) Contact your dental professional if excessive bleeding occurs after using this product or if bleeding continues to occur after 1 week of use.
- 4) If you have any questions concerning usage with a pacemaker or other implanted device, consult your physician or the implanted device manufacturer prior to use. Sonicare has been tested and is compliant with safety standards for electromagnetic devices.
- 5) The Sonicare toothbrush is a personal care appliance and is not intended for use on multiple patients in a dental office or institution.
- 6) This product is designed for cleaning teeth, gums and tongue only. Do not use it for any other purpose.

### Preparing for use

To attach the Sonicare brush head:

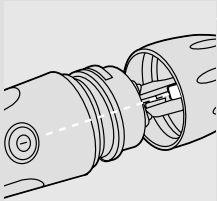
- 1** Align the front of the brush head with the front of the toothbrush handle. Match the bars on the brush head with the notches inside the top of the handle and slide the brush head onto the handle.

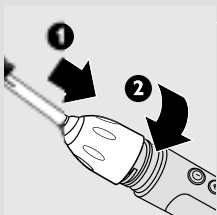
- 2** Turn the brush head nut clockwise until it is firmly tightened. If the brush head rattles during use, tighten the nut further.

## Charging the Sonicare®

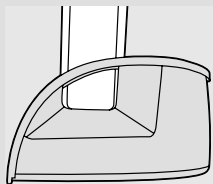Sonicare comes with the batteries completely discharged and should be charged for at least 24 hours during the first few days of use.

- 1 To charge your Sonicare, place the handle in the charger with the Sonicare logo facing forward. Make sure the bottom of the handle sits back securely in the charger, making complete contact.

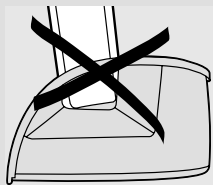- ▶ When the handle is placed incorrectly in the charger, the green charge indicator LED will still light up, but charging will not take place.

Handles of other Sonicare models do not fit in this charger and will therefore not charge effectively. Do not attempt to force a handle into this charger as this will cause damage.

It is recommended that you keep your Sonicare in the charger when not in use to maintain a full battery charge. It takes at least 24 hours to fully charge Sonicare's battery.

## Using the appliance

### Useful tips

- ▶ You can use your favourite toothpaste with your Sonicare.
- ▶ The brush head comes with a convenient travel cap. Remove the travel cap from the brush head before switching on your Sonicare.
- ▶ Place the brush head in your mouth before switching on your Sonicare, and be sure to switch it off before removing it from your mouth, otherwise toothpaste may splatter.
- ▶ Rinse the brush head and bristles after each use.
- ▶ Clean the brush head weekly, including under the nut.
- ▶ A fully charged Sonicare will provide approximately 2 weeks of typical use (2 minutes of brushing twice a day), making it easy to travel without the charger.

Remember: voltage converters and plug adapters DO NOT guarantee voltage compatibility. Always check for voltage compatibility.


### Tips for specific oral care needs

- ▶ **Stain removal**

In areas where staining occurs, a few seconds of additional brushing time may be spent to assist with stain removal.

- ▶ **Braces**

The brush should be placed at a 45-degree angle between your brackets and gumline. While applying light pressure, gently move the bristles in a small circular motion so that the long bristles reach

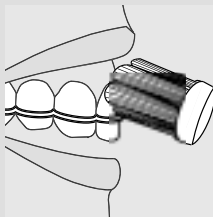between your teeth for 2 seconds. Then gently turn (roll) the brush handle so the bristles clean the bracket and the rest of the tooth before you move the bristles to the next tooth. Continue this until you have brushed all of your teeth.

#### ► Dental restorations

Use the Sonicare on dental restorations (e.g. implants, caps, crowns, bridges and fillings) as you would on your normal tooth surfaces to reduce plaque build-up and improve gum health. Work the bristles around all areas of the restoration for thorough cleaning.



#### ► Periodontal pockets

Your dental professional may have indicated that you have periodontal pockets, deep pockets where plaque bacteria collect. If this is the case, after brushing for a full 2 minutes, spend extra time gently brushing these areas to improve your gum health. You may want to take advantage of the Programmable Quadpacer® timer feature (see section 'Programmable Quadpacer timer').

### Proper brushing technique

Follow these simple steps to maximise your brushing experience:

- 1** Before switching the Sonicare on, place the bristles of the brush head on your gumline at a 45-degree angle.
- 2** Hold the handle with a light grip, switch the Sonicare on, and glide the brush across your teeth.
- 3** To maximise Sonicare's effectiveness apply only light pressure and let the Sonicare toothbrush do the brushing for you. Do not scrub.
- 4** The angled brush head and slender neck help you reach the back teeth and the tooth surfaces on the inside of the mouth. Be sure to adjust your brushing to use the angle to your advantage and keep as many bristles in contact with the teeth as possible.
- 5** To make sure you brush evenly throughout the mouth, divide your mouth into 4 sections: outside of upper teeth, inside of upper teeth, outside of lower teeth and inside of lower teeth. Begin brushing in section 1 (outside of upper teeth) and brush for 30 seconds before moving to section 2 (inside of upper teeth). Continue brushing the other two sections for 30 seconds until you have completed the 2-minute brushing cycle. If activated, the Programmable Quadpacer® timer will beep every 30 seconds to signal you to move to the next quadrant.

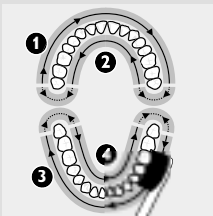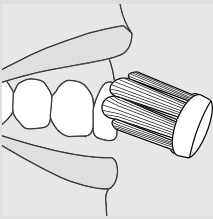- 6** After having completed the 2-minute brushing cycle, the Sonicare will switch off automatically. If you want, you can switch the toothbrush back on for another brushing cycle and spend additional time brushing the chewing surfaces of your teeth.
- 7** You may also brush your tongue, with the toothbrush switched on or off, as you prefer.
- 8** You may experience a slight tickling or tingling sensation after using an electronic toothbrush such as the Sonicare for the first time. As you become accustomed to brushing with the Sonicare this sensation will diminish.

## Special features

Features that simplify brushing:

### Easy-start® brushing power ramp-up

- ▶ All Sonicare models come with the Easy-start feature activated, which helps you get used to brushing with the Sonicare.
  - ▶ The brushing power will gently increase to full power over the first 14 brushings.
  - ▶ The Easy-start feature is active if you hear 2 beeps immediately after switching on your Sonicare. Once the appliance has reached optimal brushing power (after the 14th brushing), you will no longer hear the double beep at the beginning of the brushing cycle.
- ▶ *Note: Each of the first 14 brushings must be at least 1 minute in length to properly advance through the Easy-start ramp-up cycle.*

### To deactivate or reactivate the Easy-start feature:

- 1** Attach the brush head to the handle.
- 2** Place the handle in the plugged-in charger.

#### **3** To deactivate Easy-start:

Press and hold the on/off button for 2 seconds. You will hear 1 beep to indicate that the Easy-start feature has been deactivated.

#### **4** To reactivate Easy-start:

Press and hold the on/off button until you hear 2 beeps. This indicates that the Easy-start feature is active again.


Using the Easy-start feature beyond the initial ramp-up period is not recommended and reduces the effectiveness of the Sonicare in removing plaque.

The Easy-start feature should be deactivated when the Sonicare is used in clinical trials.

### Smartimer® 2-minute timer

- ▶ All Sonicare models come with the Smartimer feature, which automatically sets the brushing time to the 2 minutes recommended by dental professionals. After these 2 minutes, your Sonicare will automatically switch off. If you wish to brush for a longer period, simply press the on/off button and another 2-minute cycle will start.
- ▶ If you wish to pause or stop brushing during your 2-minute cycle, press the on/off button. If you press the on/off button again, the Smartimer will resume where you left off.
- ▶ The Smartimer will automatically reset to the beginning of a 2-minute cycle if:
  - brushing is paused for 30 seconds or longer;
  - the handle is returned to the charger.

### Quadpacer® 30-second interval timer

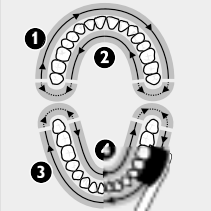- ▶ The Quadpacer feature ensures that you brush all sections of your mouth evenly and thoroughly. At 30-second intervals (at 30, 60 and 90 seconds in case of a 2-minute brushing cycle), you will hear a short beep and notice a brief pause in brushing action. This is your signal to move to the next section of your mouth.


### Programmable Quadpacer timer (specific models only)


The Programmable Quadpacer timer allows you to customise the Quadpacer performance in 1 of 3 ways:

- ▶ You may opt to keep the factory default setting of a 2-minute brushing time with Programmable Quadpacer beeps activated.
- ▶ You may choose to extend your brushing time to 2 1/2 minutes by adding an additional 30-second interval.
- ▶ You may choose a 2-minute brushing time, but with the Quadpacer beeps deactivated.

To change the Programmable Quadpacer settings:

- 1** Place the handle in the charger that is plugged into a live socket.

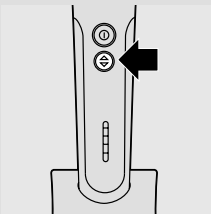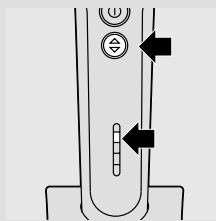
**2** Press the Dual Speed Control button briefly and note which LED lights up.

**3** Press and hold the Dual Speed Control button. Continue to hold the button while the LEDs slowly cycle through the timing interval setting options:

▶ Green LED in position 4: Programmable Quadpacer active (factory default setting).


▶ Green LED in top (5th) position: 2 1/2 minute setting - the Programmable Quadpacer active with an additional 30-second timing interval.

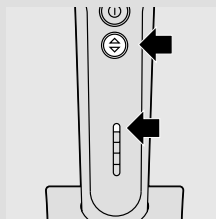▶ Yellow LED in bottom (1st) position: no timing intervals set - the Quadpacer not active.

**4** To program the desired setting, release the button when the LED corresponding to your choice lights up.

▶ *Note: Selecting 2 1/2 minutes brushing time reduces the number of brushings available from a charged battery.*
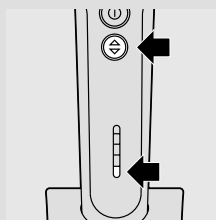
### Dual Speed Control (specific models only)


The Dual Speed Control feature allows you to choose between 2 speed settings: Normal speed for optimal brushing performance, or Reduced speed when an even more gentle brushing action is desired.

Each time you switch your Sonicare on, it will begin brushing at the Normal speed setting. To change the speed setting:

**1** From Normal to Reduced speed: press the Dual Speed Control button once during brushing.

**2** From Reduced to Normal speed: press the Dual Speed Control button a second time during brushing.


A distinctive two-tone beep indicates that the speed setting has been changed.

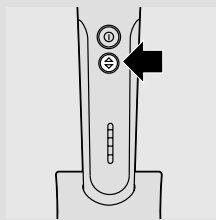It is recommended that you brush at the Normal speed setting to maintain optimal results.

The Normal speed setting should be selected when the Sonicare is used in clinical trials.

### Battery charge indicator (specific models only)

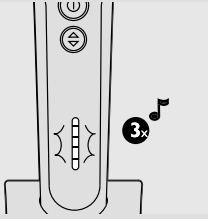- ▶ If the battery charge of your Sonicare is low, you will hear 3 beeps after you have brushed your teeth for 2 minutes and the charge indicator LED will blink yellow for 30 seconds.

A low charge indication signifies that 3 or fewer 2-minute brushings remain.

When not in use, keep the Sonicare in the charger to maintain a full battery charge.


### Multiple-LED charge indicator (specific models only)

- ▶ When the handle is placed in the charger, the green charge indicator LEDs will blink in a rising pattern as the battery charges. The number of blinking LEDs indicates the level to which the battery has been charged.

When the battery reaches full charge, the LEDs stop blinking and start to burn continuously. About once a minute, the charge indicator LEDs go out momentarily and then go back on. This indicates that the battery is being maintained at full charge. The full charge indication will persist until the handle is removed from the charger.

The blinking yellow LED indicates that fewer than 3 brushings remain.

### Single-LED charge indicator (specific models only)

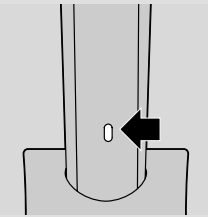- ▶ When the handle is placed in the charger, the charge indicator LED will blink green as the battery charges. When the battery reaches full charge, the charge indicator LED will burn green continuously until the handle is removed from the charger.

When your Sonicare is switched on or off, the charge indicator LED will blink yellow if the battery charge is low.

### Luxury brush holder (specific models only)


- ▶ Some Sonicare models come with a luxury brush holder.

- ▶ The luxury brush holder can be mounted to the wall with 2 screws (3.0 x 30 mm) and 2 wall plugs (not included)


## Cleaning

### Simple tips for care

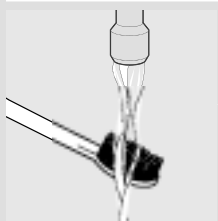
By following these simple steps, you should get years of worry-free brushing:

- 1** Rinse the brush head and bristles after each use.

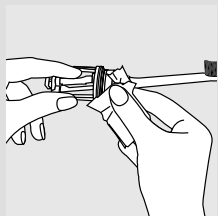
- 2** Clean the brush head weekly, including under the nut.

- 3** To remove the nut:

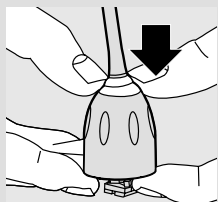
- ▶ Press down gently on the grey seal area.

Avoid using sharp objects to push on the seal as this may cause damage.

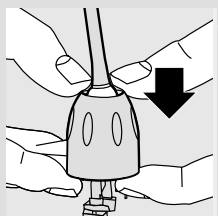- ▶ Push the nut up over the end of the brush head.

- 4** To reassemble the nut, slip the nut over the brush head shaft and insert the lower part of the brush head into the handle. Turn the nut clockwise until it is firmly tightened and the nut snaps into the seal.

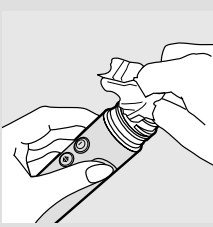- 5** Clean the handle weekly.
- 6** Clean the charger periodically.  
Unplug the charger before cleaning.
- 7** Examine the mains cord periodically. Do not use the charger if the mains cord is damaged.
- 8** Use the cord wrap facility in the bottom of the charger to store excess cord.

UNPLUG the charger before cleaning. After cleaning the charger completely dry the charger before connecting it to the mains.

Mild soap, water and a moist cloth are recommended for cleaning your Sonicare.

## Replacement

### Brush head

To maintain optimal brushing results, Sonicare brush heads should be replaced every 6 months or sooner if signs of wear appear:

Use only Sonicare Elite® series replacement brush heads.

## Environment

The handle contains a nickel cadmium battery. The battery must be recycled or disposed of in an environmentally safe way.

The battery inside your Sonicare cannot be replaced, but can be easily removed for recycling. At the end of product life and prior to disposal, please remove the battery from the handle by following the instructions below.

Please note that this process is NOT reversible.

- Do not throw the battery away with the normal household waste, but hand it in at an official collection point.**


For information about Sonicare's battery recycling programme or instructions on recycling your Sonicare, contact the Philips Customer Care Centre in your country.

### Removing the battery

Caution should be exercised when using any household tools (such as a screwdriver) to remove the battery.


- 1** Deplete the battery of any charge by removing the handle from the charger, switching the Sonicare on and allowing it to run

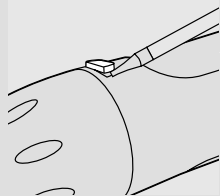
until it stops. Repeat this step until the product will no longer operate.


- 2** Remove the soft-touch rib by inserting a screwdriver under it and breaking it out.

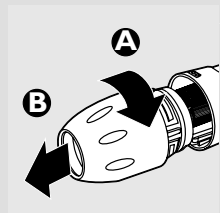- 3** Remove the nut from the brush head and place it onto the handle. Tighten the nut onto the handle (A) until the housing comes loose and then pull the nut forward (B).

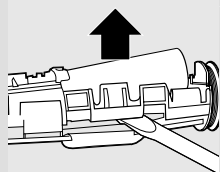- 4** Insert a screwdriver between the electronic circuit board and the plastic part close to each end of the battery; twist and break the connection between the battery and the electronic circuit board in 2 locations.

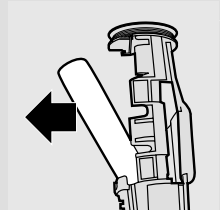- 5** Separate the rechargeable battery from the appliance by means of a screwdriver.

Please note that this process is NOT reversible

## Guarantee & service

If you need information or if you have a problem, please visit the Philips website at [www.philips.com](http://www.philips.com) or contact the Philips Customer Care Centre in your country (you will find its phone number in the worldwide guarantee leaflet). If there is no Customer Care Centre in your country, turn to your local Philips dealer or contact the Service Department of Philips Domestic Appliances and Personal Care BV.

## Guarantee restrictions

What is not covered by the terms of the international guarantee?

- Brush heads
- Damage caused by misuse, abuse, neglect or alterations
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading

## Troubleshooting

If you are unable to solve your problem by means of the troubleshooting guide below, contact your local Philips Customer Care Centre or online at [www.philips.com](http://www.philips.com).

| Problem   | Solution   |
|---|--|
| The Sonicare toothbrush does not work                                     | Recharge the Sonicare for at least 24 hours.   |
|   | Make sure the Sonicare is plugged into a live socket and not into a socket controlled by a light switch.   |
|   | You may need to re-attach or replace the brush head.   |
| I feel a tickling sensation when I use the Sonicare.                      | You may experience a slight tickling/tingling sensation when using Sonicare for the first time. As you become accustomed to brushing with Sonicare, this sensation will diminish.  |
| The brush head rattles when I use the Sonicare.                           | Sonicare should produce a gentle hum during operation. Tightening the brush head nut can eliminate excessive rattling. If the rattling does not stop, you may need to replace the brush head.  |
| The brush head chatters against my teeth.                                 | This can be avoided by keeping the brush head properly oriented, with the brush head placed on the gumline at an angle of 45 degrees.  |
| Brushing seems less powerful.   | The handle may need recharging.  |
|   | The Easy-start feature may be active and needs to be deactivated (see 'Special features').   |
|   | The Dual Speed Control button may be set to Reduced Speed.   |
|   | You may need to clean your brush head.   |
|   | You may need to replace your brush head.   |
| Unpleasant odour  | Clean the handle, brush head and charger thoroughly (see 'Cleaning').  |
|   | Replace the brush head with a new one.   |
| The LED blinks yellow. What does this mean?                               | The blinking yellow LED indicates that fewer than 3 brushings remain. Recharge the Sonicare for at least 24 hours.   |
| The appliance has an extra 30 seconds of brushing time. What is this for? | This function provides additional time for brushing troublesome areas where plaque collects and gum disease starts. It is recommended you use the extra 30 seconds of brushing time to re-brush troublesome areas, massage your gums, or brush your tongue (see 'Special features'). |

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


*Outside North America:*  
+1-425-396-2000  
**www.philips.com**

Printed in U.S.A.  
Printed on 100% recycled paper

Protected by U.S. and international patents. Other patents pending.

© 2003 Philips Oral Healthcare, Inc. All rights reserved. PHILIPS and the Philips shield are registered trademarks of Koninklijke Philips Electronics N.V. Easy-start, Quadpacer, Smartimer, Sonicare, Elite, and the Sonicare logo are registered trademarks of Philips Oral Healthcare, Inc.