
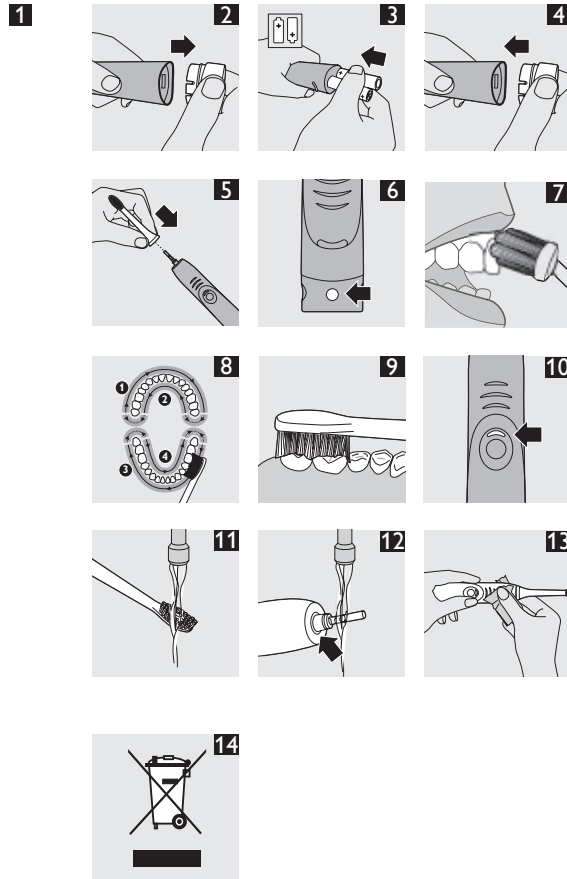
Sonicare  
 4000+ Series

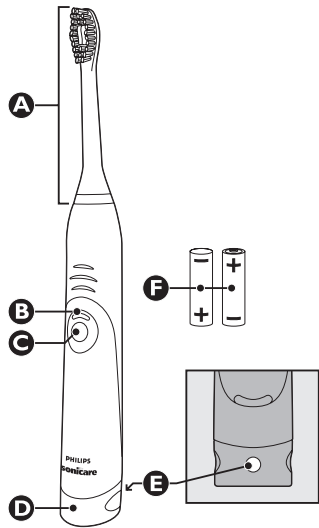**Important**

Read this user manual carefully before you use the appliance and save it for future reference.

- ▶ If you have had oral or gum surgery in the previous 2 months, consult your dentist before you use the toothbrush.
- ▶ Consult your dentist if excessive bleeding occurs after using this toothbrush or if bleeding continues to occur after 1 week of use.
- ▶ The Sonicare toothbrush complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, consult your doctor or the manufacturer of the implanted device if you have any questions about using the Sonicare.
- ▶ This appliance has only been designed for cleaning teeth, gums and tongue. Do not use it for any other purpose. Discontinue use of this product and contact a physician if discomfort or pain is experienced.
- ▶ Consult your physician prior to using the Sonicare if you have medical concerns.
- ▶ The Sonicare toothbrush is a personal care device and is not intended for use on multiple patients in a dental practice or institution.
- ▶ If this product appears damaged in any way (brush head, toothbrush handle) discontinue use. This product contains no serviceable parts. If repair is needed please refer to the Guarantee and service section of this manual and contact one of our representatives.
- ▶ Do not use the toothbrush in the bath or shower.
- ▶ Stop using a brush head with crushed or bent bristles. Replace the brush head every 3 months or sooner if signs of wear appear. Do not use other brush heads than the ones recommended by the manufacturer.
- ▶ This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- ▶ Children should be supervised to ensure that they do not play with the appliance.
- ▶ If your toothpaste contains peroxide, baking soda or bicarbonate (common in whitening toothpastes), thoroughly clean the brush head and the handle with soap and water after each use.

**Electromagnetic fields (EMF)**

This Philips appliance complies with all standards regarding electromagnetic fields (EMF). If handled properly and according to the instructions in this user manual, the appliance is safe to use based on scientific evidence available today.

**To avoid damage due to battery leakage:**

- Do not expose the appliance to high temperatures, e.g. caused by direct sunlight or heating devices.
- Do not leave empty batteries in the appliance.

**The Sonicare (Fig. 1)**

- A Replaceable HydroClean brush head
- B Battery indicator (specific types only)
- C Power on/off button
- D Battery compartment cap
- E Colour code sticker location
- F 2 AA 1.5-volt batteries

**Getting started**

**Installing or replacing batteries**

- 1 Remove the battery compartment cap from the bottom of the handle by pulling on the grips (Fig. 2).
- 2 Insert two AA 1.5-volt batteries into the handle. Make sure the + and - poles of the batteries match the indication on the inside of the handle (Fig. 3).
- 3 Firmly push the battery compartment cap back onto the handle (Fig. 4).

**Attaching the brush head**

- 1 Align the front of the brush head with the front of the toothbrush handle.

- 2 Push the brush head onto the shaft (Fig. 5).

**Colour code sticker**

- If you have multiple Sonicare battery toothbrushes in your home, you can attach the white or silver sticker to the battery cap. This helps you identify your handle (Fig. 6).

**Using the Sonicare**

- 1 Wet the bristles and apply a small amount of toothpaste.
- 2 Place the toothbrush bristles against the teeth at a slight angle towards the gumline (Fig. 7).
- 3 Hold the handle with a light grip and switch on the Sonicare. Gently move the brush head slowly across the teeth in a small back and forth motion so the longer bristles reach between your teeth. Continue this motion throughout the 2-minute brushing cycle.
- 4 To maximise effectiveness, apply only light pressure and let the toothbrush do the brushing for you. Do not scrub.
- 5 To make sure you brush evenly throughout the mouth, divide your mouth into 4 sections: outside top teeth, inside top teeth, outside bottom teeth and inside bottom teeth. Brush each section for 30 seconds for a total brushing time of 2 minutes. Start with section 1 (outside top teeth) and then brush section 2 (inside top teeth). Continue with section 3 (outside bottom teeth) and finally brush section 4 (inside bottom teeth) (Fig. 8).
- 6 After you have completed the 2-minute brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs. You may also brush your tongue, with the toothbrush switched on or off, as you prefer (Fig. 9).

The Sonicare is safe to use on:

- Braces (brush heads wear out sooner when used on braces)
- Dental restorations (fillings, crowns, veneers)
- Periodontal pockets

**Features**

**Smartimer**

- After 2 minutes, the Smartimer automatically switches off the toothbrush.
- To pause or stop brushing during the 2-minute cycle, press the on/off button. If you press the on/off button again within 30 seconds, the Smartimer picks up where you left off.

**Battery indicator (select models)**

The battery indicator shows the amount of power left in the batteries.

- Solid green light indicates there are at least 16 two-minute brushings available (Fig. 10).
- Flashing green light indicates there are less than 16 brushings available.
- Three beeps at the end of a two-minute brushing indicate that there are less than 4 two-minute brushings available.

*Note: The toothbrush switches on for one second and then immediately switches off again to indicate that you have to replace the batteries.*

*Note: The number of brushings you get from new batteries may vary depending on the type of batteries you use.*

**Easy-start**

- This Sonicare model has the Easy-start feature, which helps you get used to brushing with the Sonicare.
- The brushing power gently increases to full power over the first 14 brushings.

*Note: Each of the first 14 brushings must be at least 1 minute in length to advance through the Easy-start ramp-up cycle properly.*

*Note: Prior to a clinical test, ensure that the Easy-start process has been completed, attach a new brush head and install new AA alkaline batteries.*

**Cleaning**

Do not place brush heads and the handle in the dishwasher:

- 1 Rinse the brush head and bristles after each use. Let the brush head air-dry only (Fig. 11).
- 2 Remove the brush head once a week and clean the connection between the brush head and the handle (Fig. 12).

Do not use sharp objects to press on the rubber seal, as this may cause damage.

- 3 Clean the remainder of the handle periodically with mild soap and a moist cloth (Fig. 13).

**Storage**

If you are not going to use the toothbrush for a long time (a month or longer), remove the batteries, clean the appliance and store it in a cool and dry place.

**Replacement**

**Brush head**

- Replace Sonicare brush heads every 3 months to achieve optimal results.
- Use only Sonicare HydroClean replacement brush heads with the Sonicare B400 model.

**Environment**

- Do not throw away the appliance with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment (Fig. 14).
- Batteries contain substances that may pollute the environment. Do not throw away batteries with the normal household waste, but hand them in at an official collection point. Always remove the batteries before you discard or hand in the appliance at an official collection point.

**Guarantee and service**

If you need information or if you have a problem, please visit the Philips website at [www.philips.com](http://www.philips.com) or contact the Philips Customer Care Centre in your country (you find its phone number in the worldwide guarantee leaflet). If there is no Customer Care Centre in your country, go to your local Philips dealer or contact the Service Department of Philips Domestic Appliances and Personal Care BV.

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(C) 2006 Philips Oral Healthcare, Inc. all rights reserved.
PHILIPS and the Philips shield are registered trademarks of Koninklijke Philips Electronics N.V.,
Sonicare and the Sonicare logo are trademarks of Philips Oral Healthcare, Inc. and/or KPENV.
Philips Oral Healthcare, Inc.
U.S.A.
Printed in the U.S.A.