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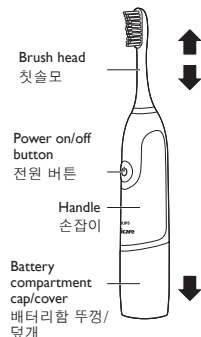
Always here to help you /

Register your product and get support at /

[www.philips.com/welcome](http://www.philips.com/welcome)

[www.philips.com/Sonicare](http://www.philips.com/Sonicare)

©2014 Koninklijke Philips N.V. (KPNV). All rights reserved. Philips and the Philips shield are trademarks of KPNV. Sonicare and the Sonicare logo are trademarks of Philips Oral Healthcare, Inc. and/or KPNV.

## ENGLISH

given supervision or instruction concerning use of the appliance by a person responsible for their safety.

-

### Caution

-

If you have medical concerns, consult your doctor before you use the Sonicare.

-

Consult your dentist if you experience significant bleeding or discomfort when you use the Sonicare.

-

Only use this appliance as intended.

-

Only use genuine Philips Sonicare replacement brush heads that are suitable for this model.

- Only use 1.5V AA alkaline batteries. Always replace both batteries; never mix old and new batteries.
- Remove the batteries from the appliance if you are not going to use it for some time.

### Electromagnetic fields (EMF)

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

### Using the appliance

- 1 Pull the battery compartment cap off the handle. Insert two 1.5V AA alkaline batteries in the handle and reattach the battery compartment cap.**

*Note: Make sure the + and - poles point in the right direction.*

- 2 To avoid splattering, place the toothbrush in your mouth before you switch it on.**

- 3 Press the on/off button to switch on the appliance.**

- 4 Brush normally with light pressure for 2 minutes.**

- 5 The appliance switches off automatically after 2 minutes.**

- 6 Replace the brush head every 3 months or sooner if needed.**

The toothbrush is safe to use on:

- Braces (brush heads wear out more quickly when used on braces)
- Dental restorations (fillings, crowns, veneers)

### Cleaning and storage

Clean the brush head with water after every use. Wipe the handle with a damp cloth regularly. Store the appliance in a cool, dry place when it is not in use.

### Ordering accessories

To buy accessories or spare parts, visit [www.shop.philips.com/service](http://www.shop.philips.com/service) or go to your Philips dealer. You can also contact the Philips Consumer Care Centre in your country (see the worldwide guarantee leaflet for contact details).

### Recycling

Do not throw away the appliance with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment.

### Guarantee and support

If you need information or support, please visit [www.philips.com/support](http://www.philips.com/support) or read the separate worldwide guarantee leaflet.

### Guarantee restrictions

The terms of the international guarantee do not cover the following:

- Brush heads.
- Replacement batteries.
- Damage caused by use of unauthorised replacement parts.
- Damage caused by misuse, abuse, neglect, alterations or unauthorised repair.
- Normal wear and tear, including chips, scratches, abrasions, discolouration or fading.