
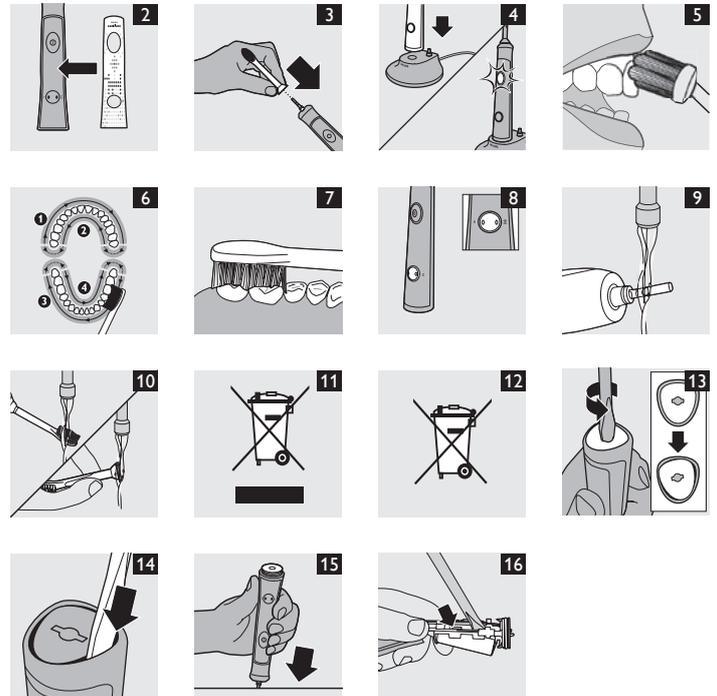
for kids

300 Series
 Rechargeable sonic toothbrush

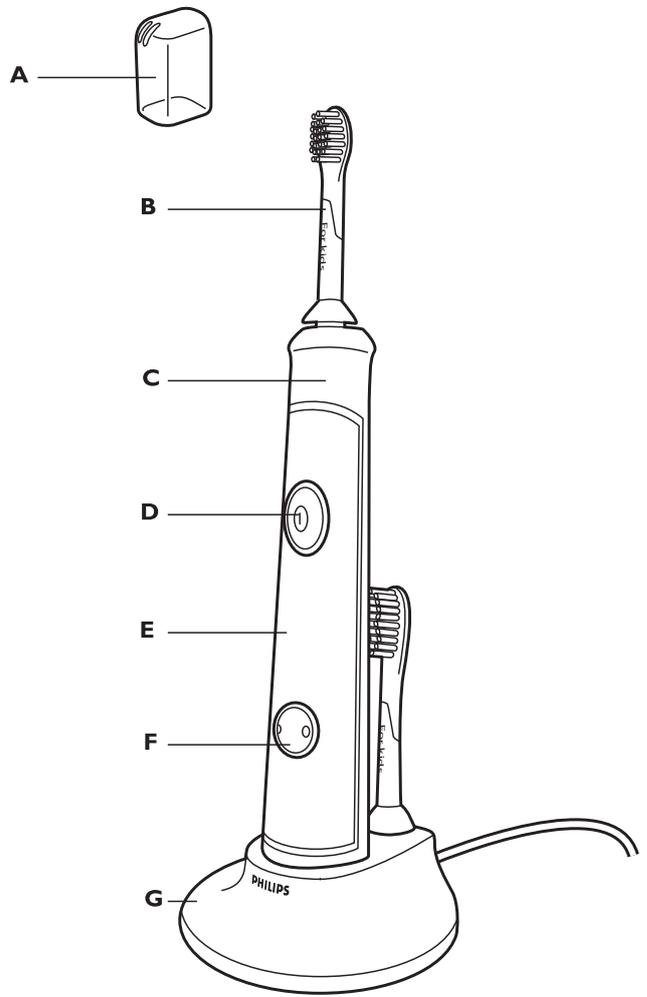www.philips.com/Sonicare

©2012 Koninklijke Philips Electronics NV (KPEENV). All rights reserved.
 PHILIPS and the Philips shield are registered trademarks of KPEENV.
 DiamondClean Easy-start, Quadpacer, Sonicare and the Sonicare logo
 are trademarks of Philips Oral Healthcare, Inc. and/or KPEENV.

Protected by U.S. and international patents. Other patents pending.
 Reg. U.S. Pat. & Tm. Off.

Introduction

Congratulations on your purchase and welcome to Philips! To fully benefit from the support that Philips offers, register your product at www.philips.com/welcome.

General description (Fig. 1)

- A Hygienic travel cap
- B Brush head
- C Handle with soft grip
- D Power on/off button with battery charge indicator
- E Interchangeable panel
- F Brushing mode button with mode lights
- G Charger with brush head holder

Important

Read this user manual carefully before you use the appliance and save it for future reference.

Danger

- Keep the charger away from water. Do not place or store it over or near water contained in a bathtub, washbasin, sink etc. Do not immerse the charger in water or any other liquid. After cleaning, make sure the charger is completely dry before you connect it to the mains.
- **This appliance is not a toy. Check brush heads regularly for cracks. Cracked brush heads may present a choke hazard.**

Warning

- The mains cord cannot be replaced. If the mains cord is damaged, discard the charger. Always have the charger replaced with one of the original type in order to avoid a hazard.
- If the appliance is damaged in any way (brush head, toothbrush handle and/or charger), stop using it.
- This appliance contains no serviceable parts. If the appliance is damaged, contact the Consumer Care Centre in your country (see chapter 'Guarantee and service').
- Do not use the charger outdoors or near heated surfaces.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.

Caution

- If you have had oral or gum surgery in the previous 2 months, consult your dentist before you use the toothbrush.
- Consult your dentist if excessive bleeding occurs after using this toothbrush or if bleeding continues to occur after 1 week of use.
- If you have medical concerns, consult your doctor before you use the Sonicare.
- The Sonicare toothbrush complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician or the device manufacturer of the implanted device prior to use.
- This appliance has only been designed for cleaning teeth, gums and tongue. Do not use it for any other purpose. Stop using the appliance and contact your doctor if you experience any discomfort or pain.
- Do not use other brush heads than the ones recommended by the manufacturer.
- If your toothpaste contains peroxide, baking soda or other bicarbonate (common in whitening toothpastes), thoroughly clean the brush head and the handle with soap and water after each use. Otherwise the plastic may crack.
- Do not clean brush heads, handle, charger, charger cover or interchangeable panels in the dishwasher or a microwave.

Electromagnetic fields (EMF)

This Philips appliance complies with all standards regarding electromagnetic fields (EMF). If handled properly and according to the instructions in this user manual, the appliance is safe to use based on scientific evidence available today.

Preparing for use

Applying the interchangeable panel

The Sonicare For Kids comes with 3 different interchangeable panels.

- 1 Remove the adhesive backing from the panel.

- 2 Apply the panel to the white section on the front of the handle (Fig. 2).

Note: Make sure you line up the cut-out areas on the panel with the buttons on the front of the handle.

Reapplying a panel that has been removed

- 1 Rinse the panel in warm water; you may also use a mild detergent to clean it.
- 2 Let the panel air dry.
- 3 Reapply the panel to the white section on the front of the handle.

Attaching the brush head

- 1 Align the brush head so the bristles point in the same direction as the front of the handle. (Fig. 3)
- 2 Firmly press the brush head onto the metal shaft until it stops.

Note: There is a small gap between the brush head and the handle.

Charging the Sonicare

- 1 Put the mains plug of the charger in the wall socket.
- 2 Place the handle on the charger (Fig. 4).
 - ▶ The battery charge indicator in the power on/off button flashes to indicate that the toothbrush is charging.
 - ▶ When the Sonicare is fully charged, the battery charge indicator stops flashing and remains on.

Note: If the battery charge of the Sonicare is low, you hear 3 beeps and the battery charge indicator flashes rapidly for 30 seconds after the brushing cycle.

Note: To keep the battery fully charged at all times, you may keep the Sonicare on the charger when not in use. When the battery is fully charged, the charger uses a minimal amount of energy. It takes at least 24 hours to fully charge the battery.

Tip: For extra convenience during travel, you can remove the travel charger and use it without the charger cover.

Using the Sonicare For Kids

Brushing instructions

For children under age 8, it is recommended that a parent or adult caregiver perform toothbrushing.

- 1 Wet the bristles and apply a small amount of toothpaste.
- 2 Place the toothbrush bristles against the teeth at a slight angle towards the gumline (Fig. 5).
- 3 Press the power on/off button to switch on the Sonicare.
- 4 Apply light pressure to maximise Sonicare's effectiveness and let the Sonicare toothbrush do the brushing for you.
- 5 Gently move the brush head slowly across the teeth in a small back and forth motion so the longer bristles reach between the teeth. Continue this motion throughout the brushing cycle.

Note: To make sure you brush evenly throughout the mouth, divide the mouth into 4 sections using the KidPacer feature (see chapter 'Features').

- 6 Begin brushing in section 1 (outside top teeth) until a tune indicates it is time to move to section 2 (inside top teeth). At the next tune, move to section 3 (outside bottom teeth). The next tune indicates it is time to move to section 4 (inside bottom teeth). The toothbrush automatically switches off at end of the brushing cycle (Fig. 6).
- 7 After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs. You may also brush your tongue, with the toothbrush switched on or off, as you prefer (Fig. 7).

The Sonicare For Kids is safe to use on:

- Braces (brush heads wear out sooner when used on braces)
- Dental restorations (fillings, crowns, veneers)

Brushing modes

The toothbrush automatically starts in the default high mode.

- 1 Press the brushing mode button to toggle between modes. (Fig. 8)
 - ▶ The green mode light indicates the selected mode.
 - High mode: Full-strength mode designed for older/independent brushers.
 - Low mode: Reduced power mode used for training and introducing younger children to the sonic experience.

Features

KidTimer

- The KidTimer plays a congratulatory tune at the end of the pre-determined brushing time. Additionally, the KidTimer automatically switches off the Sonicare when the brushing cycle is complete.
- When in high mode, the KidTimer is set at 2 minutes.
- When in low mode, the KidTimer is set at 1 minute and slowly increases to 2 minutes over time. This slow increase in brushing time helps young children reach the dentist recommended 2 minutes of brushing.

Note: Instruct your child to brush until they hear the congratulatory tune and the toothbrush switches off.

KidPacer

The KidPacer uses a short tune to alert the brusher to move to the next brushing quadrant of the mouth. This ensures that the teeth get a thorough, overall clean, see section 'Brushing instructions' in chapter 'Using the Sonicare For Kids'.

Easy-start

- The Sonicare For Kids comes with the Easy-start feature activated.
- The Easy-start feature gently increases the power over the first 14 brushings.

Note: Each of the first 14 brushings must be at least 1 minute in length to properly advance through the Easy-start ramp-up cycle.

Deactivating or activating the Easy-start feature

- 1 Attach the brush head to the handle.
- 2 Place the handle in the plugged-in charger.
 - To deactivate Easy-start:
Press and hold the power on/off button for 2 seconds. You hear 1 beep to indicate that the Easy-start feature has been deactivated.
 - To activate Easy-start:
Press and hold the power on/off button for 2 seconds. You hear 2 beeps to indicate that the Easy-start feature has been activated.

Note: Using the Easy-start feature beyond the initial ramp-up period is not recommended and reduces Sonicare's effectiveness in removing plaque.

Cleaning

Do not clean brush heads, handle, charger, charger cover or interchangeable panels in the dishwasher or a microwave.

Toothbrush handle

- 1 Remove the brush head and rinse the metal shaft area with warm water (Fig. 9).

Do not push on the rubber seal on the metal shaft with sharp objects, as this may cause damage.

- 2 Use a damp cloth to wipe the entire surface of the handle.

Brush head

- 1 Rinse the brush head and bristles after each use (Fig. 10).
- 2 Remove the brush head from the handle and rinse the brush head connection at least once a week with warm water (Fig. 10).

Charger

- 1 Unplug the charger before you clean it.
- 2 Use a damp cloth to wipe the surface of the charger.

Storage

- If you are not going to use the Sonicare for an extended period of time, unplug the charger, clean it and store it in a cool and dry place away from direct sunlight.
- You can store a brush head on the brushing head holder at the back of the charger.

Replacement

Brush head

- Replace Sonicare brush heads every 3 months to achieve optimal results.
- Use only Sonicare For Kids replacement brush heads.

Environment

- Do not throw away the appliance with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment (Fig. 11).
- The built-in rechargeable battery contains substances that may pollute the environment. Always remove the battery before you discard and hand in the appliance at an official collection point. Dispose of the battery at an official collection point for batteries. If you have trouble removing the battery, you can also take the appliance to a Philips service centre. The staff of this centre will remove the battery for you and will dispose of it in an environmentally safe way (Fig. 12).

Removing the rechargeable battery

Please note that this process is not reversible.

To remove the rechargeable battery, you need a flat-head (standard) screwdriver.

Observe basic safety precautions when you follow the procedure below. Be sure to protect your eyes, hands, fingers and the surface on which you work.

- 1 To deplete the rechargeable battery of any charge, remove the handle from the charger, switch on the Sonicare and let it run until it stops. Repeat this step until you can no longer switch on the Sonicare.
- 2 Insert a screwdriver into the slot located in the bottom of the handle and turn anticlockwise until you see a gap between the bottom cap and the handle. (Fig. 13)
- 3 Insert the screwdriver into the gap and pry the bottom cap out of the handle (Fig. 14).
- 4 Hold the handle upside down and push down on the shaft to release the internal components of the handle (Fig. 15).
- 5 Insert the screwdriver under the circuit board, next to the battery connections, and twist to break the connections. Remove the circuit board and pry the battery from the plastic carrier (Fig. 16).

Guarantee and service

If you need service or information or if you have a problem, please visit the Philips website at www.philips.com or contact the Philips Consumer Care Centre in your country (you find its phone number in the worldwide guarantee leaflet). If there is no Consumer Care Centre in your country, go to your local Philips dealer.

Guarantee restrictions

The terms of the international guarantee do not cover the following:

- Brush heads
- Interchangeable panels
- Damage caused by misuse, abuse, neglect, alterations or unauthorised repair
- Normal wear and tear; including chips, scratches, abrasions, discolouration or fading

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