# Always here to help you

Register your product and get support at  
[www.philips.com/welcome](http://www.philips.com/welcome)


Question?  
Contact  
Philips

**Philips Sonicare  
AirFloss**

**PHILIPS  
sonicare**

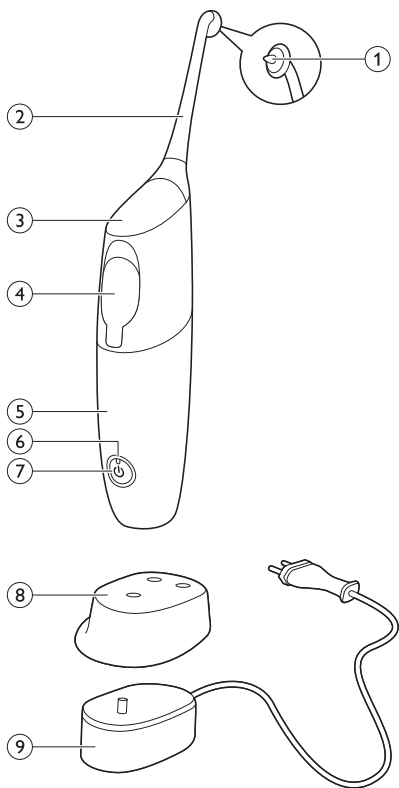**ENGLISH 6**

**DANSK 13**

**SUOMI 20**

**NORSK 27**

**SVENSKA 34**

## Introduction

Congratulations on your purchase and welcome to Philips! To fully benefit from the support that Philips offers, register your product at

**[www.philips.com/welcome](http://www.philips.com/welcome)**.

Use the AirFloss as part of your oral care routine. In addition to brushing, it helps reduce plaque.

## Important

Read this user manual carefully before you use the appliance and save it for future reference.

### Danger

- Keep the charger away from water. Do not place or store it over or near water contained in a bathtub, washbasin, etc. Do not immerse the charger in water or any other liquid. After cleaning, make sure the charger is completely dry before you connect it to the mains

### Warning

- The mains cord cannot be replaced. If the mains cord is damaged, discard the charger. Always have the charger replaced with one of the original type in order to avoid a hazard.
- If the appliance is damaged in any way (AirFloss nozzle, handle and/or charger), stop using it.
- This appliance contains no serviceable parts. If the appliance is damaged, contact the Consumer Care Centre in your country (see chapter 'Guarantee and support').
- Do not use the charger outdoors or near heated surfaces.

- This appliance can be used by children aged from 8 years and above and by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and if they understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children unless they are older than 8 and supervised.

### **Caution**

- If you have had recent oral or gum surgery, suffer from an ongoing dental condition or had extensive dental work in the previous 2 months, consult your dentist before you use this appliance.
- Consult your dentist if excessive bleeding occurs after using this appliance or if bleeding continues to occur after 1 week of use.
- If you have medical concerns, consult your doctor before you use this appliance.
- This Philips appliance complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician or the appliance manufacturer of the implanted appliance prior to use.
- Do not use attachments other than the ones recommended by the manufacturer.
- Do not exert excessive pressure on the guidance tip of the nozzle.

## 8 ENGLISH

- To avoid physical injury, do not fill the reservoir with hot water. Do not put foreign objects in the reservoir.
- Only use the AirFloss for its intended purpose as described in this user manual.
- To avoid physical injury, do not direct the spray under the tongue or into the ear, nose, eye or other sensitive area.
- Do not clean the AirFloss nozzle, the handle, the charger, the water reservoir and its cover in the dishwasher or a microwave.
- To avoid an ingestion hazard, do not put isopropyl alcohol or other cleaning fluids in the reservoir.
- The AirFloss is a personal care device and is not intended for use on multiple patients in a dental practice or institution.

### **Electromagnetic fields (EMF)**

---

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

### **General description (Fig. 1)**

- 1 Guidance tip
- 2 AirFloss nozzle
- 3 Activation button
- 4 Reservoir and cover
- 5 Handle
- 6 Charge indicator
- 7 Power on/off button
- 8 Removable nozzle holder (specific types only)
- 9 Charger

*Note: The contents of the box may vary based on the model purchased.*

## Charging

Charge the AirFloss for 24 hours before first use.

- 1** Put the plug of the charger in the wall socket (Fig. 2).
- 2** Place the handle on the charger (Fig. 3).
  - ▶ The charge indicator flashes green to indicate that the appliance is charging (Fig. 4).

*Note: When it is fully charged, the AirFloss has an operating time of up to 2 weeks (or 14 uses).*

*Note: When the charge indicator flashes yellow, the battery is low and needs to be recharged (fewer than 3 uses left).*

## Preparing for use

- 1** Align the nozzle so the tip faces the front of the handle and press the nozzle down onto the handle ('click') (Fig. 5).

*Note: Replace the nozzle every 6 months for optimal results. Also replace the nozzle if it becomes loose or no longer clicks onto the handle.*

- 2** Open the cover of the reservoir ('click') (Fig. 6).
- 3** Fill the reservoir with mouthwash or water (Fig. 7).

*Note: To operate effectively, the AirFloss has to contain either mouthwash or water.*

*Tip: For a fresher experience, use mouthwash instead of water.*

*Note: Do not use mouthwash that contains isopropyl myristate (e.g. Dentyf Active) or high concentrations of essential oils, as it may damage the appliance. This includes all mouthwash sold in glass and ceramic packaging.*

- 4** Close the cover of the reservoir ('click') (Fig. 8).
- 5** Press the on/off button (Fig. 9).
  - ▶ The green charge indicator lights up.
- 6** Press and hold the activation button until spray comes out of the tip of the nozzle (Fig. 10).

## Using the appliance

- 1 Make sure the AirFloss is switched on (Fig. 9).
- 2 Place the tip between the teeth and apply gentle pressure. Slightly close your lips over the nozzle when it is in place, to avoid splashing (Fig. 11).

*Tip: Use your finger to press the activation button.*

*Note: The AirFloss is intended to be used on the front/outside surface of the teeth only.*

- 3 Press the activation button to deliver a single burst of air and micro-droplets of mouthwash or water between the teeth. Alternatively, hold the activation button to deliver continuous bursts of mouthwash or water approximately every 1 second (select models) (Fig. 12).
- 4 Slide the guidance tip along the gum line until you feel it settle between the next two teeth.
- 5 Continue this procedure for all spaces including behind your back teeth (Fig. 13).
- 6 To switch off the AirFloss, press the on/off button again.

### Automatic shut-off function

---

The automatic shut-off function automatically switches off the AirFloss if it has not been used for 4 minutes.

## Cleaning

Do not clean the nozzle, handle or charger in the dishwasher.

Do not use cleaning agents to clean the appliance.

## Handle and nozzle

---

- 1 Remove the nozzle from the handle. Rinse the nozzle after each use (Fig. 14).
- 2 Open the reservoir and rinse it under the tap to eliminate residue (Fig. 15).

*Tip: You can also use a cotton bud to wipe the reservoir.*

*Note: Do not immerse the handle in water.*

- 3 Wipe the entire surface of the handle with a damp cloth (Fig. 16).

## Charger

---

- 1 Unplug the charger before you clean it (Fig. 17).
- 2 Wipe the surface of the charger with a damp cloth (Fig. 18).

## Storage

If you are not going to use the AirFloss for an extended period of time, follow the steps below.

- 1 Open the cover of the reservoir and empty the reservoir.
- 2 Press the activation button until no more spray comes out of the nozzle.
- 3 Unplug the charger (Fig. 17).
- 4 Clean the AirFloss nozzle, handle and charger.
- 5 Store the AirFloss in a cool, dry place away from direct sunlight.

## Ordering accessories

To buy accessories or spare parts, visit [www.shop.philips.com/service](http://www.shop.philips.com/service) or go to your Philips dealer. You can also contact the Philips Consumer Care Centre in your country (see the worldwide guarantee leaflet for contact details).

### Recycling

- This symbol on a product means that the product is covered by European Directive 2012/19/EU (Fig. 19).
- This symbol means that the product contains a built-in rechargeable battery covered by European Directive 2006/66/EC which cannot be disposed of with normal household waste. Follow the instructions in section 'Removing the rechargeable battery' to remove the battery (Fig. 20).
- Inform yourself about the local separate collection system for electrical and electronic products and rechargeable batteries. Follow local rules and never dispose of the product and rechargeable batteries with normal household waste. Correct disposal of old products and rechargeable batteries helps prevent negative consequences for the environment and human health.

### Removing the rechargeable battery

---

Please note that this process is not reversible.

- 1** To deplete the battery, repeatedly press the activation button until the AirFloss no longer produces any bursts of air (Fig. 21).
- 2** Use a screwdriver to pry open the handle (Fig. 22).
- 3** Cut the wires and remove the battery from the appliance. Observe basic safety precautions (Fig. 23).

### Guarantee and support

If you need information or support, please visit **[www.philips.com/sonicare](http://www.philips.com/sonicare)** or read the separate worldwide guarantee leaflet.

#### Guarantee restrictions

The terms of the international guarantee do not cover the following:

- AirFloss nozzle.
- Damage caused by misuse, abuse, neglect, alterations or unauthorised repair.
- Normal wear and tear, including chips, scratches, abrasions, discolouration or fading.

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©2015 Koninklijke Philips N.V. (KPNV). All rights reserved.
Philips and the Philips shield are trademarks of KPNV.
AirFloss, Sonicare and the Sonicare logo are trademarks
of Philips Oral Healthcare, Inc. and/or KPNV.