

PHILIPS

A close-up, artistic photograph of a Philips Sensiflex power toothbrush head. The brush head is shown in motion, with a significant horizontal blur effect that conveys a sense of rapid vibration. The bristles are densely packed and appear as a soft, greyish-white mass. The brush head is attached to a silver-colored handle. The background is a smooth, light blue gradient, which makes the metallic and white tones of the toothbrush stand out. The overall aesthetic is clean, modern, and emphasizes the technology and motion of the product.

sensiflex
power toothbrush

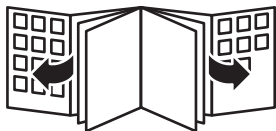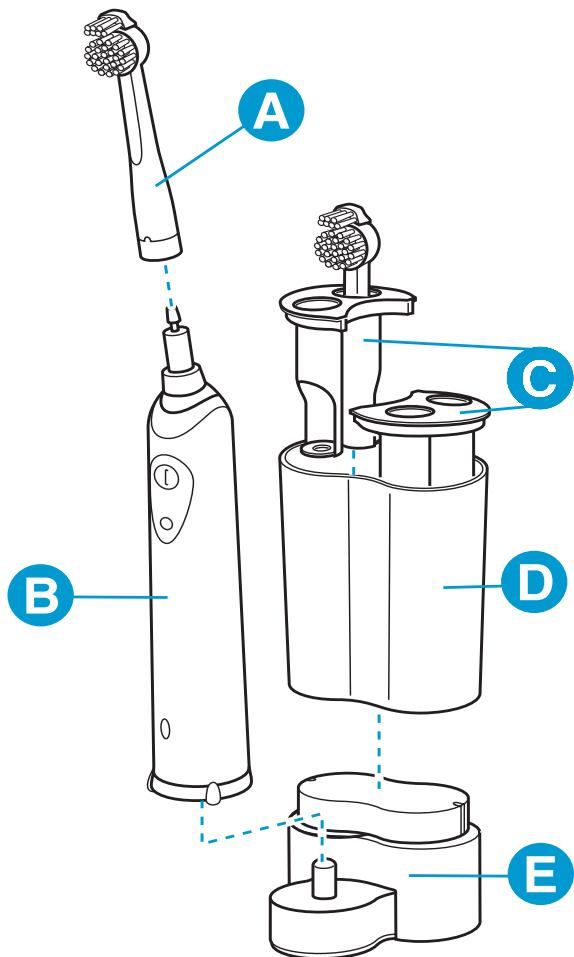**HX2585, HX2538, HX1526, HX1518, HX1515, HX1514, HX1513,
HX1511, HX1510, HX1507**

ENGLISH 6

DANSK 12

DEUTSCH 18

SUOMI 24

NORSK 30

SVENSKA 36

Introduction

This unique electric toothbrush removes plaque better than a normal, manual toothbrush. It also helps to protect your gums. When you exert too much pressure while brushing, the brush head will automatically flex back. Because of this clear signal, you will quickly adopt the right brushing force. Clinical tests have shown that, thanks to this system, brushing with the Philips Sensiflex results in a better and safer brushing technique.

General description (Fig. 1)

- A** Brush attachment
- B** Handle
- C** Holder for 2 brush attachments
- D** Storage container
- E** Charger
- F** Holder for 2 brush attachments

Important

Read these instructions for use carefully before using the appliance and save them for future reference.

- D** The cord of the charging unit cannot be replaced. If the cord is damaged, the charging unit must be discarded.
- D** Make sure the socket is live when you recharge the toothbrush. The electricity supply to sockets in mirror cabinets over washbasins or shaver sockets in bathrooms may be cut off when the light is turned off.

Preparing for use

Mounting the appliance to the wall

- 1** Mount the charging unit to the wall, using 2 screws (3.0x30mm) and 2 wall plugs (S5) (not included). (Fig. 2)

Charging the appliance

- 1** Charge the appliance for at least 16 hours before using it for the first time.

- ▶ Make sure the appliance is switched off when you start charging it. (Fig. 3)
- ▶ You can put the toothbrush in the stand and let it charge continuously.

However, when you first start using the appliance, it is recommended that you use it until the battery is empty. To do so, use the appliance in the normal way without recharging it between the brushing sessions.

Recharge the appliance when the motor almost stops running. Repeat this procedure 3 times in a row.

To extend the lifetime of the battery, repeat this procedure every 6 months.

- ▶ **Types HX2585, HX2538 and HX1526 only:** when the battery is being charged, the green pilot light is on. (Fig. 4)
- ▶ You can put the toothbrush in the stand and let it charge continuously.

Using the appliance

- 1 Mount the brush attachment onto the toothbrush by inserting it into the bayonet catch and turning it clockwise a quarter of a turn. (Fig. 5)

This will make it impossible for the brush attachment to become detached due to vibration.

- 2 Put toothpaste on the brush head.
 - 3 Place the brush head in your mouth. (Fig. 6)
 - 4 Switch the appliance on by pressing the button and start brushing your teeth. (Fig. 7)
- ▶ For the best brushing method, minimal brushing time and brushing functions of the appliance, see the 'Brushing tips' section below.

Brushing Tips

Brushing

- 1 Place the brush head against the tooth in such a way that it covers the surface of the tooth as well as the gum line. (Fig. 8)

- 2** Hold the brush head in this position for a few seconds and then move on to the next tooth. (Fig. 9)

The round brush head removes plaque from the teeth and beneath the gum line.

The individually moving Active Tip cleans between the teeth.

- 3** Guide the brush head carefully along the teeth and molars. (Fig. 10)

- ▶ Make sure you clean all surfaces of the teeth (inside, outside and chewing surface) as well as the gum line properly. A fixed brushing pattern will prevent you from skipping certain areas.

Gum Protection system

- 1** The Gum Protection system guarantees that you do not exert too much pressure on the teeth and gums. (Fig. 11)

- 2** If you exert too much pressure on the toothbrush, you will feel and hear the brush head flexing backwards. As soon as you reduce the pressure, the brush head will return to its normal position. This will protect your gums. (Fig. 12)

This will probably happen a few times when you first start using the toothbrush. After some time you will have discovered the optimal brushing pressure and the brush head will flex backwards less often.

Tip: If plaque is removed properly, the gums will have a pink colour and will firmly surround the teeth. Healthy gums do not bleed when you eat or brush your teeth (see also the section 'Optimal dental hygiene'). If, nevertheless, your gums bleed slightly the first few times you use the brush, this bleeding will generally stop after a few days. If it does not, consult your dentist or dental hygienist.

Minimum brushing time: 2 minutes

In general, 2 minutes will be enough to clean your teeth properly. Naturally, you can also brush longer than 2 minutes.

- ▶ A red pilot light will blink to indicate that 2 minutes have passed (types HX2585, HX2538 and HX1526 only). (Fig. 13)

Optimal dental hygiene

For optimal dental hygiene, it is recommended to use dental floss and toothpicks in addition to the electric toothbrush.

You will find more information about Philips Oral Healthcare products and other Philips products on our Web site: www.philips.com.

Cleaning and storage

- 1 Switch the toothbrush on and rinse the brush attachment under the tap with lukewarm water (max. 60 °C). (Fig. 14)
- 2 Switch the toothbrush off. Detach the brush attachment and rinse the inside. (Fig. 15)
- 3 Switch the toothbrush on again and rinse the drive shaft under the tap to wash away toothpaste residues. (Fig. 16)
- 4 Types HX1526, HX1515, HX1514, HX1513, HX1511, HX1510 and HX1507 only: store the brush attachment on the pin on the charging unit. Types HX2585, HX2538 and HX1526 only: store the brush attachments (max. 4) in the storage container. (Fig. 17)
 - Types HX2585 and HX2538 only: you can also remove one of the brush holders to be able to store a toothpaste tube in the storage container. (Fig. 18)

The storage container and the brush holders can be removed from the charging unit and cleaned in the dishwasher.

- 5 Clean the charging unit with a moist cloth. (Fig. 19)

Unplug the charging unit before you start cleaning it.

Replacement

- 1 Replace the brush attachment every 3 months for the best brushing results.
 - The brush attachment should also be replaced when the bristles start to bend outwards.

You can obtain new brush attachments from the shop in which you purchased the toothbrush.

All round Philips replacement brush attachments can be used.

Environment

- ▶ Do not throw the appliance away with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this you will help to preserve the environment. (Fig. 20)
 - ▶ The built-in rechargeable battery contains substances that may pollute the environment. Always remove the battery before discarding the appliance and handing it in at an official collection point. Dispose of the battery at an official collection point for batteries. If you have trouble removing the battery, you can also take the appliance to a Philips service centre, which will remove the battery for you and will dispose of it in an environmentally safe way.
- 1 The battery must be completely empty when you remove it. You should therefore let the appliance run until it stops before removing the battery.
 - 2 Pick up the charging unit and place the toothbrush on the pin. Turn the toothbrush anti-clockwise to unscrew the bottom of the toothbrush. (Fig. 21)
 - 3 Turn the toothbrush upside down and press it down on a hard surface to remove the battery holder. (Fig. 22)
 - 4 Open the lid of the battery holder by hand or by means of a screwdriver. (Fig. 23)
 - 5 Lift the battery out of the battery holder by means of a screwdriver. (Fig. 24)

After the bottom has been unscrewed, the appliance may not be watertight anymore and the guarantee will become invalid.

Guarantee & service

If you need information or if you have a problem, please visit the Philips web site at www.philips.com or contact the Philips Customer Care Centre in your country (you will find its phone number in the worldwide guarantee leaflet). If there is no Customer Care Centre in your country,

turn to your local Philips dealer or contact the Service Department of Philips Domestic Appliances and Personal Care BV.

Troubleshooting

If you are unable to solve your problem by means of the troubleshooting guide below, please contact the nearest Philips service centre, the Philips Customer Care Centre in your country or contact us online at www.philips.com.

Problem	Solution
The toothbrush does not work at all, although it has been in the charging unit overnight.	<p>Make sure that the plug of the charging unit has been properly inserted into the wall socket.</p> <p>Make sure the socket is live (see also the section 'Important').</p> <p>Make sure the toothbrush is switched off and recharge it for at least 16 hours.</p>
The motor of the toothbrush works, but the brush head doesn't.	Check if the brush attachment has been properly connected to appliance (see also the section 'Using the appliance').
The brushing power of the appliance has greatly decreased.	The toothbrush must be recharged. If you normally leave the toothbrush in the charging unit continuously, you may have to revive the battery as described in the section 'Charging'.
The teeth-cleaning performance of the toothbrush is not optimal anymore.	Remove the brush attachment and start using a new one.

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