# Always here to help you

Register your product and get support at  
[www.philips.com/welcome](http://www.philips.com/welcome)


Question?  
Contact  
Philips

**Philips Sonicare  
3 Series**

**PHILIPS**  
**sonicare**
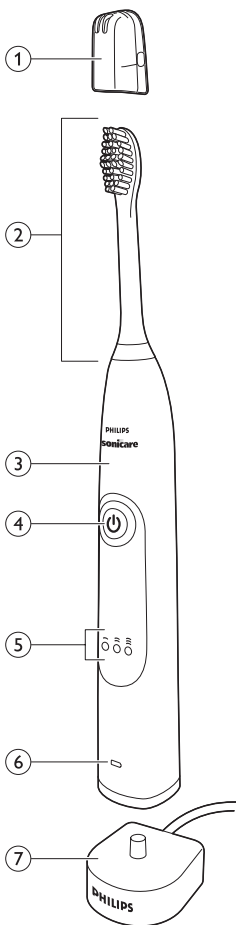
<b>ENGLISH</b>	<b>6</b>
<b>БЪЛГАРСКИ</b>	<b>16</b>
<b>ČEŠTINA</b>	<b>28</b>
<b>EESTI</b>	<b>38</b>
<b>HRVATSKI</b>	<b>48</b>
<b>MAGYAR</b>	<b>58</b>
<b>LIETUVIŠKAI</b>	<b>68</b>
<b>LATVIEŠU</b>	<b>78</b>

## Introduction

Congratulations on your purchase and welcome to Philips! To fully benefit from the support that Philips offers, register your product at [www.philips.com/welcome](http://www.philips.com/welcome).

## Important

Read this user manual carefully before you use the appliance and save it for future reference.

### **Danger**

- Keep the charger away from water. Do not place or store it over or near water contained in a bathtub, washbasin, sink etc. Do not immerse the charger in water or any other liquid. After cleaning, make sure the charger is completely dry before you connect it to the mains.

### **Warning**

- The mains cord cannot be replaced. If the mains cord is damaged, discard the charger.
- Always have the charger replaced with one of the original type in order to avoid a hazard.
- Do not use the charger outdoors or near heated surfaces.
- If the appliance is damaged in any way (brush head, toothbrush handle or charger), stop using it. This appliance contains no serviceable parts. If the appliance is damaged, contact the Consumer Care Centre in your country (see chapter 'Guarantee and support').

- This appliance can be used by children aged from 8 years and above and by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and if they understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children unless they are older than 8 and supervised.

### **Caution**

- Do not clean the brush head, the handle or the charger in the dishwasher.
- If you have had oral or gum surgery in the previous 2 months, consult your dentist before you use the toothbrush.
- Consult your dentist if excessive bleeding occurs after using this toothbrush or if bleeding continues to occur after 1 week of use. Also consult your dentist if you experience discomfort or pain when you use the Sonicare.
- The Sonicare toothbrush complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician or the device manufacturer of the implanted device prior to use.
- If you have medical concerns, consult your doctor before you use the Sonicare.

## 8 ENGLISH

- This appliance has only been designed for cleaning teeth, gums and tongue. Do not use it for any other purpose. Stop using the appliance and contact your doctor if you experience any discomfort or pain.
- The Sonicare toothbrush is a personal care device and is not intended for use on multiple patients in a dental practice or institution.
- Stop using a brush head with crushed or bent bristles. Replace the brush head every 3 months or sooner if signs of wear appear. Do not use other brush heads than the ones recommended by the manufacturer.
- If your toothpaste contains peroxide, baking soda or bicarbonate (common in whitening toothpastes), thoroughly clean the brush head with soap and water after each use. This prevents possible cracking of the plastic.

### **Electromagnetic fields (EMF)**

---

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

## General description (Fig. 1)

- 1 Hygienic travel cap
- 2 Brush head with unique identifier
- 3 Handle
- 4 Power on/off button
- 5 Intensity level LEDs
- 6 Battery level indicator
- 7 Charger

*Note: The contents of the box may vary based on the model purchased.*

## Preparing for use

### Attaching the brush head

---

- 1 Align the brush head so the bristles face the front of the handle (Fig. 2).
- 2 Firmly press the brush head down on the metal shaft until it stops.

*Note: Sonicare brush heads come with easily identifiable icons to identify your brush head.*

### Charging the appliance

---

- 1 Put the plug of the charger in a wall socket.
- 2 Place the handle on the charger (Fig. 3).
  - ▶ The battery level indicator flashes to indicate that the toothbrush is charging (Fig. 4).

*Note: It can take up to 48 hours to fully charge the battery, but the Sonicare may be used before it is fully charged.*

## Using the appliance

### Brushing instructions

---

- 1 Wet the bristles and apply a small amount of toothpaste.
- 2 Place the toothbrush bristles against the teeth at a slight angle towards the gumline (Fig. 5).
- 3 Press the power on/off button to switch on the Sonicare.
- 4 Apply light pressure to maximise Sonicare's effectiveness and let the Sonicare toothbrush do the brushing for you. Do not scrub.
- 5 Gently move the brush head slowly across the teeth in a small back and forth motion so the longer bristles reach between the teeth. Continue this motion throughout the brushing cycle.

*Note: To make sure you brush evenly throughout the mouth, divide the mouth into 4 sections using the Quadpacer feature (see chapter 'Features').*

- 6 Brush each section for 30 seconds for a total brushing time of 2 minutes. Start with section 1 (outside top teeth) and then brush section 2 (inside top teeth). Continue with section 3 (outside bottom teeth) and finally brush section 4 (inside bottom teeth) (Fig. 6).

*Note: After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs. You may also brush your tongue, with the toothbrush switched on or off, as you prefer.*

The Sonicare is safe to use on:

- Braces (brush heads may wear out more quickly when used on braces)
- Dental restorations (fillings, crowns, veneers, implants)

*Note: When the Sonicare is used in clinical studies, the handle has to be fully charged, set to the high intensity level and the Easy-start feature must be deactivated.*

## Intensity levels

---

- **Low:** 1 LED indicator light (Fig. 7)
- **Medium:** 2 LED indicator lights (Fig. 8)
- **High:** 3 LED indicator lights (Fig. 9)

*Note: To adjust the intensity level, press the power on/off button to let your Sonicare cycle through the intensity levels in the following order: High, Medium, Low, High etc. Each touch of the button allows you to select the next available intensity level.*

## Battery status (when handle is not on charger)

---

- A yellow battery level indicator indicates a low charge.
- A green battery level indicator indicates a high charge.
- The battery level indicator turns off when the battery is in between a low and a high charge.

*Note: When you hear 3 beeps and the yellow LED of the battery level indicator flashes, this means the battery charge is below 10% and the toothbrush should be recharged.*

*Note: When you hear 2 sets of 5 beeps, this means the Sonicare will not operate and needs to be recharged.*

*Note: To keep the battery fully charged at all times, you may keep the Sonicare on the charger when not in use.*

## Features

### Quadpacer

---

- The Quadpacer is an interval timer that has a short beep and pause to remind you to brush the 4 sections of your mouth. The Quadpacer feature works at all 3 intensity levels (Fig. 10).

### Smartimer

---

The Smartimer indicates that your brushing cycle is complete when it automatically switches off the toothbrush at the end of the brushing cycle. Dentists recommend brushing at least 2 minutes twice a day.

*Note: After brushing for 1.5 seconds, pressing the power on/off button causes the toothbrush to pause. After a pause of 30 seconds, the Smartimer resets.*

### Easy-start

---

The Easy-start feature gently increases power over the first 14 brushings to let you get used to brushing with the Sonicare.

The Easy-Start feature only works at the High intensity level.

This Sonicare model comes with the Easy-start feature deactivated.

#### **To activate or deactivate the Easy-start feature:**

- 1** Attach your brush head to the handle.
- 2** Place the handle in the plugged-in charger.
- 3** Activate or deactivate Easy-start:
  - ▶ To activate Easy-start: Press and hold the power on/off button for 2 seconds. You hear 2 beeps and the battery level indicator flashes green to indicate that the Easy-start feature has been activated.
  - ▶ To deactivate Easy-start: Press and hold the power on/off button for 2 seconds. You hear 1 beep and the battery level indicator flashes yellow to indicate that the Easy-start feature has been deactivated.

*Note: Each of the first 14 brushings must be at least 1 minute in length to properly move through the Easy-start ramp-up cycle.*

*Note: Using the Easy-start feature beyond the initial ramp-up period is not recommended and reduces Sonicare's effectiveness in removing plaque.*

## Cleaning

Do not clean the brush head, handle or charger in the dishwasher.

### Toothbrush handle

---

- 1 Remove the brush head and rinse the metal shaft area with warm water (Fig. 11).

Do not push on the rubber seal on the metal shaft with sharp objects, as this may cause damage.

- 2 Use a damp cloth to wipe the entire surface of the handle.

### Brush head

---

- 1 Rinse the brush head and bristles after each use (Fig. 12).
- 2 Remove the brush head from the handle and rinse the brush head connection at least once a week with warm water.

### Charger

---

- 1 Unplug the charger.
- 2 Use a damp cloth to wipe the surface of the charger.

## Storage

If you are not going to use the appliance for an extended period of time, remove the mains plug from the wall socket. Then clean the appliance and store it in a cool and dry place away from direct sunlight.

## Replacement

### Brush head

---

Replace Sonicare brush heads every 3 months to achieve optimal results. Only use genuine Philips Sonicare replacement brush heads that are suitable for this model.

### Ordering accessories

To buy accessories or spare parts, visit [www.shop.philips.com/service](http://www.shop.philips.com/service) or go to your Philips dealer. You can also contact the Philips Consumer Care Centre in your country (see the worldwide guarantee leaflet for contact details).

### Recycling

- This symbol on a product means that the product is covered by European Directive 2012/19/EU (Fig. 13).
- This symbol means that the product contains a built-in rechargeable battery covered by European Directive 2006/66/EC which cannot be disposed of with normal household waste. Follow the instructions in section 'Removing the rechargeable battery' to remove the battery (Fig. 14).
- Inform yourself about the local separate collection system for electrical and electronic products and rechargeable batteries. Follow local rules and never dispose of the product and rechargeable batteries with normal household waste. Correct disposal of old products and rechargeable batteries helps prevent negative consequences for the environment and human health.

### Removing the rechargeable battery

---

**Please note that this process is not reversible.**

To remove the rechargeable battery, you need a small standard screwdriver and tape. Observe basic safety precautions when you follow the procedure below.

- 1** To deplete the rechargeable battery of any charge, remove the handle from the charger, turn on the Sonicare and let it operate until it stops. Repeat this step until you can no longer turn on the Sonicare.
- 2** Insert a small standard screwdriver into the bottom of the handle and turn it counter-clockwise to remove the screw. Pull the end cap off of the handle (Fig. 15).

- 3** Insert the screwdriver between the handle and the internal frame. Tilt the screwdriver to loosen the handle (Fig. 16).
- 4** Remove the internal frame from the handle. (Fig. 17)
- 5** Remove the rubbers that are on top and below the battery (Fig. 18).
- 6** Insert the screwdriver between the battery and the frame. Tilt the screwdriver to pop the battery out of the frame (Fig. 19).
- 7** Place some tape over the ends of the battery to avoid any short-circuiting due to residual energy in the battery (Fig. 20).  
The rechargeable battery can now be recycled or disposed of and the rest of the product discarded appropriately.

## Guarantee and support

If you need information or support, please visit [www.philips.com/support](http://www.philips.com/support) or read the separate worldwide guarantee leaflet.

### Guarantee restrictions

---

The terms of the international guarantee do not cover the following:

- Brush heads.
- Damage caused by use of unauthorised replacement parts.
- Damage caused by misuse, abuse, neglect, alterations or unauthorised repair.
- Normal wear and tear, including chips, scratches, abrasions, discolouration or fading.

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©2015 Koninklijke Philips N.V. (KPNV). All rights reserved.
Philips and the Philips shield are trademarks of KPNV.
Sonicare and the Sonicare logo are trademarks
of Philips Oral Healthcare, Inc. and/or KPNV.