# Always here to help you

Register your product and get support at  
[www.philips.com/welcome](http://www.philips.com/welcome)


Question?  
Contact  
Philips

## Philips FlexCare Rechargeable Sonic Toothbrush

**PHILIPS**  
**sonicare**

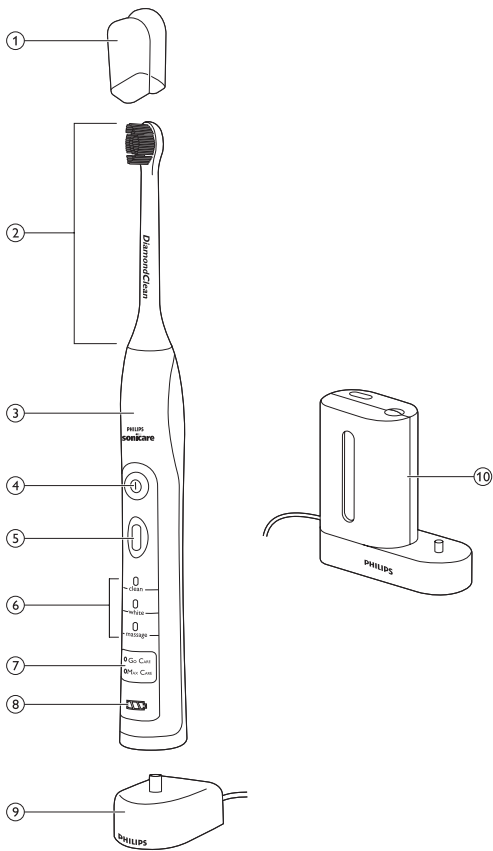**Philips FlexCare Rechargeable Sonic Toothbrush**

---

<b>ENGLISH</b>	<b>6</b>
<b>DEUTSCH</b>	<b>20</b>
<b>ΕΛΛΗΝΙΚΑ</b>	<b>36</b>
<b>ESPAÑOL</b>	<b>53</b>
<b>FRANÇAIS</b>	<b>68</b>
<b>ITALIANO</b>	<b>83</b>
<b>NEDERLANDS</b>	<b>98</b>
<b>PORTUGUÊS</b>	<b>114</b>

## Introduction

Congratulations on your purchase and welcome to Philips! To fully benefit from the support that Philips offers, register your product at [www.philips.com/welcome](http://www.philips.com/welcome).

## Important

Read this user manual carefully before you use the appliance and save it for future reference.

### **Danger**

- Keep the charger and/or sanitiser away from water. Do not place or store it over or near water contained in a bathtub, washbasin, sink etc. Do not immerse the charger and/or sanitiser in water or any other liquid. After cleaning, make sure the charger and/or sanitiser is completely dry before you connect it to the mains.

### **Warning**

- The mains cord cannot be replaced. If the mains cord is damaged, discard the charger and/or sanitiser.
- Always have the charger and/or sanitiser replaced with one of the original type in order to avoid a hazard.
- Do not use the charger and/or sanitiser outdoors or near heated surfaces. Do not operate the sanitiser where aerosol (spray) products are being used or where oxygen is being administered.

- If the appliance is damaged in any way (brush head, toothbrush handle, charger and/or sanitiser), stop using it. This appliance contains no serviceable parts. If the appliance is damaged, contact the Consumer Care Centre in your country (see chapter 'Guarantee and support').
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children unless they are older than 8 and supervised.
- Discontinue use of the sanitiser if the UV light bulb remains on when the door is open. UV light can be harmful to the human eye and skin. Keep the sanitiser out of the reach of children at all times.

### **Caution**

- Do not clean the brush head, the handle, the charger, the charger cover and/or the sanitiser in the dishwasher.
- If you have had oral or gum surgery in the previous 2 months, consult your dentist before you use the toothbrush.

## 8 ENGLISH

- Consult your dentist if excessive bleeding occurs after using this toothbrush or if bleeding continues to occur after 1 week of use. Also consult your dentist if you experience discomfort or pain when you use the Sonicare.
- The Sonicare toothbrush complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician or the device manufacturer of the implanted device prior to use.
- If you have medical concerns, consult your doctor before you use the Sonicare.
- This appliance has only been designed for cleaning teeth, gums and tongue. Do not use it for any other purpose. Stop using the appliance and contact your doctor if you experience any discomfort or pain.
- The Sonicare toothbrush is a personal care device and is not intended for use on multiple patients in a dental practice or institution.
- Stop using a brush head with crushed or bent bristles. Replace the brush head every 3 months or sooner if signs of wear appear.
- Do not use other brush heads than the ones recommended by the manufacturer.
- If your toothpaste contains peroxide, baking soda or bicarbonate (common in whitening toothpastes), thoroughly clean the brush head

with soap and water after each use. This prevents possible cracking of the plastic.

- Do not insert objects into the air openings of the sanitiser or block these air openings by placing it on a soft surface or in a position or location where objects will block these openings (examples of objects: magazine, tissue paper and hair).
- The UV light bulb is hot during and immediately after the sanitising cycle. Do not touch the UV light bulb when it is hot.
- Do not operate the sanitiser without the protective screen in place to avoid contact with a hot bulb.
- If the sanitiser has been used for 3 consecutive UV clean cycles (i.e. 3 times in a row), switch it off for at least 30 minutes before you start another UV clean cycle.

### **Electromagnetic fields (EMF)**

---

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

### **General description (Fig. 1)**

- 1** Hygienic travel cap
- 2** Brush head
- 3** Handle with soft grip
- 4** Power on/off button
- 5** Personalised brushing button
- 6** Brushing modes (specific types only)

## 10 ENGLISH

- 7** Brushing routines (specific types only)
- 8** Deluxe battery level indicator
- 9** Travel charger
- 10** UV sanitiser with integrated charger and cord wrap (specific types only)
  - Not shown: UV light bulb
  - Not shown: Sanitiser drip tray
  - Not shown: Protective screen for UV light bulb

*Note: The contents of the box will vary based on the type purchased*

### Preparing for use

#### Attaching the brush head

---

- 1** Align the brush head so the bristles face the front of the handle. (Fig. 2)
- 2** Firmly press the brush head down on the metal shaft until it stops.

*Note: When the brush head is attached, it is normal that there is a small gap between the brush head and the handle.*

#### Charging the appliance

---

- 1** Plug the charger in a wall socket.
- 2** Place the handle on the charger (Fig. 3).
  - ▶ The battery level indicator flashes to indicate that the toothbrush is charging.

#### **Deluxe battery level indicator:**

Indicates amount of charge left in battery

- 3 green LEDs: 75-100%
- 2 green LEDs: 50-74%
- 1 green LED: 25-49%
- 1 flashing yellow LED: less than 25%

*Note: If the battery charge of the Sonicare is low, you hear 3 beeps and 1 LED on the battery level indicator flashes yellow for 30 seconds after the 2-minute brushing cycle.*

*Note: To keep the battery fully charged at all times, you may keep the Sonicare on the charger when not in use. It takes at least 24 hours to fully charge the battery.*

## Using the appliance

### Brushing instructions

- 1** Wet the bristles and apply a small amount of toothpaste.
- 2** Place the toothbrush bristles against the teeth at a slight angle towards the gumline (Fig. 4).
- 3** Press the power on/off button to switch on the Sonicare.

*Note: In order to avoid toothpaste splatter, please make sure the toothbrush is in contact with your teeth before you switch on the Sonicare.*

- 4** Apply light pressure to maximise Sonicare's effectiveness and let the Sonicare toothbrush do the brushing for you.
- 5** Gently move the brush head slowly across the teeth in a small back and forth motion so the longer bristles reach between your teeth. Continue this motion throughout your brushing cycle.

*Note: To make sure you brush evenly throughout the mouth, divide your mouth into 4 sections using the Quadpacer feature (see chapter 'Features').*

- 6** Brush each section for 30 seconds for a total brushing time of 2 minutes. Start with section 1 (outside top teeth) and then brush section 2 (inside top teeth). Continue with section 3 (outside bottom teeth) and finally brush section 4 (inside bottom teeth) (Fig. 5).
- 7** After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs. You may also brush your tongue, with the toothbrush switched on or off, as you prefer (Fig. 6).

The Sonicare is safe to use on:

- Braces (brush heads wear out sooner when used on braces)
- Dental restorations (fillings, crowns, veneers)

### Personalising your brushing experience

---

Sonicare automatically starts in the default Clean mode. To personalise your brushing:

- 1 Switch on the Sonicare, press the personalised brushing button to toggle between modes and routines.
- ▶ The green LED indicates the selected mode or routine.

*Note: When the toothbrush is switched on, you can switch between modes but not between routines. Routines should be selected before you brush your teeth.*

### Brushing modes

---

#### Clean mode

Standard mode for superior teeth cleaning.

#### Sensitive mode (specific types only)

Gentle, yet thorough, cleaning of sensitive teeth and gums.

#### White mode (specific types only)

2 minutes of Clean mode, with an additional 30 seconds of White mode to focus on your visible front teeth.

#### White mode brushing instructions

- 1 Brush the first 2 minutes as explained in section 'Brushing instructions'.
- 2 After the 2 minutes of Clean mode, the White mode starts with a change in brushing sound and motion. This is your signal to start brushing the upper front teeth for 15 seconds.
- 3 At the next beep and pause, move to the bottom front teeth for the final 15 seconds of brushing.

#### Massage mode

Mode for gentle gum stimulation.

## Brushing routines

---

### Go Care routine

1-minute brushing cycle in the Clean mode, for a quick clean. You hear the Quadpacer signal at 15-second intervals.

### Max Care routine

3-minute brushing cycle that combines the Clean and Massage modes in one routine for a thorough mouth clean. There are 30 seconds of Clean mode and 15 seconds of Massage mode for each of the 4 sections of your mouth. You will hear the Quadpacer signal at 45 second intervals.

**Note: When Sonicare is used in clinical studies, the default 2-minute Clean mode must be selected. The handle should be fully charged. Deactivate the Easy-start feature and for areas where excess staining occurs, an additional 30 seconds of brushing time can be spent to assist with stain removal.**

## Features

### Easy-start

---

- This Sonicare model comes with the Easy-start feature activated.
- The Easy-start feature gently increases power over the first 14 brushings to let you get used to brushing with the Sonicare.

*Note: Each of the first 14 brushings must be at least 1 minute in length to properly move through the Easy-start ramp-up cycle properly.*

### Deactivating or activating the Easy-start feature

- 1** Attach your brush head to the handle.
- 2** Place the handle in the plugged-in charger.
  - To deactivate Easy-start:  
Press and hold the power on/off button for 5 seconds. You hear 1 beep to indicate that the Easy-start feature has been deactivated.
  - To activate Easy-start:  
Press and hold the power on/off button for 5 seconds. You hear 2 beeps to indicate that the Easy-start feature has been activated.

*Note: Using the Easy-start feature beyond the initial ramp-up period is not recommended and reduces Sonicare's effectiveness in removing plaque.*

### **Smartimer**

---

The Smartimer indicates that your brushing cycle is complete by automatically switching off the toothbrush at the end of the brushing cycle. Dental professionals recommend brushing no less than 2 minutes twice a day.

### **Quadpacer**

---

- The Quadpacer is an interval timer that has a short beep and pause to remind you to brush the 4 sections of your mouth. Depending on the brushing mode or brushing routine you have selected, the Quadpacer beeps at different intervals during the brushing cycle. See the 'Brushing instructions' and 'Personalising your brushing experience' sections (Fig. 5).

The Quadpacer feature on this model has been activated. To deactivate or reactivate the Quadpacer feature:

- 1** Place the handle with a brush head attached in the plugged-in charger.
- 2** Press and hold the personalised brushing button until:
  - A You hear 1 beep to indicate the Quadpacer has been deactivated.
  - B You hear 2 beeps to indicate the Quadpacer has been activated.

*Note: The Quadpacer feature cannot be deactivated on the pre-programmed Go Care and Max Care brushing routines.*

### **Integrated cord wrap of UV sanitiser**

---

- 1** If your Sonicare type includes the sanitiser, you can use the cord wrap feature in the bottom of the sanitiser to store excess cord (Fig. 7).

## Sanitising

### Sanitising

- With the UV sanitiser you can clean your brush head after every use (Fig. 8).

#### **Unplug the sanitiser, discontinue use and call Consumer Care if:**

- The UV light bulb remains on when the door is open.
- The window is broken or missing from the sanitiser.
- The sanitiser gives off smoke or a burning smell while it is operating.

*Note: UV light can be harmful to the human eye and skin. Keep the sanitiser out of reach of children at all times.*

- 1** After brushing, rinse the brush head and shake off excess water.

Do not place the travel cap on the brush head during sanitisation.

- 2** Place your finger in the door recess on top of the sanitiser and pull open the door (Fig. 9).

- 3** Place the brush head on one of the 2 pegs in the sanitiser with the bristles facing the light bulb. (Fig. 10)

*Note: Only clean Sonicare snap-on brush heads in the sanitiser.*

*Note: The Sonicare for Kids brush head wears and degrades in the UV sanitiser. You have to check and replace your child's brush more frequently than regular Sonicare brush heads (See chapter 'Replacement').*

- 4** Make sure the sanitiser is plugged into a wall socket.

- 5** Close the door and press the green power on/off button once to select the UV clean cycle.

*Note: You can only switch on the sanitiser if the door is properly closed.*

*Note: The sanitiser stops if you open the door during the sanitising cycle.*

*Note: The sanitiser cycle runs for 10 minutes and then automatically switches off.*

- ▶ The sanitiser is in operation when the light glows through the window.

## Cleaning

Do not clean the brush head, handle, travel charger or the UV sanitiser in the dishwasher.

### Toothbrush handle

- 1 Remove the brush head and rinse the metal shaft area with warm water (Fig. 11).

Do not push on the rubber seal on the metal shaft with sharp objects, as this may cause damage.

- 2 Use a damp cloth to wipe the entire surface of the handle.

### Brush head

- 1 Rinse the brush head and bristles after each use (Fig. 12).
- 2 Remove the brush head from the handle and rinse the brush head connection at least once a week with warm water.

### Charger

- 1 Unplug the charger.
- 2 Use a damp cloth to wipe the surface of the charger.

### Sanitiser

Do not immerse the sanitiser in water or rinse it under the faucet.

Do not clean the sanitiser when the UV light bulb is hot.

For optimal effectiveness, it is recommended to clean your sanitiser weekly.

- 1 Unplug the sanitiser.
- 2 Pull the drip tray straight out. Rinse the drip tray and wipe it clean with a damp cloth (Fig. 13).
- 3 Clean all interior surfaces with a damp cloth.

**4 Remove the protective screen in front of the UV light bulb. (Fig. 14)**

To remove the screen, grasp the edges next to the snaps, gently squeeze and pull out the protective screen.

**5 Remove the UV light bulb.**

To remove the light bulb, grasp it and pull it out of the metal clasp.

**6 Clean the protective screen and the UV light bulb with a damp cloth.****7 Reinsert the UV light bulb.**

To reinsert the light bulb, align the bottom of the light bulb with the metal clasp and push the bulb into the clasp.

**8 Reinsert the protective screen.**

To reinsert the screen, align the pegs on the screen with the slots on the reflective surface near the UV light bulb. Push the screen straight into the slots on the sanitiser.

## Storage

If you are not going to use the appliance for an extended period of time, remove the mains plug from the wall socket. Then clean the appliance and store it in a cool and dry place away from direct sunlight. You can use the cord wrap feature on the UV sanitiser to store the mains cord neatly.

## Replacement

### Brush head

---

Replace Sonicare brush heads every 3 months to achieve optimal results. Use only genuine Philips Sonicare replacement brush heads.

### UV light bulb

---

Replacement UV bulbs are available through the Philips Consumer Care Centre in your country.

## Recycling

- This symbol on a product means that the product is covered by European Directive 2012/19/EU (Fig. 15).
- This symbol means that the product contains a built-in rechargeable battery covered by European Directive 2006/66/EC which cannot be disposed of with normal household waste. Follow the instructions in section 'Removing the rechargeable battery' to remove the battery (Fig. 16).
- Inform yourself about the local separate collection system for electrical and electronic products and rechargeable batteries. Follow local rules and never dispose of the product and rechargeable batteries with normal household waste. Correct disposal of old products and rechargeable batteries helps prevent negative consequences for the environment and human health.

### Removing the rechargeable battery

---

**Only remove the rechargeable battery when you discard the appliance.  
Make sure the battery is completely empty when you remove it.**

To remove the rechargeable battery, you need a flat-head (standard) screwdriver. Observe basic safety precautions when you follow the procedure outlined below. Be sure to protect your eyes, hands, fingers, and the surface on which you work.

- 1** To deplete the rechargeable battery of any charge, remove the handle from the charger, switch on the Sonicare and let it run until it stops. Repeat this step until you can no longer switch on the Sonicare.
- 2** Insert a screwdriver into the slot located at the bottom of the handle and switch counterclockwise to release the bottom cap (Fig. 17).
- 3** Hold the handle upside down and push down on the shaft to release the internal components of the handle (Fig. 18).

- 4** Insert the screwdriver under the circuit board, next to the battery connections, and twist to break the connections. Remove the circuit board and pry the battery from the plastic carrier (Fig. 19).

The rechargeable battery can now be recycled and the rest of the product discarded appropriately.

## Guarantee and support

If you need information or support, please visit [www.philips.com/support](http://www.philips.com/support) or read the separate worldwide guarantee leaflet.

### Guarantee restrictions

---

The terms of the international guarantee do not cover the following:

- Brush heads
- Damage caused by use of unauthorised replacement parts
- Damage caused by misuse, abuse, neglect, alterations or unauthorised repair
- Normal wear and tear, including chips, scratches, abrasions, discolouration or fading
- UV light bulb

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©2015 Koninklijke Philips N.V. (KPNV). All rights reserved.
Philips and the Philips shield are trademarks of KPNV.
Sonicare and the Sonicare logo are trademarks
of Philips Oral Healthcare, LLC and/or KPNV.