# Always here to help you

Register your product and get support at  
[www.philips.com/welcome](http://www.philips.com/welcome)


Question?  
Contact  
Philips

**FlexCare+**  
**900+ series**  
**Rechargeable**  
**sonic toothbrush**

**PHILIPS**  
**sonicare**
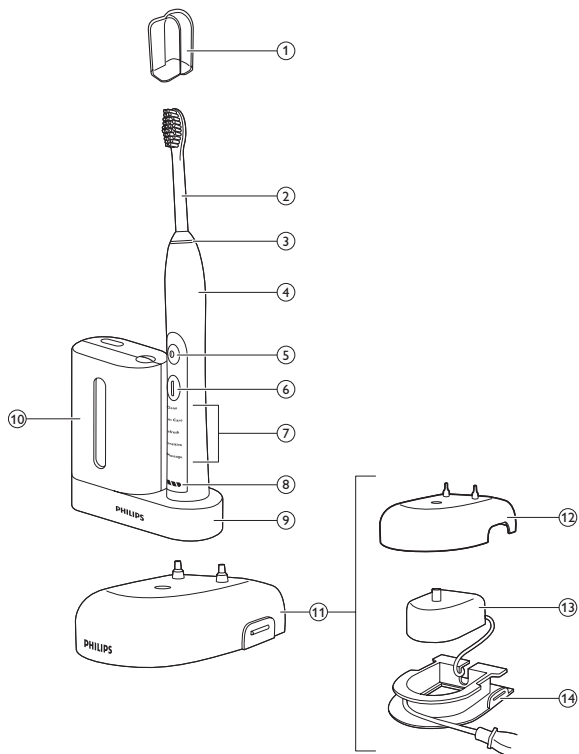
<b>ENGLISH</b>	<b>6</b>
<b>DEUTSCH</b>	<b>22</b>
<b>ΕΛΛΗΝΙΚΑ</b>	<b>40</b>
<b>ESPAÑOL</b>	<b>60</b>
<b>FRANÇAIS</b>	<b>77</b>
<b>ITALIANO</b>	<b>94</b>
<b>NEDERLANDS</b>	<b>110</b>
<b>PORTUGUÊS</b>	<b>126</b>

## Introduction

Congratulations on your purchase and welcome to Philips! To fully benefit from the support that Philips offers, register your product at [www.philips.com/welcome](http://www.philips.com/welcome).

## Important

Read this user manual carefully before you use the appliance and save it for future reference.

### **Danger**

- Keep the charger and/or sanitiser away from water. Do not place or store it over or near water contained in a bathtub, washbasin, sink etc. Do not immerse the charger and/or sanitiser in water or any other liquid. After cleaning, make sure the charger and/or sanitiser is completely dry before you connect it to the mains.

### **Warning**

- The mains cord cannot be replaced. If the mains cord is damaged, discard the charger and/or sanitiser.
- Always have the charger and/or sanitiser replaced with one of the original type in order to avoid a hazard.
- Do not use the charger and/or sanitiser outdoors or near heated surfaces.

- If the appliance is damaged in any way (brush head, toothbrush handle, charger and/or sanitiser), stop using it. This appliance contains no user-serviceable parts. If the appliance is damaged, contact the Consumer Care Centre in your country (see chapter 'Guarantee and support').
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children unless they are older than 8 and supervised.
- Discontinue use of the sanitiser if the UV light bulb remains on when the door is open. UV light can be harmful to the human eye and skin. Keep the sanitiser out of the reach of children at all times.

### **Caution**

- Do not clean the brush head, the handle, the charger, the charger cover and/or the sanitiser in the dishwasher.
- If you have had oral or gum surgery in the previous 2 months, consult your dentist before you use the toothbrush.

## 8 ENGLISH

- Consult your dentist if excessive bleeding occurs after using this toothbrush or if bleeding continues to occur after 1 week of use.
- The Sonicare toothbrush complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician or the manufacturer of the implanted device prior to use.
- If you have medical concerns, consult your doctor before you use the Sonicare.
- This appliance has only been designed for cleaning teeth, gums and tongue. Do not use it for any other purpose. Stop using the appliance and contact your doctor if you experience any discomfort or pain.
- Do not use other brush heads than the ones recommended by the manufacturer.
- If your toothpaste contains peroxide, baking soda or bicarbonate (common in whitening toothpastes), thoroughly clean the brush head with soap and water after each use. This prevents possible cracking of the plastic.
- The UV light bulb is hot during and immediately after the sanitising cycle. Do not touch the UV light bulb when it is hot.
- Do not operate the sanitiser without the protective screen in place to avoid contact with a hot bulb.

## Electromagnetic fields (EMF)

---

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

### General description (Fig. 1)

- 1 Hygienic travel cap
- 2 Brush head
- 3 Interchangeable colour code ring
- 4 Handle with soft grip
- 5 Power on/off button
- 6 Mode button
- 7 Brushing modes
- 8 Deluxe recharge gauge
- 9 UV sanitiser base with integrated charger and cord wrap (specific types only)
- 10 UV sanitiser (specific types only)
- 11 Deluxe charger (specific types only)
- 12 Charger cover with brush head holder
- 13 Travel charger
- 14 Charger base with cord wrap
  - Not shown: UV light bulb
  - Not shown: Sanitiser drip tray
  - Not shown: Protective screen for UV light bulb

## Preparing for use

### Changing the colour code ring

Sonicare brush heads come with interchangeable colour code rings to identify your brush head.

To change the colour code ring:

- 1** Pull the colour code ring from the bottom of the brush head (Fig. 2).
- 2** Install a new ring by slipping one edge of the new ring over the bottom of the brush head (Fig. 3).

### Attaching the brush head

- 1** Align the brush head so the bristles face the front of the handle. (Fig. 4)
- 2** Firmly press the brush head down on the metal shaft until it stops.

*Note: There is a small gap between the colour code ring and the handle.*

### Charging your Sonicare

- 1** Put the mains plug of the charger or sanitiser in the wall socket.
- 2** Place the handle on the charger or sanitiser (Fig. 5).
  - ▶ The flashing light of the battery gauge indicates that the toothbrush is charging.
    - 1 solid green LED indicates that the toothbrush is 34-66% charged.
    - 2 solid green LEDs indicate that the toothbrush is 67-94% charged.
    - 3 solid green LEDs indicate that the toothbrush is fully charged.

*Note: If the battery charge of your Sonicare is low, you hear 3 beeps and 1 LED of the recharge gauge flashes yellow for 30 seconds after completing the brushing cycle.*

*Note: To keep the battery fully charged at all times, you may keep your Sonicare on the charger when not in use. Once the Sonicare is fully charged, the charger draws minimal energy from the wall socket.*

*Note: It takes at least 24 hours to fully charge the battery.*

## Using the Sonicare

### Brushing instructions

- 1** Wet the bristles and apply a small amount of toothpaste.
- 2** Place the toothbrush bristles against the teeth at a slight angle towards the gumline (Fig. 6).
- 3** Press the power on/off button to switch on the Sonicare.
- 4** Apply light pressure to maximise Sonicare's effectiveness and let the Sonicare toothbrush do the brushing for you. Do not scrub.
- 5** Gently move the brush head slowly across the teeth in a small back and forth motion so the longer bristles reach between your teeth. Continue this motion throughout your brushing cycle.

*Note: To make sure you brush evenly throughout the mouth, divide your mouth into 4 sections using the Quadpacer feature (see chapter 'Features').*

- 6** Begin brushing in section 1 (outside top teeth) and brush for 30 seconds before moving to section 2 (inside top teeth). Continue brushing in section 3 (outside bottom teeth) and brush for 30 seconds before moving to section 4 (inside bottom teeth) (Fig. 7).
- 7** After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs. You may also brush your tongue, with the toothbrush switched on or off, as you prefer (Fig. 8).

Your Sonicare is safe to use on:

- Braces (brush heads wear out sooner when used on braces)
- Dental restorations (fillings, crowns, veneers)

### Brushing modes

---

Sonicare automatically starts in the default Clean mode.

- 1 Prior to switching on the Sonicare, press the mode button to toggle between modes.
  - ▶ The selected brushing mode illuminates.
  - ▶ Once a brushing mode has been selected, it remains active until another brushing mode is selected.

*Note: When the toothbrush is switched on, you cannot toggle between modes. The following brushing modes are available:*

#### **Clean**

Complete whole mouth cleaning (2 minutes).

#### **Gum Care**

Complete whole mouth cleaning (2 minutes) plus gentle cleaning for problem areas and along the gumline (1 minute).

#### **Refresh**

Touch-up for a quick clean (1 minute).

#### **Sensitive**

Gentle cleaning for sensitive gums and teeth (2 minutes).

#### **Massage**

Gentle gum stimulation (2 minutes).

*Note: When Sonicare is used in clinical studies, the default 2-minute Clean mode must be selected. The handle has to be fully charged. Deactivate the Easy-start feature. For areas where excess staining occurs, an additional 30 seconds of brushing time can be spent to help remove stains.*

## Features

### Easy-start

---

This Sonicare type comes with the Easy-start feature activated. The Easy-start feature gently increases power over the first 14 brushings to help you get used to brushing with the Sonicare.

*Note: Each of the first 14 brushings must be at least 1 minute in length to properly move through the Easy-start ramp-up cycle.*

### Deactivating or activating the Easy-start feature

- 1** Attach your brush head to the handle.
- 2** Place the handle on the plugged-in charger.
  - To deactivate Easy-start:  
Press and hold the power on/off button for 2 seconds. You hear 1 beep to indicate that the Easy-start feature has been deactivated.
  - To activate Easy-start:  
Press and hold the power on/off button for 2 seconds. You hear 2 beeps to indicate that the Easy-start feature has been activated.

*Note: Using the Easy-start feature beyond the initial ramp-up period is not recommended and reduces Sonicare's effectiveness in removing plaque.*

### Smartimer

---

The Smartimer indicates that your brushing cycle is complete when it automatically switches off the toothbrush at the end of the brushing cycle. Dentists recommend brushing at least 2 minutes twice a day.

### Quadpacer

---

- The Quadpacer is an interval timer that has a short beep and pause to remind you to brush the 4 sections of your mouth. Depending on the brushing mode you have selected, the Quadpacer beeps at different intervals during the brushing cycle. See section 'Brushing modes' in chapter 'Using the Sonicare' (Fig. 7).
- The Quadpacer feature on this type has been activated. To deactivate or reactivate the Quadpacer feature:

- 1** Place the handle with a brush head attached on the plugged-in charger.
- 2** Press and hold the mode button for 2 seconds.
  - ▶ You hear 1 beep to indicate the Quadpacer has been deactivated or you hear 2 beeps to indicate the Quadpacer has been activated.

### **Cord wrap on deluxe charger (specific types only)**

---

If your type includes the deluxe charger cover and base, the travel charger is pre-installed inside the cover. If you want to shorten the mains cord, you can store excess cord in the cord wrap feature built into the charger base.

- 1** To separate the charger cover from the charger base, press the two grey snap bars on the charger base and pull the white charger cover upwards. (Fig. 9)
- 2** Wrap excess power cord around the grey charger base as shown in the picture. Be sure to wrap the cord on the inside of the two snap levers. (Fig. 10)
- 3** When you have stored the excess cord, guide the mains cord through the small groove in the rear of the grey charger base. (Fig. 11)
- 4** To reattach the charger cover, press it down over the charger base until it snaps into place.

*Tip: For extra convenience during travel, you can remove the travel charger and use it without the charger cover and charger base.*

### **Cord wrap on sanitiser with integrated charger (specific types only)**

---

- If your type includes the sanitiser, you can store excess cord in the cord wrap feature built into the bottom of the sanitiser. (Fig. 12)

## Sanitising (specific types only)

- With the UV sanitiser you can clean your brush head after every use (Fig. 13).

Discontinue use of the sanitiser if the UV light bulb remains on when the door is open or if the blue window is broken or missing from the sanitiser. UV light can be harmful to the human eye and skin.

This appliance should be kept out of reach of children at all times.

Unplug the sanitiser and call Consumer Care if the sanitiser gives off smoke or a burning smell while it is operating.

- 1** After brushing, rinse the brush head and shake off excess water.

Do not place the travel cap on the brush head during sanitisation.

- 2** Place your finger in the door recess on top of the sanitiser and pull open the door (Fig. 14).

- 3** Place the brush head on one of the 2 pegs in the sanitiser. (Fig. 15)

- Make sure the bristles of the brush head directly face the light bulb.

*Note: Only clean Sonicare ProResults brush heads in the sanitiser.*

- 4** Make sure the sanitiser is plugged into the wall socket.

- 5** Close the door and press the green power on/off button once to select the UV clean cycle.

*Note: You can only switch on the sanitiser if the door is properly closed.*

*Note: The sanitiser stops if you open the door during the sanitising cycle.*

*Note: The sanitiser cycle runs for 10 minutes and then automatically switches off.*

- ▶ The sanitiser is in operation when the blue light glows through the window.
- ▶ When the sanitising cycle is complete, the sanitiser automatically switches off.

## Cleaning

Do not clean the brush head, handle, travel charger, charger cover and the UV sanitiser in the dishwasher.

### Toothbrush handle

---

- 1 Remove the brush head and rinse the metal shaft area with warm water (Fig. 16).

Do not push on the rubber seal on the metal shaft with sharp objects, as this may cause damage.

- 2 Use a damp cloth to wipe the entire surface of the handle.

### Brush head

---

- 1 Rinse the brush head and bristles after each use (Fig. 17).
- 2 Remove the brush head from the handle and rinse the brush head connection at least once a week with warm water.

### Charger

---

- 1 Unplug the charger.
- 2 Use a damp cloth to wipe the surface of the charger.

### Sanitiser (specific types only)

---

Do not clean the sanitiser when the UV light bulb is hot.

For optimal effectiveness, it is recommended to clean the sanitiser weekly.

- 1 Unplug the sanitiser.
- 2 Pull the drip tray straight out. Rinse the drip tray and wipe it clean with a damp cloth (Fig. 18).
- 3 Clean all reflector surfaces with a damp cloth.
- 4 Remove the protective screen in front of the UV light bulb. (Fig. 19)

To remove the screen, grasp the edges next to the snaps, gently squeeze and pull it out.

**5 Remove the UV light bulb.**

To remove the light bulb, grasp it and pull it out of the metal clasp.

**6 Clean the protective screen and the UV light bulb with a damp cloth.****7 Reinsert the UV light bulb.**

To reinsert the light bulb, align the bottom of the light bulb with the metal clasp and push the bulb into the clasp.

**8 Reinsert the protective screen.**

To reinsert the screen, align the pegs on the screen with the slots on the reflective surface near the UV light bulb. Push the screen straight into the slots on the sanitiser.

## Storage

- If you are not going to use the appliance for an extended period of time, remove the mains plug from the wall socket. Then clean the appliance and store it in a cool and dry place away from direct sunlight. You can use the cord wrap feature to store the mains cord neatly.

## Replacement

### Brush head

---

- Replace Sonicare brush heads every 3 months to achieve optimal results.
- Use only compatible Sonicare replacement brush heads.

### UV light bulb

---

- You can order replacement UV light bulbs from the Consumer Care Centre in your country or an authorised Philips service centre.

## Recycling

- This symbol on a product means that the product is covered by European Directive 2012/19/EU (Fig. 20).
- This symbol means that the product contains a built-in rechargeable battery covered by European Directive 2006/66/EC which cannot be disposed of with normal household waste. Follow the instructions in section 'Removing the rechargeable battery' to remove the battery (Fig. 21).
- Inform yourself about the local separate collection system for electrical and electronic products and rechargeable batteries. Follow local rules and never dispose of the product and rechargeable batteries with normal household waste. Correct disposal of old products and rechargeable batteries helps prevent negative consequences for the environment and human health.

## Removing the rechargeable battery

---

Please note that this process is not reversible.

To remove the rechargeable battery, you need a flat-head (standard) screwdriver:

- 1** To deplete the rechargeable battery of any charge, remove the handle from the charger, switch on the Sonicare and let it run until it stops. Repeat this step until you can no longer switch on the Sonicare.
- 2** Insert a flat-head (standard) screwdriver into the slot located at the bottom of the handle. Turn the screwdriver anticlockwise to release the bottom cap. (Fig. 22)
- 3** Hold the handle upside down and push down on the shaft to release the internal components of the handle (Fig. 23).
- 4** Insert the screwdriver under the circuit board, next to the battery connections, and twist to break the connections. Remove the circuit board and pry the battery from the plastic carrier (Fig. 24).

## Guarantee and support

If you need information or support, please visit [www.philips.com/support](http://www.philips.com/support) or read the separate worldwide guarantee leaflet.

### Guarantee restrictions

The terms of the international guarantee do not cover the following:

- Brush heads.
- Damage caused by use of unauthorised replacement parts.
- Damage caused by misuse, abuse, neglect, alterations or unauthorised repair.
- Normal wear and tear, including chips, scratches, abrasions, discolouration or fading.

## Frequently asked questions

This chapter lists the questions most frequently asked about the appliance. If you cannot find the answer to your question, contact the Consumer Care Centre in your country.

Question	Answer
Why doesn't the Sonicare toothbrush work?	You may need to recharge the Sonicare. Charge the Sonicare for at least 24 hours.
	The socket may not be live. The electricity supply to sockets in mirror cabinets over washbasins or shaver sockets in bathrooms may be cut off when the light is switched off.
	You may need to reattach or replace the brush head.

Question	Answer
Why do I feel a tickling sensation when I use the Sonicare?	You may experience a slight tickling/tingling sensation when you use Sonicare for the first time. As you become accustomed to brushing with the Sonicare, this sensation diminishes.
Why is there a brief pause in the bristle motion while I'm brushing?	This is the Quadpacer feature reminding you to move to a new section of your mouth. This feature enables you to spend equal time brushing all areas of the mouth (see chapter 'Using the Sonicare, section 'Brushing instructions').
The brush head chatters against my teeth. What can I do about this?	To avoid this, keep the brush head properly oriented, with the bristles placed on the gumline at a slight angle.
Why does the toothbrush seem less powerful?	You may need to charge the Sonicare.
	The Easy-start feature may be on. Deactivate the Easy-start feature (see chapter 'Features'). You may need to replace the brush head.
Why doesn't the mode button work during brushing?	Brushing modes should be selected prior to switching on the Sonicare.

Question	Answer
Why doesn't the UV sanitiser light go on?	You may need to plug the sanitiser into a live wall socket.
	You may need to shut the sanitiser door. You may need to replace the UV light bulb (see chapter 'Replacement').
Do I need a special charger if I travel with my Sonicare?	The travel charger, the deluxe charger and the sanitiser can be used at common worldwide voltages between 100-240V AC, at either 50 or 60Hz. All you need is a simple plug adapter.
	Travel tip: The Sonicare battery is designed to last up to 3 weeks between charges.

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©2015 Koninklijke Philips N.V. (KPNV). All rights reserved.
Philips and the Philips shield are trademarks of KPNV.
Sonicare and the Sonicare logo are trademarks
of Philips Oral Healthcare, Inc. and/or KPNV.