Philips Sonicare DiamondClean





ENGLISH 4

4 ENGLISH

1	Important Safeguards	5
2	Your Philips Sonicare	8
3	Charging your Philips Sonicare	9
4	Using your Philips Sonicare	10
5	Getting used to your Philips Sonicare	12
6	Traveling with your Philips Sonicare	13
7	Cleaning your Philips Sonicare	14
8	Storing your Philips Sonicare	15
9	Disposal	16
10	Warranty and support	18

1 IMPORTANT SAFEGUARDS

READ ALL INSTRUCTIONS BEFORE USE

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

Dangers

To reduce the risk of electrocution:

- Always place and store the charging base or charging travel case at a safe distance from water.
- Unplug the charging base or charging travel case immediately if it falls into water.
- Never use a charging base or charging travel case with a damaged cord or plug.

Warnings

To reduce the risk of burns, electrocution, fire, or physical injury:

- Do not use attachments other than those recommended by the manufacturer.
- Do not plug the charging base or charging travel case into an outlet with a voltage other than specified on the bottom of the charging base or charging travel case. This product is designed to operate within a range of 100 to 240 volts. Voltage converters do not warrant voltage compatibility.
- Never force the plug into an outlet; if the plug does not easily fit into the outlet, discontinue use.

- Discontinue use if the product appears damaged in any way (brush head, handle, charging base or charging travel case).
- Do not use the charging base or charging travel case if it has dropped into water.
- Keep the cord away from heated surfaces.
- Do not use the charging base or charging travel case outdoors.
- This product contains no user-serviceable parts. Refer to 'Warranty and support' if the appliance no longer works properly or needs repair.
- Use this appliance only for its intended use as described in this document or as recommended by your dental professional.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory, or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- Philips Sonicare is a personal care device and is not intended for use on multiple patients in a dental office or institution.
- Use UL-listed computer, USB hub, or USB duplex receptacle when you charge the toothbrush handle in the charging travel case.
- When you use a computer to recharge the toothbrush handle, the computer should be located away from sources of water (i.e., sink or tub) so as to avoid electrical shock.


Medical Warnings

- Consult your dental professional before you use this product if you have had oral or gum surgery in the previous 2 months.
- Contact your dental professional if excessive bleeding occurs after using this product or bleeding continues to occur after 1 week of use.
- Discontinue use of this product and contact a physician/ dentist if discomfort or pain is experienced.
- This product complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician or the device manufacturer prior to use.
- Consult your physician if you have other medical concerns.

Electromagnetic fields (EMF)

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

SAVE THESE INSTRUCTIONS

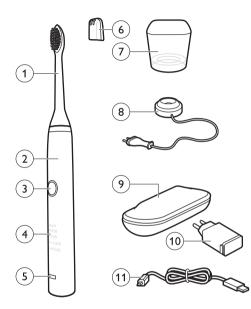

1 Brush head

2

- 2 Handle with ergonomic shape
- 3 Recessed power on/off button
- 4 Hidden brushing mode display
- 5 Charge indicator
- 6 Hygienic travel cap


- 7 Charging glass
- 8 Charging base
- 9 Charging travel case
- 10 USB wall adapter
- 11 USB cord

The contents of the box may vary based on the model purchased. See outside of box for contents included.

CHARGING YOUR PHILIPS SONICARE


3

Charge for 24 hours before first use. The charge indicator flashes green until the Philips Sonicare is fully charged. Rapid yellow flashing indicates that the battery is low. Your Philips Sonicare has been designed to last up to 3 weeks on a single charge. Read 'Traveling with your Philips Sonicare' for additional charging methods. The charging glass is designed to properly charge only one handle at a time. Do not place metallic objects in the charging glass, as this may affect the charging process and may cause the charging base to overheat.

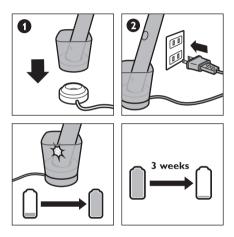

USING YOUR PHILIPS SONICARE


Brushing	 Move the brush head slowly and gently across the teeth in a small back and forth motion so the longer bristles reach between your teeth. Do not scrub like a manual toothbrush. Apply light pressure to maximize Philips Sonicare's effectiveness and let the Philips Sonicare do the brushing for you. The Quadpacer is an interval timer that beeps and briefly pauses to remind you to brush the 4 sections of your mouth. Depending on the brushing mode you have selected, the Quadpacer beeps at different intervals during the brushing cycle.
Clean	- The ultimate in plague removal (standard mode).
White	- Whitening mode of 2 minutes to remove surface stains, plus an additional 30 seconds to brighten and polish your front teeth.
Polish	- Quick 1-minute mode to brighten and polish your front teeth.
Gum care	 Complete mode that combines the Clean mode with an additional minute to gently stimulate and massage your gums to improve gum health.
Sensitive	 Extra gentle mode for sensitive teeth and gums.

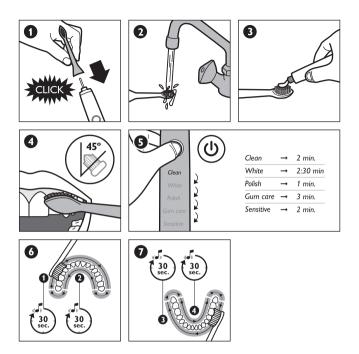

5 GETTING USED TO YOUR PHILIPS SONICARE


Easy-start - The Philips Sonicare may feel ticklish upon first use. To help you ease into the Philips Sonicare experience, the Easy-start feature gradually increases the sonic power over the first 14 brushings. Brushing must last at least 1 minute to properly advance through this ramp-up cycle.

Easy-start may be deactivated - Place the handle in the plugged-in charging base or charging travel case. Push and hold the on/off button for 2 seconds; 1 beep indicates the Easy-start feature has been deactivated; 2 beeps indicate the feature has been reactivated. Clinical studies must be performed with the Philips Sonicare fully charged and the Easy-start feature deactivated.

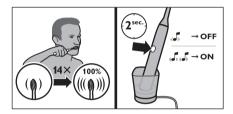


6 TRAVELING WITH YOUR PHILIPS SONICARE

Charging via wall socket - Place the handle in the charging travel case, plug the USB cord into the USB wall adapter, and plug the USB wall adapter into a wall socket. The charge indicator on the handle flashes green until the Philips Sonicare is fully charged. **Charging via USB** - Remove the USB cord from the USB wall adapter and insert the USB plug into any USB port. Only use the USB wall adapter and USB cord provided.

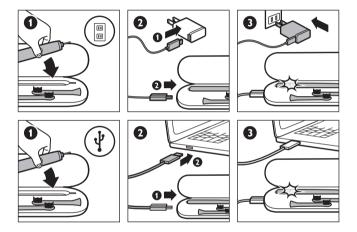

7 CLEANING YOUR PHILIPS SONICARE

Cleaning - Clean the Philips Sonicare periodically. Rinse the handle and the brush head. Wipe down all surfaces with a soft cloth and mild soap, if desired. Do not use sharp objects to clean. Only the charging glass is top rack dishwasher-safe. Do not clean the brush head, handle, charging base or charging travel case in the dishwasher.

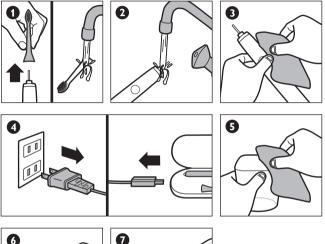





8

STORING YOUR PHILIPS SONICARE

Storing - If you are not going to use the Philips Sonicare for an extended period of time, unplug the charging base, clean all parts and store them in a cool, dry place away from direct sunlight.



9 DISPOSAL



- This appliance contains a rechargeable battery which must be disposed of properly.
- Contact your local town or city officials for battery disposal information. You can also call 1-800-8-BATTERY or visit www.rbrc.com for battery drop-off locations.
- For assistance, visit our website **www.philips.com/support** or call 1-800-243-3050 toll free.


REMOVING THE RECHARGEABLE BATTERY

This process is not reversible. To remove the rechargeable battery, you need a flat-head (standard) screwdriver. Observe basic safety precautions when you follow the procedure below. Be sure to protect your eyes, hands, fingers and the surface on which you work.








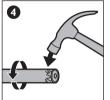



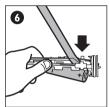



10

WARRANTY AND SUPPORT

LOCATING THE MODEL NUMBER

Look on the bottom of the unit for the model number (HX9383, HX9382, HX9362, HX9352, HX9343, HX9333 or HX9332).

Philips warrants its products for two years after the date of purchase. Defects due to faulty materials and workmanship will be repaired or replaced at Philips' expense provided that convincing proof of purchase in the qualifying period is provided. Philips recommends use of genuine Sonicare brush heads for optimal performance. Use of other brush heads may damage your Sonicare toothbrush and limit your warranty protection.

Contact our Consumer Care Center at 1-800-682-7664 (North America), outside North America contact your local Philips Consumer Care Center: Internet information: **www.sonicare.com** (North America) or **www.philips.com** (outside North America).

WARRANTY EXCLUSIONS

What is not covered under warranty:

- Brush heads.
- · Damage caused by use of unauthorized replacement parts or unauthorized brush heads.
- · Damage caused by use of unauthorized replacement parts.
- · Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.

IMPLIED WARRANTIES

ANY IMPLIED WARRANTIES, INCLUDING IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED TO THE DURATION OF THE EXPRESS WARRANTIES SET FORTH ABOVE. IN SOME STATES LIMITATIONS ON DURATION OF IMPLIED WARRANTIES DO NOT APPLY.

LIMITATION OF REMEDIES


IN NO EVENT SHALL PHILIPS OR ANY OF ITS AFFILIATED OR SUBSIDIARY COMPANIES BE LIABLE FOR ANY SPECIAL, INCIDENTAL OR CONSEQUENTIAL DAMAGES BASED UPON BREACH OF WARRANTY, BREACH OF CONTRACT, NEGLIGENCE, TORT, OR ANY OTHER LEGAL THEORY, SUCH DAMAGES INCLUDE, WITHOUT LIMITATION, LOSS OF SAVINGS OR REVENUE; LOSS OF PROFIT; LOSS OF USE; THE CLAIMS OF THIRD PARTIES INCLUDING, WITHOUT LIMITATION, DENTISTS AND DENTAL HYGIENISTS; AND COST OF ANY SUBSTITUTE EQUIPMENT OR SERVICES, SOME STATES DO NOT ALLOW THE EXCLUSION OR LIMITATION OF INCIDENTAL OR CONSEQUENTIAL DAMAGES.



www.philips.com/Sonicare

©2016 Koninklijke Philips N.V. (KPNV). All rights reserved. Philips and the Philips shield are trademarks of KPNV. Sonicare and the Sonicare logo are trademarks of Philips Oral Healthcare, LLC and/or KPNV.

Protected by U.S. and international patents. Other patents pending, Reg. U.S. Pat. & Tm. Off.

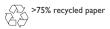

4235.020.4711.5