

6 Automatically turns off after two minutes.

7 Replace brush heads every 3 months or sooner if needed.

The toothbrush is safe to use on:

- Braces (brush heads wear out sooner when used on braces).
- Dental restorations (fillings, crowns, veneers).

CLEANING AND STORAGE

Clean head with water and wipe handle regularly with a damp cloth only. Store in a cool dry location when not used.

DISPOSAL

At product end of life - do not throw away in trash can; recycle products and batteries appropriately.

WARRANTY

This product is warranted for 1 year from the date of purchase. What is not covered under the warranty:

- Brush heads
- Replacement batteries
- Damage from misuse, abuse, neglect alterations, unauthorized repair/parts.
- Normal wear and tear including discoloration or fading.

For complete warranty details, visit www.philips.com/support and download the operator's manual with full warranty details.

Questions: Call 1-800-682-7664


Always here to help you

Register your product and get support at www.philips.com/welcome

Question?
Contact
Philips

HX3600

PHILIPS

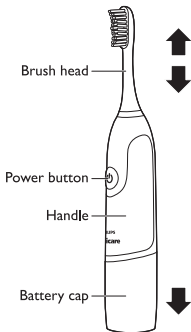

IMPOR- TANT SAFE- GUARDS

READ ALL INSTRUCTIONS
BEFORE USE

Warnings

- Use only as intended and with manufacturer attachment.
- Use only 1.5 V AA Alkaline batteries. Always replace with two new batteries; never mix old and new batteries.
- Supervise use by children and persons with reduced physical, sensory or mental capabilities.
- Discontinue use if any part of this product appears damaged in any way (brush head and handle).
- Use this product only for its intended use. This is not a toy.
- Consult dentist or physician if you experience significant bleeding, discomfort, pain or medical concerns.

SAVE THESE INSTRUC- TIONS

Electromagnetic fields (EMF)

This appliance complies with current standards regarding electromagnetic fields (EMF). If used as directed, this appliance is safe.

OPERATING INSTRUCTIONS

Note: Not intended for use on multiple patients in a dental practice or institution.

- 1** Pull battery cap down. Install 2 size 1.5 V AA Alkaline batteries. Proper battery direction is indicated under cap. Do not use the product without battery cap installed.
- 2** Install brush head by placing on metal shaft and pushing down. Remove brush head by pulling up.

Note: A gap between brush head and handle is normal.

- 3** Avoid splattering: place brush in mouth before turning on.
- 4** Press power on/off button to turn on.
- 5** Brush normally with light pressure for 2 minutes.