
FlexCare Platinum

PHILIPS
sonicare

1

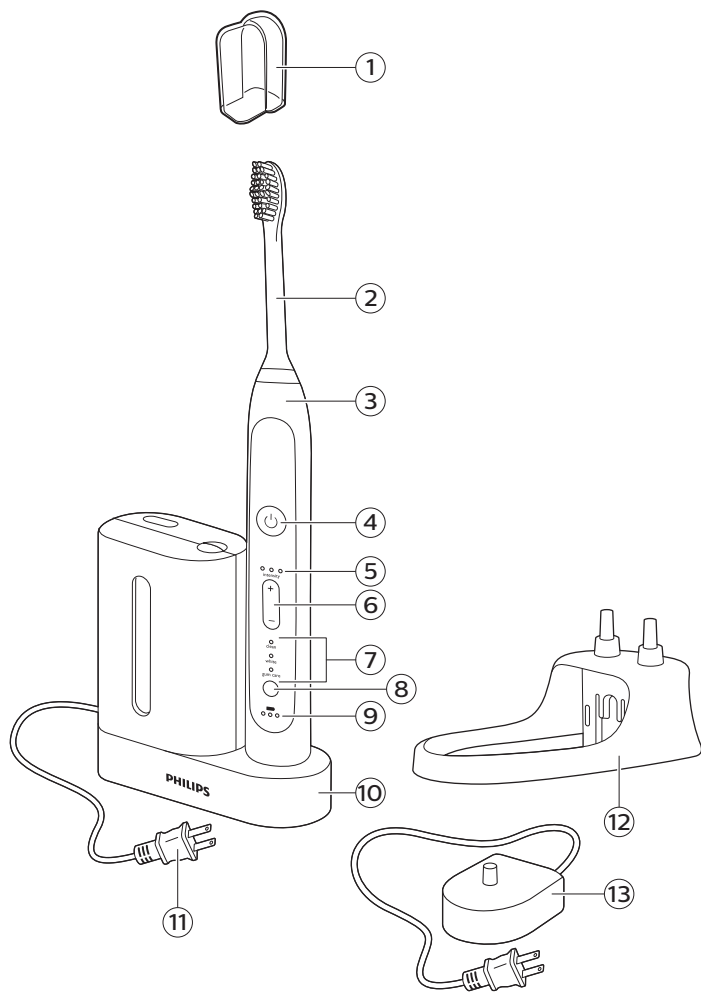— |

| —

— |

| —

English	6
Bahasa Indonesia	18
Bahasa Melayu	31
Français	44
Tiếng Việt	58
ภาษาไทย	69
العربية	79
فارسی	90

Introduction

Congratulations on your purchase and welcome to Philips! To fully benefit from the support that Philips offers, register your product at www.philips.com/welcome.

Important

Read this user manual carefully before you use the appliance and save it for future reference.

Danger

- Keep the charger away from water. Do not place or store it over or near water contained in a bathtub, washbasin, sink etc. Do not immerse the charger in water or any other liquid. After cleaning, make sure the charger is completely dry before you connect it to the mains.

Warning

- The mains cord cannot be replaced. If the mains cord is damaged, discard the charger.
- Always have the charger replaced with one of the original type in order to avoid a hazard.
- Do not use the charger outdoors or near heated surfaces.
- If the appliance is damaged in any way (brush head, toothbrush handle or charger), stop using it. This appliance contains no serviceable parts. If the appliance is damaged, contact the Consumer Care Centre in your country (see chapter 'Guarantee and support').
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.

Caution

- Do not clean the brush head, the handle, the charger, the charger cover and/or the sanitiser in the dishwasher.

- If you have had oral or gum surgery in the previous 2 months, consult your dentist before you use the toothbrush.
- Consult your dentist if excessive bleeding occurs after using this toothbrush or if bleeding continues to occur after 1 week of use. Also consult your dentist if you experience discomfort or pain when you use the Sonicare.
- The Sonicare toothbrush complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician or the device manufacturer of the implanted device prior to use.
- If you have medical concerns, consult your doctor before you use the Sonicare.
- This appliance has only been designed for cleaning teeth, gums and tongue. Do not use it for any other purpose. Stop using the appliance and contact your doctor if you experience any discomfort or pain.
- The Sonicare toothbrush is a personal care device and is not intended for use on multiple patients in a dental practice or institution.
- Stop using a brush head with crushed or bent bristles. Replace the brush head every 3 months or sooner if signs of wear appear.
- Do not use other brush heads than the ones recommended by the manufacturer.
- If your toothpaste contains peroxide, baking soda or bicarbonate (common in whitening toothpastes), thoroughly clean the brush head with soap and water after each use. This prevents possible cracking of the plastic.
- Do not insert objects into the air openings of the sanitiser or block these air openings by placing it on a soft surface or in a position or location where objects will block these openings (examples of objects:
 - magazine, tissue paper and hair).
- The UV light bulb is hot during and immediately after the sanitising cycle. Do not touch the UV light bulb when it is hot.
- Do not operate the sanitiser without the protective screen in place to avoid contact with a hot bulb.

8 English

- If the sanitiser has been used for 3 consecutive UV clean cycles (i.e. 3 times in a row), switch it off for at least 30 minutes before you start another UV clean cycle.
- Discontinue use of the sanitiser if the UV light bulb remains on when the door is open. UV light can be harmful to the human eye and skin. Keep the sanitiser out of the reach of children at all times.

Electromagnetic fields (EMF)

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

General description (Fig. 1)

- 1 Hygienic travel cap
- 2 Brush head with unique identifier
- 3 Handle
- 4 Power on/off button
- 5 Intensity level indicator
- 6 Intensity level button (+ and -)
- 7 Brushing modes
- 8 Mode button
- 9 Battery level indicator
- 10 UV sanitiser base with integrated charger and cord wrap (specific types only)
 - Not shown: UV light bulb
 - Not shown: Sanitizer drip tray
 - Not shown: Protective screen for UV light bulb
- 11 Mains cord with plug
- 12 Charger cover with brush head holders (select models)
- 13 Travel charger
 - Not shown: Travel case

The content of the box may vary based on the model purchased.

Preparing for use

Attaching the brush head

- 1 Align the brush head so the bristles point in the same direction as the front of the handle. (Fig.)
- 2 Firmly press the brush head down onto the metal shaft until it stops.

Philips Sonicare brush heads are imprinted with icons to easily identify your brush head.

Charging the appliance

- 1 Put the plug of the charger or sanitizer in a wall socket.
- 2 Place the handle on the charger.

The flashing light of the battery level indicator indicates that the toothbrush is charging. (Fig.)

Note: It takes at least 24 hours to fully charge the battery.

Using the appliance

Brushing instructions

- 1 Wet the bristles (Fig. 4) and apply a small amount of toothpaste (Fig. 5). (Fig. 4)
- 2 Place the toothbrush bristles against the teeth at a slight angle (45 degrees), pressing firmly to make the bristles reach the gumline or slightly beneath the gumline. A change in vibration (and a slight change in sound) alerts you when you apply too much pressure while brushing. (Fig. 6)

Note: Keep the center of the brush in contact with the teeth at all times.

- 3 Press the power on/off button to switch on the Philips Sonicare. (Fig. 7)
- 4 Apply light pressure to maximize Philips Sonicare's effectiveness and let the Philips Sonicare toothbrush do the brushing for you. Do not scrub.

10 English

- 5 Gently move the brush head slowly across the teeth in a small back and forth motion so the longer bristles reach between your teeth (Fig. 8). Continue this motion throughout the brushing cycle.
- 6 To clean the inside surfaces of the front teeth, tilt the brush handle semi-upright and make several vertical overlapping brushing strokes on each tooth (Fig. 9).
- 7 Begin brushing in section 1 (outside upper teeth) and brush for 30 seconds (45 seconds in Deep Clean mode) before you move to section 2 (inside upper teeth). Continue brushing in section 3 (outside lower teeth) and brush for 30 seconds (45 seconds in Deep Clean mode) before you move to section 4 (inside lower teeth (Fig. 10)).

Note: To make sure you brush evenly throughout the mouth, divide your mouth into 4 sections using the Quadpacer feature (see chapter 'Features').

Note: After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs. (Fig. 11) You can also brush your tongue, with the toothbrush turned on or off, as you prefer.

Your Philips Sonicare is safe to use on:

- Braces (brush heads wear out sooner when used on braces).
- Dental restorations (fillings, crowns, veneers).

Brushing modes

Clean mode

For an exceptional everyday clean (default mode).

White mode

2 minutes to remove surface stains and 30 seconds to brighten and polish your front teeth.

Deep Clean

Provides an invigorating 3 minutes of alternating cleaning and massaging to deliver an exceptionally clean experience.

Your Philips Sonicare comes with the Clean mode activated.

Note: When the toothbrush is switched on, you cannot toggle between modes.

- 1 Before you switch on the Philips Sonicare, press the mode button to toggle between modes.
 - The selected brushing mode lights up.
 - Once a brushing mode has been selected, it remains active until another brushing mode is selected.

Note: When the Philips Sonicare is used in clinical studies, it must be used in the default 2-minute Clean mode at high intensity or in the Deep Clean mode at high intensity.

Intensity levels per mode

- Low: 1 LED indicator light.
- Medium: 2 LED indicator lights.
- High: 3 LED indicator lights.

Note: to adjust the intensity level to the next higher or lower level, you can press the + button to increase intensity or the - button to decrease intensity at any time during the brushing cycle.

Pressure Sensor

A change in vibration (and a slight change in sound) alerts you when you apply too much pressure while brushing. See chapter 'Features'.

Battery status (when handle is not on charger)

- 1 yellow LED indicates a low charge.
- 3 green LEDs indicate a full charge.

Note: If you hear a long beep and the yellow LED of the battery level indicator flashes, your toothbrush is completely empty and switches off. Place the toothbrush on the charger to recharge it.

Note: To keep the battery fully charged at all times, you may keep your Philips Sonicare on the charger when not in use.

Features

Pressure Sensor

This Philips Sonicare model comes with the Pressure Sensor feature activated. The Pressure Sensor works in all 3 modes and at all 3 intensity levels.

12 English

Activating or deactivating the Pressure Sensor feature

- 1 Attach the brush head to the handle.
- 2 Place the handle in the plugged-in charger or sanitizer.
 - To deactivate Pressure Sensor:
Press and hold the power on/off and mode button for 2 seconds. You hear 1 beep to indicate that the Pressure Sensor feature has been deactivated.
 - To activate Pressure Sensor:
Press and hold the power on/off and mode button for 2 seconds. You hear 2 beeps to indicate that the Pressure Sensor feature has been activated.

Quadpacer

- The Quadpacer is an interval timer that has a short beep and pause to remind you to brush the 4 sections of your mouth. Depending on the brushing mode or brushing routine you have selected, the Quadpacer beeps at different intervals during the brushing cycle. The Quadpacer works in all 3 modes and at all 3 intensity levels. (Fig. 10)

The Quadpacer feature on this model has been activated. To deactivate or reactivate the Quadpacer feature:

- 1 Place the handle in the plugged-in charger. or sanitizer.
- 2 Press and hold the mode button until:
 - a You hear 1 beep to indicate the Quadpacer has been deactivated.
 - b You hear 2 beeps to indicate the Quadpacer has been activated.

Smartimer

The Smartimer indicates that your brushing cycle is complete by automatically turning off the toothbrush at the end of the brushing cycle.

Dental professionals recommend brushing no less than 2 minutes twice a day.

Easy-start

This Sonicare comes with the Easy-start deactivated.

The Easy-start feature gently increases the power over the first 14 brushings to help you get used to the brushing with the Philips Sonicare.

Note: Each of the first 14 brushings must be at least 1 minute in length to properly advance through the Easy-start ramp-up cycle.

Activating or deactivating the Easy-start feature

- 1 Attach the brush head to the handle.
- 2 Place the handle in the plugged-in charger or sanitizer.
 - To activate Easy-start:
Press and hold the power on/off and mode button for 2 seconds. You hear 2 beeps to indicate that the Easy-start feature has been activated.
 - To deactivate Easy-start:
Press and hold the power on/off and mode button for 2 seconds. You hear 1 beep to indicate that the Easy-start feature has been deactivated.

Note: Using the Easy-start feature beyond the initial ramp-up period is not recommended and reduces Philips Sonicare effectiveness in removing plaque.

Note: Each of the first 14 brushings must be at least 1 minute in length to properly move through the Easy-start ramp-up cycle.

Sanitizing

- With the UV sanitizer you can clean your brush head after every use (Fig. 12).
- Unplug the sanitizer, discontinue use and call Consumer Care if:
- The UV light bulb remains on when the door is open.
 - The window is broken or missing from the sanitizer.
 - The sanitizer gives off smoke or a burning smell while it is operating.

Note: UV light can be harmful to the human eye and skin. This device should be kept out of reach of children at all times.

- 1 After brushing, rinse the brush head and shake off excess water.

Warning: Do not place the travel cap on the brush head during sanitization.

- 2 Place your finger in the door recess on top of the sanitizer and pull open the door (Fig. 13).
- 3 Place the brush head on one of the 2 pegs in the sanitizer with the bristles facing the light bulb (Fig. 14).

Note: Only clean Philips Sonicare snap-on brush heads in the sanitizer.

Note: Do not clean Philips Sonicare for Kids brush heads in the sanitiser.

- 4 Make sure the sanitizer is plugged into a live outlet of appropriate voltage.
- 5 Close the door and press the green power on/off button once to select the UV clean cycle.

Note: You can only turn on the sanitizer if the door is properly closed.

Note: The sanitizer stops if you open the door during the sanitizing cycle.

Note: The sanitizer cycle runs for 10 minutes and then automatically shuts off.

- The sanitizer is in operation when the light glows through the window.

Cleaning

Warning: Do not clean the brush head, handle, travel charger or the UV sanitizer in the dishwasher.

Toothbrush handle

- 1 Remove the brush head and rinse the metal shaft area with warm water. Make sure you remove any residual toothpaste (Fig. 15).

Caution: Do not push on the rubber seal on the metal shaft with sharp objects, as this may cause damage.

- 2 Wipe the entire surface of the handle with a damp cloth.

Brush head

- 1 Rinse the brush head and bristles after each use (Fig. 16).
- 2 Remove the brush head from the handle and rinse the brush head connection with warm water at least once a week.

Sanitizer

Warning: Do not immerse the sanitizer in water or rinse it under the faucet.

Warning: Do not clean the sanitizer when the UV light bulb is hot.

For optimal effectiveness, it is recommended to clean your sanitizer weekly.

- 1 Unplug the sanitizer.
- 2 Pull the drip tray straight out. Rinse the drip tray and wipe it clean with a damp cloth (Fig. 17).
- 3 Clean all interior surfaces with a damp cloth.
- 4 Remove the protective screen in front of the UV light bulb (Fig. 18).

To remove the screen, grasp the edges next to the snaps, gently squeeze and pull out the protective screen.

- 5 Remove the UV light bulb.

To remove the light bulb, grasp it and pull it out of the metal clasp.

- 6 Clean the protective screen and the UV light bulb with a damp cloth.
- 7 Reinsert the UV light bulb.

To reinsert the light bulb, align the bottom of the light bulb with the metal clasp and push the bulb into the clasp.

- 8 Reinsert the protective screen.

To reinsert the screen, align the pegs on the screen with the slots on the reflective surface near the UV light bulb. Push the screen straight into the slots on the sanitizer.

Storage

If you are not going to use the product for a long time, unplug it, clean it (see chapter 'Cleaning') and store it in a cool and dry place away from direct sunlight. You can use the integrated cord wrap to store excess cord of the UV sanitizer in an organized manner.

Replacement

Brush head

- Replace Philips Sonicare brush heads every 3 months to achieve optimal results.
- Use only Philips Sonicare replacement brush heads.

Ordering accessories

To buy accessories or spare parts, visit

www.shop.philips.com/service or go to your Philips dealer. You can also contact the Philips Consumer Care Center in your country (see the worldwide guarantee leaflet for contact details).

Recycling

- Do not throw away the product with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment.

Removing the rechargeable battery

Warning: PLEASE NOTE THAT THIS PROCESS IS NOT REVERSIBLE.

- To remove the rechargeable battery, you need a flat-head (standard) screwdriver. Observe basic safety precautions when you follow the procedure outlined below. Be sure to protect your eyes, hands, fingers, and the surface on which you work.
- 1 To deplete the rechargeable battery of any charge, remove the handle from the charger, turn on the Philips Sonicare and let it run until it stops. Repeat this step until you can no longer switch on the Philips Sonicare.

- 2 Insert a screwdriver into the slot located at the bottom of the handle and turn counterclockwise to release the bottom cap (Fig. 19).
- 3 Hold the handle upside down and push down on the shaft to release the internal components of the handle (Fig. 20).
- 4 Insert the screwdriver under the circuit board, next to the battery connections, and twist to break the connections. Remove the circuit board and pry the battery from the plastic carrier (Fig. 21).

The rechargeable battery can now be recycled and the rest of the product discarded appropriately.

Guarantee and support

If you need information or support, please visit

www.philips.com/support or read the separate worldwide guarantee leaflet.

Guarantee restrictions

The terms of the international guarantee do not cover the following:

- Brush heads.
- Damage caused by use of unauthorised replacement parts.
- Damage caused by misuse, abuse, neglect, alterations or unauthorised repair.
- Normal wear and tear, including chips, scratches, abrasions, discolouration or fading.

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©2016 Koninklijke Philips N.V. (KPNV). All rights reserved.
Philips and the Philips shield are trademarks of KPNV.
Sonicare and the Sonicare logo are trademarks
of Philips Oral Healthcare, LLC and/or KPNV.