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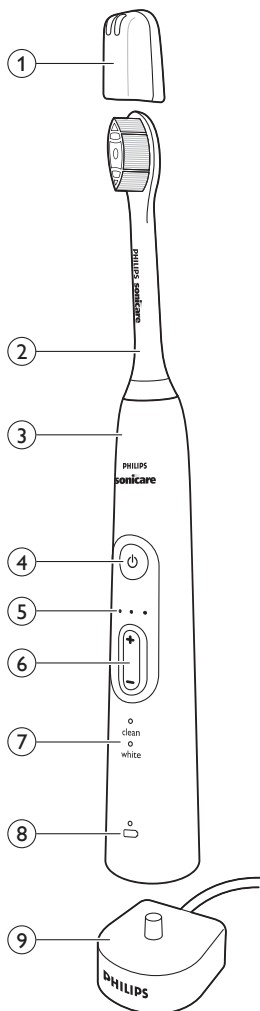
Introduction

Congratulations on your purchase and welcome to Philips! To fully benefit from the support that Philips offers, register your product at www.philips.com/welcome.

IMPORTANT SAFEGUARDS

Warning: READ ALL INSTRUCTIONS BEFORE USE

DANGERS

To reduce the risk of electrocution:

- 1 Always place and store the charger at a safe distance from water.
- 2 Do not reach for a charger that has fallen into water. Unplug immediately.
- 3 Never use a charger with a damaged cord or plug.

WARNINGS

To reduce the risk of burns, electrocution, fire or physical injury:

- 1 Do not use attachments other than those recommended by the manufacturer.
- 2 Do not plug the charger into an outlet with a voltage other than specified on the bottom of the charger. This product is designed to operate within a range of 100 to 240 volts. Voltage converters do not warrant voltage compatibility.

- 3** Never force the plug into an outlet; if the plug does not easily fit into the outlet, discontinue use.
- 4** Discontinue use if the product appears damaged in any way (brush head, handle or charger).
- 5** The mains cord cannot be replaced. If the mains cord is damaged, discard the charger.
- 6** Keep the cord away from heated surfaces.
- 7** Do not use the charger outdoors.
- 8** This product contains no user-serviceable parts. Refer to 'Warranty and support' if the product no longer works properly or needs repair.
- 9** Use this product only for its intended use as described in this document or as recommended by your dental professional.
- 10** This product is not intended for use by persons (including children) with reduced physical, sensory, or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- 11** Children should be supervised to ensure that they do not play with the product.
- 12** Philips Sonicare is a personal care device and is not intended for use on multiple patients in a dental office or institution.

MEDICAL WARNINGS

- 1 Consult your dental professional before you use this product if you have had oral or gum surgery in the previous 2 months.
- 2 Contact your dental professional if excessive bleeding occurs after using this product or bleeding continues to occur after 1 week of use.
- 3 Discontinue use of this product and contact a physician/dentist if discomfort or pain is experienced.
- 4 This product complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician or the device manufacturer prior to use.
- 5 Consult your physician if you have other medical concerns.

Electromagnetic fields (EMF)

This Philips device complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

SAVE THESE INSTRUCTIONS

Your Philips Sonicare (Fig. 1)


- 1 Hygienic travel cap
- 2 Brush head
- 3 Handle
- 4 Power on/off and mode button

- 5 Intensity level indicator
- 6 Intensity level button (+ and -)
- 7 Mode indicator
- 8 Battery level indicator
- 9 Charger



Note: The content of the box may vary based on the model purchased.

Getting started

Attaching the brush head

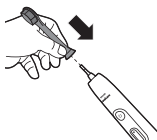
- 1 Align the brush head so the bristles point in the same direction as the front of the handle.
- 2 Firmly press the brush head down onto the metal shaft until it stops.

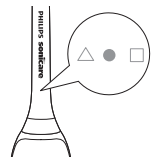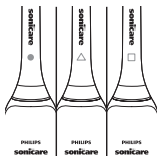Note: Philips Sonicare brush heads are imprinted with icons to easily identify your brush head.

Charging your Philips Sonicare

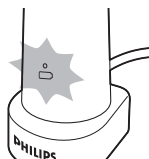
- 1 Plug the charger into an electrical outlet.
- 2 Place the handle on the charger.
 - The flashing light of the battery level indicator indicates that the toothbrush is charging.
 - When the battery is fully charged, the battery level indicator lights up solid green for 10 seconds before it turns off.

Note: It can take up to 48 hours to fully charge the battery, but the Philips Sonicare may be used before it is fully charged.

Battery status (when handle is not on charger)

- A flashing yellow battery level indicator indicates a low charge.
- A flashing green battery level indicator indicates a moderate charge.
- A solid green battery level indicator indicates a full charge.

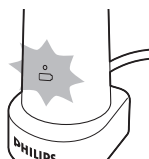
Note: If you hear 5 beeps and the battery level indicator flashes yellow, this means that you have to recharge the battery. When the battery is completely empty, the Philips Sonicare turns off. Place the Philips Sonicare on the charger to charge it.

To keep the battery fully charged at all times, you may keep your Sonicare on the charger when not in use.

Using your Philips Sonicare

Brushing instructions

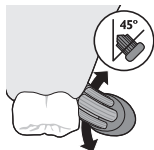- 1 Wet the bristles and apply a small amount of toothpaste.
- 2 Place the toothbrush bristles against the teeth at a slight angle (45 degrees), pressing firmly to make the bristles reach the gumline or slightly beneath the gumline.

Note: Keep the center of the brush in contact with the teeth at all times.

- 3 Press the power on/off and mode button to turn on the Philips Sonicare.
- 4 Apply light pressure to maximize Philips Sonicare's effectiveness and let the Philips Sonicare toothbrush do the brushing for you. Do not scrub.
- 5 Gently move the brush head slowly across the teeth in a small back and forth motion so the longer bristles reach between your teeth. Continue this motion throughout the brushing cycle.
- 6 To clean the inside surfaces of the front teeth, tilt the brush handle semi-upright and make several vertical overlapping brushing strokes on each tooth.


Note: To make sure you brush evenly throughout the mouth, divide your mouth into 4 sections using the Quadpacer feature (see chapter 'Features').

- 7 Begin brushing in section 1 (outside upper teeth) and brush for 30 seconds before you move to section 2 (inside upper teeth). Continue brushing in section 3 (outside lower teeth) and brush for 30 seconds before you move to section 4 (inside lower teeth).

Note: When using the White mode, two additional 15 second intervals can be used for polishing your front teeth.

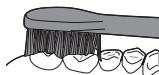- 8 After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs. You can also brush your tongue, with the toothbrush turned on or off, as you prefer.

Your Philips Sonicare is safe to use on:

- Braces (brush heads wear out sooner when used on braces).
- Dental restorations (fillings, crowns, veneers).

Note: Ensure that braces or dental restorations are properly adhered and not compromised. If not, braces or dental restorations may be damaged when they come into contact with the brush head.

Note: When the Philips Sonicare is used in clinical studies, the default 2-minute Clean mode with high intensity must be selected. The handle should be fully charged and the Easy-start feature deactivated.

Brushing modes

Clean mode

The ultimate in plaque removal (default mode).
Clean mode lasts for 2 minutes and includes four 30 second Quadpacer intervals.

White mode

To brighten and polish your front teeth.
White mode lasts for 2.5 minutes and includes four 30 second Quadpacer intervals and two and two additional 15 second Quadpacer intervals to focus on your visible front teeth.

- 1 Press the power on/off and mode button to turn on your Philips Sonicare.
- 2 Press the power on/off and mode button again within 2.5 seconds to switch to White mode.
- 3 Press the power on/off and mode button after 2.5 seconds to pause your Philips Sonicare.

Intensity levels

- **Low:** 1 LED indicator light
- **Medium:** 2 LED indicator lights
- **High:** 3 LED indicator lights

Note: To adjust the intensity level to the next higher or lower level, press the + intensity level button to increase intensity or the - intensity level button to decrease intensity at any time during the brushing cycle

Features

Smartimer

The Smartimer indicates that your brushing cycle is complete by automatically turning off the toothbrush at the end of the brushing cycle. Dental professionals recommend brushing no less than 2 minutes twice a day.

Easy-start

The Easy-start feature gently increases the power over the first 14 brushings to help you get used to the brushing with the Philips Sonicare. This Philips Sonicare comes with the Easy-start feature deactivated.

Activating or deactivating the Easy-start feature

- 1 Attach the brush head to the handle.
- 2 Place the handle in the plugged-in charger.

To activate Easy-start:

Press and hold the power on/off and mode button for 2 seconds. You hear 2 beeps and the battery level indicator lights up green to indicate that the Easy-start feature has been activated.

To deactivate Easy-start:

Press and hold the power on/off and mode button for 2 seconds. You hear 1 beep and the battery level indicator lights up yellow to indicate that the Easy-start feature has been deactivated.

Note: Using the Easy-start feature beyond the initial ramp-up period is not recommended and reduces Philips Sonicare effectiveness in removing plaque.

Quadpacer

The Quadpacer is an interval timer that has a short beep and pause to remind you to brush the 4 sections of your mouth evenly and thoroughly.

Depending on the brushing mode you have selected, the Quadpacer beeps at different intervals during the brushing cycle. See the "brushing modes" section for more details.

Activating or deactivating the Quadpacer feature:

Place the handle in the plugged-in charger.

- To activate Quadpacer:
Press and hold the intensity level button down (-) for 2 seconds. You hear 2 beeps and the battery level indicator lights up green to indicate that the Quadpacer feature has been activated.
- To deactivate Quadpacer:
Press and hold the intensity level button down (-) for 2 seconds. You hear 1 beep and the battery level indicator lights up yellow to indicate that the Quadpacer feature has been deactivated.

Charger confirmation

When the charger confirmation feature is activated, the handle beeps and vibrates briefly when it is placed on the charger base. This Philips Sonicare model comes with the charger confirmation feature deactivated.

Activating or deactivating the charger confirmation feature:

Place the handle on the charger.

- To activate charger confirmation:

Press and hold the intensity level button up (+) for 2 seconds. You hear 2 beeps and the battery level indicator lights up green twice to indicate that the charger confirmation feature has been activated.

- To deactivate charger confirmation:

Press and hold the intensity level button up (+) for 2 seconds. You hear 1 beep and the battery level indicator lights up yellow once to indicate that the charger confirmation feature has been deactivated.

Cleaning

Toothbrush handle


- 1 Remove the brush head and rinse the metal shaft area with warm water. Make sure you remove any residual toothpaste.

Caution: Do not push on the rubber seal on the metal shaft with sharp objects, as this may cause damage.

- 2 Wipe the entire surface of the handle with a damp cloth.

Brush head

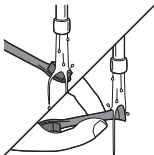- 1 Rinse the brush head and bristles after each use.
- 2 Remove the brush head from the handle and rinse the brush head connection with warm water at least once a week.

Charger

- 1 Unplug the charger before you clean it.
- 2 Wipe the surface of the charger with a damp cloth. Make sure you remove any toothpaste or other residue from the charger surfaces.

Storage

If you are not going to use the Sonicare for an extended period of time, unplug the charger from the electrical outlet, clean it and store it in a cool and dry place away from direct sunlight.

Replacement

Brush head

Replace Philips Sonicare brush heads every 3 months to achieve optimal results. Use only Philips Sonicare replacement brush heads.

Disposal

This appliance contains a rechargeable battery which must be disposed of properly. Contact your local town or city officials for battery disposal information. You can also call 1-800-8-BATTERY or visit www.rbrc.com for battery drop-off locations.


For assistance, visit our website www.philips.com/support or call 1-800-243-3050 toll free.

Removing the rechargeable battery

Warning: Only remove the rechargeable battery when you discard the appliance. Make sure the battery is completely empty when you remove it.


To remove the rechargeable battery, you need a towel or cloth, a hammer and a flat-head (standard) screwdriver. **Observe basic safety precautions when you follow the procedure outlined below.** Be sure to protect your eyes, hands, fingers, and the surface on which you work.

- 1 To deplete the rechargeable battery of any charge, remove the handle from the charger, turn on the Philips Sonicare and let it run until it stops. Repeat this step until you can no longer turn on the Philips Sonicare.
- 2 Remove and discard the brush head. Cover the entire handle with a towel or cloth.

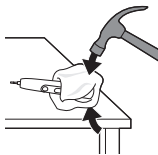
- 3 Hold the top of the handle with one hand and strike the handle housing 0.5 inch above the bottom end. Strike firmly with a hammer on all 4 sides to eject the end cap.


Note: You may have to hit on the end several times to break the internal snap connections.

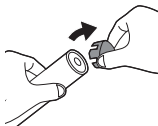
- 4 Remove the end cap from the toothbrush handle. If the end cap does not release easily from the housing, repeat step 3 until the end cap is released.

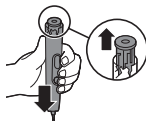
- 5** Holding the handle upside down, press the shaft down on a hard surface. If the internal components do not easily release from the housing, repeat step 3 until the internal components are released.

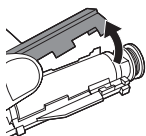- 6** Remove the rubber battery cover.

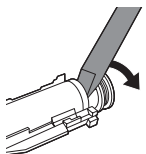
- 7** Wedge the screwdriver between the battery and the black frame at the bottom of the internal components. Then pry the screwdriver away from the battery to break the bottom of the black frame.

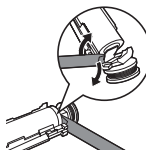- 8** Insert the screwdriver between the bottom of the battery and the black frame to break the metal tab connecting the battery to the green printed circuit board. This will release the bottom end of the battery from the frame.

18 English

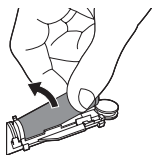- 9 Grab the battery and pull it away from the internal components to break the second metal battery tab.

Caution: Be aware of the sharp edges of the battery tabs so as to avoid injury to your fingers.

- 10 Cover the battery contacts with tape to prevent any electrical short from residual battery charge. The rechargeable battery can now be recycled and the rest of the product discarded appropriately.

Warranty and support

Philips warrants its products for two years after the date of purchase. Defects due to faulty materials and workmanship will be repaired or replaced at Philips' expense provided that convincing proof of purchase in the qualifying period is provided. Philips recommends use of genuine Sonicare brush heads for optimal performance. Use of other brush heads may damage your Sonicare toothbrush and limit your warranty protection. Contact our Consumer Care Center at 1-800-682-7664 (North America), outside North America contact your local Philips Consumer Care Center. Internet information: **www.sonicare.com** (North America) or **www.philips.com/support** (outside North America).

WARRANTY EXCLUSIONS

What is not covered under warranty:

- Brush heads.
- Damage caused by use of unauthorized replacement parts or unauthorized brush heads.
- Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.

IMPLIED WARRANTIES

ANY IMPLIED WARRANTIES, INCLUDING IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED TO THE DURATION OF THE EXPRESS WARRANTIES SET FORTH ABOVE. IN SOME STATES LIMITATIONS ON DURATION OF IMPLIED WARRANTIES DO NOT APPLY.

LIMITATION OF REMEDIES

IN NO EVENT SHALL PHILIPS OR ANY OF ITS AFFILIATED OR SUBSIDIARY COMPANIES BE LIABLE FOR ANY SPECIAL, INCIDENTAL OR CONSEQUENTIAL DAMAGES BASED UPON BREACH OF WARRANTY, BREACH OF CONTRACT, NEGLIGENCE, TORT, OR ANY OTHER LEGAL THEORY. SUCH DAMAGES INCLUDE, WITHOUT LIMITATION, LOSS OF SAVINGS OR REVENUE; LOSS OF PROFIT; LOSS OF USE; THE CLAIMS OF THIRD PARTIES INCLUDING, WITHOUT LIMITATION, DENTISTS AND DENTAL HYGIENISTS; AND COST OF ANY SUBSTITUTE EQUIPMENT OR SERVICES. SOME STATES DO NOT ALLOW THE EXCLUSION OR LIMITATION OF INCIDENTAL OR CONSEQUENTIAL DAMAGES.

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©2015 Koninklijke Philips N.V. (KPNV). All rights reserved.
Philips and the Philips shield are trademarks of KPNV.
Sonicare and the Sonicare logo are trademarks
of Philips Oral Healthcare, LLC. and/or KPNV.
de Philips Oral Healthcare, LLC. et (ou) KPNV.