
AirFloss Pro

PHILIPS
sonicare

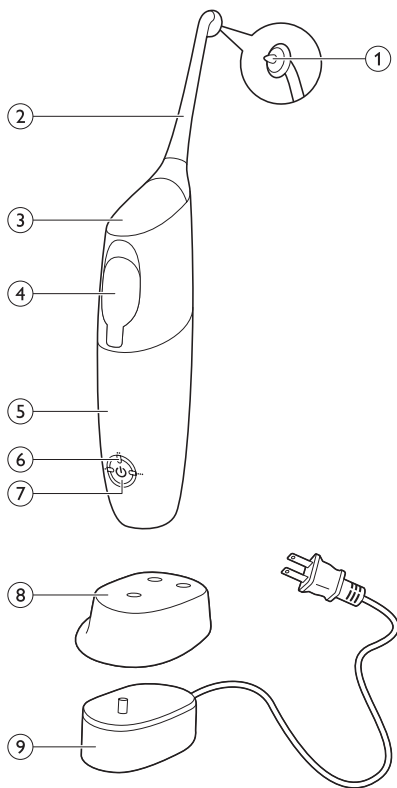Philips Sonicare AirFloss Pro

| | |
|-------------------|-----------|
| ENGLISH | 6 |
| DEUTSCH | 14 |
| ΕΛΛΗΝΙΚΑ | 24 |
| ESPAÑOL | 34 |
| FRANÇAIS | 43 |
| ITALIANO | 52 |
| NEDERLANDS | 61 |
| PORTUGUÊS | 70 |

Introduction

Congratulations on your purchase and welcome to Philips!

To fully benefit from the support that Philips offers, register your appliance at **www.philips.com/welcome**.

Use the AirFloss Pro as part of your oral care routine. In addition to brushing, it helps reduce plaque.

Use the AirFloss Pro with your favourite mouthwash as part of your oral care routine. In addition to brushing, the AirFloss Pro helps reduce the plaque in between teeth in order to improve your gum health.

Important

Read this user manual carefully before you use the appliance and save it for future reference.

Danger

- Keep the charger away from water. Do not place or store it over or near water contained in a bathtub, washbasin, etc. Do not immerse the charger in water or any other liquid. After cleaning, make sure the charger is completely dry before you connect it to the mains.

Warning

- The mains cord cannot be replaced. If the mains cord is damaged, discard the charger. Always have the charger replaced with one of the original type in order to avoid a hazard.
- Use only charger type HX6100 provided with the AirFloss Pro to recharge the battery.
- If the appliance is damaged in any way (AirFloss Pro nozzle, handle and/or charger), stop using it.

- This appliance contains no serviceable parts. If the appliance is damaged, contact the Consumer Care Centre in your country (see chapter 'Guarantee and support').
- Do not use the charger outdoors or near heated surfaces.
- This appliance can be used by children aged from 8 years and above and by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and if they understand the hazards involved. Cleaning and user maintenance shall not be made by children unless they are older than 8 and supervised. Keep the appliance and its cord out of reach of children aged less than 8 years.
- Children shall not play with the appliance.

Caution

- Consult your dentist before you use this appliance if you have had recent oral or gum surgery, suffer from an ongoing dental condition, had extensive dental work in the previous 2 months, or have suspected issues with your dental work (filings, crowns, etc.).
- Consult your dentist if excessive bleeding occurs after using this appliance or if bleeding continues to occur after 1 week of use.
- If you have medical concerns, consult your doctor before you use this appliance.

8 ENGLISH

- This Philips appliance complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician or the appliance manufacturer of the implanted appliance prior to use.
- Do not use attachments other than the ones recommended by the manufacturer.
- Do not exert excessive pressure on the guidance tip of the nozzle.
- To avoid physical injury, do not fill the reservoir with hot water. Do not put foreign objects in the reservoir.
- Only use the AirFloss Pro for its intended purpose as described in this user manual.
- To avoid physical injury, do not direct the spray under the tongue or into the ear, nose, eye or other sensitive area.
- Do not clean the AirFloss Pro nozzle, the handle, the charger, the water reservoir and its cover in the dishwasher or a microwave.
- To avoid an ingestion hazard, do not put isopropyl alcohol or other cleaning fluids in the reservoir.
- The AirFloss Pro is a personal care appliance and is not intended for use on multiple patients in a dental practice or institution.

Electromagnetic fields (EMF)

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

General description (Fig. 1)

- 1 Guidance tip
- 2 AirFloss Pro nozzle
- 3 Activation button
- 4 Reservoir and reservoir cover
- 5 Handle
- 6 Charge & burst mode indicators
 - Single burst: 1 LED
 - Double burst: 2 LEDs
 - Triple burst: 3 LEDs
- 7 Power on/off button
- 8 Removable nozzle holder (specific types only)
- 9 Charger

Note: The contents of the box may vary based on the type purchased.

Charging

Charge the AirFloss Pro for 24 hours before first use.

- 1** Put the plug of the charger in the wall socket (Fig. 2).
- 2** Place the handle on the charger (Fig. 3).
 - ▶ The charge indicators flash green to indicate that your AirFloss Pro is charging (Fig. 4).

Note: When it is fully charged, the AirFloss Pro has an operating time of up to 11 days when it is used in triple burst mode, and up to 33 days when it is used in single burst mode.

Note: When the charge indicators flash yellow, the battery is low and needs to be recharged (fewer than 3 uses left).

Preparing for use

- 1** Align the nozzle so the tip faces the front of the handle and press the nozzle down onto the handle ('click') (Fig. 5).

Note: Replace the nozzle every 6 months for optimal results. Also replace the nozzle if it becomes loose or no longer 'clicks' onto the handle.

10 ENGLISH

Note: The AirFloss Pro only supports AirFloss Pro nozzles. Do not attempt to use non-AirFloss Pro nozzles on the AirFloss Pro handle.

2 Open the cover of the reservoir ('click') (Fig. 6).

3 Fill the reservoir with mouthwash or water (Fig. 7).

Tip: For the best result and a fresher experience, we advise you to use mouthwash with your AirFloss Pro.

Note: Do not use mouthwash that contains isopropyl myristate or high concentrations of essential oils, as it may damage the appliance.

4 Close the cover of the reservoir by pressing ('click') (Fig. 8).

5 Press the power on/off button to switch on the appliance. (Fig. 9)
▶ The LEDs light up.

6 Press and hold the activation button until spray comes out of the tip of the nozzle (Fig. 10).

Using the appliance

- 1** Make sure the AirFloss Pro is switched on. If it is off, press and release the power on/off button to switch on the appliance. (Fig. 9)
- ▶ You can change the modes to have the desired amounts of bursts each time you press the activation button. For optimal results, use the triple burst mode (the default setting).
 - ▶ To change modes, press the power on/off button to cycle through the modes until you reach your desired mode. The LEDs show the selected mode:
 - Single burst: 1 LED (Fig. 11)
 - Double burst: 2 LEDs (Fig. 12)
 - Triple burst: 3 LEDs (Fig. 13)

Note: The AirFloss Pro remembers what mode was selected last and it is ready for the next use. If your AirFloss Pro becomes fully drained of power, it goes back to the triple burst default mode once it has been recharged.

2 Place the guidance tip between two teeth at your gum line. Slightly close your lips over the nozzle when it is in place, to avoid splashing (Fig. 14).

- 3 Press the activation button to deliver 1-3 bursts (depending on the selected mode) of air and micro-droplets of mouthwash or water between the teeth. (Fig. 15)
- 4 Slide the guidance tip along the gum line until you feel it settle between the next two teeth.
- 5 Continue this procedure for all spaces between your teeth, including behind your back teeth (Fig. 16).

Note: Refill the reservoir if needed.

Auto-burst function

Hold the activation button to deliver bursts continuously, approximately one burst every second.

- 1 To switch off AirFloss Pro, press and hold the power on/off button for one second.

Automatic shut-off function

The AirFloss Pro automatically switches off if it has not been used for one minute.

Cleaning

Do not clean the nozzle, handle or charger in the dishwasher.

Do not use cleaning agents to clean the appliance.

Handle and nozzle

- 1 Remove the nozzle from the handle. Rinse the nozzle after each use (Fig. 17).
- 2 Open the reservoir and rinse it under the tap to eliminate residue (Fig. 18).

Tip: You can also use a cotton bud to wipe the reservoir.

Note: Do not immerse the handle in water.

- 3 Wipe the entire surface of the handle with a damp cloth (Fig. 19).

12 ENGLISH

Charger

- 1 Unplug the charger before you clean it (Fig. 20).
- 2 Wipe the surface of the charger with a damp cloth (Fig. 21).

Storage

If you are not going to use the AirFloss Pro for an extended period of time, follow the steps below.

- 1 Open the cover of the reservoir and empty it.
- 2 Press the activation button until no more spray comes out of the nozzle.
- 3 Unplug the charger (Fig. 20).
- 4 Clean the AirFloss Pro nozzle, handle and charger. See chapter 'Cleaning'.
- 5 Store the AirFloss Pro in a cool, dry place away from direct sunlight.

Environment

- This symbol on a product means that the product is covered by European Directive 2012/19/EU (Fig. 22).
- This symbol means that the product contains a built-in rechargeable battery covered by European Directive 2006/66/EC which cannot be disposed of with normal household waste. Follow the instructions in section 'Removing the rechargeable battery' to remove the battery (Fig. 23).
- Inform yourself about the local separate collection system for electrical and electronic products and rechargeable batteries. Follow local rules and never dispose of the product and rechargeable batteries with normal household waste. Correct disposal of old products and rechargeable batteries helps prevent negative consequences for the environment and human health.

Removing the rechargeable battery

Please note that this process is not reversible.

- 1** To deplete the battery, repeatedly press the activation button until the AirFloss no longer produces any bursts of air (Fig. 24).
- 2** Remove the nozzle from the handle. (Fig. 25)
- 3** Firmly grab the handle with one hand at the top and one hand on the bottom (Fig. 26).
- 4** Snap the handle into two parts. (Fig. 27)

Note: This step requires a good amount of physical force. The top and bottom half of the handle should be separated once you snap the handle.

- 5** Remove the bottom half of the handle (Fig. 28).
- 6** Cut all 6 wires on the internal component with scissors (Fig. 29).
- 7** Completely pull apart the internal component from the top part of the handle (Fig. 30).
- 8** The battery is located inside the top half of the handle. Remove the battery and dispose of it properly (not in household waste) (Fig. 31).

Guarantee and support

If you need information or support, please visit www.philips.com/support or read the separate worldwide guarantee leaflet.

Guarantee restrictions

The terms of the international guarantee do not cover the following:

- AirFloss Pro nozzle.
- Damage caused by misuse, abuse, neglect, alterations or unauthorised repair.
- Normal wear and tear, including chips, scratches, abrasions, discolouration or fading.

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©2016 Koninklijke Philips N.V. (KPNV). All rights reserved.
Philips and the Philips shield are trademarks of KPNV.
AirFloss, Sonicare and the Sonicare logo are trademarks
of Philips Oral Healthcare, LLC and/or KPNV.
de Philips Oral Healthcare, LLC et (ou) KPNV.