

11) 保証とアフターサービス

保証書と修理サービスについて
(必ずお読みください)

- 保証は、本製品を最初に購入された方のみ適用されます。
- 保証書は、必ず「お買い上げ日・販売店名」等の記入をお確かめのうえ、販売店から受け取っていただき、内容よくお読みになったあと大切に保管してください。

保証期間	お買い上げ日から2年間
対象部分	ハンドル、充電器(ブラシヘッドは除く)

[補修用性能部品の保有期間]

- 補修用性能部品の保有期間は製造打ち切り後6年です。
- 性能部品とは、その製品の機能を維持するために必要な部品です。

[ご不明な点や修理に関するご相談]

- 修理に関するご相談などございましたら、フィリップスサポートセンターにお問い合わせください。
- 修理を依頼されるときは

- 修理をご依頼されまことにこの取扱説明書をよくお読みいただき、再度点検していただき。尚、異常のあるときは使用を中止し、フィリップスサポートセンターにお問い合わせください。

[保証期間中は]

- 製品と保証書をご用意いただき、フィリップスサポートセンターにお問い合わせください。保証書の記載内容により無料修理いたします。

[保証期間が過ぎているときは]

- 修理によって製品の機能が維持できる場合は、補修用性能部品の保有期間内であれば、ご希望により有料で修理させていただきます。

[海外での本製品の保証及びアフターサービスについて]

- 本製品の保証は海外においても効力です。(同シリーズ商品の取り扱いがある国へ限ります)

- 日本国外以外のフィリップスサービス部門においても保証期間内及び保証期間の経過後のアフターサービスを受けることができますが、この場合多少日数を要するものがあります。

- 海外にてアフターサービスを受けられる場合は、現地のフィリップスサービス部門にお問い合わせください。尚、お困りの点がございましょう下記まで連絡ください。

Philips Oral Healthcare, LLC. 22100 Bothell Everett Hwy Bothell, WA 98021 U.S.A.

Phone: +1(800) 662-7664










[お客様の個人情報のお取り扱いについて]

- お受けしましたお客様の個人情報は、株式会社フィリップス・ジャパンのホームページ

<http://www.philips.co.jp/>に掲載されています

お客様の個人情報の取り扱いについての当社の方針に基づき適切に管理いたします。

ENGLISH HX3293/80

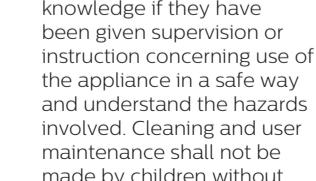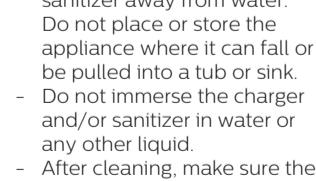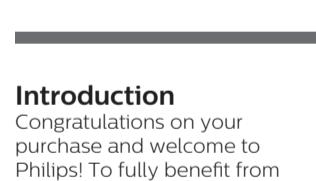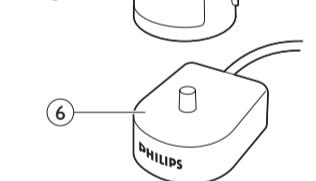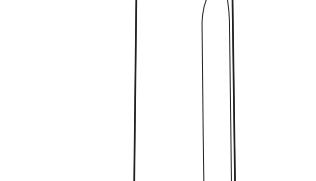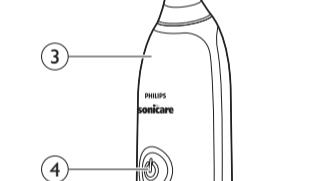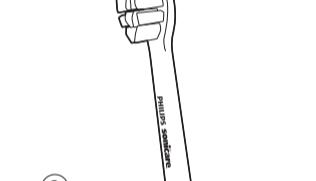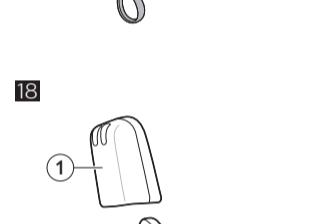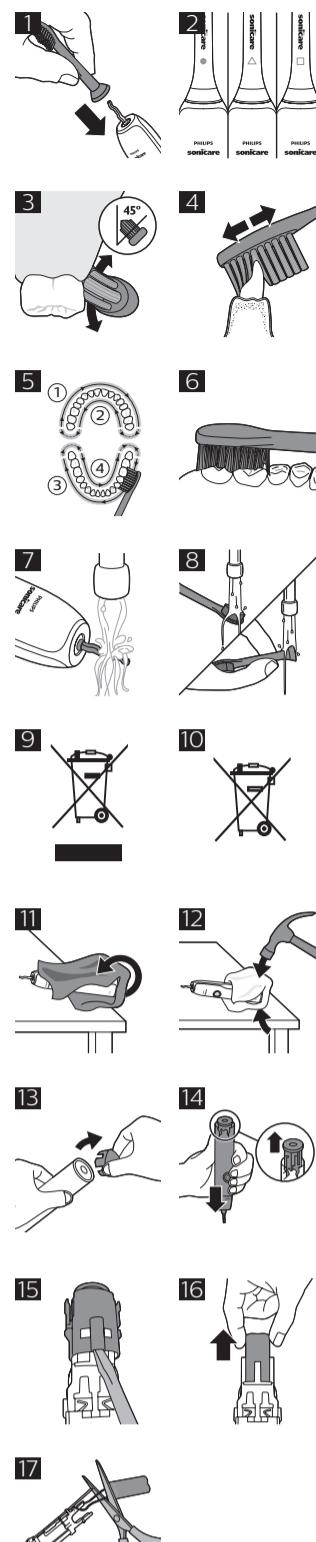

toothbrush handle, charger and/or sanitizer), stop using it.

- The mains cord cannot be replaced. If the mains cord is damaged, discard the charger and/or sanitizer.
- Always replace the charger and/or sanitizer with one of the original type in order to avoid a hazard.
- To recharge the battery with the travel case, only use the detachable USB cable and wall adapter if provided with this appliance.
- This appliance contains no user-serviceable parts. If the appliance is damaged, contact the Consumer Care Center in your country (see 'Warranty and support').
- Keep the cord away from heated surfaces.
- Do not use the charger and/or sanitizer outdoors or near heated surfaces.

CAUTION

- Do not clean the brush head, handle, charger, charger cover or sanitizer in the dishwasher.
- If you have had oral or gum surgery in the previous 2 months, consult your dentist before you use this appliance.
- Consult your dentist if excessive bleeding occurs after using this appliance or if bleeding continues to occur after 1 week of use. Also consult your dentist if you experience discomfort or pain when you use this appliance.
- This Philips appliance complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your doctor or the appliance manufacturer of the implanted appliance prior to use.
- If you have medical concerns, consult your doctor before you use this appliance.
- This appliance has only been designed for cleaning teeth, gums and tongue. Do not use it for any other purpose. Stop using the appliance and contact your doctor if you experience any discomfort or pain.
- This appliance is a personal care device and is not intended for use on multiple patients in a dental practice or institution.
- Do not use other brush heads than the ones recommended by Philips Sonicare.
- Stop using a brush head with crushed or bent bristles. Replace the brush head every 3 months or sooner if signs of wear appear.
- If your toothpaste contains peroxide, baking soda or bicarbonate (common in whitening toothpastes), thoroughly clean the brush head with soap and water after each use. This prevents possible cracking of the plastic.
- Avoid direct contact with products that contain essential oils or coconut oil. Contact may result in bristle tufts dislodging.

Electromagnetic fields (EMF)

- This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.
- Your Philips Sonicare toothbrush**
- 1 Hygienic travel cap
 - 2 Brush head
 - 3 Handle
 - 4 Power on/off button
 - 5 Battery level indicator
 - 6 Charger
- Note: The content of the box may vary based on the model purchased.

Intensities (specific types only)

Your power toothbrush gives you the option to choose between low and high intensity.

Note: When you use the toothbrush for the first time, the default setting is the low intensity.

- Press the power button once to turn on the toothbrush.

After 2 seconds of brushing, pressing the power button again will also pause the toothbrush.

Smartimer

The Smartimer indicates that your brushing cycle is complete by automatically turning off the toothbrush at the end of the brushing cycle.

Dental professionals recommend brushing no less than 2 minutes twice a day.

Note: If you press the power on/off button after you started the brushing cycle, the toothbrush will pause. After a pause of 30 seconds, the Smartimer resets.

Easy-start

The Easy-start feature gently increases power over the first 14 brushings to help you get used to brushing with the Philips Sonicare toothbrush. This Philips Sonicare toothbrush comes with the Easy-start feature activated.

Activating or deactivating Easy-start

- 1 Put the handle on the plugged-in charger.
- 2 Press and hold the power on/off button while the handle remains on the charger.

- The light of the battery level indicator flashes and the appliance beeps twice.

This indicates that the toothbrush is charging.

Note: It can take up to 48 hours to fully charge the battery, but you can use the Philips Sonicare toothbrush before it is fully charged.

Using your Philips Sonicare toothbrush

Brushing instructions

- 1 Wet the bristles and apply a small amount of toothpaste.

2 Place the toothbrush bristles against the teeth at a slight angle (45 degrees), pressing firmly to make the bristles reach the gumline or slightly beneath the gumline (Fig. 3).

Note