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Philips Oral Healthcare, LLC.  
22100 Bothell Everett Hwy  
Bothell, WA 98021 U.S.A.  
Phone +1 (800) 682-7664

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## ENGLISH

## ■Introduction

Congratulations on your purchase and welcome to Philips! To fully benefit from the support that Philips offers, register your product at [www.philips.com/welcome](http://www.philips.com/welcome).

## IMPORTANT SAFEGUARDS

Read this user manual carefully before you use the appliance and save it for future reference.

## READ ALL INSTRUCTIONS BEFORE USE

## DANGERS

To reduce the risk of electrocution:

- 1 Do not place or store the appliance where it can fall or be pulled into a tub or sink.
- 2 Do not place the product in or drop it into water or any other liquid.
- 3 Do not reach for a product that has fallen into water. Unplug immediately.
- 4 Do not use the product while bathing.

## WARNINGS


To reduce the risk of electrocution, fire, or injury to persons:

- 1 Discontinue use if this product appears damaged in any way (brush head, handle, charger). This product contains no serviceable parts. Refer to chapter 'Guarantee and service' on your Sonicare no longer works properly or needs repair.
- 2 Never use the charger if it has a damaged cord or plug.
- 3 Keep the cord away from heated surfaces.
- 4 Do not use the charger outdoors.
- 5 Do not use the charger if dropped into water.
- 6 This product is designed to clean your teeth, gums and tongue only. Use this product only for its intended use as described in this booklet. Discontinue use of this product and contact a physician/dentist if discomfort or pain is experienced.
- 7 This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- 8 Children should be supervised to ensure that they do not play with the appliance.
- 9 Do not use attachments other than those recommended by the manufacturer.
- 10 If your toothpaste includes peroxide, baking soda or other bicarbonate (common in whitening toothpastes), thoroughly clean your brush head and handle with soap and water after each use. Otherwise plastic cracking may occur.

## MEDICAL WARNINGS

1 Consult your dental professional before you use this product if you have had oral or gum surgery in the previous 2 months.

- 2 Contact your dental professional if excessive bleeding occurs after using this product or bleeding continues to occur after 1 week of use.
- 3 If you have a pacemaker or other implanted device contact your physician or the device manufacturer prior to use.
- 4 Sonicare has been tested and is compliant with safety standards for electromagnetic devices.
- 5 Consult your physician prior to using the Sonicare if you have medical concerns.

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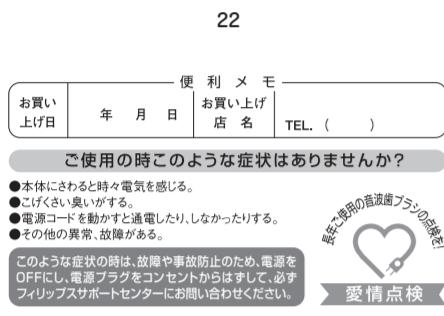## SAVE THESE INSTRUCTIONS

## ■Your Sonicare (Fig. 1)

- A Replaceable brush head  
B Battery charge indicator  
C Power on/off button  
D Handle  
E Charger  
F Hygienic travel cap

## ■Getting started

## Attaching the brush head

- 1 Align the front of the brush head with the front of the toothbrush handle.
- 2 Push the brush head onto the metal shaft (Fig. 2).

Note: There is a small gap between the brush head and the handle.

## Charging the Sonicare

- 1 Put the plug of the charger in the electrical outlet and place the handle on charger.

- Flashing light indicates that the toothbrush is charging
- Solid light indicates that the charging process is completed.

Note: In order to save energy the solid light will disappear about 30 seconds after the charging process is completed.

- 2 During brushing or within 30 seconds after brushing, the charge indicator shows the approximate remaining battery charge (Fig. 3).
  - Flashing light indicates that the battery is low (20% or less)
  - Flashing light and 3 beeps indicate that the battery is extremely low (10% or less) and toothbrush has to be recharged

Note: To keep the battery fully charged at all times, you may keep your Sonicare on the charger when not in use. It takes at least 24 hours to fully charge the battery.

## ■Using your Sonicare

- 1 Wet the bristles and apply a small amount of toothpaste.

- 2 Place the toothbrush bristles against the teeth at a slight angle towards the gumline (Fig. 4).

- 3 Press the power on/off button to switch on the Sonicare.

- 4 Apply light pressure to maximize Sonicare's effectiveness and let the Sonicare toothbrush do the brushing for you. Do not scrub.

- 5 Gently move the brush head slowly across the teeth in a small back and forth motion so the longer bristles reach between your teeth. Continue this motion throughout the brushing cycle.

Note: To make sure you brush evenly throughout the mouth, divide your mouth into 4 sections using the Quadpacer feature (see chapter 'Features').

- 6 Brush each section for 30 seconds. Begin brushing in section 1 (outside top teeth) and then brush section 2 (inside top teeth). Continue with section 3 (outside bottom teeth) and finally brush section 4 (inside bottom teeth) (Fig. 5).

Note: To make sure you brush evenly throughout the mouth, divide your mouth into 4 sections using the Quadpacer feature (see chapter 'Features').

## ■WARRANTY EXCLUSIONS

- Brush not covered under warranty:
- Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.

If you need service or information or if you have a question, please visit the Philips website at [www.philips.com](http://www.philips.com) or contact the Philips Consumer Care Center in your country (you find its phone number in the worldwide guarantee leaflet). If there is no Consumer Care Center in your country, go to your local Philips dealer.


## ■Guarantee and service

## Consumer Service Information

If you need service or information or if you have a question, please visit the Philips website at [www.philips.com](http://www.philips.com) or contact the Philips Consumer Care Center in your country (you find its phone number in the worldwide guarantee leaflet). If there is no Consumer Care Center in your country, go to your local Philips dealer.


## I

PHILIPS  
sonicare
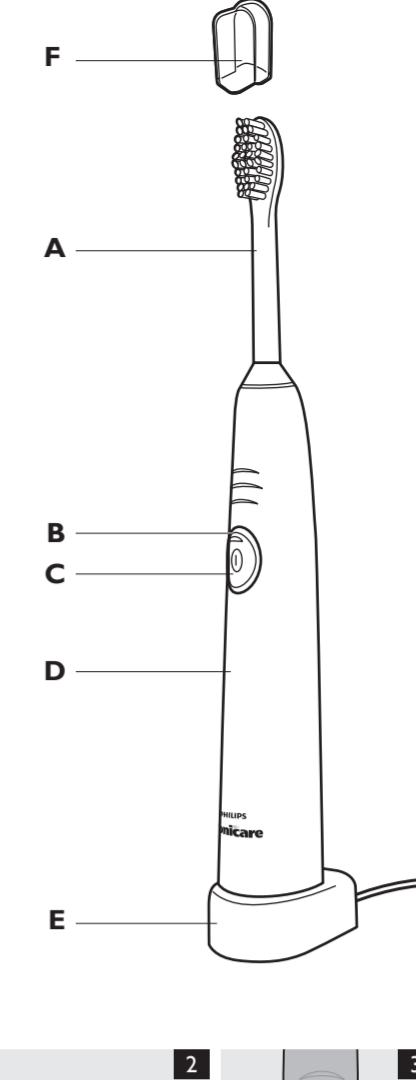
## EasyClean

Rechargeable sonic toothbrush

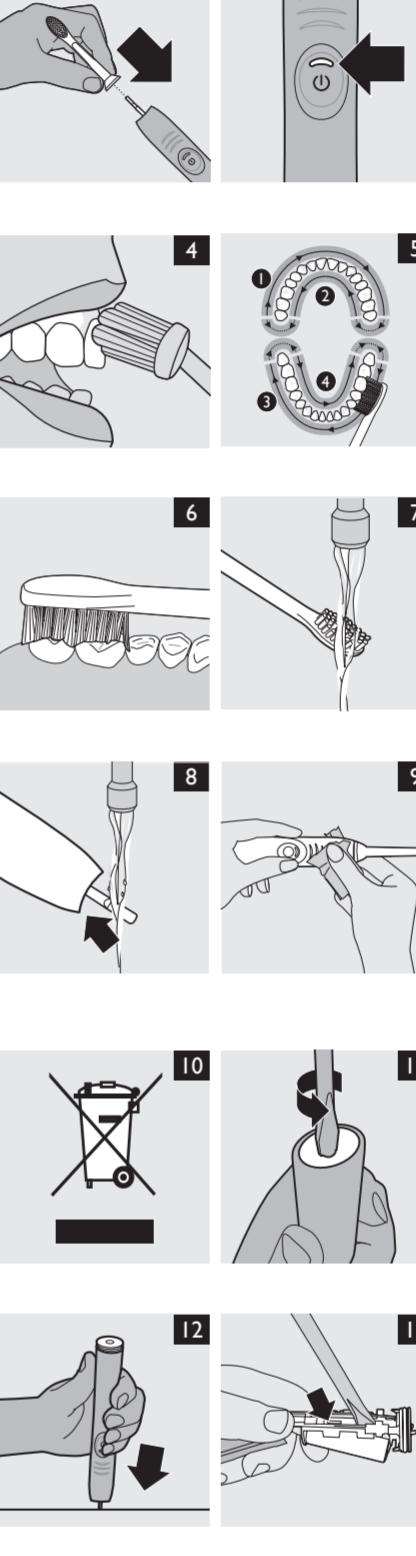[www.philips.com](http://www.philips.com)

Philips Oral Healthcare, LLC.  
22100 Bothell Everett Hwy  
Bothell, WA 98021 U.S.A.

Outside North America  
+1-425-396-2000

Protected by U.S. and international patents. Other patents pending.  
©2010 Koninklijke Philips N.V. (KPNV). All rights reserved.  
Philips and the Philips shield are trademarks of KPNV.  
Sonicare and the Sonicare logo are trademarks of Philips Oral Healthcare, LLC. and/or KPNV.  
Reg. U.S. Pat. & Tm. Off.

7 After you have completed the 2-minute brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs. You may also brush your tongue, with the toothbrush turned on or off, as you prefer (Fig. 6).

Your Sonicare is safe to use on:
 

- Braces (brush heads wear out sooner when used on braces)
- Dental restorations (fillings, crowns, veneers)

## ■Features

## Easy-start

This Sonicare model comes with the Easy-start feature activated. The Easy-start feature gently increases the power over the first 14 brushings.

Note: Each of the first 14 brushings must be at least 1 minute in length to advance through the Easy-start ramp-up cycle properly.

## Deactivating or activating the Easy-start feature

1 Attach the brush head to the handle.

2 Place the handle in the plugged-in charger.

- To deactivate Easy-start: Press and hold the power on/off button for 2 seconds. You hear 1 beep to indicate that the Easy-start feature has been deactivated.

- To activate Easy-start: Press and hold the power on/off button for 2 seconds. You hear 2 beeps to indicate that the Easy-start feature has been activated.

Note: Using the Easy-start feature beyond the initial ramp-up period is not recommended and reduces Sonicare's effectiveness in removing plaque.

## Smartimer

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©2010 Koninklijke Philips N.V. KPNV All rights reserved.