PHILIPS

Sonicare
for Kids

PHILIPS
sonicare

www.philips.com/Sonicare

©2019 Koninklijke Philips N.V. (KPNV). All rights reserved. Philips and the Philips shield are trademarks of KPNV. Sonicare and the Sonicare logo are trademarks of Philips Oral Healthcare, LLC and/or KPNV.


©2019 Koninklijke Philips N.V. (KPNV). 保留所有权利。飞利浦及飞利浦盾牌标志是 KPNV 的商标。Sonicare 和 Sonicare 徽标是 Philips Oral Healthcare, LLC 和/或 KPNV 的商标。

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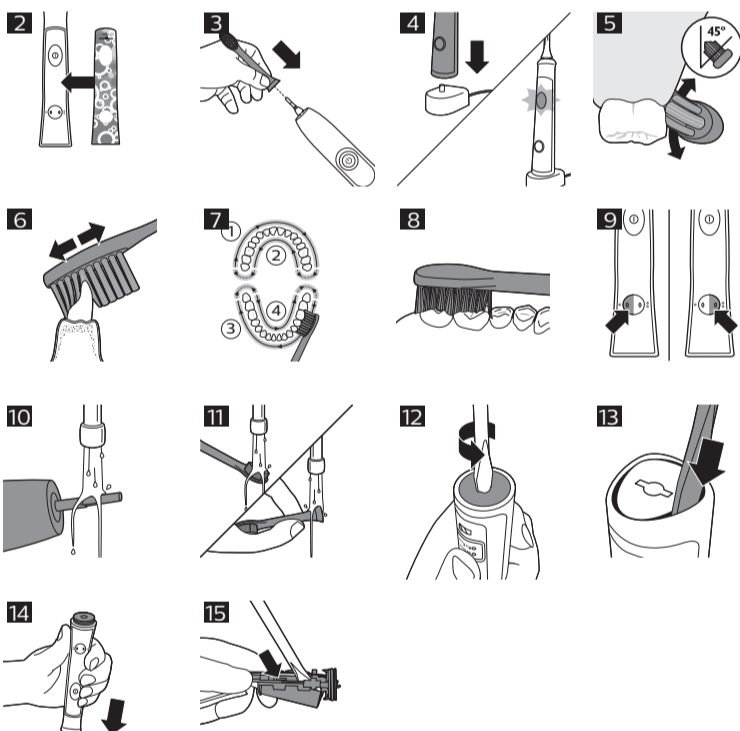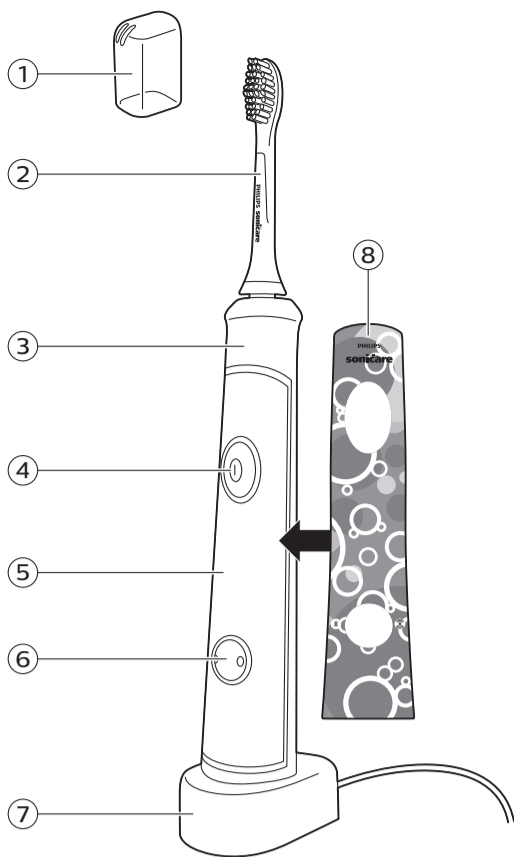
English

Important safety information

Read this important information carefully before you use the appliance and save it for future reference.

Note: The accessories and features of different models vary. For example, some models come with a sanitizer or are equipped with Bluetooth.

Danger

- Keep the charger and/or sanitizer away from water. Do not place or store the appliance where it can fall or be pulled into a tub or sink.
- Do not immerse the charger and/or sanitizer in water or any other liquid.
- After cleaning, make sure the charger and/or sanitizer is completely dry before you connect it to the wall socket.

Warning

- This appliance can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.
- Children shall not play with the appliance.
- If the appliance is damaged in any way (brush head, toothbrush handle, charger and/or sanitizer), stop using it.
- The mains cord cannot be replaced. If the mains cord is damaged, discard the charger and/or sanitizer.
- Always replace the charger and/or sanitizer with one of the original type in order to avoid a hazard.
- This appliance contains no user-serviceable parts. If the appliance is damaged, contact the Consumer Care Center in your country (see "Warranty and support").
- Keep the cord away from heated surfaces.

Caution

- Do not clean the brush head, handle, charger, charger cover or sanitizer in the dishwasher.
- If you have had oral or gum surgery in the previous 2 months, consult your dentist before you use this appliance.
- Consult your dentist if excessive bleeding occurs after using this appliance or if bleeding continues to occur after 1 week of use. Also consult your dentist if you experience discomfort or pain when you use this appliance.
- This Philips appliance complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your doctor or the appliance manufacturer of the implanted appliance prior to use.
- If you have medical concerns, consult your doctor before you use this appliance.
- This appliance has only been designed for cleaning teeth, gums and tongue. Do not use it for any other purpose. Stop using the appliance and contact your doctor if you experience any discomfort or pain.
- This appliance is a personal care device and is not intended for use on multiple patients in a dental practice or institution.
- Do not use other brush heads than the ones recommended by Philips Sonicare.
- Stop using a brush head with crushed or bent bristles. Replace the brush head every 3 months or sooner if signs of wear appear.
- If your toothpaste contains peroxide, baking soda or bicarbonate (common in whitening toothpastes), thoroughly clean the brush head with soap and water after each use. This prevents possible cracking of the plastic.
- Avoid direct contact with products that contain essential oils or coconut oil. Contact may result in bristle tufts dislodging.

Electromagnetic fields (EMF)

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

Radio Equipment Directive

Hereby, Philips declares that electric toothbrushes with radio frequency interface (Bluetooth) are in compliance with Directive 2014/53/EU. The full text of the EU declaration of conformity is available at the following internet address: www.philips.com/support

Introduction

Congratulations on your purchase and welcome to Philips! To fully benefit from the support that Philips offers, register your product at www.philips.com/welcome.

General description (Fig. 1)

- 1 Hygienic travel cap
- 2 Brush head
- 3 Handle with soft grip
- 4 Power on/off button with battery charge indicator
- 5 Panel for your favorite sticker
- 6 Brushing mode button with mode lights
- 7 Charger
- 8 Sticker

Preparing for use

Applying your favorite sticker

The Sonicare For Kids comes with replaceable stickers.

- 1 Apply your favorite sticker on the cut-out areas on the sticker (Fig. 2).

Note: Make sure your sticker lines up on the handle areas on the sticker with the buttons on the front of the handle.

Attaching the brush head

- 1 Align the brush head so the bristles point in the same direction as the front of the handle (Fig. 3).
- 2 Firmly press the brush head down onto the metal shaft until it stops.

Note: There is a small gap between the brush head and the handle.

Charging your Sonicare

- 1 Plug the charger into an electrical outlet.
- 2 Place the handle on the charger (Fig. 4).
 - The battery charge indicator in the power on/off button flashes and the handle beeps twice to indicate that the toothbrush is charging.
 - When the Philips Sonicare is fully charged, the battery charge indicator stops flashing and remains on.

Note: If the battery charge of the Philips Sonicare is low, you hear 3 beeps and the battery charge indicator flashes rapidly for 30 seconds after the brushing cycle.

Note: To keep the battery fully charged at all times, you may keep the Philips Sonicare on the charger when not in use. The charging indicator flashes until your toothbrush is fully charged.

Bluetooth connectivity

Before you start the connectivity procedure, make sure your mobile device runs on Bluetooth 4.0 or higher.

Note: For the latest information on device compatibility, please visit the App Store or Google Play.

Connecting your toothbrush

- 1 Download the Philips Sonicare For Kids app from the App Store or Google Play.
- 2 Confirm Bluetooth is switched on on your mobile device before you open the app. Check the settings on your mobile device to turn on Bluetooth.

- Open the Philips Sonicare For Kids app and follow the app instructions displayed on your mobile device.
- Connect the toothbrush to the app, turn on the toothbrush by pressing either the power on/off button or the brushing mode button.

Note: Make sure the toothbrush is not on the charger when you turn it on.

- The app automatically connects with the toothbrush. If the toothbrush does not connect, follow the troubleshooting steps provided in the app.

Note: Keep your mobile device close to your toothbrush when you use your toothbrush. Make sure your mobile device is safe in a dry place.

Using your Philips Sonicare

For children under age 8, it is recommended that a parent or adult caregiver perform toothbrushing.

Brushing instructions

- Wet the bristles and apply a small amount of toothpaste.
- Place the toothbrush bristles against the teeth at a slight angle (45 degrees), pressing firmly to make the bristles reach the gumline or slightly beneath the gumline. (Fig. 5)

Note: Keep the center of the brush in contact with the teeth at all times.

- Press the power on/off button to turn on the Philips Sonicare.
- Apply light pressure to maximize Philips Sonicare's effectiveness and let the Philips Sonicare toothbrush do the brushing for you. Do not scrub.
- Gently move the brush head slowly across the teeth in a small back and forth motion so the longer bristles reach between your teeth. Continue this motion throughout the brushing cycle.
- To clean the inside surfaces of the front teeth, tilt the brush handle semi-upright and make several vertical overlapping brushing strokes on each tooth (Fig. 6).

Note: To make sure you brush evenly throughout the mouth, divide your mouth into 4 sections using the Quadpacer feature.

- Begin brushing in section 1 (outside upper teeth) and brush for 30 seconds before you move to section 2 (inside upper teeth). Continue brushing in section 3 (outside lower teeth) and brush for 30 seconds before you move to section 4 (inside lower teeth) (Fig. 7).
- After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs (Fig. 15). You can also brush your tongue, with the toothbrush turned on or off, as you prefer.

Your Philips Sonicare is safe to use on:

- Braces (brush heads wear out sooner when used on braces).
- Dental restorations (fillings, crowns, veneers).

Note: When the Sonicare is used in clinical studies, the handle should be fully changed and the Easy-start feature deactivated

Brushing modes

The Philips Sonicare automatically starts in the default high mode.

- Press the brushing mode button to toggle between modes (Fig. 9). The green mode light indicates the selected mode.
 - High mode: Full-strength mode designed for older/independent brushers.
 - Low mode: Reduced power mode used for training and introducing younger children to the sonic experience.

Features

KidTimer

- The KidTimer plays a congratulatory sound sequence at the end of the pre-determined brushing time. Additionally, the KidTimer automatically turns off the handle when the brushing cycle is complete.
- When in high mode, the KidTimer is set at 2 minutes.
- When in low mode the KidTimer is set at 1 minute and slowly increases over time. The progression of time helps young children reach the dentist recommended 2 minutes of brushing.

Note: Instruct your child to brush until they hear the congratulatory tune and the toothbrush turns off.

KidPacer

The Sonicare For Kids uses a short series of tunes to alert the brusher to move to the next brushing quadrant of the mouth and concludes the brushing cycle with a congratulatory tune before it automatically turns off. This ensures that the child's teeth get a thorough, overall clean (see "Brushing instructions").

Note: Instruct your child to brush until they hear the congratulatory tune and the toothbrush turns off.

Easy-start

The Easy-start feature gently increases the power over the first 14 brushings to help you get used to the brushing with the Philips Sonicare. The Sonicare For Kids comes with the Easy-start feature activated.

Note: Each of the first 14 brushings must be at least 1 minute in length to properly advance through the Easy-start ramp-up cycle.

Activating or deactivating the Easy-start feature

- Place the handle on the charger.
- To activate Easy-start: Press and hold the power on/off button for 2 seconds. You hear 2 beeps to indicate that the Easy-start feature has been activated. To deactivate Easy-start: Press and hold the power on/off button for 2 seconds. You hear 1 beep to indicate that the Easy-start feature has been deactivated.

Note: Using the Easy-start feature beyond the initial ramp-up period is not recommended and reduces Philips Sonicare effectiveness in removing plaque.

Cleaning

Warning: Do not clean brush heads, handle, charger in the dishwasher.

Toothbrush handle

- Remove the brush head and rinse the metal shaft area with warm water. Make sure you remove any residual toothpaste (Fig. 10).

Caution: Do not push on the rubber seal on the metal shaft with sharp objects, as this may cause damage.

- Wipe the entire surface of the handle with a damp cloth.

Brush head

- Rinse the brush head and bristles after each use (Fig. 11).

Charger

- Unplug the charger before you clean it.
- Wipe the surface of the charger with a damp cloth. Make sure you remove any toothpaste or other residue from the charger surfaces.

Storage

If you are not going to use the Sonicare for an extended period of time, unplug the charger from the electrical outlet, clean it and store it in a cool and dry place away from direct sunlight.

Replacement

Brush head

- Replace Sonicare brush heads every 3 months to achieve optimal results.
- Use only Sonicare For Kids replacement brush heads.
- Brush head can be thrown away with the normal household waste.

Troubleshooting

This chapter summarizes the most common problems you could encounter with the appliance. If you are unable to solve the problem with the information below, visit www.philips.com/support for a list of frequently asked questions or contact the Consumer Care Center in your country.

Problem	Possible cause	Solution
I cannot fully attach the brush head. There is a gap between the brush head and handle.		This is normal. A small gap is necessary to allow movement in order to deliver the vibrations needed to clean your teeth.
My Philips Sonicare toothbrush vibration is less powerful than before	The Easy-start function is activated, the intensity level has been lowered or the brush head needs to be replaced.	<ul style="list-style-type: none"> Turn off the Easy-start function and/or raise the intensity level. If this does not resolve the problem, replace the brush head. We recommend you to replace your Philips Sonicare brush head after 3 months of use (based on brushing twice a day for 2 minutes each session). Activating or deactivating the Easy-start feature: <ol style="list-style-type: none"> Place the handle on the charger. <p>To activate Easy-start: Press and hold the power on/off button for two seconds. You will hear two beeps to indicate that the Easy-start feature has been activated.</p> <p>To deactivate Easy-start: Press and hold the power on/off button for two seconds. You will hear one beep to indicate that the Easy-start feature has been deactivated.</p>

Problem	Possible cause	Solution
I cannot connect the toothbrush to my smartphone.	You have not downloaded the app	Download the Philips Sonicare for Kids app (in the Google play store or App Store)
	Your smartphone does not meet the app requirements.	Check the device for these requirements: <ul style="list-style-type: none"> - Make sure the your smartphone is equipped with Bluetooth 4.0 or higher - Android 5.0 OS or higher. - iOS 10.0 OS or higher - Check the Apple App Store or Android Play Store for the latest OS version supported by the App.
	The Bluetooth cache on your smartphone needs to be cleared	<p>Android</p> Go to Settings Select 'Application Manager' Press 'More' or the 3 dots in the upper right corner Select 'Show system apps' Select 'Bluetooth' or 'Bluetooth share' Select Storage Tap 'Clear Cache' Restart your phone Try to connect again following the steps on the app
		<p>iOS</p> Go to 'Settings' Select 'Bluetooth' Make sure that Bluetooth is on Under 'my devices' tap the 'info icon' next to the Sonicare device Select 'Forget this device' Try to connect again following the steps on the app

Recycling

- Do not throw away the product with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment.
- This product contains a built-in rechargeable battery which shall not be disposed of with normal household waste. Please take your product to an official collection point or a Philips service center to have a professional remove the rechargeable battery.
- Follow your country's rules for the separate collection of electrical and electronic products and rechargeable batteries. Correct disposal helps prevent negative consequences for the environment and human health.

Removing the rechargeable battery

Warning: Only remove the rechargeable battery when you discard the appliance. Make sure the battery is completely empty when you remove it.

To remove the rechargeable battery, you need a flat-head (standard) screwdriver. Observe basic safety precautions when you follow the procedure outlined below. Be sure to protect your eyes, hands, fingers, and the surface on which you work.

- 1 To deplete the rechargeable battery of any charge, remove the handle from the charger, turn on the Philips Sonicare and let it run until it stops. Repeat this step until you can no longer turn on the Philips Sonicare.
- 2 Insert a screwdriver into the slot located in the bottom of the handle and turn counterclockwise until you see a gap between the bottom cap and the handle (Fig. 12).
- 3 Insert the screwdriver into the gap and pry the bottom cap out of the handle (Fig. 13).
- 4 Hold the handle upside down and push down on the shaft to release the internal components of the handle (Fig. 14).
- 5 Insert the screwdriver under the circuit board, next to the battery connections, and twist to break the connections. Remove the circuit board and pry the battery from the plastic carrier (Fig. 15).

The rechargeable battery can now be recycled or disposed of and the rest of the product discarded appropriately.

Warranty and support

If you need information or support, please visit www.philips.com/support or read the international warranty leaflet.

Warranty restrictions

The terms of the international warranty do not cover the following:

- Brush heads.
- Damage caused by use of unauthorized replacement parts.
- Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.

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