
PHILIPS


Lumea

BRI947

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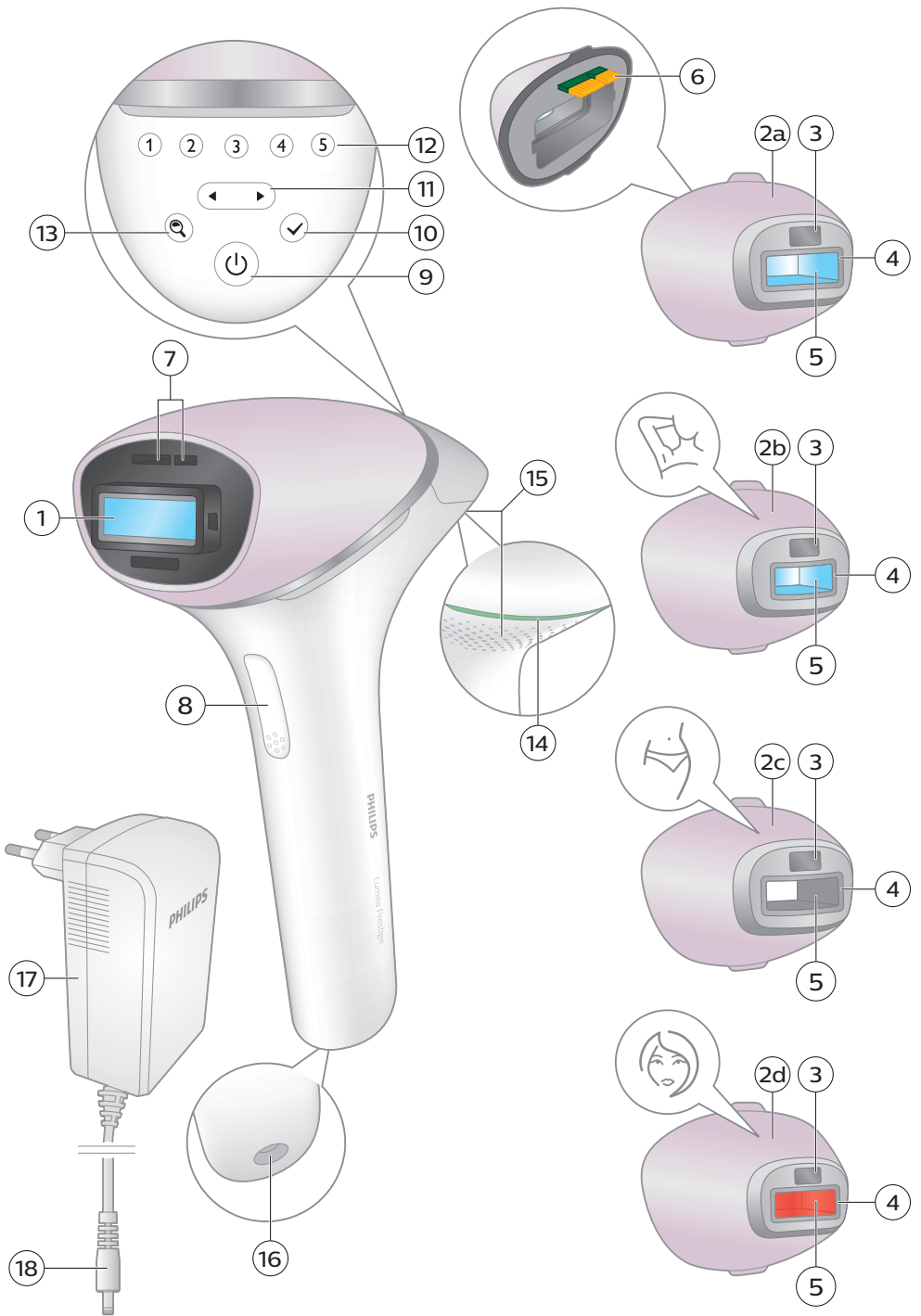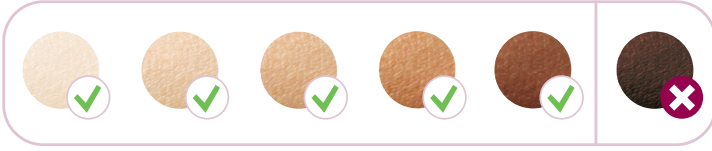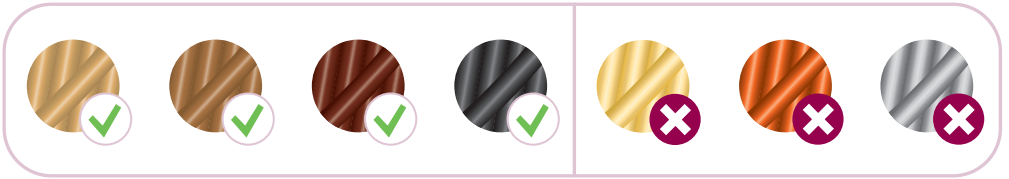
English 6
Español 30

Contents

Welcome	6
Device overview	7
Contraindications	7
Important safety instructions	9
How IPL works	11
Recommended treatment schedule	12
What to expect	12
How to use your Lumea Prestige	13
Learning more about your Lumea Prestige	18
Attachments	18
SmartSkin sensor	19
Selecting the light intensity	20
Two treatment modes: Stamp & Flash and Slide & Flash	21
Treatment time per area	22
Using your Lumea Prestige before and after tanning	23
Traveling with the device	23
After use	23
Troubleshooting	25
Warranty and support	27
Recycling	27
Technical specifications	27
Electromagnetic compatibility - Compliance information	28
Explanation of symbols	28

Welcome

Welcome to the beauty world of Lumea Prestige! You are only a few weeks away from silky-smooth skin.

Philips Lumea Prestige uses Intense Pulsed Light (IPL) technology, known as one of the most effective methods to continuously prevent hair regrowth. In close collaboration with skin experts we adapted this light-based technology, originally used in professional beauty salons, for easy and effective use in the comfort of your home. Philips Lumea Prestige is gentle and offers convenient and effective treatment at a light intensity that you find comfortable. Unwanted hairs are finally a thing of the past.

Philips Lumea Prestige has been designed and developed for women, but may be suitable for use by men.

To fully benefit from the support that Philips offers, register your product and find further information about your product at

www.philips.com/support to find our experts' advice, tutorial videos and FAQs and make the most of your Lumea Prestige.

Note: Keep these instructions with your product at all times.

Device overview

- 1 Light exit window with integrated UV filter
- 2 Attachments
 - a Body attachment
 - b Underarm attachment 
 - c Bikini attachment 
 - d Facial attachment 
- 3 SmartSkin sensor
- 4 Integrated safety system
- 5 Reflector inside the attachment
- 6 Electronic contacts
- 7 Opening for electronic contacts
- 8 Flash button
- 9 On/off button
- 10 Confirmation button
- 11 Toggle buttons
- 12 Intensity lights
- 13 SmartSkin button
- 14 'READY to flash' light
- 15 Air vents
- 16 Device socket
- 17 Adapter
- 18 Small plug
- 19 Pouch (not shown)
- 20 Cleaning cloth (not shown)

Contraindications

Conditions

Never use the device:

- If you have skin type VI (you rarely to never sunburn, very dark tanning). In this case you run a high risk of developing skin reactions, such as hyperpigmentation and hypopigmentation, strong redness or burns.
- If you are pregnant or breastfeeding as the device was not tested on pregnant or breastfeeding women.
- If you have any active implants such as a pacemaker, neurostimulator, insulin pump etc.

Medications/History

Never use the device if you take any of the medications listed below:

- If your skin is currently being treated with or has recently been treated in the past week with Alpha-Hydroxy Acids (AHAs), Beta-Hydroxy Acids (BHAs), topical isotretinoin and azelaic acid.
- If you have taken any form of isotretinoin Accutane or Roaccutane in the last six months. This treatment can make skin more susceptible to tears, wounds and irritations.

- If you are taking photosensitizing agents or medications, check the package insert of your medicine and never use the device if it is stated that it can cause photo-allergic reactions, photo-toxic reactions or if you have to avoid sun when taking this medicine.
- If you take anticoagulation medications, including heavy use of aspirin, in a manner which does not allow for a minimum 1-week washout period prior to each treatment.

Never use the device:

- If you have received radiation therapy or chemotherapy within the past 3 months.
- If you are on painkillers which reduce the skin's sensitivity to heat.
- If you take immunosuppressive medications.

Pathologies/Disorders

Never use the device:

- If you have diabetes or other systemic or metabolic diseases.
- If you have congestive heart disease.
- If you have a disease related to photosensitivity, such as polymorphic light eruption (PMLE), solar urticaria, porphyria etc.
- If you have a history of collagen disorder, including a history of keloid scar formation or a history of poor wound healing.
- If you have epilepsy with flashlight sensitivity.
- If your skin is sensitive to light and easily develops a rash or an allergic reaction.
- If you have a skin disease such as active skin cancer, you have a history of skin cancer or any other localized cancer in the areas to be treated.
- If you have a history of vascular disorder, such as the presence of varicose veins or vascular ectasia in the areas to be treated.
- If you have any bleeding disorder.
- If you have a history of immunosuppressive disease (including HIV infection or AIDS).

Skin conditions

Never use the device:

- If you have infections, eczema, burns, inflammation of hair follicles, open lacerations, abrasions, herpes simplex (cold sores), wounds or lesions and haematomas (bruises) in the areas to be treated.
- On irritated (red or cut), sunburned, recently tanned or fake-tanned skin.
- On the following areas: On moles, freckles, large veins, darker pigmented areas, scars, skin anomalies without consulting your doctor. This can result in a burn and a change in skin color, which makes it potentially harder to identify skin-related diseases.
- On the following areas: warts, tattoos or permanent make-up.

Location/areas

Never use the device on the following areas:

- Around the eyes and on or near the eyebrows.
- On the lips, nipples, areolas, labia minora, vagina, anus and the inside of the nostrils and ears.
- If you have had surgery in the areas to be treated in the last three weeks.
- Men must not use the device on the face and neck including all beard-growing areas, nor on the whole genital area.
- On areas where you use long-lasting deodorants. This can result in skin reactions.
- Over or near anything artificial like silicone implants, subcutaneous injection ports (insulin dispenser) or piercings.

Note: This list is not exhaustive. If you are not sure whether you can use the device, we advise you to consult your doctor.

Important safety instructions

Warnings


- Please read and follow these instructions to optimize the use of the device and prevent/ minimize the risk of injury, skin reactions and side effects.
- Before you use Lumea Prestige, you should clean your skin and make sure it is hair-free, entirely dry and free from any lotion or gel.
- Keep the device and the adapter dry.
- This device is not washable. Never immerse the device in water and do not rinse it under the tap.
- Water and electricity are a dangerous combination. Do not use this appliance in wet surroundings (e.g. near a filled bath, a running shower or a filled swimming pool). Moisture can damage the device and lead to a potential safety hazard.
- Always check the device before you use it. Do not use the device or adapter if it is damaged. Always replace a damaged part with one of the original type.
- If the device or adapter is broken or damaged, do not touch any inner part to avoid electric shock.
- Do not use the device if the UV filter of the light exit window and/or attachment is broken, as safe operation cannot be guaranteed under these circumstances.
- Do not attempt to open or repair your device. Opening Philips Lumea Prestige may expose you to dangerous electrical components and to pulsed light energy, either of which may cause bodily and/or eye injury.
- Always return the device to a service center authorized by Philips for examination or repair. Repair by unqualified people could cause an extremely hazardous situation for the user.


- No modification of this equipment is allowed. Do not modify, override or cut off any part of the device (e.g. adapter or the cord), as this causes a hazardous situation.
- Do not look at the flash while using the device. Use the device in a well-lit room so that the light is less glaring to your eyes. While not necessary for safety, you may wear sunglasses or goggles for your convenience if you find the light uncomfortable or too bright.
- Do not treat the same skin area more than once during a session. This may increase the risk of skin reactions.
- For hygienic reasons, the device should only be used by one person.
- Do not use hair removal creams in the area (to be) treated with the device, as chemicals may cause skin reactions.
- Do not use any pencil or pen to mark the areas to be treated. This may cause burns on your skin.
- The device can be stored at a temperature between -25 °C and 70 °C. Use the device at a temperature between 5 °C and 35 °C. The device will not work properly outside of this temperature range. This may cause burns on your skin.
- Hair removal by intense pulsed light sources can cause increased hair growth in some individuals. Based upon currently available data, the highest risk groups for this response are females of Mediterranean, Middle Eastern and South Asian heritage treated on the face and neck.
- If you experience skin irritation or skin redness after treatment, wait until it disappears before applying any product to your skin. If you experience skin irritation after applying a product to your skin, wash it off with water. If you still get skin reactions, stop using the device and consult your doctor.
- If shaving, epilating or waxing causes skin irritation, we advise you not to use the device until skin irritation has disappeared.
- This device is not intended for use by persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the device by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the device. Supervision is required to avoid potential risks such as exposure to the light output, electrocution, strangulation from the cables or choking on the attachments.
- Tanning with natural or artificial sunlight might influence the sensitivity and color of your skin. Perform a skin test to determine the appropriate light intensity setting.
- After each treatment, wait at least 48 hours before tanning to ensure any skin reaction has been resolved.
- If you notice a skin tone change since the last treatment, we advise you to perform a skin test and to wait 24 hours before your next treatment.
- Always check the device before you use it. Make sure that the SmartSkin sensor window is clean when using the device, to avoid improper operation (see 'D Clean & store after treatment').

- The adapter, light exit window and the filter of the attachments can become very hot (>210 °C/>410 °F) during usage. Always let the adapter, the light exit window as well as the filter and inner parts of the attachments cool down before you touch them.
- This device has a detachable power supply unit. Only use the adapter provided with this device. Do not use an extension cord. The reference number (AD2069x20020HF) can be found on your device. The 'x' in this number refers to the plug type for your country.


How IPL works

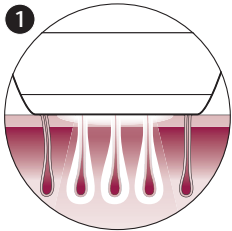With Intense Pulsed Light technology, gentle pulses of light are applied to the skin and absorbed by the hair root. The lighter the skin and the darker the hair, the better the pulses of light are absorbed.

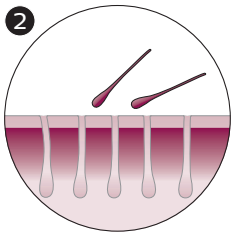The pulses of light stimulate the hair follicle to go into a resting phase. As a consequence, the hair sheds naturally and hair regrowth is prevented.

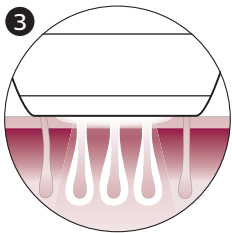The cycle of hair growth consists of different phases. IPL technology is only effective when the hair is in its growing phase. Not all hairs are in the growing phase at the same time. This is why we recommend you to follow the initial treatment phase (4 treatments, every treatment 2 weeks apart) and then the follow-up treatment phase (touch-ups every 4 weeks) to make sure all hairs are effectively treated in the growing phase.

Tip: To assure long lasting hair removal, touch-ups every 4 weeks are recommended.

Suitable body hair colors

Treatment with Lumea Prestige is not effective if you have light blond, grey, red or white hairs as hairs with these colors do not absorb enough light.

Note: To check if your body hair color allows usage of the device, consult the hair color table on the foldout page.

Recommended treatment schedule

Initial phase

To achieve effective hair reduction, follow the initial treatment schedule.

Initial treatment schedule

For the first 4 treatments, use Lumea Prestige once every two weeks to ensure that all hairs are treated.

- Treat yourself within 3 days before or after the planned treatment date.
- If you have not treated yourself within 3 days before or after the planned treatment date, restart the initial phase. Plan in 4 treatments, one treatment every two weeks.

Note: If you do not restart the initial phase, you may not see effective hair reduction.

Note: Replacing one of the IPL treatments with another hair removal method (waxing, epilating, etc.) will not help to reach the desired hair reduction.

If you want to remove hairs in between the Lumea Prestige treatments, you can use your normal hair removal method (excluding hair removal creams).

Touch-up phase

After the initial phase (4 treatments), we recommend touch-ups every 4 weeks.

Touch-up treatment schedule

Treat yourself every 4 weeks. Repeat this 8 times to achieve effective hair reduction. This is to maintain results and enjoy smooth skin for months.

- We advise you to treat yourself within 4 days before or after the planned treatment date.
- The results may vary based on your individual hair regrowth and also across different body areas.
- If you have completed the touch-up phase and you see the need to continue treatment (hair starts coming back), we advise you to restart the initial phase.

Tip: Mark the treatment schedule on your calendar to remind yourself of the planned treatments.

Note: Using the device more often does not enhance the effectiveness.

What to expect

After initial phase

- After the first treatment, it can take 1 to 2 weeks for the hairs to fall out. In the first weeks following the initial treatments, you still see some hairs growing. These are likely to be hairs that were not in their growing phase during the first treatments.


- After 2-3 treatments, you should see a noticeable reduction in hair growth. However, to effectively treat all hairs, it is important to keep on treating according to the recommended treatment schedule.
- After 4 treatments, you should see a significant reduction of hair growth in the areas that you treated with Lumea Prestige. A reduction of hair density should be visible as well.

During touch-up phase


- Keep on treating with frequent touch-ups (every 4 weeks) to maintain the result.

How to use your Lumea Prestige

A Remove hair & clean your skin

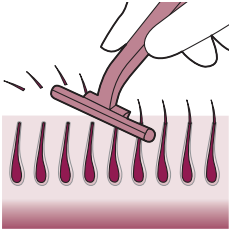
- 1 Remove hair before every treatment as long as hair is still visible. You can either shave, epilate or wax. Make sure you start to treat before new hair is visible on your skin. If you choose to wax, wait 24 hours before you use Lumea Prestige.

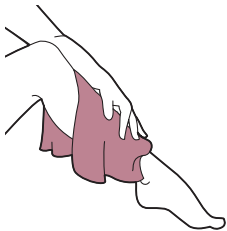
- 2 Clean and dry your skin. Make sure that it is free from any lotion or gel.

B Test skin 24 hrs before treatment

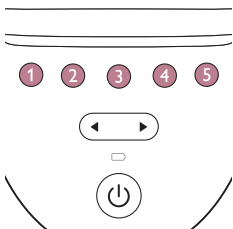
- 1 Test your skin to find the correct light setting:
 - For every new body area.
 - After recent tanning.

2 Select the correct attachment to test your skin (see 'Attachments').

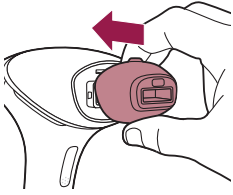3 Press the on/off button  to turn on the device.

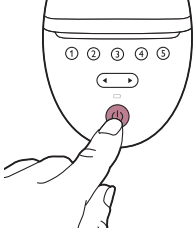
4 Choose an area for your skin test.


5 Press the device firmly onto your skin and hold it in place.


The integrated safety system prevents unintentional flashing without full skin contact.

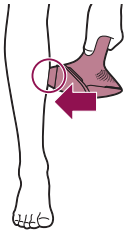6 Confirm that the 'READY to flash' light lights up white to ensure that there is full contact with your skin.

- If the 'READY to flash' light is orange, your skin is too dark to treat. Try Lumea Prestige on a lighter body area.

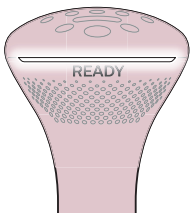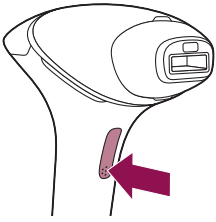7 Try each setting on your skin if comfortable:


- When the 'READY to flash' light is white, press the flash button on the handle once, to treat one spot with setting ①

The device makes a soft popping sound. The flash gives you a warm sensation.


• You will see a flash of light •

Note: You can toggle through the light intensity settings by pressing the toggle button.

- Move to next area. Treat one spot with setting ②

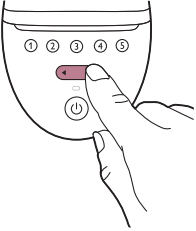- Move to next area. Treat one spot with setting ③
- Move to next area. Treat one spot with setting ④
- Move to next area. Treat one spot with setting ⑤

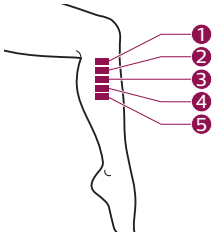8 Wait 24 hours.

24 hrs

9 Check your skin. The highest setting that...

- still felt comfortable and
- did not show a skin reaction (e.g. redness, irritation, rash)


... is the setting to start the treatment with.


C Treatment

IMPORTANT: First finish section A (Remove hair & clean your skin) and B (Test skin 24 hrs before treatment).


1 Connect the power adapter to the device and plug into a power outlet.

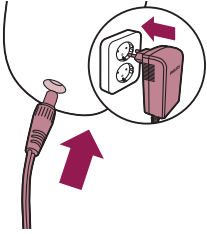
2 Press the on/off button  to turn on the device.

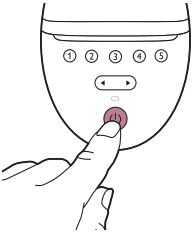3 Choose the setting you selected after the skin test (section B, step 9).

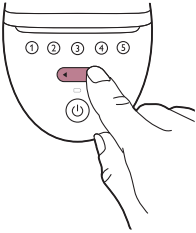
4 Press the device firmly onto your skin and hold it in place. The 'READY to flash' light should light up white.

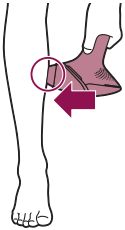
5 Press the flash button to release a flash on your skin. The flash treats the skin in that spot.

The device makes a soft popping sound. The flash gives you a warm sensation.


After each flash, it takes up to 3.5 seconds until the device is ready to flash again. You can release a flash when the 'READY to flash' light lights up white.

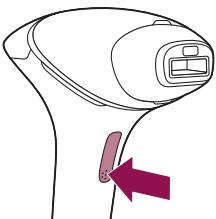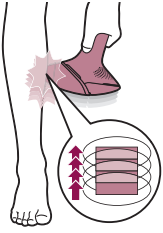
- 6 Move to an adjacent spot and press the flash button again. The surface of the attachment that is in contact with your skin is slightly larger than the light exit window. The surface of the attachment should overlap the previous spot to make sure you flash the complete skin area. **DO NOT** flash the same spot more than 3 times.
 - 7 Continue until you have treated the full area. There are two treatment modes (see 'Two treatment modes: Stamp & Flash and Slide & Flash').
 - 8 Turn off the device when you finished the treatment.
 - 9 Remove hair before every next treatment, as long as hair is still visible.
- Note: Do not try the device on difficult or sensitive areas (ankle and bony area).


D Clean & store after treatment

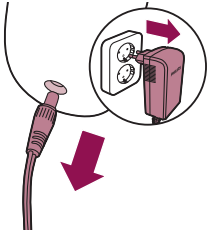
- 1 Unplug the device from the power outlet.
- 2 Let the device cool down.

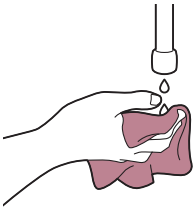
- 3 Moisten the soft cloth supplied with a few drops of water.

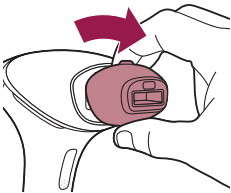- 4 Take the attachment off the device.

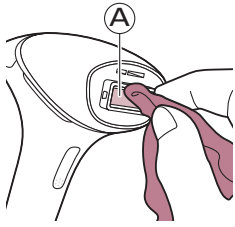5 Clean the light exit window of the device (A).

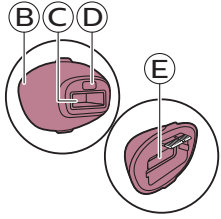6 Clean all surfaces of the attachment, including the metal reflector and window of the attachment (parts B, C, D, E).

7 Let all parts dry thoroughly.

8 Store in a dust-free place at room temperature.


Learning more about your Lumea Prestige

Attachments


For optimal results, it is important to use the correct attachment intended for the specific body area you want to treat. Lumea Prestige offers full-body treatment with several attachments designed for use on specific body areas. SenseIQ technology adapts the treatment program when a different attachment is connected.

Note: The device may not work anymore and show an error when there is dirt on the attachment connector. Clean the contact leads when this occurs.

To place the attachment, simply snap it onto the light exit window.

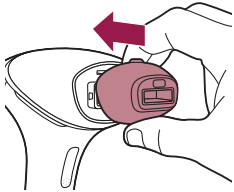To remove the attachment, pull it off the light exit window.

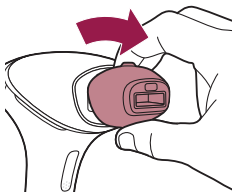Body attachment

The body attachment has the largest treatment window and a curved-in design to effectively cover and treat areas below the neckline, especially large areas such as legs, arms and stomach.

Facial attachment

The facial attachment has a precise flat design with extra integrated filter for safe and precise treatment on the sensitive skin above the upper lip, chin and sideburns. The device is only intended for removing unwanted body hair from areas below the cheekbones. When you use the device on areas with sharp contours like the jawline and chin, it may be difficult to establish full skin contact and release a flash.

Tip: Place your tongue between your upper lip and your teeth or puff out your cheeks to make the treatment easier.

Underarm attachment

The underarm attachment has a specifically curved-out design to treat underarm hairs that might be hard to reach.

Bikini attachment

The bikini attachment has a specialized design for effective treatment of the bikini area. It has a curved-out design and has a higher light intensity compared to body attachment. Hair in this area tends to be thicker and stronger than leg hair.

SmartSkin sensor

Lumea Prestige offers SenseIQ technology for a personalized hair removal treatment. The integrated SmartSkin sensor measures the skin tone at the start of each session and during the session. This provides two features:

- 1 The SmartSkin sensor prevents you to treat areas of your skin that are too dark: if it detects a skin tone that is too dark for treatment with Lumea Prestige, the 'READY to flash' light starts blinking orange. It does not flash when you press the flash button. The device automatically disables. This prevents you from developing skin reactions.
- 2 The SmartSkin sensor will help you to select the light intensity setting that is comfortable for your skin. For this you can use the SmartSkin button.

Selecting the light intensity

Lumea Prestige provides 5 different light intensities, you will always be able to change to a light intensity setting you find comfortable.

- 1 Press the on/off button to turn on the device. The device starts to operate at light intensity setting 1.

Note: To adjust the light intensity setting manually without using the SmartSkin button, use the toggle button until you have reached the required setting. The corresponding intensity light level lights up white.

- 2 Press the SmartSkin button  and put the device to your skin.
- 3 After the device has scanned your skin, the light intensity indicators will blink white to indicate which intensities are most comfortable based on your skin tone.
- 4 Press the confirmation button  to confirm that you want to use the indicated setting.

Lumea Prestige automatically uses the highest indicated intensity, which is shown by a blinking intensity indicator.

Lumea Prestige gives you the freedom to select the light intensity that you find most comfortable. Skin sensitivity can differ per person.


Therefore, the skin test is the most important guide for selecting your light intensity setting (see 'B Test skin 24 hrs before treatment').

- 5 Using Lumea Prestige should never be uncomfortably painful. If you experience discomfort, reduce the light intensity setting. You can do this by using the toggle buttons.

Note: The device is automatically disabled when your skin tone is too dark (brownish black or darker), to prevent you from developing skin reactions. The 'READY to flash' light will blink orange to indicate if the skin tone is too dark.

- 6 After moving to another body area or after recent tanning (see 'Using your Lumea Prestige before and after tanning'), perform a skin test to determine a comfortable light intensity setting. To re-activate the setting indication feature, press the SmartSkin button.


Note: Perform a skin test before you start to treat another body area or after a skin tone change. Determine the light intensity setting for each body area separately.

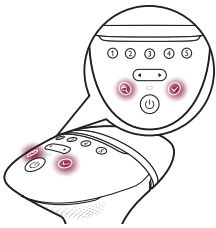
Two treatment modes: Stamp & Flash and Slide & Flash

Your Philips Lumea Prestige has two treatment modes for more convenient use on different body areas:


- The Stamp & Flash mode is ideal to treat small or curvy areas like knees and underarms. Simply press and release the flash button to release a single flash.

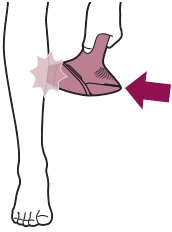- The Slide & Flash mode offers convenient use on larger areas like legs. Keep the flash button pressed while you slide the device over your skin to release several flashes in a row.

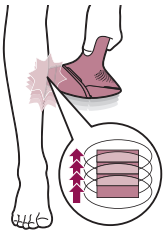
Treatment time per area

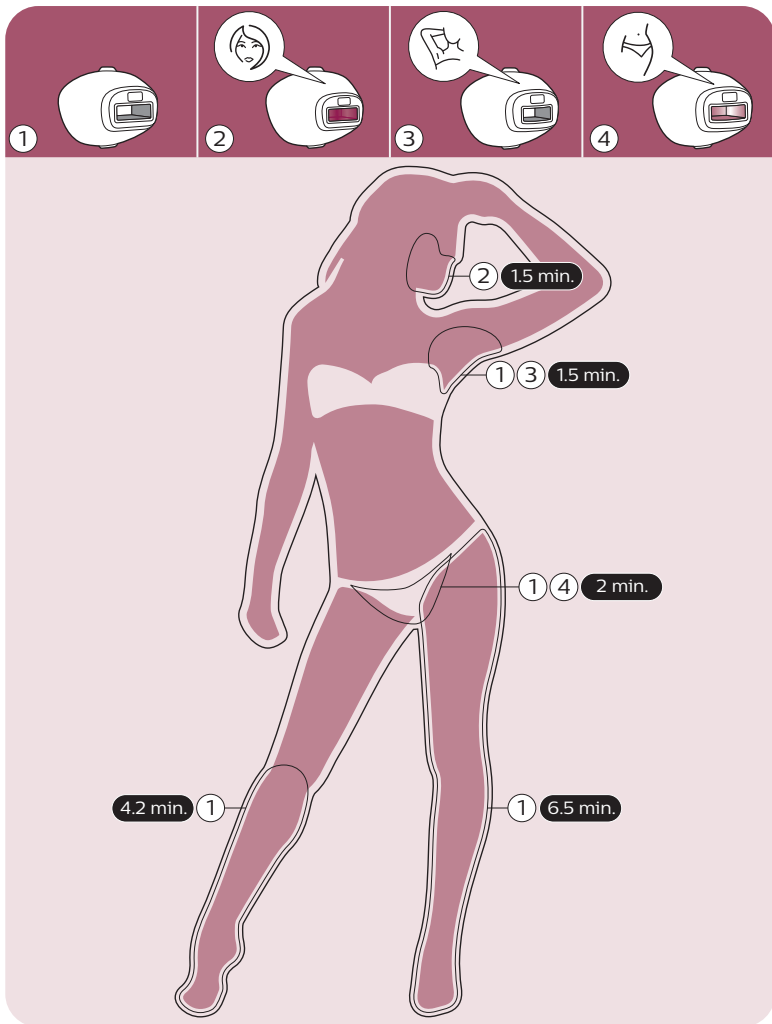Using your Lumea Prestige before and after tanning

Tanning with natural or artificial sunlight

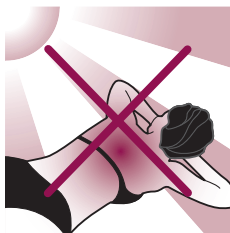Intentionally exposing your skin to natural or artificial sunlight with the aim of developing a tan influences the sensitivity and color of your skin.

Therefore the following is important:

- After each treatment, wait at least 48 hours before tanning. Even after 48 hours, make sure that the treated skin does not show any redness from the treatment anymore.
- If you expose your skin to the sun (without tanning intentionally) in the 48 hours after treatment, use a sunblock SPF 50+ on the treated areas. After this period, you can use a sunblock SPF 30+ for two weeks.
- After tanning, wait at least 2 weeks before you use Lumea Prestige and perform a skin test. This is to determine the appropriate light intensity setting (see 'B Test skin 24 hrs before treatment').
- Do not use Lumea Prestige on sunburned body areas.

Note: Occasional and indirect sun exposure does not qualify as tanning.

Tanning with creams

If you have used an artificial tanning lotion, wait until the artificial tan is completely gone before you use the device.

Traveling with the device

When you travel with the device, check with the airplane carrier to confirm that the device can be carried and/or used on the airplane.

After use

Aftercare

After use, you can apply lotions, creams, deodorant, moisturizer or cosmetics to the treated areas.

Warning: If you experience skin irritation or skin redness after treatment, wait until it disappears before you apply any product to your skin. If you experience skin irritation after applying a product to your skin, wash it off with water. If you still get skin reactions, stop using the device and consult your doctor.

Common skin reactions

Your skin may show slight redness and/or may prickle, tingle or feel warm. This reaction disappears quickly.

Dry skin and itching may occur due to shaving or a combination of shaving and light treatment. You can cool the area with an ice pack or a wet face

cloth. If dryness persists, you can apply a non-scented moisturizer on the treated area.

Rare side effects

- Burns, excessive redness (e.g. around hair follicles) and swelling: these reactions occur rarely. They are the result of using a light intensity that is too high for your skin tone. If these reactions do not disappear within 3 days, we advise you to consult a doctor. Wait with the next treatment until the skin has healed completely and make sure you use a lower light intensity.
- Skin discoloration: this occurs very rarely. Skin discoloration manifests itself as either a darker patch (hyperpigmentation) or a lighter patch (hypopigmentation) than the surrounding area. This is the result of using a light intensity that is too high for your skin tone. If the discoloration does not disappear within 2 weeks, we advise you to consult a doctor. Do not treat discolored areas until the discoloration has disappeared and your skin has regained its normal skin tone.
- Skin infection is very rare but is a possible risk following a (micro)wound, a skin burn, skin irritation etc.
- Epidermal heating (a sharply defined brownish area which often occurs with darker skin tones and is not accompanied with skin dryness): This reaction occurs very rarely. In case this reaction does not disappear within 1 week, we advise you to consult a doctor. Wait with the next treatment until the skin has healed completely and make sure you use a lower light intensity.
- Blistering (looks like small bubbles on the surface of the skin): this occurs very rarely. In case this reaction does not disappear within 1 month or when the skin gets infected, we advise you to consult a doctor. Wait with the next treatment until the skin has healed completely and make sure you use a lower light intensity.
- Scarring: often the secondary effect of a burn, which can take longer than a month to heal.
- Folliculitis (swelling around hair follicles combined with pustule formation): this reaction occurs very rarely and is the result of bacteria penetrating the damaged skin. In case this reaction occurs, we advise you to consult a doctor as folliculitis may need antibiotic ointment.
- Pain: this can occur during or after treatment if you use the device on skin that is not hair-free, if you use the device at a light intensity that is too high for your skin tone, if you flash the same area more than once and if you use the device on open wounds, inflammations, infections, tattoos, burns, etc.

Troubleshooting

Problem	Possible cause	Solution
The device/adapter becomes warm during use.	It is normal for the device and adapter to become warm (but not too hot to touch) during use.	Use the device in a slightly cooler environment and/ or let it cool down before continuing use.
	The fan is not working.	Check if the attachment is well connected. Clean the connectors on the attachment if needed. In case the attachment is well connected and it is possible to produce flashes with the device, contact the Consumer Care Center in your country, your Philips dealer or a Philips service center.
	The cooling airflow of the fan is blocked by hands or a towel.	Make sure the cooling airflow of the fan is unblocked.
When I place the device on the skin, it does not release a flash. The 'ready to flash' light blinks orange.	Your skin tone in the area to be treated is too dark.	Treat other body areas with lighter skin tones with Lumea Prestige.
The 'ready to flash' light blinks orange and all 5 intensity lights blink as well.	The device needs to be reset.	To reset the device, take the plug out of the socket, wait for 30 minutes to let the device cool down. The device should function normally again. In case it does not work again, contact the Consumer Care Center in your country.
The 'ready to flash' light does not light up white.	The device is not completely in contact with your skin.	Place the device at a 90 degree angle on the skin so that the integrated safety system is in contact with your skin.
The device produces a strange smell.	The light exit window or the SmartSkin sensor is dirty.	Clean the light exit window and the SmartSkin sensor carefully.
	You have not removed the hairs on the area to be treated properly. These hairs may get burned and can cause the smell.	Pretreat your skin before you use Lumea Prestige.
The device does not flash, the fan is not switched on and all 5 intensity lights blink.	The attachment is not attached properly.	Make sure you attach the attachment completely. If necessary, clean the electronic contacts on the attachment.

Problem	Possible cause	Solution
The skin feels more sensitive than usual during treatment. I experience discomfort when I use the device.	<p>The light intensity setting you use is too high.</p> <p>You did not remove the hairs on the areas to be treated.</p> <p>The UV filter of the light exit window is broken.</p>	<p>Check if you have selected a comfortable light intensity setting. If necessary, select a lower setting.</p> <p>Pretreat your skin before you use Lumea Prestige.</p> <p>If the UV filter is broken, do not use the device anymore. Contact the Consumer Care Center in your country, your Philips dealer or a Philips service center.</p>
The skin feels more sensitive than usual during treatment. I experience discomfort when I use the device.	You treated an area for which the device is not intended.	<p>Never use the device on the following areas: on moles, freckles, large veins, darker pigmented areas, scars, skin anomalies, warts, tattoos or permanent make-up, around the eyes and on or near the eyebrows, on the lips, nipples, areolas, labia minora, vagina, anus and the inside of the nostrils and ears. Men must not use the device on the face and neck including all beard-growing areas, nor on the whole genital area.</p> <p>Never use the device on moles, freckles, large veins, darker pigmented areas, scars, skin anomalies without consulting your doctor.</p> <p>Never use the device on warts, tattoos, permanent make-up.</p> <p>Never use the device irritated (red or cut), sunburned, recently tanned or fake-tanned skin.</p> <p>Never use the device if you have infections, eczema, burns, inflammation of hair follicles, open lacerations, abrasions, herpes simplex (cold sores), wounds or lesions and haematomas (bruises) in the areas to be treated.</p>
There is no glass filter in my attachment.	This is normal.	No action required: there is no filter in the body and underarm attachment. Only the facial, precision and bikini attachment have a specialized filter.
The skin reaction after treatment lasts longer than usual.	You have used a light intensity setting which is too high for you.	Select a lower intensity next time (see 'Selecting the light intensity').
The flash is very bright to my eyes. Do I need to wear goggles?	No, Philips Lumea Prestige does not hurt your eyes.	The scattered light produced by the device is harmless to your eyes. Do not look at the flash while using the device. It is not necessary to wear goggles during use. Use the device in a well-lit room so that the light is less glaring to your eyes. Be sure to make good skin contact to avoid scattered light.

Problem	Possible cause	Solution
The results of the treatment are not satisfactory.	You have used a light intensity setting which is too low for you.	Select a higher setting next time.
	You did not flash an area adjacent to an area you treated before.	You have to release the flashes close to, but not overlapping, the previous spot.
	The device is not effective on your body hair color.	If you have light blond, grey, red or white hair, the treatment is not effective.
	You do not use the device as often as recommended.	To remove all hairs successfully, we advise you to follow the recommended treatment schedule. You can reduce the time between treatments, but do not treat more often than once every two weeks.
	You respond more slowly to IPL treatment.	Continue using the device for at least 6 months, as hair regrowth can still decrease over the course of this period.

Warranty and support

If you need information or support, please visit www.philips.com/support or read the international warranty leaflet.

Recycling

- Do not throw away the product with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment.
- Follow your country's rules for the separate collection of electrical and electronic products. Correct disposal helps prevent negative consequences for the environment and human health.

Technical specifications

Transport and storage conditions	Lumea Prestige remains operational in NORMAL USE within its specification after transport or storage in the following environmental range: -25 °C to +70 °C up to 90 % RH, non-condensing, 700 hPa to 1060 hPa
Operating temperature conditions	Lumea Prestige complies with its specifications when operated in NORMAL USE under the following environmental operating conditions: +5 °C to +35 °C up to 90 % RH, non-condensing, 700 hPa to 1060 hPa
Emitted wavelengths	565 to 1400 nm
Pulse duration	1 - 2.5 ms, for all settings
Flash interval	Between 1 - 3.5 s

Optical exposure	Body attachment: 2.5 – 5.6 J/cm ² Underarm attachment: 2.5 – 6.5 J/cm ² Precision/Bikini attachment: 2.5 – 6.5 J/cm ² Facial attachment: 2.4 – 5.9 J/cm ²
------------------	--

Adapter input	100 - 240 V ~ 50/60 Hz 1.5 A
---------------	------------------------------

Adapter output	19.5 V  3.33 A 65.0 W
----------------	--

Average active efficiency	≥88.0 %
---------------------------	---------

Efficiency at low load (10 %)	≥79.0 %
----------------------------------	---------

No-load power consumption	<0,21 W
------------------------------	---------

Electromagnetic compatibility - Compliance information

This device complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

Explanation of symbols

The warning signs and symbols are essential to ensure that you use this product safely and correctly and to protect you and others from injury. Below you find the meaning of the warning signs and symbols on the label and in the user manual.

This symbol identifies the manufacturer of the product: Philips Consumer Lifestyle B.V. Tussendiepen 4, 9206 AD Drachten. The Netherlands

This symbol means that you have to read the user manual carefully before you use the device. Save the user manual for future reference.

This symbol means that the user should consult the instructions for use for important safety information such as cautions that cannot, for a variety of reasons, be presented on the device itself.

This symbol means that the user should consult the instructions for use for important safety information such as contraindications and warnings that cannot, for a variety of reasons, be presented on the device itself.

To identify the manufacturer's serial number.

To identify the manufacturer's catalog number.

Symbol to identify equipment meeting the safety requirements specified for Class II equipment.

This symbol indicates that the equipment is suitable for direct current only.

Symbol indicates that the device should not be used around the eyes and near the eyebrows.

This symbol means: 'Do not use this device in a bathtub, shower or water-filled reservoir.'

This symbol means that the transformer is designed to be short-circuit proof.

This symbol means that the transformer is considered a Switch Mode Power Supply (SMPS) unit.

This symbol means that the device is designed for indoor use.

t_a

This symbol indicates the maximum temperature at which the transformer may be operated continuously under normal conditions of use.

A recycling symbol to identify the material of which an item is made, to facilitate recycling or other reprocessing. The symbol may include a number and/or abbreviation.

This symbol identifies the switch or switch position by means of which part of the equipment is switched on in order to bring it into the stand-by condition, and to identify the control to shift to or to indicate the state of low power consumption.

Symbol indicates the energy efficiency of the adapter.
Energy efficiency complies with NOM-029-ENER-2017.

Certification mark for Mexico that shows that this product has been evaluated and determined to meet the relevant safety requirements

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