

PHILIPS

sonicare

Sonic electric
toothbrush with app

DiamondClean
9000

Connected brushing made easy

Built-in pressure sensor

Smart brush head recognition

4 modes, 3 intensities

HX9917/88



Whiter, healthier teeth for life

Remove up to 100% more stains***

Philips Sonicare's best whitening in our most elegant sonic electric toothbrush.
Switch to Philips Sonicare.

Enjoy a stylish and user-friendly design

- Simple to charge and style to boot

DiamondClean users experience improved oral health

- All-in-One brush head for complete care
- Up to 7x healthier gums**** in just 2 weeks

Let the DiamondClean 9000 adapt to your needs

- Personalized brushing experience
- Lets you know if you're pressing too hard
- BrushSync automatically selects the best mode for you
- Always know when to replace your brush heads

Stay on track with a personalized progress report

- Start and keep up healthy habits


Highlights

Premium All-in-One brush head




Clean without compromise with the A3 Premium All-in-One brush head. Our all-time best has angled bristles to help remove up to 20x more plaque even in hard-to-reach areas*. Triangular tips work to remove up to 100% more stains in less than two days***. And longer bristles clean deep for gums that are up to 15x healthier in just two weeks**. All that without ever switching brush heads. Did you know: brush heads become less effective after 3 months of use? Our BrushSync™ feature can remind you when it's almost time for a replacement.

Care for your gums

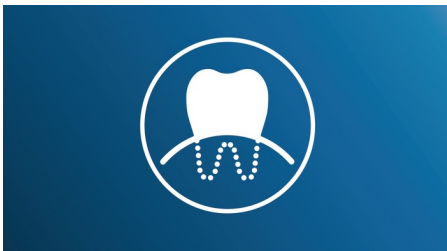

Click on the G3 Premium Gum Care brush head to improve your gum health. Its smaller size and targeted gum line bristles deliver gentle yet effective cleaning along the gum line, where gum disease starts. It's clinically proven to deliver up to 100% less gum inflammation**** and up to 7x healthier gums in just two weeks.****

Progress Report

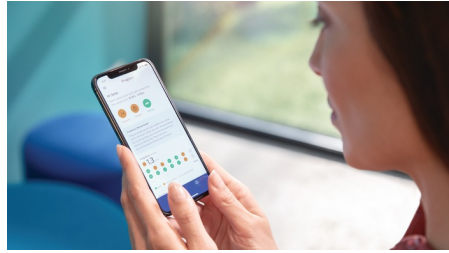


The DiamondClean 9000 gives you the guidance you need to improve and maintain healthy brushing habits in between dentist checkups. Built-in smart sensors let you know when you're using too much pressure, and by connecting your brushing experience to the Sonicare app, a personalized Progress Report helps you stay on track, to see how much you've improved over time.

4 modes, 3 intensities



DiamondClean 9000 comes with Clean, White+, Gum Health and Deep Clean+ to take care of your brushing needs. Clean mode is for exceptional everyday cleaning, White+ is for removing stains, Gum Health provides a gentle yet effective clean for gums and Deep Clean+ gives you an invigorating deep clean. Three intensities allow you to choose between a higher setting to boost your clean and a lower one for more-sensitive mouths.

Built-in pressure sensor

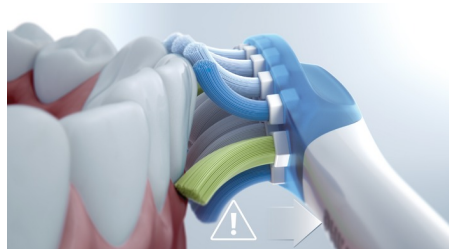

You might not notice if you're brushing too hard, but your DiamondClean 9000 will. If you need to ease up, the toothbrush will make a

pulsing sound. It's a heads up to let your brush head do the work. 7 out of 10 people found this feature helped them become a better brusher.

BrushSync mode pairing



Smart brush heads ensure you're using the right mode and intensity for the best possible clean. For an example, say you're using the W3 Premium White brush head. Your DiamondClean 9000's BrushSync technology will automatically sync your brush head with the White+ mode to help whiten your teeth.

BrushSync replacement reminder



All brush heads wear out over time, so you'll want to keep an eye on yours to make sure you're still getting a great clean. Our BrushSync technology tracks how long you've been using your brush head for, and how hard you've been brushing. The BrushSync replacement reminder on your handle and a short beep gives you a heads up when it's time to replace it.

Long-life battery



Your DiamondClean 9000 will last for up to two weeks on a single charge with regular use.

Specifications

Modes

Clean: For exceptional everyday clean

White+: To remove surface stains

Deep Clean+: For an invigorating deep clean

Gum Health: Special attention to molars

3 intensities: Low, Medium, High

Items included

Handle: 1 DiamondClean 9000

Brush heads: 1 A3 Premium All-in-One, 1 G3 Premium Gum Care

Glass charger: 1

Travel case: USB charging travel case

BrushSync Mode Pairing

C3 Premium Plaque Control: Pairs with Clean mode

G3 Premium Gum Care brush head: Pairs with Gum Health mode

W3 Premium White brush head: Pairs with White+ mode

Connectivity

Bluetooth® wireless technology: Connected brushing app

Compatibility

Android compatibility: Android phones, Bluetooth 4.0 enabled tablets

iOS compatibility: iPad 3rd Gen or higher, iPhone 4S or higher, with iOS7 or higher

Design and finishing

Color: White

Cleaning performance

Speed: Up to 62000 brush movement/min

Performance: Removes up to 20x more plaque*

Health benefits: Up to 15x healthier gums**

Whitening benefits: Up to 100% less stains***

Pressure feedback: Ring lights up in purple, Vibrates handle to alert user

TouchUp: Ensures 100% coverage

Timer: BrushPacer and SmarTimer

Smart sensor technology

Pressure sensor: Alerts when brushing too hard

Location sensor: Tracks & improves coverage

Scrubbing sensor: Guides to reduce scrubbing

3D mouth map: Highlights trouble areas

TouchUp: To address missed spots

Ease of use

Handle compatibility: Easy click-on brush heads

Brush head recognition: Syncs with optimal mode

Replacement reminder: To always ensure best results, reminder icon lights up

Battery indicator: Illuminated icon indicates battery life

Handle: Slim ergonomic design

Display: Illuminated display

Technical specifications

Battery: Rechargeable

Battery type: Lithium ION

Operating time (full to empty): 14 days*****

Power

Voltage: 110-220 V

Service

Warranty: 2-year limited warranty



* vs. a manual toothbrush

** vs. a manual toothbrush in 6 weeks

*** in lab test vs. a manual toothbrush

**** in Gum Health mode vs a manual toothbrush at 2 weeks

***** based on two periods of two-minute brushings per day on standard mode