
**PHILIPS**

*Colgate*

**SonicPro 30**

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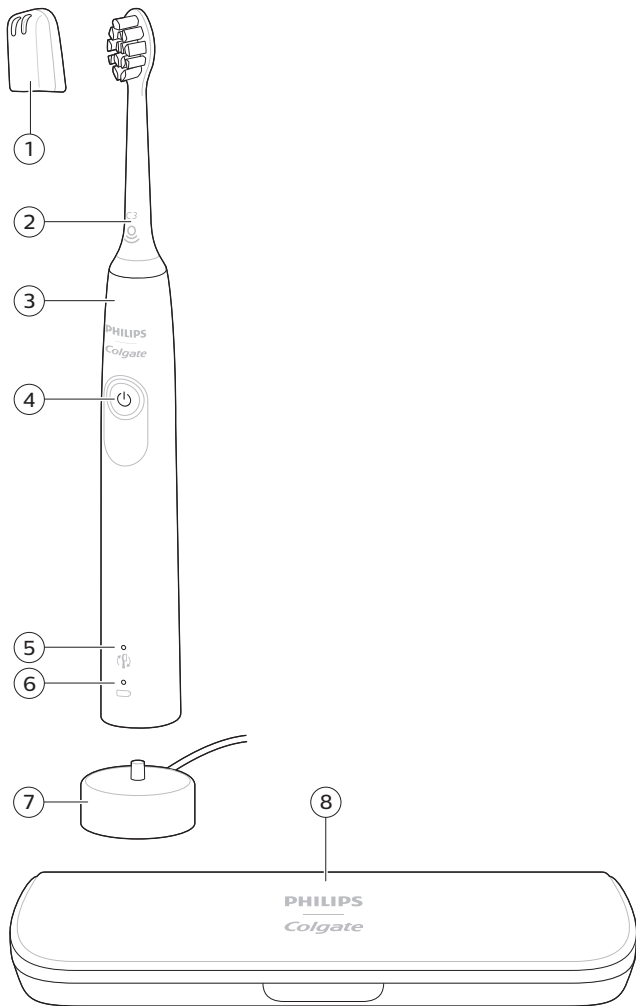## Introduction

Congratulations on your purchase and welcome to Philips Colgate! To fully benefit from the support that Philips Colgate offers, register your product at [www.philips.com/welcome](http://www.philips.com/welcome).

## Important safety information

Only use the product for its intended purpose. Read this important information carefully before you use the product and its batteries and accessories, and save it for future reference. Misuse can lead to hazards or serious injuries.

## Warnings

- Charge only with Philips WAA1001, WAA2001, UL or equivalent safety approved adapters with input ratings (100 - 240V~; 50/60Hz; 3.5W) and output ratings (5Vdc; 2.5W). For the list of adapters, please see [www.philips.com/support](http://www.philips.com/support)
- Keep the charger away from water.
- Do not place or store the appliance where it can fall or be pulled into a tub or sink.
- Do not immerse the charger in water or any other liquid.
- Make sure the handle is completely dry before connecting to the charger.
- This appliance can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.
- Children shall not play with the appliance.
- Only use original Philips Colgate accessories or consumables. Only use the charger provided with the product.

- This appliance contains no user-serviceable parts. If the appliance is damaged, stop using it and contact the Consumer Care Center in your country (see 'Warranty and support').
- Do not charge the product outdoors or near heated surfaces.
- Do not clean any part of the product in the dishwasher.
- This appliance has only been designed for cleaning teeth, gums and tongue.
- Stop using the appliance and consult your dentist/doctor if excessive bleeding occurs after use, if bleeding continues to occur after 1 week of use or if you experience discomfort or pain.
- If you have had oral or gum surgery in the previous 2 months, consult your dentist before you use this appliance.
- If you have a pacemaker or other implanted device, contact your doctor or the manufacturer of the implanted device prior to use.
- If you have medical concerns, consult your doctor before you use this appliance.
- This appliance is a personal care device and is not intended for use on multiple patients in a dental practice or institution.
- Stop using a brush head with crushed or bent bristles. Replace the brush head every 3 months or sooner if signs of wear appear.
- If your toothpaste contains peroxide, baking soda or bicarbonate (common in whitening toothpastes), thoroughly clean the brush head with soap and water after each use. This prevents possible cracking of the plastic.
- Avoid direct contact with products that contain essential oils or coconut oil. Contact may result in bristles dislodging.
- Charge, use and store the product at a temperature between 32 °F / 0 °C and 104 °F / 40 °C.

- Keep product and batteries away from fire and do not expose them to direct sunlight or high temperatures.
- If the product becomes abnormally hot or smelly, changes color or if charging takes longer than usual, stop using and charging the product and contact Philips Colgate.
- Do not place products and their batteries in microwave ovens or on induction cookers.
- Do not open, modify, pierce, damage or dismantle the product or battery to prevent batteries from heating up or releasing toxic or hazardous substances. Do not short-circuit, overcharge or reverse charge batteries.
- If batteries are damaged or leaking, avoid contact with the skin or eyes. If this occurs, immediately rinse well with water and seek medical care.

### **Electromagnetic fields (EMF)**

This Philips Colgate appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

### **Intended use**

SonicPro toothbrushes are intended to remove adherent plaque and food debris from the teeth to reduce tooth decay and improve and maintain oral health. SonicPro power toothbrushes are intended for consumer home use. Use by children should be with adult supervision.

### **Your SonicPro toothbrush**

- 1 Hygienic travel cap
- 2 Brush head
- 3 Handle
- 4 Power on/off button
- 5 Brush head replacement reminder
- 6 Battery level indicator
- 7 USB charger (wall adapter not included)
- 8 Travel case (Specific models only)


Note: Wall adapter not included. Only use 5V IPX4 adapter (Type A).

Note: Brush heads may vary.

## Getting started

### Attaching the brush head

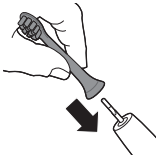- 1 Align the brush head so the bristles point in the same direction as the front of the handle.

- 2 Firmly press the brush head down onto the metal shaft until it stops.

Note: It is normal to see a slight gap between the brush head and the handle. This allows the brush head to vibrate properly.

### Charging your SonicPro

- 1 Plug the USB charging cord into a wall adapter that is connected to a electrical outlet.
- 2 Place the handle on the charger.
  - The light of the battery level indicator illuminates and the handle beeps twice. This indicates that the toothbrush is charging.

Note: It can take up to 24 hours to fully charge the battery, but you can use the SonicPro toothbrush before it is fully charged.


### Using your SonicPro

If you are using your SonicPro toothbrush for the first time, it is normal to feel more vibration than when using a non-electric toothbrush. It is common for first time users to apply too much pressure. Only apply gentle pressure and let the toothbrush do the brushing for you. Please follow brushing instructions below for the best experience. To help you transition

to your new SonicPro power toothbrush, it comes with the EasyStart feature turned on. This feature gently increases the power over the first 14 brushings to help you adapt to brushing with a Philips Colgate SonicPro toothbrush.

### Brushing instructions

- 1 Wet the bristles and apply a small amount of toothpaste.
- 2 Place the toothbrush bristles against the teeth at a slight angle (45 degrees), pressing gently to make the bristles reach the gumline or slightly beneath the gumline.

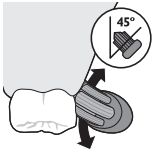Note: Keep the center of the brush in contact with the teeth at all times.

- 3 Press the power on/off button to turn on the SonicPro.
- 4 Gently keep the bristles placed on the teeth and in the gumline. Brush your teeth with small back and forth motion so the bristles reach between the teeth. Continue this motion throughout your brushing cycle.

Note: The bristles should slightly flare. It is not recommended to scrub as you would with a manual toothbrush.


- 5 To clean the inside surfaces of the front teeth, tilt the brush handle semi-upright and make several vertical overlapping brushing strokes on each tooth.

Note: To make sure you brush evenly throughout the mouth, divide your mouth into 4 sections using the Quadpacer feature (see chapter 'Features').

- 6 Begin brushing in section 1 (outside upper teeth) and brush for 30 seconds before you move to section 2 (inside upper teeth). Continue brushing in section 3 (outside lower teeth) and brush for 30 seconds before you move to section 4 (inside lower teeth).

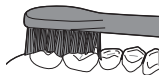- 7 After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs. You can also brush your tongue, with the toothbrush turned on or off, as you prefer. Your SonicPro toothbrush should be safe to use on braces (brush heads wear out sooner when used on braces) and dental restorations (fillings, crowns, veneers) if they are properly adhered and not compromised. If a problem occurs, please follow up with a dental professional.

Note: When the SonicPro toothbrush is used in clinical studies, the handle should be fully charged and the EasyStart feature deactivated.

## Charging and battery status

- 1 Attach the USB charging cord to a wall adapter, plug the wall adapter into an electrical wall outlet.
- 2 Place the handle on the charger.

Note: The handle will emit 2 short beeps to confirm that the handle is placed properly on the charger and is charging.

Note: This SonicPro toothbrush is designed to provide at least 28 brushing sessions, each session lasting 2 minutes (14 days if used twice per day). The battery level indicator shows the battery status when you complete the 2-minute brushing session, when you pause the handle, or while it is charging.

### Battery status (when handle is on a working charger)

Battery status	Battery indicator
Full	Solid green (turns off after 30 seconds)
Charging	Flashing green

## Battery status (When handle is not placed on charger)

Battery status	Battery indicator	Sound
Full	Solid green	-
Partially full	Flashing green	-
Low	Flashing orange	3 beeps
Empty	Flashing orange	2 sets of 5 beeps

Note: To save energy, the battery light will turn off when not in use.

Note: When the battery is completely empty, the SonicPro toothbrush turns off. Place the SonicPro toothbrush on the charger to charge it.

Note: To keep the battery fully charged at all times, you may keep your SonicPro toothbrush on the charger when not in use.

## Features

- Intensities
- EasyStart
- Brush head replacement reminder
- Pressure sensor
- Quadpacer
- SmarTimer

## Intensities

Your power toothbrush gives you the option to choose between low and high intensity.

Note: When you use the toothbrush for the first time, the default setting is the low intensity.

- Press the power button once to turn on the toothbrush.
- Press a second time **within 2 seconds** to change the intensity.
- Press a third time **within 2 seconds** to pause.

**After 2 seconds** of brushing, pressing the power button again will also pause the toothbrush.

## EasyStart


This SonicPro model comes with the EasyStart feature activated. The EasyStart feature gently increases the power over the first 14 brushings to help you get used to the brushing with the SonicPro.

To deactivate EasyStart see 'Activating or deactivating features'.

## Brush head replacement reminder

Your SonicPro is equipped with BrushSync technology that tracks the wear of your brush head.

- 1 When attaching a new smart brush head for the first time the handle recognizes that you have a Philips Colgate smart brush head with BrushSync technology and begins tracking brush head wear.
- 2 Over time, based on the pressure you apply and the amount of time used, the handle will track brush head wear in order to determine the optimal time to change your brush head. This feature gives you the guarantee for the best cleaning and care of your teeth.
- 3 When the brush head replacement reminder light lights up amber, you should replace your brush head.

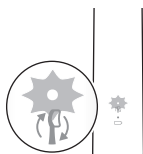Note: All smart brush heads come with the brush head replacement reminder feature activated. To deactivate the brush head replacement reminder (see 'Activating or deactivating features').

## Pressure Sensor

Your SonicPro is equipped with an advanced sensor that measures the pressure you apply while brushing. If you apply too much pressure, the toothbrush will provide immediate feedback to indicate that you need to reduce the pressure. This feedback is given by a change in vibration and therefore brushing sensation and sound. Additionally, the brush head replacement reminder light will flash amber.

To deactivate the pressure sensor see 'Activating or deactivating features'.

## QuadPacer

The QuadPacer is an interval timer that has a short beep and pause to remind you to brush the 4 sections of your mouth evenly and thoroughly. This SonicPro toothbrush comes with the QuadPacer feature activated.

## SmarTimer

The SmarTimer indicates that your brushing cycle is complete by automatically turning off the toothbrush at the end of the brushing cycle. Dental professionals recommend brushing no less than 2 minutes twice a day.

Note: If you press the power on/off button after you started the brushing cycle, the toothbrush will pause. After a pause of 30 seconds, the SmarTimer resets.

## Activating or deactivating features

You can activate or deactivate the following features of your toothbrush:

- EasyStart
- Brush head replacement reminder
- Pressure sensor

To activate or deactivate these features, follow the instructions below:

### EasyStart

- 1 Put the handle on the plugged-in charger.
- 2 Press and hold the power on/off button while the handle remains on the charger.
- 3 Keep the power on/off button pressed until you hear a single short beep (after 3 seconds).
- 4 Release the power on/off button.
  - Triple tone of low-med-high means the EasyStart feature has been activated. The battery level indicator light will also blink green 2 times to confirm activation.

- Triple tone of high-med-low means the EasyStart feature has been deactivated. The battery level indicator light will blink amber 1 time to confirm deactivation.

Note: To achieve clinical efficacy, EasyStart needs to be deactivated.

## Brush head replacement reminder

Note: You can only activate or deactivate the brush head replacement reminder when a smart brush head is on the handle.

- 1 Put the handle on the plugged-in charger.
- 2 Press and hold the power on/off button while the handle remains on the charger.
- 3 Keep the power on/off button pressed until you hear a series of two short beeps (after 5 seconds).
- 4 Release the power on/off button.
  - Triple tone of low-med-high means the brush head replacement reminder feature has been activated. The battery level indicator light will also blink green 2 times to confirm activation.
  - Triple tone of high-med-low means the brush head replacement reminder feature has been deactivated. The battery level indicator light will blink amber 1 time to confirm deactivation.

## Pressure sensor

- 1 Put the handle on the plugged-in charger.
- 2 Press and hold the power on/off button while the handle remains on the charger.
- 3 Keep the power on/off button pressed until you hear a series of three short beeps (after 7 seconds).
- 4 Release the power on/off button.
  - Triple tone of low-med-high means the pressure sensor feature has been activated. The battery level indicator light will also blink green 2 times to confirm activation.

- Triple tone of high-med-low means the pressure sensor feature has been deactivated. The battery level indicator light will blink amber 1 time to confirm deactivation.

If you continue holding the power on/off button after the three short beeps, the activate/deactivate sequence repeats.

## Cleaning

You should clean your SonicPro toothbrush regularly to remove toothpaste and other residue. Failure to clean your SonicPro toothbrush may result in an unhygienic product and damage to it may occur.

**Warning: Do not clean product or accessories with sharp objects or in dishwasher, microwave, with chemicals or in boiling hot water**

### Toothbrush handle


- 1 Remove the brush head and rinse the metal shaft area with warm water. Make sure you remove any residual toothpaste.

**Caution: Do not push on the rubber seal on the metal shaft with sharp objects, as this may cause damage.**

- 2 Wipe the entire surface of the handle with a damp cloth.

**Note: Do not tap the handle on the sink to remove excess water.**

### Brush head

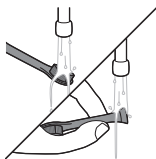- 1 Rinse the brush head and bristles after each use.
- 2 Remove the brush head from the handle and rinse the brush head connection with warm water at least once a week. Rinse the travel cap as often as needed.

### Charger

- 1 Unplug the charger before you clean it.
- 2 Wipe the surface of the charger with a damp cloth.

## Travel case

Rinse with warm water and use a damp cloth to clean inside of travel case

## Storage

If you are not going to use the product for a long time, unplug it, clean it (see chapter 'Cleaning') and store it in a cool and dry place away from direct sunlight.

## Replacement

### Brush head

Replace SonicPro brush heads every 3 months to achieve optimal results. Use only SonicPro replacement brush heads.

## Locating the model number

Look on the bottom of the SonicPro toothbrush handle for the model number (PC08XX).

## Warranty and support

If you need information or support, please visit **[www.philips.com/support](http://www.philips.com/support)** or read the international warranty leaflet.

## Warranty restrictions

The terms of the international warranty do not cover the following:

- Brush heads.
- Damage caused by use of unauthorized replacement parts.
- Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.

## Recycling

- This symbol means that electrical products and batteries shall not be disposed of with normal household waste.
- Follow your country's rules for the separate collection of electrical products and batteries.

## Removal of built-in rechargeable battery

The built-in rechargeable battery must only be removed by a qualified professional when the appliance is discarded. Before removing the battery, make sure that the battery is completely empty.

**Take any necessary safety precautions when you handle tools to open the appliance and when you dispose of the rechargeable battery. Be sure to protect your eyes, hands, fingers, and the surface on which you work.**

**When you handle batteries, make sure that your hands, the product and the batteries are dry.**

**To avoid accidental short-circuiting of batteries after removal, do not let battery terminals come into contact with metal objects (e.g. coins, hairpins, rings). Do not wrap batteries in aluminum foil. Tape battery terminals or put batteries in a plastic bag before you discard them.**

**Incorrectly discarding batteries may bring risk/damages to human health and environment.**


## Removing the rechargeable battery

To remove the rechargeable battery, you need a towel or cloth, a hammer and a flat-head (standard) screwdriver.

- 1 To deplete the rechargeable battery of any charge, remove the handle from the charger, turn on the SonicPro and let it run until it stops. Repeat this step until you can no longer turn on the SonicPro.
- 2 Remove and discard the brush head. Cover the entire handle with a towel or cloth.


- 3 Hold the top of the handle with one hand and strike the handle housing 0.5 inch above the bottom end. Strike firmly with a hammer on all 4 sides to eject the end cap.

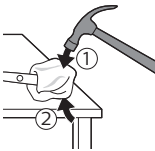
Note: You may have to hit on the end several times to break the internal snap connections.


- 4 Remove the end cap from the toothbrush handle. If the end cap does not release easily from the housing, repeat step 3 until the end cap is released.

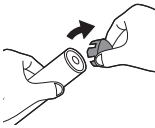
- 5 Holding the handle upside down, press the shaft down on a hard surface. If the internal components do not easily release from the housing, repeat step 3 until the internal components are released.

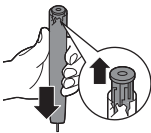- 6 Wedge the screwdriver between the battery and the white frame at the bottom of the internal components. Then pry the screwdriver away from the battery to break the bottom of the white frame.

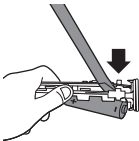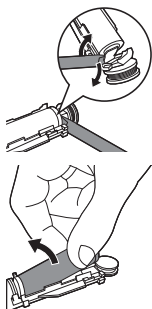
- 7 Insert the screwdriver between the bottom of the battery and the white frame to break the metal tab connecting the battery to the green printed circuit board. This will release the bottom end of the battery from the frame.

- 8 Grab the battery and pull it away from the internal components to break the second metal battery tab.

**Caution: Be aware of the sharp edges of the battery tabs so as to avoid injury to your fingers.**

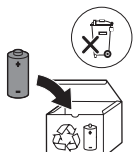- 9 Cover the battery contacts with tape to prevent any electrical short from residual battery charge. The rechargeable battery can now be recycled and the rest of the product discarded appropriately.

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© 2021 Koninklijke Philips N.V. (KPNV) All rights reserved. Philips and the Philips shield are trademarks of KPNV.
Colgate and the Colgate logo are the registered trademarks of Colgate-Palmolive Company.