
PHILIPS





Colgate

SonicPro

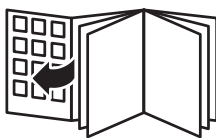1

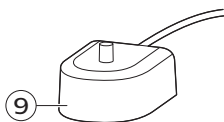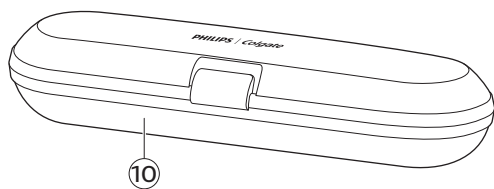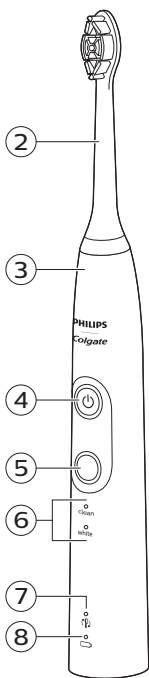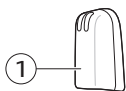English 6
Español 26

Introduction

Congratulations on your new Philips Colgate SonicPro power toothbrush! Superior plaque removal, whiter teeth and healthier gums are at your fingertips. Using a combination of gentle sonic technology and clinically developed and proven features, you can be confident that you're getting the very best clean, every time. Further support and product registration are available to you at: www.philips.com/welcome.

Welcome to the Philips Colgate SonicPro family!

Intended Use

Philips Colgate SonicPro power toothbrushes are intended to remove adherent plaque and food debris from the teeth to reduce tooth decay and improve and maintain oral health. They are intended for consumer home use. Use by children should be with adult supervision.

Important safety information

Only use the product for its intended purpose. Read this important information carefully before you use the product and its batteries and accessories, and save it for future reference. Misuse can lead to hazards or serious injuries.

Danger

- Keep the charger away from water. Do not place or store the appliance where it can fall or be pulled into a tub or sink.
- Do not immerse the charger in water or any other liquid.
- After cleaning, make sure the charger is completely dry before you connect it to the wall socket.

Warning

- This appliance can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.
- Children shall not play with the appliance.
- If the appliance is damaged in any way (brush head, toothbrush handle and charger), stop using it.
- The mains cord cannot be replaced. If the mains cord is damaged, discard the charger.
- Always replace the charger with one of the original type in order to avoid a hazard.
- This appliance contains no user-serviceable parts. If the appliance is damaged, contact the Consumer Care Center in your country (see 'Warranty and support').
- Keep the cord away from heated surfaces.

Caution

- Do not use the charger outdoors or near heated surfaces.
- Do not clean the brush head, handle, or charger in the dishwasher.
- If you have had oral or gum surgery in the previous 2 months, consult your dentist before you use this appliance.
- Consult your dentist if excessive bleeding occurs after using this appliance or if bleeding continues to occur after 1 week of use. Also consult your dentist if you experience discomfort or pain when you use this appliance.
- This Philips Colgate SonicPro toothbrush complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician or the device manufacturer of the implanted device prior to use.
- If you have medical concerns, consult your doctor before you use this appliance.
- This appliance has only been designed for cleaning teeth, gums and tongue. Do not use it for any other purpose. Stop using the appliance and contact your doctor if you experience any discomfort or pain.
- This appliance is a personal care device and is not intended for use on multiple patients in a dental practice or institution.

- Do not use other brush heads than the ones recommended by Philips Colgate SonicPro.
- Stop using a brush head with crushed or bent bristles. Replace the brush head every 3 months or sooner if signs of wear appear.
- Avoid direct contact with products that contain essential oils or coconut oil. Contact may result in bristles dislodging.

Battery safety instructions

- Only use original Philips Colgate SonicPro accessories or consumables. Only use detachable supply unit HX6100.
- Charge, use and store the product at a temperature between 0 °C and 40 °C.
- Keep product and batteries away from fire and do not expose them to direct sunlight or high temperatures.
- If the product becomes abnormally hot or smelly, changes color or if charging takes longer than usual, stop using and charging the product and contact Philips.
- If the product becomes abnormally hot, gives off an abnormal smell, changes color or if charging takes much longer than usual, stop using and charging the product and contact Philips directly via a service center or warranty line at www.philips.com/support, or toll free at 800-504-6200.
- Do not place products and their batteries in microwave ovens or on induction cookers.

- Do not open, modify, pierce, damage or dismantle the product or battery to prevent batteries from heating up or releasing toxic or hazardous substances. Do not short-circuit, overcharge or reverse charge batteries.
- If batteries are damaged or leaking, avoid contact with the skin or eyes. If this occurs, immediately rinse well with water and seek medical care.

Electromagnetic fields (EMF)

This Philips Colgate SonicPro appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

Radio Equipment Directive

Hereby, Philips declares that electric toothbrushes with radio frequency interface (Bluetooth) are in compliance with Directive 2014/53/EU.

The full text of the EU declaration of conformity is available at the following internet address: **www.philips.com/support**

The radio frequency interface in this product operates at 13.56 MHz.

The maximum RF power transmitted by the appliance is 30.16 dBm.

Your SonicPro (Fig. 1)

- 1 Hygienic travel cap
- 2 Smart brush head(s)*
- 3 Handle
- 4 Power on/off button
- 5 Mode/intensity button*
- 6 Mode light
- 7 Brush head replacement reminder light
- 8 Battery indication
- 9 Charging base
- 10 Travel case (specific types only)*

*Note: The content of the box may vary based on the model purchased.

Brush heads

Your Philips Colgate SonicPro comes with one or more brush heads which are designed to deliver superior results for your oral care needs.

BrushSync Technology

BrushSync technology enables your brush heads to communicate with your handle using a microchip. This enables the brush head replacement reminder feature.

Brushing modes

Your power toothbrush is equipped with two brushing modes. To toggle between modes press the mode button located below the on/off button before switching the toothbrush on.

Clean mode

Clean is a 2-minute mode. It gives you a thorough and complete clean and is recommended to be used with the Deep Clean brush head.

White mode

White is a 2-minute 30-seconds mode, recommended to be used with the Whitening brush head.

Intensities

Your power toothbrush gives you the option to choose between low and high intensity.

Note: When you use the toothbrush for the first time, the default setting is the low intensity.

The intensity can be changed by pressing the mode/intensity button while brushing. The intensity cannot be changed as long as the handle is powered off or paused.

Battery status and charging

Charging

- 1 Put the plug of the charger in an electrical outlet.
- 2 Place the handle on the charger.
 - The charger will emit 2 short beeps to confirm that the handle is placed properly.

- The flashing light of the battery level indicator shows that the toothbrush is charging.
- When handle is fully charged on the charger, the battery light shows solid green for 30 seconds and switches off.

Note: Your toothbrush comes pre-charged for first use. After first use, charge for at least 24 hours.


Battery status (when handle is not on charger)

When removing the SonicPro from the charger, the battery light at the bottom of the toothbrush will indicate the status of the battery life.

- Solid green LED: full battery
- Flashing green LED: medium battery
- Flashing amber LED and three beeps: low battery
- Flashing amber LED and two sets of five beeps: No brushing sessions left (charge toothbrush)

Using your SonicPro


Brushing instructions

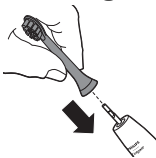- 1 Firmly press the brush head down onto the metal shaft of the handle until it stops. The brush head bristles should be facing the same direction as the front of the handle.

Note: It is normal to see a slight gap between the brush head and the handle.


14 English

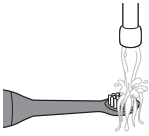2 Wet the bristles.


3 Apply a small amount of toothpaste on the bristles.

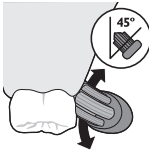4 Place the toothbrush bristles against the teeth at a slight angle (45 degrees), pressing firmly to make the bristles reach the gumline or slightly beneath the gumline.

Note: Keep the center of the brush in contact with the teeth at all times.

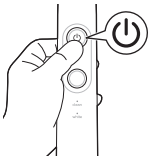5 Press the power on/off button to turn on the toothbrush.



6 Gently keep the bristles placed on the teeth and in the gumline. Brush your teeth with small back and forth motion so the bristles reach between the teeth.

Note: The bristles should slightly flare. Do not scrub.

- 7** To clean the inside surfaces of the front teeth, tilt the brush handle semi-upright and make several vertical overlapping brushing strokes on each tooth.

Note: A change in vibration (and a slight change in sound) alerts you when you apply too much pressure while brushing.

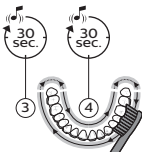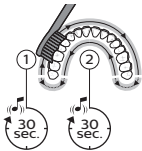- 8** To make sure you brush evenly throughout the mouth, divide the mouth into 4 sections using the Quadpacer feature. The Quadpacer feature beeps every 30 seconds over 2 minutes.

- 9** Brush for the full 2 minutes following the Quadpacer.
- In Clean mode, after 2 minutes the handle will automatically shut off.
 - In White mode, after 2 minutes the handle will continue to brush for an additional 30 seconds so you can polish and brighten your front teeth.
 - In Gum Care mode, after 2 minutes the handle will continue to brush for an additional minute so you can gently stimulate and massage your gums.

After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs.

Note: Your SonicPro is safe to use on:

- Braces (brush heads wear out sooner when used on braces)
- Dental restorations (fillings, crowns, veneers)

Note: When the Philips Colgate SonicPro toothbrush is used in clinical studies, it must be used in the Clean mode (for plaque and gum health claims) or in White mode (for whitening claims) at high intensity with EasyStart turned off.

Features

- EasyStart
- Brush head replacement reminder
- Pressure sensor
- Quadpacer


EasyStart

The EasyStart feature gently increases the power over the first 14 brushings to help you get used to brushing with the SonicPro. The EasyStart feature works in all brushing modes.

Note: This SonicPro comes with the EasyStart feature deactivated. To activate EasyStart see 'Activating or deactivating features'.

Note: To achieve clinical efficacy, EasyStart needs to be deactivated.

Brush head replacement reminder

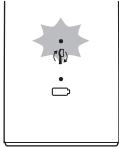Your SonicPro is equipped with BrushSync technology that tracks the wear of your brush head.

- 1** When attaching a new smart brush head for the first time the brush head replacement reminder light blinks green three times. This confirms you have a Philips brush head with BrushSync technology.
- 2** Over time, based on the pressure you apply and the amount of time used, the handle will track brush head wear in order to determine the optimal time to change your brush head. This feature gives you the guarantee for the best cleaning and care of your teeth.
- 3** When the brush head replacement reminder light lights up amber, you should replace your brush head.

Note: This SonicPro comes with the brush head replacement reminder activated. To deactivate the brush head replacement reminder see 'Activating or deactivating features'.

Note: The brush head replacement reminder only works in combination with SonicPro smart brush heads (see 'Brush heads').

Note: You can only activate or deactivate the brush head replacement reminder when a smart brush head is on the handle.

Note: New smart brush heads come with the brush head replacement reminder activated.

Pressure sensor

Your SonicPro is equipped with an advanced sensor that measures the pressure you apply while brushing. If you apply too much pressure, the toothbrush will provide immediate feedback to indicate that you need to reduce the pressure. This feedback is given by a change in vibration and therefore brushing sensation.

Note: This SonicPro comes with the pressure sensor activated. To deactivate the pressure sensor see 'Activating or deactivating features'.

Quadpacer

The Quadpacer is an interval timer that has a short beep and pause to remind you to brush the different sections of your mouth. Depending on the brushing mode you have selected (see 'Brushing modes'), the Quadpacer beeps at different intervals during the brushing cycle.


Activating or deactivating features

You can activate or deactivate the following features of your toothbrush:

- EasyStart
- Brush head replacement reminder
- Pressure sensor

To activate or deactivate these features, follow the instructions below:

Step 1: Place the handle on the charger.

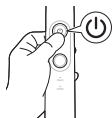
Step 2: Press and hold the power button for:

EasyStart	Brush head replacement reminder	Pressure sensor
Up to 3 seconds	Up to 5 seconds	Up to 7 seconds
↓	↓	↓


Step 3: Let go power button when you hear

1 beep	2 beeps	3 beeps
↓	↓	↓

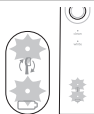If you see the brush head replacement light and battery indication light blink green 2 times, and hear 3 tones from low to high, then the feature has been activated.

OR

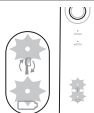
If you see the brush head replacement light and battery indication light blink amber once, and hear 3 tones from high to low, then the feature has been deactivated.

Cleaning

Warning: Do not clean the brush head, travel cap, handle, charger or travel case in the dishwasher.

Note: Do not use isopropyl rubbing alcohol, vinegar or bleach to clean the handle, brush head or charger, as this may cause discoloration.

Toothbrush handle

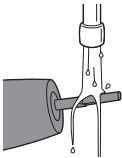
- 1 Remove the brush head from the handle and rinse the metal shaft area with warm water. Make sure you remove any residual toothpaste.

Caution: Do not push on the rubber seal at the top of the handle with sharp objects, as this may cause damage.

- 2 Wipe the entire surface of the handle with a damp cloth.


Brush head

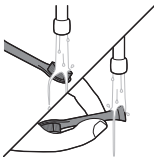- 1 Remove the brush head from the handle and rinse it thoroughly.
- 2 Rinse the brush head connection with warm water at least once a week.

Charger

- 1 Unplug the charger before you clean it.

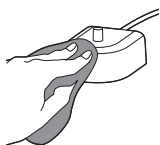- 2 Wipe the surface of the charger with a damp cloth. Make sure you remove any toothpaste or other residue from the charger surfaces.

Storage

If you are not going to use the product for an extended period of time, unplug the charger from the electrical outlet, clean it and store it in a cool and dry place away from direct sunlight.

Locating the model number

Look on the bottom of the SonicPro toothbrush handle for the model number (HX68xx).

Recycling


- This symbol means that electrical products and batteries shall not be disposed of with normal household waste.
- Follow your country's rules for the separate collection of electrical products and batteries.

Removing the rechargeable battery

Warning: Only remove the rechargeable battery when you discard the appliance. Make sure the battery is completely empty when you remove it. Any attempt to remove the battery will void the warranty.

- To remove the rechargeable battery, you need a towel or cloth, a hammer and a flat-head (standard) screwdriver. Observe basic safety precautions when you follow the procedure outlined below. Be sure to protect your eyes, hands, fingers, and the surface on which you work.


- 1 To deplete the rechargeable battery of any charge, remove the handle from the charger, turn on the toothbrush and let it run until it stops. Repeat this step until you can no longer turn on the toothbrush.
- 2 Remove and discard the brush head. Cover the entire handle with a towel or cloth.

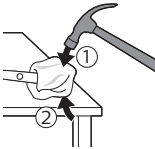
- 3 Hold the top of the handle with one hand and strike the handle housing 0.5 inch above the bottom end. Strike firmly with a hammer on all 4 sides to eject the end cap.


Note: You may have to hit on the end several times to break the internal snap connections.

- 4 Remove the end cap from the toothbrush handle. If the end cap does not release easily from the housing, repeat step 3 until the end cap is released.

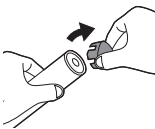
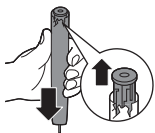
5 Holding the handle upside down, press the shaft down on a hard surface. If the internal components do not easily release from the housing, repeat step 3 until the internal components are released.

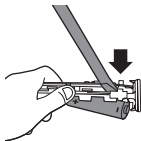6 Wedge the screwdriver between the battery and the black frame at the bottom of the internal components. Then pry the screwdriver away from the battery to break the bottom of the black frame.

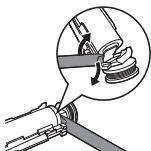
7 Insert the screwdriver between the bottom of the battery and the black frame to break the metal tab connecting the battery to the green printed circuit board. This will release the bottom end of the battery from the frame.

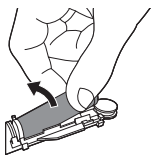8 Grab the battery and pull it away from the internal components to break the second metal battery tab.

Caution: Be aware of the sharp edges of the battery tabs so as to avoid injury to your fingers.

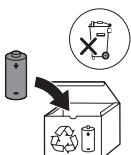9 Cover the battery contacts with tape to prevent any electrical short from residual battery charge. The rechargeable battery can now be recycled and the rest of the product discarded appropriately.

Warranty and support

If you need information or support, contact our Consumer Care Center at 800-504-6200, visit **www.philips.com/support**, or read the international warranty leaflet.

Warranty restrictions

The terms of the international warranty do not cover the following:

- Brush heads.
- Damage caused by use of unauthorized replacement parts.
- Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.

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© 2021 Koninklijke Philips N.V. (KPNV) All rights reserved. Philips and the Philips shield are trademarks of KPNV.
Colgate and the Colgate logo are the registered trademarks of Colgate-Palmolive Company.