PHILIPS



One
HY1100

300005273224 (04/26/2022)
Philips Oral Healthcare LLC
22100 Bothell Everett Highway Bothell,
WA 98021
www.philips.com/one

©2022 Koninklijke Philips N.V. (KPNV).
All rights reserved.
Philips and the Philips shield are trademarks
of KPNV.
Sonicare and the Sonicare logo are trademarks
of Philips Oral Healthcare, LLC and/or KPNV.

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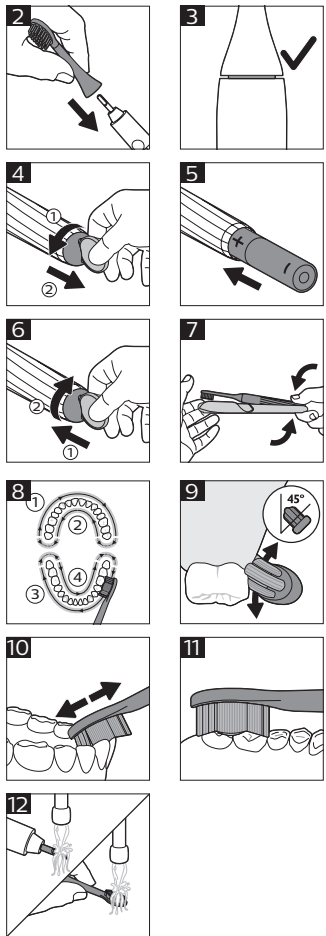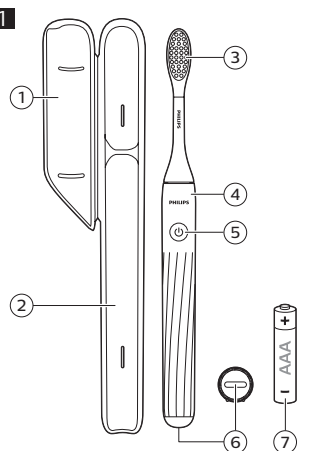IMPORTANT SAFEGUARDS

READ ALL INSTRUCTIONS BEFORE USE

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

DANGERS

- To reduce the risk of electrocution:
- Do not place or store product where it can fall or be pulled into a tub or sink.
- Do not place in or drop into water or other liquid.
- Do not use while bathing.

WARNINGS

- To reduce the risk of burns, electrocution, fire or injury to persons:
- Discontinue use if this product appears damaged in any way (brush head and/or handle). This appliance contains no user-serviceable parts. Refer to chapter 'Warranty and support' for your Philips One no longer works properly or needs repair.

- This appliance has only been designed for cleaning teeth, gums and tongue only. Use this product only for its intended use as described in this user manual.
- Stop using the appliance and consult your dentist/doctor if excessive bleeding occurs after use, if bleeding continues to occur after 1 week of use or if you experience discomfort or pain.
- This appliance can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.
- Children shall not play with the appliance. Only use original Philips accessories or consumables.
- Do not clean any part of the product in the dishwasher.
- Stop using a brush head with crushed or bent bristles. Replace the brush head every 3 months or sooner if signs of wear appear.
- If your toothpaste contains peroxide, baking soda or bicarbonate (common in whitening toothpastes), thoroughly clean the brush head with soap and water after each use. This prevents possible cracking of the plastic.
- This appliance is a personal care device and is not intended for use on multiple patients in a dental practice or institution.
- Avoid direct contact with products that contain essential oils or coconut oil. Contact may result in bristles dislodging.
- Accessories supplied may vary for different products.

MEDICAL WARNINGS

- Consult your dental professional before you use this product if you have had oral or gum surgery in the previous 2 months.
- Contact your dental professional if excessive bleeding occurs after using this product or bleeding continues to occur after 1 week of use.
- If you have a pacemaker or other implanted device contact your physician or the device manufacturer prior to use.
- This product has been tested and is compliant with safety standards for electromagnetic devices.
- Consult your physician prior to using the product if you have medical concerns.

Battery safety instructions

- Only use the product for its intended purpose. Read this information carefully before you use the product and its batteries and accessories, and save it for future reference. Misuse can lead to hazards or serious injuries.
- Only operate this product with type of disposable battery specified in the user manual.
- Always keep coin/button cell batteries out of reach of children. Swallowing these batteries can result in serious injury or death.
- Do not recharge non-rechargeable batteries.
- Do not burn products and their batteries. Keep product and batteries away from fire and do not expose them to direct sunlight or high temperatures. Batteries may explode if overheated.
- If the product becomes abnormally hot or smelly, or changes color, stop using the product and contact Philips support directly or Philips service center at www.philips.com/support or toll free at 1-800-243-3050.
- Do not place products and their batteries in microwave ovens or on induction cookers.
- Do not modify, pierce, damage or dismantle the product or battery to prevent batteries from heating up or releasing toxic or hazardous substances. Do not short-circuit, overcharge or reverse charge batteries.
- If batteries are damaged or leaking, avoid contact with the skin or eye. If this occurs, immediately rinse well with water and seek medical care.
- When you handle batteries, make sure that your hands, the product and the batteries are dry.

- Insert batteries with the + and - poles pointing in the direction indicated in the battery compartment or holder.
- Do not mix different brands and types of batteries, do not mix new and used batteries and do not use batteries with different date codes.
- Remove batteries from the product if you are not going to use it for some time.
- To avoid accidental short-circuiting of terminals come into contact with metal objects (e.g. coins, hairpins, rings). Do not wrap batteries in aluminium foil. Tape battery terminals or put batteries in a plastic bag before you discard them.
- Do not leave empty disposable batteries inside the product.

SAVE THESE INSTRUCTIONS

Electromagnetic fields (EMF)

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

Introduction

Congratulations on your purchase and welcome to Philips! To fully benefit from the support that Philips offers, register your product at www.philips.com/welcome.

Intended use

Philips One power toothbrushes are intended to remove adherent plaque and food debris from the teeth to reduce tooth decay and improve and maintain oral health. Philips One power toothbrushes are intended for consumer home use. Use by children should be with adult supervision.

Your Philips One (Fig. 1)

- 1 Travel case door
- 2 Travel case
- 3 Brush head
- 4 Handle
- 5 Power button
- 6 Removable end cap with coin slot
- 7 Replaceable AAA-battery

Note: The content of the box may vary based on the model purchased.

Getting started

Attaching the brush head

- 1 Align the brush head so the bristles point in the same direction as the front of the handle. Press the brush head down onto the metal shaft (Fig. 2).

Note: It is normal to see a slight gap between the brush head and the handle (Fig. 3). This allows the brush head to vibrate properly.

Battery replacement

Your power toothbrush is powered by a replaceable 1.5V Alkaline AAA battery. A single AAA battery should provide enough power for two 2-minute brushing sessions every day for 3 months. If used more often, the battery will drain before 3 months.

Note: Make sure your hands and the product are dry when you insert the battery.

- 1 Insert a coin into the coin slot at the bottom of the handle and turn counterclockwise (Fig. 4). Gently shake the battery out of the handle.
- 2 Insert a fresh 1.5V Alkaline AAA battery into the battery compartment (Fig. 5). Make sure the + and - poles of the battery point in the right direction.
- 3 Place the removable end cap back on the handle and use a coin to turn clockwise until it is secure and the lock indicator markings on the handle and end cap line up (Fig. 6).

Note: If the indication on the handle and end cap do not line up, it may result in your toothbrush not working

Travel case

Your Philips One comes with a travel case to store your power toothbrush when traveling.

- 1 Place your Philips One toothbrush with the bristles facing downwards in the travel case and snap closed.
- 2 To remove, open and press down on the bottom of the handle (Fig. 7).

Using your Philips One toothbrush

If you are using this toothbrush for the first time, it is normal to feel more vibration than when using a manual toothbrush. The Philips One uses microvibrations to help gently clean your teeth. Please follow the brushing instructions below for the best experience.

Philips One comes with a SmarTimer and QuadPacer to help you brush evenly throughout your mouth.

SmarTimer

The SmarTimer indicates that your 2-minute brushing cycle is complete by automatically turning off the toothbrush at the end of 2 minutes.

Note: Dental professionals recommend brushing no less than 2 minutes twice a day.

QuadPacer

The QuadPacer is an interval timer that beeps every 30 seconds to remind you to brush all teeth evenly by dividing your mouth into 4 sections (Fig. 8).

Begin brushing in section 1 (outside upper teeth) and brush for 30 seconds before you move to section 2 (inside upper teeth). Continue brushing in section 3 (outside lower teeth) and brush for 30 seconds before you move to section 4 (inside lower teeth).

Note: If you press the power button after you started the brushing cycle, the toothbrush will pause. After a pause of 30 seconds, the SmarTimer and QuadPacer resets.

Brushing instructions

Please follow the brushing instructions for the best experience.

- 1 Wet the bristles and apply a small amount of toothpaste.
- 2 Place the toothbrush bristles against the teeth at a slight angle (45 degrees) (Fig. 9), pressing gently to make the bristles reach the gumline or slightly beneath the gumline.
- 3 Press the power button to turn brush on.
- 4 Brush your teeth with back and forth strokes so that the bristles reach all surfaces of the teeth. Change to another section when the QuadPacer beeps. Continue this motion throughout the 2-minute brushing cycle.

Note: The bristles should slightly flare.

- 5 To clean the inside surfaces of the front teeth, tilt the brush vertically and make several up-and-down strokes (Fig. 10).
- 6 After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs (Fig. 11). You can also brush your tongue, with the toothbrush turned on or off, as you prefer.

Your power toothbrush should be safe to use on braces (brush heads wear out sooner when used on braces) and dental restorations (fillings, crowns, veneers) if they are properly adhered and not compromised. If a problem occurs, please follow up with a dental professional.

Battery status

This product does not come with a battery status indicator. As the battery gets weaker, you may experience a drop in vibration power. This is normal. While the microvibrations are still effective, you can increase the power by replacing your battery.

Cleaning

You should clean your toothbrush and accessories regularly to remove toothpaste and other residue. Failure to clean your product may result in an unhygienic product and damage to it may occur.

Warning: Do not clean product or accessories in dishwasher, microwave, with chemicals or with boiling hot water.

Toothbrush handle

- 1 Remove the brush head and rinse the metal shaft area with warm water (Fig. 12). Make sure you remove any residual toothpaste.
- 2 Wipe the entire surface of the handle with a damp cloth.

Note: Do not tap the handle on the sink to remove excess water.

Brush head and travel case

- 1 Rinse the brush head and bristles after each use.

- 2 Remove the brush head from the handle and rinse the brush head connection with warm water at least once a week. Rinse the travel case weekly with warm water and wipe clean with a damp cloth.

Storage

If you are not going to use the product for a long time, clean it (see 'Cleaning') and store it in a cool and dry place away from direct sunlight.

Warning: If product will not be used on a regular basis, remove the battery to prevent possible damage due to chemical leakage. If battery does leak, remove it carefully. Do not allow bare skin to touch leaking fluid.

Replacement

Brush head

Replace the brush heads every 3 months to achieve optimal results. Only Philips One brush heads are compatible with Philips One.

Locating the model number

Remove brush head from handle. Look at the top of the handle near the metal shaft for the model number (HY110X).

Warranty and support

Philips warrants this product for one year after the date of purchase. For questions, call 1-800-682-7664 or visit www.philips.com/support. Philips recommends use of genuine Sonicare brush heads for optimal performance. Use of other brush heads may damage your Sonicare toothbrush and limit your warranty protection.

Warranty restrictions

The terms of the international warranty do not cover the following:

- Brush heads.
- Replacement batteries.
- Damage caused by use of unauthorized replacement parts.
- Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.

Disposal

Your product is designed and manufactured with high quality materials and components, which can be recycled and reused. For recycling information, please contact your local waste management facilities or visit www.recycle.philips.com

- This product contains batteries:
- Dispose of batteries properly. Do not incinerate. Batteries may explode if overheated.
- Do not wrap in metal or aluminum foil. Tape the waste battery terminals before discarding.
- It is suggested that you contact your local town or city to determine proper battery redemption site(s) in your area.
- To remove the battery, follow the instructions in the battery replacement chapter.

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