
PHILIPS


ExpertClean

7300, 7500

PHILIPS
sonicare

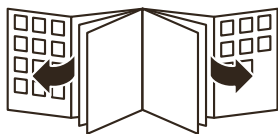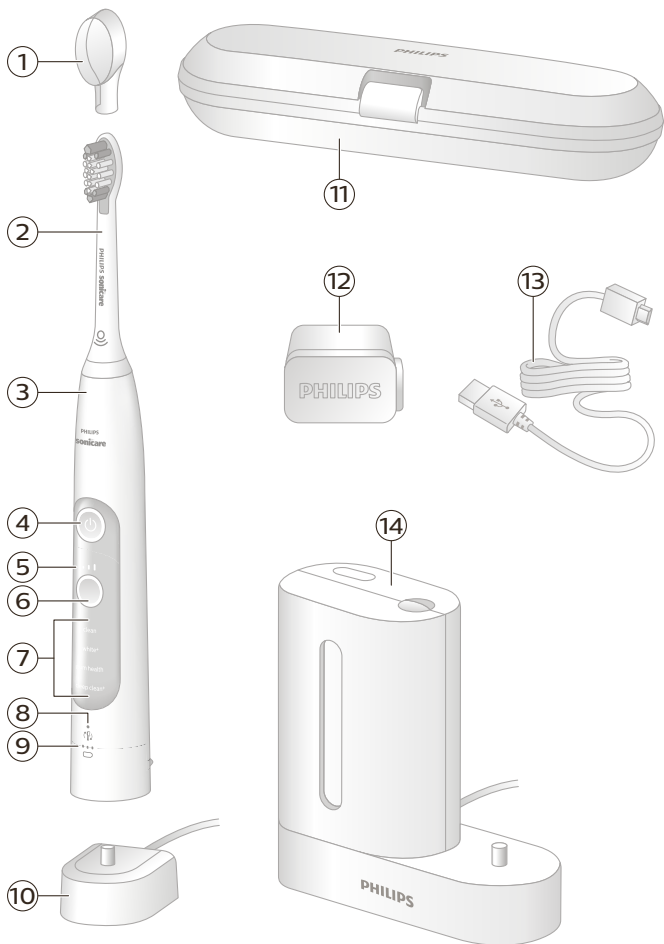English	6
Čeština	30
Eesti	56
Hrvatski	81
Latviešu	108
Lietuviškai	133
Magyar	160
Polski	187
Română	216
Slovenščina	243
Slovensky	269
Srpski	296
Български	318
Українська	348

Introduction

Congratulations on your new Philips Sonicare power toothbrush! Superior plaque removal, whiter teeth and healthier gums are at your fingertips. Using Sonicare's combination of gentle sonic technology and clinically developed and proven features, you can be confident that you're getting the very best clean, every time. Further support and product registration are available to you at: <http://www.philips.com/welcome>. Welcome to the Philips Sonicare family!

Intended Use

ExpertClean series power toothbrushes are intended to remove adherent plaque and food debris from the teeth to reduce tooth decay and improve and maintain oral health. ExpertClean series power toothbrushes are intended for consumer home use. Use by children should be with adult supervision.

Important safety information

Read this important information carefully before you use the appliance and save it for future reference.

Note: The accessories and features of different models vary. For example, some models come with a sanitizer or are equipped with Bluetooth.

Danger

- Keep the charger and/or sanitizer away from water. Do not place or store the appliance where it can fall or be pulled into a tub or sink.
- Do not immerse the charger and/or sanitizer in water or any other liquid.
- After cleaning, make sure the charger and/or sanitizer is completely dry before you connect it to the wall socket.

Warning

- This appliance can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.
- Children shall not play with the appliance.
- Keep the sanitizer out of the reach of children at all times.
- Stop using the sanitizer if the UV light bulb remains on when the door is open. UV light can be harmful to the human eye and skin.
- If the appliance is damaged in any way (brush head, toothbrush handle, charger and/or sanitizer), stop using it.
- The mains cord cannot be replaced. If the mains cord is damaged, discard the charger and/or sanitizer.

- Always replace the charger and/or sanitizer with one of the original type in order to avoid a hazard.
- This appliance contains no user-serviceable parts. If the appliance is damaged, contact the Consumer Care Center in your country (see 'Warranty and support').
- Keep the cord away from heated surfaces.
- Do not use the charger and/or sanitizer outdoors or near heated surfaces.

Caution

- Do not clean the brush head, handle, charger, charger cover or sanitizer in the dishwasher.
- If you have had oral or gum surgery in the previous 2 months, consult your dentist before you use this appliance.
- Consult your dentist if excessive bleeding occurs after using this appliance or if bleeding continues to occur after 1 week of use. Also consult your dentist if you experience discomfort or pain when you use this appliance.

- This Philips appliance complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your doctor or the appliance manufacturer of the implanted appliance prior to use.
- If you have medical concerns, consult your doctor before you use this appliance.
- This appliance has only been designed for cleaning teeth, gums and tongue. Do not use it for any other purpose. Stop using the appliance and contact your doctor if you experience any discomfort or pain.
- This appliance is a personal care device and is not intended for use on multiple patients in a dental practice or institution.
- Do not use other brush heads than the ones recommended by Philips Sonicare.
- Stop using a brush head with crushed or bent bristles. Replace the brush head every 3 months or sooner if signs of wear appear.

- Avoid direct contact with products that contain essential oils or coconut oil. Contact may result in bristles dislodging.
- The UV light bulb of the sanitizer is hot during and immediately after the sanitizing cycle. Do not touch the UV light bulb when it is hot.
- Do not operate the sanitizer without the protective screen in place to avoid contact with a hot UV light bulb.
- If you have used the sanitizer for 3 UV clean cycles in a row, switch it off for at least 30 minutes before you start another UV clean cycle.

Electromagnetic fields (EMF)

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

Radio Equipment Directive

Hereby, Philips declares that electric toothbrushes with radio frequency interface (Bluetooth) are in compliance with Directive 2014/53/EU.

The full text of the EU declaration of conformity is available at the following internet address: **www.philips.com/support**

The radio frequency interface in this product operates at 13.56 MHz.

The maximum RF power transmitted by the appliance is 30.16 dBm.

Your Philips Sonicare (Fig. 1)

- 1 Hygienic travel cap
- 2 Smart brush head(s)


- 3 Handle
- 4 Power on/off button
- 5 Intensity light
- 6 Mode/intensity button
- 7 Modes
- 8 Brush head replacement reminder light
- 9 Battery indication
- 10 Charger
- 11 Travel case *
- 12 USB wall adapter*
- 13 Micro USB cable *
- 14 UV Sanitizer*

* Note: The content of the box may vary based on the model purchased.

Using your Philips Sonicare

Brushing instructions

- 1 Push the brush head firmly onto the handle. Wet bristles and apply toothpaste.


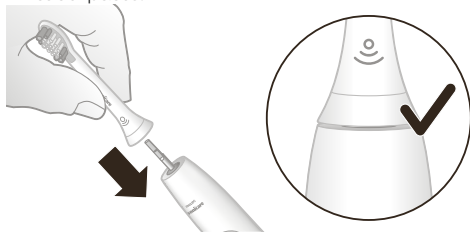
Note: It is normal to see a slight gap between the brush head and the handle. This allows the brush head to vibrate properly.


- 2 Place the toothbrush bristles against the teeth at a slight angle (45 degrees), pressing firmly to make the bristles reach the gumline or slightly beneath the gumline.

Note: Keep the center of the brush in contact with the teeth at all times.

- 3 Press the power on/off button to turn on the Philips Sonicare.

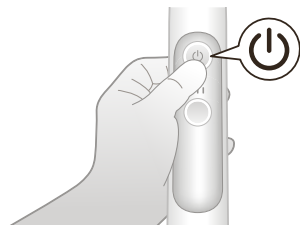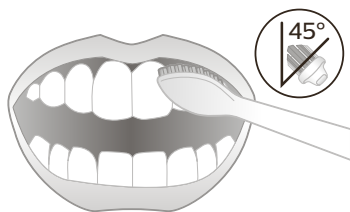
- 4 Apply light pressure to maximize Philips Sonicare's effectiveness and let the Philips Sonicare toothbrush do the brushing for you. Gently move the brush head slowly across the teeth in a small back and forth motion so the longer bristles reach between your teeth. Continue this motion throughout the brushing cycle.

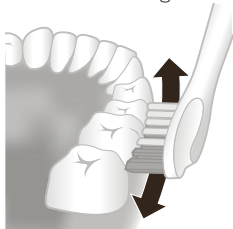Note: To learn how to brush with specific modes, see 'Modes'.

- 5 To clean the inside surfaces of the front teeth, tilt the brush handle semi-upright and make several vertical overlapping brushing strokes on each tooth.

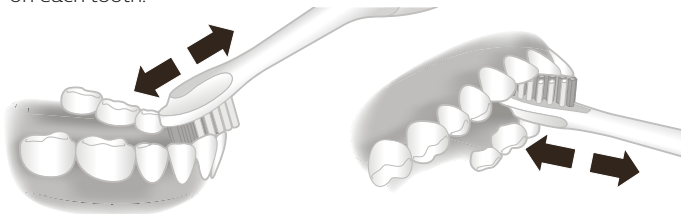Note:

- A change in vibration from the handle alerts you when you apply too much pressure (see Pressure sensor)
- Your Philips Sonicare toothbrush is safe to use on braces ,dental

moderations and dental restorations (fillings, crowns, veneers)
 - Do not use your power toothbrush in the shower.

Brush heads

To see what type of brush head you have, check the label at the bottom of the brush head. The symbol at the bottom of the brush head indicates that it is a smart brush head, equipped with the BrushSync technology.

Brush head name	Premium Plaque Control	Premium White	Premium Gum Care
Brush head label	C3	W3	G3
Benefit	Plaque removal	Plaque and stain removal	Healthier Gums

To explore our full range of brush heads, go to www.philips.com/toothbrush-heads for more information.

*Note: The type of brush heads vary based on the model purchased.

BrushSync Technology

BrushSync technology enables your brush heads to communicate with your handle using a microchip. The symbol at the bottom of the brush head indicates that it is equipped with the BrushSync technology. This enables the following features:


- Brush head mode pairing
- Brush head replacement reminder

Brushing modes

Your power toothbrush is equipped with up to 4 different modes*: Clean, Deep Clean+, White+* and Gum health*.

Based on the smart brush head you use, the mode will be selected automatically.

To manually change the mode, press the mode/intensity button while the toothbrush is off.

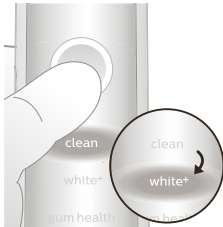* Note:







- The type and number of modes vary based on the model purchased.
- All brush heads work with each mode.
- When the Philips Sonicare toothbrush is used in clinical studies, it must be used at high intensity, with the EasyStart feature turned off.

Brush head mode pairing

Brush heads with BrushSync technology are equipped with a microchip which communicates to the handle and automatically pairs it with the recommended mode and intensity.

Note: If you change modes/intensity, the brush head will remember the last used setting. If you wish to brush with the recommended setting, you can find the recommended mode/intensity in 'Brushing Modes'.

	Clean and Deep Clean+	White+	Gum Health
Benefit	Plaque removal	Plaque and surface stains removal	Plaque removal and gum massage
Recommended brush head	Premium Plaque Control* (C3)	Premium White* (W3)	Premium Gum Care* (G3)
Recommended intensity level	3	3	3
Total brushing time	Clean: 2 minutes Deep Clean: 3 minutes*	2 minutes and 40 seconds	3 minutes and 20 seconds

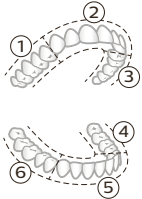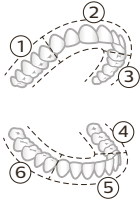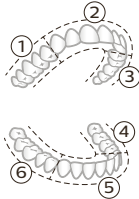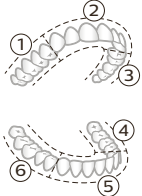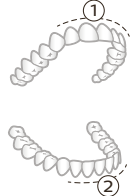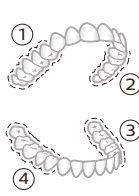	Clean and Deep Clean+	White+	Gum Health
How to brush	Clean: Brush each segment for 20 seconds	Brush each segment for 20 seconds,	Brush each segment for 20 seconds,
			
	Deep Clean+: Brush each segment for 30 seconds	then brush segments 1 and 2 for 20 seconds each.	then brush segments 1, 2, 3 and 4 for 20 seconds each.
			

***Note:**

- The modes and type brush head vary based on the model purchased.
- When brushing with the Sonicare App in 'Deep Clean' mode, the brushing time will be 2 minutes.

Intensity settings

Your power toothbrush comes with 3 different intensity settings: low, medium and high. When attached to the handle, the smart brush head will automatically select the recommended intensity. To manually change the intensity, turn on the handle and press the mode/intensity button.

While toothbrush is ON, press this button to change intensity setting

Features

- EasyStart
- Sonicare App - Connected brushing experience
- Brush head replacement reminder
- Pressure sensor
- BrushPacer

EasyStart

The EasyStart feature gently increases the power of the handle over the first 14 brushings to help you get used to the brushing power of the Philips Sonicare.

The EasyStart feature works with all modes.

Note: The EasyStart feature comes deactivated with your product. To activate EasyStart, see 'Activating or deactivating features'.

Sonicare App – Getting Started

The Sonicare App pairs with your toothbrush to provide you with interactive features. By connecting your toothbrush to your Sonicare App account, you will be able to:

- Easily review your weekly Progress Report for ways to improve your brushing routine
- Receive personalized tips and recommendations for managing your oral health

To start using the Sonicare App:

- Download the Sonicare App to your phone
- Open the App and follow the guided steps
- Pair your toothbrush with the App
- Create your account
- Personalize your App experience by answering questions about your brushing habits and interests

- Brush regularly
- Receive weekly updates to improve your oral care

Note: Make sure your phone's Bluetooth is turned on. Your toothbrush uses the Bluetooth connection to transfer your brushing data to the App. If you have questions about why your data is collected, be sure to review the Sonicare Privacy Statement, available throughout the App setup process.

Sonicare App – Your oral healthcare dashboard

Your Sonicare power toothbrush is equipped with smart sensors to provide you with feedback about your brushing habits, including:

- Daily brushing sessions
- Time spent during each brushing session
- Applied brushing pressure
- Brush head replacement reminder (based on its actual usage and effectiveness)

The Sonicare App collects the data provided by the smart sensors and tracks your brushing habits over time.

Each of the interactive features available in the App are designed to help you understand your brushing habits and provide you with simple, dentist-approved recommendations for improving and maintaining your oral health. You will receive personalized brushing information in the Sonicare App.

Brush Head Replacement Reminder

Over time, the handle will track the wear of your smart brush head by measuring:

- The overall pressure you apply while brushing
- The total time you have brushed with your brush head

Once your brush head is no longer effective, the Brush Head Replacement Reminder will blink amber and audibly beep.

Note: The Brush Head Replacement Reminder comes activated with your product. To deactivate this feature, see 'Activating or deactivating features'.

Pressure sensor

Your Philips Sonicare measures the pressure you apply while brushing to protect your gums and teeth from damage.

If you apply too much pressure, the handle will change its vibration until you reduce the pressure.

Note: The Pressure Sensor comes activated with your product. To deactivate this feature, see 'Activating or deactivating features'.

BrushPacer

The BrushPacer features indicates when you are supposed to move onto the next section of your teeth via a brief change in vibration.

The toothbrush will automatically turn off at the end of the brushing session.

Charging and battery status

Charging

- 1 Plug the charger into an electrical outlet.
- 2 Place the handle on the charger.
- 3 The charging indicator will flash green until your toothbrush is fully charged.

Note: Your handle comes pre-charged for first use. After first use, charge for at least 24 hours.

Charging on charger or sanitizer

- 1 Plug the charger or sanitizer into an electrical outlet.
- 2 Place the handle on the charger or sanitizer.
- 3 The charger will charge your toothbrush until it is fully charged.

Note: Your handle comes pre-charged for first use. After first use, charge for at least 24 hours.

Charging with the travel case*

* Note: The charging travel case may not be included depending on the model purchased.

- 1 Connect the charging travel case to 5V DC input power either from:
 - The USB adapter provided

- a UL-listed or certified to ANSI/UL 60950-1 personal computer's standard USB connector or,
- a UL-listed or certified to ANSI/UL 60950-1 LPS power supply with USB output port or,
- a UL-listed or certified 1310 Class 2 power supply.

Important: The wall adapter must be a safe extra low voltage source, compatible with the electrical ratings marked on the bottom of the device.

2 The travel case will charge your toothbrush until your toothbrush is fully charged.

Note: Charging the handle with the charging travel case may take longer than on the charger.

Battery Status (handle on plugged in charger)

When the handle is placed on the charger, sanitizer or travel case, the battery light will flash green based on the current battery level.

Number of Flashing LEDs	Battery charging status
1	Low
2	Medium
3	High

Note: When the handle is fully charged, the three battery indicator lights light up green. After 30 seconds the battery indicator lights switch off.

Battery Status (when handle is not placed on Charger)

When removing the handle from the charger, sanitizer or travel case, the battery light at the bottom of the handle will indicate the status of the battery.


LED Color	Number of LEDs	Battery status
Solid Green	3	High
Solid Green	2	Medium
Solid Green	1	Low
Flashing Amber	1 (with beeps)	Empty or near empty

UV Sanitizer

- With the UV sanitizer you can clean your brush head after every use.
- Unplug the sanitizer, discontinue use and call Consumer Care if:
 - The UV light bulb remains on when the door is open.
 - The window is broken or missing from the sanitizer.
 - The sanitizer gives off smoke or a burning smell while it is operating.

Note: UV light can be harmful to the human eye and skin. This device should be kept out of reach of children at all times.

- 1 After brushing, rinse the brush head and shake off excess water.

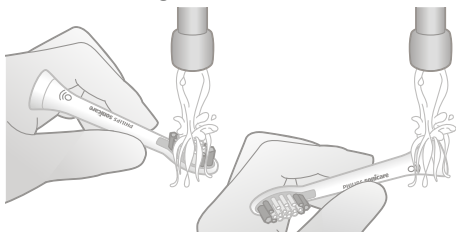Warning: Do not place the travel cap on the brush head during sanitization.

- 2 Place your finger in the door recess on top of the sanitizer and pull open the door.


- 3 Place the brush head on one of the 2 pegs in the sanitizer with the bristles facing the light bulb.

Note: Only clean Philips Sonicare snap-on brush heads in the sanitizer.

Note: Do not clean Philips Sonicare for Kids brush heads in the sanitizer.

- 4 Make sure the sanitizer is plugged into a live outlet of appropriate voltage.

- 5 Close the door and press the green power on/off button once to select the UV clean cycle.

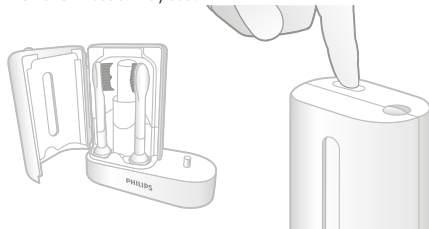Note: You can only turn on the sanitizer if the door is properly closed.

Note: The sanitizer stops if you open the door during the sanitizing cycle.

Note: The sanitizer cycle runs for 10 minutes and then automatically shuts off.

The sanitizer is in operation when the light glows through the window.

Activating or deactivating features

You can activate or deactivate the following features of your toothbrush:

- EasyStart
- Brush Head Replacement Reminder
- Pressure Sensor

Step 1: Place the handle on the charger

Step 2: Press and hold power button for:


	EasyStart	Brush Head Replacement Reminder	Pressure Sensor
	3 seconds	Up to 5 seconds	Up to 7 seconds
	↓	↓	↓
Step 3: Let go power button when you hear			
	1 beep	1 beep and then 2 beeps	1 beep, 2 beeps and then 3 beeps
	↓	↓	↓
 	If you see battery indicator lights flashing green 2 times and hear Low to High triple tones, then the feature has been activated.		
OR			
 	If you see battery indicator lights flashing amber once and hear High to Low triple tones, then the feature has been deactivated.		

Cleaning


- All parts can be cleaned by damp cloth or rinsed with warm water.


- Remove the brush head from the handle and rinse it thoroughly.

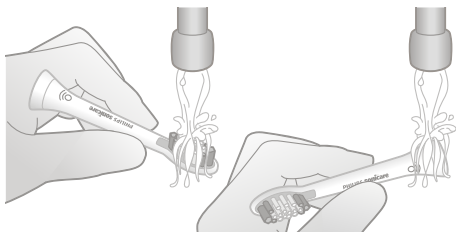- Rinse the entire handle, especially the brush head connection, at least once a week.

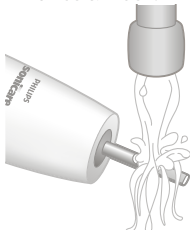Caution: Do not push on the rubber seal at the top of the handle. Gently clean around the rubber seal.

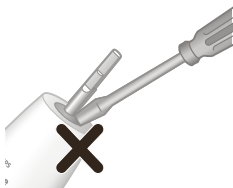Cautions:

- Do not clean the brush head, travel cap, handle, charger, UV sanitizer or travel case in the dishwasher.
- Do not use isopropyl alcohol, vinegar, bleach, or any other house hold cleaning products, to clean the handle as this may cause discoloration.
- Unplug the charger before you clean it.
- Do not use essential oils to clean the brush heads as this can cause damage.


Cleaning the UV Sanitizer

Warning: Do not immerse the sanitizer in water or rinse it under the faucet.

Warning: Do not clean the sanitizer when the UV light bulb is hot.

For optimal effectiveness, it is recommended to clean your sanitizer weekly.

- 1 Unplug the sanitizer.
- 2 Pull the drip tray straight out. Clean the interior of the sanitizer with a damp cloth.
 - Rinse the drip tray and wipe it clean with a damp cloth.

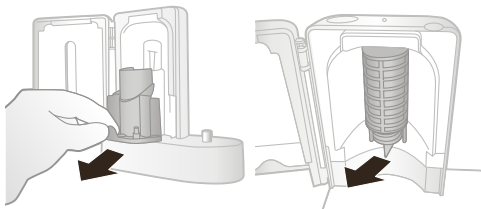- 3 Remove the protective screen in front of the UV light bulb.
To remove the screen, grasp the edges next to the snaps, gently squeeze and pull out the protective screen.
- 4 Remove the UV light bulb.
To remove the light bulb, grasp it and pull it out of the metal clasp.
- 5 Clean the protective screen and the UV light bulb with a damp cloth.
- 6 Reinsert the UV light bulb.
To reinsert the light bulb, align the bottom of the light bulb with the metal clasp and push the bulb into the clasp.
- 7 Reinsert the protective screen.
To reinsert the screen, align the pegs on the screen with the slots on the reflective surface near the UV light bulb. Push the screen straight into the slots on the sanitizer.

Storage

If you are not going to use the product for an extended period of time, unplug it from the electrical outlet, clean it and store it in a cool and dry place away from direct sunlight.

Frequently asked questions

Question	Answer
Can I leave my Sonicare on the charger between brushings?	Yes, you can always keep your Philips Sonicare toothbrush on a plugged-in charger. This will not affect the battery lifetime.
Will other brush heads fit my Philips Sonicare toothbrush?	All Philips Sonicare electric toothbrush brush heads that snap on will fit your toothbrush. However, not all can be used for BrushSync technology. You can identify the brush heads with BrushSync technology by the icon at the base of the brush head.
How often do I need to replace my brush head?	We recommend you to replace your Philips Sonicare brush head after 3 months of use (based on brushing twice a day for 2 minutes each session). The Brush Head Replacement Reminder will also notify you when it is time to change your brush head.
Is my Philips Sonicare toothbrush waterproof?	Yes, your Philips Sonicare toothbrush is waterproof. However, never put your handle in a dishwasher to clean it. And do not use it under the shower.
What phones are compatible with the Philips Sonicare App?	The Philips Sonicare app support Android OS, and Apple iOS. For the most up to date version of the Philips Sonicare App check out the Android Play Store or Apple App Store. Mobile devices must have Bluetooth 4.0 and newer.

Question	Answer
What other Philips Sonicare toothbrushes can I use with the Philips Sonicare App?	You can only use Philips Sonicare Flexcare Platinum Connected, DiamondClean Smart, expert results and ExpertClean toothbrushes with the Sonicare app. Look for the Bluetooth symbol on the front panel of the package.

Recycling

- This symbol means that this product shall not be disposed of with normal household waste (2012/19/EU).
- This symbol means that this product contains a built-in rechargeable battery which shall not be disposed of with normal household waste (2006/66/EC). Please take your product to an official collection point or a Philips service center to have a professional remove the rechargeable battery.
- The UV light bulb of the sanitizer contains mercury. Do not throw the UV light bulb away with the normal household waste at the end of its life. Hand it in at an official collection point for recycling.
- Follow your country's rules for the separate collection of electrical and electronic products and rechargeable batteries. Correct disposal helps prevent negative consequences for the environment and human health.


Removing the rechargeable battery

Warning: Only remove the rechargeable battery when you discard the appliance. Make sure the battery is completely empty when you remove it. To remove the rechargeable battery, you need a towel or cloth, a hammer and a flat-head (standard) screwdriver. Observe basic safety precautions when you follow the procedure outlined below. Be sure to protect your eyes, hands, fingers, and the surface on which you work.

- 1 To deplete the rechargeable battery of any charge, remove the handle from the charger, turn on the Philips Sonicare and let it run until it stops. Repeat this step until you can no longer turn on the Philips Sonicare.


- 2** Remove and discard the brush head. Cover the entire handle with a towel or cloth.

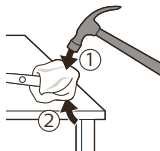
- 3** Hold the top of the handle with one hand and strike the handle housing 0.5 inch above the bottom end. Strike firmly with a hammer on all 4 sides to eject the end cap.


Note: You may have to hit on the end several times to break the internal snap connections.

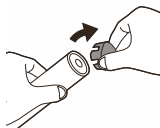
- 4** Remove the end cap from the toothbrush handle. If the end cap does not release easily from the housing, repeat step 3 until the end cap is released.

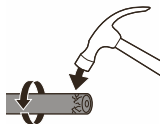
- 5** Holding the handle upside down, press the shaft down on a hard surface. If the internal components do not easily release from the housing, repeat step 3 until the internal components are released.

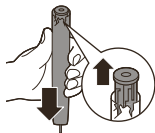- 6** Wedge the screwdriver between the battery and the black frame at the bottom of the internal components. Then pry the screwdriver away from the battery to break the bottom of the black frame.

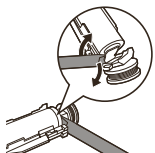- 7** Insert the screwdriver between the bottom of the battery and the black frame to break the metal tab connecting the battery to the green printed circuit board. This will release the bottom end of the battery from the frame.

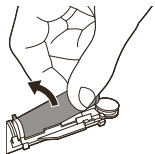- 8** Grab the battery and pull it away from the internal components to break the second metal battery tab.

Caution: Be aware of the sharp edges of the battery tabs so as to avoid injury to your fingers.

- 9** Cover the battery contacts with tape to prevent any electrical short from residual battery charge. The rechargeable battery can now be recycled and the rest of the product discarded appropriately.

Note: Do not throw away the appliance with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment.

Warranty and support

If you need information or support, please visit www.philips.com/support or read the international warranty leaflet.

Warranty restrictions

The terms of the international warranty do not cover the following:

- Brush heads.
- Damage caused by use of unauthorized replacement parts.
- Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.
- UV light bulb.

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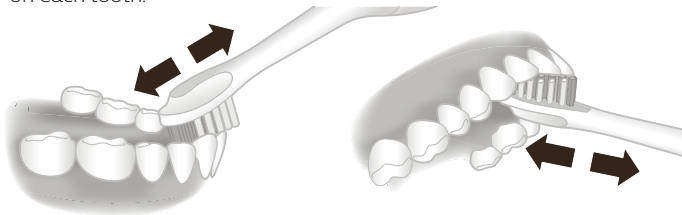
Note:

- A change in vibration from the handle alerts you when you apply too much pressure (see Pressure sensor)
- Your Philips Sonicare toothbrush is safe to use on braces, dental moderations and dental restorations (fillings, crowns, veneers)
- Do not use your power toothbrush in the shower.

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To explore our full range of brush heads, go to **www.philips.com/toothbrush-heads** for more information.

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BrushPacer

The BrushPacer features indicates when you are supposed to move onto the next section of your teeth via a brief change in vibration. The toothbrush will automatically turn off at the end of the brushing session.

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- 1 Connect the charging travel case to 5V DC input power either from:

- The USB adapter provided
- a UL-listed or certified to ANSI/UL 60950-1 personal computer's standard USB connector or,
- a UL-listed or certified to ANSI/UL 60950-1 LPS power supply with USB output port or,
- a UL-listed or certified 1310 Class 2 power supply.

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Up to 5 seconds

Up to 7 seconds

Step 3: Let go power button when you hear

1 beep

1 beep and then
2 beeps

1 beep, 2 beeps and
then 3 beeps

If you see battery indicator lights flashing green 2 times and hear Low to High triple tones, then the feature has been activated.

OR


If you see battery indicator lights flashing amber once and hear High to Low triple tones, then the feature has been deactivated.

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
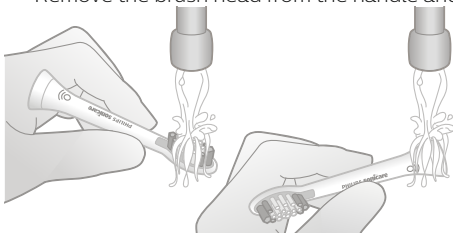
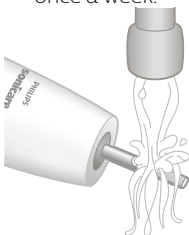
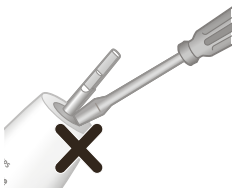
- 1 Isključite uređaj za sanitizaciju.
- 2 Pull the drip tray straight out. Clean the interior of the sanitizer with a damp cloth.
 - Rinse the drip tray and wipe it clean with a damp cloth.

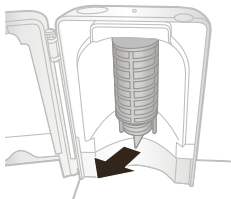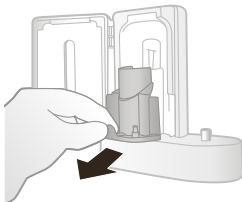- 3 Remove the protective screen in front of the UV light bulb.
Da biste uklonili foliju, uhvatite iverice pored jezičaka, nežno stisnite i izvucite zaštitnu foliju.
- 4 Uklonite UV lampu.
Da biste uklonili lampu, uhvatite je i izvucite iz metalne kopče.
- 5 Očistite zaštitnu foliju i UV lampu vlažnom krpom.
- 6 Vratite UV lampu na mesto.
To reinsert the light bulb, align the bottom of the light bulb with the metal clasp and push the bulb into the clasp.
- 7 Vratite zaštitnu foliju.
To reinsert the screen, align the pegs on the screen with the slots on the reflective surface near the UV light bulb. Push the screen straight into the slots on the sanitizer.

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Ako proizvod nećete koristiti duže vreme, isključite ga iz zidne utičnice, očistite ga i odložite na hladno i suvo mesto, dalje od direktne sunčeve svetlosti.

Najčešća pitanja

Question	Answer
Can I leave my Sonicare on the charger between brushings?	Yes, you can always keep your Philips Sonicare toothbrush on a plugged-in charger. This will not affect the battery lifetime.
Will other brush heads fit my Philips Sonicare toothbrush?	All Philips Sonicare electric toothbrush brush heads that snap on will fit your toothbrush. However, not all can be used for BrushSync technology. You can identify the brush heads with BrushSync technology by the icon at the base of the brush head.

Question	Answer
How often do I need to replace my brush head?	We recommend you to replace your Philips Sonicare brush head after 3 months of use (based on brushing twice a day for 2 minutes each session). The Brush Head Replacement Reminder will also notify you when it is time to change your brush head.
Is my Philips Sonicare toothbrush waterproof?	Yes, your Philips Sonicare toothbrush is waterproof. However, never put your handle in a dishwasher to clean it. And do not use it under the shower.
What phones are compatible with the Philips Sonicare App?	The Philips Sonicare app support Android OS, and Apple iOS. For the most up to date version of the Philips Sonicare App check out the Android Play Store or Apple App Store. Mobile devices must have Bluetooth 4.0 and newer.
What other Philips Sonicare toothbrushes can I use with the Philips Sonicare App?	You can only use Philips Sonicare Flexcare Platinum Connected, DiamondClean Smart, expert results and ExpertClean toothbrushes with the Sonicare app. Look for the Bluetooth symbol on the front panel of the package.

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

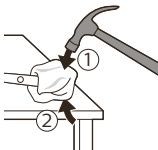
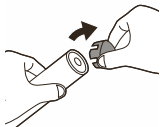
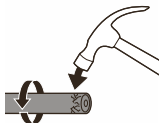
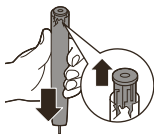
- 6 Wedge the screwdriver between the battery and the black frame at the bottom of the internal components. Then pry the screwdriver away from the battery to break the bottom of the black frame.

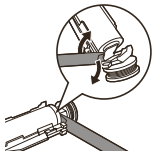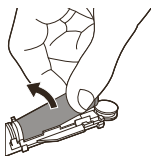- 7 Insert the screwdriver between the bottom of the battery and the black frame to break the metal tab connecting the battery to the green printed circuit board. This will release the bottom end of the battery from the frame.

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©2021 Koninklijke Philips N.V. (KPNV). All rights reserved. Philips and the Philips shield are trademarks of KPNV. Sonicare and the Sonicare logo are trademarks of Philips Oral Healthcare, LLC and/or KPNV.