www.philips.com/one

Philips Oral Healthcare LLC.  
22100 Bothell Everett Highway Bothell,  
WA 98021





©2020 Koninklijke Philips N.V. (KPNV). All rights reserved.

Philips and the Philips shield are trademarks of KPNV.

Sonicare and the Sonicare logo are trademarks of Philips Oral Healthcare, LLC and/or KPNV.

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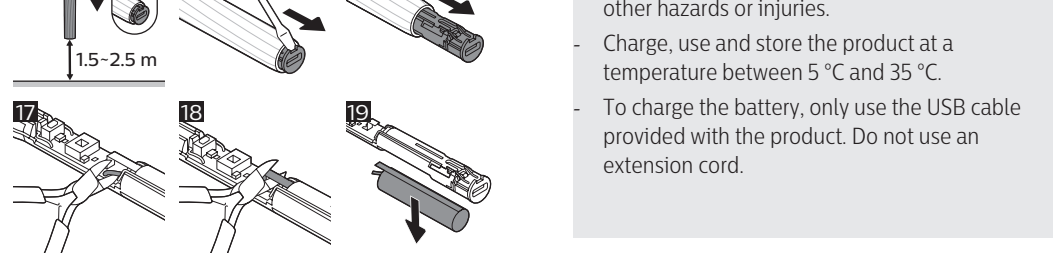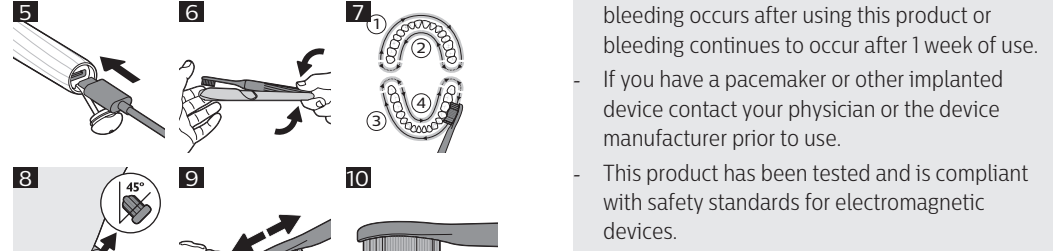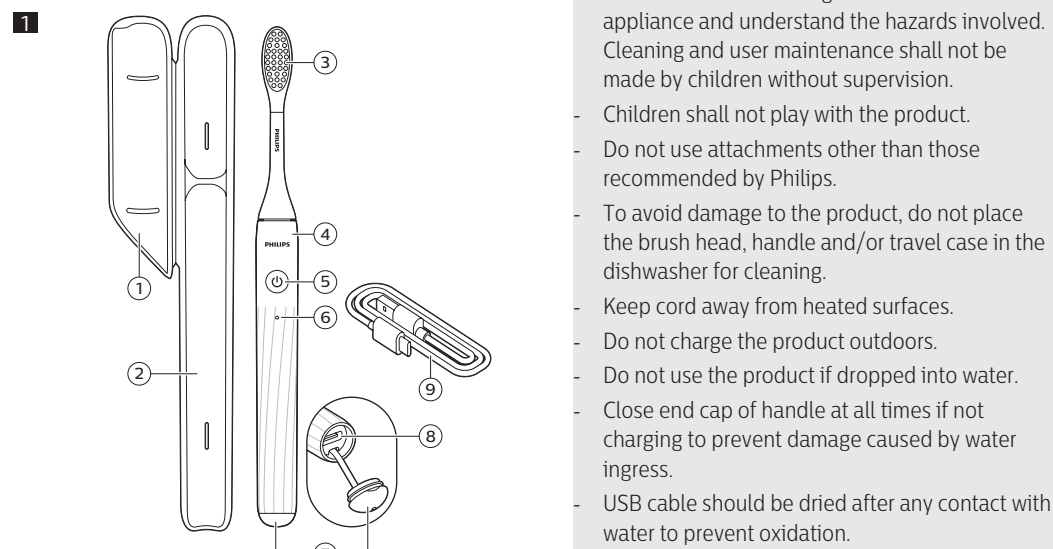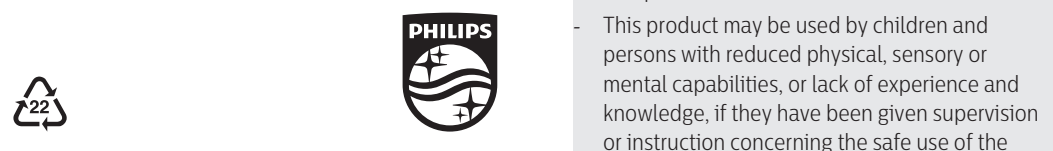## IMPORTANT SAFEGUARDS

### READ ALL INSTRUCTIONS BEFORE USE

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

## DANGERS

To reduce the risk of electrocution:

- Do not place or store product where it can fall or be pulled into a tub or sink.
- Do not place in or drop into water or other liquid. Do not use while bathing.
- Do not reach for a product that has fallen into water. Unplug immediately.
- After cleaning, make sure the toothbrush handle is completely dry before you connect it to the mains.

## WARNINGS

To reduce the risk of burns, electrocution, fire or injury to persons:

- Charge only with Philips HX9200 or any UL approved adapters with input ratings (100-240 V -, 50/60 Hz; 3.5W) and output ratings (5Vdc; 2.5W). For the list of adapters, please see [Philips.com/support](http://Philips.com/support).
- Discontinue use if this product appears damaged in any way (brush head, handle, travel case, and/or USB cable). This product contains no serviceable parts. Refer to chapter 'Warranty and support' if your Philips One no longer works properly or needs repair.
- This product is designed to clean your teeth, gums and tongue only. Use this product only for its intended use as described in this user manual. Discontinue use of this product and contact a physician/dentist if discomfort or pain is experienced.
- This product may be used by children and persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, if they have been given supervision or instruction concerning the safe use of the appliance and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision. Children shall not play with the product. Do not use attachments other than those recommended by Philips.
- To avoid damage to the product, do not place the brush head, handle and/or travel case in the dishwasher for cleaning.
- Keep cord away from heated surfaces.
- Do not charge the product outdoors.
- Do not use the product if dropped into water.
- Close end cap of handle at all times if not charging to prevent damage caused by water ingress.
- USB cable should be dried after any contact with water to prevent oxidation.

## MEDICAL WARNINGS

- Consult your dental professional before you use this product if you have had oral or gum surgery in the previous 2 months.
- Contact your dental professional if excessive bleeding occurs after using this product or bleeding continues to occur after 1 week of use.
- If you have a pacemaker or other implanted device contact your physician or the device manufacturer prior to use.
- This product has been tested and is compliant with safety standards for electromagnetic devices.
- Consult your physician prior to using the product if you have medical concerns.

## BATTERY SAFETY INSTRUCTIONS

- Only use this toothbrush for its intended purpose and follow the general and battery safety instructions as described in this user manual. Any misuse can cause electric shock, burns, fire and other hazards or injuries.
- Charge, use and store the product at a temperature between 5 °C and 35 °C.
- To charge the battery, only use the USB cable provided with the product. Do not use an extension cord.

- Do not burn the toothbrush handle and its battery and do not expose them to direct sunlight or to high temperatures (e.g. in hot cars or near hot stoves). Batteries may explode if overheated.
- If the product becomes abnormally hot, gives off an abnormal smell, changes color or if charging takes much longer than usual, stop using and charging the product and contact your local Philips dealer.

- Do not place the toothbrush handle and its battery inside microwave ovens or on induction cookers.
- This product contains a rechargeable battery that is non-replaceable. Do not open the toothbrush handle to replace the rechargeable battery.
- When you handle the battery to remove them during final product disposal, make sure that your hands, the product and the battery are dry.
- To prevent the battery from heating up or releasing toxic or hazardous substances, do not modify, pierce or damage the toothbrush handle and battery, and do not disassemble, short circuit, overcharge or reverse charge the battery.
- To avoid accidental short-circuiting of the battery after removal, do not let battery terminals come into contact with metal object (e.g. coins, hairpins, rings). Do not wrap the battery in aluminum foil. Tape battery terminals or put the battery in a plastic bag before you discard it.
- If the battery is damaged or leaking, avoid contact with the skin or eyes. If this occurs, immediately rinse well with water and seek medical care.

## SAVE THESE INSTRUCTIONS

**Electromagnetic fields (EMF)**  
This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

**Introduction**  
Congratulations on your purchase and welcome to Philips! To fully benefit from the support that Philips offers, register your product at [www.philips.com/welcome](http://www.philips.com/welcome).

**Intended use**  
Philips One power toothbrushes are intended to remove adherent plaque and food debris from the teeth to reduce tooth decay and improve and maintain oral health. Philips One power toothbrushes are intended for consumer home use. Use by children should be with adult supervision.

**Your Philips One (Fig. 1)**  
1 Travel case lid  
2 Travel case  
3 Brush head  
4 Handle  
5 Power Button  
6 Battery light indicator  
7 End cap  
8 USB-C charging port  
9 USB-C charging cable

**Getting started**  
**Attaching the brush head**  
1 Align the brush head so the bristles point in the same direction as the front of the handle. Press the brush head down onto the metal shaft (Fig. 2).

**Note:** It is normal to see a slight gap between the brush head and the handle (Fig. 3). This allows the brush head to vibrate properly.

**Charging your Philips One**  
1 Open the end cap at the bottom of the toothbrush handle (Fig. 4) and plug the USB-C charging cable into the USB-C charging port (Fig. 5).  
2 Connect the USB cable into a 5V USB wall adapter (Type A).

**Note:** Wall adapter not included. Only use 5V IPX-4 adapter (Type A).

3 Your Philips One beeps twice to indicate that the toothbrush is charging.  
4 Your toothbrush continues to charge until it is fully charged.  
- When the toothbrush is charging, the battery indicator fades up and down in white.  
- When fully charged, the battery indicator lights up solid white and turns off after 30 seconds of being fully charged.

**Note:** It can take up to 8 hours to fully charge the battery, but you can use the toothbrush before it is fully charged.

When not charging your Philips One, be sure the end cap is securely closed.

**Travel case**  
Your Philips One comes with a travel case to store your power toothbrush when traveling or between uses. Put your Philips One toothbrush in the travel case with the bristles facing downwards. Snap close the lid of the travel case.

To remove your Philips One toothbrush in the travel case, follow these steps:  
1 Flip open the travel case lid.  
2 While holding both ends of the travel case, pinch the base of the toothbrush (Fig. 6).  
3 The toothbrush will raise. Lift the toothbrush.

**Using your Philips One toothbrush**  
If you are using this toothbrush for the first time, it is normal to feel more vibration than when using a manual toothbrush. The Philips One uses microvibrations to help gently clean your teeth. Please follow the brushing instructions below for the best experience.

Philips One comes with a SmarTimer and QuadPacer to help you brush evenly throughout your mouth.  
**SmarTimer**  
The SmarTimer indicates that your 2-minute brushing cycle is complete by beeping 3 times and automatically turning off.

**Note:** Dental professionals recommend brushing no less than 2 minutes twice a day.

**QuadPacer**  
The QuadPacer is an interval timer that beeps every 30 seconds to remind you to brush all teeth evenly by dividing your mouth into 4 sections (Fig. 7). Begin brushing in section 1 (outside upper teeth) and brush for 30 seconds before you move to section 2 (inside upper teeth). Continue brushing in section 3 (outside lower teeth) and brush for 30 seconds before you move to section 4 (inside lower teeth).

**Note:** If you press the power button after you started the brushing cycle, the toothbrush will pause. After a pause of 30 seconds, the SmarTimer and QuadPacer reset.

**Brushing instructions**  
Please follow the brushing instructions for the best experience.  
1 Wet the bristles and apply a small amount of toothpaste.  
2 Place the toothbrush bristles against the teeth at a slight angle (45 degrees) (Fig. 8), pressing gently to make the bristles reach the gumline or slightly beneath the gumline.

**Note:** The bristles should slightly flare.

- 3 Press the power button to turn on the toothbrush.
- 4 Brush your teeth with back and forth strokes so that the bristles reach all surfaces of the teeth. Change to another section when the QuadPacer beeps. Continue this motion throughout the 2-minute brushing cycle.
- 5 To clean the inside surfaces of the front teeth, tilt the brush vertically and make several up-and-down strokes (Fig. 9).
- 6 After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs (Fig. 10). You can also brush your tongue, with the toothbrush turned on or off, as you prefer.

Your power toothbrush should be safe to use on braces (brush heads wear out sooner when used on braces) and dental restorations (fillings, crowns, veneers) if they are properly adhered and not compromised. If a problem occurs, please follow up with a dental professional.

**Battery status**  
This Philips One toothbrush is designed to provide at least 60 brushing sessions after it has been fully charged for 8 hours. This allows sufficient charge for 1 month's usage if you brush 2 times a day (each session lasting 2 minutes).

When the battery light indicator flashes amber it is time to charge your Philips One (see chapter Charging your Philips One). The table below specifies the battery status:

Number of brushings remaining	Battery indicator	Beep sound
40-60 brushings	White- solid	-
7-39 brushings	White- fades up and down	-
≤6 brushings	Amber- flashing	5 beeps after brushing
0 brushings remain	Amber- flashing	5 beeps x 2

- This product contains batteries:  
- Dispose of batteries properly. Do not incinerate. Batteries may explode if overheated.

**Note:** When the battery is completely empty, the toothbrush turns off. The toothbrush will not turn on until battery is recharged. Toothbrush does not turn on while charging.

**Note:** Battery usage time is expected to reduce over product lifetime.

**Cleaning**  
You should clean your toothbrush and accessories regularly to remove toothpaste and other residue. Failure to clean your product may result in an unhygienic product and damage to it may occur.

**Warning: Do not clean product or accessories in dishwasher, microwave, with chemicals or with boiling hot water.**

**Toothbrush handle**  
1 Remove the brush head and rinse the metal shaft area with warm water (Fig. 11). Make sure you remove any residual toothpaste.  
2 Wipe the entire surface of the handle with a damp cloth.

**Note:** Do not tap the handle on the sink to remove excess water.

**Cleaning brush head**  
1 Rinse the brush head and bristles after each use.  
2 Remove the brush head from the handle and rinse the brush head connection with warm water at least once a week.

**Cleaning travel case**  
Rinse the travel case weekly with warm water and wipe clean with a damp cloth.

**Cleaning USB charging cable and port**  
1 Use compressed air to clean the USB port on the toothbrush handle and connector on the USB charging cable.  
2 Use a slightly dampened cloth to wipe the cable

**Note:** Avoid contact of the USB port and connections with water. If contact occurs, dry immediately.

**Storage**  
If you are not going to use the product for a long time, clean it (see 'Cleaning'), charge it fully, and store it in a cool and dry place away from direct sunlight.

**Replacement**  
**Brush head**  
Replace the brush heads every 3 months to achieve optimal results. Only Philips One brush heads are compatible with Philips One.

**Locating the model number**  
Remove brush head from handle. Look at the top of the handle near the metal shaft for the model number (HY120X).

**Warranty and support**  
Philips warrants this product for one year after the date of purchase. For questions, call 1-800-682-7664 or visit [www.philips.com/support](http://www.philips.com/support). Philips recommends use of genuine Sonicare toothbrush heads for optimal performance. Use of other brush heads may damage your Sonicare toothbrush and limit your warranty protection.

**Warranty restrictions**  
The terms of the international warranty do not cover the following:  
- Brush heads.  
- Replacement batteries.  
- Damage caused by use of unauthorized replacement parts.  
- Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.  
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.

**Disposal**  
- This product contains a rechargeable nickel metal hydride battery which must be disposed of properly (Fig. 12).  
- Contact your local town or city officials for battery disposal information You can also call 1-800-822-8837 or visit [www.call2recycle.org](http://www.call2recycle.org) for battery drop off locations.  
- Your product is designed and manufactured with high quality materials and components, which can be recycled and reused. For recycling information, please contact your local waste management facilities or visit [www.recycle.philips.com](http://www.recycle.philips.com).  
- This product contains batteries:  
- Dispose of batteries properly. Do not incinerate. Batteries may explode if overheated.

- Do not wrap in metal or aluminum foil. Tape the waste battery terminals before discarding.
- It is suggested that you contact your local town or city to determine proper battery redemption site(s) in your area. To remove the battery, follow the instructions in Removing the rechargeable battery chapter.

**Removing the rechargeable battery**  
The rechargeable battery inside your Philips One cannot be replaced, but it is easy to remove for recycling. At the end of product life and prior to disposal, please remove the rechargeable battery from the handle by following the instructions below.

**Warning: Only remove the rechargeable battery when you discard the appliance. Make sure the battery is completely empty when you remove it.**

**Warning: Please note that this process is not reversible.**

To remove the rechargeable battery, you need wire cutters/scissors and a screwdriver. Observe basic safety precautions when you follow the procedure outlined below. Be sure to protect your eyes, hands, fingers, and the surface on which you work.  
1 To deplete the rechargeable battery of any charge, remove the handle from the USB cable, turn on the Philips One toothbrush and let it run until it stops. Repeat this step until you can no longer turn on the Philips One toothbrush.  
2 Remove and discard the brush head.  
3 Open the end cap and cut the tether with scissors to separate the end cap (Fig. 13). You can discard this piece.  
4 Drop the handle vertically (with the USB port facing downwards) onto a hard surface from approximately 6.5ft (2m) (Fig. 14).  
5 The internal structure will become dislodged (Fig. 14). If the internals have not been exposed, repeat step 5.  
6 Use a screwdriver to push the internal bracket out of the handle (Fig. 15, Fig. 16).  
7 Once the internal components have been removed, locate the wire and steel plate that are in contact with the battery. Cut the steel plate and battery wire from the circuit board close to the battery (Fig. 17, Fig. 18).

**Caution: Be aware of any sharp edges to avoid injury to your fingers**

8 Grab the battery and pull it out of the internal structure (Fig. 19).  
9 Cover the battery contacts with tape to prevent any electrical short from residual battery charge. The rechargeable battery can now be recycled and the rest of the product discarded appropriately.

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