**PHILIPS**








ProtectiveClean

6100

**PHILIPS**  
**sonicare**

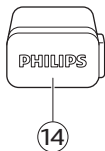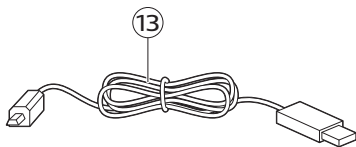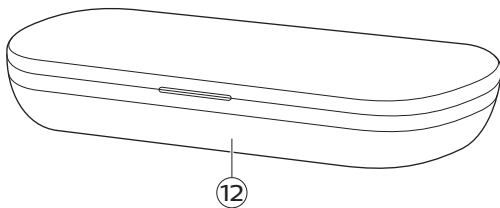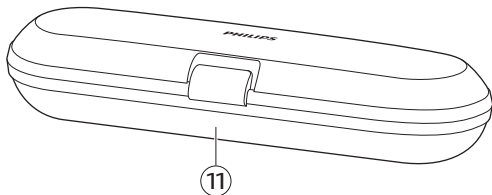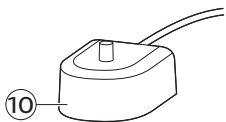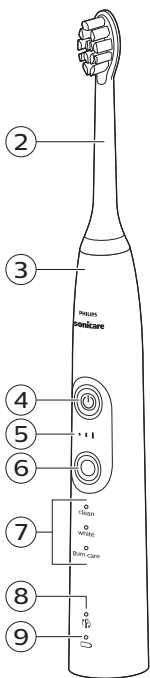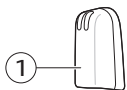**English 6**  
**Español 28**  
**Français (Canada) 52**

# IMPORTANT SAFEGUARDS

## READ ALL INSTRUCTIONS BEFORE USE

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

## DANGERS

To reduce the risk of electrocution:

- Do not place or store product where it can fall or be pulled into a tub or sink.
- Do not place in or drop into water or other liquid.
- Do not use while bathing.
- Do not immerse the charger in water or any other liquid.
- After cleaning, make sure the charger is completely dry before you connect it to a power outlet.
- Do not reach for a product that has fallen into water. Unplug immediately.

## WARNINGS

To reduce the risk of burns, electrocution, fire or physical injury:

- Do not clean the brush head, the handle or the charger in the dishwasher.
- Do not use attachments other than those recommended by Philips Sonicare.
- Never use a charger with a damaged cord or plug.
- Always replace the charger with one of the original type in order to avoid a hazard.

- Do not plug into an outlet with a voltage other than specified on the bottom of the product. This product is designed to operate within a range of 100 to 240 volts. Voltage converters do not warrant voltage compatibility.
- Never force the plug into an outlet; if the plug does not easily fit into the outlet, discontinue use.
- Discontinue use if the product appears damaged in any way (brush head, handle or charger, etc.).
- The power cord cannot be replaced. If the power cord is damaged, discard the charger.
- Keep the power cord away from heated surfaces.
- This product contains no user-serviceable parts. Refer to 'Warranty and support' if the product no longer works properly or needs repair.
- Do not use the charger outdoors.
- Do not use the charger if dropped into water
- This product is designed to clean your teeth, gums and tongue only. Use this product only for its intended use as described in this booklet. Discontinue use of this product and contact a physician/dentist if discomfort or pain is experienced.
- This product is not intended for use by persons (including children) with reduced physical, sensory, or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the product.
- Philips Sonicare is a personal care device and is not intended for use on multiple patients in a dental office or institution.

## 8 English

- This appliance contains no serviceable parts. If the appliance is damaged, contact the Consumer Care Center in your country (see chapter 'Warranty and support').
- For connection to a supply not in the U.S.A., use a charger of the proper configuration for the power outlet.
- Use a UL-listed computer, USB Hub, or USB duplex receptacle when charging the handle in the travel case.
- Avoid direct contact with products that contain essential oils or coconut oil. Contact may result in bristle tufts dislodging.
- If your toothpaste includes peroxide, baking soda or other bicarbonate (common in whitening toothpastes), thoroughly clean your brush head and handle with soap and water after each use. Otherwise plastic cracking may occur.

## MEDICAL WARNINGS

- 1 Consult your dental professional before you use this product if you have had oral or gum surgery in the previous 2 months.
- 2 Contact your dental professional if excessive bleeding occurs after using this product or bleeding continues to occur after 1 week of use.
- 3 Discontinue use of this product and contact a physician/dentist if discomfort or pain is experienced.
- 4 This product complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician or the device manufacturer prior to use.
- 5 Consult your physician if you have other medical concerns.

## SAVE THESE INSTRUCTIONS

## Electromagnetic fields (EMF)

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

## Radio Equipment Directive

- Radio Equipment in this product operates at 13.56 MHz
- Maximum RF power transmitted by the Radio Equipment is 30.16dBm

## Introduction

Congratulations on your purchase and welcome to Philips! To fully benefit from the support that Philips offers, register your product at [www.philips.com/welcome](http://www.philips.com/welcome).

## Intended use

ProtectiveClean series power toothbrushes are intended to remove adherent plaque and food debris from the teeth to reduce tooth decay and improve and maintain oral health.

## Your Philips Sonicare (Fig. 1)

- 1 Hygienic travel cap
- 2 Smart brush head(s)\*
- 3 Handle
- 4 Power on/off button
- 5 Intensity light
- 6 Mode/intensity button
- 7 Mode light
- 8 Brush head replacement reminder light
- 9 Battery indication
- 10 Charging base

- 11 Travel case (specific types only)\*
- 12 Charging travel case (specific types only)\*
- 13 USB cable (specific types only)\*
- 14 USB wall adapter (specific types only)\*

\*Note: The content of the box may vary based on the model purchased.


## Brush heads

Your Philips Sonicare comes with one or more brush heads which are designed to deliver superior results for your oral care needs.

### Philips Sonicare BrushSync Technology


Your smart brush heads use a microchip to sync with your handle, enabling brush head replacement reminder. For more information regarding the brush head replacement reminder, see 'Features'.

This symbol indicates the brush head is equipped with BrushSync technology.
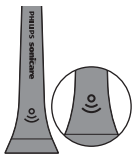
## Getting started

### Attaching the brush head

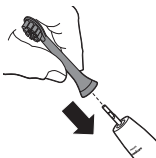- 1 Align the brush head so the bristles point in the same direction as the front of the handle.
- 2 Firmly press the brush head down onto the metal shaft until it stops.

Note: It is normal to see a slight gap between the brush head and the handle.

## Brushing modes

Your power toothbrush is equipped with 3 different modes to accommodate all your oral care needs.

To change modes, press the mode/intensity button while the toothbrush is off.

The following modes are available:

### Clean mode

Clean is a 2-minute mode, recommended to be used with intensity level 3. It gives you a thorough and complete clean and is recommended to be used with the C2 Optimal Plaque Control brush head.


### White mode

White is a 2-minute 30-seconds mode, recommended to be used with intensity level 3 and the W2 Optimal White brush head.

### Gum Care mode

Gum Care is a 3-minute mode, recommended to be used with intensity level 3 and the G2 Optimal Gum Care brush head.


## Intensity settings

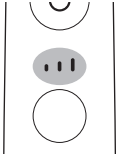Your power toothbrush gives you the option to choose between 3 different intensity levels: low, medium and high. Based on the brush head you attach, the intensity will be selected. The intensity can be changed by pressing the mode/intensity button while brushing. For optimal results it is recommended to use the automatically selected intensity. The intensity cannot be changed as long as the handle is powered off or paused.


## Using your Philips Sonicare

### Brushing instructions

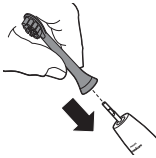- 1 Press the brush head onto the handle in such a way that the bristles will face the front of the handle.
- 2 Firmly press the brush head down onto the metal shaft until it stops.


Note: It is normal to see a slight gap between the brush head and the handle.

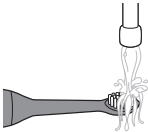- 3 Wet the bristles.


- 4** Apply a small amount of toothpaste on the bristles.

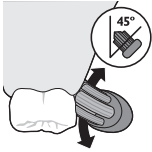- 5** Place the toothbrush bristles against the teeth at a slight angle (45 degrees), pressing firmly to make the bristles reach the gumline or slightly beneath the gumline. A change in vibration (and a slight change in sound) alerts you when you apply too much pressure while brushing.

Note: Keep the center of the brush in contact with the teeth at all times.

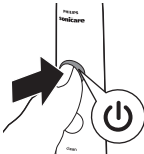- 6** Press the power on/off button to turn on the Philips Sonicare.

- 7** Gently keep the bristles placed on the teeth and in the gumline. Brush your teeth with small back and forth motion so the bristles reach between the teeth.

Note: The bristles should slightly flare. Do not scrub.


- 8** To clean the inside surfaces of the front teeth, tilt the brush handle semi-upright and make several vertical overlapping brushing strokes on each tooth.

## 14 English


- 9** To make sure you brush evenly throughout the mouth, divide the mouth into 4 sections using the Quadpacer feature. The quadpacer feature beeps after 30 seconds.

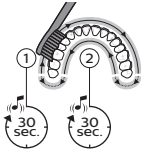- 10** Brush for the full 2 minutes following the Quadpacer. After 2 minutes, the toothbrush will automatically shut off.

- In white mode, after 2 minutes the handle will continue to brush for 30 seconds so that you can brighten and polish your front teeth.
- In gum care mode, after 2 minutes the handle will continue to brush for an additional minute so you can gently stimulate and massage your gums.
- After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs. You may also brush your tongue, with the toothbrush turned on or off, as you prefer.

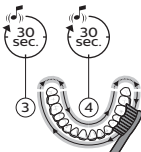Your Philips Sonicare toothbrush is safe to use on:

- Braces (brush heads wear out sooner when used on braces)
- Dental restorations (fillings, crowns, veneers)

Note: When the Philips Sonicare toothbrush is used in clinical studies, it must be used in the Clean mode (for plaque and gum health claims) or in White mode (for whitening claims) at high intensity with EasyStart turned off.

## Features

- EasyStart
- Brush head replacement reminder
- Pressure sensor
- Brush head mode pairing
- Quadpacer

## Easy-start


This Philips Sonicare model comes with the EasyStart feature deactivated. The EasyStart feature gently increases the power over the first 14 brushings to help you get used to the brushing with the Philips Sonicare. The EasyStart feature works in all modes.

To activate EasyStart see 'Activating or deactivating features'.

## Brush head replacement reminder

Your Philips Sonicare is equipped with BrushSync technology that tracks the wear of your brush head.

## 16 English

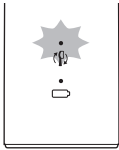- 1 When attaching a new smart brush head for the first time the brush head replacement reminder light blinks green three times. This confirms you have a Philips brush head with BrushSync technology.
- 2 Over time, based on the pressure you apply and the amount of time used, the handle will track brush head wear in order to determine the optimal time to change your brush head. This feature gives you the guarantee for the best cleaning and care of your teeth.
- 3 When the brush head replacement reminder light lights up amber, you should replace your brush head.

To deactivate the brush head replacement reminder see 'Activating or deactivating features'.

## Pressure sensor

Your Philips Sonicare is equipped with an advanced sensor that measures the pressure you apply while brushing. If you apply too much pressure, the toothbrush will provide immediate feedback to indicate that you need to reduce the pressure. This feedback is given by a change in vibration and therefore brushing sensation.

To deactivate the pressure sensor see 'Activating or deactivating features'.

## Brush head mode pairing

An RFID chip inside the brush head communicates with the toothbrush handle. When you attach a brush head, the handle automatically recognizes the brush head and selects the right mode and intensity level to optimize your brushing performance.

## Quadpacer

The Quadpacer is an interval timer that has a short beep and pause to remind you to brush the different sections of your mouth. Depending on the brushing mode you have selected (see 'Brushing modes'), the Quadpacer beeps at different intervals during the brushing cycle.

## Activating or deactivating features

You can activate or deactivate the following features of your toothbrush:

- EasyStart
- Brush head replacement reminder
- Pressure sensor

To activate or deactivate these features, follow the instructions below:

## EasyStart

- 1** Put the handle on the plugged-in charger.
- 2** Press and hold the power on/off button while the handle remains on the charger.
- 3** Keep the power on/off button pressed until you hear a single short beep (after 2 seconds).

- 4 Release the power on/off button.
  - Triple tone of low-med-high means the EasyStart feature has been activated. The brush head replacement light and battery light will also blink green 2 times in unison to confirm activation.
  - Triple tone of high-med-low means the EasyStart feature has been deactivated. The brush head replacement light and battery light will also blink amber 2 times in unison to confirm deactivation.

Note: To achieve clinical efficacy, EasyStart needs to be deactivated.

### Brush head replacement reminder

- 1 Put the handle on the plugged-in charger.
- 2 Press and hold the power on/off button while the handle remains on the charger.
- 3 Keep the power on/off button pressed until you hear a series of two short beeps (after 4-5 seconds).
- 4 Release the power on/off button.
  - Triple tone of low-med-high means the Brush head replacement reminder feature has been activated. The brush head replacement light and battery light will also blink green 2 times in unison to confirm activation.

- Triple tone of high-med-low means the Brush head replacement reminder feature has been deactivated. The brush head replacement light and battery light will also blink amber 2 times in unison to confirm deactivation.

## Pressure sensor

- 1 Put the handle on the plugged-in charger.
- 2 Press and hold the power on/off button while the handle remains on the charger.
- 3 Keep the power on/off button pressed until you hear a series of three short beeps (after 6-7 seconds).
- 4 Release the power on/off button.
  - Triple tone of low-med-high means the pressure sensor feature has been activated. The brush head replacement light and battery light will also blink green 2 times in unison to confirm activation.
  - Triple tone of high-med-low means the pressure sensor feature has been deactivated. The brush head replacement light and battery light will also blink amber 1 time in unison to confirm deactivation.

If you continue holding the power on/off button after the three short beeps, the activate/deactivate sequence repeats.

## Battery status and charging

- 1 Put the plug of the charger in an electrical outlet.
- 2 Place the handle on the charger.
  - The charger will emit 2 short beeps to confirm that the handle is placed properly.
  - The flashing light of the battery level indicator shows that the toothbrush is charging.
  - When handle is fully charged on the charger, the battery light shows solid green for 30 seconds and switches off.

Note: Your toothbrush comes pre-charged for first use. After first use, charge for at least 24 hours.

## Charging with the charging travel case \*

Note: The charging travel case may not be included depending on the model purchased.

- 1 Take the USB cable from the bottom of the charging travel case.
- 2 Connect the charging travel case to 5V DC input power either from:
  - The USB adapter provided
  - a UL-listed or certified to ANSI/UL 60950-1 personal computer's standard USB connector or,

- a UL-listed or certified to ANSI/UL 60950-1 LPS power supply with USB output port or,
- a UL-listed or certified 1310 Class 2 power supply.

**Important: The wall adapter must be a safe extra low voltage source, compatible with the electrical ratings marked on the bottom of the device.**

- 3** The travel case charges your toothbrush until your toothbrush is fully charged.
- The charger will emit 2 short beeps to confirm that the handle is placed properly.
  - The flashing light of the battery level indicator shows that the toothbrush is charging.
  - When handle is fully charged on the charger, the battery light shows solid green for 30 seconds and switches off.

### **Battery status (when handle is not on charger)**

When removing the Philips Sonicare from the charger, the battery light at the bottom of the toothbrush will indicate the status of the battery life.


- Solid green LED: full battery
- Flashing green LED: medium battery
- Flashing amber LED and three beeps: low battery

- Flashing amber LED and two sets of five beeps: No brushing sessions left (charge toothbrush)

## Cleaning

Note: Do not clean the brush head, handle or travel case in the dishwasher.

### Toothbrush handle

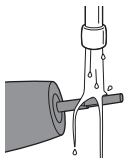
- 1 Remove the brush head and rinse the metal shaft area with warm water. Make sure you remove any residual toothpaste.

Caution: Do not push on the rubber seal on the metal shaft with sharp objects, as this may cause damage.

- 2 Wipe the entire surface of the handle with a damp cloth.

Note: Do not use isopropyl rubbing alcohol, vinegar or bleach to clean handle as this may cause discoloration.

### Brush head

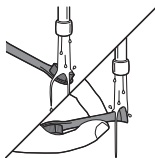- 1 Rinse the brush head and bristles after each use.
- 2 Remove the brush head from the handle and rinse the brush head connection with warm water at least once a week.

### Charger

- 1 Unplug the charger before you clean it.
- 2 Wipe the surface of the charger with a damp cloth.

## Storage

If you are not going to use your Philips Sonicare for an extended period of time, unplug the charger from the wall socket, clean it and store it in a cool and dry place away from direct sunlight.

## Locating the model number

Look on the bottom of the Philips Sonicare toothbrush handle for the model number (HX684x/HX685x).

## Disposal

This appliance contains a rechargeable battery which must be disposed of properly. Contact your local town or city officials for battery disposal information. You can also call 1-800-8-BATTERY or visit [www.rbrc.com](http://www.rbrc.com) for battery drop-off locations.


For assistance, visit our website [www.philips.com/support](http://www.philips.com/support) or call 1-800-243-3050 toll free.

## Removing the rechargeable battery

**Warning: Only remove the rechargeable battery when you discard the appliance. Make sure the battery is completely empty when you remove it. Any attempt to remove the battery will void the warranty.**

- To remove the rechargeable battery, you need a towel or cloth, a hammer and a flat-head (standard) screwdriver. Observe basic safety precautions when you follow the procedure outlined below. Be sure to protect your eyes, hands, fingers, and the surface on which you work.


- 1** To deplete the rechargeable battery of any charge, remove the handle from the charger, turn on the Philips Sonicare and let it run until it stops. Repeat this step until you can no longer turn on the Philips Sonicare.
- 2** Remove and discard the brush head. Cover the entire handle with a towel or cloth.

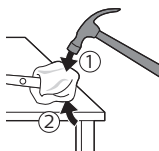
- 3** Hold the top of the handle with one hand and strike the handle housing 0.5 inch above the bottom end. Strike firmly with a hammer on all 4 sides to eject the end cap.

Note: You may have to hit on the end several times to break the internal snap connections.


- 4** Remove the end cap from the toothbrush handle. If the end cap does not release easily from the housing, repeat step 3 until the end cap is released.

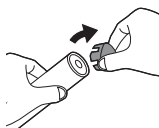
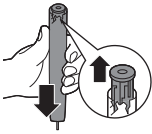
- 5 Holding the handle upside down, press the shaft down on a hard surface. If the internal components do not easily release from the housing, repeat step 3 until the internal components are released.

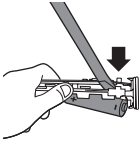- 6 Wedge the screwdriver between the battery and the black frame at the bottom of the internal components. Then pry the screwdriver away from the battery to break the bottom of the black frame.

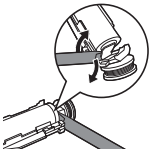- 7 Insert the screwdriver between the bottom of the battery and the black frame to break the metal tab connecting the battery to the green printed circuit board. This will release the bottom end of the battery from the frame.

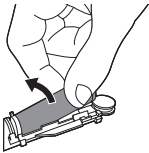- 8 Grab the battery and pull it away from the internal components to break the second metal battery tab.

**Caution: Be aware of the sharp edges of the battery tabs so as to avoid injury to your fingers.**

- 9 Cover the battery contacts with tape to prevent any electrical short from residual battery charge. The rechargeable battery can now be recycled and the rest of the product discarded appropriately.

## Warranty and support

Philips warrants its products for two years after the date of purchase. Defects due to faulty materials and workmanship will be repaired or replaced at Philips' expense provided that convincing proof of purchase in the qualifying period is provided. Philips recommends use of genuine Sonicare brush heads for optimal performance. Use of other brush heads may damage your Sonicare toothbrush and limit your warranty protection. Contact our Consumer Care Center at 1-800-682-7664 (North America), outside North America contact your local Philips Consumer Care Center. Internet information: **[www.sonicare.com](http://www.sonicare.com)** (North America) or **[www.philips.com/support](http://www.philips.com/support)** (outside North America).

### WARRANTY EXCLUSIONS

What is not covered under warranty:

- Brush heads.
- Damage caused by use of unauthorized replacement parts or unauthorized brush heads.
- Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.

### IMPLIED WARRANTIES

ANY IMPLIED WARRANTIES, INCLUDING IMPLIED WARRANTIES OF MERCHANTABILITY

AND FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED TO THE DURATION OF THE EXPRESS WARRANTIES SET FORTH ABOVE. IN SOME STATES LIMITATIONS ON DURATION OF IMPLIED WARRANTIES DO NOT APPLY.

## **LIMITATION OF REMEDIES**

IN NO EVENT SHALL PHILIPS OR ANY OF ITS AFFILIATED OR SUBSIDIARY COMPANIES BE LIABLE FOR ANY SPECIAL, INCIDENTAL OR CONSEQUENTIAL DAMAGES BASED UPON BREACH OF WARRANTY, BREACH OF CONTRACT, NEGLIGENCE, TORT, OR ANY OTHER LEGAL THEORY. SUCH DAMAGES INCLUDE, WITHOUT LIMITATION, LOSS OF SAVINGS OR REVENUE; LOSS OF PROFIT; LOSS OF USE; THE CLAIMS OF THIRD PARTIES INCLUDING, WITHOUT LIMITATION, DENTISTS AND DENTAL HYGIENISTS; AND COST OF ANY SUBSTITUTE EQUIPMENT OR SERVICES. SOME STATES DO NOT ALLOW THE EXCLUSION OR LIMITATION OF INCIDENTAL OR CONSEQUENTIAL DAMAGES.

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©2018 Koninklijke Philips N.V. (KPNV). All rights reserved. Philips and the Philips shield are trademarks of KPNV. Sonicare and the Sonicare logo are trademarks of Philips Oral Healthcare, LLC and/or KPNV.