

# PHILIPS

## sonicare

Sonic electric  
toothbrush with app

DiamondClean  
Smart

Personalized coaching  
Smart brush head sensors  
Smart brush head recognition  
5 modes, 3 intensities



HX9924/62

## Complete system for a healthier mouth

Our best toothbrush with smart sensor technology

Get the ultimate clean with DiamondClean Smart. Eliminate guesswork with real time feedback and personalized coaching for superior teeth, gum, whitening & fresh breath. It's more than just a toothbrush - it's a complete oral care solution.

### Proven to improve oral health

- Enjoy long-lasting fresh breath
- Up to 7x healthier gums\*\*\*\* in just 2 weeks
- All-in-One brush head for complete care

### Personalized in-app coaching to perfect brushing

- Guides you to brush and highlights missed areas
- Takes you back to any missed spots for 100% coverage
- Start taking better care of your trouble areas
- Alerts you if you're pressing too hard
- Guides you to reduce your scrubbing

### Optimize your brushing

- Brush heads automatically select the optimal settings
- Make sure you're getting the most from your brush head
- Choose from 5 modes and 3 intensity settings
- Premium charging glass and travel case


# DiamondClean Smart

Sonic electric toothbrush with app

HX9924/62


## Highlights

### Premium All-in-One brush head

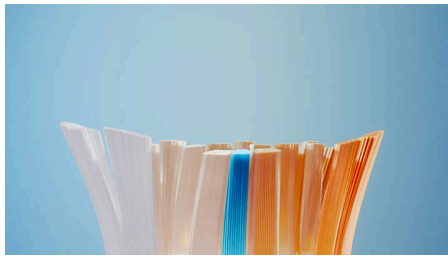

Clean without compromise with the A3 Premium All-in-One brush head. Our all-time best has angled bristles to help remove up to 20x more plaque even in hard-to-reach areas\*. Triangular tips work to remove up to 100% more stains in less than two days\*\*\*. And longer bristles clean deep for gums that are up to 15x healthier in just two weeks\*\*. All that without ever switching brush heads. Did you know: brush heads become less effective after 3 months of use? Our BrushSync™ feature can remind you when it's almost time for a replacement.

### Care for your gums

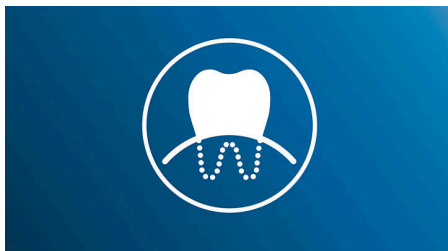

Click on the G3 Premium Gum Care brush head to improve your gum health. Its smaller size and targeted gum line bristles deliver gentle yet effective cleaning along the gum line, where gum disease starts. It's clinically proven to deliver up to 100% less gum inflammation\*\*\*\* and up to 7x healthier gums in just two weeks.\*\*\*\*


### TongueCare+ tongue brush



Click on the TongueCare+ tongue brush to gently remove odor causing bacteria from the pores of your tongue. Its 240 specially designed


MicroBristles get between all your tongue's ridges and grooves to remove bacteria and debris that lead to bad breath. Team with our antibacterial BreathRx tongue spray for superior cleaning and super fresh breath.

### Focus areas and goal-setting

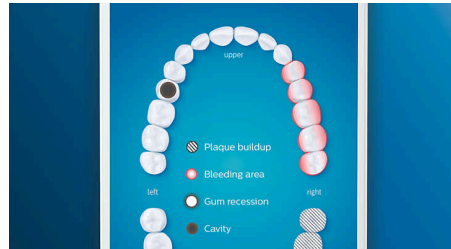


Any trouble areas your dentist has pointed out? We'll highlight them on your in-app 3D mouth map, and remind you to pay extra attention to these areas.

### Never miss a spot

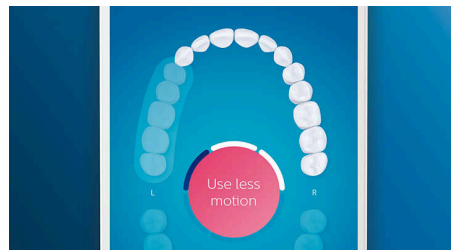


The Philips Sonicare app lets you know when you've achieved a thorough clean, and coaches you to be a better, more mindful brusher.

### Too much pressure?

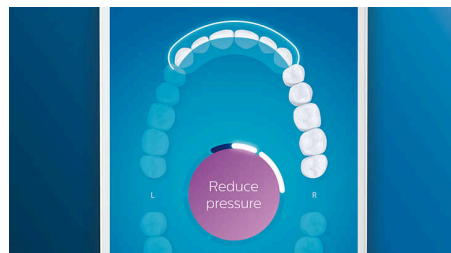


You might not notice when you brush too hard, but your DiamondClean Smart will. If you use too much pressure, the light ring on the end of your handle will flash. This is a gentle reminder to ease off the pressure, and let your brush head do the work. 7 out of 10 people found this feature helped them become a better brusher.

### Easy does it

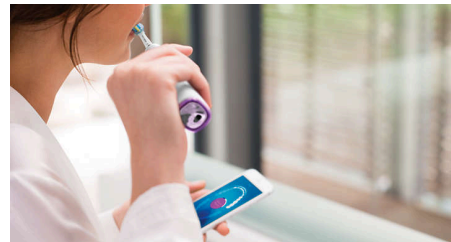


With an electric toothbrush, you let the brush do the work. We've built a scrubbing sensor into your handle as a gentle reminder to reduce scrubbing. This way, you can improve your technique and get a gentler, more effective clean.

### TouchUp feature

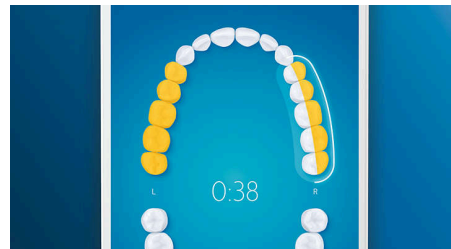

If you happen to miss a spot while brushing your teeth, your app's TouchUp feature will show you. You can then go back for a second pass, and be confident you are getting a complete clean, every time.

### BrushSync mode pairing

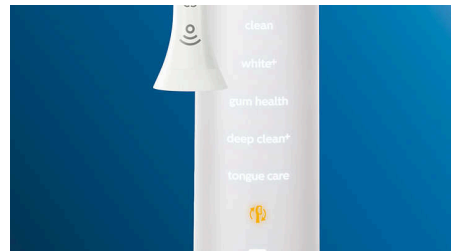

Wondering which mode and intensity you should use? Wonder no longer. A microchip in your brush head tells your DiamondClean Smart which one you're using. So if you click on a gum-care brush head, for example, your toothbrush knows to select the optimal mode and intensity to provide a gentle yet effective care for your gums. All you need to do is to press the power button.


# DiamondClean Smart

Sonic electric toothbrush with app

HX9924/62

## Highlights

### Smart replacement reminders

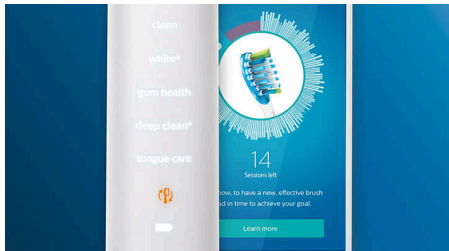

Brush heads wear-out over time and become less effective at removing plaque. With smart brush head recognition you'll never lose track of when to replace your brush head. Your toothbrush tracks how long and hard you've brushed with each brush head, and will alert you when it's time to replace. You can also monitor the lifetime of your brush heads in the Philips Sonicare app, and even order replacements so you're never caught short.

### 5 modes, 3 intensity settings



DiamondClean Smart comes with five modes which cover a range of brushing needs: Clean mode is for exceptional daily cleaning, White+ is for removing stains, Deep Clean+ gives you an invigorating deep clean, Gum Health provides a gentle yet effective clean for gums, and Tongue Care is the perfect mode for cleaning your tongue. It breaks up tongue debris and sweeps away bad-breath bacteria. Three intensities allow you to choose between a higher setting to boost your clean and a lower one for more-sensitive mouths.

### Design with storage in mind



Charging a toothbrush has never been so effortless. Simply place your toothbrush into the sleek glass holder after brushing, and it instantly starts charging. Not only does your glass charge fit stylishly into your bathroom, it also doubles as a rinsing glass. Enjoy two weeks of regular use from a single full charge. Our premium travel case lets you store your toothbrush hygienically when you're on the go.

# DiamondClean Smart

Sonic electric toothbrush with app

HX9924/62

## Specifications

### Connectivity

**Bluetooth® wireless technology:** Connected brushing app

### Power

**Voltage:** 110-220 V

### Technical specifications

**Operating time (full to empty):** 14 days\*\*\*\*\*

**Battery:** Rechargeable

**Battery type:** Lithium ION

### Design and finishing

**Color:** Rose gold

### Service

**Warranty:** 2-year limited warranty

### Compatibility

**Android compatibility:** Android phones,

Bluetooth 4.0 enabled tablets

**iOS compatibility:** iPad 3rd Gen or higher, iPhone

4S or higher, with iOS7 or higher

### Ease of use

**Display:** Illuminated display

**Battery indicator:** Illuminated icon indicates battery life

**Handle:** Slim ergonomic design

**Handle compatibility:** Easy click-on brush heads  
**Brush head recognition:** Syncs with optimal mode

**Replacement reminder:** To always ensure best results, reminder icon lights up

### Items included

**Handle:** 1 DiamondClean Smart

**Travel case:** Premium travel case

**Brush heads:** 2 A3 Premium All-in-One, 1 G3 Premium Gum Care, 1 TongueCare+ tongue brush

**Glass charger:** 1

### Cleaning performance

**Speed:** Up to 62000 brush movement/min

**Health benefits:** Up to 15x healthier gums\*\*

**Performance:** Removes up to 20x more plaque\*

**Timer:** BrushPacer and SmarTimer

**Whitening benefits:** Up to 100% less stains\*\*\*

**Pressure feedback:** Ring lights up in purple, Vibrates handle to alert user

**TouchUp:** Ensures 100% coverage

### Modes

**3 intensities:** Low, Medium, High

**Clean:** For exceptional everyday clean

**Deep Clean+:** For an invigorating deep clean

**Gum Health:** Special attention to molars

**TongueCare+:** For long lasting clean breath

**White+:** To remove surface stains

### BrushSync Mode Pairing

**C3 Premium Plaque Control:** Pairs with Clean mode

**G3 Premium Gum Care brush head:** Pairs with Gum Health mode

**W3 Premium White brush head:** Pairs with White+ mode

**TongueCare+ tongue brush:** Pairs with TongueCare mode

### Smart sensor technology

**Pressure sensor:** Alerts when brushing too hard

**Location sensor:** Tracks & improves coverage

**Scrubbing sensor:** Guides to reduce scrubbing

**3D mouth map:** Highlights trouble areas

**TouchUp:** To address missed spots

