
PHILIPS


Sonicare
for Kids

PHILIPS
sonicare

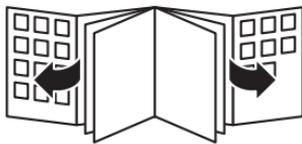1

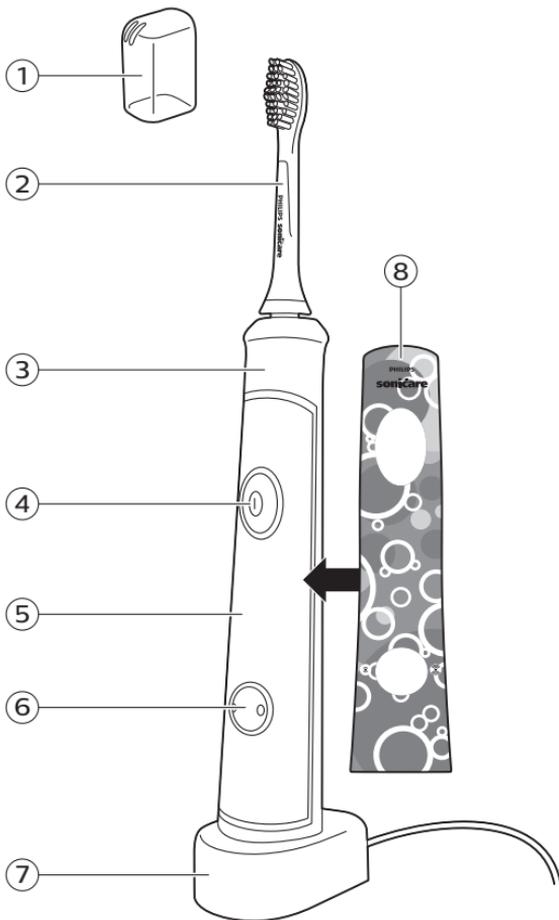— |

| —

— |

| —

English 6
Français (Canada) 15

IMPORTANT SAFEGUARDS

READ ALL INSTRUCTIONS BEFORE USE

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

DANGERS

To reduce the risk of electrocution:

- Always place and store the charger at a safe distance from water.
- Do not reach for a charging pad that has fallen into water. Unplug immediately.
- Never use a charging pad with a damaged cord.
- Do not place or store product where it can fall or be pulled into a tub or sink
- Do not place in or drop into water or other liquid.
- Do not reach for a product that has fallen into water. Unplug immediately.
- Do not use while bathing.
- After cleaning, make sure the charger and/or sanitiser is completely dry before you connect it to the mains.

WARNINGS

To reduce the risk of burns, electrocution, fire or physical injury:

- Do not clean the brush head, the handle or the charger in the dishwasher.
- Do not use attachments other than those recommended by the manufacturer.
- Do not plug into an outlet with a voltage other than specified on the bottom of the product. This product is designed to operate within a range of 100 to 240 volts. Voltage converters do not warrant voltage compatibility.
- Never force the plug into an outlet; if the plug does not easily fit into the outlet, discontinue use.
- Discontinue use if the product appears damaged in any way (brush head, handle or charger, etc.).

- The power cord cannot be replaced. If the mains cord is damaged, discard the charger.
- Keep the power cord away from heated surfaces.
- Do not use the base charger outdoors or near heated surfaces.
- This product contains no user-serviceable parts. Refer to 'Warranty and support' if the product no longer works properly or needs repair.
- Keep the Philips Sonicare out of direct sunlight.
- Use this product only for its intended use as described in this document or as recommended by your dental professional.
- This product is not intended for use by persons (including children) with reduced physical, sensory, or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the product.
- Philips Sonicare is a personal care device and is not intended for use on multiple patients in a dental office or institution.
- The product is not intended to be repaired by service personnel in case of failure or component defect (unit can be thrown away).
- For connection to a supply not in the U.S.A., use a charger of the proper configuration for the power outlet.
- Use a UL listed computer, USB Hub, or USB duplex receptacle when charging the handle.
This appliance has only been designed for cleaning teeth, gums and tongue. Do not use it for any other purpose. Stop using the appliance and contact your doctor if you experience any discomfort or pain.
- Avoid direct contact with products that contain essential oils or coconut oil. Contact may result in bristle tufts dislodging.
- Never use the charger if it has a damaged cord or plug.
- Do not use the charger if dropped into water.
- If your toothpaste includes peroxide, baking soda or other bicarbonate (common in whitening toothpastes), thoroughly clean your brush head and handle with soap and water after each use. Otherwise plastic cracking may occur.

MEDICAL WARNINGS

- Consult your dental professional before you use this product if you have had oral or gum surgery in the previous 2 months.
- Contact your dental professional if excessive bleeding occurs after using this product or bleeding continues to occur after 1 week of use.
- This product complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician or the device manufacturer prior to use.
- Consult your physician if you have other medical concerns.
- Sonicare has been tested and is compliant with safety standards for electromagnetic devices.

SAVE THESE INSTRUCTIONS

Electromagnetic fields (EMF)

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

Introduction

Congratulations on your purchase and welcome to Philips! To fully benefit from the support that Philips offers, register your product at www.philips.com/welcome.

You and your child now own a Sonicare connected toothbrush. The Bluetooth connected toothbrush interacts with a fun brushing mobile device application to keep your child engaged while learning proper brushing techniques.

General description (Fig. 1)

- 1 Hygienic travel cap
- 2 Brush head
- 3 Handle with soft grip
- 4 Power on/off button with battery charge indicator
- 5 Panel for your favorite sticker
- 6 Brushing mode button with mode lights
- 7 Charger

8 Sticker

Preparing for use

Applying your favorite sticker.

The Sonicare For Kids comes with replaceable stickers.

- 1 Apply your favorite sticker on the handle (Fig. 2).

Note: Make sure you line up the cut-out areas on the sticker with the buttons on the front of the handle.

Attaching the brush head

- 1 Align the brush head so the bristles point in the same direction as the front of the handle (Fig. 3).
- 2 Firmly press the brush head down onto the metal shaft until it stops.

Note: There is a small gap between the brush head and the handle.

Charging your Sonicare

- 1 Plug the charger into an electrical outlet.
- 2 Place the handle on the charger (Fig. 4).
 - The battery charge indicator in the power on/off button flashes and the handle beeps twice to indicate that the toothbrush is charging.
 - When the Philips Sonicare is fully charged, the battery charge indicator stops flashing and remains on.

Note: If the battery charge of the Philips Sonicare is low, you hear 3 beeps and the battery charge indicator flashes rapidly for 30 seconds after the brushing cycle.

Note: To keep the battery fully charged at all times, you may keep the Philips Sonicare on the charger when not in use. When the battery is fully charged, the charger uses a minimal amount of energy. It takes at least 24 hours to fully charge the battery.

Using your Philips Sonicare

For children under age 8, it is recommended that a parent or adult caregiver perform toothbrushing.

Brushing instructions

- 1 Wet the bristles and apply a small amount of toothpaste.
- 2 Place the toothbrush bristles against the teeth at a slight angle (45 degrees), pressing firmly to make the bristles reach the gumline or slightly beneath the gumline. (Fig. 5)

Note: Keep the center of the brush in contact with the teeth at all times.

- 3 Press the power on/off button to turn on the Philips Sonicare.
- 4 Apply light pressure to maximize Philips Sonicare's effectiveness and let the Philips Sonicare toothbrush do the brushing for you. Do not scrub.
- 5 Gently move the brush head slowly across the teeth in a small back and forth motion so the longer bristles reach between your teeth. Continue this motion throughout the brushing cycle.
- 6 To clean the inside surfaces of the front teeth, tilt the brush handle semi-upright and make several vertical overlapping brushing strokes on each tooth (Fig. 6).

Note: To make sure you brush evenly throughout the mouth, divide your mouth into 4 sections using the Quadpacer feature.

- 7 Begin brushing in section 1 (outside upper teeth) and brush for 30 seconds before you move to section 2 (inside upper teeth). Continue brushing in section 3 (outside lower teeth) and brush for 30 seconds before you move to section 4 (inside lower teeth) (Fig. 7).
- 8 After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs. You can also brush your tongue, with the toothbrush turned on or off, as you prefer. (Fig. 8)

Your Philips Sonicare is safe to use on:

- Braces (brush heads wear out sooner when used on braces).
- Dental restorations (fillings, crowns, veneers).

Note: When the Sonicare is used in clinical studies, the handle should be fully charged. Deactivate the Easy-start feature and for areas where excess staining occurs, an additional 30 seconds of brushing time can be spent to assist with stain removal.

Brushing modes

The Philips Sonicare automatically starts in the default high mode.

- 1 Press the brushing mode button to toggle between modes (Fig. 9).

The green mode light indicates the selected mode.

- High mode: Full-strength mode designed for older/independent brushers.
- Low mode: Reduced power mode used for training and introducing younger children to the sonic experience.

Features

KidTimer

- The KidTimer plays a congratulatory sound sequence at the end of the pre-determined brushing time. Additionally, the KidTimer automatically turns off the handle when the brushing cycle is complete.
- When in high mode, the KidTimer is set at 2 minutes.
- When in low mode the KidTimer is set at 1 minute and slowly increases over time. The progression of time helps young children reach the dentist recommended 2 minutes of brushing.

Note: Instruct your child to brush until they hear the congratulatory tune and the toothbrush turns off.

KidPacer

The Sonicare For Kids uses a short series of tunes to alert the brusher to move to the next brushing quadrant of the mouth and concludes the brushing cycle with a congratulatory tune before it automatically turns off. This ensures that the child's teeth get a thorough, overall clean (see "Brushing instructions").

Note: Instruct your child to brush until they hear the congratulatory tune and the toothbrush turns off.

Easy-start

The Easy-start feature gently increases the power over the first 14 brushings to help you get used to the brushing with the Philips Sonicare. The Sonicare For Kids comes with the Easy-start feature activated.

12 English

Note: Each of the first 14 brushings must be at least 1 minute in length to properly advance through the Easy-start ramp-up cycle.

Activating or deactivating the Easy-start feature

- 1 Place the handle on the charger.

To activate Easy-start: Press and hold the power on/off button for 2 seconds. You hear 2 beeps to indicate that the Easy-start feature has been activated. To deactivate Easy-start: Press and hold the power on/off button for 2 seconds. You hear 1 beep to indicate that the Easy-start feature has been deactivated.

Note: Using the Easy-start feature beyond the initial ramp-up period is not recommended and reduces Philips Sonicare effectiveness in removing plaque.

Cleaning

Warning: Do not clean brush heads, handle, charger, charger cover or replaceable panel in the dishwasher.

Toothbrush handle

- 1 Remove the brush head and rinse the metal shaft area with warm water. Make sure you remove any residual toothpaste (Fig. 10).

Caution: Do not push on the rubber seal on the metal shaft with sharp objects, as this may cause damage.

- 2 Wipe the entire surface of the handle with a damp cloth.

Brush head

- 1 Rinse the brush head and bristles after each use (Fig. 11).
- 2 Remove the brush head from the handle and rinse the brush head connection with warm water at least once a week.

Charger

- 1 Unplug the charger before you clean it.
- 2 Wipe the surface of the charger with a damp cloth. Make sure you remove any toothpaste or other residue from the charger surfaces.

Storage

If you are not going to use the Sonicare for an extended period of time, unplug the charger from the electrical outlet, clean it and store it in a cool and dry place away from direct sunlight.

Replacement

Brush head

- Replace Sonicare brush heads every 3 months to achieve optimal results.
- Use only Sonicare For Kids replacement brush heads.
- Brush head can be thrown away with the normal household waste.

Recycling

- This symbol means that this product shall not be disposed of with normal household waste (2012/19/EU) (Fig. 16).
- This symbol means that this product contains a built-in rechargeable battery which shall not be disposed of with normal household waste (Fig. 17) (2006/66/EC). Please take your product to an official collection point or a Philips service center to have a professional remove the rechargeable battery.
- The UV light bulb of the sanitizer contains mercury. Do not throw the UV light bulb away with the normal household waste at the end of its life. Hand it in at an official collection point for recycling.
- Follow your country's rules for the separate collection of electrical and electronic products and rechargeable batteries. Correct disposal helps prevent negative consequences for the environment and human health.

Removing the rechargeable battery

Warning: Only remove the rechargeable battery when you discard the appliance. Make sure the battery is completely empty when you remove it.

To remove the rechargeable battery, you need a flat-head (standard) screwdriver. Observe basic safety precautions when you follow the procedure outlined below. Be sure to protect your eyes, hands, fingers, and the surface on which you work.

14 English

- 1 To deplete the rechargeable battery of any charge, remove the handle from the charger, turn on the Philips Sonicare and let it run until it stops. Repeat this step until you can no longer turn on the Philips Sonicare.
- 2 Insert a screwdriver into the slot located in the bottom of the handle and turn counterclockwise until you see a gap between the bottom cap and the handle (Fig. 12).
- 3 Insert the screwdriver into the gap and pry the bottom cap out of the handle (Fig. 13).
- 4 Hold the handle upside down and push down on the shaft to release the internal components of the handle (Fig. 14).
- 5 Insert the screwdriver under the circuit board, next to the battery connections, and twist to break the connections. Remove the circuit board and pry the battery from the plastic carrier (Fig. 15).

The rechargeable battery can now be recycled or disposed of and the rest of the product discarded appropriately.

Warranty and support

If you need information or support, please visit

www.philips.com/support or read the international warranty leaflet.

Warranty restrictions

The terms of the international warranty do not cover the following:

- Brush heads.
- Replaceable stickers.
- Damage caused by use of unauthorized replacement parts.
- Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.

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©2018 Koninklijke Philips N.V. (KPNV). All rights reserved. Philips and the Philips shield are trademarks of KPNV. Sonicare and the Sonicare logo are trademarks of Philips Oral Healthcare, LLC and/or KPNV.