

**PHILIPS**

VisaCare

SC6240

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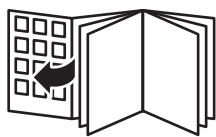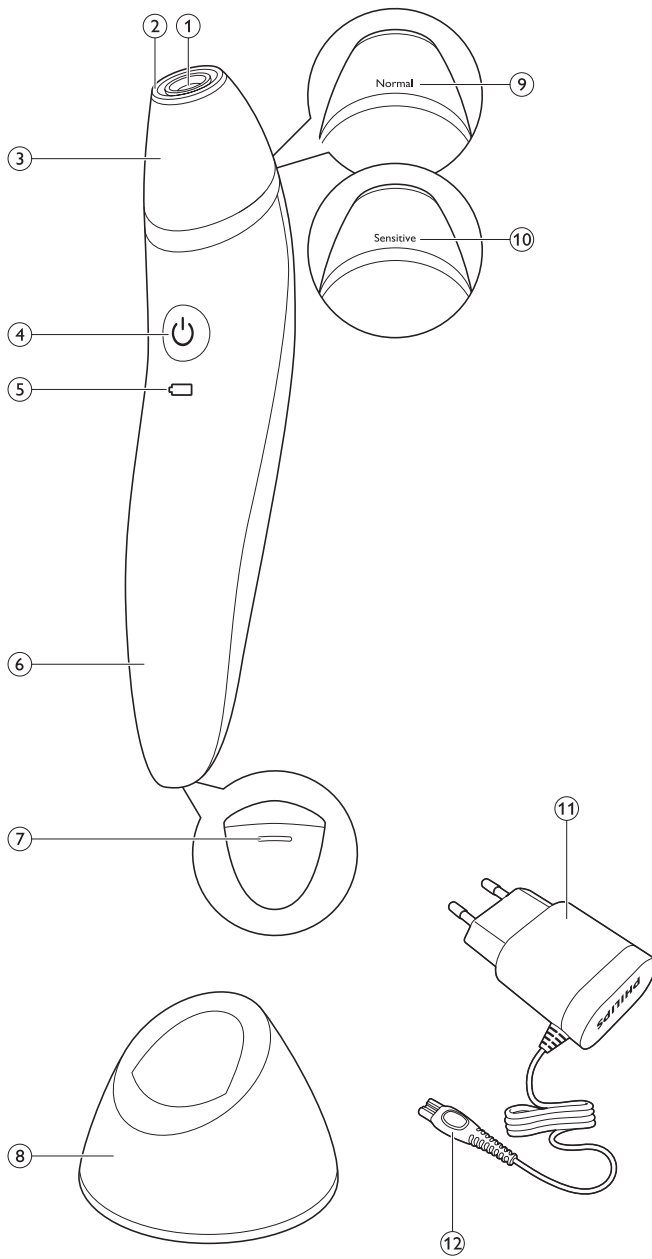## Introduction

Congratulations on your purchase and welcome to Philips! To fully benefit from the support that Philips offers, register your product at [www.philips.com/welcome](http://www.philips.com/welcome).

What is VisaCare and how can it help your skin? Age, lack of sleep, emotional stress or environmental elements can take their toll on your skin. When you start to see the first signs of skin ageing and your skin does not seem to recover so quickly anymore, your skin requires extra care and attention. That is why we have developed VisaCare, our exciting new approach to address the first signs of ageing, for a youthful looking skin today and tomorrow. Depending on your skin, you can experience immediate effects such as smoother and revitalised skin after one session or a few sessions. Effects such as skin firmness and a more even skin tone can be expected within 4-6 weeks.

You can use VisaCare in the comfort of your own home to enjoy smoother, firmer and revitalised skin. VisaCare comes with a dual-action system that exfoliates and massages at the same time. This system is more effective and controlled than exfoliating your skin with your hands only. Designed to enhance your current routine, use VisaCare after cleansing your skin and before applying a nourishing mask, day cream or night cream. After use, your skin is prepared for better absorption of the mask or cream. Ready to get started? Use VisaCare at home twice a week to enjoy real and visible results.

## How Philips Microdermabrasion technology works

Philips VisaCare offers a DualAction Air Lift & Exfoliation system, which massages and exfoliates at the same time. This stimulates the resurfacing process of the skin for firmer, smoother and visibly radiant skin.

- 1 The exfoliating tip renews and unveils skin that is smoother, with more refined texture, and firmer by stimulating your natural cell renewal process. It gently nudges away roughness, flakiness to refine the skin and to prepare the skin for better absorption of your favorite cream after use.
- 2 The Air Lift System™ creates a gentle, stimulating massaging action encouraging blood circulation while removing dead skin cells and revitalising your skin.

## General description (Fig. 1)

- 1 Filter
- 2 Exfoliation ring
- 3 Microdermabrasion tip
- 4 On/off button
- 5 Charging light and 'battery low' indication
- 6 Handle
- 7 Air outlet
- 8 Charging stand
- 9 'Normal' tip for regular use
- 10 'Sensitive' tip for mild use
- 11 Adapter
- 12 Small plug

- Pouch (not shown)
- Cleaning brush (not shown)

## Important

Read this user manual carefully before you use the appliance and save it for future reference.

## Danger

- Keep the adapter dry.
- Do not charge this appliance in wet surroundings (e.g. near a filled bath, a running shower or near a filled swimming pool).

## Warning

- To charge the battery, only use the detachable supply unit (SSW-2600) provided with the appliance.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- Do not use the appliance or the supply unit if it is damaged or broken to avoid injuries.
- The adapter contains a transformer. Do not cut off the adapter to replace it with another plug, as this causes a hazardous situation.
- If the adapter is damaged, always have it replaced with one of the original type in order to avoid a hazard.

- Do not use the appliance if the tip is broken, cracked, damaged or dirty. Always inspect the tip before use. When in doubt, replace the tip.
- Do not insert metal-containing material into the socket for the small plug to avoid short circuiting.
- This appliance contains a magnet. Keep magnets away from cardiac pacemakers. Persons with cardiac pacemakers or other devices implanted in their bodies should not come into close proximity with magnets. Specialist medical opinion must be sought before such persons handle or come into contact with magnets.

### Caution

- Do not use the appliance for any other purpose than to massage and gently remove the dead upper skin layer to speed up normal skin turnover.
- Do not use the appliance more than twice a week, as this may cause skin irritation.
- Do not use the appliance in combination with (scrub) creams or with chemical or mechanical peelings, as this may cause damage to the appliance and your skin.
- Use, charge and store the appliance at a temperature between 5°C/40°F and 35°C/95°F.
- Fully charge the appliance at least every 3 months to sustain battery lifetime.

- VisaCare enhances the absorption of your creams and facial masks. Due to this enhanced uptake, a negative skin reaction can occur. The advice is to not change between creams during the first 2-4 weeks of using VisaCare. When you start to use a new cream after this period, first test this cream on a small area. Do not continue to use this cream if there is a negative skin reaction when you apply it on the skin on which you used VisaCare. When you start to use the device, or when you restart after discontinuing for more than 2 weeks, begin again with the sensitive tip. This will allow your skin to once again become used to the device.

## Contraindications and other warnings

This device has been developed to address the first signs of facial ageing. As the epidermal skin layer becomes thinner over the years due to extrinsic and intrinsic factors, rough skin and signs of photodamage become more prominent. The focus of this device is facial skin (excluding the area underneath and above the eyes).

### **Do not use this appliance:**

- If you have skin infections, open wounds, rashes, burns, sunburn, inflammation, eczema, psoriasis, rosacea, telangiectasia (couperose or spider veins), sensitive skin, active infections like herpes or active acne breakouts in the areas to be treated.
- If you have large moles, tattoos, port wine stains, permanent make-up or any skin irritation in the areas to be treated.
- If your skin is still recovering from a cosmetic procedure, such as a chemical or mechanical peeling, laser resurfacing, fractional laser treatment, or mesotherapy treatment or any other treatment with dermal fillers in the past two weeks.
- If you are taking medication that may affect your skin, such as steroids or acne medication (Accutane, Isotretinoine or Roaccutane). Consult your doctor.
- If you have warts, do not try to erase them with this device. It can spread the warts.
- If you have scars, consult your doctor or dermatologist before you start using the device.
- In combination with any other home skincare system containing an active peeling ingredient.

### **General advice**

- After each session it is advised to use a sun protection factor (SPF) of 30 or above on the areas on which you used the device, for a minimum of 3 days.

- If you have a significant skin condition (including but not limited to melanoma, basal cell carcinoma, actinic keratosis, melasma, systemic lupus erythematosus, porphyria, coagulation abnormality and collagen abnormalities such as keloid wounds, excessive scar formation or a history of bad wound healing), consult your doctor before you start using the device.
- If you plan to enjoy the sun (e.g. summer holidays) we recommend that you stop using the Firming tip and only use the Brightening tip one week before the holidays and as long as you are extensively exposed to the sun.

## Electromagnetic fields (EMF)

This Philips device complies with all applicable standards and regulations regarding exposure to electromagnetic fields.


## Charging

### Charging the appliance


Note: You can use the appliance before the battery is fully charged. However, its performance and operating time will not be optimal in this case.


Note: Charging normally takes approximately 8 hours. When the appliance is fully charged, it has a cordless operating time of up to 30 minutes.

- 1 Make sure the appliance is switched off.
- 2 Insert the small plug in the socket of the charging stand and insert the adapter in the wall socket.

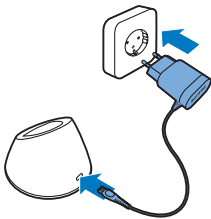- 3 Put the appliance in the charging stand.

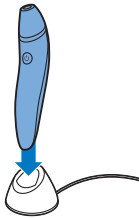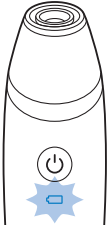- The charging light flashes white to indicate that the appliance is charging.

- When the battery is fully charged, the charging light lights up white continuously. After approx. 30 minutes, the charging light goes out to save energy.

## 'Battery low' indication

When the battery is almost empty, the charging light flashes amber for 4 seconds during use and after you have switched off the appliance. The battery still contains enough energy for a complete 5-minute session. When the battery is completely empty, the device switches off automatically and the charging light flashes amber very quickly.

## Preparing for use

Make sure you have charged the appliance for 8 hours before you use it for the first time.

### Preparation

Always clean your face and gently dry it before you use the appliance.

**Note:** The appliance is intended for dry use only. Do not use it on a wet face or in the shower or bathtub.

### Application tips

VisaCare comes with two different tips:

- 'Normal' tip for regular use.
- 'Sensitive' tip for mild use.

We advise you to use the 'Sensitive' tip during the first sessions to see how your skin reacts. If the skin reacts well to the 'Sensitive' tip, you can try the 'Normal' tip. Normal and adverse skin reactions are described elsewhere in this user manual (see 'Troubleshooting').

**Note:** Use the QR-code on the front cover or follow the link to watch the online video for instructions on correct use of the appliance.

## Using the appliance

Note: Use the appliance twice a week. Do not use the appliance more often, as this may cause skin irritation.

1 Press the on/off button to switch on the appliance.

2 When you use the appliance for the first time, test the tip on your forearm before you use it on your face. Move the appliance over your skin once and then wait for approximately one hour. If a possible skin reaction diminishes, you can safely use the tip on your face.

Note: Always visually check the tip before use to see if it is damaged. A damaged tip scratches your skin. If the tip is damaged, replace it.

Note: Do not use the device without the tip as this does not benefit your skin and may cause an unwanted skin reaction.

3 Put the tip on your face in such a way that you feel the suction. Make sure there is always full skin contact between the tip and the area to be treated. However, do not press the tip hard onto the skin's surface.

4 Gently move the tip in straight lines. Start on the chin, and then move to the cheeks, forehead and nose. make sure you treat all areas.

Tip: Use your free hand to stretch and hold the skin. Move the device away from your hand. This makes it easier to achieve proper suction, especially on the temples and cheeks.

See facial map in the back of the user manual for detailed instructions on how to move the appliance across your face.

Note: Do not move the appliance over large moles, tattoos, port wine stains or permanent makeup.

Note: Do not use the appliance on the skin below or above the eyes, as it is thin and delicate. When you use the appliance on the fine lines on the side of the eye, put a finger of the other hand in the corner of the eyelid and move the device away from the finger towards the temple.

Note: Do not hold the treatment tip in one spot when you have placed it on the skin.

5 Depending on the type of tip and the condition of your skin, start to treat each area 2-4 times. When you skin is used to the treatment, increase to 4-6 times or start to use the normal treatment tip. When you use the normal treatment tip, start to treat each area 3-4 times. When your skin is used to the treatment, increase to 4-6 times.

- 6 When you have treated all areas, press the on/off button to switch off the appliance. A complete session takes approximately 5 minutes.

Depending on your skin type, you may experience a warm and tingling sensation and your skin may become slightly red. This is only a temporary effect that should diminish within 5-10 minutes.

**Tip:** In areas where full skin contact is harder to achieve, such as the forehead, try different ways of holding the appliance. Try vertical and horizontal movements to find the best way to establish proper contact with your skin. You can also hold your hand still and move your head instead.


**Caution:** If you suspect damage to the appliance (such as after a drop or bump), examine the tip carefully. If the tip is damaged, replace it. Do not use a damaged tip, as this may cause scratches.

## After treatment


After the treatment, apply a nourishing mask, day cream or night cream and make sure you use sun protection when you go outside. Do not use any other method of mechanical or chemical exfoliation for 3 days. Do not undergo electrolysis, facial waxing or use depilatories within 3 days after you have used the appliance. You can apply make-up immediately after the treatment.

**Caution:** Do not immerse the appliance in water. Do not use chemicals to clean the appliance.


**Warning:** Never use scouring pads, abrasive cleaning agents or aggressive liquids such as petrol or acetone to clean the handle or tips.

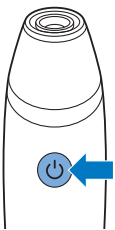- 1 Make sure the appliance is switched off before you start cleaning it.

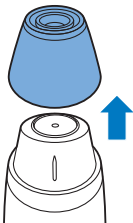- 2 Remove the tip from the appliance.

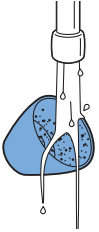- 3 Clean the tip under the tap with hot water (max. 70°C) and some washing-up liquid. If the exfoliation ring is dirty, let the tip soak in hot water with washing-up liquid or in alcohol solution for 5-10 minutes. You can use the brush to clean it. Make sure you clean the filter of the tip upside down to rinse off any dirt left in the filter. Dry the tip with a lint-free towel or leave it to air dry. Place it back onto the handle when it is dry.

- 4 Clean the handle with a damp cloth.
- 5 Place the dry treatment tip back onto the appliance.
- 6 Wipe the charging stand with a damp cloth.

**Danger:** Keep the adapter and the charging stand dry. Never immerse them in water or rinse them under the tap.

Clean the treatment tip after each use for hygienic reasons and to prevent clogging of the filter. Only use washing up-liquid or an alcohol solution (max 70%) to clean the treatment tip.

## Storage

Store the appliance with the tip attached in the charging stand or pouch.

**Note:** Do not store a wet treatment tip on the appliance or in the pouch.

## Replacement

For optimal results, it is recommended to replace the tips every six months (or earlier if the abrasive structure is deformed or damaged). Replacement tips are available on our website. The tip can be identified by the product code. You can find this code on the inside of the tip.

## Ordering accessories

To buy accessories or spare parts, visit [www.shop.philips.com/service](http://www.shop.philips.com/service) or go to your Philips dealer. You can also contact the Philips Consumer Care Center in your country (see the international warranty leaflet for contact details).

## Recycling

- This symbol means that this product shall not be disposed of with normal household waste (2012/19/EU).


- This symbol means that this product contains a built-in rechargeable battery which shall not be disposed of with normal household waste (2006/66/EC). Please take your product to an official collection point or a Philips service center to have a professional remove the rechargeable battery.
- Follow your country's rules for the separate collection of electrical and electronic products and rechargeable batteries. Correct disposal helps prevent negative consequences for the environment and human health.


## Removing the rechargeable battery

**Only remove the rechargeable batteries when you discard the appliance. Before you remove the batteries, make sure that the appliance is disconnected from the wall socket and that the batteries are completely empty.**

**Take any necessary safety precautions when you handle tools to open the appliance and when you dispose of the rechargeable battery.**

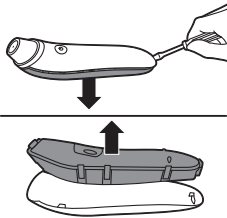- 1 Remove the back and front panel of the appliance with a screwdriver and take out the inner housing of the appliance.

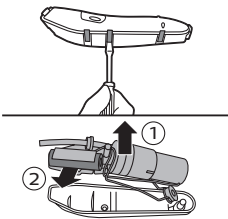- 2 Remove the clips from the inner body using a screwdriver.
- 3 Remove the electronics from the inner housing (1) and detach the rechargeable battery (2).

## Warranty and support

If you need information or support, please visit [www.philips.com/support](http://www.philips.com/support) or read the international warranty leaflet.

## Troubleshooting

This chapter summarizes the most common problems you could encounter with the appliance. If you are unable to solve the problem with the

information below, visit [www.philips.com/support](http://www.philips.com/support) for a list of frequently asked questions or contact the Consumer Care Center in your country.

Problem	Solution
I do not know if the appliance is suitable for use on my skin.	Please check if the appliance is suitable for your skin.
The appliance does not charge.	Make sure the socket to which you connect the appliance is live. Check if the voltage indicated on the appliance corresponds to the local mains voltage. Check if you have placed the appliance in the stand properly. If you use a socket in a bathroom cabinet, you may need to switch on the light to activate the socket. If the light on the appliance still does not light up or if the appliance still does not charge, contact the Philips Consumer Care Center in your country.
The appliance does not work when it is charging.	This is a safety measure. The appliance is not intended to work when it is charging. Let the appliance charge approximately 1 hour for a 5-minute session.
I don't get optimal results when I use the appliance on wet skin.	Only use the appliance on dry skin. The appliance only works optimally on dry skin.
My skin is reacting to usage of the appliance.	Using this appliance may lead to mild redness directly after use. You may experience a warm, tingling sensation during use. These effects usually disappear within several minutes and are generally experienced as pleasant. It is possible that you experience dry or flaky skin, prolonged redness or swelling. In that case you may have pressed the appliance too hard on your face or have passed it over the same area too many times. Try to exert less pressure onto the appliance during the next session or pass the appliance over the area fewer times. In rare cases the appliance may cause a bruise or blister, scratch or pigmentation change. Almost all skin reactions will disappear within a week, pigmentation change will decrease within several weeks or a few months. To avoid possible side effects, use a day cream with sun protection factor or a sun protection cream after the session. This helps to avoid sunburn. The appliance has been tested and confirmed safe for use when used as described in the user manual. Medicines or products and their possible side effects can affect individual people in different ways. Since there is always the possibility that a rare or previously unknown side effect may occur with the use of this appliance or any other product, please contact the Philips Consumer Care Center in case you experience a side effect which is not mentioned above.
It is difficult to slide the appliance over my skin/I do not feel the suction/the effect seems to be less than after previous sessions.	Check if the filter is clogged. Clean the tip by rinsing it under the tap. If the exfoliation ring is dirty, let the tip soak in hot water with washing-up liquid for 5-10 minutes. You can use the cleaning brush to clean it. Make sure you clean the filter with the tip upside down to rinse off any dirt left in the filter. Dry the tip with a lint-free towel or let it air-dry. Check if you can feel the suction of the appliance when the tip is taken off: you remove the tip and then switch on the appliance. When you put your finger on the hole, you should be able to feel the suction. If you do not feel the suction, the appliance might be broken. Please contact your Philips Consumer Care Center in your country.

<b>Problem</b>	<b>Solution</b>
I have trouble maintaining the suction on my forehead.	Try different ways of holding the appliance in order to find the right angle of the tip with the forehead for proper skin contact. Make sure to move the appliance slowly over the forehead. You can try to hold your hand still and move your head instead.
I do not know when to start using the ' Normal' tip.	As soon as you do not experience side effects from using the 'Sensitive' tip, you may try the ' Normal' tip. Please note that the first time you use a new tip, you may have a different skin reaction. Give your skin some time to get used to the new tip.

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© 2017 Koninklijke Philips N.V. All rights reserved. Specifications are subject to change without notice. Trademarks are the property of Koninklijke Philips N.V. (Royal Philips) or their respective owners.