
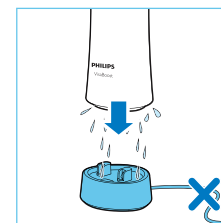
Register your product and get support at  
在以下网站注册并获得支持  
[www.philips.com.cn/welcome](http://www.philips.com.cn/welcome)

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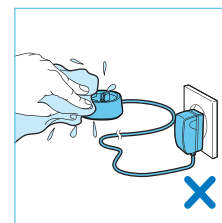3140 035 42662

EN User manual

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PHILIPS


**EN** Apply preferred serum, cream or gel to your cleansed face, or moisten the probe with your skincare products before using.

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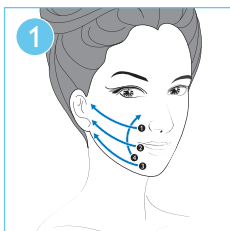
**EN** Start from the corner of your mouth and move toward your ear, and glide the probe from chin to the lower part of your cheekbone. Repeat each step 2-3 times with each motion lasting for 5 sec.

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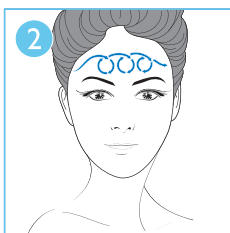
**EN** Make circular movements on your forehead from left temple to the right. Repeat it 2-3 times.

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
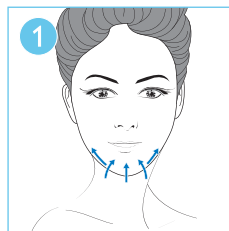
**EN** Glide the probe over your chin by following upward directions. Repeat steps 3 times. Each motion should last for 5 sec.

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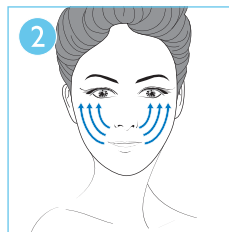
**EN** Move the probe with upward directions around your cheekbone. Repeat this step 2-3 times with each lasting for 5 sec.

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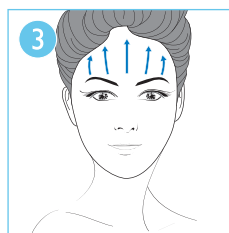**EN** Put the probe above your brow, and then glide to your hairline. Repeat this step 2-3 times with each motion lasting for 5 sec.

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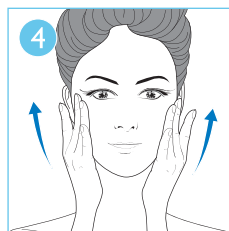**EN** Massage your skin with palms when finished, preferably using an upward direction.

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Congratulations on your purchase, and welcome to Philips!  
To fully benefit from the support that Philips offers, register your product at [www.philips.com.cn/welcome](http://www.philips.com.cn/welcome).

## Important

Read this user manual carefully before you use the appliance and keep it for future reference.

### Danger

- Keep the adapter dry. Never immerse the handle into water.
- Do not charge this appliance in wet surroundings (e.g. near a filled bath, a running shower or near a filled swimming pool).
- Never throw the appliance into fire.

### Warning

- Only use the appliance in combination with the adapter supplied.
- Do not use the appliance, the attachment or the adapter if it is damaged or broken to avoid injuries.
- If the adapter is damaged, always have it replaced with one of the original type in order to avoid a hazard.
- Do not cut off the adapter to replace it with another plug, as this causes a hazardous situation.
- Fully charge the appliance at least every 3 months to sustain battery lifetime.

- This appliance is not intended for use by persons younger than 12 years and persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- This appliance is NOT a toy. Keep the appliance out of the reach of children younger than 12 years old.
- Children should be supervised to ensure that they do not play with the appliance.
- Do not insert metal-containing material into the socket of the appliance plug to avoid short circuiting.
- Women who are pregnant or lactating are not recommended to use this appliance.

### **Caution**

- Do not use the appliance for any other purpose than described in this manual.
- Charge, use and store the appliance at a temperature between 5°C/40°F and 35°C/95°F.
- Do not use the appliance if you are taking steroid based medication.
- Do not use this appliance if your skin is still recovering from a cosmetic procedure,

### 3 ENGLISH

such as a chemical or mechanical peeling, laser resurfacing, fractional laser treatment, or mesotherapy treatment or any other treatment with dermal fillers in the past two weeks.

#### **General**

- Regular cleaning of the appliance ensures optimal results and a longer life of the appliance.
- Do not use when sleeping or drowsy.

#### **Contraindications and other warnings**

- Do not apply this device directly to your ear as this may cause temporary hearing loss.
- Do not use this device if you have any ear problems (Swimmer's ear, ear infection, any anomaly).
- Do not use the appliance on chapped skin, open wounds, or if you suffer from skin diseases or skin irritations, such as severe acne, sunburn, skin infection, etc.
- Do not use this device if you have a pacemaker, implanted device, or have a metal prosthesis.
- Do not use the appliance on the sensitive skin area around the eyes, as it is thin and delicate.
- Do not use the appliance on any mucosal area (e.g. inside nose, mouth, genital area) as it may cause local irritation.

## Electromagnetic fields (EMF)

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

## Environment

Do not throw away the appliance with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment.

## Introduction

Do you often feel that skincare products can not penetrate well into the skin and stay on the surface? After using loads of creams, lotions, serums everyday, it still does not give you the ideal benefits as you want.

The new Philips VisaBoost Absorption Booster Device introduces a unique innovative ultrasound technology working together with vibrating massage which encourages better absorption of skincare products in the deep layer of the skin.

The technology has been especially designed to help skincare products be better absorbed into the deep layer of skin, where they can achieve maximum effect for the appearance of healthy, moist, dewy skin today and tomorrow.


The ultrasound waves instantly dissolve skincare products into micro-molecules which are more easily transported into your skin's epidermis. It can stimulate your skin's metabolism and encourages natural cell renewal process. Combined with vibration, it gently massages the skin and stimulates micro-circulation, which further boosts absorption of skincare products and delivers a healthy glow continuously.

The device is innovated to complement your current skincare routine. It is recommended to use after applying a serum, lotion or cream. Before application, please clean and gently dry your face until it is completely dry; during application, please apply your regular skincare products on the probe and also your face as step instruction, in order to better absorb skincare essence. Once you start to use the VisaBoost, you may notice an immediate improvement in your skin's appearance. After long-term usage, your skin becomes more healthy and radiant.

We hope you and your skin will enjoy using this VisaBoost Absorption Booster Device. For maximum results, we recommend to complement your regular facial treatment routine by using the VisaPure Cleansing device and VisaCare Firming device.

## General description

- ① Probe
- ② Power-on indicator
- ③ Power on/off button
- ④ Charging indicator
- ⑤ Charging stand
- ⑥ Appliance socket
- ⑦ Charging pins
- ⑧ Handle
- ⑨ Protective cap


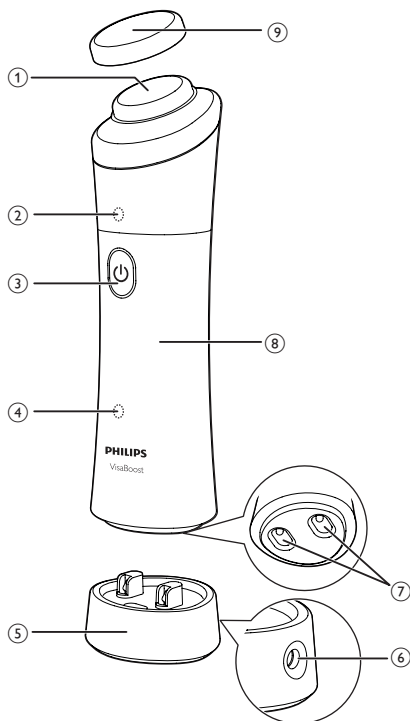
### Auto Shut-off

The appliance will automatically shut off after 10 minutes of continuous use.

## Preparing for use

### Charging the appliance

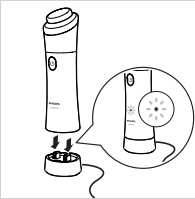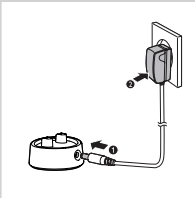**Note:** Make sure the appliance is fully charged before you use it for the first time.

Charging the appliance takes less than 2.5 hours. A fully charged appliance has a cordless operating time of up to 20 minutes.

- 1** Make sure the appliance is switched off.
- 2** Insert the small plug in the socket of the charging stand, and put the adapter in the socket.
- 3** Put the appliance in the charging stand.

▶ The charging indicator will light up when charging. When the battery is fully charged, the indicator goes out.

**Note:** You can use the appliance before the battery is fully charged, but its performance and operating time will not be optimal in this case.

### Cleansing your face

- 1** Cleanse your face with a suitable facial cleanser. Make sure all dirt and makeup are removed.
- 2** Rinse your face and pat it dry gently with a clean towel.

## Using the appliance


We suggest to use this appliance three times per week for best results.

### Notes:

- Do not put the dry probe on your face when the appliance is on, as it may irritate your skin.
- Make sure to take off the protection cap before applying serum, gel or cream onto the probe.


- 1** Moisten the probe with one or two drops of your usual serum, cream or gel. Evenly distribute the drops with the help of your fingertips.

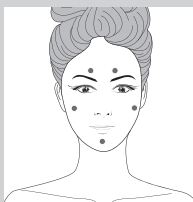- 2** Apply a few drops of your serum, cream or gel around your face as shown, and gently tap it into your skin with the help of your fingertips.

*Tips:*

- For optimal results, do not use oil based skincare products.

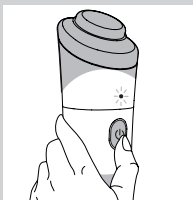
- 3** Press the on/off button once to switch on the appliance.

- 4** Put the probe on your face. Make sure it is fully contacting your face. Gently move the appliance from your chin all the way up to your ears. This motion should take about 5 seconds. To emphasize the result, repeat this step 2-3 times on each side of your cheek. (Refer to “Quick Treatment” and “Emphasized Treatment” in inner leaflet)

- *Tip: To treat the areas where full skin contact is hard to achieve, such as the corner on either side of your nose, you can try different ways of holding the appliance for instance, vertical and horizontal ways, or with the help of the other hand to stretch and hold the skin to be treated. You can apply this method to your chin and forehead areas.*

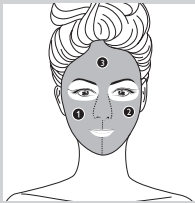
## Notes:

- Do not treat any zone longer than 3 minutes (Refer to the zone separation ① ② ③ as shown).
- It is normal for the probe to feel warmer during use due to continuous vibration.

- 5 To treat your forehead, make overlapping circular movements with either inwards or outwards directions.

## After use:

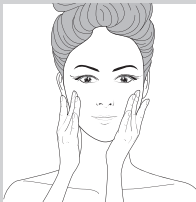- 1 Press the on/off button to turn off the appliance. The appliance is equipped with an auto-shutoff function, it will be shut off automatically after ten minutes of continuous use.
- 2 Tap extra serum or gel on your face with your fingertips. Alternatively use your palms to gently lift both sides of your cheek two to three times.



## Cleaning and maintenance

### Cleaning the appliance

## Notes:

- Make sure the appliance is switched off before you start cleaning.
- Make sure the adapter is unplugged when you clean it. Only wipe the adapter with a dry cloth.
- Always detach the charging stand from the adapter before you clean it.
- Never use scouring pads, abrasive cleaning agents or aggressive liquids such as alcohol or acetone to clean the appliance.
- Clean the appliance and the probe thoroughly after each use to ensure optimum performance.

## 9 ENGLISH

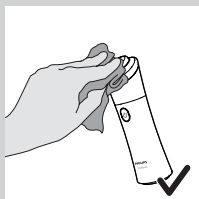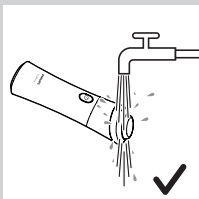- The probe and front end of the appliance is waterproof and can be rinsed under tap water, but never rinse the whole appliance under water.
- Make sure all parts are dry before you reconnect to the adapter.

- 1** Rinse only the probe under a lukewarm tap. You can also use a mild facial cleanser to clean the probe.
- 2** Wipe these parts with a clean and dry towel.
  - *Tip: For easy cleaning, you can also use a damp cloth to remove residue from the probe.*

### Storage


- 1** Store the appliance with the protective cap in a dry and dust free place. You can also keep the packaging for future storage.
- 2** Ensure the protective cap is in place to protect the appliance from damage when storing or traveling.

### Discard the rechargeable battery


#### Warning:

- Do not attempt to replace the rechargeable battery.
- Do not reconnect the appliance to the mains cord after you have removed the rechargeable battery.


Remove the rechargeable battery before discarding the appliance to help preserve the environment.

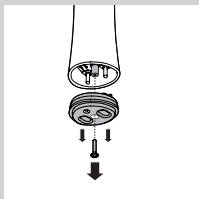
- 1** Remove the appliance from the charging stand, switch on the appliance and let it operate until it switches off automatically. Repeat this action until the rechargeable battery is empty.
- 2** Loosen and release the screw and detach the bottom part off the appliance.

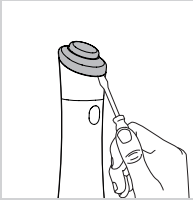- 3** Insert a small flathead screwdriver between the golden color ring and the appliance. Pry loose the top part off the appliance.

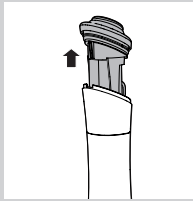- 4** Carefully pull the top part with the battery holder out of the housing.

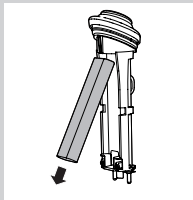- 5** Remove the rechargeable battery from the battery holder, and dispose of it at an official collection point for batteries.

## Guarantee and service

If you need service or information or if you have a problem, please visit the Philips website at [www.philips.com.cn/support](http://www.philips.com.cn/support) or contact the Philips Consumer Care Centre in your country. You find its phone number in the worldwide guarantee leaflet. If there is no Consumer Care Centre in your country, go to your local Philips dealer.

## Troubleshooting

This chapter summarises the most common problems you could encounter with the appliance. If you are unable to solve the problem with the information below, visit [www.philips.com.cn/support](http://www.philips.com.cn/support) for a list of frequently asked questions or contact the Consumer Care Centre in your country.

Problem	Solution
I do not know if the appliance is suitable for use on my skin.	Do not use the appliance on dry, chapped skin, open wounds, or when you suffer from skin diseases or skin irritation, such as severe acne, sunburn, skin infection, etc. Do not use the appliance if you are taking steroid-based medication.
The appliance does not charge.	<p>Make sure the socket to which you connect the appliance is live. Check whether the voltage indicated on the appliance corresponds to the local mains voltage. Check whether contact pins on the appliance are in proper contact with the contact points in the stand.</p> <p>When the indicator on the appliance still does not light up or when the appliance still does not charge, take it to your Philips dealer or the nearest Philips Service Centre.</p>
The appliance does not work anymore.	Make sure that you have pressed the on/off button properly. Charge the appliance according to the instructions in this user manual. Verify whether there is a power failure and if the wall socket is live. Check if the charging light on the appliance lights up to make sure the appliance is charging. If it does not light up or if the appliance still does not work, take it to your Philips dealer or the nearest Philips Service Centre.

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© 2016 Koninklijke Philips N.V.