
PHILIPS

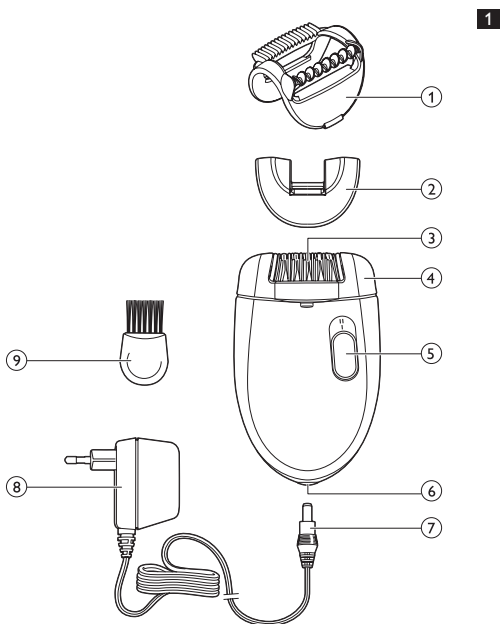ENGLISH

Introduction

Congratulations on your purchase and welcome to Philips! To fully benefit from the support that Philips offers, register your product at www.philips.com/welcome.

General description (Fig. 1)

- 1 Hair-lifting attachment with massage element
- 2 Sensitive area cap
- 3 Epilating discs
- 4 Epilating head
- 5 On/off slide
 - O = off
 - I = normal speed
 - II = high speed (only to be used for epilating head)
- 6 Socket for appliance plug
- 7 Small plug
- 8 Adapter
- 9 Cleaning brush
- Storage pouch (not shown)

Important

Read this user manual carefully before you use the appliance and save it for future reference.

Danger

- Keep the appliance and the adapter dry.
- Do not use the appliance near or over a washbasin or bath filled with water (Fig. 2).
- Do not use the appliance in the bath or in the shower (Fig. 3).
- If you use the appliance in the bathroom, do not use an extension cord (Fig. 4).

Warning

- Only use the appliance in combination with the adapter supplied.
- Do not use the appliance or the adapter if it is damaged.
- If the adapter is damaged, always have it replaced with one of the original type in order to avoid a hazard.
- The adapter contains a transformer. Do not cut off the adapter to replace it with another plug, as this causes a hazardous situation.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.

Caution

- Only epilate the underarms and bikini line with the sensitive area cap attached.
- This appliance is only intended for removing women's body hair on areas below the neck. Do not use it for any other purpose.
- To prevent damage and injuries, keep the **operating** appliance (with or without attachment) away from scalp hair, eyebrows, eyelashes, clothes, threads, cords, brushes etc.
- Do not use the appliance on irritated skin or skin with varicose veins, rashes, spots, moles (with hairs) or wounds without consulting your doctor first. People with a reduced immune response or people who suffer from diabetes mellitus, haemophilia or immunodeficiency should also consult their doctor first.
- Your skin may become a little red and irritated the first few times you use the epilator. This phenomenon is absolutely normal and quickly disappears. As you use the appliance more often, your skin gets used to epilation, skin irritation decreases and regrowth becomes thinner and softer. If the irritation has not disappeared within three days, we advise you to consult a doctor.
- Do not use the epilating head and attachments if they are damaged or broken, as this may cause injury.

Electromagnetic fields (EMF)

This Philips appliance complies with all standards regarding electromagnetic fields (EMF).

Using the appliance

Note: Clean the epilating head before first use.

Epilation tips

- Epilation is easier just after taking a bath or shower. Make sure your skin is completely dry when you start epilating.
- Make sure your skin is clean, entirely dry and free from grease. Do not use any cream before you start epilating.
- Epilation is easier and more comfortable when the hair has the optimum length of 3–4 mm. If the hairs are longer, we recommend that you either shave first and epilate the shorter regrowing hairs after 1 or 2 weeks, or that you use a trimmer with comb to pre-cut hairs to this ideal length.
- When you epilate for the first time, we advise you to epilate in the evening, so that any reddening that develops can disappear overnight.
- On some areas, e.g. the underarms, hair may grow in different directions. In this case, it may be helpful to move the appliance in different directions to achieve optimal results.

Connecting the appliance

- 1 Insert the appliance plug into the appliance (Fig. 5).
- 2 Insert the adapter into the wall socket.

Epilating the legs

We advise you to start epilating with the hair-lifting attachment with massage element. This attachment ensures optimal hair lifting, relaxes the skin and minimises the pulling sensation of epilation.

Only use the hair-lifting attachment on your legs.

- 1 Put the hair-lifting attachment on the epilating head (Fig. 6).

The procedure for placing the attachments onto the appliance is the same for all attachments.

- 2 Select the desired speed to switch on the appliance (Fig. 7).

Preferably use speed II. Switch to speed I for areas with little hair growth, for hard-to-reach areas and for areas where bones are located directly beneath the skin, such as knees and ankles.

- 3 Stretch the skin with your free hand to make the hairs stand upright. Then place the epilator on the skin at an angle of 90° with the on/off slide pointing in the direction in which you are going to move the appliance (Fig. 8).

- 4 Move the appliance slowly over the skin, against the direction of hair growth.

Press the appliance lightly onto the skin.

Epilating the underarms and bikini line

Your epilator comes with a sensitive area cap that reduces the number of active epilating discs.

- 1 Thoroughly clean the area to be epilated to remove any residues of deodorant, cream etc. Then dab the area dry.
 - 2 Put the sensitive area cap on the epilating head (Fig. 9).
 - 3 Switch on the appliance.
 - 4 Stretch the skin with your free hand to make the hairs stand upright. Then place the epilator on the skin at an angle of 90° with the on/off slide pointing in the direction in which you are going to move the appliance.
 - 5 Move the appliance slowly over the skin, against the direction of the hair growth.
- Epilate your bikini line as shown in the figure (Fig. 10).
 - Epilate your underarms as shown in the figure. Move the appliance in different directions (Fig. 11).

Cleaning and maintenance

Note: Make sure the appliance is switched off and unplugged before you start to clean it.

Never use scouring pads, abrasive cleaning agents or aggressive liquids such as petrol or acetone to clean the appliance.

- Never rinse the appliance or the adapter under the tap (Fig. 12).
- Keep the appliance and the adapter dry (Fig. 13).

Cleaning the epilating head and attachments

- 1 Remove any attachment from the epilating head (Fig. 14).
- 2 Push the epilating head in the direction of the arrows on the back of the appliance (1) and remove it (2) (Fig. 15).
- 3 Remove loose hairs from the epilating head and attachments with the cleaning brush (Fig. 16).
- 4 Rinse the epilating head and attachments under the tap while turning them (Fig. 17).
- 5 Shake the epilating head and attachments firmly. Then dry them thoroughly with a tea towel (Fig. 18).

Note: Make sure the epilating head and attachments are completely dry before you place them back onto the appliance.

Ordering accessories

To purchase accessories for this appliance, please visit our online shop at www.shop.philips.com/service. If the online shop is not available in your country, go to your Philips dealer or a Philips service centre. If you have any difficulties obtaining accessories for your appliance, please contact the Philips Consumer Care Centre in your country. You find its contact details in the worldwide guarantee leaflet.

Environment

- Do not throw away the appliance with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment (Fig. 19).

Guarantee and service

Guarantee

- Philips Consumer Lifestyle guarantees its products for a period of two years after the date of purchase. If any defect due to faulty materials and workmanship occurs within this two-year guarantee period, Philips Consumer Lifestyle will repair or replace the product at its expense.

- Philips Consumer Lifestyle will only pay for repair or replacement if convincing proof is provided, for instance by means of a receipt, that the day on which service is claimed is within the guarantee period.
- The guarantee does not cover products and/or product parts that are subject to wear, that can be considered as consumable parts by their nature or that are made of glass.
- The guarantee is not valid if a defect is due to damage caused by incorrect use, poor maintenance (e.g. parts blocked by scale) or if alterations or repairs have been carried out by persons not authorised by Philips Consumer Lifestyle.
- For the appliance to be used correctly, the user should strictly adhere to all instructions included in this user manual and should abstain from any actions or uses that are described as undesired or which are warned against in this user manual.
- For some products, the guarantee may be subject to certain restrictions, for instance as a result of the type of material used.
- These guarantee restrictions do not affect your statutory rights.

Service

Service during and after the guarantee period can be obtained in all countries where the product is officially distributed by Philips. In countries where the product is not distributed by Philips, service is provided by the national Philips organisation. In this case there may be some delay if the required spare parts are not readily available. If you need service or information or if you have a problem, please visit the Philips website at www.philips.com/support or contact the Philips Consumer Care Centre in your country. You find its contact details below. If there is no Consumer Care Centre in your country, go to your local Philips dealer or contact an authorised Philips service centre. Addresses and phone numbers of these service centres can be found in the separate repair-shop leaflet.

Contact details Philips Consumer Care Centre

- 0800 701 0203 (Demais localidades)
- (11) 2121 0203 (Grande Sao Paulo)

Troubleshooting

This chapter summarises the most common problems you could encounter with the appliance. If you are unable to solve the problem with the information below, contact the Consumer Care Centre in your country.

Problem	Cause	Solution
Poor epilating results.	You move the appliance in the wrong direction.	Move the appliance against the direction of hair growth, with the on/off slide pointing in the direction in which you move the appliance.

Problem	Cause	Solution
	You place the epilating head on the skin at the wrong angle.	Make sure that you place the appliance on the skin at an angle of 90° (see chapter 'Using the appliance').
	The hairs are too short.	Sometimes hairs are too short to be caught. The optimal length for easy removal of the hairs is 3-4mm. The hairs that are not caught during this session will be long enough to be removed when you epilate the next time. You can also try to catch the hairs by moving the appliance across the skin a few more times.
The appliance does not work.	The socket to which the appliance is connected is not live.	Make sure the socket is live. If you have connected the epilator to the socket in a bathroom cabinet, you may need to switch on the bathroom light to activate the socket.

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