ENGLISH 4

FRANÇAIS 12

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Introduction

The Cellesse Vitesse is the latest product in the Philips Cellesse range. Clinical tests show that the deep, pleasant massage by its improved massage system can help you reduce or prevent unattractive cellulite more effectively than ever before. More importantly, with this new product the treatment time for each leg has been reduced by 50%.

This booklet explains what causes cellulite and describes how the Cellesse Vitesse can help you solve this typically feminine problem.

Important: Before you start using the appliance, please take the time to read this booklet thoroughly (especially the section headed 'Read carefully before using Cellesse Vitesse').

Cellulite - a typically feminine problem.

Cellulite is the name given to bumpy, pitted skin that arises when lumpy fat deposits situated just under the skin become visible. Because of its appearance, it's often called "orange peel skin".

Cellulite only affects certain parts of the body, particularly the hips, thighs and buttocks.

Cellulite is the result of two things: the mattress-like structure of women's skin and the tendency to retain fluid in the body.


► The natural mattress-like structure of women's skin

Men and women have different skin structures. Women's skin naturally has a mattress-like structure. Fat is deposited in 'pockets' that stand at right angles to the surface, giving a structure that resembles a mattress with springs. This mattress-like structure becomes clearly visible when the skin is squeezed. In men's skin, by contrast, fat is deposited in 'pockets' that tend to run horizontally under the skin. As a result of their skin structure, all women run the risk of developing cellulite.

► Fluid retention

In most cases, the mattress-like structure becomes visible as a result of excess fluid in the fat layer. Women tend to retain more fluid than men, especially in their legs. Excess fluid may accumulate around the fat cells and cause the cellulite layer to swell. If this fluid retention is allowed to continue for months or even years, the fat deposits may thicken and harden. Eventually they'll become hard lumps which can be seen on the surface of the skin even when it's not squeezed.

When it reaches an advanced stage, cellulite may be painful when the skin is squeezed or touched. At that point, cellulite is more difficult to reverse. At earlier stages, however, its effects can be considerably reduced.

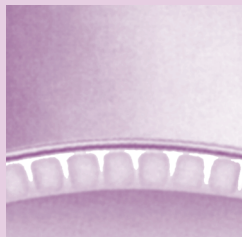Cellulite - what can you do about it?

Once developed, cellulite won't disappear or diminish on its own. On the contrary, it's likely to increase. The only solution is to make the fight against cellulite part of your regular personal care routine and to begin treatment right away, before the problem gets worse.

The good news is that cellulite can be reduced. It will never disappear completely, however, and will return if you stop doing something about it.

The Cellesse Vitesse: the cellulite massage system that's been proven to work!

Philips has therefore specially developed the compact and convenient the Cellesse Vitesse to fight cellulite in an even more effective and comfortable way. The Cellesse Vitesse gets right down to the cellulite layer within the skin and can easily be used in the privacy of your own home.

Cellesse Vitesse massages the cellulite layer in your skin in an intensive, pleasant way. This massage mobilises the excess fluid, which is then removed naturally by the body through the normal circulation and lymphatic drainage systems. The result is a reduction in the bumpiness of the cellulite layer and therefore a smoother skin.

The Cellesse Vitesse in action

The Cellesse Vitesse has been specially developed for those parts of the body that are prone to cellulite (thighs, hips and buttocks). The appliance works through a massage system that's been shown in tests to be highly effective - Vacumotion™.

► The 3-in-1 Vacumotion™



The 3-in-1 Vacumotion™ system works by providing a continuous lifting and rolling movement, which massages the cellulite layer in an intensive and pleasant way. The 3-in-1 Vacumotion™ system consists of four rollers with a suction opening in between. Because the elements are electronically coordinated, the 3-in-1 action is integrated into a smooth, continuous massage.




► 1. Roll-up

The front roller, made of soft-touch material and with small raised pads, gently picks up the skin and guides it into the suction bell.

► 2. Suction

The vacuum lifts the skin even more to form a fold, enabling the massage to reach the cellulite layer effectively.

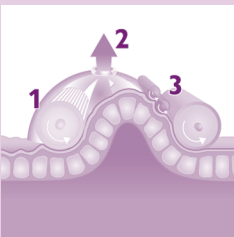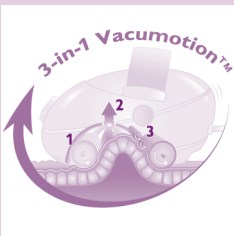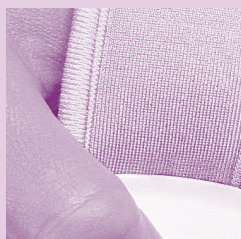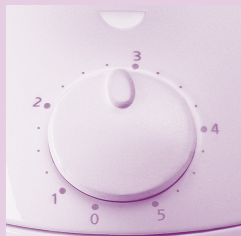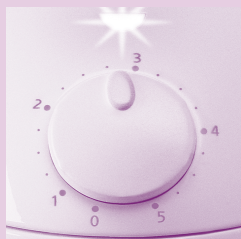► 3. Massage

Four rollers, working together, give the cellulite layer a pleasant and deep massage, removing the main cause of cellulite - the excess fluid surrounding the fat cells. The rear rollers then gently guide the skin away.

► Sensor control and indicator light

To make the treatment as effective as possible, the Vacumotion system will only initiate its suction action if the Cellesse Vitesse is in proper contact with the skin. A special sensor 'feels' whether this is the case, and makes a small indicator light on the front of the appliance go on to show that proper contact has been realised. In this way you can always be sure that you're using the Cellesse Vitesse correctly.

► Intensity settings

You can determine the intensity of the massage yourself, using the continuously variable control. Remember that the massage should always feel pleasant. If you experience the treatment as painful, reduce the intensity to a comfortable level.

► Vacumotion Stop

Vacumotion Stop is a simple and safe way to interrupt or stop your treatment at any time. Press the Vacumotion Stop button and the suction action will stop. You can then lift the appliance off your skin without any effort.

Following an effective treatment programme for the best results

We recommend you to use the Cellesse Vitesse three times a week, spending only 7.5 minutes on each leg each time. Clinical tests have shown that this type of regular treatment can significantly smoothen the cellulite layer in the skin within one month. Naturally, results will vary from individual to individual, depending on how bad the cellulite is to start with and on the individual's lifestyle.

Don't stop the treatment! If you do, the cellulite layer will start to swell again after a while, as a result of the natural tendency of women's bodies to retain fluid. We recommend you to use the appliance three times a week during the first two or three months. After this time you can reduce the frequency and follow a maintenance programme of 15 minutes once a week or 5 minutes three times a week. This should be enough to keep your fluid retention under control.

Proven effective in clinical tests

Although the undesirable consequences of cellulite are visible on the surface, its causes lie deep inside the skin. Independent experts at the

University of Brussels carried out two types of clinical tests to establish the effectiveness of the 3-in-1 Vacumotion™ system.

- ▶ **One type of test showed that using the 3-in-1 Vacumotion™ system only three times a week smoothenes the cellulite layer inside the skin by up to 49% in just one month.**

These tests used ultrasound scanning to measure the bumpiness of the cellulite layer - the cause of cellulite.

The results achieved by women using the Cellesse Vitesse were compared with those of a control group of non-users.

Proven to visibly reduce cellulite

Additional tests showed that the Cellesse Vitesse is also highly effective on the outside of the skin, producing visibly smoother skin.

The extensive tests carried out by the University of Brussels employed a technique known as fringe projection topography. This uses a computer to visualise the surface of the skin and measures the reduction of cellulite.

READ CAREFULLY BEFORE USING THE CELLESSE VITESSE

Do not use the Cellesse Vitesse on hypersensitive skin, problem skin or areas with skin ailments, irritated skin, varicose veins, scars, spots or wounds. Please consult your doctor before using the Cellesse Vitesse if you're pregnant or have a heart condition, diabetes or circulation problems.

Only use the appliance to massage the outside of your thighs, hips and buttocks. Don't use it on any other parts of the body!

Side-effects during and after use

After using the Cellesse Vitesse, your skin may become slightly red and feel rather dry. Don't worry: these are the normal side-effects of a massage as the blood circulation improves. They should soon disappear. If necessary, you can apply a body lotion or cream after using the Cellesse Vitesse to moisturise the skin. If these effects do not disappear within 24 hours or if any other undesirable skin reaction occurs, stop using the appliance and consult your doctor.


If you are prone to bruising or have small red veins clearly visible, you should only use higher settings with care, especially during the first few weeks of use. Whenever necessary, you can always press the Vacumotion Stop button to immediately lift the appliance off your skin.

For best results, the appliance should only be used as described on the following pages.


Important

- ▶ Before using the appliance for the first time, read the section 'Read carefully before using the Cellesse Vitesse'.
- ▶ This appliance is meant for household use only.
- ▶ Check if the voltage indicated on the mains plug corresponds to the mains voltage in your home before you connect the appliance.
- ▶ Use the appliance only with the cord supplied.
- ▶ Do not use the appliance on a wet skin.
- ▶ The appliance has been specially designed for treatment of thighs, hips and buttocks. Do not use the appliance on the upper arms, the inside of the thighs, the belly or the breasts.
- ▶ Keep the appliance away from water!
- ▶ Keep the appliance away from children.
- ▶ Regular cleaning ensures optimum results and a longer life for the appliance.
- ▶ If the mains cord of this appliance is damaged, it must always be replaced by a cord of the original type. Please contact your Philips dealer or the Philips organisation in your country for a replacement.


Preparing the Cellesse Vitesse for use

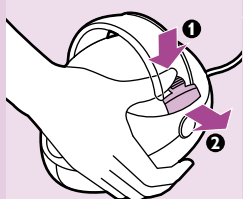
- 1** You can adjust the strap to the size of your hand. Press the adjustment button on the inside of the strap with your thumb and push the clasp out of the appliance.

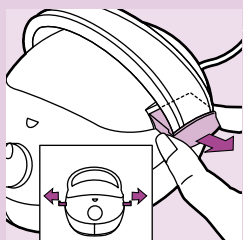- 2** Pull the lug at the end of the hand-strap and remove the pin of the strap from its slots.

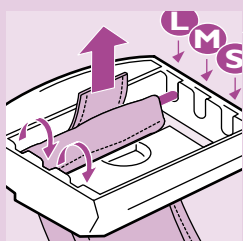
- 3** Press the pin into one of the other two slots. There are three settings: L for large hands, M for normal hands and S for small hands.

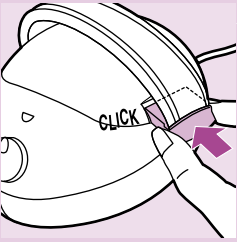- 4** Slide the clasp back into the appliance until it snaps home. You can do the same with the other end of the strap to tighten or loosen the strap even more. This way you can select the strap size that fits your hand best.

Using the Cellesse Vitesse

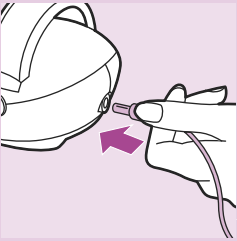
- 1** Make sure your skin is clean, dry and free of grease or other skincare products.

Do not apply any creams or lotions to the skin prior to the massage.

- 2** Make sure that clothes or long hair do not get caught between the rollers when using the Cellesse Vitesse. Always tie back long hair before you start using the appliance.


- 3** Insert the appliance plug into the appliance.

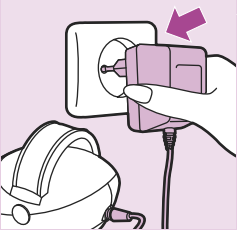- 4** Put the transformer plug in the wall socket.


- 5** Slide your hand under the hand strap and hold the Cellesse Vitesse by clasping your fingers around the back of it (with your fingertips pointing towards the socket for the appliance plug).

- 6** Use the Cellesse Vitesse while you are sitting down, standing up or lying down.


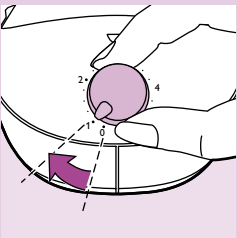
- 7** Before starting, set the on/off switch to the preferred setting. Higher settings will provide more suction power and therefore intensify the massage.

We advise you to start the first treatment on setting 1. Later, either during the first treatment or when your skin has become used to the appliance, you can gradually increase the intensity. If treatment becomes painful, decrease the intensity. You will achieve the best results when you experience the massage as pleasant. Since everyone's skin is different, you will need to discover which intensity setting suits you best.

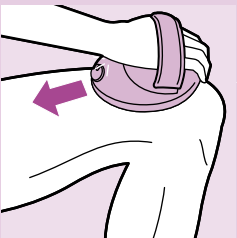- 8** Gently position the appliance on your skin, a little above the knee. Make sure the on/off switch is pointing in the direction of the upper part of your body.

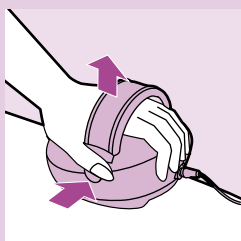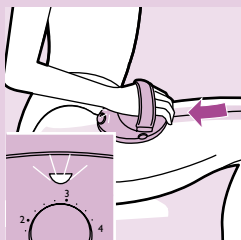- 9** Wait until the indicator light on the front of the appliance comes on. This is the signal that the Cellesse Vitesse is in proper contact with the skin and is ready for use.
- 10** Always move the appliance in one direction only: upwards towards the heart.
- 11** Move the appliance slowly and gradually from a little above the knee to the upper part of your thigh.
The Vacumotion massage system will lift the skin, massaging the cellulite layer in an intensive, effective and pleasant way.
- 12** Let your hand follow the movement of the appliance.
A complete stroke from just above the knee to the upper part of the thigh will take approximately five seconds.
- 13** When you have finished a stroke, press the Vacumotion Stop button. The indicator light goes out and the Vacumotion action is temporarily switched off.
The Vacumotion Stop system is particularly convenient when you are using the appliance at higher intensities.
- 14** Lift the appliance off your skin.
- 15** Place the Cellesse Vitesse back onto your skin, right next to where you began the first stroke, and start a new stroke.
- 16** Repeat this procedure until you have treated your entire thigh, except the inside, with a series of slow, upward, adjacent strokes.


- ▮ If you wish, you may apply a cream, lotion, talcum powder or other type of skincare product to your skin after using the Cellesse Vitesse.
- ▮ The recommended frequency for initial treatment is three times a week, 7.5 minutes on each leg. Read the section headed 'Following an effective treatment programme'.

Cleaning and maintenance

Never immerse the appliance in water or rinse it!

- ▮ For the best results, clean the appliance after every treatment with a moist cloth and, if necessary, a drop of washing-up liquid. Make sure the filter of the suction pump does not get damp.


Do not use alcohol, acetone, petrol, abrasives, scourers, etc., to clean the appliance.

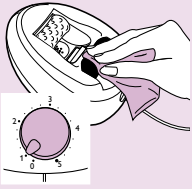- ▶ When cleaning the appliance, pay special attention to the rollers. Select setting I to make the rollers start turning and then clean them with a moist cloth as indicated.
- ▶ The hand-strap (with clasps) can be washed in a washing machine at max. 40°C.

Washing may cause the colour of the strap to change slightly.

Storage

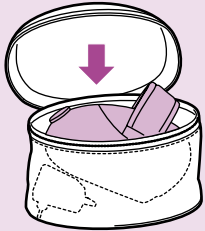- ▶ After use, remove the transformer plug from the wall socket and the small appliance plug from the appliance and store the Cellesse Vitesse together with the cord in the pouch supplied.

Guarantee & service

If you need information or if you have a problem, please visit the Philips Web site at www.philips.com or contact the Philips Customer Care Centre in your country (you will find its phone number in the worldwide guarantee leaflet). If there is no Customer Care Centre in your country, turn to your local Philips dealer or contact the Service Department of Philips Domestic Appliances and Personal Care BV.

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